



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2981	1 CardioGuide Patient App Graphical User Interface		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2982	Document Number: 91033813		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI7	2 Purpose		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2695	This document provides a visual reference to the CardioGuide Patient Application Graphical User Interface (GUI). This document defines the visual aspects of the display including front size, format, color, and placement. The specific text content and behavior of Components, Panels, and Screens is provided in this module.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI8	3 Scope		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2696	This document applies to the CardioGuide Patient Application as a part of the CardioMEMS HF system interacts with the Cloud, and Merlin.net Portal.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI9	4 Reference Documents		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2697	91016560 DRx - myCardioMEMS Patient Application Software Requirements		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2700	91033817 Heart Failure Patient App Translatable Strings		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI10	5 Document Conventions		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2701	This section defines the terms and format used for documenting GUI related information and design. This format is designed to make it easy to document common GUI elements without duplication. GUI elements often appear in multiple locations. If separate requirements are used to document the same element, these specifications eventually get out of sync which may lead to defects.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2703	GUI specifications need to provide enough information to determine what is displayed. This informs software engineers on what to implement and verification engineers on what to test. At a high level, text specifications perform this function. Text specifications are used to indicate which GUI elements are shown in each display and the content of each GUI element. Text is not used to convey exact visual characteristics such as font, placement, shape and color. These details are documented in visual representations provided in the UI Wireframe objects.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3664	5.1 English (US) Strings Attribute		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3665	Some line items in the GUI Specification display strings as part of their expected behavior. The intent of this attribute is to display the English (US) string that is detailed within the GUI element description as well as linked with a GUI element. This allows for easy review of expected behavior for the development team. If the development team needs to review other languages that are supported by the Patient Application, they will need to review the listed string in the Heart Failure Patient App Translatable Strings (91033817) document.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3660	5.2 Object Type Attribute		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3662	Each line item in the GUI Specification includes a 'Object Type' attribute. The intent of this attribute is to demonstrate what kind of information is displayed within the line item. The attribute can take on values including the following: Heading - no information, used for document organization Info - information related to reading the GUI Specification UI Info - information detailing GUI element behavior UI Design - images of GUI elements		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3661	5.3 Status Attribtute		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3663	Each line item in the GUI Specification includes a 'Status' attribute. The intent of this attribute is to track the documented cross-functional review status for each line item. The attribute can take on values including the following: New - item has never been reviewed by the applicable cross-functional team Pending - item has been reviewed by the applicable cross-functional team and awaits further modification before approval Modified - item has been reviewed by the applicable cross-functional team, has since been changed, and awaits a review to accept the changes Approved - item has been reviewd by the applicable cross-fucntional team and has been approved as expected design behavior Deleted - item has been reviewed by the applicable cross-functional team and has been removed as it is no longer part of the GUI Specification		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3650	5.4 Project Attribute		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3651	Each line item in the GUI Specification includes a 'Project' attribute. The intent of this attribute is to demonstrate which version of the Patient Application the line item applies to. The attribute can take on values of CloudDRx 1.0, ClouDRxRTR, and Future or a combination of these values. Note: Future allows for any design work that is currently out of scope for CloudDRx to be documented to reduce the need for redesign at a later date.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3652	6 Defined Terms		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3653	Element - an element is the generic term to describe any part of the GUI display.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3654	Screen - representation of the complete GUI display. Only one Screen is active at any one time. A collection of Screens is used to define the complete GUI Specification. This allows the functionality of the GUI to be divided into logically collected areas. Behaviors can cause the system to transition from one Screeen to another. The text specifications for a Screen list the contained Elements and can also describe behaviors for the Screen (these are under a 'Display' header). The details of spacial placement within the Screen are shown in the UI Design objects (under a 'Layout' header). Screens are assigned a unique name that is in all capital letters ends with the word SCREEN.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3655	Panel - section of related elements. Panels may contain other Panels. This allows for cases where the heirarchy makes the organization clearer. In general, however, Panels contain Components. The purpose of Panels is best explained by considering cases where they can be used to simplify requirements. The details of spacial placement within the SCREEN are shown in the UI Design objects under a 'Layout' header for a Screen or for a Panel. Panels are assigned a unique name that is in all capital letters and ends with the word PANEL.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3656	Component - an individual GUI element. A Component can be as simple as an icon or as complex as a dynamic graph. For dynamic elements, this includes the source of the information being displayed and any state related changes. The spatial placement and visual details fo Components are shown in the UI Design objects with a <Picture> label within the Component details. Components are allocated to Screens in the text specifications under a 'Elements' header. Components are assiged unique name that is in all capital letters and ends with the word COMPONENT.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3657	Shared Element - an element displayed in at least two locations such as two different Panels, two different Screens, or a Screen and a Panel displays under a different Screen. The shared elements are listed within a given Screen or Panel and point to the Shared Element section at the end of the document. The Shared Elements are organized by Panel or Component and then listed in alphabetical order.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI11	8 User Interface		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3436	8.1 Application Widget		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3437	8.1.1 Application Widget (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3440	The APPLICATION WIDGET is displayed in the user's mobile device home screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3438	8.1.2 Application Widget (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4083	<Picture> APPLICATION WIDGET (example: in the mobile device home screen, iOS) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3439	8.1.3 Application Widget (Elements)		Heading	CloudDRx 1.0

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI4082	Component: APPLICATION WIDGET State: N/A Display: Always displayed in the user's mobile device home screen. Label: <label.gb.cardio_guide> Behavior: Icon, text, and trigger. Upon activation, the APPLICATION WIDGET displays the STARTUP SCREEN.	PString693 {0} N/A - controlled by local OS label.gb.cardio_guide	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4032	<Picture> APPLICATION WIDGET 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI12	8.2 Startup Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI13	8.2.1 Startup Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1159	The STARTUP SCREEN is displayed when the application opens.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1160	The STARTUP SCREEN consists of: OS STATUS PANEL PATIENT APP LOGO ICON COMPONENT PATIENT APP TRADEMARK COMPONENT START REGISTRATION ACTION TRIGGER COMPONENT USER SAME DEVICE ERROR PANEL (See Shared Elements: Panels) LOGIN ACTION TRIGGER COMPONENT ABBOTT LOGO COMPONENT (See Shared Elements: Components) CMEMS INTRODUCTION ACTION TRIGGER COMPONENT VIDEO PLAYER PANEL (See Shared Elements: Components) RECOMMENDED SOFTWARE UPDATE DIALOG PANEL REQUIRED SOFTWARE UPDATE DIALOG PANEL MFA MAXIMUM ATTEMPTS DIALOG PANEL (See Shared Elements:Panels)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

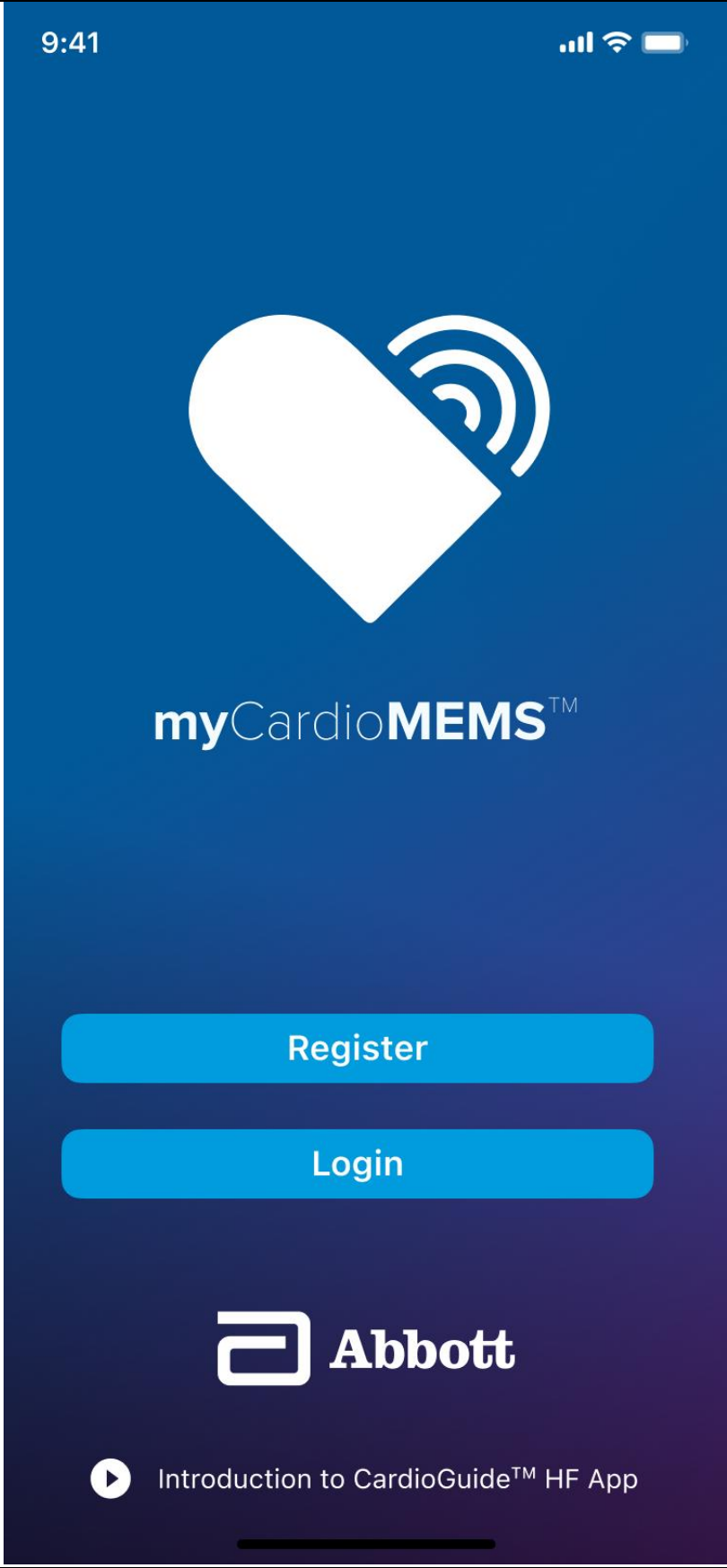
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI14	8.2.2 Startup Screen (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4081	<Picture> STARTUP SCREEN		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

	 The image shows the splash screen of the myCardioMEMS patient application. The background is a dark blue gradient. At the top left, the time is 9:41. At the top right, there are icons for cellular signal, Wi-Fi, and battery. In the center, there is a white heart icon with a Wi-Fi signal symbol inside it. Below the heart, the text 'myCardioMEMS™' is displayed in white. Further down, there are two blue buttons with white text: 'Register' and 'Login'. At the bottom, the Abbott logo is shown in white. Below the logo, there is a play button icon followed by the text 'Introduction to CardioGuide™ HF App'.				
--	---	--	--	--	--

of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI15	8.2.3 Startup Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1163	8.2.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1180	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1164	8.2.3.2 PATIENT APP LOGO ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1533	Component: PATIENT APP LOGO ICON COMPONENT State: N/A Display: Always displayed as a component of the STARTUP SCREEN. Label: N/A Behavior: Icon.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1534	<Picture> PATIENT APP LOGO ICON COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR
PAGUI1165	8.2.3.3 PATIENT APP TRADEMARK COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4079	Component: PATIENT APP TRADEMARK COMPONENT State: N/A Display: Always displayed as a component of the STARTUP SCREEN. Label: <label.gb.cardio_guide> Behavior: Text.	PString693 {0} N/A - controlled by local OS label.gb.cardio_guide	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4078	<Picture> PATIENT APP TRADEMARK COMPONENT 		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1166	8.2.3.4 START REGISTRATION ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1184	Component: START REGISTRATION ACTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the STARTUP SCREEN. Label: <label.gb.register> Behavior: Text and trigger. Upon activation, the START REGISTRATION ACTION TRIGGER COMPONENT validates the user's registration request. If a required software update is available, display REQUIRED SOFTWARE UPDATE PANEL. If validation is not successful, display USER SAME DEVICE ERROR PANEL. If validation is successful, display POLICY SCREEN.	PString30 Register register label.gb.register	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1185	<Picture> START REGISTRATION ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1597	8.2.3.5 USER SAME DEVICE ERROR PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1598	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1173	8.2.3.6 LOGIN ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1206	Component: LOGIN ACTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the STARTUP SCREEN. Label: <label.gb.login> Behavior: Text and trigger. Upon activation, the the LOGIN ACTION TRIGGER COMPONENT displays the SIGN IN POPUP DIALOG PANEL. If a required software update is available, display REQUIRED SOFTWARE UPDATE PANEL.	PString458 Login login label.gb.login	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1207	<Picture> LOGIN ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1541	8.2.3.7 ABBOTT LOGO COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1542	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3389	8.2.3.8 CMEMS INTRODUCTION ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI3390	Component: CMEMS INTRODUCTION ACTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the STARTUP SCREEN. Label: <label.gb.cmems_intro> Behavior: Text and trigger. Upon activation, the CMEMS INTRODUCTION ACTION TRIGGER COMPONENT displays the VIDEO PLAYER PANEL and the associated video.	PString663 Introduction to {0} introductionText label.gb.cmems_intro	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3391	<Picture> CMEMS INTRODUCTION ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3392	8.2.3.9 VIDEO PLAYER PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3393	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3396	8.2.3.10 RECOMMENDED SOFTWARE UPDATE PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3397	8.2.3.10.1 RECOMMENDED SOFTWARE UPDATE PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3400	The RECOMMENDED SOFTWARE UPDATE PANEL is displayed when the user opens the Patient Application and there is a recommended software update available. The RECOMMENDED SOFTWARE UPDATE PANEL is displayed over the screen, does not displace any elements, and any part of the screen that is not covered by the RECOMMENDED SOFTWARE UPDATE PANEL is greyed out.			CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3401	The RECOMMENDED SOFTWARE UPDATE PANEL is composed of: RECOMMENDED SOFTWARE UPDATE LABEL COMPONENT SKIP UPDATE ACTION TRIGGER COMPONENT UPDATE ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3398	8.2.3.10.2 RECOMMENDED SOFTWARE UPDATE PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3402	<Picture> RECOMMENDED SOFTWARE UPDATE PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	There is an update available for your CardioGuide HF App. Please update the app at your earliest convenience. If you need assistance contact technical support at 1 8776963754. Skip Update			
PAGUI3399	8.2.3.10.3 RECOMMENDED SOFTWARE UPDATE PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3403	8.2.3.10.3.1 RECOMMENDED SOFTWARE UPDATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3404	Component: RECOMMENDED SOFTWARE UPDATE LABEL COMPONENT State: N/A Display: Always displayed as a component of the RECOMMENDED SOFTWARE UPDATE PANEL. Label: <msg.gb.sw_update_recommended> Behavior: Text.	PString664 There is a software update available for your device. Please download the new software at your earliest convenience. soft_update_message msg.gb.sw_update_recommended PString1023 There is an update available for your CardioGuide HF App. Please update the app at your earliest convenience. If you need assistance contact technical support at {16}. soft_update_message msg.gb.sw_update_recommended	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3405	<Picture> RECOMMENDED SOFTWARE UPDATE LABEL COMPONENT		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

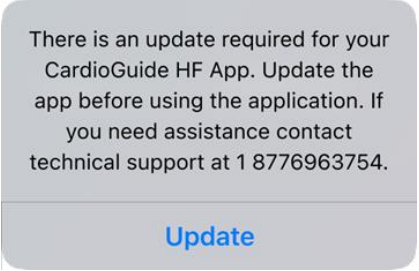


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	There is an update available for your CardioGuide HF App. Please update the app at your earliest convenience. If you need assistance contact technical support at 1 8776963754.			
PAGUI3524	8.2.3.10.3.2 SKIP UPDATE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3525	Component: SKIP UPDATE ACTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the RECOMMENDED SOFTWARE UPDATE PANEL. Label: <msg.gb.skip> Behavior: Text and trigger. Upon activation, the SKIP UPDATE ACTION TRIGGER COMPONENT closes the RECOMMENDED SOFTWARE UPDATE PANEL.	PString294 Skip globalSkip label.gb.skip	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3526	<Picture> SKIP UPADATE ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3406	8.2.3.10.3.3 UPDATE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3407	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3408	8.2.3.11 REQUIRED SOFTWARE UPDATE PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3409	8.2.3.11.1 REQUIRED SOFTWARE UPDATE PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3410	The REQUIRED SOFTWARE UPDATE PANEL is displayed when the user opens the Patient Application and there is a required software update available. The REQUIRED SOFTWARE UPDATE PANEL may also be displayed after the user selects the START REGISTRATION ACTION TRIGGER COMPONENT or the LOGIN ACTION TRIGGER COMPONENT. The REQUIRED SOFTWARE UPDATE PANEL is displayed over the screen, does not displace any elements, and any part of the screen that is not covered by the REQUIRED SOFTWARE UPDATE PANEL is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3411	The REQUIRED SOFTWARE UPDATE PANEL is composed of: REQUIRED SOFTWARE UPDATE LABEL COMPONENT UPDATE ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3412	8.2.3.11.2 REQUIRED SOFTWARE UPDATE PANEL (Layout)		Heading	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3413	<Picture> REQUIRED SOFTWARE UPDATE PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3414	8.2.3.11.3 REQUIRED SOFTWARE UPDATE PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3415	8.2.3.11.3.1 REQUIRED SOFTWARE UPDATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3416	Component: REQUIRED SOFTWARE UPDATE LABEL COMPONENT State: N/A Display: Always displayed as a component of the REQUIRED SOFTWARE UPDATE PANEL. Label: <msg.gb.sw_update_required> Behavior: Text.	PString665 There is a software update for your device. Please download the new software before using the application. hard_update_message msg.gb.sw_update_required PString1024 There is an update required for your CardioGuide HF App. Update the app before using the application. If you need assistance contact technical support at {16}. hard_update_message msg.gb.sw_update_required	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3417	<Picture> REQUIRED SOFTWARE UPDATE LABEL COMPONENT		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



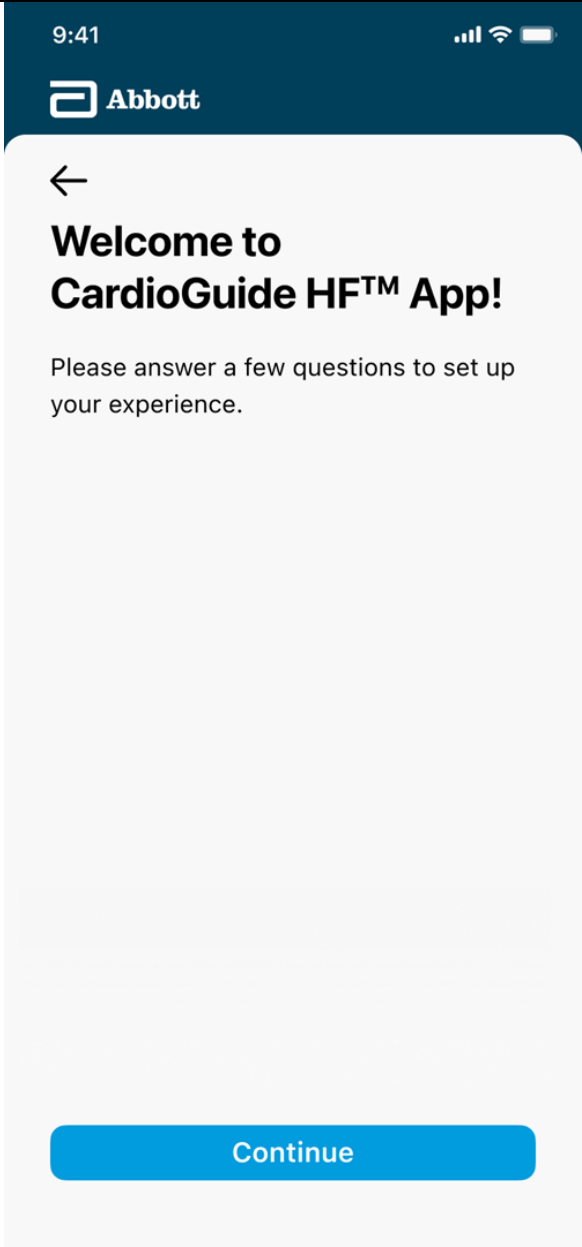
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	There is an update required for your CardioGuide HF App. Update the app before using the application. If you need assistance contact technical support at 1 8776963754.			
PAGUI3418	8.2.3.11.3.2 UPDATE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3419	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4700	8.2.3.11.3.2.1 MFA MAXIMUM ATTEMPTS DIALOG PANEL		Heading	CloudDRxRTR
PAGUI4701	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI1467	8.3 Policy Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1468	8.3.1 Policy Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1469	The POLICY SCREEN is a series of panels that are displayed after the user selects START REGISTRATION ACTION TRIGGER COMPONENT, or after the user selects the LOGIN ACTION TRIGGER COMPONENT, if the user has not reviewed the most up to date Privacy Policy or End User License Agreement.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4077	The POLICY SCREEN consists of: POLICY WELCOME PANEL END USER LICENSE AGREEMENT PANEL SHORT PRIVACY POLICY PANEL FULL PRIVACY POLICY PANEL		UI Info	CloudDRxRTR
PAGUI1471	8.3.2 Policy Screen (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4076	Not applicable since this is a series of screens.		UI Design	CloudDRxRTR CloudDRxPS
PAGUI1473	8.3.3 Policy Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4494	8.3.3.1 POLICY WELCOME PANEL		Heading	CloudDRxRTR
PAGUI4505	8.3.3.1.1 POLICY WELCOME PANEL (Display)		Heading	CloudDRxRTR
PAGUI4499	The POLICY WELCOME PANEL is displayed after the user selects START REGISTRATION ACTION TRIGGER COMPONENT or LOGIN ACTION TRIGGER COMPONENT for the first time after installing the app.		UI Info	CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
	The POLICY WELCOME PANEL is displayed over the entire screen except the OS STATUS PANEL and GLOBAL HEADING PANEL.			
PAGUI4500	The POLICY WELCOME PANEL consists of: OS STATUS PANEL GLOBAL HEADING PANEL BACK ACTION TRIGGER COMPONENT (Shared Elements: Components) POLICY WELCOME HEADER COMPONENT POLICY WELCOME LABEL COMPONENT POLICY WELCOME CONTINUE ACTION TRIGGER COMPONENT		UI Info	CloudDRxRTR
PAGUI4506	8.3.3.1.2 POLICY WELCOME PANEL (Layout)		Heading	CloudDRxRTR
PAGUI4508	<Picture>		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1476	8.3.3.1.3.2 GLOBAL HEADING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1477	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1539	8.3.3.1.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1540	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1478	8.3.3.1.3.4 POLICY WELCOME HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4074	Component: POLICY WELCOME HEADER COMPONENT State: N/A Display: Always displayed as a component of the POLICY SCREEN. Label: <label.gb.welcome_message_cardio_guide> Behavior: Text.	PString694 Welcome to {0}! registerTitleIntro label.gb.welcome_message_cardio_guide	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4075	<Picture> POLICY WELCOME HEADER COMPONENT Welcome to CardioGuide HF™ App!	PString694 Welcome to {0}! registerTitleIntro label.gb.welcome_message_cardio_guide	UI Design	CloudDRxRTR CloudDRxPS
PAGUI1481	8.3.3.1.3.5 POLICY WELCOME LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1482	Component: POLICY WELCOME LABEL COMPONENT State: N/A Display: Always displayed as a component of the POLICY WELCOME PANEL. Label: <label.gb.answer_question> Behavior: Text.	PString459 Please answer a few questions to help us set up your experience. registerIntroPrimaryText label.gb.answer_question PString1014 Please answer a few questions to set up your experience. registerIntroPrimaryText label.gb.answer_question	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1483	<Picture> POLICY WELCOME LABEL COMPONENT Please answer a few questions to set up your experience.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

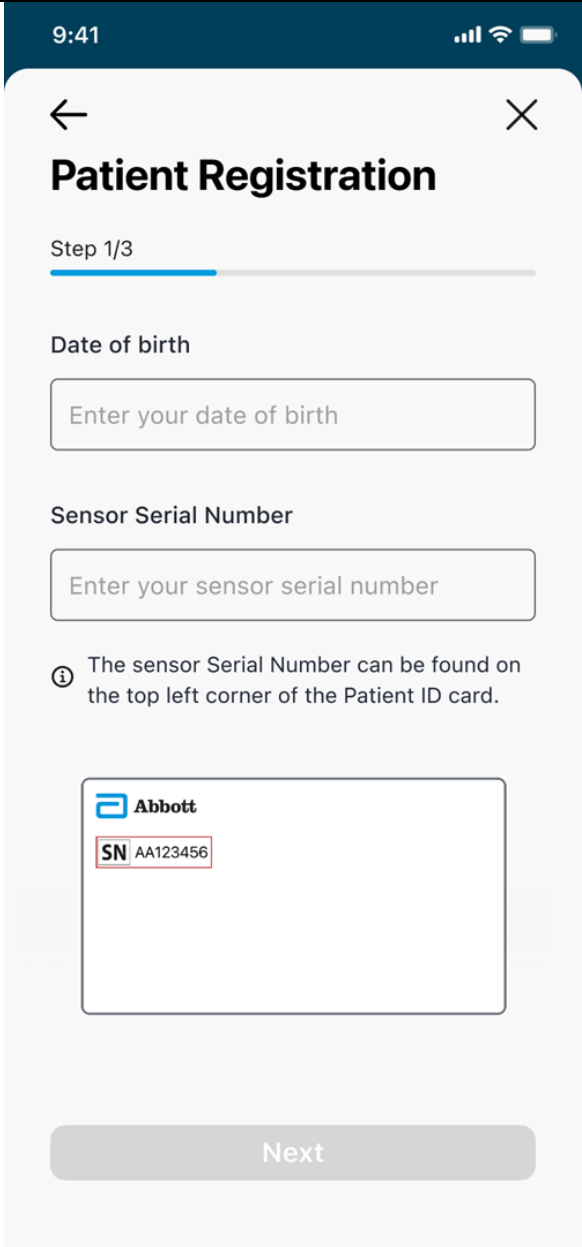
ID		English US Strings	Object Type	Project
PAGUI1487	8.3.3.1.3.6 POLICY WELCOME CONTINUE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1488	Component: POLICY WELCOME CONTINUE ACCTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the POLICY SCREEN. Label: <label.gb.continue> Behavior: Text and trigger. Upon activation, the the POLICY CONTINUE ACTION TRIGGER COMPONENT displays the END USER LICENSE AGREEMENT PANEL.	PString268 Continue continue label.gb.continue	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1489	<Picture> POLICY WELCOME CONTINUE ACTION TRIGGER COMPONENT		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1510	8.3.3.2 END USER LICENSE AGREEMENT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1511	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3900	8.3.3.4 SHORT PRIVACY POLICY PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI4203	See Shared Elements: Panels.			CloudDRxRTR CloudDRxPS
PAGUI3792	8.3.3.6 FULL PRIVACY POLICY PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI3793	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS
PAGUI16	8.4 Registration Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI17	8.4.1 Registration Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1225	The REGISTRATION SCREEN is a series of panels that are used to identify the user and register the application. The REGISTRATION SCREEN is displayed after the user selects the START REGISTRATION ACTION TRIGGER COMPONENT and finishes the POLICY SCREEN.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1226	The REGISTRATION SCREEN consists of: PATIENT IDENTIFICATION PANEL PATIENT INFORMATION PANEL CONNECT HEALTH APP PANEL (Shared Elements: Panels)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4473	The number of steps displayed in these panels will differ between iOS and Android, due to the fact that Androids do not have the AppleConnect step (CONNECT HEALTH APP PANEL). This is described in the REGISTRATION PROGRESS BAR COMPONENT in Shared Elements:Components.		UI Info	CloudDRxRTR
PAGUI18	8.4.2 Registration Screen (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1227	Not applicable - the REGISTRATION SCREEN is a series of panels; layout details for each panel can be seen in Registration Screen (Elements).		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI19	8.4.3 Registration Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1228	8.4.3.1 PATIENT IDENTIFICATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1229	8.4.3.1.1 PATIENT IDENTIFICATION PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1232	The PATIENT IDENTIFICATION PANEL is displayed after the user selects the CREATE NEW ACCOUNT ACTION TRIGGER COMPONENT. The PATIENT IDENTIFICATION PANEL is displayed over the entire screen except the OS STATUS PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1233	The PATIENT IDENTIFICATION PANEL consists of: BACK ACTION TRIGGER COMPONENT (Shared Elements: Components) CLOSE PANEL ACTION TRIGGER COMPONENT (Shared Elements: Components) PATIENT IDENTIFICATION HEADER COMPONENT REGISTRATION PROGRESS BAR COMPONENT (See Shared Elements: Components) DATE OF BIRTH LABEL COMPONENT DATE OF BIRTH INPUT COMPONENT CALENDAR COMPONENT SERIAL NUMBER LABEL COMPONENT SERIAL NUMBER INPUT COMPONENT ALPHABETICAL KEYBOARD COMPONENT (See Shared Elements: Components) SERIAL NUMBER INFORMATION LABEL COMPONENT PATIENT IDENTIFICATION CARD DIAGRAM COMPONENT PATIENT IDENTIFCATION ERROR PANEL NEXT ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1230	8.4.3.1.2 PATIENT IDENTIFICATION PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1234	<Picture> PATIENT IDENTIFICATION PANEL (default: no input, iOS)		UI Design	CloudDRx 1.0

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1235	<Picture> PATIENT IDENTIFICATION PANEL (example: all required input, iOS)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	9:41 ← × Patient Registration Step 1/3 Date of birth 01-01-2000 Sensor Serial Number 112923 ⓘ The sensor Serial Number can be found on the top left corner of the Patient ID card. Abbott SN AA123456 Next			
PAGUI1231	8.4.3.1.3 PATIENT IDENTIFICATION PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

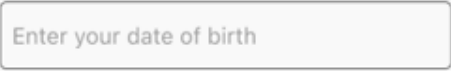


myCardioMEMS Patient Application Graphical User Interface Specification

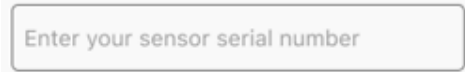
ID		English US Strings	Object Type	Project
PAGUI1236	8.4.3.1.3.1 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1247	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1237	8.4.3.1.3.2 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1248	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1238	8.4.3.1.3.3 PATIENT IDENTIFICATION HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1249	Component: PATIENT IDENTIFICATION HEADER COMPONENT State: N/A Display: Always displayed in the PATIENT IDENTIFICATION PANEL. Label: <label.gb.patient_registration> Behavior: Text.	PString460 Patient Registration registerStep1_title label.gb.patient_registratio n	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1250	<Picture> PATIENT IDENTIFICATION HEADER COMPONENT Patient Registration		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1239	8.4.3.1.3.4 REGISTRATION PROGRESS BAR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1251	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1240	8.4.3.1.3.5 DATE OF BIRTH LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1252	Component: DATE OF BIRTH LABEL COMPONENT State: N/A Display: Always displayed in the PATIENT IDENTIFICATION PANEL. Label: <label.gb.date_birth> Behavior: Text.	PString461 Date of birth dateOfBirthLabel label.gb.date_birth	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1253	<Picture> DATE OF BIRTH LABEL COMPONENT Date of birth		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1241	8.4.3.1.3.6 DATE OF BIRTH INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1254	Component: DATE OF BIRTH INPUT COMPONENT State: Enabled/ Disabled Display: Always displayed in the PATIENT IDENTIFICATION PANEL. Default state is disabled. Label: <label.gb.ender_dob> Behavior: Text and input display. When disabled, the DATE OF BIRTH INPUT COMPONENT displays label.When enabled, the DATE OF BIRTH INPUT COMPONENT displays the input from the CALENDAR COMPONENT.	PString462 Enter your date of birth dateOfBirthplaceholder label.gb.ender_dob	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1255	<Picture> DATE OF BIRTH INPUT COMPONENT (Default: disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1256	<Picture> DATE OF BIRTH INPUT COMPONENT (Example: enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1262	8.4.3.1.3.7 CALENDAR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1553	Component: CALENDAR COMPONENT State: N/A Display: Conditionally displayed in the PATIENT IDENTIFICATION PANEL when the user slects the DATE OF BIRTH INPUT COMPONENT and continues to be displayed until the user selects part of the screen that is not CALENDAR COMPONENT or by selecting the CLOSE CALENDAR ACTION TRIGGER COMPONENT. Label: N/A - determined by local third party vendor package Behavior: Calendar and trigger. The user can select a date and the date selected is displayed in the DATE OF BIRTH INPUT COMPONENT by scrolling through day, month, and year options.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1554	<Picture> CALENDAR COMPONENT (example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3644	8.4.3.1.3.7.1 CLOSE CALENDAR ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

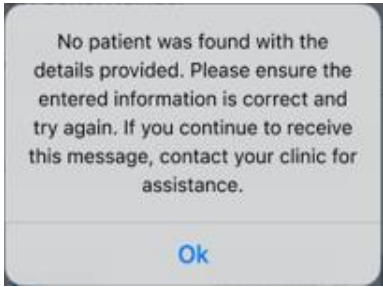
ID		English US Strings	Object Type	Project
PAGUI3645	Component: CLOSE CALENDAR ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the CALENDAR COMPONENT. Label: N/A - determined by local Operating System Behavior: Text and trigger. Upon activation, the CLOSE CALENDAR ACTION TRIGGER COMPONENT closes the CALENDAR COMPONENT and saves the selected date within the DATE OF BIRTH INPUT COMPONENT.	PString204 Done done label.gb.done	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3646	<Picture> CLOSE CALENDAR ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1242	8.4.3.1.3.8 SERIAL NUMBER LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1257	Component: SERIAL NUMBER LABEL COMPONENT State: N/A Display: Always displayed in the PATIENT IDENTIFICATION PANEL. Label: <label.gb.sensor_serial> Behavior: Text.	PString463 Sensor Serial Number ssnLabel label.gb.sensor_serial	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1258	<Picture> SERIAL NUMBER LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1243	8.4.3.1.3.9 SERIAL NUMBER INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1259	Component: SERIAL NUMBER INPUT COMPONENT State: Enabled/ Disabled Display: Always displayed in the PATIENT IDENTIFICATION PANEL. Default state is disabled. Label: <label.gb.enter_sensor_serial> Behavior: Text and input display. When disabled, the SERIAL NUMBER INPUT COMPONENT displays label. When enabled, the SERIAL NUMBER INPUT COMPONENT displays the input from the NUMERIC KEYBOARD COMPONENT and a cursor. Input is automatically capitalized as it is entered.	PString464 Enter your sensor serial number ssnPlaceholder label.gb.enter_sensor_seria l	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1260	<Picture> SERIAL NUMBER INPUT COMPONENT (Default: disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1261	<Picture> SERIAL NUMBER INPUT COMPONENT (Example: enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1244	8.4.3.1.3.10 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

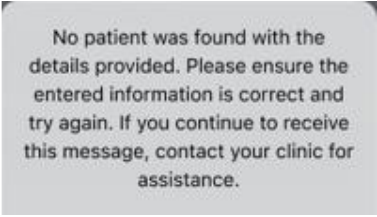
ID		English US Strings	Object Type	Project
PAGUI1264	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1245	8.4.3.1.3.11 SERIAL NUMBER INFORMATION LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1265	Component: SERIAL NUMBER INFORMATION LABEL COMPONENT State: N/A Display: Always displayed in the PATIENT IDENTIFICATION PANEL. Label: <msg.gb.serial_identify> Behavior: Text.	PString43 The sensor Serial Number can be found on the top left corner of the Patient ID card. ssnHelpText msg.gb.serial_identify	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1266	<Picture> SERIAL NUMBER INFORMATION LABEL COMPONENT  The sensor Serial Number can be found on the top left corner of the Patient ID card.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1427	8.4.3.1.3.12 PATIENT IDENTIFICATION CARD DIAGRAM COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1428	Component: PATIENT IDENTIFICATION CARD DIAGRAM COMPONENT State: N/A Display: Always displayed in the PATIENT IDENTIFICATION PANEL. Label: N/A Behavior: Icon.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1429	<Picture> PATIENT IDENTIFICATION CARD DIAGRAM COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1246	8.4.3.1.3.13 NEXT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1267	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1430	8.4.3.1.3.14		Heading	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1431	8.4.3.1.3.14.1 PATIENT IDENTIFICATION ERROR PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1434	The PATIENT IDENTIFICATION ERROR PANEL is displayed after the user selects the NEXT ACTION TRIGGER COMPONENT and the input is not valid. This could be because the serial number is incorrect, the patient is not enrolled on the clinican portal, or another situation in which input is not valid. The PATIENT IDENTIFICATION ERROR PANEL is displayed over the screen and the background is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1435	The PATIENT IDENTIFICATION ERROR PANEL consists of: PATIENT IDENTIFICATION ERROR LABEL COMPONENT OK ACTION TRIGGER COMPONENTET (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1432	8.4.3.1.3.14.2 PATIENT IDENTIFICATION ERROR PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1436	<Picture> PATIENT IDENTIFICATION ERROR PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1433	8.4.3.1.3.14.3 PATIENT IDENTIFICATION ERROR PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1437	8.4.3.1.3.14.3.1 PATIENT IDENTIFICATION ERROR LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1438	Component: PATIENT IDENTIFICATION ERROR LABEL COMPONENT State: N/A Display: Always displayed in the PATIENT IDENTIFICATION ERROR PANEL. Label: <msg.bg.no_patient_found> Behavior: Text.	PString331 No patient was found with the details provided. Please ensure the entered information is correct and try again. If you continue to receive this message, contact your clinic for assistance. invalidText msg.gb.no_patient_found	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1439	<Picture> PATIENT IDENTIFICATION ERROR LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1440	8.4.3.1.3.14.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1441	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1268	8.4.3.2 PATIENT INFORMATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1269	8.4.3.2.1 PATIENT INFORMATION PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1272	The PATIENT INFORMATION PANEL is displayed after the user selects the NEXT ACTION TRIGGER COMPONENT. The PATIENT IDENTIFICATION PANEL is displayed over the entire screen except the OS STATUS PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1273	The PATIENT INFORMATION PANEL consists of: BACK ACTION TRIGGER COMPONENT (Shared Elements: Components) CLOSE PANEL ACTION TRIGGER COMPONENT (Shared Elements: Components) PATIENT INFORMATION HEADER COMPONENT REGISTRATION PROGRESS BAR COMPONENT (See Shared Elements: Components) FULL NAME LABEL COMPONENT FULL NAME INPUT COMPONENT FULL NAME ERROR COMPONENT EMAIL LABEL COMPONENT EMAIL INPUT COMPONENT EMAIL FORMAT ERROR COMPONENT PHONE LABEL COMPONENT PHONE INPUT COMPONENT PHONE ERROR COMPONENT PASSWORD LABEL COMPONENT (See Shared Elements: Components) PASSWORD INPUT COMPONENT PASSWORD ERROR COMPONENT PASSWORD REQUIREMENT LABEL COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR

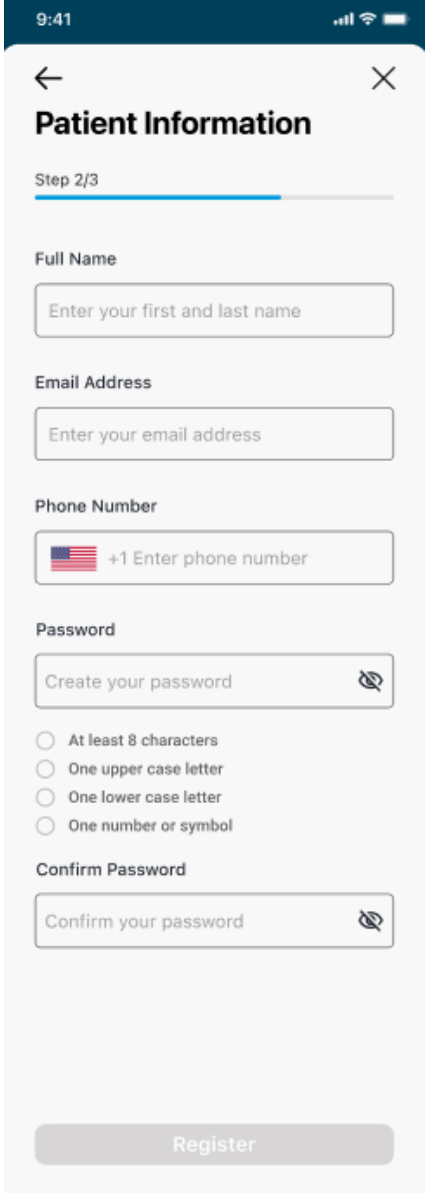
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	CONFIRM PASSWORD LABEL COMPONENT (See Shared Elements: Components) CONFIRM PASSWORD INPUT COMPONENT CONFIRM PASSWORD ERROR COMPONENT ALPHABETIC KEYBOARD COMPONENT (See Shared Elements: Components) NUMERIC KEYBOARD COMPONENT (See Shared Elements: Components) REGISTRATION FAILED ERROR PANEL REGISTRATION ACTION TRIGGER COMPONENT EMAIL IN USE POPUP ERROR COMPONENT USER AUTHENTICATION SCREEN (See User Authentication Screen)			
PAGUI1270	8.4.3.2.2 PATIENT INFORMATION PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4129	<Picture> PATIENT INFORMATION PANEL (default: no input, iOS)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1275	<Picture> PATIENT INFORMATION PANEL (example: completed input, iOS)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

	9:41 ← × Patient Information Step 3/3 Full Name Alexander Stallone Email Address Astallone@gmail.com Phone Number  +1 321-457-5743 Password 🙁 ✓ At least 8 characters ✓ One upper case letter ✓ One lower case letter ✓ One number or symbol Confirm Password 🙁 Register			
--	---	--	--	--

entia document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1271	8.4.3.2.3 PATIENT INFORMATION PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1276	8.4.3.2.3.1 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1277	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1278	8.4.3.2.3.2 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1279	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1280	8.4.3.2.3.3 PATIENT INFORMATION HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1281	Component: PATIENT INFORMATION HEADER COMPONENT State: N/A Display: Always displayed in the PATIENT INFORMATION PANEL Label: <label.gb.patient_information> Behavior: Text.	PString465 Patient Information registerStep2Title label.gb.patient_informatio n	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1282	<Picture> PATIENT INFORMATION HEADER COMPONENT Patient Information		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1283	8.4.3.2.3.4 REGISTRATION PROGRESS BAR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1284	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1285	8.4.3.2.3.5 FULL NAME LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1286	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1288	8.4.3.2.3.6 FULL NAME INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1641	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1556	8.4.3.2.3.7 FULL NAME ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1557	Component: FULL NAME ERROR COMPONENT State: N/A Display: Conditionally displayed in the PATIENT INFORMATION PANEL if the FULL NAME INPUT COMPONENT user input fails validation. Label: <msg.gb.full_name> Behavior: Text.	PString466 Please enter your full name nameFieldRequired msg.gb.full_name	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1558	<Picture> FULL NAME ERROR COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1293	8.4.3.2.3.8 EMAIL LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1294	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1296	8.4.3.2.3.9 EMAIL INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1652	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1559	8.4.3.2.3.11 EMAIL FORMAT ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1560	Component: EMAIL FORMAT ERROR COMPONENT State: N/A Display: Conditionally displayed in the PATIENT INFORMATION PANEL if the EMAIL INPUT COMPONENT user input fails email format validation. Label: <msg.gb.email_format> Behavior: Text.	PString468 Please enter your email in format: example@domain.com emailFieldInvalid msg.gb.email_format	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1561	<Picture> EMAIL FORMAT ERROR COMPONENT		UI Design	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

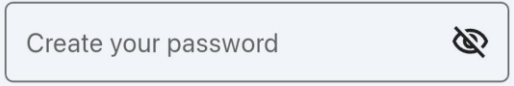
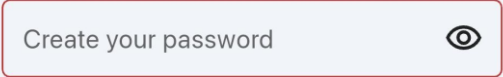
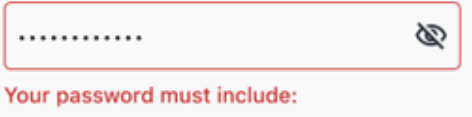


myCardioMEMS Patient Application Graphical User Interface Specification

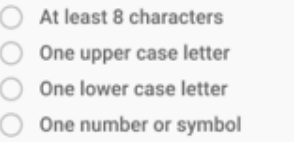
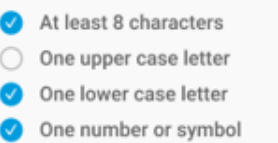
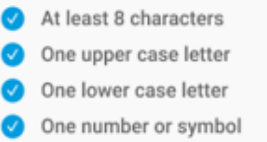
ID		English US Strings	Object Type	Project
	Astallone@gmail.com Please enter your email in format: example@domain.com			CloudDRxPS
PAGUI1301	8.4.3.2.3.12 PHONE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1649	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1304	8.4.3.2.3.13 PHONE INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1650	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1563	8.4.3.2.3.15 PHONE ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1564	Component: PHONE ERROR COMPONENT State: N/A Display: Conditionally displayed in the PATIENT INFORMATION PANEL if the PHONE INPUT COMPONENT user input fails validationi. Label: Behavior: Text.	PString469 Please enter a valid phone number phoneFieldInvalid msg.gb.valid_phone	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1565	<Picture> PHONE ERROR COMPONENT 12 Please enter a valid phone number		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1308	8.4.3.2.3.16 PASSWORD LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1309	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1311	8.4.3.2.3.17 PASSWORD INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

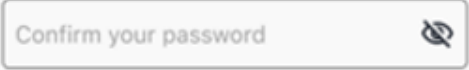
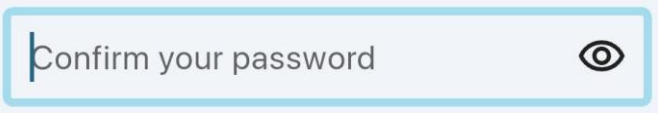
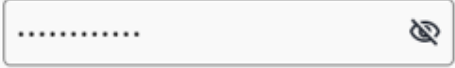
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4031	Component: PASSWORD INPUT COMPONENT State: No Input/ Input and Hide Input/ Show Input Display: Always displayed in the PATIENT INFORMATION PANEL. The default status is no input and hide input. Label: <label.gb.create_your_password> Behavior: Input box. In the no input state, the PASSWORD INPUT COMPONENT displays default text. In the input state, the PASSWORD INPUT COMPONENT displays the input given by the user through the ALPHABETICAL KEYBOARD COMPONENT. The PASSWORD INPUT COMPONENT changes to show input state when the user selects the icon on the right side of the input box. Input is validated as the user provides input to the PASSWORD INPUT COMPONENT.	PString696 Create your password passwordFieldPlaceholder label.gb.create_your_passw ord	UI Info	CloudDRxRTR CloudDRxPS
PAGUI3692	<Picture> PASSWORD INPUT COMPONENT (default: no input, default: hide input) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3693	<Picture> PASSWORD INPUT COMPONENT (default: no input, example: show input) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3694	<Picture> PASSWORD INPUT COMPONENT (example: input, default: hide input) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3695	<Picture> PASSWORD INPUT COMPONENT (example: input, example: show input) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI1448	8.4.3.2.3.18 PASSWORD ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1449	Component: PASSWORD ERROR COMPONENT State: N/A Display: Conditionally displayed in the PATIENT INFORMATION PANEL if PASSWORD INPUT COMPONENT fails validation. Label: <msg.gb.password_format> Behavior: Text.	PString470 Your password must include: invalid msg.gb.password_format	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1450	<Picture> PASSWORD ERROR COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1317	8.4.3.2.3.19 PASSWORD REQUIREMENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1318	Component: PASSWORD REQUIREMENT LABEL COMPONENT	PString471	UI Info	CloudDRx 1.0

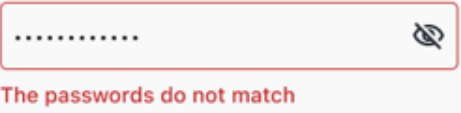
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	State: Met/ Not Met Display: Always displayed in the PATIENT INFORMATION PANEL. The default state is not met. Label: <label.gb.password_characters>, <label.gb.password_cap_characters>, <label.gb.password_lower_characters>, <label.gb.password_special_characters> Behavior: Icon and text. The checkmark icons to the left of each requirement listed changes color as the specific requirement changes state.	At least 8 characters rule1 label.gb.password_characters PString472 One upper case letter rule2 label.gb.password_cap_characters PString473 One lower case letter rule3 label.gb.password_lower_characters PString474 One number or symbol rule4 label.gb.password_special_characters		CloudDRxRTR CloudDRxPS
PAGUI1319	<Picture> PASSWORD REQUIREMENT LABEL COMPONENT (default: not met) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1451	<Picture> PASSWORD REQUIREMENT LABEL COMPONENT (example: partially met) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1320	<Picture> PASSWORD REQUIREMENT LABEL COMPONENT (example: met) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1321	8.4.3.2.3.20 CONFIRM PASSWORD LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1322	See Shared Elements: Components.		Info	CloudDRx 1.0

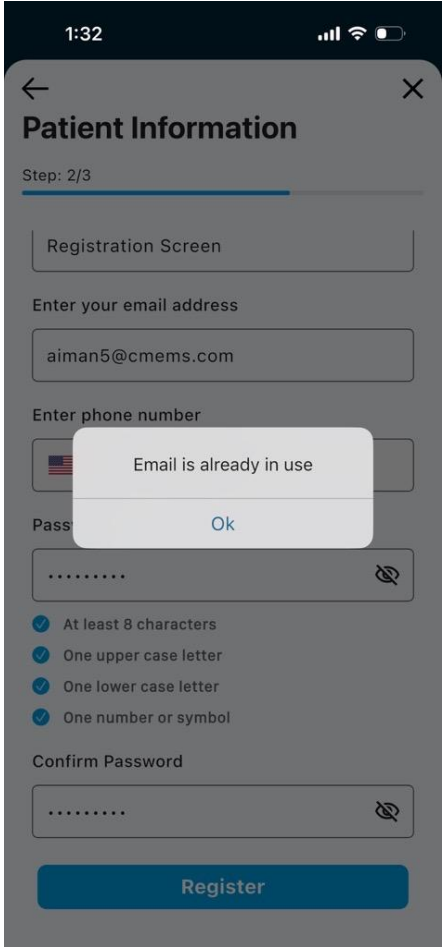
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1324	8.4.3.2.3.21 CONFIRM PASSWORD INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1325	Component: CONFIRM PASSWORD INPUT COMPONENT State: No Input/ Input and Hide Input/ Show Input Display: Always displayed in the PATIENT INFORMATION PANEL. The default status is no input and hide input. Label: <label.gb.confirm_your_password> Behavior: Input box. In the no input state, the CONFIRM PASSWORD INPUT COMPONENT displays default text. In the input state, the CONFIRM PASSWORD INPUT COMPONENT displays the input given by the user through the ALPHABETICAL KEYBOARD COMPONENT. The CONFIRM PASSWORD INPUT COMPONENT changes to show input state when the user selects the icon on the right side of the input box. Input is validated while the user enters input into CONFIRM PASSWORD INPUT COMPONENT.	PString695 Confirm your password confirmPasswordFieldPlace holder label.gb.confirm_your_pass word	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1326	<Picture> CONFIRM PASSWORD INPUT COMPONENT (default: no input, default: hide input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1327	<Picture> CONFIRM PASSWORD INPUT COMPONENT (default: no input, example: show input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1328	<Picture> CONFIRM PASSWORD INPUT COMPONENT (example: input, default: hide input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1329	<Picture> CONFIRM PASSWORD INPUT COMPONENT (example: input, show input, invalid input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1452	8.4.3.2.3.22 CONFIRM PASSWORD ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1453	Component: CONFIRM PASSWORD ERROR COMPONENT State: N/A Display: Conditionally displayed in the PATIENT INFORMATION PANEL. The CONFIRM PASSWORD ERROR COMPONENT is displayed when there is an error associated with the CONFIRM PASSWORD INPUT. Label: <msg.gb.password_not_match> Behavior: Text.	PString476 The passwords do not match confirmPasswordFieldInvali d msg.gb.password_not_mat ch	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

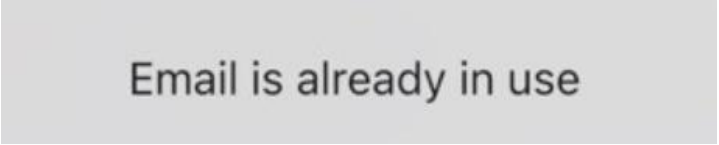
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1454	<Picture> CONFIRM PASSWORD ERROR COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1330	8.4.3.2.3.23 ALPHABETIC KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1331	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1332	8.4.3.2.3.24 NUMERIC KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1333	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1336	8.4.3.2.3.25 REGISTRATION ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1353	Component: REGISTRATION ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Always displayed in the PATIENT INFORMATION PANEL. The default state is disabled. Label: <label.gb.register> Behavior: Text and trigger. The REGISTRATION ACTION TRIGGER COMPONENT changes state when all required input in the PATIENT INFORMATION PANEL has been completed and formatting is valid and no error components are displayed. Upon activation in the enabled state, the REGISTRATION ACTION TRIGGER COMPONENT sends the patient information input in the PATIENT INFORMATION PANEL to the backend and displays the USER AUTHENTICATION SCREEN.	PString30 Register register label.gb.register	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1354	<Picture> REGISTRATION ACTION TRIGGER COMPONENT (default: disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1355	<Picture> REGISTRATION ACTION TRIGGER COMPONENT (enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4832	8.4.3.2.3.26 EMAIL IN USE POPUP ERROR PANEL		Heading	CloudDRxRTR
PAGUI4833	8.4.3.2.3.26.1 EMAIL IN USE POPUP ERROR PANEL (Display)		Heading	CloudDRxRTR
PAGUI4834	The EMAIL IN USE POPUP ERROR PANEL is displayed as a compoenent of the PATIENT INFORMATION PANEL if the email the user is attempting to register with is already in use with another CardioMEMS sensor. The EMAIL IN USE POPUP ERROR PANEL is displayed over the screen and the background is greyed out.		UI Info	CloudDRxRTR
PAGUI4840	The EMAIL IN USE POPUP ERROR PANEL consists of:		UI Info	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	EMAIL IN USE POPUP ERROR MESSAGE COMPONENT OKAY ACTION TRIGGER COMPONENT (Shared Elements: Components)			
PAGUI4835	8.4.3.2.3.26.2 EMAIL IN USE POPUP ERROR PANEL (Layout)		Heading	CloudDRxRTR
PAGUI4836	<Picture> EMAIL IN USE POPUP ERROR PANEL 		UI Design	CloudDRxRTR
PAGUI4837	8.4.3.2.3.26.3 EMAIL IN USE POPUP ERROR PANEL (Elements)		Heading	CloudDRxRTR
PAGUI4838	8.4.3.2.3.26.3.1 EMAIL IN USE POPUP ERROR MESSAGE COMPONENT		Heading	CloudDRxRTR
PAGUI4839	<Picture> EMAIL IN USE POPUP ERROR MESSAGE COMPONENT		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4841	8.4.3.2.3.26.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4842	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI1455	8.4.3.2.3.28 REGISTRATION FAILED PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1456	8.4.3.2.3.28.1 REGISTRATION FAILED PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1459	The REGISTRATION FAILED PANEL is displayed as a component of the PATIENT INFORMATION PANEL if there are unknown registration errors such as re-registering a patient who has already been registered with their serial number and date of birth. The REGISTRATION FAILED PANEL is displayed over the screen and the background is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1460	The REGISTRATION FAILED PANEL consists of: REGISTRATION FAILED LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1457	8.4.3.2.3.28.2 REGISTRATION FAILED PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1461	<Picture> REGISTRATION FAILED PANEL 		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1458	8.4.3.2.3.28.3 REGISTRATION FAILED PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1462	8.4.3.2.3.28.3.1 REGISTRATION FAILED LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1463	Component: REGISTRATION FAILED LABEL COMPONENT State: N/A Display: Always displayed in the REGISTRATION FAILED PANEL	PString477 Registration failed. If in the US, dial {16}. If outside the	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

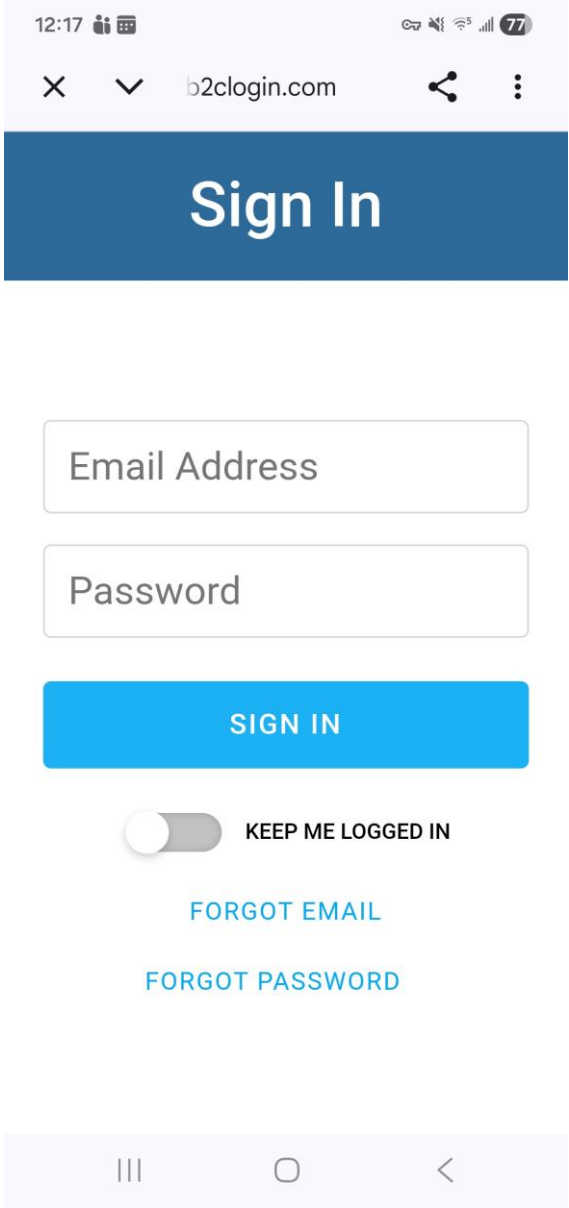
ID		English US Strings	Object Type	Project
	Label: <msg.gb.reg_fail>, <msg.gb.registration_fail> Behavior: Text.	US, contact your Abbott Representative. registrationFailed msg.gb.registration_failed		
PAGUI1464	<Picture> REGISTRATION FAILED LABEL COMPONENT Registration failed. If in the US, dial 1 8776963754. If outside the US, contact your Abbott Representative.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1465	8.4.3.2.3.28.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1466	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4481	8.4.3.2.3.30 USER AUTHENTICATION SCREEN		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4482	See User Authentication Screen		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1338	8.4.3.3 CONNECT HEALTH APP PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1339	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI620	8.5 User Authentication Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI621	8.5.1 User Authentication Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI624	The USER AUTHENTICATION SCREEN is the a series of screens that display upon activation of the LOGIN ACTION TRIGGER COMPONENT or REGISTER ACTION TRIGGER COMPONENT (when registration is successful). The USER AUTHENTICATION SCREEN is a website outside of the application and the user is automatically brought back to the application once login is successful. This section is for reference only.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI625	The USER AUTHENTICATION SCREEN is composed of: LOGIN PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR



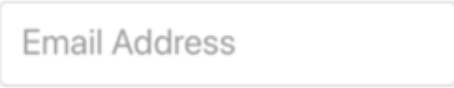
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	TWO FACTOR AUTHENTICATION PANEL (See Shared Elements: Panels)			
PAGUI622	8.5.2 User Authentication Screen (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4569	Not Applicable. This Screen is a series of panels.		Info	CloudDRxRTR
PAGUI623	8.5.3 User Authentication Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4565	8.5.3.1 LOGIN PANEL		Heading	CloudDRxRTR
PAGUI4568	8.5.3.1.1 LOGIN PANEL (Display)		Heading	CloudDRxRTR
PAGUI4571	The LOGIN PANEL displays upon activation of the LOGIN ACTION TRIGGER COMPONENT or REGISTER ACTION TRIGGER COMPONENT when registration is successful. The LOGIN PANEL is a website outside of the application and the user is automatically brought back to the application once login is successful. This section is for reference only.		UI Info	CloudDRxRTR
PAGUI4570	The LOGIN PANEL is composed of: LOGIN HEADER COMPONENT CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) LOGIN EMAIL INPUT COMPONENT LOGIN PASSWORD INPUT COMPONENT ALPHABETICAL KEYBOARD COMPONENT (See Shared Elements: Components) LOGIN ACTION TRIGGER COMPONENT KEEP ME LOGGED IN TRIGGER COMPONENT FORGOT USER ID ACTION TRIGGER COMPONENT FORGOT USER ID PANEL FORGOT PASSWORD ACTION TRIGGER COMPONENT FORGOT PASSWORD PANEL (See Shared Elements: Panels) ENTER CREDENTIALS ERROR PANEL INVALID USER INLINE ERROR USER SAME DEVICE ERROR PANEL (See Shared Elements: Panels) KEEP ME LOGGED IN PANEL		UI Info	CloudDRxRTR
PAGUI4567	8.5.3.1.2 LOGIN PANEL (Layout)		Heading	CloudDRxRTR
PAGUI626	<Picture> LOGIN PANEL (default)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

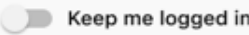
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 A mobile application login screen mockup. At the top, a status bar shows the time 12:17, signal strength, and battery level at 77%. Below the status bar is a browser-like address bar with a close button (X), a dropdown arrow (v), the URL 'b2clogin.com', a share icon, and a menu icon (three dots). The main content area features a large blue 'Sign In' button. Below this are two white input fields with rounded corners, labeled 'Email Address' and 'Password'. Under the password field is a blue 'SIGN IN' button. Below the button is a toggle switch labeled 'KEEP ME LOGGED IN', which is currently turned off. Below the toggle are two links: 'FORGOT EMAIL' and 'FORGOT PASSWORD'. At the bottom of the screen is a white navigation bar with three icons: a hamburger menu (three horizontal lines), a home icon (a circle), and a back arrow (<).			
PAGUI4566	8.5.3.1.3 LOGIN PANEL (Elements)		Heading	CloudDRxRTR
PAGUI634	8.5.3.1.3.1 LOGIN HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI635	Component: LOGIN HEADER COMPONENT State: N/A Display: Always displayed as a component of the LOGIN PANEL Label: N/A - controlled by a third party Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI649	<Picture> LOGIN HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4572	8.5.3.1.3.3 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4573	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI646	8.5.3.1.3.4 LOGIN EMAIL INPUT COMPONENT		Heading	CloudDRxRTR
PAGUI648	Component: LOGIN EMAIL LABEL COMPONENT State: No Input/ User Input Display: Always displayed as a component of the LOGIN PANEL. The default state is no input. Label: N/A - controlled by a third party Behavior: Input box. The EMAIL LABEL COMPONENT displays generic text in the no input state. Upon activation in the no input state, the EMAIL LABEL COMPONENT displays the ALPHABETICAL KEYBOARD COMPONENT. The EMAIL LABEL COMPONENT displays input given by the user through the ALPHABETICAL KEYBOARD COMPONENT.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI651	<Picture> LOGIN EMAIL INPUT COMPONENT (default: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI655	8.5.3.1.3.5 LOGIN PASSWORD INPUT COMPONENT		Heading	CloudDRxRTR
PAGUI656	Component: LOGIN PASSWORD LABEL COMPONENT State: No Input/ User Input Display: Always displayed as a component of the LOGIN PANEL. The default state is no input. Label: N/A - controlled by a third party Behavior: Input box. The PASSWORD LABEL COMPONENT displays generic text in the no input state. Upon activation in the no input state, the PASSWORD LABEL COMPONENT displays the ALPHABETICAL KEYBOARD COMPONENT. The PASSWORD LABEL COMPONENT displays input given by the user through the ALPHABETICAL KEYBOARD COMPONENT.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI657	<Picture> LOGIN PASSWORD INPUT COMPONENT (default: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1581	8.5.3.1.3.6 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

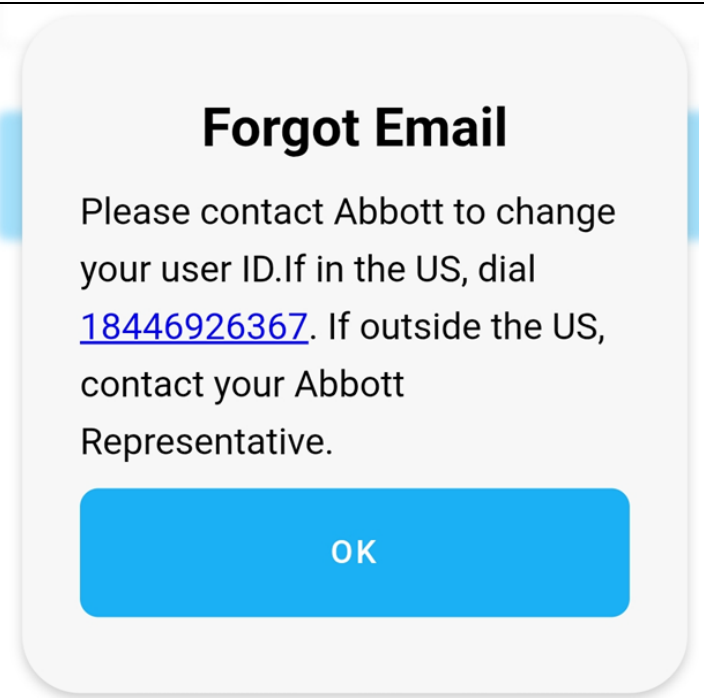
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1582	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI640	8.5.3.1.3.7 LOGIN ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI641	Component: LOGIN ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the LOGIN PANEL. Label: N/A - controlled by a third party Behavior: Text and trigger. Upon activation, the LOGIN ACTION TRIGGER COMPONENT begins the user validation process. Display INVALID USER PANE if the user input does not match an existing user. Display KEEP ME LOGGED IN PANEL if the user input matches an existing user and the KEEP ME LOGGED IN TOGGLE COMPONENT is in the on state. Upon completing the initial validation, the LOGIN ACTION TRIGGER COMPONENT displays the TWO FACTOR AUTHENTICATION PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI684	<Picture> LOGIN ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1365	8.5.3.1.3.8 KEEP ME LOGGED IN TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI639	Component: KEEP ME LOGGED IN TRIGGER COMPONENT State: On/ Off Display: Always displayed in the LOGIN PANEL. The default state is off. Label: <label.gb.keep_me_logged_in> Behavior: Text and trigger. Upon activation in the off state, the KEEP ME LOGGED IN TRIGGLR COMPONENT enables keeping the user logged in and displays the KEEP ME LOGGED IN PANEL. Upon activation in the on state, the KEEP ME LOGGED IN TRIGGER COMPONENT disables keeping the user logged in. The KEEP ME LOGGED IN TRIGGER COMPONENT changes from the on state to the off state if the user selects the CANCEL ACTION TRIGGER COMPONENT in the KEEP ME LOGGED IN PANEL.	PString267 Keep me logged in remember_me label.gb.keep_me_logged_i n	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1367	<Picture> KEEP ME LOGGED IN TRIGGER COMPONENT (example: on) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1368	<Picture> KEEP ME LOGGED IN TRIGGER COMPONENT (default: off) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3291	8.5.3.1.3.9 FORGOT USER ID ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI3292	Component: FORGOT USER ID ACTION TRIGGER COMPONENT State: N/A	PString688 Forgot Email	UI Info	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Display: Always displayed as a component of the LOGIN PANEL. Label: <label.gb.forgot_email> Behavior: Text and trigger. Upon activation, the FORGOT USER ID ACTION TRIGGER COMPONENT displays the FORGOT USER ID PANEL.	label.gb.forgot_email		CloudDRxPS
PAGUI3293	<Picture> FORGOT USER ID ACTION TRIGGER COMPONENT FORGOT EMAIL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3294	8.5.3.1.3.10 FORGOT USER ID PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3295	8.5.3.1.3.10.1 FORGOT USER ID PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3298	The FORGOT USER ID PANEL is conditionally displayed in the LOGIN PANEL after the user selects the FORGOT USER ID ACTION TRIGGER COMPONENT. The panel is displayed as a pop up and any portion of the screen that is not covered by the panel is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3299	The FORGOT USER ID PANEL consists of: FORGOT USER ID HEADER LABEL COMPONENT FORGOT USER ID INSTRUCTIONS COMPONENT OKAY ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3296	8.5.3.1.3.10.2 FORGOT USER ID PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3300	<Picture> FORGOT USER ID PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 A light gray dialog box with rounded corners. It has a title "Forgot Email" in bold black text. Below the title is a paragraph of text: "Please contact Abbott to change your user ID.If in the US, dial 18446926367 . If outside the US, contact your Abbott Representative." The text "18446926367" is a blue hyperlink. At the bottom of the dialog is a blue button with the text "OK" in white.			
PAGUI3297	8.5.3.1.3.10.3 FORGOT USER ID PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4584	8.5.3.1.3.10.3.1 FORGOT USER ID HEADER LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI4585	Component: FORGOT USER ID HEADER LABEL COMPONENT State: N/A Display: Always displayed as a component of the FORGOT USER ID PANEL. Label: <msg.gb.forgot_email> Behavior: Text	PString688 Forgot Email label.gb.forgot_email	UI Info	CloudDRxRTR
PAGUI4586	<Picture> FORGOT USER ID HEADER LABEL COMPONENT  A light gray rectangular box containing the text "Forgot Email" in bold black font.		UI Design	CloudDRxRTR
PAGUI3301	8.5.3.1.3.10.3.2 FORGOT USER ID INSTRUCTIONS COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3302	Component: FORGOT USER ID INSTRUCTIONS COMPONENT State: N/A Display: Always displayed as a component of the FORGOT USER ID PANEL.	PString657 Please contact {22} to reset your user ID. If in the US,	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Label: <msg.gb.forgot_user_id> Behavior: Text.	dial {16}. If outside the US, contact your Abbott Representative. forgotUserIDPasswordBody msg.gb.forgot_user_id		
PAGUI3303	<Picture> FORGOT USER ID INSTRUCTIONS COMPONENT Please contact Abbott to change your user ID.If in the US, dial 18446926367 . If outside the US, contact your Abbott Representative.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3304	8.5.3.1.3.10.3.3 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3305	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1104	8.5.3.1.3.11 FORGOT PASSWORD ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1105	Component: FORGOT PASSWORD ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the LOGIN PANEL. Label: <label.gb.forgot_password> Behavior: Text and trigger. Upon activation, the FORGOT PASSWORD ACTION TRIGGER COMPONENT displays the FORGOT PASSWORD PANEL.	PString34 Forgot Password forgotPassword label.gb.forgot_password	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1106	<Picture> FORGOT PASSWORD ACTION TRIGGER COMPONENT FORGOT PASSWORD		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI664	8.5.3.1.3.12 FORGOT PASSWORD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI665	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1583	8.5.3.1.3.14 ENTER CREDENTIALS IN LINE ERROR		Heading	CloudDRx 1.0 CloudDRxRTR

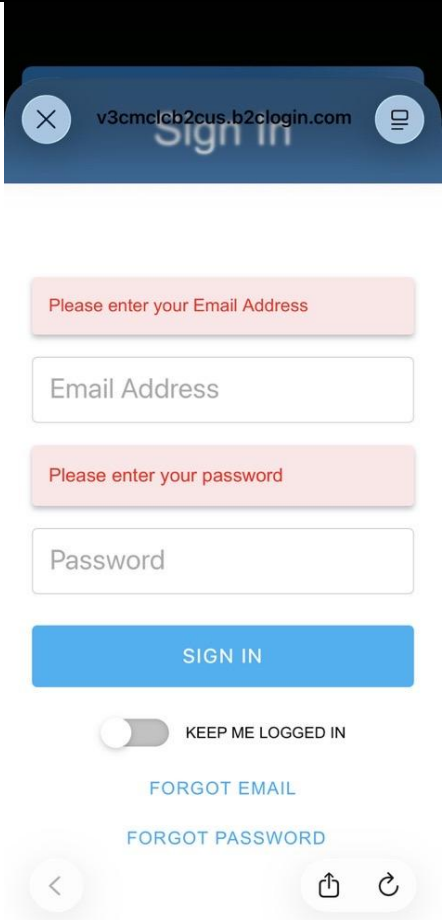
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

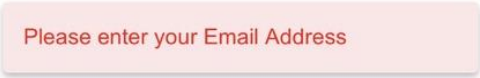
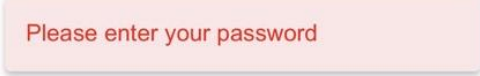
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1584	8.5.3.1.3.14.1 ENTER CREDENTIALS IN LINE ERROR (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1585	The ENTER CREDENTIALS IN LINE ERROR is conditionally displayed in the LOGIN PANEL if the user selects the LOGIN ACTION TRIGGER COMPONENT and at least one of the input fields is empty. The ENTER CREDENTIALS IN LINE ERROR PANEL is displayed within the LOGIN PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1586	The ENTER CREDENTIALS IN LINE ERROR is composed of: ENTER EMAIL ADDRESS ERROR LABEL COMPONENT ENTER PASSWORD ERROR LABEL COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1587	8.5.3.1.3.14.2 ENTER CREDENTIALS IN LINE ERROR (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1588	<Picture> ENTER CREDENTIALS IN LINE ERROR		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

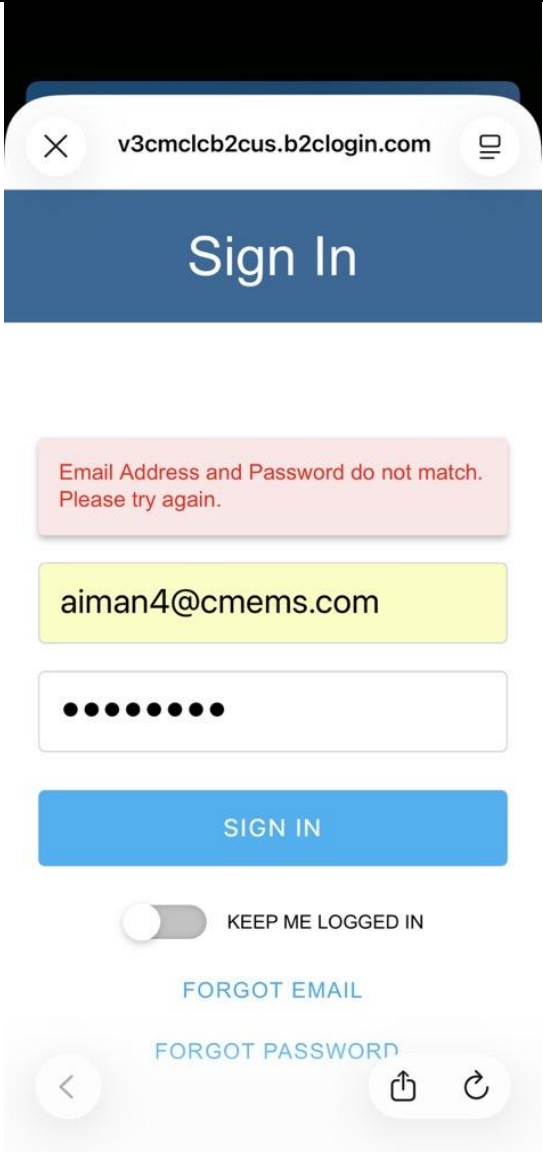
ID		English US Strings	Object Type	Project
				
PAGUI1589	8.5.3.1.3.14.3 ENTER CREDENTIALS IN LINE ERROR (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1590	8.5.3.1.3.14.3.1 ENTER EMAIL ADDRESS IN LINE ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1591	Component: ENTER EMAIL ADDRESS ERROR LABEL COMPONENT State: N/A Display: Displayed when the email address field is empty Label: ots.control.enter_email_address Behavior: Text.	PString78 Please enter both email address and password. requiredField_password requiredField_generic msg.gb.empty_login	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



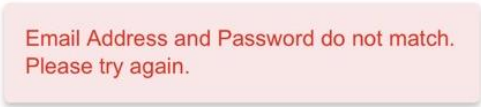
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		PString957 Please enter your Email Address ots.control.enter_email_address		
PAGUI1592	<Picture> ENTER EMAIL ADDRESS ERROR LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1593	8.5.3.1.3.14.3.2 ENTER PASSWORD IN LINE ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1594	Component: ENTER PASSWORD ERROR LABEL COMPONENT State: N/A Display: Displayed when the password field is empty Label: ots.control.enter_password Behavior: Text.	PString958 Please enter your password ots.control.enter_password	Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4216	<Picture> ENTER PASSWORD ERROR LABEL COMPONENT 		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI685	8.5.3.1.3.15 INVALID USER IN LINE ERROR		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI686	8.5.3.1.3.15.1 INVALID USER IN LINE ERROR (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI687	The INVALID USER IN LINE ERROR is conditionally displayed in the LOGIN PANEL if user validation is not sucessful.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI688	The INVALID USER PANEL is composed of: INVALID USER IN LINE ERROR LABEL COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI689	8.5.3.1.3.15.2 INVALID USER IN LINE ERROR (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI690	<Picture> INVALID USER IN LINE ERROR		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

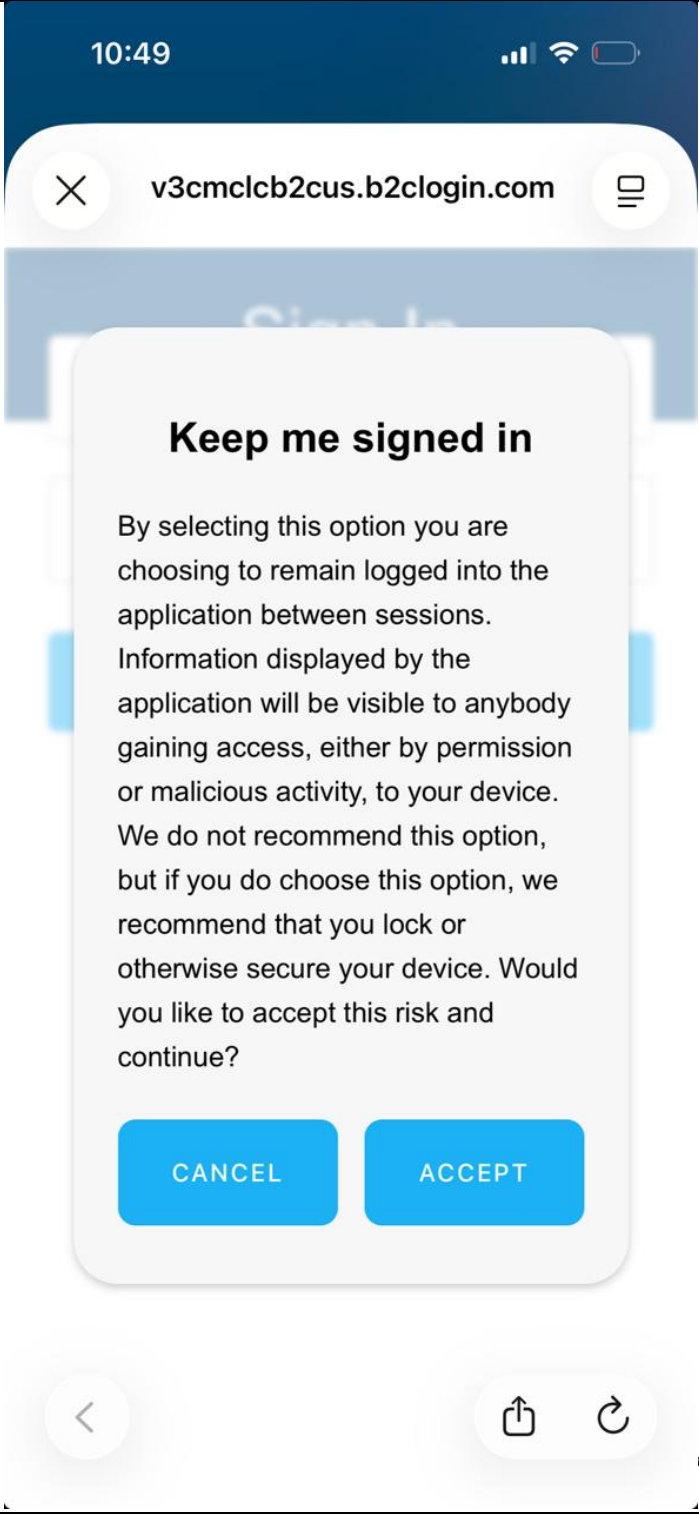
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI691	8.5.3.1.3.15.3 INVALID USER IN LINE ERROR (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI692	8.5.3.1.3.15.3.1 INVALID USER IN LINE ERROR LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI693	Component: INVALID USER LABEL COMPONENT State: N/A Display: Always displayed in the INVALID USER INLINE ERROR. Label: <msg.gb.invalid_credentials> Behavior: Text.	PString490 Email address and Password do not match. Please try again. emailAndPasswordDoNotMatch msg.gb.invalid_credentials	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI694	<Picture> INVALID USER LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1595	8.5.3.1.3.16 USER SAME DEVICE ERROR PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1596	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI697	8.5.3.1.3.17 KEEP ME LOGGED IN PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI698	8.5.3.1.3.17.1 KEEP ME LOGGED IN PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI699	The KEEP ME LOGGED IN PANEL is conditionally displayed in the LOGIN PANEL when the KEEP ME LOGGED IN TOGGLE COMPONENT transitions to the on state. The KEEP ME LOGGED IN PANEL is displayed over the LOGIN PANEL, does not displace any elements, and any part of the LOGIN PANEL that is not covered by the KEEP ME LOGGED IN PANEL is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI700	The KEEP ME LOGGED IN PANEL is composed of: KEEP ME LOGGED IN HEADER LABEL COMPONENT KEEP ME LOGGED IN LABEL COMPONENT CANCEL ACTION TRIGGER COMPONENT ACCEPT ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI701	8.5.3.1.3.17.2 KEEP ME LOGGED IN PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4135	<Picture> KEEP ME LOGGED IN PANEL		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

				
--	--	--	--	--

t is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI703	8.5.3.1.3.17.3 KEEP ME LOGGED IN PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4587	8.5.3.1.3.17.3.1 KEEP ME LOGGED IN HEADER LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI4588	Component: KEEP ME LOGGED IN HEADER LABEL COMPONENT State: N/A Display: Always displayed in the KEEP ME LOGGED IN PANEL. Label: <label.gb.keep_me_logged_in> Behavior: Text.	PString267 Keep me logged in remember_me label.gb.keep_me_logged_i n	UI Info	CloudDRxRTR
PAGUI4589	<Picture> KEEP ME LOGGED IN HEADER LABEL COMPONENT Keep me logged in		UI Design	CloudDRxRTR
PAGUI704	8.5.3.1.3.17.3.2 KEEP ME LOGGED IN LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4136	Component: KEEP ME LOGGED IN LABEL COMPONENT State: N/A Display: Always displayed in the KEEP ME LOGGED IN PANEL. Label: <label.gb.keep_me_logged_in_legal> Behavior: Text.	PString991 By selecting this option you are choosing to remain logged into the application between sessions. Information displayed by the application will be visible to anybody gaining access, either by permission or malicious activity, to your device. We do not recommend this option, but if you do choose this option, we recommend that you lock or otherwise secure your device. Would you like to accept this risk and continue? ots.authentication.keep_me_logged_in_instruction PString998 By selecting this option you are choosing to remain	UI Info	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		logged into the {0} between sessions. Information displayed by the application will be visible to anybody gaining access to your device. keepMeLoggedInBody label.gb.keep_me_logged_in_legal		
PAGUI706	<Picture> KEEP ME LOGGED IN LABEL COMPONENT By selecting this option you are choosing to remain logged into the application between sessions. Information displayed by the application will be visible to anybody gaining access, either by permission or malicious activity, to your device. We do not recommend this option, but if you do choose this option, we recommend that you lock or otherwise secure your device. Would you like to accept this risk and continue?		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI707	8.5.3.1.3.17.3.3 CANCEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI708	Component: CANCEL ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the KEEP ME LOGGED IN PANEL. Label: <label.gb.cancel> Behavior: Text and trigger. Upon activation, the CANCEL ACTION TRIGGER COMPONENT closes the KEEP ME LOGGED IN PANEL and changes the KEEP ME LOGGED IN TRIGGER COMPONENT from the on state to the off state.	PString104 Cancel cancel label.gb.cancel	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI709	<Picture> CANCEL ACTION TRIGGER COMPONENT		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

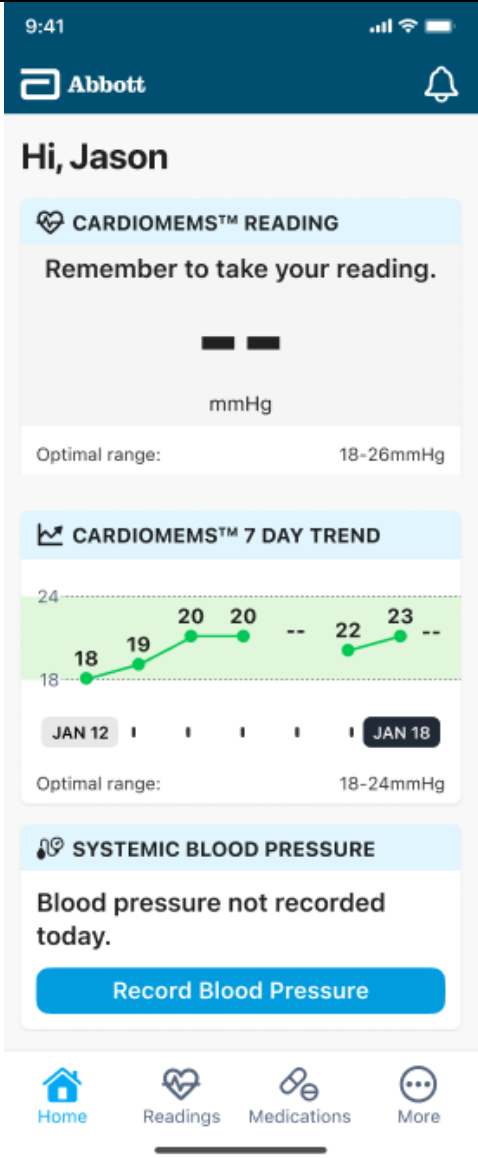
ID		English US Strings	Object Type	Project
				
PAGUI710	8.5.3.1.3.17.3.4 ACCEPT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI711	Component: ACCEPT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the KEEP ME LOGGED IN PANEL. Label: <label.gb.accept> Behavior: Text and trigger. Upon activation, the ACCEPT ACTION TRIGGER COMPONENT displays the LOGIN PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI712	<Picture> ACCEPT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1599	8.5.3.2 TWO FACTOR AUTHENTICATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1600	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI20	8.6 Home Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI21	8.6.1 Home Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI65	The HOME SCREEN is the default screen for the application. Users are sent to this screen once user authentication is completed successfully. The screen is optimized to give key status information for input required from the user.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3968	The HOME SCREEN is composed of: OS STATUS PANEL (See Shared Elements: Panels)		UI Info	CloudDRxRTR CloudDRxPS



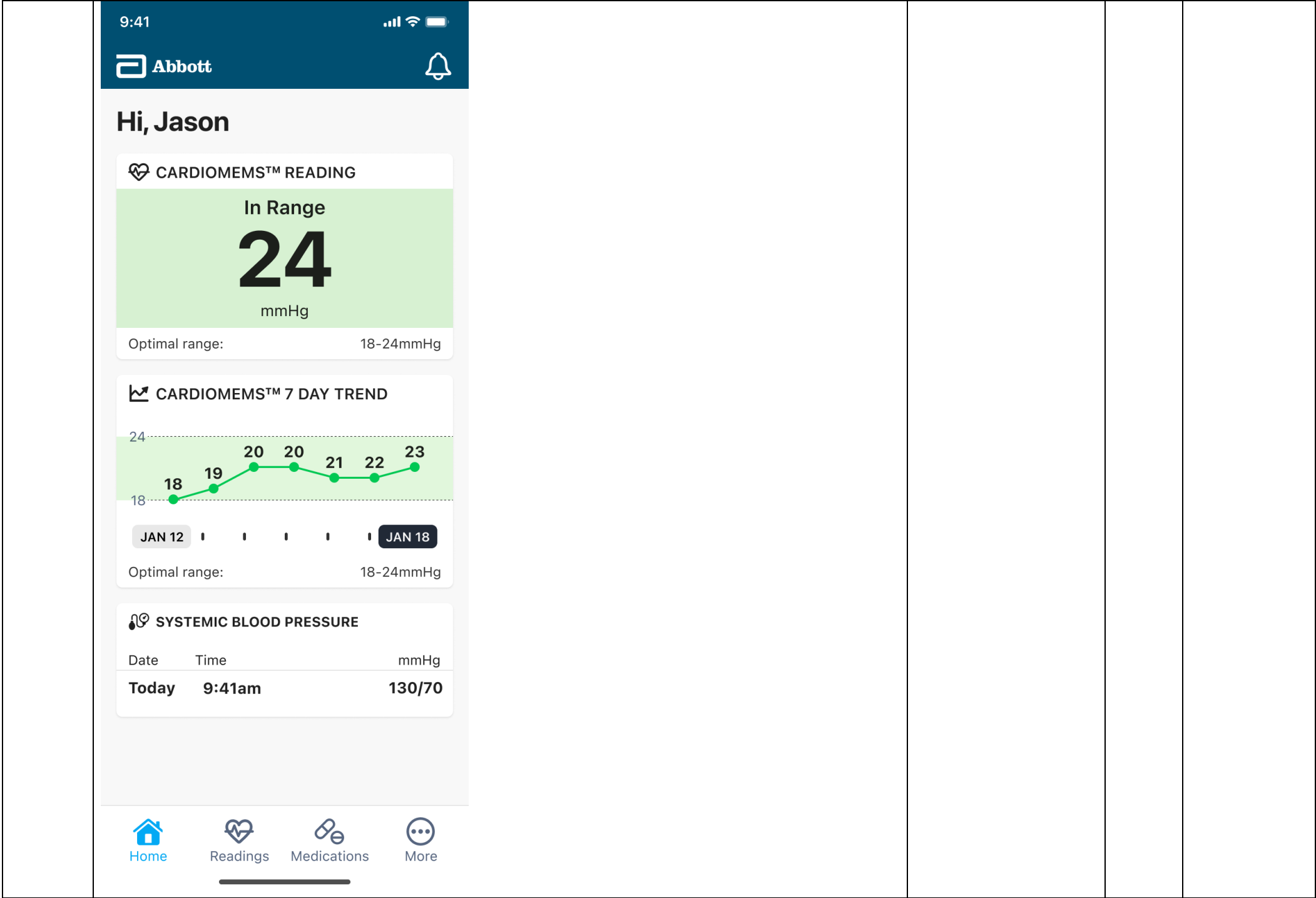
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	GLOBAL HEADING PANEL (See Shared Elements: Panels) HOME LABEL COMPONENT ELEVATION CHANGE MESSAGE CARD PANEL (See Shared Elements: Panels) MEDICATION MESSAGE CARD PANEL (See Shared Elements: Panels) DRX MEDICATION MESSAGE CARD PANEL (See Shared Elements: Panels) CMEMS PAP READING CARD PANEL (See Shared Elements: Panels) CMEMS PAP TREND CARD PANEL (See Shared Elements: Panels) BP READING CARD PANEL (See Shared Elements: Panels) RECORD BP READING PANEL (See Shared Elements: Panels) MEDICATION NOTIFICATION DIALOG PANEL (See Shared Elements: Panels) ELEVATION NOTIFICATION DIALOG PANEL (See Shared Elements: Panels) GLOBAL NAVIGATION PANEL (See Shared Elements: Panels) PA PRESSURE ENABLED PANEL PA PRESSURE DISABLED PANEL			
PAGUI22	8.6.2 Home Screen (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI67	<Picture> HOME SCREEN (default: no readings taken, PA Visibility enabled)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI68	<Picture> HOME SCREEN (readings taken, PA Visibility enabled, Reading in optimal range)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification



This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification					
ID		English US Strings	Object Type	Project	
PAGUI4333	<Picture> HOME SCREEN (PA Visibility disabled, PA reading taken, BP not taken) 9:41  Abbott  Hi, Jason CARDIOMEMS™ READING Reading taken  CARDIOMEMS™ 7 DAY TREND Feb 29 30 31 1 2 3 4        SYSTEMIC BLOOD PRESSURE Blood pressure not recorded today. Record Blood Pressure Home Readings Medications More				CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI23	8.6.3 Home Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS	
PAGUI71	8.6.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS	



myCardioMEMS Patient Application Graphical User Interface Specification

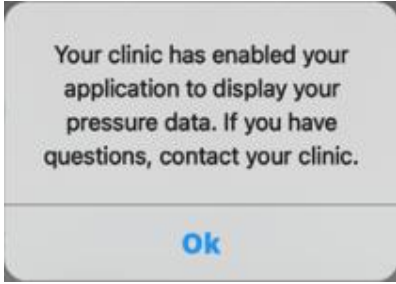
ID		English US Strings	Object Type	Project
PAGUI72	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI73	8.6.3.2 GLOBAL HEADING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI74	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI75	8.6.3.3 HOME LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI76	Component: HOME LABEL COMPONENT State: N/A Display: Always displayed as a component of the HOME SCREEN. Default state is Hi. Label: <label.gb.hi> Behavior: Text. The HOME LABEL COMPONENT displays generic text followed by a name that changes depending on the user profile.	PString421 Good Morning, {5} goodMorning label.gb.good_morning PString422 Good Afternoon, {5} goodAfternoon label.gb.good_afternoon PString423 Good Evening, {5} goodEvening label.gb.good_evening PString1041 Hi, {5} home_title label.gb.hi	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3973	8.6.3.7 ELEVATION CHANGE MESSAGE CARD PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI3974	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS
PAGUI1608	8.6.3.8 MEDICATION MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1609	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3208	8.6.3.9 DRX MEDICATION MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

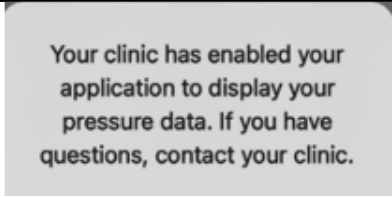
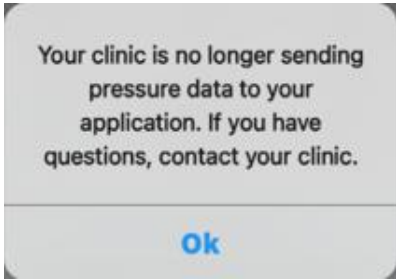
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI3209	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI83	8.6.3.10 CMEMS PAP READING CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI84	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4299	8.6.3.11 CMEMS PAP TREND CARD PANEL		Heading	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4300	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI85	8.6.3.12 BP READING CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI86	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3203	8.6.3.13 RECORD BP READING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3204	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4185	8.6.3.15 MEDICATION NOTIFICATION DIALOG PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI4186	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS
PAGUI4447	8.6.3.16 ELEVATION NOTIFICATION DIALOG PANEL		Heading	CloudDRxRTR
PAGUI4448	See: Shared Elements: Panels		Info	CloudDRxRTR
PAGUI89	8.6.3.18 GLOBAL NAVIGATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI90	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3576	8.6.3.19 PA PRESSURE ENABLED PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3577	8.6.3.19.1 PA PRESSURE ENABLED PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3580	The PA PRESSURE ENABLED PANEL is displayed the first time a user opens the Patient Application after the clinician has enabled the patient to view the PA Pressure values. The PA PRESSURE ENABLED PANEL does not displace any elements and any part of the screen that is not covered by the PA PRESSURE ENABLED PANEL is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3581	The PA PRESSURE ENABLED PANEL consists of: PA PRESSURE ENABLED LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3578	8.6.3.19.2 PA PRESSURE ENABLED PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3582	<Picture> PA PRESSURE ENABLED PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3579	8.6.3.19.3 PA PRESSURE ENABLED PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3583	8.6.3.19.3.1 PA PRESSURE ENABLED LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3584	Component: PA PRESSURE ENABLED LABEL COMPONENT State: N/A Display: Always displayed in the PA PRESSURE ENABLED PANEL. Label: <msg.gb.send_reading_data_switch_on> Behavior: Text.	PString405 Your clinic has enabled your application to display your pressure data. If you have questions, contact your clinic. showPapReadingsOn msg.gb.send_reading_data_switch_on	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3585	<Picture> PA PRESSURE ENABLED LABEL COMPONENT			CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

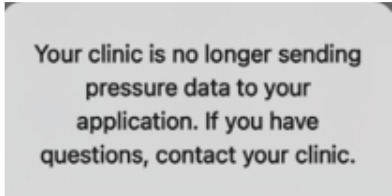
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3586	8.6.3.19.4 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3587	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3588	8.6.3.20 PA PRESSURE DISABLED PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3589	8.6.3.20.1 PA PRESSURE DISABLED PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3590	The PA PRESSURE DISABLED PANEL is displayed the first time a user opens the Patient Application after the clinician has disabled the patient to view the PA Pressure values. The PA PRESSURE DISABLED PANEL does not displace any elements and any part of the screen that is not covered by the PA PRESSURE DISABLED PANEL is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3591	The PA PRESSURE DISABLED PANEL consists of: PA PRESSURE DISABLED LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3592	8.6.3.20.2 PA PRESSURE DISABLED PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3593	<Picture> PA PRESSURE DISABLED PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3594	8.6.3.20.3 PA PRESSURE DISABLED PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3595	8.6.3.20.3.1 PA PRESSURE DISABLED LABEL COMPONENT		Heading	CloudDRx 1.0

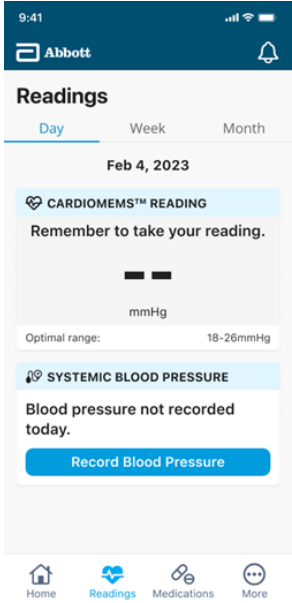
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3596	Component: PA PRESSURE DISABLED LABEL COMPONENT State: N/A Display: Always displayed in the PA PRESSURE DISABLED PANEL. Label: <msg.gb.send_reading_data_switch> Behavior: Text.	PString389 Your clinic is no longer sending pressure data to your application. If you have questions, contact your clinic. showPapReadingsOff msg.gb.send_reading_data_switch	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3597	<Picture> PA PRESSURE DISABLED LABEL COMPONENT 			CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3598	8.6.3.20.4 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3599	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI24	8.7 Readings Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI25	8.7.1 Readings Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI43	The READINGS SCREEN is the application view used to view PAP and BP readings data for the patient. USERS can access the screen by selecting the READINGS SCREEN ICON in the navigation tab. The view is optimized to show historical readings, trend data, and reminders for the patient to take their daily readings.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI44	The READINGS SCREEN consists of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) READINGS LABEL COMPONENT INTERVAL TAB COMPONENT (See Shared Elements: Components) DATE LABEL COMPONENT (See Shared Elements: Components) DATE SELECTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	HISTORICAL DATE SELECTION CALENDAR COMPONENT (See Shared Elements: Components) CMEMS PAP READING CARD PANEL (See Shared Elements: Components) CMEMS PAP TREND CARD PANEL (See Shared Elements: Components) BP READING CARD PANEL (See Shared Elements: Components) GLOBAL NAVIGATION PANEL (See Shared Elements: Panels) BP ERROR PANEL (See Shared Elements: Panels)			
PAGUI26	8.7.2 Readings Screen (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI45	<Picture> READINGS SCREEN (no readings taken today, day view, PA Visibility enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI91	<Picture> READINGS SCREEN (default: no readings taken today, week view, PA Visibility enabled)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

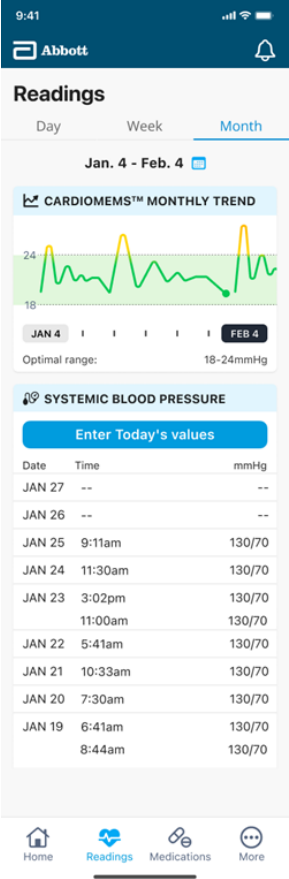


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI92	<Picture> READINGS SCREEN (no readings taken today, month view, PA Visibility enabled)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI27	8.7.3 Readings Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI94	8.7.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI95	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI96	8.7.3.2 GLOBAL HEADING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI97	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI98	8.7.3.3 READINGS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI99	Component: READINGS LABEL COMPONENT State: N/A Display: Always displayed as a component of the READINGS SCREEN. Label: <label.gb.readings_header> Behavior: Text.	PString427 Readings readings label.gb.readings_header	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1610	<Picture> READINGS LABEL COMPONENT Readings		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI146	8.7.3.4 INTERVAL TAB COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI147	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI148	8.7.3.5 DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI149	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI150	8.7.3.6 DATE SELECTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI151	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI152	8.7.3.7 HISTORICAL DATE SELECTION CALENDAR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI153	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI100	8.7.3.8 CMEMS PAP READING CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI101	See Shared Elements: Panels.		Info	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI4301	8.7.3.9 CMEMS PAP TREND CARD PANEL		Heading	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4302	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI102	8.7.3.10 BP READING CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI103	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4376	8.7.3.11 BP ERROR PANEL		Heading	CloudDRxRTR
PAGUI4377	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI110	8.7.3.12 GLOBAL NAVIGATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI111	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI28	8.8 Medications Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI29	8.8.1 Medications Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI353	The MEDICATIONS SCREEN is the application view used to view medications. The view is optimized to show medication dosage and changes in medication.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI354	The MEDICATIONS SCREEN consists of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) MEDICATIONS LABEL COMPONENT MEDICATIONS PANEL GLOBAL NAVIGATION PANEL (See Shared Elements: Panels) MEDICATIONS ENABLED PANEL MEDICATIONS DISABLED PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI30	8.8.2 Medications Screen (Layout)		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI356	<Picture> MEDICATIONS SCREEN (example: before scrolling) 9:41 Abbott Medications Clorthalidone Trade name 30mg 1 Tablet Once daily Furosemide Trade name 30mg 1 Tablet Once daily Telmisartan Zebeta 30mg 2 Tablets Once daily Rivaroxaban Zebeta 30mg 4 Tablets Once daily Home Readings Medications More		UI Design	CloudDRx 1.0 CloudDRxRTR
PAGUI31	8.8.3 Medications Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI359	8.8.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI360	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI361	8.8.3.2 GLOBAL HEADING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI362	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI363	8.8.3.3 MEDICATIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI364	Component: MEDICATIONS LABEL COMPONENT State: N/A Display: Always displayed in the MEDICATIONS SCREEN. Label: <label.gb.medications> Behavior: Text.	PString152 Medications medications label.gb.medications	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1920	<Picture> MEDICATIONS LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI371	8.8.3.6 MEDICATIONS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI395	8.8.3.6.1 MEDICATIONS PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI398	The MEDICATIONS PANEL consists of: MEDICATION INFORMATION SCROLLABLE ACTION TRIGGER COMPONENT SCROLL DOWN ICON COMPONENT (Shared Elements: Components) MEDICATION INFORMATION LABEL COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI396	8.8.3.6.2 MEDICATIONS PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI399	<Picture> MEDICATIONS PANEL (example: top of panel)	PString728 List fullList	UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Clorthalidone Trade name 30mg 1 Tablet Once daily Furosemide Trade name 30mg 1 Tablet Once daily Telmisartan Zebeta 30mg 2 Tablets Once daily Rivaroxaban Zebeta 30mg 4 Tablets Once daily	label.gb.list		
PAGUI4889	<Picture> MEDICATIONS PANEL (example: bottom of panel) Clorthalidone Trade name 30mg 1 Tablet Once daily Furosemide Trade name 30mg 1 Tablet Once daily Telmisartan Zebeta 30mg 2 Tablets Once daily Rivaroxaban Zebeta 30mg 4 Tablets Once daily Med Name Trade name 3mg 2 Tablets Once daily This may not be your full medication list. Do not stop any medications without your clinician's direction.		UI Design	CloudDRxRTR
PAGUI397	8.8.3.6.3 MEDICATIONS PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI400	8.8.3.6.3.1 MEDICATION INFORMATION SCROLLABLE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

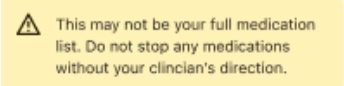
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI401	Component: MEDICATION INFORMATION SCROLLABLE ACTION TRIGGER COMPONENT State: Medication Management Enabled / Medication Management Disabled Display: Display the medication name, trade name, dosage amount, medication format, and number of times the medication is taken a day. Header: N/A Label: <data.Cod.Med_Frequency_Cd.Code_Sesc.QD>, <data.Code.Med_Frequency_Cd.Code_Desc.BID>, <data.Code.Med_Frequency_Cd.Code_Desc.TID>, <data.Code.Med_Frequency_Cd.Code_Desc.QID>, <data.Code.Med_Frequency_Cd.Code_Desc.QHS>, <data.Code.Med_Frequency_Cd.Code_Desc.QAM>, <data.Code.Med_Frequency_Cd.Code_Desc.QPM>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.TABLET>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.CAPSULE>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.IV>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.INJECTION>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.PATCH>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.OTHER>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.SPRAY> Medication names, formats, frequency, and dosage amounts are data collected from Merlin.net and simply displayed in the Patient Application. Behavior: Text and scroll. The MEDICATION INFORMATION SCROLLABLE ACTION TRIGGER COMPONENT may be scrolled vertically if the user's medication list requires more space than what is visible on the screen.	PString140 Once daily qd data.Code.Med_Frequency_Cd.Code_Desc.QD PString141 Twice daily bid data.Code.Med_Frequency_Cd.Code_Desc.BID PString142 3 times daily tid data.Code.Med_Frequency_Cd.Code_Desc.TID PString143 4 times daily qid data.Code.Med_Frequency_Cd.Code_Desc.QID PString146 Every night at bedtime qhs data.Code.Med_Frequency_Cd.Code_Desc.QHS PString147 Daily before noon qam data.Code.Med_Frequency_Cd.Code_Desc.QAM PString148 Daily after noon qpm data.Code.Med_Frequency_Cd.Code_Desc.QPM PString154 Tablet tablet data.Code.Dosage_Form_Unit_Cd.Code_Desc.TABLET	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

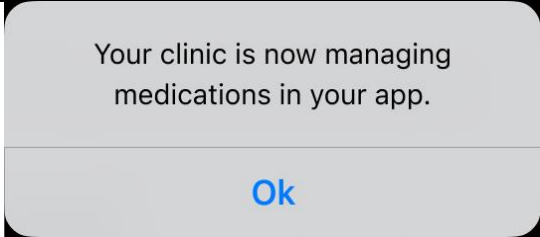
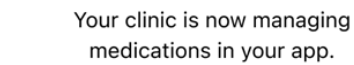
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		PString155 Capsule capsule data.Code.Dosage_Form_U nit_Cd.Code_Desc.CAPSULE PString156 IV iv data.Code.Dosage_Form_U nit_Cd.Code_Desc.IV PString157 Injection injection data.Code.Dosage_Form_U nit_Cd.Code_Desc.INJECTIO N PString158 Patch patch data.Code.Dosage_Form_U nit_Cd.Code_Desc.PATCH PString159 Other other data.Code.Dosage_Form_U nit_Cd.Code_Desc.OTHER PString161 Spray spray data.Code.Dosage_Form_U nit_Cd.Code_Desc.SPRAY		
PAGUI1927	<Picture> MEDICATIONS INFORMATION SCROLLABLE ACTION TRIGGER COMPONENT (default: standard) Rivaroxaban Zebeta 30mg 4 Tablets Once daily		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4860	8.8.3.6.3.2 SCROLL DOWN ICON COMPONENT (confirm cross functional alignment)		Heading	CloudDRxRTR
PAGUI4861	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI4378	8.8.3.6.3.3 MEDICATION INFORMATION LABEL COMPONENT		Heading	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

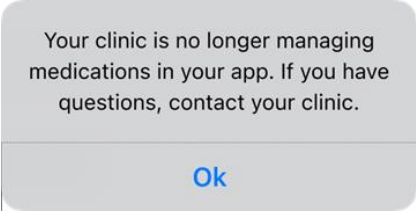
ID		English US Strings	Object Type	Project
PAGUI4379	<Picture> MEDICATION INFORMATION LABEL COMPONENT (example) 	PString1054 Your tailored list of medications to support your heart and blood pressure management label.gb.full_list_header PString1093 This may not be your full medication list. Do not stop any medications without your clinician's direction. label.gb.medList_warningBox	UI Design	CloudDRxRTR
PAGUI373	8.8.3.8 GLOBAL NAVIGATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI374	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1962	8.8.3.9 MEDICATIONS ENABLED PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1963	8.8.3.9.1 MEDICATIONS ENABLED PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1966	The MEDICATIONS ENABLED PANEL is displayed when the user selects MEDICATIONS SCREEN for the first time after the clinician enables medication tracking in the backend. The MEDICATIONS ENABLED PANEL is displayed over the entire screen and any part of the screen that is not covered is grayed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1967	The MEDICATIONS ENABLED PANEL consists of: MEDICATIONS ENABLED LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1964	8.8.3.9.2 MEDICATIONS ENABLED PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1968	<Picture> MEDICATIONS ENABLED PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1965	8.8.3.9.3 MEDICATIONS ENABLED PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1969	8.8.3.9.3.1 MEDICATIONS ENABLED LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1970	Component: MEDICATIONS ENABLED LABEL COMPONENT State: N/A Display: Always displayed in the MEDICATIONS ENABLED PANEL. Label: <msg.gb.med_manage_on> Behavior: Text.	PString512 Medication information has been enabled. If you have questions, contact your clinic. medicationsManageMedOn msg.gb.med_manage_on PString1020 Your clinic is now managing medications in your app. medicationsManageMedOn msg.gb.med_manage_on	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1971	<Picture> MEDICATIONS ENABLED LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1972	8.8.3.9.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1973	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1974	8.8.3.10 MEDICATIONS DISABLED PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1975	8.8.3.10.1 MEDICATIONS DISABLED PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

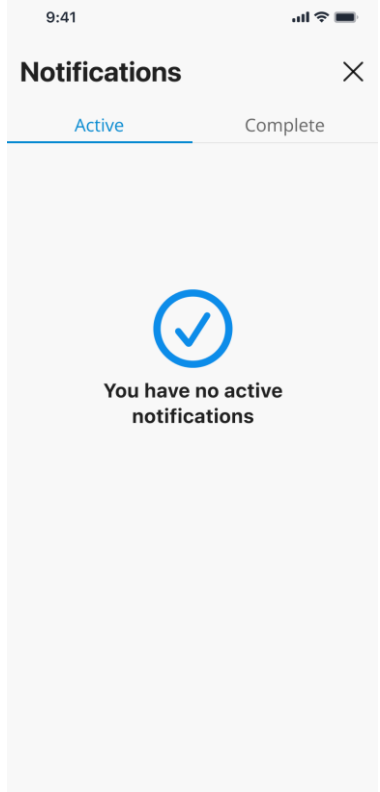
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1976	The MEDICATIONS DISABLED PANEL is displayed when the user selects MEDICATIONS SCREEN for the first time aftet the clinician disables medication tracking in the backend. The MEDICATIONS DISABLED PANEL is displayed over the entire screen and any part of the screen that is not covered is grayed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1977	The MEDICATIONS DISABLED PANEL consists of: MEDICATIONS DISABLED LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1978	8.8.3.10.2 MEDICATIONS DISABLED PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1979	<Picture> MEDICATIONS DISABLED PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1980	8.8.3.10.3 MEDICATIONS DISABLED PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1981	8.8.3.10.3.1 MEDICATIONS DISABLED LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1982	Component: MEDICATIONS DISABLED LABEL COMPONENT State: N/A Display: Always displayed in the MEDICATIONS DISABLED PANEL. Label: <msg.gb.med_manage_off> Behavior: Text.	PString324 Medication information is currently not available. Please contact your clinic if you have questions. medicationsManageMedOff msg.gb.med_manage_off PString1019 Your clinic is no longer managing medications in your app. If you have questions, contact your clinic. medicationsManageMedOff msg.gb.med_manage_off	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1983	<Picture> MEDICATIONS DISABLED LABEL COMPONENT		UI Design	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Your clinic is no longer managing medications in your app. If you have questions, contact your clinic.			CloudDRxPS
PAGUI1984	8.8.3.10.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1985	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI36	8.10 Notifications Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI37	8.10.1 Notifications Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI583	The NOTIFICATIONS SCREEN is the application view used to view any notifications sent to the user. The view is optimized to show unread notifications in priority order and historical notifications.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI584	The NOTIFICATIONS SCREEN consists of: OS STATUS PANEL (See Shared Elements: Panels) NOTIFICATIONS HEADER COMPONENT CLOSE PANEL ACTION TRIGGER COMPONENT NOTIFICATIONS VIEW TAB COMPONENT ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI38	8.10.2 Notifications Screen (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI585	<Picture> NOTIFICATIONS SCREEN (default: active, example: no messages)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI735	<Picture> NOTIFICATIONS SCREEN (example: complete, example: completed messages)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	9:41 Notifications × Active Complete ✓ MEDICATION ADJUSTMENT FRIDAY, MAR 10 Effective: 2022-06-24 Increase Clorthalidone Clorthalidone Trade name 1mg Tablet Once daily 2mg 1 Tablet Once daily ✓ CONTACT CLINIC FRIDAY, MAR 10 Contact clinic You had 3 high pressure readings in a row. Please contact your clinician for information about your treatment.			
PAGUI39	8.10.3 Notifications Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI586	8.10.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI587	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI590	8.10.3.3 NOTIFICATIONS HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

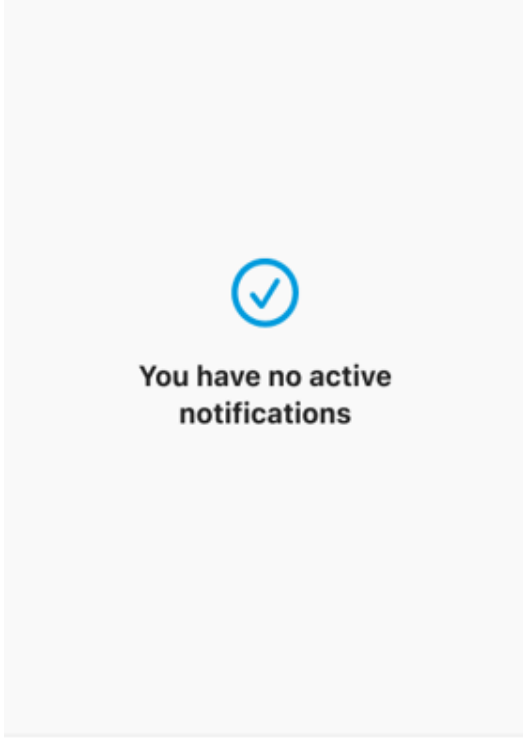
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI591	Component: NOTIFICATIONS HEADER COMPONENT State: N/A Display: Always displayed in the NOTIFICATIONS SCREEN Label: <label.gb.notifications> Behavior: Text.	PString191 Notifications notifications label.gb.notifications	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI592	<Picture> NOTIFICATIONS HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4492	8.10.3.4 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4491	See Shared Elements:Components		Info	CloudDRxRTR
PAGUI593	8.10.3.5 NOTIFICATIONS VIEW TAB COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI594	Component: NOTIFICATIONS VIEW TAB COMPONENT State: Active/ Complete Display: Always displayed in the NOTIFICATIONS SCREEN. The default state is active. Label: <label.gb.med_adj_complete>, <label.gb.adj_active> Behavior: Text and trigger. In the active state, display the ACTIVE NOTIFICATIONS PANEL. In the complete state, display the COMPLETE NOTIFICATIONS PANEL. Upon activation, the NOTIFICATIONS VIEW TAB COMPONENT will display different panels according to the state selected by the user. The tab corresponding to the active state is a different color from the other states.	PString115 Complete complete label.gb.med_adj_complete PString550 Active active label.gb.adj_active	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI595	<Picture> NOTIFICATIONS VIEW TAB COMPONENT (default: active) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1986	<Picture> NOTIFICATIONS VIEW TAB COMPONENT (example: complete) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1987	8.10.3.6 ACTIVE NOTIFICATIONS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1988	8.10.3.6.1 ACTIVE NOTIFICATIONS PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1991	The ACTIVE NOTIFICATIONS PANEL is displayed in the NOTIFICATIONS PANEL as default and after the user elects the active tab in the NOTIFICATIONS VIEW TAB COMPONENT. The ACTIVE NOTIFICATIONS PANEL is displayed under the NOTIFICATIONS VIEW TAB COMPONENT and allows the user to scroll down if the notifications can't all fit in the screen at one time. Notifications are displayed in priority order; if two notifications are the same priority level, they are displayed in chronological order.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3969	The ACTIVE NOTIFICATIONS PANEL consists of: NO NOTIFICATIONS LABEL COMPONENT		UI Info	CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	High Priority Notifications ELEVATION CHANGE MESSAGE CARD PANEL (See Shared Elements: Panels) DRX MEDICATION MESSAGE CARD PANEL (See Shared Elements: Panels) MEDICATION MESSAGE CARD PANEL (See Shared Elements: Panels) LAB REQUEST MESSAGE CARD PANEL (See Shared Elements: Panels) ELEVATED PRESSURE MESSAGE CARD PANEL (See Shared Elements: Panels) CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (See Shared Elements: Panels) DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (See Shared Elemtns: Panels) FREEFORM MESSAGE CARD PANEL (See Shared Elements: Panels) DRX FREEFORM MESSAGE CARD PANEL (See Shared Elements: Panels) Low Priority Notifications READINGS REVIEWED MESSAGE CARD PANEL (See Shared Elements: Panels) READING REMINDER MESSAGE CARD PANEL (See Shared Elements: Panels)			
PAGUI1989	8.10.3.6.2 ACTIVE NOTIFICATIONS PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1993	<Picture> ACTIVE NOTIFICATIONS PANEL (example: no notifications) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



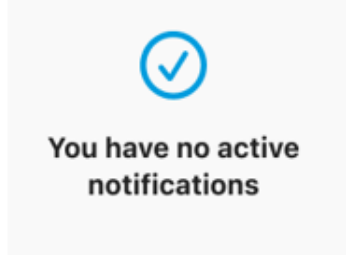
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1994	<Picture> ACTIVE NOTIFICATIONS PANEL (example: active notifications)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

	MEDICATION ADJUSTMENT FRIDAY, MAR 10 You had 3 high CardioMEMS™ readings in a row. Your clinician recommends an adjustment. 19 21 23 23 25 24 25 24 18 JAN 12 JAN 18 Optimal range: 18-26mmHg See Adjustment CONTACT CLINIC FRIDAY, MAR 10 You had 3 high pressure readings in a row. Please contact your clinician for information about your treatment. Call Clinic Mark as Read			
PAGUI1990	8.10.3.6.3 ACTIVE NOTIFICATIONS PANEL (Elements)		Heading	CloudDRx 1.0

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1995	8.10.3.6.3.1 NO NOTIFICATIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1996	Component: NO NOTIFICATIONS LABEL COMPONENT State: N/A Display: Conditionally displayed in the ACTIVE NOTIFICATIONS PANEL. The NO NOTIFICATIONS LABEL COMPONENT is displayed when there are no active notifications. Label: Behavior: Icon and text.	PString551 You have no active notifications noUnreadMessages msg.no_notifications	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1997	<Picture> NO NOTIFICATIONS LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3970	8.10.3.6.3.2 ELEVATION CHANGE MESSAGE CARD PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI3971	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS
PAGUI1999	8.10.3.6.3.3 DRX MEDICATION MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2000	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2001	8.10.3.6.3.4 MEDICATION MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2002	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2003	8.10.3.6.3.5 LAB REQUEST MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2004	See Shared Elements: Panels.		Info	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

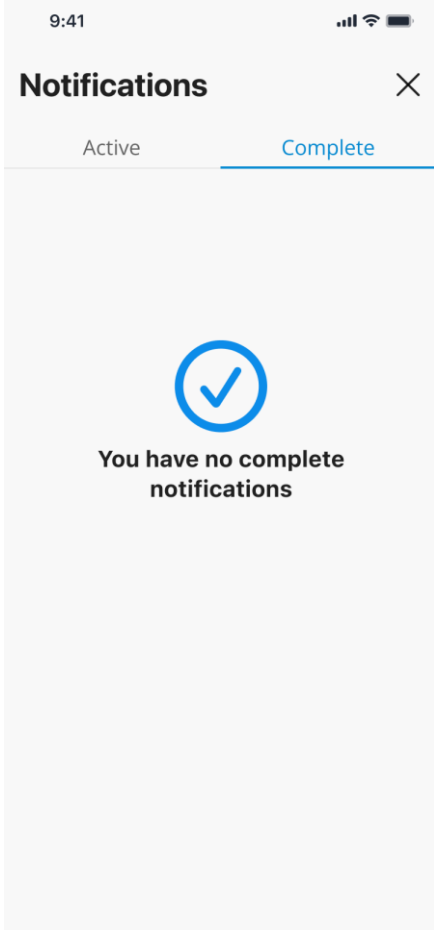
ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI2005	8.10.3.6.3.6 ELEVATED PRESSURE MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2006	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2007	8.10.3.6.3.7 CALL CLINIC DIRECTCALL MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2008	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2009	8.10.3.6.3.8 DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2010	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2011	8.10.3.6.3.9 FREEFORM MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2012	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2013	8.10.3.6.3.10 DRX FREEFORM MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2014	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2015	8.10.3.6.3.11 READINGS REVIEWED MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2016	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2017	8.10.3.6.3.12 READING REMINDER MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI2018	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2021	8.10.3.7 COMPLETE NOTIFICATIONS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2022	8.10.3.7.1 COMPLETE NOTIFICATIONS PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2023	The COMPLETE NOTIFICATIONS PANEL is conditionally displayed in the NOTIFICATIONS PANEL after the user selects the complete tab in the NOTIFICATIONS VIEW TAB COMPONENT. The COMPLETE NOTIFICATIONS PANEL is displayed under the NOTIFICATIONS VIEW TAB COMPONENT and allows the user to scroll down if the notifications can't all fit in the screen at one time. Notifications are displayed in chronological order.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3977	The COMPLETE NOTIFICATIONS PANEL consists of: NO NOTIFICATIONS LABEL COMPONENT ELEVATION CHANGE MESSAGE CARD PANEL (See Shared Elements: Panels) DRX MEDICATION MESSAGE CARD PANEL (See Shared Elements: Panels) MEDICATION MESSAGE CARD PANEL (See Shared Elements: Panels) LAB REQUEST MESSAGE CARD PANEL (See Shared Elements: Panels) ELEVATED PRESSURE MESSAGE CARD PANEL (See Shared Elements: Panels) CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (See Shared Elements: Panels) DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (See Shared Elements: Panels) FREEFORM MESSAGE CARD PANEL (See Shared Elements: Panels) DRX FREEFORM MESSAGE CARD PANEL (See Shared Elements: Panels) READINGS REVIEWED MESSAGE CARD PANEL (See Shared Elements: Panels) READING REMINDER MESSAGE CARD PANEL (See Shared Elements: Panels)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI2025	8.10.3.7.2 COMPLETE NOTIFICATIONS PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2026	<Picture> COMPLETE NOTIFICATIONS PANEL (example: no notifications)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI2027	<Picture> COMPLETE NOTIFICATIONS PANEL (example: complete notifications)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	✓ MEDICATION ADJUSTMENT FRIDAY, MAR 10 Effective: 2022-06-24 Increase Clorthalidone Clorthalidone Trade name 1mg Tablet Once daily 2mg 1 Tablet Once daily ✓ CONTACT CLINIC FRIDAY, MAR 10 Contact clinic You had 3 high pressure readings in a row. Please contact your clinician for information about your treatment.			
PAGUI2028	8.10.3.7.3 COMPLETE NOTIFICATIONS PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3707	8.10.3.7.3.1 NO NOTIFICATIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3708	Component: NO NOTIFICATIONS LABEL COMPONENT State: N/A Display: Conditionally displayed in the COMPLETE NOTIFICATIONS PANEL. The NO NOTIFICATIONS LABEL COMPONENT is displayed when there are no unread notifications available after the Earliest History Calendar Date. Label: Behavior: Icon and text.	PString708 You have no read messages. noReadMessages msg.gb.no_read_notifications PString1021 You have no complete notifications	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		noReadMessages msg.gb.no_read_notifications		
PAGUI3709	<Picture> NO NOTIFICATIONS LABEL COMPONENT You have no complete notifications		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3975	8.10.3.7.3.2 ELEVATION CHANGE MESSAGE CARD PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI3976	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS
PAGUI2032	8.10.3.7.3.3 DRX MEDICATION MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2033	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2034	8.10.3.7.3.4 MEDICATION MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2035	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2036	8.10.3.7.3.5 LAB REQUEST MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2037	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2038	8.10.3.7.3.6 ELEVATED PRESSURE MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2039	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2040	8.10.3.7.3.7 CALL CLINIC DIRECTCALL MESSAGE CARD PANEL		Heading	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

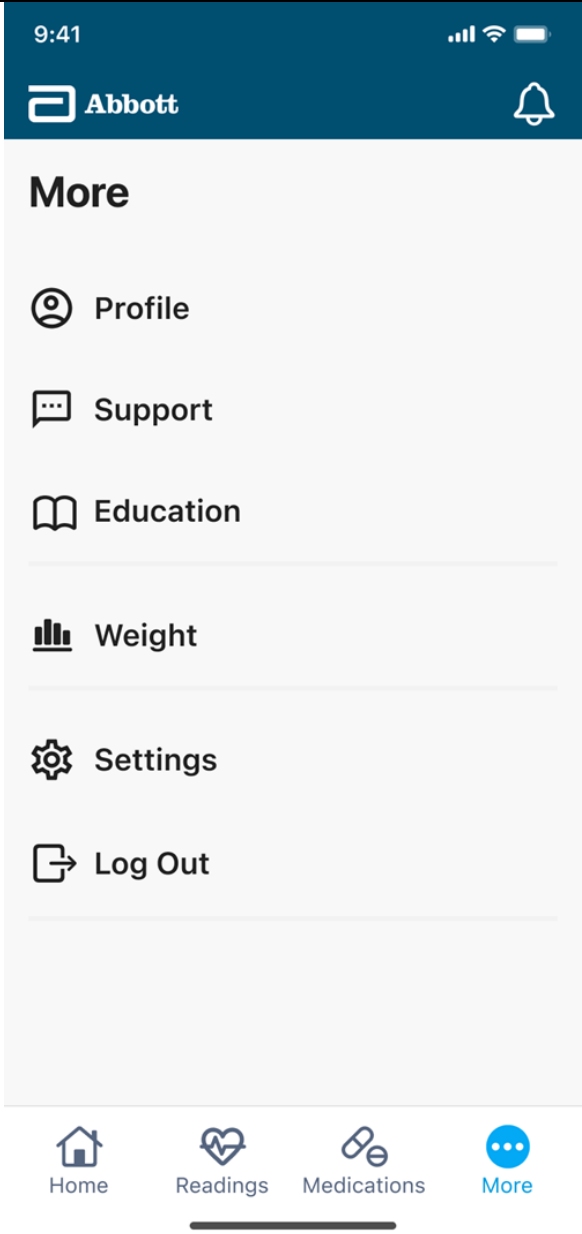
ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI2041	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2042	8.10.3.7.3.8 DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2043	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2044	8.10.3.7.3.9 FREEFORM MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2045	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2046	8.10.3.7.3.10 DRX FREEFORM MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2047	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2048	8.10.3.7.3.11 READINGS REVIEWED MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2049	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2050	8.10.3.7.3.12 READING REMINDER MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2051	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI551	8.11 More Screen		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI552	8.11.1 More Screen (Display)		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI555	The MORE SCREEN is displayed over the whole screen except for the OS STATUS PANEL and GLOBAL HEADING PANEL and is conditionally displayed as a component of the GLOBAL NAVIGATION PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI556	The MORE SCREEN consists of: MORE HEADER COMPONENT PROFILE ACTION TRIGGER COMPONENT SETTINGS ACTION TRIGGER COMPONENT EDUCATION ACTION TRIGGER COMPONENT SUPPORT ACTION TRIGGER COMPONENT WEIGHT ACTION TRIGGER COMPONENT LOG OUT ACTION TRIGGER COMPONENT LOG OUT PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI553	8.11.2 More Screen (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI557	<Picture> MORE SCREEN		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI554	8.11.3 More Screen (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI559	8.11.3.1 MORE HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI560	Component: MORE HEADER COMPONENT State: N/A Display: Always displayed in the MORE SCREEN. Label: <label.gb.menu> Behavior: Text.	PString582 Menu menu label.gb.menu PString1049 More more_label label.gb.more	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI561	<Picture> MORE HEADER COMPONENT More		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI564	8.11.3.3 PROFILE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI565	Component: PROFILE ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the MORE SCREEN. Label: <label.gb.profile> Behavior: Icon, text, and trigger. Upon activation, the PROFILE ACTION TRIGGER COMPONENT displays the PROFILE PANEL.	PString205 Profile profile label.gb.profile	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI566	<Picture> PROFILE ACTION TRIGGER COMPONENT  Profile		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI573	8.11.3.3.1 SUPPORT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI574	Component: SUPPORT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the MORE SCREEN. Label: <label.gb.support> Behavior: Icon, text, and trigger. Upon activation, the SUPPORT ACTION TRIGGER COMPONENT displays the SUPPORT PANEL.	PString585 Support support label.gb.support	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI575	<Picture> SUPPORT ACTION TRIGGER COMPONENT  Support		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI570	8.11.3.4 EDUCATION ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

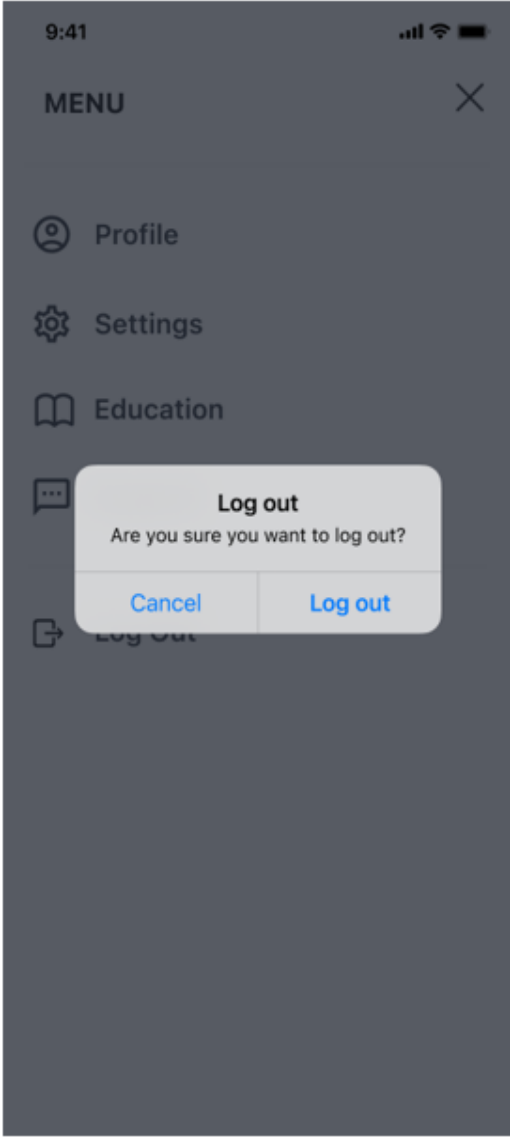
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI571	Component: EDUCATION ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the MORE SCREEN. Label: <label.gb.education> Behavior: Icon, text, and trigger. Upon activation, the EDUCATION ACTION TRIGGER COMPONENT displays the EDUCATION PANEL.	PString584 Education education label.gb.education	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI572	<Picture> EDUCATION ACTION TRIGGER COMPONENT  Education		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4330	8.11.3.5 WEIGHT ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4331	Component: WEIGHT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the MORE SCREEN. Label: <label.gb.education> Behavior: Icon, text, and trigger. Upon activation, the WEIGHT ACTION TRIGGER COMPONENT displays the WEIGHT PANEL.		UI Info	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4332	<Picture> WEIGHT ACTION TRIGGER COMPONENT  Weight		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI567	8.11.3.6 SETTINGS ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI568	Component: SETTINGS ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the MORE SCREEN. Label: <label.gb.settings> Behavior: Icon, text, and trigger. Upon activation, the SETTINGS ACTION TRIGGER COMPONENT displays the SETTINGS PANEL.	PString497 Settings settingsHome_title label.gb.settings	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI569	<Picture> SETTINGS ACTION TRIGGER COMPONENT  Settings		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI576	8.11.3.7 LOG OUT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI577	Component: LOG OUT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the MORE SCREEN. Label: <label.gb.logout>	PString187 Log out logout label.gb.logout	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

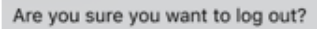
ID		English US Strings	Object Type	Project
	Behavior: Icon, text, and trigger. Upon activation, the LOG OUT ACTION TRIGGER COMPONENT displays the LOG OUT PANEL.			
PAGUI578	<Picture> LOG OUT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3029	8.11.3.8 LOG OUT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3030	8.11.3.8.1 LOG OUT PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3033	The LOG OUT PANEL is conditionally displayed as an element of the MORE SCREEN after the user selects the LOG OUT ACTION TRIGGER COMPONENT.The LOG OUT PANEL is displayed over the MORE SCREEN and the background is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3034	The LOG OUT PANEL consists of: LOG OUT HEADER COMPONENT LOG OUT LABEL COMPONENT CANCEL LOG OUT ACTION TRIGGER COMPONENT CONFIRM LOG OUT ACTION TRIGGER COMPONENT		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3031	8.11.3.8.2 LOG OUT PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3035	<Picture> LOG OUT PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

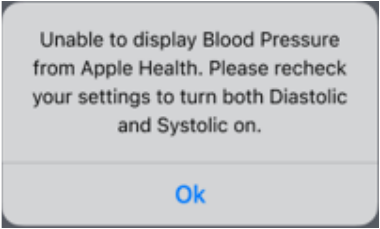
ID		English US Strings	Object Type	Project
				
PAGUI3032	8.11.3.8.3 LOG OUT PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3036	8.11.3.8.3.1 LOG OUT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3037	Component: LOG OUT HEADER COMPONENT State: N/A Display: Always displayed in the LOG OUT PANEL. Label: <label.gb.logout> Behavior: Text.	PString187 Log out logout label.gb.logout	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3038	<Picture> LOG OUT HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3039	8.11.3.8.3.2 LOG OUT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3040	Component: LOG OUT LABEL COMPONENT State: N/A Display: Always displayed in the LOG OUT PANEL. Label: <msg.gb.confirm_logout> Behavior: Text.	PString222 Are you sure you want to log out? logoutBody msg.gb.confirm_logout	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3041	<Picture> LOG OUT LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3042	8.11.3.8.3.3 CANCEL LOG OUT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3043	Component: CANCEL LOG OUT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the LOG OUT PANEL. Label: <label.gb.cancel> Behavior: Icon, text, and trigger. Upon activation, the CANCEL LOG OUT ACTION TRIGGER COMPONENT displays the MORE SCREEN.	PString104 Cancel cancel label.gb.cancel	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3044	<Picture> CANCEL LOG OUT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3045	8.11.3.8.3.4 CONFIRM LOG OUT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3046	Component: CONFIRM LOG OUT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the LOG OUT PANEL. Label: <label.gb.logout> Behavior: Icon, text, and trigger. Upon activation, the CONFIRM LOG OUT ACTION TRIGGER COMPONENT logs the user out of the Patient Application.	PString187 Log out logout label.gb.logout	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3047	<Picture> CONFIRM LOG OUT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI40	8.12 Shared Elements		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI41	8.12.1 Panels		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1619	8.12.1.3 BP ERROR PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1620	8.12.1.3.1 BP ERROR PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1623	The BP ERROR PANEL is conditionally displayed as a component of READINGS PANEL and SETTINGS PANEL if only Diastolic Blood Pressure or Systolic Blood Pressure is shared from the connected health app in Apple iOS 17 and 18. The BP ERROR PANEL is displayed over the screen, does not displace any elements, and any part of the screen that is not covered by the BP ERROR PANEL is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1625	The BP ERROR PANEL is a component of: READINGS SCREEN SETTINGS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1624	The BP ERROR PANEL is composed of: BP ERROR LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1621	8.12.1.3.2 BP ERROR PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1626	<Picture> BP ERROR PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1622	8.12.1.3.3 BP ERROR PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1627	8.12.1.3.3.1 BP ERROR LABEL COMPONENT		Heading	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1628	Component: BP ERROR LABEL COMPONENT State: N/A Display: Always displayed in the BP ERROR PANEL in Apple iOS 17 and 18. Label: <msg.gb.bp_permissions> Behavior: Text.	PString554 Unable to display Blood Pressure from {19}. Please recheck your settings to turn both Diastolic and Systolic on. missingSystolicAndDiastolicPermission msg.gb.bp_permissions	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1629	<Picture> BP ERROR LABEL COMPONENT Unable to display Blood Pressure from Apple Health. Please recheck your settings to turn both Diastolic and Systolic on.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1630	8.12.1.3.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1631	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI47	8.12.1.4 BP READING CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI112	8.12.1.4.1 BP READING CARD PANEL (DISPLAY)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI118	The BP READING CARD PANEL is displayed in the READINGS SCREEN and HOME SCREEN. If the user has not taken a BP Reading for the day, display the BP READING HEADER COMPONENT in enabled state and provide a reminder to take daily reading as described in BP READING REMINDER LABEL COMPONENT. Also, provide RECORD BP READING ACTION TRIGGER COMPONENT so user has a way of inputting BP data. If the user has taken a BP Reading for the day, display the BP READING HEADER COMPONENT in the disabled state and show BP READING DISPLAY COMPONENT.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI326	The BP READING CARD PANEL is an element of: HOME SCREEN READINGS SCREEN		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI122	The BP READING CARD PANEL consists of: BP READING HEADER COMPONENT BP READING REMINDER LABEL COMPONENT RECORD BP READING ACTION TRIGGER COMPONENT BP READING DISPLAY COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

Commented [RC1]: Example: BP Reading Card Panel

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project																											
PAGUI113	8.12.1.4.2 BP READING CARD PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS																											
PAGUI2989	<Picture> BP READING CARD (example: Day, default: no reading for the day)  SYSTEMIC BLOOD PRESSURE Blood pressure not recorded today. Record Blood Pressure		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS																											
PAGUI2990	<Picture> BP READING CARD (default: week, example: reading complete, Screen: Readings)  SYSTEMIC BLOOD PRESSURE <table><tr><th>Date</th><th>Time</th><th>mmHg</th></tr><tr><td>Today</td><td>9:41am</td><td>130/70</td></tr><tr><td>Jan 29</td><td>7:45am</td><td>115/77</td></tr><tr><td></td><td>7:40am</td><td>110/90</td></tr><tr><td>Jan 28</td><td>8:02am</td><td>110/90</td></tr><tr><td>Jan 27</td><td>8:14am</td><td>121/87</td></tr><tr><td>Jan 26</td><td>7:56am</td><td>116/85</td></tr><tr><td>Jan 25</td><td>7:40am</td><td>115/68</td></tr><tr><td>Jan 24</td><td>8:05am</td><td>119/85</td></tr></table>	Date	Time	mmHg	Today	9:41am	130/70	Jan 29	7:45am	115/77		7:40am	110/90	Jan 28	8:02am	110/90	Jan 27	8:14am	121/87	Jan 26	7:56am	116/85	Jan 25	7:40am	115/68	Jan 24	8:05am	119/85		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
Date	Time	mmHg																													
Today	9:41am	130/70																													
Jan 29	7:45am	115/77																													
	7:40am	110/90																													
Jan 28	8:02am	110/90																													
Jan 27	8:14am	121/87																													
Jan 26	7:56am	116/85																													
Jan 25	7:40am	115/68																													
Jan 24	8:05am	119/85																													
PAGUI2991	<Picture> BP READING CARD (example: Month, default: no reading for the day, Screen: Readings)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS																											

myCardioMEMS Patient Application Graphical User Interface Specification

SYSTEMIC BLOOD PRESSURE

Enter today's values

Date	Time	mmHg
Today	--	--
Jan 29	7:45am	115/77
	7:40am	110/90
Jan 28	8:02am	110/90
Jan 27	8:14am	121/87
Jan 26	7:56am	116/85
Jan 25	7:40am	115/68
Jan 24	8:05am	119/85
Jan 23	7:45am	130/70
Jan 22	7:45am	115/77
Jan 21	8:02am	110/90
Jan 20	8:14am	121/87
Jan 19	7:56am	116/85
Jan 18	7:40am	115/68
Jan 17	8:05am	119/85
Jan 16	7:45am	130/70
Jan 15	7:45am	115/77
Jan 14	8:02am	110/90
Jan 13	8:14am	121/87
Jan 12	7:56am	116/85
Jan 11	7:40am	115/68
Jan 10	8:05am	119/85
Jan 09	7:45am	130/70
Jan 08	7:45am	115/77
Jan 07	8:02am	110/90
Jan 06	8:14am	121/87
Jan 05	7:56am	116/85
Jan 04	7:40am	115/68
Jan 03	8:05am	119/85

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

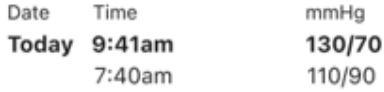


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI114	8.12.1.4.3 BP READING CARD PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI115	8.12.1.4.3.1 BP READING HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI120	Component: BP READING HEADER COMPONENT State: Enabled/Disabled Display: Always displayed in the BP READING CARD PANEL. The default state is enabled. If daily reading has been taken, display the header text with white background as no action is needed (disabled). If daily reading has not been taken, display header text with blue background (enabled). Label: <label.gb.bp> Behavior: Icon and text.	PString447 Systemic Blood Pressure bpReadingTitle label.gb.bp	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2992	<Picture> BP READING CARD (default: Enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2993	<Picture> BP READING CARD (example: Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI116	8.12.1.4.3.2 BP READING REMINDER LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI121	Component: BP READING REMINDER LABEL COMPONENT State: No Reading/ Reading Display: Conditionally displayed in the BP READING CARD PANEL in the no reading state on the HOME SCREEN. In READINGS SCREEN, the BP READING REMINDER LABEL COMPONENT is displayed in the Day tab. In reading state, the BP READING REMINDER LABEL COMPONENT is not present. Label: <label.gb.bp_reminder> Behavior: Text. The BP READING REMINDER LABEL COMPONENT changes from the no reading state to the reading state when a BP Reading is recorded in the Patient Application or in the Connected Health Application.	PString448 It's time to record today's blood pressure. bpReadingSubtitle label.gb.bp_reminder PString1012 Blood pressure not recorded today. bpReadingSubtitle label.gb.bp_reminder	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2994	<Picture> BP READING REMINDER LABEL COMPONENT (reading not taken) Blood pressure not recorded today.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2995	8.12.1.4.3.3 RECORD BP READING ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2996	Component: RECORD BP READING ACTION TRIGGER COMPONENT State: No Reading/ Reading	PString449 Record Blood Pressure	UI Info	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Display: Conditionally displayed in the BP READING CARD PANEL in the no reading state. In the READINGS SCREEN, the RECORD BP READING ACTION TRIGGER COMPONENT displays different text in the 'Day' view than it does in the 'Week' and 'Month' view. Label: <label.gb.bp_record>, <label.gb.bp_enter> Behavior: Text and trigger. Upon activation in the no reading state, the RECORD BP READING ACTION TRIGGER COMPONENT displays the RECORD BP READING PANEL. The RECORD BP READING ACTION TRIGGER COMPONENT changes from the no reading state to the reading state when a BP Reading is recorded in the Patient Application or in the Connected Health Application.	bpReadingButtonText label.gb.bp_record PString450 Enter Today's Values enterTodaysValues label.gb.bp_enter		CloudDRxPS
PAGUI2997	<Picture> RECORD BP READING ACTION TRIGGER COMPONENT (default: HOME SCREEN and READINGS SCREEN in the 'Day' view) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2998	<Picture> RECORD BP READING ACTION TRIGGER COMPONENT (example: READINGS SCREEN in the 'Week' and 'Month' view) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI117	8.12.1.4.3.4 BP READING DISPLAY COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI51	Component: BP READING DISPLAY COMPONENT State: N/A Display: Conditionally displayed in the BP READING CARD PANEL if a reading has been taken for the day or if the BP READING CARD PANEL is in the READINGS SCREEN under the 'Week' or 'Month' view. Display the number of available readings in table format per the interval specified in the READINGS INTERVAL SELECTION COMPONENT and READINGS DATE SELECTION COMPONENT. If no readings are available for a displayed date, display an indication that no values are available. Multiple readings can be sent in one day by using a Connected Health Application. Label: N/A - determined by a third party vendor package Behavior: Text.	PString94 Today globalToday label.gb.today PString224 Time time label.gb.time PString451 Date date label.gb.date PString452 mmHg mmHg label.gb.mmhg	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2999	<Picture> BP READING DISPLAY COMPONENT (default for HOME SCREEN: Day view, example: multiple readings via Connected Health Application) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project																											
PAGUI3000	<Picture> BP READING DISPLAY COMPONENT (default for READINGS SCREEN: Week view, default: no reading for today) <table><tr><th>Date</th><th>Time</th><th>mmHg</th></tr><tr><td>Today</td><td>--</td><td>--</td></tr><tr><td>Jan 29</td><td>7:45am</td><td>115/77</td></tr><tr><td></td><td>7:40am</td><td>110/90</td></tr><tr><td>Jan 28</td><td>8:02am</td><td>110/90</td></tr><tr><td>Jan 27</td><td>8:14am</td><td>121/87</td></tr><tr><td>Jan 26</td><td>7:56am</td><td>116/85</td></tr><tr><td>Jan 25</td><td>7:40am</td><td>115/68</td></tr><tr><td>Jan 24</td><td>8:05am</td><td>119/85</td></tr></table>	Date	Time	mmHg	Today	--	--	Jan 29	7:45am	115/77		7:40am	110/90	Jan 28	8:02am	110/90	Jan 27	8:14am	121/87	Jan 26	7:56am	116/85	Jan 25	7:40am	115/68	Jan 24	8:05am	119/85		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
Date	Time	mmHg																													
Today	--	--																													
Jan 29	7:45am	115/77																													
	7:40am	110/90																													
Jan 28	8:02am	110/90																													
Jan 27	8:14am	121/87																													
Jan 26	7:56am	116/85																													
Jan 25	7:40am	115/68																													
Jan 24	8:05am	119/85																													
PAGUI3001	<Picture> BP READING DISPLAY COMPONENT (default for READINGS SCREEN: Week view, example: reading for the day) <table><tr><th>Date</th><th>Time</th><th>mmHg</th></tr><tr><td>Today</td><td>9:41am</td><td>130/70</td></tr><tr><td>Jan 29</td><td>7:45am</td><td>115/77</td></tr><tr><td></td><td>7:40am</td><td>110/90</td></tr><tr><td>Jan 28</td><td>8:02am</td><td>110/90</td></tr><tr><td>Jan 27</td><td>8:14am</td><td>121/87</td></tr><tr><td>Jan 26</td><td>7:56am</td><td>116/85</td></tr><tr><td>Jan 25</td><td>7:40am</td><td>115/68</td></tr><tr><td>Jan 24</td><td>8:05am</td><td>119/85</td></tr></table>	Date	Time	mmHg	Today	9:41am	130/70	Jan 29	7:45am	115/77		7:40am	110/90	Jan 28	8:02am	110/90	Jan 27	8:14am	121/87	Jan 26	7:56am	116/85	Jan 25	7:40am	115/68	Jan 24	8:05am	119/85		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
Date	Time	mmHg																													
Today	9:41am	130/70																													
Jan 29	7:45am	115/77																													
	7:40am	110/90																													
Jan 28	8:02am	110/90																													
Jan 27	8:14am	121/87																													
Jan 26	7:56am	116/85																													
Jan 25	7:40am	115/68																													
Jan 24	8:05am	119/85																													
PAGUI3002	<Picture> BP READING DISPLAY COMPONENT (default for READINGS SCREEN: Month view, default: no reading today)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS																											



myCardioMEMS Patient Application Graphical User Interface Specification

Date	Time	mmHg			
Today	--	--			
Jan 29	7:45am	115/77			
	7:40am	110/90			
Jan 28	8:02am	110/90			
Jan 27	8:14am	121/87			
Jan 26	7:56am	116/85			
Jan 25	7:40am	115/68			
Jan 24	8:05am	119/85			
Jan 23	7:45am	130/70			
Jan 22	7:45am	115/77			
Jan 21	8:02am	110/90			
Jan 20	8:14am	121/87			
Jan 19	7:56am	116/85			
Jan 18	7:40am	115/68			
Jan 17	8:05am	119/85			
Jan 16	7:45am	130/70			
Jan 15	7:45am	115/77			
Jan 14	8:02am	110/90			
Jan 13	8:14am	121/87			
Jan 12	7:56am	116/85			
Jan 11	7:40am	115/68			
Jan 10	8:05am	119/85			
Jan 09	7:45am	130/70			
Jan 08	7:45am	115/77			
Jan 07	8:02am	110/90			
Jan 06	8:14am	121/87			
Jan 05	7:56am	116/85			
Jan 04	7:40am	115/68			
Jan 03	8:05am	119/85			

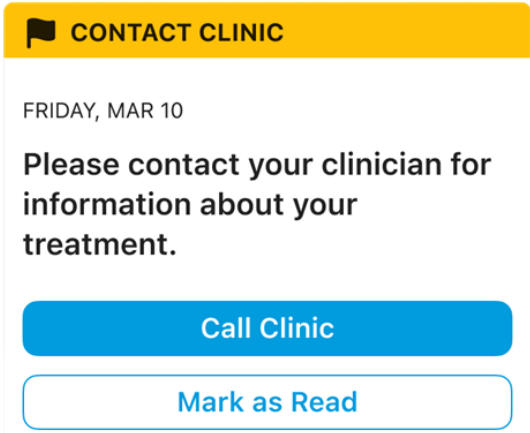
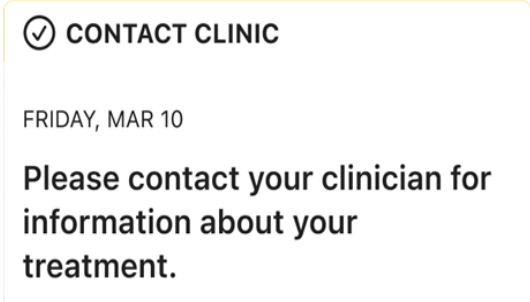
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project																																																									
PAGUI3003	<Picture> BP READING DISPLAY COMPONENT (default for READINGS SCREEN: Month view, example: reading for the day) <table><tr><th>Date</th><th>Time</th><th>mmHg</th></tr><tr><td>Today</td><td>9:41am</td><td>130/70</td></tr><tr><td>Jan 29</td><td>7:45am</td><td>115/77</td></tr><tr><td></td><td>7:40am</td><td>110/90</td></tr><tr><td>Jan 28</td><td>8:02am</td><td>110/90</td></tr><tr><td>Jan 27</td><td>8:14am</td><td>121/87</td></tr><tr><td>Jan 26</td><td>7:56am</td><td>116/85</td></tr><tr><td>Jan 25</td><td>7:40am</td><td>115/68</td></tr><tr><td>Jan 24</td><td>8:05am</td><td>119/85</td></tr><tr><td>Jan 23</td><td>7:45am</td><td>130/70</td></tr><tr><td>Jan 22</td><td>7:45am</td><td>115/77</td></tr><tr><td>Jan 21</td><td>8:02am</td><td>110/90</td></tr><tr><td>Jan 20</td><td>8:14am</td><td>121/87</td></tr><tr><td>Jan 19</td><td>7:56am</td><td>116/85</td></tr><tr><td>Jan 18</td><td>7:40am</td><td>115/68</td></tr><tr><td>Jan 17</td><td>8:05am</td><td>119/85</td></tr><tr><td>Jan 16</td><td>7:45am</td><td>130/70</td></tr><tr><td>Jan 15</td><td>7:45am</td><td>115/77</td></tr><tr><td>Jan 14</td><td>8:02am</td><td>110/90</td></tr></table>	Date	Time	mmHg	Today	9:41am	130/70	Jan 29	7:45am	115/77		7:40am	110/90	Jan 28	8:02am	110/90	Jan 27	8:14am	121/87	Jan 26	7:56am	116/85	Jan 25	7:40am	115/68	Jan 24	8:05am	119/85	Jan 23	7:45am	130/70	Jan 22	7:45am	115/77	Jan 21	8:02am	110/90	Jan 20	8:14am	121/87	Jan 19	7:56am	116/85	Jan 18	7:40am	115/68	Jan 17	8:05am	119/85	Jan 16	7:45am	130/70	Jan 15	7:45am	115/77	Jan 14	8:02am	110/90		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
Date	Time	mmHg																																																											
Today	9:41am	130/70																																																											
Jan 29	7:45am	115/77																																																											
	7:40am	110/90																																																											
Jan 28	8:02am	110/90																																																											
Jan 27	8:14am	121/87																																																											
Jan 26	7:56am	116/85																																																											
Jan 25	7:40am	115/68																																																											
Jan 24	8:05am	119/85																																																											
Jan 23	7:45am	130/70																																																											
Jan 22	7:45am	115/77																																																											
Jan 21	8:02am	110/90																																																											
Jan 20	8:14am	121/87																																																											
Jan 19	7:56am	116/85																																																											
Jan 18	7:40am	115/68																																																											
Jan 17	8:05am	119/85																																																											
Jan 16	7:45am	130/70																																																											
Jan 15	7:45am	115/77																																																											
Jan 14	8:02am	110/90																																																											
PAGUI2301	8.12.1.5 CALL CLINIC DIRECTCALL MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS																																																									
PAGUI2302	8.12.1.5.1 CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS																																																									
PAGUI2303	The CALL CLINIC DIRECTCALL MESSAGE CARD PANEL is conditionally displayed if a Call Clinic DirectCall Message is received by the Patient Application at least once since the Earliest History Calendar Date. A unique CALL CLINIC DIRECTCALL MESSAGE CARD PANEL is displayed for every Call Clinic DirectCall Message that is received by the Patient Application. Multiple CALL CLINIC DIRECTCALL MESSAGE CARD PANELS may be displayed in the same screen or panel. The CALL CLINIC DIRECTCALL MESSAGE CARD PANEL is in the active state by default. The CALL CLINIC DIRECTCALL MESSAGE CARD PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS																																																									

myCardioMEMS Patient Application Graphical User Interface Specification

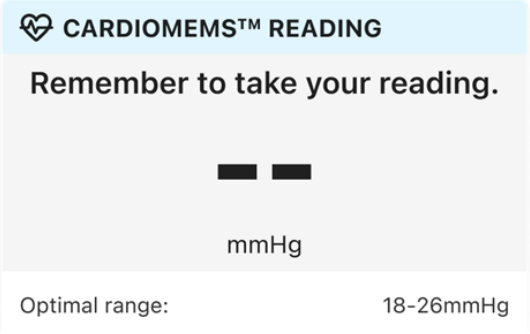
ID		English US Strings	Object Type	Project
	CARD PANEL changes state from active to complete when the user completes all required actions in the card. The elements in the CALL CLINIC DIRECTCALL MESSAGE CARD PANEL listed below change behavior along these two states.			
PAGUI2304	The CALL CLINIC DIRECTCALL MESSAGE CARD PANEL is an element of: ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2305	The CALL CLINIC DIRECTCALL MESSAGE CARD PANEL consists of: CONTACT CLINIC BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) CALL CLINIC CONTENT LABEL COMPONENT (See: Shared Elements: Components) CALL CLINIC ACTION TRIGGER COMPONENT (See Shared Elements: Components) MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2306	8.12.1.5.2 CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2307	<Picture> CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (default: active, example, DTP OFF) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2308	<Picture> CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (complete, example, DTP OFF) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2309	8.12.1.5.3 CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

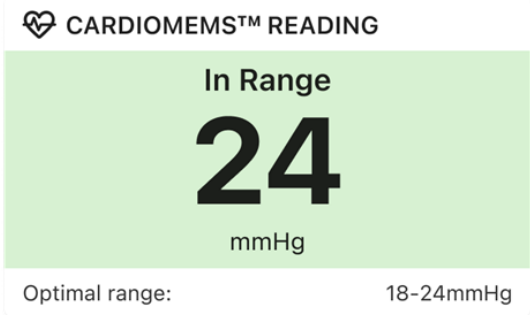
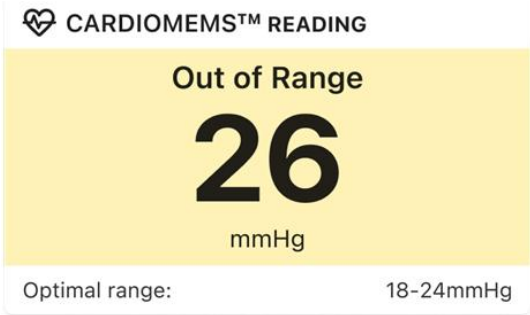
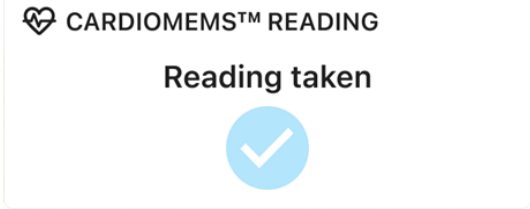
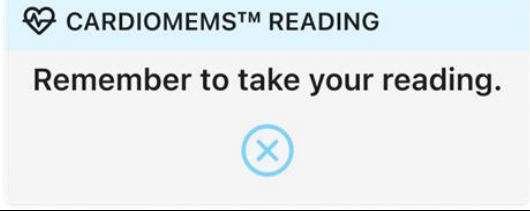
ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI2324	8.12.1.5.3.1 CONTACT CLINIC BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2325	Component: CONTACT CLINIC BANNER COMPONENT State: Active/ Complete Display: Always displayed in the CALL CLINIC DIRECTCALL MESSAGE CARD PANEL. The default state is active. In the active state, the CALL CLINIC BANNER COMPONENT displays the active icon and the active background color. In the complete state, the CALL CLINIC BANNER COMPONENT displays the complete icon and complete background color. Label: <label.gb.contact_clinic>; Text is all capitalized. Behavior: Icon and text. The CALL CLINIC BANNER COMPONENT changes state with the message card panel it is an element of.	PString604 Contact Clinic contactClinic label.gb.contact_clinic	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2326	<Picture> CALL CLINIC BANNER COMPONENT (default: Active) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2327	<Picture> CALL CLINIC BANNER COMPONENT (Complete) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2312	8.12.1.5.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2313	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2317	8.12.1.5.3.3 CALL CLINIC CONTENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2318	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2322	8.12.1.5.3.4 CALL CLINIC ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2323	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2320	8.12.1.5.3.5 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

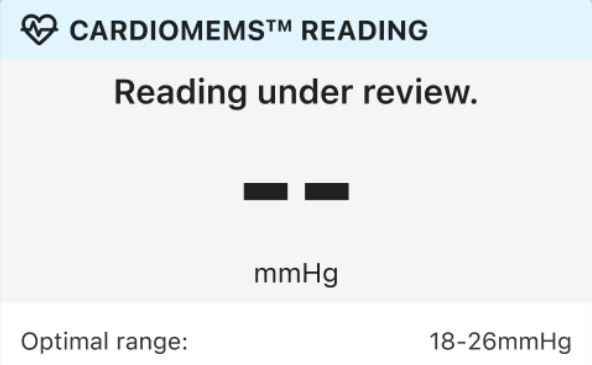
ID		English US Strings	Object Type	Project
PAGUI2321	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI46	8.12.1.6 CMEMS PAP READING CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI55	8.12.1.6.1 CMEMS PAP READING CARD PANEL (DISPLAY)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI57	The CMEMS PAP READING CARD PANEL displays PAP related information and reminders. The CMEMS PAP READING CARD PANEL displays different PAP information depending on if the patient has completed their PAP Reading for the day, and if PA Visibility is enabled. Upon activation, display the CMEMS PAP READING INFORMATION PANEL. If the USER has not taken a CMEMS PAP Reading for the day, enable panel activation. If the USER has taken a CMEMS PAP Reading for the day, disable panel activation (display only). The CMEMS PAP READING CARD PANEL is only displayed in the Day tab of the READINGS SCREEN.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI206	The CMEMS PAP READING CARD PANEL is an element of: HOME SCREEN READINGS SCREEN		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI58	The CMEMS PAP READING CARD PANEL consists of: CMEMS PAP READING HEADER COMPONENT CMEMS PAP READING DISPLAY COMPONENT OPTIMAL RANGE LABEL COMPONENT CMEMS PAP READING INSTRUCTIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI54	8.12.1.6.2 CMEMS PAP READING CARD PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2975	<Picture> CMEMS PAP READING CARD PANEL (default for HOME SCREEN: Day, default: No Reading, PA Visibility enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3194	<Picture> CMEMS PAP READING CARD PANEL (default for HOME SCREEN: Day, example: Reading, example: Reading In range, PA Visibility enabled)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

Commented [RC2]: Example: CMEMS PAP READING CARD PANEL

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4291	<Picture> CMEMS PAP READING CARD PANEL (default for HOME SCREEN: Day, example: Reading, example: Reading Out of Range, PA Visibility enabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI3195	<Picture> CMEMS PAP READING CARD PANEL (default for HOME SCREEN: Day, example: Reading taken, example: PA Visibility Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4488	<Picture> CMEMS PAP READING CARD PANEL (default for HOME SCREEN: Day, example: no reading, example: PA Visibility Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4888	<Picture> CMEMS PAP READING CARD PANEL (example: reading, example: Out of Home Elevation)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI56	8.12.1.6.3 CMEMS PAP READING CARD PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI59	8.12.1.6.3.1 CMEMS PAP READING HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI62	Component: CMEMS PAP READING HEADER COMPONENT State: Enabled/Disabled Display: Always displayed in the CMEMS PAP READING CARD PANEL. The default state is enabled. If daily reading has been taken, display the header text with white background as read-only (disabled). If daily reading has not been taken, display header text with blue background and allow component activation upon click (enabled). Label: <label.gb.cmems_readings_header> Behavior: Text and trigger. Upon activation in the enabled state, display the CMEMS PAP READING INSTRUCTIONS PANEL.	PString1013 CardioMEMS™ Reading papReadingTitle label.gb.cmems_readings_header	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2978	<Picture> CMEMS PAP READING HEADER COMPONENT (default: Enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2979	<Picture> CMEMS PAP READING HEADER COMPONENT (example: Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4395	8.12.1.6.3.3 CMEMS PAP READING DISPLAY COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4405	The CMEMS PAP READING DISPLAY COMPONENT consists of: READING STATUS COMPONENT READING VALUE COMPONENT READING UNIT COMPONENT		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4406	8.12.1.6.3.3.1 READING STATUS COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4400	Component: READING STATUS State: PA Visibility enabled/disabled, reading taken/not taken, READING VALUE COMPONENT	PString454	UI Info	CloudDRxRTR CloudDRxPS

Commented [RC3]: Reading Status Component

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Display: Always displayed in the Home Screen and Readings Screen in Day View Label: <msg.gb.readings_reminder> <label.gb.in_range> <label.gb.out_of_range> <msg.gb.reading_taken> <msg.gb.readingUnderReview> Behavior: Text and color will change depending on state.	Remember to take your reading. papReadingSubtitle msg.gb.readings_reminder PString1050 In Range pap_reading_in label.gb.in_range PString1051 Out of Range pap_reading_out label.gb.out_of_range PString1066 Reading taken readingTaken msg.gb.reading_taken PString1095 Reading under review.		
PAGUI4401	<Picture> READING STATUS COMPONENT (example: reading taken, in range, PA Visibility enabled)		UI Design	CloudDRxRTR CloudDRxPS
	In Range			
PAGUI4402	<Picture> READING STATUS COMPONENT (example: reading taken, out of range, PA Visibility enabled)		UI Design	CloudDRxRTR CloudDRxPS
	Out of Range			
PAGUI4403	<Picture> READING STATUS COMPONENT (reading not taken, PA visibility enabled/disabled)		UI Design	CloudDRxRTR CloudDRxPS
	Remember to take your reading.			
PAGUI4404	<Picture> READING STATUS COMPONENT (reading taken, PA Visibility disabled)		UI Design	CloudDRxRTR CloudDRxPS
	Reading taken			
PAGUI4887	<Picture> READING STATUS COMPONENT (reading taken, Out of Home Elevation)		UI Design	CloudDRxRTR
	Reading under review.			
PAGUI4407	8.12.1.6.3.3.2 READING VALUE COMPONENT		Heading	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI4408	Component: READING VALUE COMPONENT State: PA Visibility enabled/diabled, reading taken/not taken Display: Always displayed in the Home Screen and Readings Screen in Day View Label: N/A Behavior: Text or image will change depending on state and the reading value. Values are pulled from the backend.		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4409	<Picture> READING VALUE COMPONENT (example: reading taken, in range, PA Visibility enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4410	<Picture> READING VALUE COMPONENT (example: reading taken, out of range, PA Visibility enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4411	<Picture> READING VALUE COMPONENT (reading not taken, PA visibility enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4412	<Picture> READING VALUE COMPONENT (reading taken, PA visibility disabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4489	<Picture> READING VALUE COMPONENT (reading taken, PA visibility enabled) 		UI Design	CloudDRxRTR
PAGUI4413	8.12.1.6.3.3.3 READING UNIT COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4414	Component: READING UNIT COMPONENT State: PA Visibility enabled/disabled Display: Always displayed in the Home Screen and Readings Screen in Day View when PA Visibility is enabled. Label: <label.gb.mmhg> Behavior: Text will be present when PA Visibility is enabled. It will not present when PA Visibility is disabled.	PString452 mmHg mmHg label.gb.mmhg	UI Info	CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4415	<Picture> READING UNIT COMPONENT (example: PA Visibility enabled) mmHg		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4416	8.12.1.6.3.4 OPTIMAL RANGE LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4417	See Shared Elements: Components.		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3466	8.12.1.6.3.9 PAP READING INSTRUCTIONS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3467	8.12.1.6.3.9.1 PAP READING INSTRUCTIONS PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3470	The PAP READING INSTRUCTIONS PANEL is displayed when the user selects the CMEMS PAP READING HEADER COMPONENT in the enabled state. The PAP READING INSTRUCTIONS PANEL is displayed over the entire screen except the OS STATUS PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3471	The PAP READING INSTRUCTIONS PANEL consists of: CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) INSTRUCTIONS HEADER COMPONENT VIDEO ARTICLE ACTION TRIGGER COMPONENT (See Shared Elements: Components) VIDEO PLAYER PANEL (See Shared Elements: Panels) PAP INSTRUCTIONS LABEL COMPONENT GOT IT ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3468	8.12.1.6.3.9.2 PAP READING INSTRUCTIONS PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3472	<Picture> PAP READING INSTRUCTIONS PANEL (example: top of scrollable panel)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
	9:41 Instructions  1. Place device Place CardioMEMS onto the surface where you will lie down to take your readings. 2. Press 'Clear' button To take a reading begin by pressing the clear button on the back of your system device to turn it on. The button lights up blue when the system is on. 3. Remove handheld from storage Place the handheld unit from the storage cradle by the plastic handle and place near the system. Make sure to completely unwind the			
PAGUI3473	<Picture> PAP READING INSTRUCTIONS PANEL (exapmle: bottom of scrollable panel)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	4. Lie down on the pillow Lie down and a large green circle that reads start will appear on the screen of the handheld unit. Press and release the green button once avoid placing the handheld unit directly on your chest during a reading. Adjust your position on the pillow until the bar on the screen turns green and you hear good position on pillow. 5. Stay still If you hear shift slightly on pillow, change your body position by making slight movements towards the shoulder opposite of where your censor is located. Please be as still as possible during this 18 second reading period when the music finishes. 6. Reading complete! You may get up at this point, step away from the electronics system but do not push the power button or unplug your system. The reading is automatically stored and sent to your doctor. Got it			
PAGUI3469	8.12.1.6.3.9.3 PAP READING INSTRUCTIONS PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3474	8.12.1.6.3.9.3.1 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3475	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3476	8.12.1.6.3.9.3.2 INSTRUCTIONS HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3477	Component: INSTRUCTIONS HEADER COMPONENT State: N/A Display: Always displayed in the PAP READING INSTRUCTIONS PANEL. Label: <label.gb.instructions> Behavior: Text.	PString670 Instructions instructions label.gb.instructions	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3478	<Picture> INSTRUCTIONS HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3479	8.12.1.6.3.9.3.3 VIDEO ARTICLE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3480	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3481	8.12.1.6.3.9.3.4 VIDEO PLAYER PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3482	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3483	8.12.1.6.3.9.3.5 PAP INSTRUCTIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3484	Component: PAP INSTRUCTIONS LABEL COMPONENT State: N/A Display: Conditionally displayed in the PAP READING INSTRUCTIONS PANEL. The PAP INSTRUCTIONS LABEL COMPONENT is displayed when there are PAP reading instructions available in the content management location for the particular panel. The number of PAP INSTRUCTIONS LABEL COMPONENTS displayed in a panel is determined by how many unique instruction steps there are stored in the content management location. Label: <label.gb.reading_instructions> Behavior: Text.	PString922 1. Place device {n}Place CardioMEMS™ onto the surface where you will lie down to take your readings.2. Press ‘Clear’ button	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		{n}To take a reading begin by pressing the clear button on the back of your system device to turn it on. {n}The button lights up blue when the system is on.3. Remove handheld from storage {n}Place the handheld unit from the storage cradle by the plastic handle and place near the system. {n}Make sure to completely unwind the cable from its grove.4. Lie down on the pillow {n}Lie down and a large green circle that reads start will appear on the screen of the handheld unit. Press and release the green button once avoid placing the handheld unit directly on your chest during a reading. {n}Adjust your position on the pillow until the bar on the screen turns green and you hear good position on pillow.5. Stay still {n}If you hear shift slightly on pillow, change your body position by making slight movements towards the shoulder opposite of where your sensor is located. {n}Please be as still as possible during this 18 second reading period when the music finishes.6. Reading complete!		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

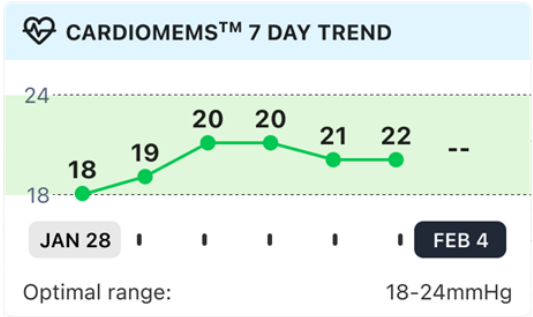
myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
		{n}You may get up at this point, step away from the electronics system but do not push the power button or unplug your system. The reading is automatically stored and sent to your doctor. papReadingsInstructionsText label.gb.reading_instruction		
PAGUI3485	<Picture> PAP INSTRUCTIONS LABEL COMPONENT (example: step 1) 1. Place device Place CardioMEMS onto the surface where you will lie down to take your readings.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3486	<Picture> PAP INSTRUCTIONS LABEL COMPONENT (example: step 2) 2. Press 'Clear' button To take a reading begin by pressing the clear button on the back of your system device to turn it on. The button lights up blue when the system is on.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3487	<Picture> PAP INSTRUCTIONS LABEL COMPONENT (example: step 3) 3. Remove handheld from storage Place the handheld unit from the storage cradle by the plastic handle and place near the system. Make sure to completely unwind the cable from its grove.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3488	<Picture> PAP INSTRUCTIONS LABEL COMPONENT (example: step 4)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

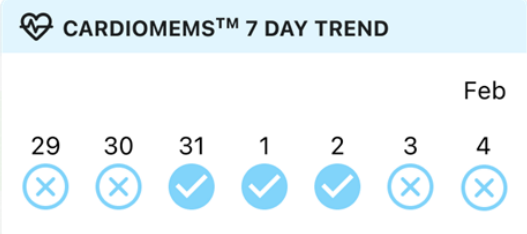
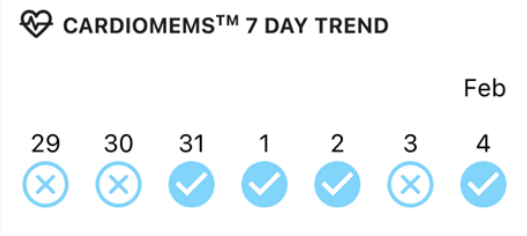
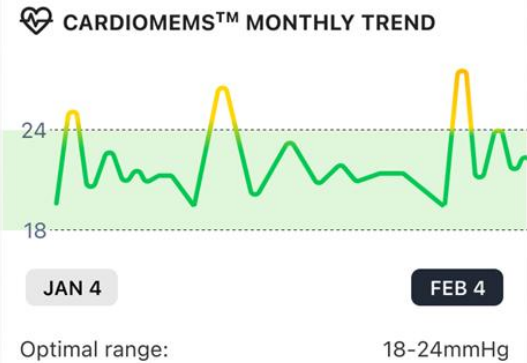
ID		English US Strings	Object Type	Project
	4. Lie down on the pillow Lie down and a large green circle that reads start will appear on the screen of the handheld unit. Press and release the green button once avoid placing the handheld unit directly on your chest during a reading. Adjust your position on the pillow until the bar on the screen turns green and you hear good position on pillow.			
PAGUI3489	<Picture> PAP INSTRUCTIONS LABEL COMPONENT (example: step 5) 5. Stay still If you hear shift slightly on pillow, change your body position by making slight movements towards the shoulder opposite of where your censor is located. Please be as still as possible during this 18 second reading period when the music finishes.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3490	<Picture> PAP INSTRUCTIONS LABEL COMPONENT (example: step 6) 6. Reading complete! You may get up at this point, step away from the electronics system but do not push the power button or unplug your system. The reading is automatically stored and sent to your doctor.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3491	8.12.1.6.3.9.3.6 GOT IT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3492	Component: GOT IT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the PAP READING INSTRUCTIONS PANEL. Label: <label.gb.got_it> Behavior: Text and trigger. Upon activation, close the PAP READING INSTRUCTIONS PANEL and display the most recently displayed screen.	PString671 Got it papReadingInstructionsButt onText label.gb.got_it	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

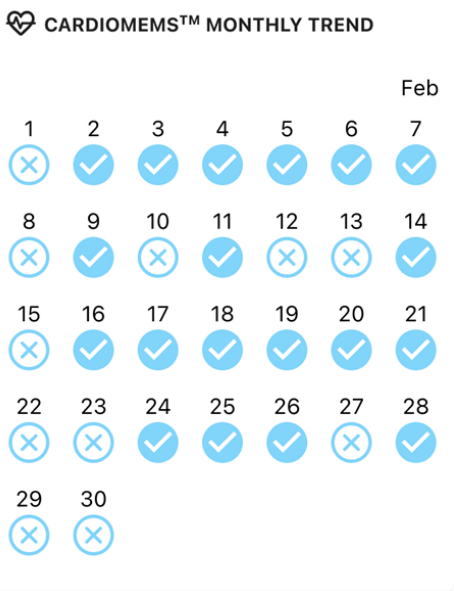
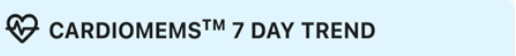
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3493	<Picture> GOT IT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4292	8.12.1.7 CMEMS PAP TREND CARD PANEL		Heading	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4295	8.12.1.7.1 CMEMS PAP TREND CARD PANEL (Display)		Heading	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4296	The CMEMS PAP TREND CARD PANEL displays PAP trend information. The CMEMS PAP TREND CARD PANEL displays different PAP trend information depending on the selected range of dates (week or month). If PA Visibility is enabled, trend information is shown in graphical format. If PA Visibility is disabled, trend information displays an indicator for if the patient has taken a reading for each day in the trend.		UI Info	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4304	The CMEMS PAP TREND CARD PANEL consists of: CMEMS PAP TREND BANNER COMPONENT CMEMS PAP TREND DISPLAY PANEL		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4303	The CMEMS PAP TREND CARD PANEL is an element of: HOME SCREEN READINGS SCREEN		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4294	8.12.1.7.2 CMEMS PAP TREND CARD PANEL (Layout)		Heading	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4305	<Picture> CMEMS PAP TREND CARD PANEL (default for HOME SCREEN: Week, default: No Reading, example: PA Visibility enabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4306	<Picture> CMEMS PAP TREND CARD PANEL (for READING SCREEN: Week, default: No Reading, example: PA Visibility disabled)		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR

myCardioMEMS Patient Application Graphical User Interface Specification

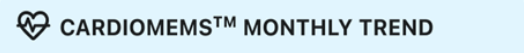
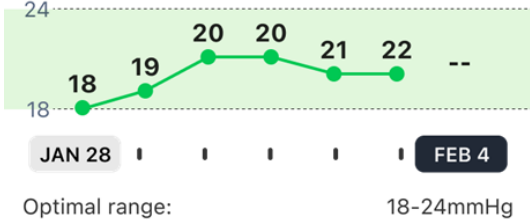
ID		English US Strings	Object Type	Project
				
PAGUI4310	<Picture> CMEMS PAP TREND CARD PANEL (for READING SCREEN: Week, Reading taken, example: PA Visibility disabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4307	<Picture> CMEMS PAP TREND CARD PANEL (for READING SCREEN: Month: Reading taken, example: PA Visibility enabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4309	<Picture> CMEMS PAP TREND CARD PANEL (for READING SCREEN: Month: Reading not taken, example: PA Visibility disabled)		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4293	8.12.1.7.3 CMEMS PAP TREND CARD PANEL (Elements)		UI Info	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4311	8.12.1.7.3.1 CMEMS PAP TREND BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4312	Component: CMEMS PAP TREND BANNER COMPONENT State: Enabled/Disabled Display: Always displayed in the CMEMS PAP TREND CARD PANEL. The default state is enabled. If daily reading has been taken, display the header text with white background (disabled). If daily reading has not been taken, display header text with blue background (enabled). The default state is enabled. Label: <label.gb.cardiomems_trend_title_week>; Text is all capitalized. Behavior: Text.	PString1053 CardioMEMS™ 7 Day Trend home_weekTrendTitle label.gb.cardiomems_trend_title_week PString1086 CardioMEMS™ Monthly Trend home_monthTrendTitle label.gb.cardiomems_trend_title_month	UI Info	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4313	<Picture> CMEMS PAP TREND BANNER COMPONENT (default: week, enabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4314	<Picture> CMEMS PAP TREND BANNER COMPONENT (example: week, disabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR

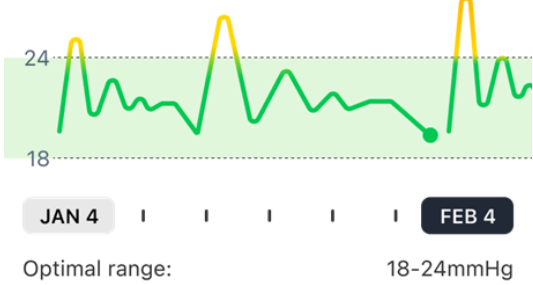
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

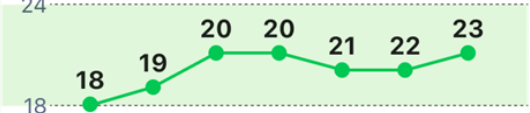
ID		English US Strings	Object Type	Project
PAGUI4315	<Picture> CMEMS PAP TREND BANNER COMPONENT (example: month, enabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4316	<Picture> CMEMS PAP TREND BANNER COMPONENT (default: month, disabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4453	8.12.1.7.3.2 CMEMS PAP TREND DISPLAY PANEL		Heading	CloudDRxRTR
PAGUI4454	See: Shared Elements: Panels		Info	CloudDRxRTR
PAGUI4317	8.12.1.8 CMEMS PAP TREND DISPLAY PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI4455	8.12.1.8.1 CMEMS PAP TREND DISPLAY PANEL (DISPLAY)		Heading	CloudDRxRTR
PAGUI4458	The CMEMS PAP TREND DISPLAY PANEL displays PAP trend information. The CMEMS PAP TREND DISPLAY PANEL displays different PAP trend information depending on the selected range of dates (week or month). If PA Visibility is enabled, trend information is shown in graphical format. If PA Visibility is disabled, trend information displays an indicator for if the patient has taken a reading for each day in the trend. In the DRX MEDICATION MESSAGE CARD PANEL, the CMEMS PAP TREND DISPLAY PANEL is displayed only in the Week state with PA Visibility enabled.		UI Info	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4459	The CMEMS PAP TREND DISPLAY PANEL consists of: CMEMS PAP TREND VIEW COMPONENT TREND DATE LABEL COMPONENT OPTIMAL RANGE LABEL COMPONENT		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4461	The CMEMS PAP TREND DISPLAY PANEL is an element of: CMEMS PAP TREND CARD PANEL DRX MEDICATION MESSAGE CARD PANEL		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4457	8.12.1.8.2 CMEMS PAP TREND DISPLAY PANEL (LAYOUT)		Heading	CloudDRxRTR
PAGUI4465	<Picture> CMEMS PAP TREND DISPLAY PANEL (example, Week, PA Visibility enabled) 		UI Design	CloudDRxRTR
PAGUI4466	<Picture> CMEMS PAP TREND DISPLAY PANEL (example: Month, PA Visibility enabled)		UI Design	CloudDRxRTR



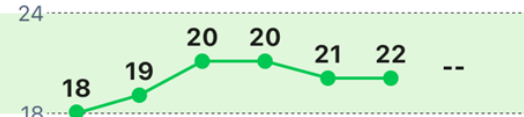
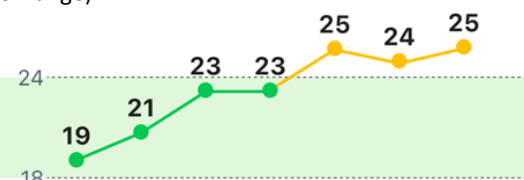
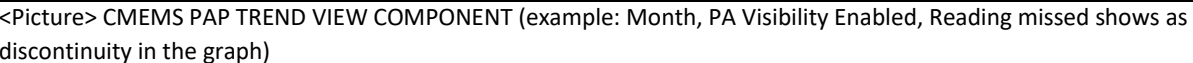
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4467	<Picture> CMEMS PAP TREND DISPLAY PANEL (example: Week, PA Visibility disabled) Feb 29 30 31 1 2 3 4 		UI Design	CloudDRxRTR
PAGUI4472	<Picture> CMEMS PAP TREND DISPLAY PANEL (example: Month, PA Visibility disabled) Feb 1 2 3 4 5 6 7 		UI Design	CloudDRxRTR
PAGUI4456	8.12.1.8.3 CMEMS PAP TREND DISPLAY PANEL (ELEMENTS)		Heading	CloudDRxRTR
PAGUI4460	8.12.1.8.3.1 CMEMS PAP TREND VIEW COMPONENT		Heading	CloudDRxRTR
PAGUI4318	Component: CMEMS PAP TREND VIEW COMPONENT State: Week/Month and Reading/No Reading and PA Visibility Enabled/Disabled		UI Info	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Display: Displayed in the CMEMS PAP TREND CARD PANEL. Default state is No Reading and PA Visibility enabled.The CMEMS PAP TREND VIEW COMPONENT is displayed by default in the week state. Conditionally displayed in the DRX MEDICATION MESSAGE PANEL if PA Visibility is enabled. The CMEMS PAP TREND VIEW COMPONENT will be displayed in the week state. Label: N/A Behavior: Text and graph or indicators. In the DRX MEDICATION MESSAGE CARD PANEL, display the available readings in the week state, when PA Visibility is enabled. In the READINGS SCREEN, display the number of available readings per the interval specified in the READINGS INTERVAL SELECTION COMPONENT and READINGS DATE SELECTION COMPONENT. Display available readings with graph or indicators depending on PA visibility enabled/disabled. In the HOME SCREEN, display the available readings in the week state with graph or indicators depending on PA visibility enabled/disabled. PA Visibility Disabled The CMEMS PAP TREND VIEW COMPONENT displays an indication of if the reading has been completed for each date within the specified interval. PA Visibility Enabled In the no reading state, the CMEMS PAP TREND VIEW COMPONENT displays dashes for no readings in the week state, and displays nothing for missing readings in the month state. In the week state, the CMEMS PAP TREND VIEW COMPONENT displays a trend graphic with all available readings within the specified interval. In the month state, the CMEMS PAP TREND VIEW COMPONENT displays a trend graphic with all available readings within the specified interval. The graphic displays a different background color for the 'optimal range' values in the week state and month state with PA Visibility enabled. The trend line in the display changes color if the trend line leaves the optimal range.			
PAGUI4319	<Picture> CMEMS PAP TREND VIEW COMPONENT (default: Week, Reading taken, PA Visibility Enabled) 		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4321	<Picture> CMEMS PAP TREND VIEW COMPONENT (default: Week, Reading taken, PA Visibility Disabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4322	<Picture> CMEMS PAP TREND VIEW COMPONENT (example: Week, No Reading, PA Visibility Enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4323	<Picture> CMEMS PAP TREND VIEW COMPONENT (example: Week, Reading taken, PA Visibility Enabled, some readings out of range) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4324	<Picture> CMEMS PAP TREND VIEW COMPONENT (example: Month, PA Visibility Disabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4325	<Picture> CMEMS PAP TREND VIEW COMPONENT (example: Month, PA Visibility Enabled, Reading missed shows as discontinuity in the graph) 		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4334	8.12.1.8.3.2 TREND DATE LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4335	8.12.1.8.3.2.1 <i>Component: TREND DATE LABEL COMPONENT</i> State: N/A Display: Always displayed in the CMEMS PAP TREND DISPLAY PANEL. Label: <label.gb.month_label>, <label.gb.month_day_start>, <label.gb.month_day_end> Behavior: Text and values When PA Visibility is enabled, the TREND DATE LABEL COMPONENT displays the start month and day and the end month and day of the trend. When PA Visibility is disabled, the TREND DATE LABEL COMPONENT displays only the month label.	PString1069 {0} {1} label.gb.month_day_start PString1070 {0} {1} label.gb.month_day_end PString1071 {0} label.gb.month_label	UI Info	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4339	<Picture> TREND DATE LABEL COMPONENT (example: PA Visibility enabled) 		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4340	<Picture> TREND DATE LABEL COMPONENT (example: PA Visibility disabled) Feb		UI Design	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4327	8.12.1.8.3.3 OPTIMAL RANGE LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS CloudDRxLMR
PAGUI4478	See Shared Elements: Components		Info	CloudDRxRTR
PAGUI2852	8.12.1.9 CONFIRM INPUT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2853	8.12.1.9.1 CONFIRM INPUT PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2854	The CONFIRM INPUT PANEL allows users to review input before the value is saved; the user can then choose to return to the original input panel or save the value.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2875	The CONFIRM INPUT PANEL is an element of:		UI Info	CloudDRx 1.0

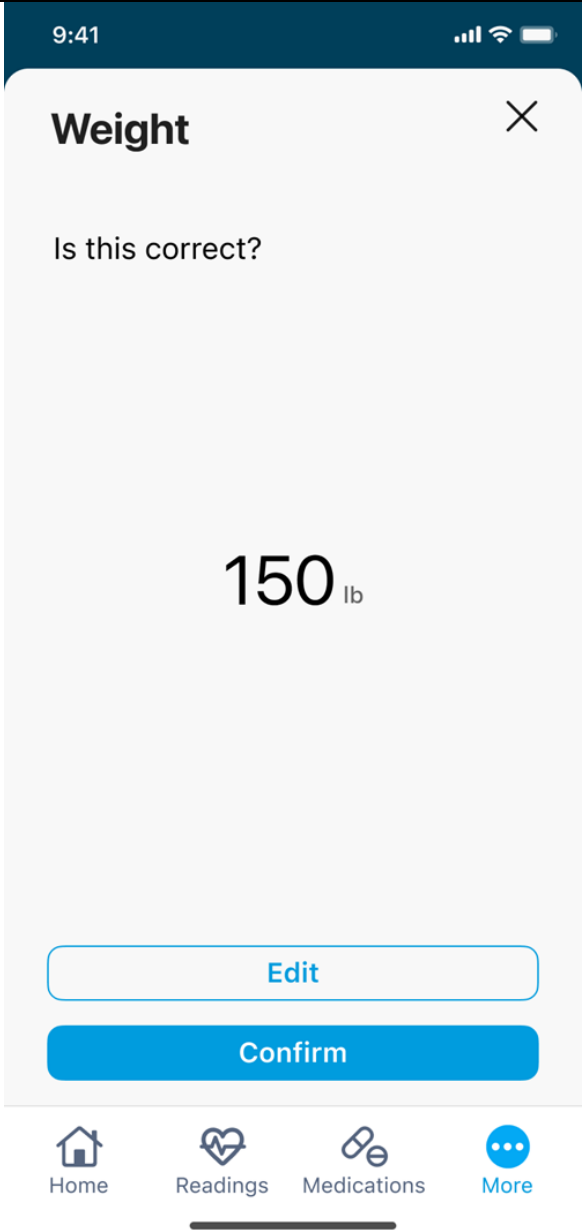
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



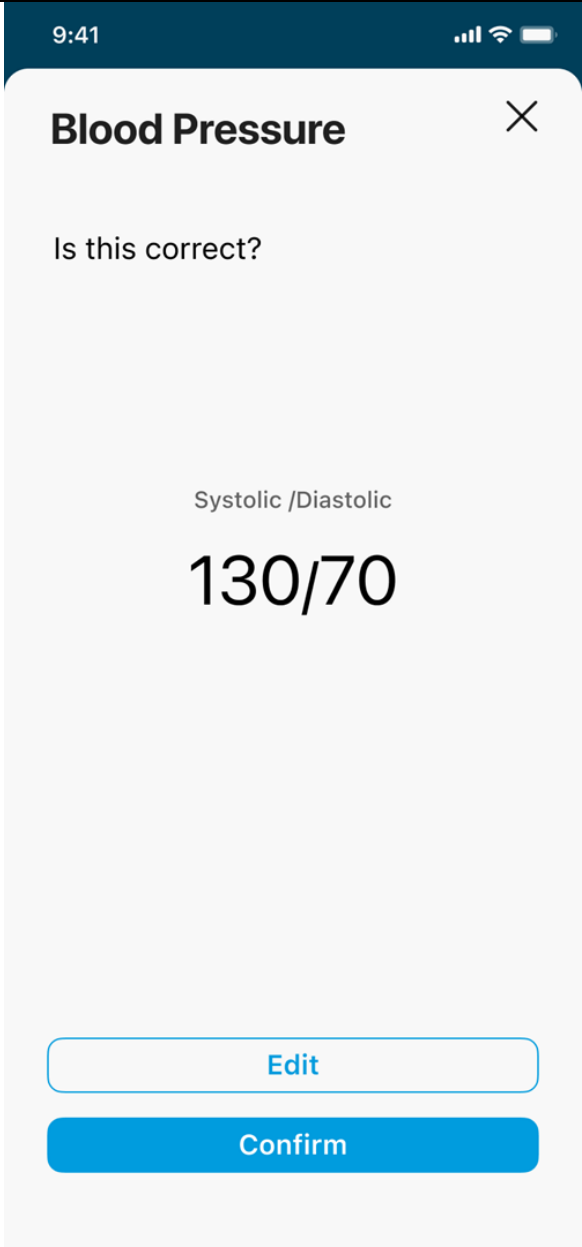
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	RECORD WEIGHT PANEL RECORD BP READING PANEL			CloudDRxRTR CloudDRxPS
PAGUI2855	The CONFIRM INPUT PANEL consists of: CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) CONFIRM INPUT HEADER COMPONENT CORRECT LABEL COMPONENT CURRENT INPUT LABEL COMPONENT EDIT ACTION TRIGGER COMPONENT (See Shared Elements: Components) CONFIRM ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2856	8.12.1.9.2 CONFIRM INPUT PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2857	<Picture> CONFIRM INPUT PANEL (example: Weight)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4341	<Picture> CONFIRM INPUT PANEL (example: Blood Pressure)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI2858	8.12.1.9.3 CONFIRM INPUT PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2859	8.12.1.9.3.1 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2860	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2861	8.12.1.9.3.2 CONFIRM INPUT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2862	Component: CONFIRM INPUT HEADER COMPONENT State: Weight / Blood Pressure Display: Always displayed as a component of the CONFIRM INPUT PANEL. Displays the weight state when navigated to from the RECORD WEIGHT PANEL. Displays the Blood pressure reading when navigated to from the RECORD BP READING PANEL. Label: <label.gb.weight> in Weight state or <label.gb.blood_pressure> in Blood Pressure State Behavior: Text.	PString521 Weight weight label.gb.weight PString616 Blood Pressure bloodPressure label.gb.blood_pressure	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2863	<Picture> CONFIRM INPUT HEADER COMPONENT (example: Weight) Weight		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2864	8.12.1.9.3.3 CORRECT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2865	Component: CORRECT LABEL COMPONENT State: N/A Display: Always displayed as a component of the CONFIRM INPUT PANEL. Label: <label.gb.is_correct> Behavior: Text.	PString556 Is this correct? isThisCorrect label.gb.is_correct	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2866	<Picture> CORRECT LABEL COMPONENT Is this correct?		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2867	8.12.1.9.3.4 CURRENT INPUT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2868	Component: CURRENT INPUT LABEL COMPONENT State: Weight or Blood Pressure Display: Always displayed as a component of the CONFIRM INPUT PANEL. Displays the weight state when navigated to from the RECORD WEIGHT PANEL. Displays the blood pressure reading when navigated to from the RECORD BP READING PANEL Label: <data.Code.Weight_Unit_LB>, <data.Code.Med_Unit_Cd.Code_Desc.KG>	PString342 kg kg data.Code.Med_Unit_Cd.Code_Desc.KG PString520	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

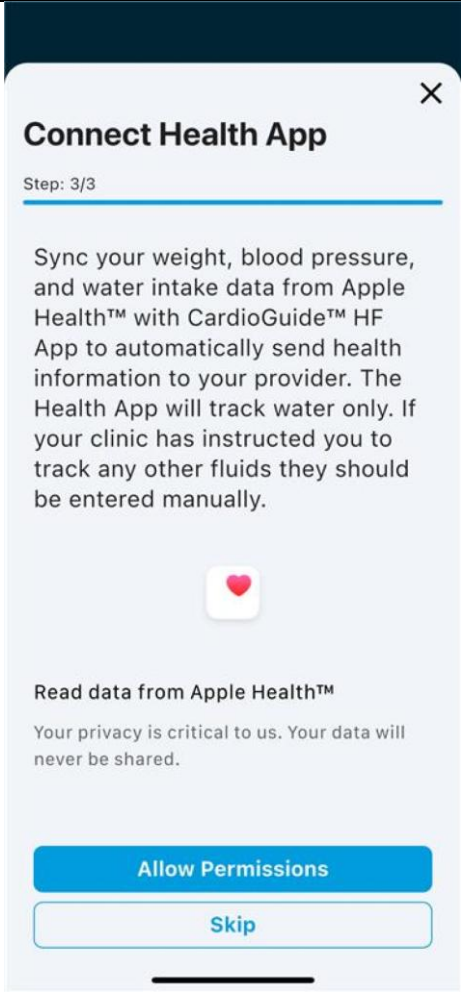
ID		English US Strings	Object Type	Project
	<label.gb.bp_systolic>, <label.gb.bp_diastolic> Behavior: Text. The CURRENT INPUT LABEL COMPONENT displays the input value from the previous panel and the label or associated with that input.	lb lb data.Code.Weight_Unit_LB		
PAGUI2869	<Picture> CONFIRM INPUT INPUT COMPONENT (example: Weight) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4381	<Picture> CONFIRM INPUT INPUT COMPONENT (example: Blood Pressure) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2871	8.12.1.9.3.5 EDIT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2878	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2880	8.12.1.9.3.6 CONFIRM ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2881	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1036	8.12.1.10 CONNECT HEALTH APP PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1037	8.12.1.10.1 CONNECT HEALTH APP PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1123	The CONNECT HEALTH APP PANEL is displayed in REGISTRATION SCREEN and after the user navigates to it from the SETTINGS PANEL. The CONNECT HEALTH APP PANEL is displayed over the entire screen except the OS STATUS PANEL and GLOBAL HEADING PANEL. The default state is not connected. From registration, if the user chooses to connect a supported health application and complete the connection process, the HOME SCREEN will be displayed when the user returns to the Patient Application. From the SETTINGS PANEL, the user returns to the SETTINGS PANEL after the user chooses to connect a supported health application. The CONNECT HEALTH APP PANEL is only displayed when the user is on a device wherein the supported health application(s) is(are) available.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1124	The CONNECT HEALTH APP PANEL consists of: CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) CONNECT HEALTH APP HEADER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



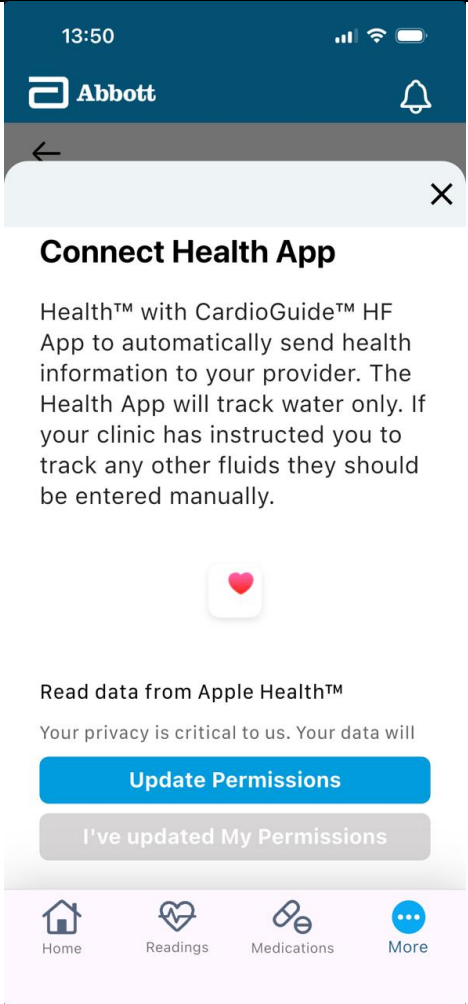
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	REGISTRATION PROGRESS BAR COMPONENT (See Shared Elements: Components) SYNC HEALTH DATA INSTRUCTION LABEL COMPONENT HEALTH APP ICON COMPONENT DATA PRIVACY LABEL COMPONENT ALLOW PERMISSIONS ACTION TRIGGER COMPONENT UPDATE PERMISSIONS INFORMATION PANEL SKIP ACTION TRIGGER COMPONENT UPDATE PERMISSIONS ACTION TRIGGER COMPONENT CONFIRM UPDATE PERMISSIONS ACTION TRIGGER COMPONENT UPDATE PERMISSIONS INSTRUCTION PANEL			
PAGUI1038	8.12.1.10.2 CONNECT HEALTH APP PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1125	<Picture> CONNECT HEALTH APP PANEL (default: not connected in REGISTRATION SCREEN; example: Apple device)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1341	<Picture> CONNECT HEALTH APP PANEL (default: not connected in SETTINGS PANEL, example: Apple device)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

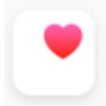
ID		English US Strings	Object Type	Project
				
PAGUI1039	8.12.1.10.3 CONNECT HEALTH APP PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1128	8.12.1.10.3.1 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1143	See Shared Elements: Components		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1129	8.12.1.10.3.2 CONNECT HEALTH APP LABEL COMPONENT		Heading	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1144	Component: CONNECT HEALTH APP LABEL COMPONENT State: N/A Display: Always displayed in the CONNECT HEALTH APP PANEL. Label: <label.gb.connect_health_app> Behavior: Text.	PString557 Connect Health App connectHealthApp label.gb.connect_health_app	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1145	<Picture> CONNECT HEALTH APP LABEL COMPONENT (example: Apple) Connect Health App		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1345	8.12.1.10.3.3 REGISTRATION PROGRESS BAR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1346	See Shared Elements: Components		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1130	8.12.1.10.3.4 SYNC HEALTH DATA INSTRUCTION LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1146	Component: SYNC HEALTH DATA INSTRUCTION LABEL COMPONENT State: N/A Display: Always displayed in the CONNECT HEALTH APP PANEL. Label: <msg.gb.connect_health_app> Behavior: Text.	PString1030 Sync your weight, blood pressure, and water intake data from {19 or 20} with {0} to automatically send health information to your provider. The Health App will track water only. If your clinic has instructed you to track any other fluids they should be entered manually. syncFromGoogleFit syncFromAppleHealth msg.gb.connect_health_app	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1147	<Picture> SYNC HEALTH DATA INSTRUCTION LABEL COMPONENT (example: Apple)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Sync your weight, blood pressure, and water intake data from Apple Health™ with CardioGuide™ HF App to automatically send health information to your provider. The Health App will track water only. If your clinic has instructed you to track any other fluids they should be entered manually.			
PAGUI1131	8.12.1.10.3.5 HEALTH APP ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1149	Component: HEALTH APP ICON COMPONENT State: N/A Display: Always displayed in the CONNECT HEALTH APP PANEL. Label: N/A Behavior: Icon.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1150	<Picture> HEALTH APP ICON COMPONENT (example: Apple) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1132	8.12.1.10.3.6 DATA PRIVACY LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1152	Component: DATA PRIVACY LABEL COMPONENT State: N/A Display: Always displayed in the CONNECT HEALTH APP PANEL. Label: <label.gb.read_data>, <msg.gb.privacy> Behavior: Text.	PString559 Read data from {19 or 20}. readDataFromAppleHealth readDataFromGoogleFit label.gb.read_data PString560 Your privacy is critical to us. Your data will never be shared. yourPrivacylsCriticalToUs msg.gb.privacy	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1153	<Picture> DATA PRIVACY LABEL COMPONENT (example: Apple)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Read data from Apple Health™ Your privacy is critical to us. Your data will never be shared.			
PAGUI1133	8.12.1.10.3.7 ALLOW PERMISSIONS ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1155	Component: ALLOW PERMISSIONS ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the CONNECT HEALTH APP PANEL if in the REGISTRATION SCREEN. The ALLOW PERMISSIONS TRIGGER COMPONENT is not displayed when the CONNECT HEALTH APP PANEL is accessed from the SETTINGS PANEL. Label: <label.gb.allow_permissions> Behavior: Text and trigger. Upon activation, the ALLOW PERMISSIONS ACTION TRIGGER COMPONENT displays a settings panel directly connected with the health application. When the user completes the SOUP permissions flow, the HOME SCREEN is displayed.	PString561 Allow Permissions allowPermissions label.gb.allow_permissions	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1156	<Picture> ALLOW PERMISSIONS ACTION TRIGGER COMPONENT (example: Apple) Allow Permissions		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1348	8.12.1.10.3.8 SKIP ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1349	Component: SKIP ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the CONNECT HEALTH APP PANEL if in the REGISTRATION SCREEN. The SKIP ACTION TRIGGER COMPONENT is not displayed when the CONNECT HEALTH APP PANEL is accessed from the SETTINGS PANEL. Label: <label.gb.skip> Behavior: Text and trigger. Upon activation, the SKIP ACTION TRIGGER COMPONENT closes the REGISTRATION SCREEN and displays the HOME SCREEN.	PString294 Skip globalSkip label.gb.skip	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1350	<Picture> SKIP ACTION TRIGGER COMPONENT (example: Apple) Skip		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1566	8.12.1.10.3.9 UPDATE PERMISSIONS ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1567	Component: UPDATE PERMISSIONS ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Conditionally displayed in the CONNECT HEALTH APP PANEL if in SETTINGS PANEL. The default state is enabled. Label: <label.gb.update_permissions>	PString562 Update Permissions updatePermissions label.gb.update_permissions	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

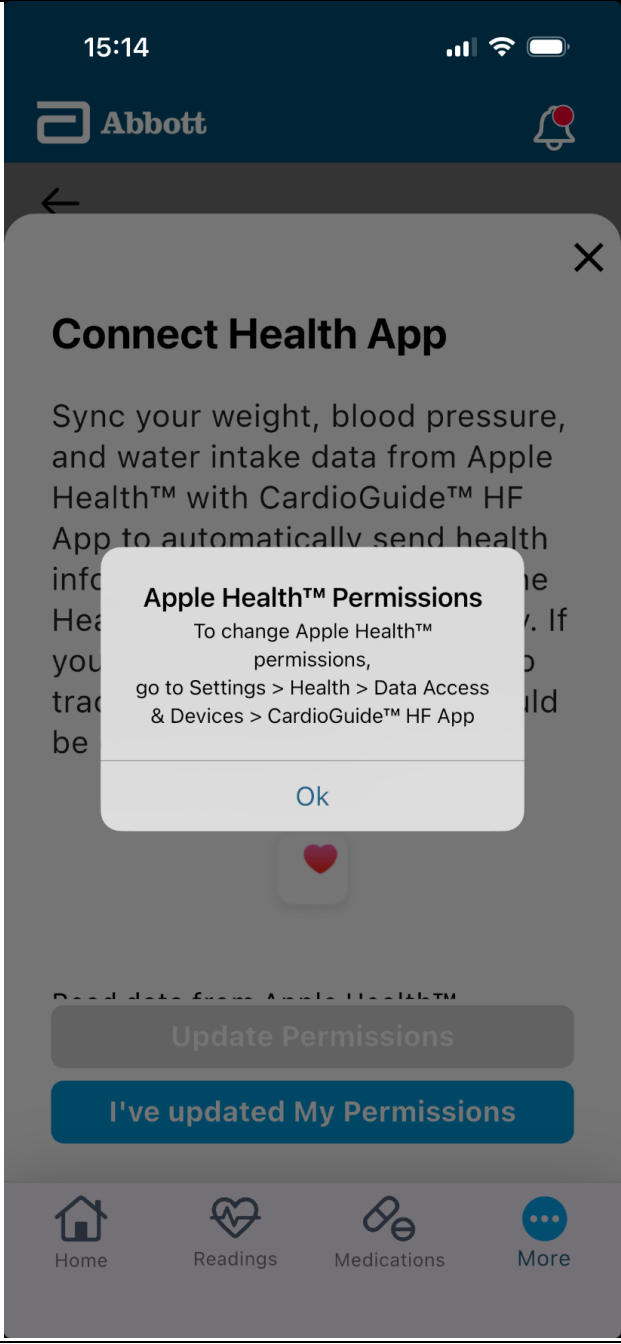
ID		English US Strings	Object Type	Project
	Behavior: Text and trigger. Upon activation in the enabled state and before the health app is connected, the UPDATE PERMISSIONS ACTION TRIGGER COMPONENT displays a settings panel directly connected with the health application. Upon activation in the enabled state and after the health app is connected, the UPDATE PERMISSIONS ACTION TRIGGER COMPONENT is disabled and displays the UPDATE PERMISSIONS INSTRUCTIONS PANEL. The UPDATE PERMISSIONS ACTION TRIGGER COMPONENT changes from the the disabled state to the enabled state after the user selects the CONFIRM UPDATE PERMISSIONS ACTION TRIGGER COMPONENT.			
PAGUI1568	<Picture> UPDATE PERMISSIONS ACTION TRIGGER COMPONENT (default: Enabled, example: Apple) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3388	<Picture> UPDATE PERMISSIONS ACTION TRIGGER COMPONENT (example: Disabled, example: Apple) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3384	8.12.1.10.3.10 CONFIRM UPDATE PERMISSIONS ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3385	Component: UPDATE PERMISSIONS ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Conditionally displayed in the CONNECT HEALTH APP PANEL if in SETTINGS PANEL. The default state is disabled. Label: <label.gb.confirm_update_permissions> Behavior: Text and trigger. The UPDATE PERMISSIONS ACTION TRIGGER COMPONENT changes from the disabled state to the enabled state after the user selects the UPDATE PERMISSIONS ACTION TRIGGER COMPONENT.	PString679 I've updated My Permissions permission_update_acknow ledgement label.gb.confirm_update_p ermissons	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3386	<Picture> CONFIRM UPDATE PERMISSIONS ACTION TRIGGER COMPONENT (default: Disabled, example: Apple) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3387	<Picture> CONFIRM UPDATE PERMISSIONS ACTION TRIGGER COMPONENT (example: Enabled, example: Apple) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1569	8.12.1.10.3.11 UPDATE PERMISSIONS INSTRUCTION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1570	8.12.1.10.3.11.1 UPDATE PERMISSIONS INSTRUCTION PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1573	The UPDATE PERMISSIONS INSTRUCTION PANEL is displayed as a component of the CONNECT HEALTH APP PANEL when the user selects the UPDATE PERMISSIONS ACTION TRIGGER COMPONENT and the health app has been connected in the past. The UPDATE PERMISSIONS INSTRUCTION PANEL is displayed over the CONNECT HEALTH APP PANEL and the background is greyed out. The UPDATE PERMISSIONS INSTRUCTION PANEL is specific to Apple.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1575	The UPDATE PERMISSIONS INSTRUCTION PANEL consists of: UPDATE PERMISSIONS HEADER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	UPDATE PERMISSIONS LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (Shared Elements: Components)			CloudDRxPS
PAGUI1571	8.12.1.10.3.11.2 UPDATE PERMISSIONS INSTRUCTION PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1574	<Picture> UPDATE PERMISSIONS INSTRUCTION PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1572	8.12.1.10.3.11.3 UPDATE PERMISSIONS INSTRUCTION PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4603	8.12.1.10.3.11.3.1 UPDATE PERMISSIONS HEADER COMPONENT		Heading	CloudDRxRTR
PAGUI4604	Component: UPDATE PERMISSIONS HEADER COMPONENT State: N/A Display: Always displayed in the UPDATE PERMISSIONS INSTRUCTION PANEL. Label: <msg.gb.health_permissions> Behavior: Text.	PString1082 {0} Permissions healthPermissionsUpdate msg.gb.health_permissions	UI Info	CloudDRxRTR
PAGUI4605	<Picture> UPDATE PERMISSIONS HEADER COMPONENT (example: Apple) Apple Health™ Permissions		UI Design	CloudDRxRTR
PAGUI1576	8.12.1.10.3.11.3.2 UPDATE PERMISSIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1577	Component: UPDATE PERMISSIONS LABEL COMPONENT State: N/A Display: Always displayed in the UPDATE PERMISSIONS INSTRUCTION PANEL. Label: <msg.gb.aapple_permissions_header>, <msg.gb.update_apple_permissions> Behavior: Text.	PString762 To change {19} permissions, go to Settings > Health > Data Access & Devices > {0} toChangePermissions msg.gb.update_apple_perm issions	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1578	<Picture> UPDATE PERMISSIONS LABEL COMPONENT To change Apple Health™ permissions, go to Settings > Health > Data Access & Devices > CardioGuide™ HF App		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1579	8.12.1.10.3.11.3.3 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1580	See Shared Elements: Components		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2178	8.12.1.14 DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2179	8.12.1.14.1 DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL (Display)		Heading	CloudDRx 1.0

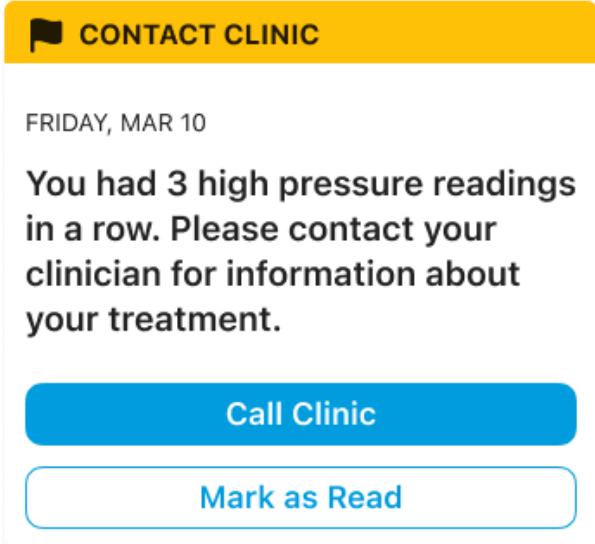
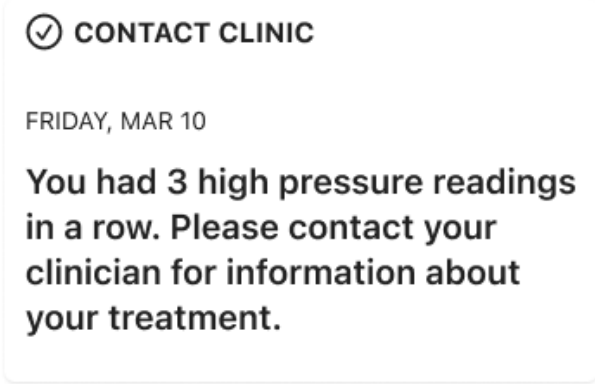
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI2180	The DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL is conditionally displayed if a DRx Call Clinic DirectCall Message is received by the Patient Application at least once since the Earliest History Calendar Date. Patients with DTP enabled will be able to recieve a DRx Call Clinic message. A unique DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL is displayed for every DRx Call Clinic DirectCall Message that is received by the Patient Application. Multiple DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANELS may be displayed in the same screen or panel. The DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL is in the active state by default. The DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL CARD PANEL changes state from active to complete when the user completes all required actions in the card. The elements in the DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL listed below change behavior along these two states.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2181	The DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL is an element of: ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4168	The DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL consists of: DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) DRX TREND MESSAGE LABEL COMPONENT (See Shared Elements: Components) CALL CLINIC CONTENT LABEL COMPONENT (See: Shared Elements: Components) CALL CLINIC ACTION TRIGGER COMPONENT (See Shared Elements: Components) MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI2183	8.12.1.14.2 DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4169	<Picture> DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL (default: active)		UI Design	CloudDRxRTR CloudDRxPS

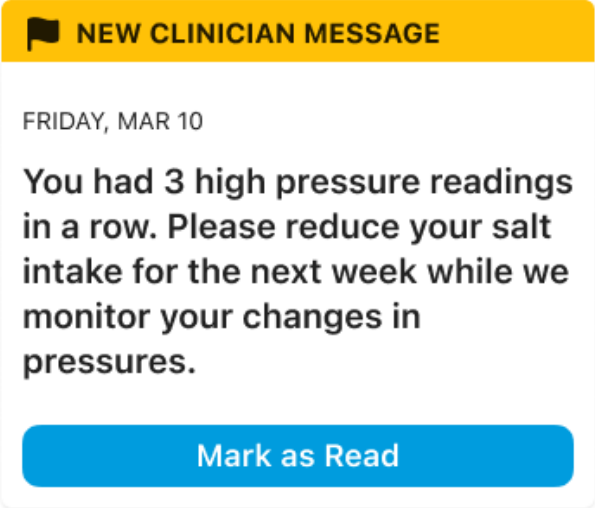
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4170	<Picture> DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL (complete) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2187	8.12.1.14.3 DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2198	8.12.1.14.3.1 DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4427	Component: DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL BANNER COMPONENT State: Active/ Complete Display: Always displayed in the DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL. The default state is active. In the active state, the DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL BANNER COMPONENT displays the active icon and	PString604 Contact Clinic contactClinic label.gb.contact_clinic	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

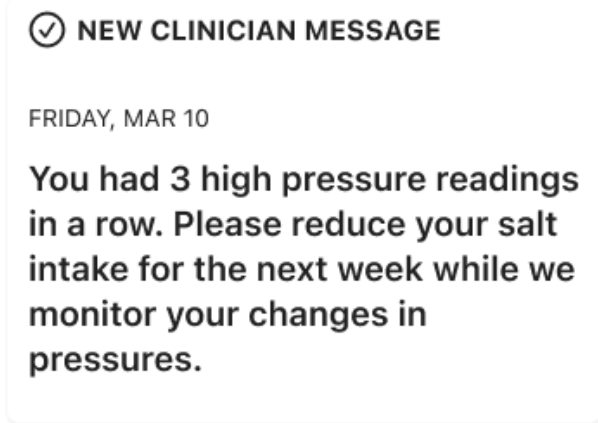
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	the active background color. In the complete state, the DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL BANNER COMPONENT displays the complete icon and complete background color. Label: <label.gb.contact_clinic>; Text is all capitalized. Behavior: Icon and text. The DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL BANNER COMPONENT changes state with the message card panel it is an element of.			
PAGUI2199	<Picture> DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL BANNER COMPONENT (active) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4372	<Picture> DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL BANNER COMPONENT (complete) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2200	8.12.1.14.3.3 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2201	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4171	8.12.1.14.3.4 DRX TREND MESSAGE LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4172	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI2205	8.12.1.14.3.5 CALL CLINIC CONTENT LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4385	See Shared Elements: Components		Info	CloudDRxRTR CloudDRxPS
PAGUI2210	8.12.1.14.3.6 CALL CLINIC ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2209	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2208	8.12.1.14.3.7 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2211	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2218	8.12.1.15 DRX FREEFORM MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2219	8.12.1.15.1 DRX FREEFORM MESSAGE CARD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2220	The DRX FREEFORM MESSAGE CARD PANEL is conditionally displayed if a DRx Freeform Message is received by the Patient Application at least once since the Earliest History Calendar Date. Patients with DTP enabled will be able to recieve a DRx Freeform message. A unique DRX FREEFORM MESSAGE CARD PANEL is displayed for every DRx Freeform Message that is received by the Patient Application. Multiple DRX FREEFORM MESSAGE CARD PANELS may be displayed in the same screen or panel. The DRX FREEFORM MESSAGE CARD PANEL is in the active state by default. The DRX FREEFORM MESSAGE CARD PANEL CARD PANEL changes state from active to complete when the user completes all required actions in the card. The elements in the DRX FREEFORM MESSAGE CARD PANEL listed below changes behavior along these two states.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2221	The DRX FREEFORM MESSAGE CARD PANEL is an element of: ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4175	The DRX FREEFORM MESSAGE CARD PANEL consists of: DRX FREEFORM MESSAGE BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) DRX TREND MESSAGE LABEL COMPONENT (See Shared Elements: Components) FREEFORM CONTENT LABEL COMPONENT MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI2223	8.12.1.15.2 DRX FREEFORM MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4176	<Picture> DRX FREEFORM MESSAGE CARD PANEL (default: active, example, DTP On) 		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4177	<Picture> DRX FREEFORM MESSAGE CARD PANEL (complete, example) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2226	8.12.1.15.3 DRX FREEFORM MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2227	8.12.1.15.3.1 DRX FREEFORM MESSAGE BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4426	Component: DRX FREEFORM BANNER COMPONENT State: Active/ Complete Display: Always displayed in the DRX FREEFORM MESSAGE CARD PANEL. The default state is active. In the active state, the DRX FREEFORM BANNER COMPONENT displays the active icon and the active background color. In the complete state, the DRX FREEFORM BANNER COMPONENT displays the complete icon and complete background color. Label: <label.gb.new_message>; Text is all capitalized. Behavior: Icon and text. The DRX FREEFORM BANNER COMPONENT changes state with the message card panel it is an element of.	PString179 New Clinician Message HCM label.gb.new_message	UI Info	CloudDRxRTR
PAGUI2228	<Picture> DRX FREEFORM MESSAGE BANNER COMPONENT (Default: active, DTP on/Off) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4362	<Picture> DRX FREEFORM MESSAGE BANNER COMPONENT (complete, DTP On/Off) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2229	8.12.1.15.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2230	See Shared Elements: Components.		Info	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



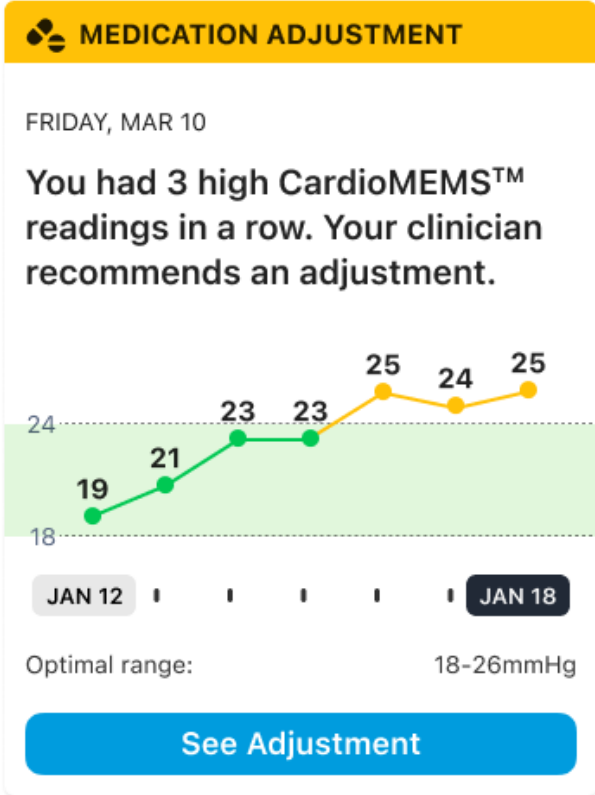
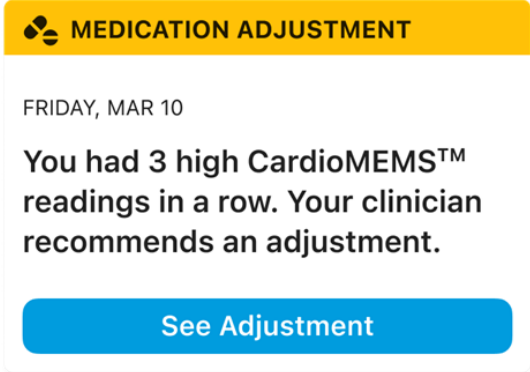
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI4173	8.12.1.15.3.4 DRX TREND MESSAGE LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4174	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI2234	8.12.1.15.3.5 DRX FREEFORM CONTENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2235	Component: DRX FREEFORM CONTENT LABEL COMPONENT State: N/A Display: Always displayed in the DRX FREEFORM MESSAGE CARD PANEL. Label: N/A Behavior: Text. DRX The FREEFORM CONTENT LABEL COMPONENT is determined by clinician input in the portal.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2236	<Picture> FREEFORM CONTENT LABEL COMPONENT (example) We've noticed a trend of 3 days with high pressure readings.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2237	8.12.1.15.3.6 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2238	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2059	8.12.1.16 DRX MEDICATION MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2060	8.12.1.16.1 DRX MEDICATION MESSAGE CARD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2063	The DRX MEDICATION MESSAGE CARD PANEL is conditionally displayed if a DRx Medication Message is received by the Patient Application at least once since the Earliest History Calendar Date. Patients with DTP enabled will be able to recieve a DRx Medication message. A unique DRX MEDICATION MESSAGE CARD PANEL is displayed for every DRx Medication Message that is received by the Patient Application. Multiple DRX MEDICATION MESSAGE CARD PANELS may be displayed in the same screen or panel. State: Active/Complete/Out of Date/Pre-acknowledged Additional Display Behavior: Active		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Only the most recent DRX MEDICATION MESSAGE CARD in the active state is displayed in the HOME SCREEN. The DRX MEDICATION MESSAGE CARD PANEL is in the active state by default. Complete The DRX MEDICATION MESSAGE CARD PANEL CARD PANEL changes state from active to complete when the user completes all required actions in the card. Out of Date The DRX MEDICATION MESSAGE CARD PANEL changes state from active to out of date if a more recent DRX MEDICATION MESSAGE CARD PANEL is received before the user is able to complete the required actions. Pre-acknowledged The DRX MEDICATION MESSAGE CARD PANEL changes state from active to pre-acknowledged if the clinician acknowledges the medication change on behalf of the patient (for example, when the clinician has discussed the medication change directly with the patient). The elements in the DRX MEDICATION MESSAGE CARD PANEL listed below change behavior along these three states.			
PAGUI2064	The DRX MEDICATION MESSAGE CARD PANEL is an element of: HOME SCREEN ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2065	The DRX MEDICATION MESSAGE CARD PANEL consists of: DRX MEDICATION ADJUSMENT BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) DRX MEDICATION TREND LABEL COMPONENT DRX MEDICATION CONTENT LABEL COMPONENT CMEMS PAP TREND DISPLAY PANEL SEE ADJUSTMENT ACTION TRIGGER COMPONENT (See Shared Elements: Components) MEDICATION UPDATE CARD PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2061	8.12.1.16.2 DRX MEDICATION MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2066	<Picture> DRX MEDICATION MESSAGE CARD PANEL (default: active, DTP ON, PA Visibility ON, Screen: Notifications)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 The image shows a medication adjustment message card. At the top is a yellow header with a pill icon and the text 'MEDICATION ADJUSTMENT'. Below this, the date 'FRIDAY, MAR 10' is displayed. The main text reads: 'You had 3 high CardioMEMS™ readings in a row. Your clinician recommends an adjustment.' Below the text is a line graph showing blood pressure readings over time. The y-axis ranges from 18 to 24 mmHg. The x-axis shows dates from JAN 12 to JAN 18. The graph shows a green line for readings below 24 mmHg and a yellow line for readings at or above 24 mmHg. The readings are: 19, 21, 23, 23, 25, 24, 25. Below the graph is a date range selector with 'JAN 12' and 'JAN 18' and four vertical tick marks. Below this is the text 'Optimal range: 18-26mmHg'. At the bottom is a blue button with the text 'See Adjustment'.			
PAGUI4443	<Picture> DRX MEDICATION MESSAGE CARD PANEL (default: active, DTP ON, PA Visibility ON, Screen: Home)  The image shows a medication adjustment message card. At the top is a yellow header with a pill icon and the text 'MEDICATION ADJUSTMENT'. Below this, the date 'FRIDAY, MAR 10' is displayed. The main text reads: 'You had 3 high CardioMEMS™ readings in a row. Your clinician recommends an adjustment.' Below the text is a blue button with the text 'See Adjustment'.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4363	<Picture> DRX MEDICATION MESSAGE CARD PANEL (default: active, DTP ON, PA Visibility OFF, Screen: Notifications, Home) (example: Return to Optimal)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
	 MEDICATION ADJUSTMENT FRIDAY, MAR 10 Based on your recent CardioMEMS™ readings, your clinician is prescribing a change to your medications. See Adjustment			
PAGUI2067	<Picture> DRX MEDICATION MESSAGE CARD PANEL (complete, example, Notifications Complete Screen)  MEDICATION ADJUSTMENT FRIDAY, MAR 10 Effective: 2022-06-24 Increase Clorthalidone Clorthalidone Trade name 1mg Tablet Once daily 2mg 1 Tablet Once daily		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2068	<Picture> DRX MEDICATION MESSAGE CARD PANEL (example, out of date, Notifications Complete Screen)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 MEDICATION ADJUSTMENT FRIDAY, MAR 10 Based on your recent CardioMEMS™ readings, your clinician is prescribing a change to your medications. See Adjustment			
PAGUI2062	8.12.1.16.3 DRX MEDICATION MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4386	8.12.1.16.3.1 DRX MEDICATION ADJUSTMENT BANNER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4444	Component: DRX MEDICATION ADJUSTMENT BANNER COMPONENT State: Active/ Complete/ Out of Date / Pre-acknowledged Display: Always displayed in the DRX MEDICATION MESSAGE CARD PANEL. The default state is active. In the active state, the DRX MEDICATION ADJUSTMENT BANNER COMPONENT displays the active icon and the active background color. In the complete state, the DRX MEDICATION ADJUSTMENT BANNER COMPONENT displays the complete icon and complete background color. In the out-of-date and pre-acknowledged states, the DRX MEDICATION ADJUSTMENT BANNER COMPONENT displays the out-of-date icon and the complete background color. Label: <label.gb.medication_adj>; Text is all capitalized. Behavior: Icon and text. The DRX MEDICATION ADJUSTMENT BANNER COMPONENT changes state with the message card panel it is an element of the DRX MEDICATION MESSAGE CARD PANEL.	PString264 Medication Adjustment notificationTitleMEDADJ label.gb.medication_adj	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4387	<Picture> DRX MEDICATION ADJUSTMENT BANNER COMPONENT (active) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4388	<Picture> DRX MEDICATION ADJUSTMENT BANNER (complete) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4442	<Picture> DRX MEDICATION ADJUSTMENT BANNER (example, out of date) (example, pre-acknowledged) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2140	8.12.1.16.3.4 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI2141	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4449	8.12.1.16.3.5 DRX MEDICATION TREND LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI4450	Component: DRX MEDICATION TREND LABEL COMPONENT State: Active/ Complete/ Out of Date / Pre-acknowledged Display: Conditionally displayed in the DRX MEDICATION MESSAGE CARD PANEL when PA Visibility is ON and DTP is ON at the same time, and only when the state is Active. Not displayed when the message from Merlin returns 0 days high/low (this happens when the patient returns to their optimal range and medications reset to daily titrations). Label: <notif.dtp.1_high>, <notif.dtp.many_high>, <notif.dtp.1_low>, <notif.dtp.many_low>, <notif.med_change> Behavior: Text.	PString1003 You had {26} high CardioMEMS™ reading. highPressureReadings notif.dtp.1_high PString1004 You had {26} high CardioMEMS™ readings in a row. highPressureReadingsInARow notif.dtp.many_high PString1005 You had {27} low CardioMEMS™ reading. lowPressureReadings notif.dtp.1_low PString1006 You had {27} low CardioMEMS™ readings in a row. lowPressureReadingsInARow notif.dtp.many_low PString1068 Your clinician recommends an adjustment. medicationAdjustmentBody1 notif.med_change	UI Info	CloudDRxRTR
PAGUI2147	<Picture> DRX MEDICATION TREND LABEL COMPONENT (example 1: Active, DTP ON, PA Visibility On, Screen: Notifications, Nonzero days outside of optimal range) (example 2: Active, DTP ON, PA Visibility On, Screen: Home, Nonzero days outside of optimal range)		UI Design	CloudDRx 1.0 CloudDRxRTR

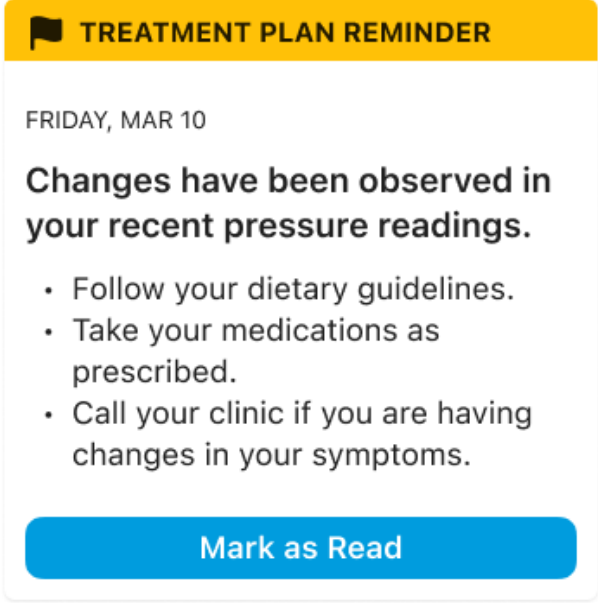


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	You had 3 high CardioMEMS™ readings in a row. Your clinician recommends an adjustment.			
PAGUI2145	8.12.1.16.3.7 DRX MEDICATION CONTENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2146	Component: DRX MEDICATION CONTENT LABEL COMPONENT State: Active/ Complete/ Out of Date / Pre-acknowledged Display: Conditionally displayed in the DRX MEDICATION MESSAGE CARD PANEL. DRX MEDICATION CONTENT LABEL COMPONENT is displayed in all Active states except the cases when PA Visibility is ON and DTP is ON and the message from Merlin returns a nonzero number of days high/low. Default state is active. Label: <notif.manual.med_change> Behavior: Text.	PString1008 Based on your recent CardioMEMS™ readings, your clinician is prescribing a change to your medications. medicationAdjustmentBody 1 notif.manual.med_change	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4366	<Picture> DRX MEDICATION CONTENT LABEL COMPONENT (default: Active, DTP ON, PA Visibility Off, Screen: Notifications) Based on your recent CardioMEMS™ readings, your clinician is prescribing a change to your medications.		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2163	8.12.1.16.3.8 CMEMS PAP TREND DISPLAY PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2164	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2149	8.12.1.16.3.9 SEE ADJUSTMENT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2150	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4470	8.12.1.16.3.10 MEDICATION UPDATE CARD PANEL		Heading	CloudDRxRTR
PAGUI4471	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI2265	8.12.1.17 ELEVATED PRESSURE MESSAGE CARD PANEL		Heading	CloudDRxRTR
PAGUI2266	8.12.1.17.1 ELEVATED PRESSURE MESSAGE CARD PANEL (Display)		Heading	CloudDRxRTR
PAGUI2267	The ELEVATED PRESSURE MESSAGE CARD PANEL is conditionally displayed if a Elevated Pressure Message is received by the Patient Application at least once since the Earliest History Calendar Date. A unique ELEVATED PRESSURE MESSAGE CARD		UI Info	CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	PANEL is displayed for every Elevated Pressure Message that is received by the Patient Application. Multiple ELEVATED PRESSURE MESSAGE CARD PANELS may be displayed in the same screen or panel. The ELEVATED PRESSURE MESSAGE CARD PANEL is in the active state by default. The ELEVATED PRESSURE MESSAGE CARD PANEL changes state from active to complete when the user completes all required actions in the card. The elements in the ELEVATED PRESSURE MESSAGE CARD PANEL listed below change behavior along these two states.			
PAGUI2268	The ELEVATED PRESSURE MESSAGE CARD PANEL is an element of: ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRxRTR
PAGUI2269	The ELEVATED PRESSURE MESSAGE CARD PANEL consists of: TREATMENT PLAN BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) PRESSURE CHANGE CONTENT LABEL COMPONENT PRESSURE CHANGE INSTRUCTIONS LABEL COMPONENT MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRxRTR
PAGUI2270	8.12.1.17.2 ELEVATED PRESSURE MESSAGE CARD PANEL (Layout)		Heading	CloudDRxRTR
PAGUI2271	<Picture> ELEVATED PRESSURE MESSAGE CARD PANEL (default: active, example) 		UI Design	CloudDRxRTR
PAGUI2272	<Picture> ELEVATED PRESSURE MESSAGE CARD PANEL (complete, example)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 TREATMENT PLAN REMINDER FRIDAY, MAR 10 Changes have been observed in your recent pressure readings. - Follow your dietary guidelines. - Take your medications as prescribed. - Call your clinic if you are having changes in your symptoms.			
PAGUI2274	8.12.1.17.3 ELEVATED PRESSURE MESSAGE CARD PANEL (Elements)		Heading	CloudDRxRTR
PAGUI2275	8.12.1.17.3.1 TREATMENT PLAN BANNER COMPONENT		Heading	CloudDRxRTR
PAGUI2298	Component: TREATMENT PLAN BANNER COMPONENT State: Active/ Complete Display: Always displayed in the ELEVATED PRESSURE MESSAGE CARD PANEL. The default state is active. In the active state, the TREATMENT PLAN BANNER COMPONENT displays the active icon and the active background color. In the complete state, the TREATMENT PLAN BANNER COMPONENT displays the complete icon and complete background color. Label: <label.gb.elevated_pressure>; Text is all capitalized. Behavior: Icon and text. The TREATMENT PLAN BANNER COMPONENT changes state with the notification panel (active/complete) it is an element of.	PString340 Treatment Plan Reminder HEP label.gb.elevated_pressure	UI Info	CloudDRxRTR
PAGUI2299	<Picture> TREATMENT PLAN BANNER COMPONENT (default: Active) 		UI Design	CloudDRxRTR
PAGUI2300	<Picture> TREATMENT PLAN BANNER COMPONENT (Complete) 		UI Design	CloudDRxRTR
PAGUI2279	8.12.1.17.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI2280	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI2281	8.12.1.17.3.3 PRESSURE CHANGE CONTENT LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI2282	Component: PRESSURE CHANGE CONTENT LABEL COMPONENT State: N/A Display: Always displayed in ELEVATED PRESSURE MESSAGE CARD PANEL. Label: <notif.standard.pressure_change> Behavior: Text.	PString572 Changes have been observed in your recent pressure readings. HEPBody1	UI Info	CloudDRxRTR

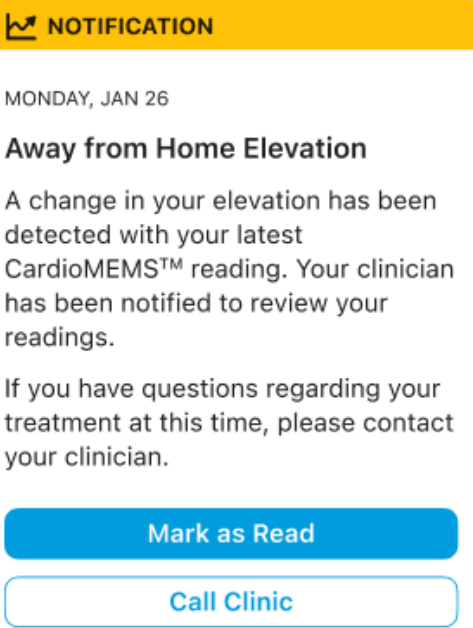
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		notif.standard.pressure_change		
PAGUI2283	<Picture> PRESSURE CHANGE CONTENT LABEL COMPONENT Changes have been observed in your recent pressure readings.		UI Design	CloudDRxRTR
PAGUI2284	8.12.1.17.3.4 PRESSURE CHANGE INSTRUCTIONS LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI2285	Component: PRESSURE CHANGE INSTRUCTIONS LABEL COMPONENT State: N/A Display: Always displayed in the ELEVATED PRESSURE MESSAGE CARD PANEL Label: <notif.standard.pressure_change_instruction> Behavior: Text.	PString746 •Follow your dietary guidelines. •Take your medications as prescribed. Call your clinic if you are having changes in your symptoms. HEPBody2 notif.standard.pressure_change_instruction	UI Info	CloudDRxRTR
PAGUI2286	<Picture> PRESSURE CHANGE INSTRUCTIONS LABEL COMPONENT (example) - Follow your dietary guidelines. - Take your medications as prescribed. - Call your clinic if you are having changes in your symptoms.		UI Design	CloudDRxRTR
PAGUI2288	8.12.1.17.3.5 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI2289	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI3978	8.12.1.18 ELEVATION CHANGE MESSAGE CARD PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI3979	8.12.1.18.1 ELEVATION CHANGE MESSAGE CARD PANEL (Display)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3982	The ELEVATION CHANGE MESSAGE CARD PANEL is conditionally displayed if an Elevation Change Message is received by the Patient Application at least once since the Earliest History Calendar Date. A unique ELEVATION CHANGE MESSAGE CARD PANEL is dispalyed displayed for every Elevation Change Message that is received by the Patient Application. The ELEVATION CHANGE MESSAGE CARD PANEL is displayed in the active state in the HOME SCREEN and the ACTIVE NOTIFICATIONS PANEL. If there is at least one ELEVATION CHANGE MESSAGE CARD PANEL, the user is not able to navigate to the MEDICATIONS SCREEN in any way. The ELEVATION CHANGE MESSAGE CARD PANEL is displayed in the complete state in the COMPLETE NOTIFICATIONS PANEL. The ACTIVE NOTIFICATIONS PANEL changes state from active to complete when the user selects the MARK AS READ ACTION TRIGGER COMPONENT.		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3983	The ELEVATION CHANGE MESSAGE CARD PANEL is an element of:		UI Info	CloudDRxRTR

Commented [RC4]: Example: elevation

Commented [RC5]: Fix typo

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	HOME SCREEN ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL			CloudDRxPS
PAGUI3984	The ELEVATION CHANGE MESSAGE CARD PANEL consists of: NOTIFICATION BANNER COMPONENT (See Shared Elements: Components) NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) ELEVATION CHANGE HEADER COMPONENT ELEVATION CHANGE CONTENT COMPONENT ELEVATION CHANGE INSTRUCTIONS COMPONENT MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components) CALL CLINIC ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3980	8.12.1.18.2 ELEVATION CHANGE MESSAGE CARD PANEL (Layout)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3985	<Picture> ELEVATION CHANGE MESSAGE CARD PANEL (default: active) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3986	<Picture> ELEVATION CHANGE MESSAGE CARD PANEL (complete)		UI Design	CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	✔ NOTIFICATION MONDAY, JAN 26 Away from Home Elevation A change in your elevation has been detected with your latest CardioMEMS™ reading. Your clinician has been notified to review your readings. If you have questions regarding your treatment at this time, please contact your clinician.			
PAGUI3981	8.12.1.18.3 ELEVATION CHANGE MESSAGE CARD PANEL (Elements)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3987	8.12.1.18.3.1 NOTIFICATION BANNER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3988	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI3989	8.12.1.18.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3990	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI3991	8.12.1.18.3.3 ELEVATION CHANGE HEADER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3992	Component: ELEVATION CHANGE HEADER COMPONENT State: N/A Display: Always displayed as a component of the ELEVATION CHANGE MESSAGE CARD PANEL. Label: <notif.all.elevation_change_header> Behavior: Text.	PString793 Out of Home Elevation altitudeChangeHeader notif.all.altitude_change_header PString1087 Away from Home Elevation elevationChangeHeader notif.all.elevation_change_header	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4851	<Picture> ELEVATION CHANGE HEADER COMPONENT Away from Home Elevation		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

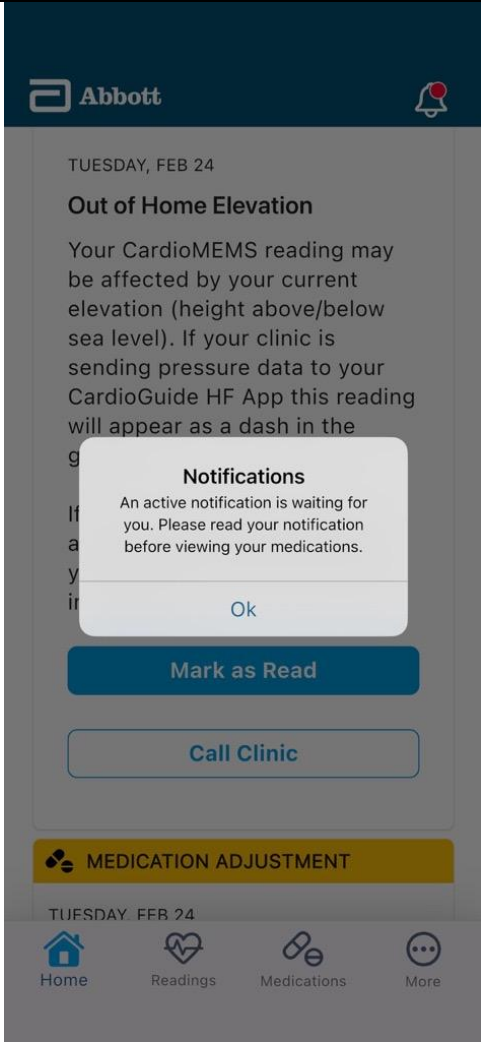
ID		English US Strings	Object Type	Project
PAGUI3994	8.12.1.18.3.4 ELEVATION CHANGE CONTENT COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3995	Component: ELEVATION CHANGE CONTENT COMPONENT State: N/A Display: Always displayed as a component of the ELEVATION CHANGE MESSAGE CARD PANEL. Label: <notif.altitude_change_content> Behavior: Text.	PString794 Your CardioMEMS™ reading may be affected by your current elevation (height above/below sea level). If your clinic is sending pressure data to your CardioGuide HF App this reading will appear as a dash in the graph. altitudeChangeActiveContent notif.all.altitude_change_active_content PString1088 A change in your elevation has been detected with your latest CardioMEMS™ reading. Your clinician has been notified to review your readings. elevationChangeContent notif.altitude_change_content	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4852	<Picture> ELEVATION CHANGE CONTENT COMPONENT A change in your elevation has been detected with your latest CardioMEMS™ reading. Your clinician has been notified to review your readings.			CloudDRxRTR
PAGUI3998	8.12.1.18.3.5 ELEVATION CHANGE INSTRUCTIONS COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3999	Component: ELEVATION CHANGE INSTRUCTIONS COMPONENT State: active/complete Display: Always displayed as a component of the ELEVATION CHANGE MESSAGE CARD PANEL when the ELEVATION CHANGE MESSAGE CARD PANEL. Label: <notif.all.altitude_elevation_instructions> Behavior: Text.	PString796 If you receive a medication adjustment in the app, contact your clinic to confirm treatment instructions.	UI Info	CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		altitudeChangeInstructions notif.all.altitude_change_instructions PString1089 If you have questions regarding your treatment at this time, please contact your clinician. elevationInstructions notif.all.altitude_elevation_instructions		
PAGUI4853	<Picture> ELEVATION CHANGE INSTRUCTIONS COMPONENT (default: active) If you have questions regarding your treatment at this time, please contact your clinician.			CloudDRxRTR
PAGUI4001	8.12.1.18.3.6 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4002	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI4003	8.12.1.18.3.7 CALL CLINIC ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4004	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI4204	8.12.1.19 ELEVATION NOTIFICATION DIALOG PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI4205	8.12.1.19.1 ELEVATION NOTIFICATION DIALOG PANEL (Display)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4206	The ELEVATION NOTIFICATION DIALOG PANEL is displayed after the user selects the MEDICATIONS ACTION TRIGGER COMPONENT when there is an active ELEVATION CHANGE MESSAGE CARD. The ELEVATION NOTIFICATION DIALOG PANEL is displayed over the screen, does not displace any elements, and any part of the screen that is not covered by the ELEVATION DIALOG PANEL is greyed out.			CloudDRxRTR CloudDRxPS
PAGUI4207	The ELEVATION NOTIFICATION DIALOG PANEL is composed of: ELEVATION NOTIFICATION LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4208	8.12.1.19.2 ELEVATION NOTIFICATION DIALOG PANEL (Layout)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4209	<Picture> ELEVATION NOTIFICATION DIALOG PANEL		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

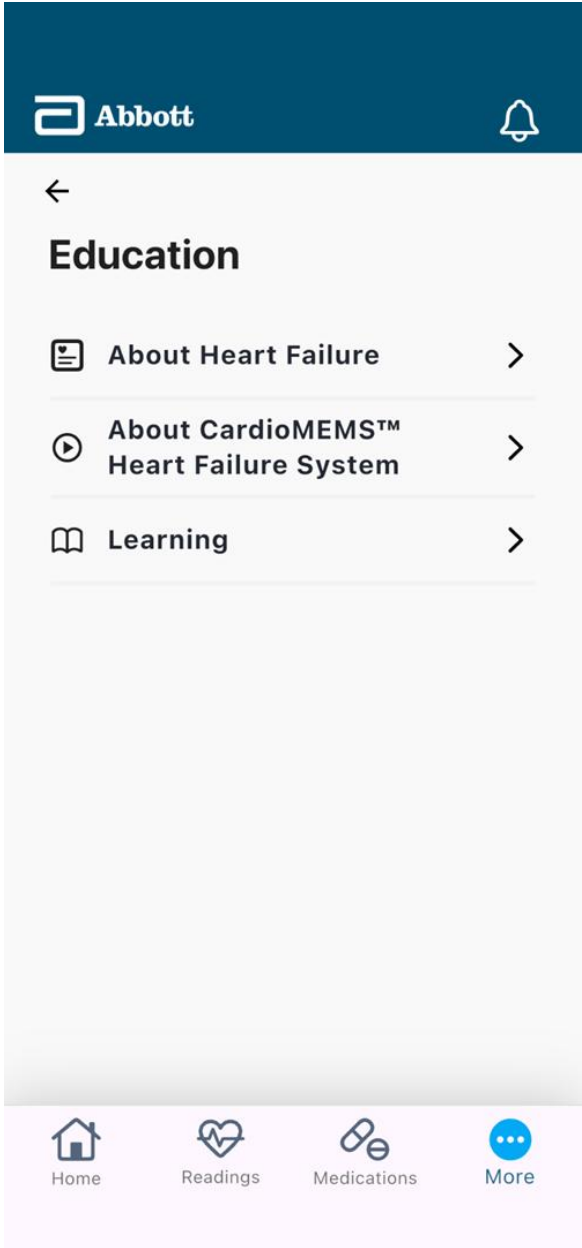
ID		English US Strings	Object Type	Project
				
PAGUI4210	8.12.1.19.3 ELEVATION NOTIFICATION DIALOG PANEL (Elements)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4211	8.12.1.19.3.1 ELEVATION NOTIFICATION LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4212	Component: ELEVATION NOTIFICATION LABEL COMPONENT State: N/A Display: Always displayed as a component of the ELEVATION NOTIFICATION DIALOG PANEL. Label: <label.gb.notifications> <msg.gb.elevation_change> Behavior: Text.	PString191 Notifications notifications label.gb.notifications PString1025	UI Info	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		An active notification is waiting for you. Please read your notification before viewing your medications. msg.gb.elevation_change		
PAGUI4213	<Picture> ELEVATION NOTIFICATION LABEL COMPONENT 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4214	8.12.1.19.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4215	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI579	8.12.1.20 EDUCATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI860	8.12.1.20.1 EDUCATION PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI863	The EDUCATION PANEL is displayed after the user navigates to it from the MORE SCREEN. The EDUCATION PANEL is displayed over the entire screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI864	The EDUCATION PANEL is composed of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) EDUCATION HEADER COMPONENT ABOUT HEART FAILURE ACTION TRIGGER COMPONENT ABOUT HEART FAILURE PANEL ABOUT HEART FAILURE SYSTEM ACTION TRIGGER COMPONENT ABOUT HEART FAILURE SYSTEM PANEL LEARNING ACTION TRIGGER COMPONENT LEARNING PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI861	8.12.1.20.2 EDUCATION PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI865	<Picture> EDUCATION PANEL 		UI Design	CloudDRxPS CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI862	8.12.1.20.3 EDUCATION PANEL (Elements)		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1725	8.12.1.20.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1726	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4594	8.12.1.20.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4595	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI866	8.12.1.20.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI867	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI870	8.12.1.20.3.4 EDUCATION HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI871	Component: EDUCATION HEADER COMPONENT State: N/A Display: Always displayed in the EDUCATION PANEL. Label: <label.gb.education> Behavior: Text.	PString574 Education education label.gb.education	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI872	<Picture> EDUCATION HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI873	8.12.1.20.3.5 ABOUT HEART FAILURE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI874	Component: ABOUT HEART FAILURE ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the EDUCATION PANEL. Label: <label.gb.about_hf> Behavior: Icon, text, and trigger. Upon activation, the ABOUT HEART FAILURE ACTION TRIGGER COMPONENT displays the ABOUT HEART FAILURE PANEL.	PString575 About Heart Failure aboutHeartFailureTitle label.gb.about_hf	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI875	<Picture> ABOUT HEART FAILURE ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI876	8.12.1.20.3.6 ABOUT HEART FAILURE PANEL		Heading	CloudDRx 1.0

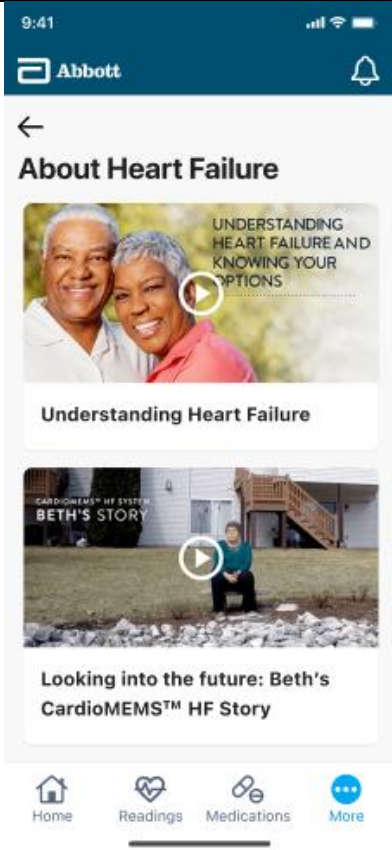
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

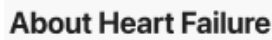
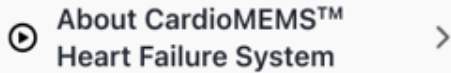
ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI877	8.12.1.20.3.6.1 ABOUT HEART FAILURE PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI880	The ABOUT HEART FAILURE PANEL is conditionally displayed in the EDUCATION PANEL when the user selects the ABOUT HEART FAILURE ACTION TRIGGER COMPONENT.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI881	The ABOUT HEART FAILURE PANEL is composed of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) ABOUT HEART FAILURE HEADER COMPONENT VIDEO ARTICLE ACTION TRIGGER COMPONENT (See Shared Elements: Components) VIDEO PLAYER PANEL (See Shared Elements: Panels)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI878	8.12.1.20.3.6.2 ABOUT HEART FAILURE PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI882	<Picture> ABOUT HEART FAILURE PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI879	8.12.1.20.3.6.3 ABOUT HEART FAILURE PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1727	8.12.1.20.3.6.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1728	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4596	8.12.1.20.3.6.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4597	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI883	8.12.1.20.3.6.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI884	See Shared Elements: Components.		Info	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

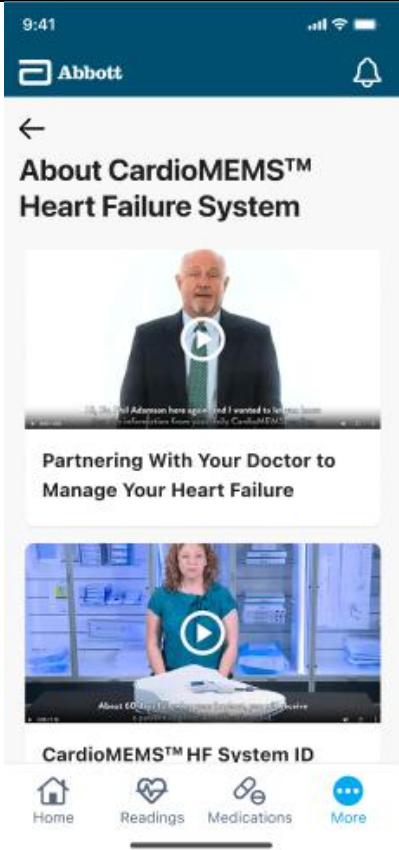
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI885	8.12.1.20.3.6.3.4 ABOUT HEART FAILURE HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI886	Component: ABOUT HEART FAILURE HEADER COMPONENT State: N/A Display: Always displayed in the ABOUT HEART FAILURE PANEL. Label: <label.gb.about_hf> Behavior: Text.	PString575 About Heart Failure aboutHeartFailureTitle label.gb.about_hf	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI887	<Picture> ABOUT HEART FAILURE HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI888	8.12.1.20.3.6.3.5 VIDEO ARTICLE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI889	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1733	8.12.1.20.3.6.3.6 VIDEO PLAYER PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1734	See Shared Elements: Panel.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1749	8.12.1.20.3.7 ABOUT HEART FAILURE SYSTEM ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1750	Component: ABOUT HEART FAILURE SYSTEM ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the EDUCATION PANEL. Label: label.gb.about_hf_system Behavior: Icon, text, and trigger. Upon activation, the ABOUT HEART FAILURE SYSTEM ACTION TRIGGER COMPONENT displays the ABOUT HEART FAILURE SYSTEM PANEL.	PString576 About {21} aboutCardioMEMSTMHeart FailureSystemTitle label.gb.about_hf_system	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1751	<Picture> ABOUT HEART FAILURE SYSTEM ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1752	8.12.1.20.3.8 ABOUT HEART FAILURE SYSTEM PANEL		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

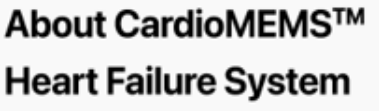
myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1753	8.12.1.20.3.8.1 ABOUT HEART FAILURE SYSTEM PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1754	The ABOUT HEART FAILURE SYSTEM PANEL is conditionally displayed in the EDUCATION PANEL when the user selects the ABOUT HEART FAILURE SYSTEM ACTION TRIGGER COMPONENT.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1755	The ABOUT HEART FAILURE SYSTEM PANEL is composed of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) ABOUT HEART FAILURE SYSTEM HEADER COMPONENT VIDEO ARTICLE ACTION TRIGGER COMPONENT (See Shared Elements: Components) VIDEO PLAYER PANEL (See Shared Elements: Panels)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1756	8.12.1.20.3.8.2 ABOUT HEART FAILURE SYSTEM PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1757	<Picture> ABOUT HEART FAILURE SYSTEM PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1758	8.12.1.20.3.8.3 ABOUT HEART FAILURE SYSTEM PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1759	8.12.1.20.3.8.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1760	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4598	8.12.1.20.3.8.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4599	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI1761	8.12.1.20.3.8.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1762	See Shared Elements: Components.		Info	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1763	8.12.1.20.3.8.3.4 ABOUT HEART FAILURE SYSTEM HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1764	Component: ABOUT HEART FAILURE SYSTEM HEADER COMPONENT State: N/A Display: Always displayed in the ABOUT HEART FAILURE SYSTEM PANEL. Label: <label.gb.about_hf_system> Behavior: Text.	PString576 About {21} aboutCardioMEMSTMHeart FailureSystemTitle label.gb.about_hf_system	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1765	<Picture> ABOUT HEART FAILURE SYSTEM HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1766	8.12.1.20.3.8.3.5 VIDEO ARTICLE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1767	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1768	8.12.1.20.3.8.3.6 VIDEO PLAYER PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1769	See Shared Elements: Panel.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI891	8.12.1.20.3.9 LEARNING ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI892	Component: LEARNING ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the EDUCATION PANEL. Label: <label.gb.learning> Behavior: Icon, text, and trigger. Upon activation, the LEARNING ACTION TRIGGER COPONENT displays the LEARNING PANEL.	PString199 Learning learning label.gb.learning	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI893	<Picture> LEARNING ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI894	8.12.1.20.3.10 LEARNING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI895	8.12.1.20.3.10.1 LEARNING PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI898	The LEARNING PANEL is conditionally displayed in the EDUCATION PANEL when the user selects the LEARNING ACTION TRIGGER COMPONENT.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI899	The LEARNING PANEL is composed of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) LEARNING HEADER COMPONENT IN-APP TRAINING ACTION TRIGGER COMPONENT IN-APP TRAINING PANEL FAQ TOPIC PANEL (See Shared Elements: Panels)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI896	8.12.1.20.3.10.2 LEARNING PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI900	<Picture> LEARNING PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

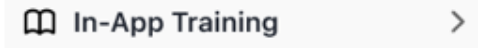
myCardioMEMS Patient Application Graphical User Interface Specification

	9:41    Learning  In-App Training  About the CardioGuide™ HF App What does the CardioGuide HF App do?  How does the CardioGuide HF App communicate with the clinic?  Who can use the CardioGuide HF App?  Getting Started LANGUAGE How do I select the language of the CardioGuide HF App?  STAYING CONNECTED How can I stay connected to my clinic?  Do I always need to be connected to the Internet?   Home  Readings  Medications  More				
--	--	--	--	--	--

entia document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



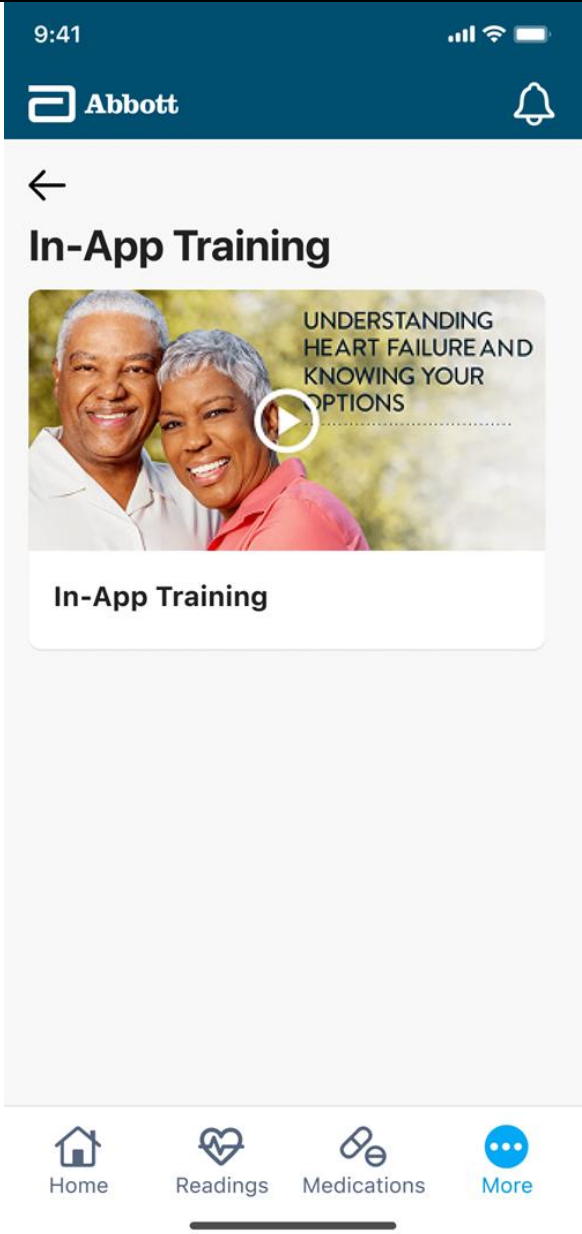
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI897	8.12.1.20.3.10.3 LEARNING PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1770	8.12.1.20.3.10.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1771	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4712	8.12.1.20.3.10.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4713	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI901	8.12.1.20.3.10.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI902	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI903	8.12.1.20.3.10.3.4 LEARNING HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI904	Component: LEARNING HEADER COMPONENT State: N/A Display: Always displayed in the LEARNING PANEL. Label: <label.gb.learning> Behavior: Text.	PString199 Learning learning label.gb.learning	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI905	<Picture> LEARNING HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI909	8.12.1.20.3.10.3.5 IN-APP TRAINING ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI910	Component: IN-APP TRAINING ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the LEARNING PANEL. Label: Behavior: Icon, text, and trigger. Upon activation, the IN-APP TRAINING ACTION TRIGGER COMPONENT displays the IN-APP TRAINING PANEL.	PString577 In-App Training inAppTrainingTitle label.gb.in_app_training	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI911	<Picture> IN-APP TRAINING ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3004	8.12.1.20.3.10.3.6 IN-APP TRAINING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3005	8.12.1.20.3.10.3.6.1 IN-APP TRAINING PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3008	The IN-APP TRAINING PANEL is displayed as an element of the LEARNING PANEL after the user selects the IN-APP TRAINING ACTION TRIGGER COMPONENT.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3432	The IN-APP TRAINING PANEL consists of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) IN-APP TRAINING HEADER COMPONENT VIDEO ARTICLE ACTION TRIGGER COMPONENT (See Shared Elements: Components) VIDEO PLAYER PANEL (See Shared Elements: Panels)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3006	8.12.1.20.3.10.3.6.2 IN-APP TRAINING PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3433	<Picture> IN-APP TRAINING PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
				
PAGUI3007	8.12.1.20.3.10.3.6.3 IN-APP TRAINING PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

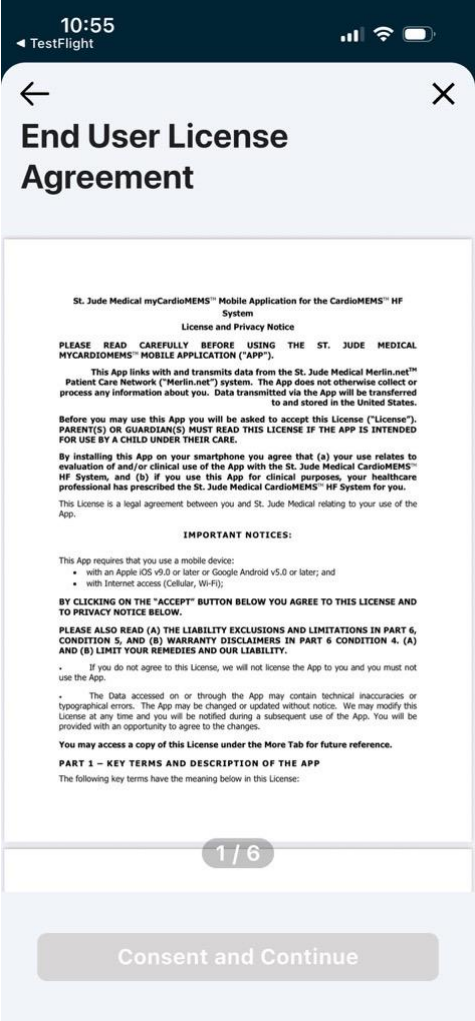
ID		English US Strings	Object Type	Project
PAGUI3011	8.12.1.20.3.10.3.6.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3012	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4600	8.12.1.20.3.10.3.6.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4601	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI3013			Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3014	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3015	8.12.1.20.3.10.3.6.3.4 IN-APP TRAINING HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3016	Component: IN-APP TRAINING HEADER COMPONENT State: N/A Display: Always displayed in the IN-APP TRAINING PANEL. Label: <label.gb.in_app_training> Behavior: Text.	PString577 In-App Training inAppTrainingTitle label.gb.in_app_training	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3017	<Picture> IN-APP TRAINING HEADER COMPONENT In-App Training		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3434	8.12.1.20.3.10.3.6.3.8 VIDEO ARTICLE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3435	See Shared Elements: Components		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3027	8.12.1.20.3.10.3.6.3.9 VIDEO PLAYER PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3028	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1772	8.12.1.20.3.10.3.7 FAQ TOPIC PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



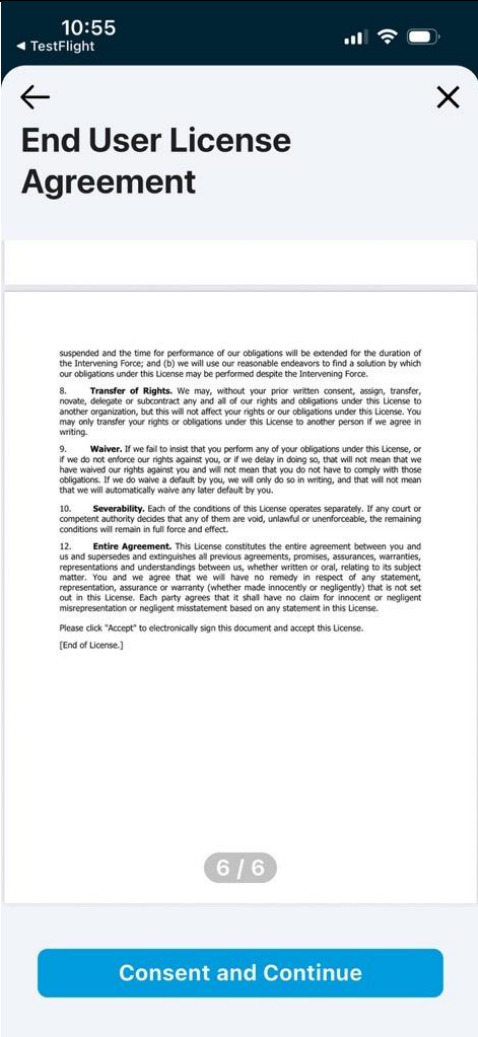
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1773	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1851	8.12.1.21 END USER LICENSE AGREEMENT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1852	8.12.1.21.1 END USER LICENSE AGREEMENT PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1853	The END USER LICENSE AGREEMENT PANEL is displayed after the user selects the POLICY WELCOME CONTINUE ACTION TRIGGER COMPONENT on the POLICY WELCOME PANEL or the END USER LICENSE AGREEMENT ACTION TRIGGER COMPONENT on the ABOUT PANEL. The END USER LICENSE AGREEMENT PANEL is conditionally displayed after the user selects the LOGIN ACTION TRIGGER COMPONENT if there is an update to the End User License Agreement. The END USER LICENSE AGREEMENT PANEL is displayed over the entire screen except the OS STATUS PANEL. In the ABOUT PANEL, the GLOBAL HEADING PANEL and DONE ACTION TRIGGER COMPONENT is displayed. In the POLICY SCREEN, the GLOBAL HEADING PANEL is not displayed and the CONSENT ACTION TRIGGER COMPONENT is displayed in lieu of the DONE ACTION TRIGGER COMPONENT.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1890	The END USER LICENSE AGREEMENT PANEL is an element of: POLICY SCREEN ABOUT PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4010	The END USER LICENSE AGREEMENT PANEL consists of: OS STATUS PANEL (Shared Elements: Panels) GLOBAL HEADING PANEL (Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (Shared Elements: Components) CLOSE PANEL ACTION TRIGGER COMPONENT (Shared Elements: Components) END USER LICENSE AGREEMENT HEADER COMPONENT END USER LICENSE AGREEMENT LABEL COMPONENT CONSENT ACTION TRIGGER COMPONENT DONE ACTION TRIGGER COMPONENT (Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI1855	8.12.1.21.2 END USER LICENSE AGREEMENT PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1856	<Picture> END USER LICENSE AGREEMENT PANEL (example: from POLICY SCREEN, top of PDF)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

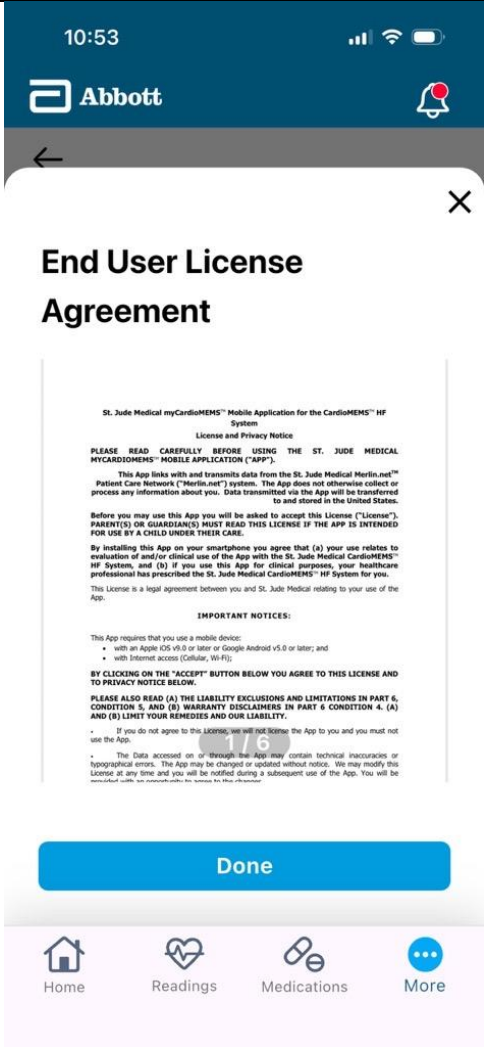
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4684	<Picture> END USER LICENSE AGREEMENT PANEL (example: from POLICY SCREEN, bottom of PDF)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

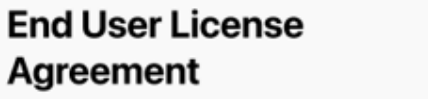
ID		English US Strings	Object Type	Project
				
PAGUI1888	<Picture> END USER LICENSE AGREEMENT PANEL (example: from ABOUT PANEL)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1857	8.12.1.21.3 END USER LICENSE AGREEMENT PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4634	8.12.1.21.3.1 OS STATUS PANEL		Heading	CloudDRxRTR
PAGUI4635	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI4685	8.12.1.21.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4686	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI1858	8.12.1.21.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1859	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4636	8.12.1.21.3.4 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4637	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI1860	8.12.1.21.3.5 END USER LICENSE AGREEMENT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
	Component: END USER LICENSE AGREEMENT HEADER COMPONENT State: N/A Display: Always displayed in the END USER LICENSE AGREEMENT PANEL. Label: <label.gb.EULA> Behavior: Text.	PString581 End User License Agreement termsOfUseTitle label.gb.EULA	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1862	<Picture> END USER LICENSE AGREEMENT HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1863	8.12.1.21.3.6 END USER LICENSE AGREEMENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1864	Component: END USER LICENSE AGREEMENT LABEL COMPONENT State: N/A Display: Always displayed in the END USER LICENSE AGREEMENT PANEL. Label: <other.pdf.EULA> - the Patient Application displays a PDF document located with the associated website link. Behavior: Text. The content for the END USER LICENSE AGREEMENT LABEL COMPONENT is managed through content management in Merlin.net to enable updates as needed. The PDF document can be multiple pages with the numbr of pages displayed on the bottom and allows scrolling.	PString702 Display a PDF of the applicable document. pdfTermsOfUse other.pdf.EULA	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1865	<Picture> END USER LICENSE AGREEMENT LABEL COMPONENT (example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Lorem ipsum dolor sit amet, consectetur adipiscing elit. Sed enim nunc, euismod vel porta tincidunt, condimentum eget risus. Nulla sodales, leo sit amet sagittis lobortis, sem risus rutrum nulla, sodales iaculis mi tellus pellentesque nulla. Curabitur sit amet felis et orci scelerisque semper. Pellentesque in dui viverra, cursus urna non, fermentum eros. Sed sed lectus vestibulum, semper augue non, auctor lorem. Aliquam ex nisl, finibus non aliquam sit amet, scelerisque venenatis dui. Quisque at sapien accumsan, finibus dolor vel, gravida magna. Nulla eleifend nibh ligula, quis sollicitudin magna dignissim at. Nunc imperdiet nunc vitae quam volutpat, at rhoncus eros dignissim. Ut convallis dui tortor, quis luctus erat egestas sed. Donec et scelerisque dolor, et posuere neque. Sed vitae justo ultrices, lobortis ipsum quis, interdum tortor.			
PAGUI4016	8.12.1.21.3.8 CONSENT ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4017	Component: CONSENT ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Displayed in the END USER LICENSE AGREEMENT PANEL when in the POLICY SCREEN. Default state is disabled. Label: <label.gb.consent> Behavior: Text and trigger. The CONSENT ACTION TRIGGER COMPONENT is enabled when the user scrolls to the bottom of the END USER LICENSE AGREEMENT LABEL COMPONENT. Upon activation of the CONSENT ACTION TRIGGER COMPONENT in the enabled state the SHORT PRIVACY PANEL is displayed.	PString611 Consent and Continue consentAndContinue label.gb.consent	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4018	<Picture> CONSENT ACTION TRIGGER COMPONENT (default: disabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4019	<Picture> CONSENT ACTION TRIGGER COMPONENT (enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4843	8.12.1.21.3.9 DONE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4844	See Shared Elements: Components.		Info	CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1774	8.12.1.22 FAQ TOPIC PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1775	8.12.1.22.1 FAQ TOPIC PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1778	The FAQ TOPIC PANEL is conditionally displayed when there is a general topic in the content mangement location. The number of FAQ TOPIC PANELS displayed is determined by how many unique general topics there are stored in the content management location.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1779	The FAQ TOPIC PANEL is an element of: LEARNING PANEL FAQ PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1785	The FAQ TOPIC PANELS include: About myCardioMEMS App What does it do? How does it communicate? Who can use it? Getting Started What does it do? How does it communicate with clinic? Who can use it? App Features Passwords		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4084	The FAQ TOPIC PANEL that describes general information about the Patient Application is displayed in the LEARNING PANEL and the FAQ PANEL. The order of questions and answers is listed below. Header: <label.about_header> Question 1: <label.about_q1> Answer 1: <label.about_a1> Question 2: <label.about_q2> Answer 2: <label.about_a2> Question 3: <label.about_q3> Answer 3: <label.about_a3> Question 4: <label.about_q4> Answer 4: <label.about_a4>	PString806 About the {0} faqCategoryAboutTheCardi oApp label.about_header PString807 What does the {0} do? faqCategoryAboutTheCardi oAppFaq1_question label.about_q1 PString809 The {0} is part of a remote care system used to monitor a patient's heart failure. The app is used to view and track PA pressure sensor readings, receive messages from the clinics,	UI Info	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		provide medication reminders, and record and track Systemic Blood Pressure and Weight. {n}{n}Note: PA pressure visibility and medication reminders are only available if your clinician provides permission to access the information via the clinician's software platform. faqCategoryAboutTheCardioAppFaq1_answer label.about_a1 PString811 How does the {0} communicate with the clinic? faqCategoryAboutTheCardioAppFaq2_question label.about_q2 PString814 The {0} uses Wi-Fi or cellular data to receive information from the clinician's software platform. Information will be sent automatically by the clinician's software platform or when the clinician uses the platform to send a message. The {0} must be running in the background at all times. Do not force close the {0}. If the {0} is force closed, it will not be able to perform its activities as programmed by the clinic. Ensure		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		notifications are enabled for the {0} within your phone settings to know when your clinician has sent you a message. faqCategoryAboutTheCardioAppFaq2_answer label.about_a2 PString815 Who can use the {0}? faqCategoryAboutTheCardioAppFaq3_question label.about_q3 PString817 As part of the {21}, the {0} is intended to allow patients/caregivers to input health metric information and communicate with clinicians regarding their heart failure care. Depending on clinician-defined settings, the {0} may also allow users to view {21} data and receive clinician-directed automated treatment plan instructions. faqCategoryAboutTheCardioAppFaq3_answer label.about_a3 PString818 What is {34} (PCN) Heart Failure (HF) Portal? TBD label.about_q4 PString819 {34} (PCN) Heart Failure (HF) Portal is the software that your clinicians use to monitor your {21} readings.		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		TBD label.about_a4		
PAGUI4085	The FAQ TOPIC PANEL that describes basic functions the user may need to start using the Patient Application is displayed in the LEARNING PANEL and the FAQ PANEL. The display order is described below. Header: <label.start_header> Sub-Header: <label.start_language> Question 1: <label.start_language_q1> Answer 1: <label.start_language_a1> Sub-Header: <label.start_connected> Question 1: <label.start_connected_q1> Answer 1: <label.start_connected_a1> Question 2: <label.start_connected_q2> Answer 2: <label.start_connected_a2> Question 3: <label.start_connected_q3> Answer 3: <label.start_connected_a3> Question 4: <label.start_connected_q4> Answer 4: <label.start_connected_a4> Sub-Header: <label.start_use> Question 1: <label.start_use_q1> Answer 1: <label.start_use_a1> Question 2: <label.start_use_q2> Answer 2: <label.start_use_a2> Question 3: <label.start_use_q3> Answer 3: <label.start_use_a3> Question 4: <label.start_use_q4> Answer 4: <label.start_use_a4> Question 5: <label.start_use_q5> Answer 5: <label.start_use_a5> Question 6: <label.start_use_q6> Answer 6: <label.start_use_a6>	PString821 Getting Started faqCategoryGettingStarted label.start_header PString822 Language language label.start_language PString823 How do I select the language of the {0}? faqCategoryLanguageFaq1_question label.start_language_q1 PString824 The {0} supports multiple languages. The {0} uses the language of your mobile device to determine which language to display. If the language of your mobile device is supported by the {0}, the app will use this language when displaying information. If the language of your mobile device is not supported, the {0} will use English as the default language. faqCategoryLanguageFaq1_answer label.start_language_a1 PString825 Staying Connected stayingConnected label.start_connected PString826 How can I stay connected to my clinic?	UI Info	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		faqCategoryStayingConnectedFaq1_question label.start_connected_q1 PString828 The {0} must be running in the background to be able to complete its functions (tasks). Please make sure that your mobile device is connected to the Internet at all times. In addition, make sure that notifications on your phone are turned on. Do not force close the {0}. If the app is force closed, it will not be able to perform its activities as programmed by the clinic. {n}If you would like to contact your clinic, you can find your clinic's contact number in the Menu. Otherwise, you may continue to contact your clinic with your standard method of communication. faqCategoryStayingConnectedFaq1_answer label.start_connected_a1 PString829 Do I always need to be connected to the Internet? faqCategoryStayingConnectedFaq2_question label.start_connected_q2 PString831 The {0} uses Wi-Fi or cellular data to receive data from your clinic. Thus, your mobile device should		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		always be connected to the Internet to make sure that messages from your clinic are received in a timely manner. If your device is not connected to the internet, you can open the {0} with limited functions. Any information sent by your clinic will be received in the {0} when you are able to reconnect with Wi-Fi or cellular data. faqCategoryStayingConnectedFaq2_answer label.start_connected_a2 PString832 How do I turn on Wi-Fi? faqCategoryStayingConnectedFaq3_question label.start_connected_q3 PString833 You can turn on and set up Wi-Fi from your mobile device settings. For further help, you can consult your mobile device manual. faqCategoryStayingConnectedFaq3_answer label.start_connected_a3 PString834 How do I know if I am connected to a cellular data network? faqCategoryStayingConnectedFaq4_question label.start_connected_q4 PString835		

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		Your cellular Internet connection is indicated on the top of your mobile device screen. Please consult your mobile device manual for details. The {0} will also indicate when network connectivity is lost. faqCategoryStayingConnectedFaq4_answer label.start_connected_a4 PString836 Using the {0} usingTheMyCardioMEMApp label.start_use PString837 Can I use other apps on my mobile device while using the {0}? faqCategoryUsingTheCardioAppFaq1_question label.start_use_q1 PString838 Yes, but if you are using an Apple mobile digital device do not kill or force quit the {0}. The {0} must be running in the background to be able to complete its functions (tasks). faqCategoryUsingTheCardioAppFaq1_answer label.start_use_a1 PString839 Can I terminate the {0} after it has finished its task to save my phone battery? faqCategoryUsingTheCardioAppFaq2_question label.start_use_q2 PString840		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		If you have an Apple device, terminating the app will prevent the {0} from receiving notifications. If you turn off your phone, please remember to open the {0} when you power your phone back on and leave it running in the background. faqCategoryUsingTheCardio AppFaq2_answer labe.start_use_a2 PString841 How does the {0} use notifications? faqCategoryUsingTheCardio AppFaq3_question label.start_use_q3 PString843 The {0} uses notifications to let you know if there are new messages from the clinic. {n}Messages within the notification screen in the {0} have been sent from your clinic and may provide recommendations for your ongoing treatment. {n}Please ensure the notifications on your mobile device are enabled for the {0}. faqCategoryUsingTheCardio AppFaq3_answer label.start_use_a3 PString844		

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		I usually turn off notifications for all my apps. Will that affect the {0}'s ability to receive messages and data? faqCategoryUsingTheCardio AppFaq4_question label.start_use_q4 PString846 It is recommended that you leave notifications on so that the app can function optimally. The {0} uses notifications to let you know if there are new messages from the clinic. If you turn off notifications, you may not realize that there are new messages from the clinic until you reopen the app. faqCategoryUsingTheCardio AppFaq4_answer label.start_use_a4 PString847 How do I turn on notifications? faqCategoryUsingTheCardio AppFaq5_question label.start_use_q5 PString848 If you opted out of notifications during the app setup you can turn them back on from your mobile device settings. You can also consult your mobile device manual. faqCategoryUsingTheCardio AppFaq5_answer label.start_use_a5		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		PString849 Does the {0} collect my location information? Does the {0} track me? faqCategoryUsingTheCardioAppFaq6_question label.start_use_q6 PString850 No, {22} and the {0} do not track your location or keep your location information. faqCategoryUsingTheCardioAppFaq6_answer label.start_use_a6		
PAGUI4086	The FAQ TOPIC PANEL that describes key features of the Patient Application is displayed in the LEARNING PANEL and the FAQ PANEL. The display order is described below. Header: <label.feature_header> Sub-Header: <label.feature_reading> Question 1: <label.feature_reading_q1> Answer 1: <label.feature_reading_a1> Question 2: <label.feature_reading_q2> Answer 2: <label.feature_reading_a2> Sub-Header: <label.feature_history> Question 1: <label.feature_history_q1> Answer 1: <label.feature_history_a1> Question 2: <label.feature_history_q2> Answer 2: <label.feature_hisotry_a2> Question 3: <label.feature_history_q3> Answer 3: <label.feature_history_a3> Question 4: <label.feature_history_q4> Answer 4: <label.feature_history_a4> Sub-Header: <label.feature_message> Question 1: <label.feature_messages_q1> Answer 1: <label.feature_messages_a1> Question 2: <label.feature_messages_q2> Answer 2: <label.feature_messages_a2> Question 3: <label.feature_messages_q3> Answer 3: <label.feature_messages_a3> Sub-Header: <label.feature_med_messages> Information: <label.feature_med_messages_info>	PString852 App Features faqCategoryAppFeatures label.feature_header PString854 Readings Page in the CardioGuide™ HF App readingComplianceStatus label.feature_reading PString855 What is Reading Compliance Status? faqCategoryReadingComplianceStatusFaq1_question label.feature_reading_q1 PString857 Patients with a {32} are encouraged to use their home electronics (\pillow\) to take and send a device reading daily or as prescribed by their clinic. The {0} provides an indication if the {35} has successfully received a reading each day.	UI Info	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Question 1: <label.feature_med_messages_q1> Answer 1: <label.feature_med_messages_a1> Questions 2: <label.feature_med_messages_q2> Answer 2: <label.feature_med_messages_a2> Question 3: <label.feature_med_messages_q3> Answer 3: <label.feature_med_messages_a3> Question 4: <label.feature_med_reminder_q1> Answer 4: <label.feature_med_reminder_a1> Question 5: <label.feature_med_reminder_q2> Answer 5: <label.feature_med_reminder_a2> Question 6: <label.feature_med_reminder_q3> Answer 6: <label.feature_med_reminder_a3> Sub-Header: <label.feature_passwords> Question 1: <label.feature_passwords_q1> Answer 1: <label.feature_passwords_a1> Question 2: <label.feature_passwords_q2> Answer 2: <label.feature_passwords_a2> Question 3: <label.feature_passwords_q3> Answer 3: <label.feature_passwords_a3>	{n}The {0} will provide a reminder until the {35} has successfully received a reading for the day. faqCategoryReadingComplianceStatusFaq1_answer label.feature_reading_a1 PString859 I sent a {21} Reading, but it isn't showing up in the {0}. faqCategoryReadingComplianceStatusFaq2_question label.feature_reading_q2 PString861 It may take some time for the information to be communicated to the clinic. Contact your clinic, if you have any questions. faqCategoryReadingComplianceStatusFaq2_answer label.feature_reading_a2 PString862 History Calendar glabal_historyCalendar label.feature_history PString863 What is the History Calendar? faqCategoryHistoryCalendarFaq1_question label.feature_history_q1 PString864 The History Calendar shows a calendar month view that indicates the days where a PA Pressure Sensor reading was received and the number of messages received on each day.		

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		faqCategoryHistoryCalendarFaq1_answer label.feature_history_a1 PString865 How much history is available? faqCategoryHistoryCalendarFaq2_question label.feature_history_q2 PString866 The History Calendar displays up to 180 days of data. Previous messages received will display starting on the day in which you registered the {0}. faqCategoryHistoryCalendarFaq2_answer label.feature_history_a2 PString867 What do the symbols on each day mean? faqCategoryHistoryCalendarFaq3_question label.feature_history_q3 PString868 Reading Indicators are shown directly above the numerical date value: {n}{n}- A date with a green circle and a checkmark indicates a reading was received on that day. {n}- A date with a green circle and a question mark indicates a reading has not yet been received on that day. {n}- A date with a gray X indicates a reading was not		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		received on that day.\n- An empty circle indicates a future date. {n}{n}Message Indicators are shown in a small circle at the upper-right section of the reading indicator: {n}{n}- The lack of a message indicator indicates that no messages were received on that day. {n}- A number in a message indicator indicates the number of messages received on that day. {n}- When a number is present in the message indicator, the color of the message indicator represents the highest priority of the message(s) received that day. The colors are: Red (high priority), Yellow (medium priority), Green (low priority), Blue (lowest priority). faqCategoryHistoryCalendarFaq3_answer label.feature_history_a3 PString869 How do I see more details of a given day? faqCategoryHistoryCalendarFaq4_question label.feature_history_q4 PString870 Upon selecting a date on the displayed calendar, the		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		{0} will show a summary of the reading and message events for that day below the calendar. faqCategoryHistoryCalendarFaq4_answer label.feature_history_a4 PString871 Messages glabal_messages label.feature_messages PString872 What are messages? faqCategoryMessagesFaq1_question label.feature_messages_q1 PString873 The clinic may send messages to the app based on the clinic’s review of the patient’s{32} reading. All messages are initiated by the clinician. {n}Message types include:\n- Reading Reviewed: This message may be sent by the clinician after they have reviewed the {32} reading received on the {35}. {n}- Treatment Plan Reminder: This message may be sent by the clinician to remind the patient to follow their dietary guidelines and to take their prescribed medications. This message serves as a reminder to continue with instructions previously		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		provided by the clinic and requires no other immediate action by the patient. However, if the patient is experiencing changes in their heart failure symptoms, they should call the heart failure clinic. {n}- Reading Reminder: This message may be sent by the clinician when they notice that the PA pressure sensor reading has not been received. Upon receiving this message, the patient should use their home electronics to send a reading. {n}- Call Clinic: This message may be sent by the clinician when they want the patient to call the clinic to discuss their treatment. Upon receiving this message, the patient should call their heart failure clinic. {n}- Lab Visit: This message may be sent by the clinician to request that the patient goes to their lab to have a blood test performed. Upon receiving this message, the patient should schedule an appointment with their lab. {n}- Custom Message: This message contains clinician entered information or instructions. Upon receiving this message, the patient should follow any		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		instructions provided in the message. faqCategoryMessagesFaq1_answer label.feature_messages_a1 PString874 Can I reply to messages using the app? faqCategoryMessagesFaq2_question label.feature_messages_q2 PString875 No, the app does not provide the capability to reply to messages. If the you have questions about a message or need to reply to a message, you should call your heart failure clinic directly. faqCategoryMessagesFaq2_answer label.feature_messages_a2 PString876 How do I clear messages? faqCategoryMessagesFaq3_question label.feature_messages_q3 PString877 Messages can be cleared from the Home page by selecting the “Mark as Read” button at the bottom of the message. Once cleared, messages can be viewed in the Calendar history view and are never deleted. Messages will not be displayed to “Other” user types until the “Patient/Primary Caregiver”		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		user marks the messages as read. faqCategoryMessagesFaq3_answer label.feature_messages_a3 PString878 Medication Messages medicationMessages label.feature_med_messages PString879 The patient will only receive Medication messages if: {n}- The clinic uses the {34} to communicate medication information to the patient. (The clinic may also use other methods to provide medication instructions.) {n}- The clinic has enabled the patient's {0} to receive medication information and the clinic uses the {35} in the way described above. {n}- The patient is a Patient/Primary Caregiver user type. Note: Other user types will see medication messages on the Calendar history page after they have been viewed and acknowledged by the Patient/Primary Caregiver user. {n}IF YOU EVER HAVE QUESTIONS ABOUT YOUR MEDICATIONS OR RECEIVED MEDICATION INSTRUCTIONS, PLEASE CALL YOUR CLINIC.		



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		faqMedicationMessagesSupportingText label.feature_med_messages_info PString880 What are Medication Messages? faqCategoryMedicationMessagesFaq1_question label.feature_med_messages_q1 PString881 Medication Messages provide medication instructions from your clinic. They are sent by the clinic when they want the patient to start, change, or stop a medication. All messages are initiated by the clinician. faqCategoryMedicationMessagesFaq1_answer label.feature_med_messages_a1 PString882 What do I do when I receive a medication message? faqCategoryMedicationMessagesFaq2_question label.feature_med_messages_q2 PString883 When you receive a Medication Message, select the “Begin” button at the bottom of the message card. If you are taking multiple medications, you will be asked to confirm that you are taking your		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		other current Heart Failure medications correctly. {n}While you are performing this confirmation step, it is important to answer the questions accurately in order to help ensure that you are receiving appropriate medication information. {n}If you confirm that you are taking your current Heart Failure medications correctly, you will be provided the new medication information. Review the new medication information closely. If you have any questions about the new medication information, call your clinic. {n}If you are unable to confirm that you are taking your current Heart Failure medications correctly, you will be instructed to call your clinic. faqCategoryMedicationMessagesFaq2_answer label.feature_med_messages_a2 PString884 Can I take a Medication Message to my pharmacy to get a prescription filled? faqCategoryMedicationMessagesFaq3_question label.feature_med_messages_q3 PString885		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		No. Medication messages are NOT considered formal prescriptions. Medication messages should only contain information about medications that have already been prescribed by your clinician. If you do not have the medication available, please call your clinic. faqCategoryMedicationMessagesFaq3_answer label.feature_med_messages_a3 PString886 What are Medication Reminders? faqCategoryMedicationRemindersFaq1_question label.feature_med_reminder_q1 PString887 If the patient receives medication information on the {0}, the patient will be able to use the {0} to provide reminders when it is time to take the medication as prescribed. faqCategoryMedicationRemindersFaq1_answer label.feature_med_reminder_a1 PString888 How do I setup the medication reminders? faqCategoryMedicationRemindersFaq2_question label.feature_med_reminder_q2		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		PString889 Medication reminders are specific to each medication. The number of reminders sent each day depends on the frequency of the medication. To set up reminder for a medication, select the medication from the Medications page. The patient can change the default reminder times by adjusting the time(s) when to receive the reminder. faqCategoryMedicationRemindersFaq2_answer label.feature_med_reminder_a2 PString890 How do I change the times in which reminders are sent each day? faqCategoryMedicationRemindersFaq3_question label.feature_med_reminder_q3 PString891 To change the times in which medication reminders are sent, select the reminder time from the Account Preference page under "Medication Times". faqCategoryMedicationRemindersFaq3_answer label.feature_med_reminder_a3 PString892 Passwords passwords label.feature_passwords		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		PString893 How do I change my password? faqCategoryPasswordsFaq1_question label.feature_passwords_q1 PString894 To change your password, select the Account button from the Home page. From there, select “Change Password” under My Account. faqCategoryPasswordsFaq1_answer label.feature_passwords_a1 PString895 What should I do if I forget my password? faqCategoryPasswordsFaq2_question label.feature_passwords_q2 PString896 If you forget your password, call {22} Remote Care Technical Support at +1-877-myCMEMS ({16}) Monday – Friday 8:00 am to 8:00 pm (EST). Outside of the US, please contact your local Abbott representative or your clinic.", faqCategoryPasswordsFaq2_answer label.feature_passwords_a2 PString897		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		Can I access the app without entering a password each time? faqCategoryPasswordsFaq3_question label.feature_passwords_q3 PString898 By default, you will be required to enter your password each time you log into the app. You may select the “Keep me logged in” feature to stay logged in and help ensure that notifications are received in a timely manner. However, for security purposes, it is important to note that if you choose to remain continuously logged into the app, your device should be password protected. faqCategoryPasswordsFaq3_answer label.feature_passwords_a3		
PAGUI4088	The FAQ TOPIC PANEL that describes common issues that may be seen in the Patient Application and their respective fixes is displayed in the FAQ PANEL. The display order is described below. Header: <label.troubleshooting_header> Sub-Header: <label.troubleshooting_help> Question 1: <label.troubleshooting_help_q1> Answer 1: <label.troubleshooting_help_a1> Sub-Header: <label.troubleshooting_travel> Question 1: <label.troubleshooting_travel_q1> Answer 1: <label.troubleshooting_travel_a1> Question 2: <label.troubleshooting_travel_q2> Answer 2: <label.troubleshooting_travel_a2> Question 2: <label.troubleshooting_travel_q3> Answer 2: <label.troubleshooting_travel_a3>	PString900 Troubleshooting faqCategoryTroubleshooting label.troubleshooting_header PString901 General Help generalHelp label.troubleshooting_help PString902 Who should I contact if I need help?	UI Info	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Sub-Header: <label.troubleshooting_internet> Question 1: <label.troubleshooting_internet_q1> Answer 1: <label.troubleshooting_internet_a1> Sub-Header: <label.troubleshooting_app> Question 1: <label.troubleshooting_app_q1> Answer 1: <label.troubleshooting_app_a1> Sub-Header: <label.troubleshooting_security> Question 1: <label.troubleshooting_security_q1> Answer 1: <label.troubleshooting_security_a1>	faqCategoryGeneralHelpFa q1_question label.troubleshooting_help _q1 PString903 Travel travel label.troubleshooting_trave l PString904 In the US, you can call {22} Remote Care Technical Support at +1-877- myCMEMS ({16}) Monday – Friday 8:00 am to 8:00 pm (EST). Outside of the US, please contact your local Abbott representative or your clinic. faqCategoryGeneralHelpFa q1_answer label.troubleshooting_help _a1 PString905 I am traveling. Do I need to take my mobile device with me? faqCategoryTravelFaq1_que stion label.troubleshooting_trave l_q1 PString906 Yes. You should have your mobile device with you and your app running so that you can receive messages from your clinic. faqCategoryTravelFaq1_ans wer label.troubleshooting_trave l_a1		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		PString907 I am traveling overseas and might not have access to internet for a few days. What should I do? faqCategoryTravelFaq2_question label.troubleshooting_travel_q2 PString908 We recommend that you contact your clinic to make arrangements in case the clinic needs to communicate with you during your travel. faqCategoryTravelFaq2_answer label.troubleshooting_travel_a2 PString909 Internet Connection Problems internetConnectionProblems label.troubleshooting_internet PString910 I turned on Wi-Fi. Now what? faqCategoryInternetConnectionProblemsFaq1_question label.troubleshooting_internet_q1 PString911 If Wi-Fi and/or cellular data is enabled on your phone, the {0} will continue to receive data as long as there is sufficient network		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		connectivity. The {0} must be running in the background to be able to complete its functions (tasks). Do not kill or force close the {0}. If the app is force closed, it will not be able to perform its activities as programmed by the clinic. faqCategoryInternetConnectionProblemsFaq1_answer label.troubleshooting_internet_a1 PString912 Troubleshooting App troubleshootingApp label.troubleshooting_app PString913 The {0} displays a notification indicating that the {0} or the operating system needs to be updated. What should I do? faqCategoryTroubleshootingAppFaq1_question label.troubleshooting_app_q1 PString914 To help ensure the {0} is operating correctly, we recommend that you keep the {0} and operating system updated to the latest versions. Abbott suggests leaving app auto-updates enabled so that the app can update itself as new updates are available.		

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		{n}- If the {0} or operating system are out of date but are still functional, you will be able to dismiss the notification and proceed to use the app. However, in this case, be aware that the full functionality of the {0} may not be supported. {n}- If the {0} or operating system are out of date to the extent that functionality is no longer supported, you will not be able to dismiss the notification. In this case, you will need to update the {0}, your operating system, or both prior to using the {0}. faqCategoryTroubleshootingAppFaq1_answer label.troubleshooting_app_a1 PString915 Mobile Device Security mobileDeviceSecurity label.troubleshooting_security PString916 Is there any sensitive data displayed on the {0}? faqCategoryMobileDeviceSecurityFaq1_question label.troubleshooting_security_q1 PString917 Depending upon how the {0} is configured and utilized by your clinic, the {0} may display private health information such as		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		private messages from your clinician, your medication information, and data related to your PA Sensor. You can keep your information on your mobile device more secure by following the steps below: {n}- Use a secure screen lock: Set up a passcode, pattern, Apple {n}- Touch ID fingerprint identity, sensor, Apple Face ID secure authentication, etc. to unlock your phone. {n}- Set a short screen lock time (for example, 30 seconds): If your mobile device is unlocked and you are not using it, it will automatically lock after a short amount of time. {n}- Update the {0} and your mobile device operating system whenever updates are available. {n}Alterations to your mobile device operating system, such as rooting, will make your mobile device more unstable or insecure. faqCategoryMobileDeviceSecurityFaq1_answer label.troubleshooting_security_a1 PString923 I am traveling and received a notification about elevation. What is it and what should I do?		



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		faqCategoryTravelFaq3_question label.troubleshooting_travel_q3 PString924 The {21} may be affected when measurements are taken at a significantly different elevation. For example, if you normally take readings at sea level and you take a reading in a mountainous area, the difference in elevation may impact the {21} reading. When a {21} reading indicates that you are at a significantly different elevation than normal, a notification will appear in the {0}. When this notification is received, please contact your clinic to confirm your treatment instructions. faqCategoryTravelFaq3_answer label.troubleshooting_travel_a3		
PAGUI1780	The FAQ TOPIC PANEL consists of: TOPIC HEADER COMPONENT TOPIC SUB-SECTION HEADER COMPONENT ARTICLE DROP DOWN ACTION TRIGGER COMPONENT ARTICLE CONTENT LABEL COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1776	8.12.1.22.2 FAQ TOPIC PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1781	<Picture> FAQ TOPIC PANEL (example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
	About the CardioGuide™ HF App What does the CardioGuide™ HF App do? ▾ How does the CardioGuide™ HF App communicate with the clinic? ▾ Who can use the CardioGuide™ HF App? ▾			
PAGUI1900	<Picture> FAQ TOPIC PANEL (example: with sub-sections)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
	Getting Started LANGUAGE How do I select the language of the myCardioMEMS™ app? ▾ STAYING CONNECTED How do I stay connected to my clinic? ▾ Do I always need to be connected to the Internet? ▾ How do I turn on Wi-Fi? ▾			
PAGUI1777	8.12.1.22.3 FAQ TOPIC PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1782	8.12.1.22.3.1 TOPIC HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1783	Component: TOPIC HEADER COMPONENT State: N/A Display: Always displayed in the FAQ TOPIC PANEL. The TOPIC HEADER COMPONENT displays the topic from the content managment location.	PString806 About the {0} faqCategoryAboutTheCardioApp	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

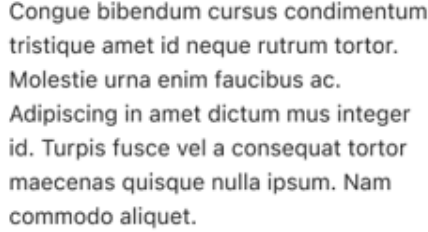
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Label: N/A Behavior: Text.	label.about_header PString821 Getting Started faqCategoryGettingStarted label.start_header PString852 App Features faqCategoryAppFeatures label.feature_header		
PAGUI1784	<Picture> TOPIC HEADER COMPONENT (example) About the CardioGuide™ HF App		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1786	<Picture> TOPIC HEADER COMPONENT (example) Getting Started		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1787	<Picture> TOPIC HEADER COMPONENT (example) App Features		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1901	8.12.1.22.3.2 TOPIC SUB-SECTION HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1902	Component: TOPIC SUB-SECTION HEADER COMPONENT State: N/A Display: Conditionally displayed in the FAQ TOPIC PANEL. The TOPIC HEADER COMPONENT displays the topic sub-section if applicable from the content mangement location. Label: N/A Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1903	<Picture> TOPIC SUB-SECTION HEADER COMPONENT (example) LANGUAGE		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1904	<Picture> TOPIC SUB-SECTION HEADER COMPONENT (example) STAYING CONNECTED		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1788	8.12.1.22.3.3 ARTICLE DROP DOWN ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1789	Component: ARTICLE DROP DOWN ACTION TRIGGER COMPONENT State: Open/ Closed Display: Conditionally displayed in the FAQ TOPIC PANEL. The ARTICLE DROP DOWN ACTION TRIGGER COMPONENT is displayed when there is an article available in the content management location for the particular FAQ TOPIC PANEL. The		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	default state is closed. The number of ARTICLE DROP DOWN ACTION TRIGGER COMPONENTS displayed in a panel is determined by how many unique articles there are stored in the content management location. The ARTICLE DROP DOWN ACTION TRIGGER COMPONENT displays a title for the associated article and an icon to indicate if the drop dow is open or closed. Label: N/A Behavior: Text, icon, and trigger. Upon activation in the Closed state, the ARTICLE DROP DOWN ACTION TRIGGER COMPONENT expands the FAQ TOPIC PANEL down to display the associated ARTICLE CONTENT LABEL COMPONENT and changes the orientation of the icon on the right side. Upon activation in the Open state, the ARTICLE DROP DOWN ACTION TRIGGER COMPONENT stops displaying the associated ARTICLE CONTENT LABEL COMPONENT to shrink the FAQ TOPIC PANEL and changes the orientation of the icon on the right side.			
PAGUI1791	<Picture> ARTICLE DROP DOWN ACTION TRIGGER COMPONENT (default: closed) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1792	<Picture> ARTICLE DROP DOWN ACTION TRIGGER COMPONENT (example: open) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1793	8.12.1.22.3.4 ARTICLE CONTENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1794	Component: ARTICLE CONTENT LABEL COMPONENT State: N/A Display: Conditionally displayed in the FAQ TOPIC PANEL. The ARTICLE CONTENT LABEL COMPONENT is displayed after the ARTICLE DROP DOWN ACTION TRIGGER COMPONENT is activated in the Closed state. The text displayed is pulled from the associated content managment location. Label: N/A Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1795	<Picture> ARTICLE CONTENT LABEL COMPONENT (example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1107	8.12.1.23 FORGOT PASSWORD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1108	8.12.1.23.1 FORGOT PASSWORD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1109	The FORGOT PASSWORD PANEL is A2B2C series of screens that display once the user presses FORGOT PASSWORD ACTION TRIGGER. The FORGOT PASSWORD PANEL is controlled by a third-party vendor and this section serves as examples of the screens.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1110	The FORGOT PASSWORD PANEL is an element of: LOGIN PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1119	The FORGOT PASSWORD PANEL consists of: USER DETAILS EMAIL ID PANEL USER DETAILS CONFIRMATION PANEL USER DETAILS NEW PASSWORD PANEL (See Shared Elements: Panels)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1111	8.12.1.23.2 FORGOT PASSWORD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1112	The FORGOT PASSWORD PANEL is a series of panels, and their details are described below.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1113	8.12.1.23.3 FORGOT PASSWORD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1114	8.12.1.23.3.1 USER DETAILS EMAIL ID PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2070	8.12.1.23.3.1.1 USER DETAILS EMAIL ID PANEL (DISPLAY)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2077	The USER DETAILS EMAIL ID PANEL is displayed after the user selects FORGOT PASSWORD ACTION TRIGGER COMPONENT. The USER DETAILS EMAIL ID PANEL is displayed over the entire screen except for the OS STATUS PANEL. The USER DETAILS EMAIL ID PANEL is controlled by a third party and configurable for the Patient Application; the third party may change design details like strings, font, or icons. This section is for reference only.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2078	The USER DETAILS EMAIL ID PANEL consists of: CANCEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) USER DETAILS ICON COMPONENT (See Shared Elements: Components) USER DETAILS HEADER COMPONENT (See Shared Elements: Components) REQUIRED FIELD ERROR COMPONENT (See Shared Elements: Components) CLAIM NOT VERIFIED ERROR COMPONENT VERIFICATION INSTRUCTIONS LABEL COMPONENT INFORMATION REQUIRED ERROR COMPONENT (See Shared Elements: Components) INVALID EMAIL ERROR COMPONENT EMAIL ID INPUT COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

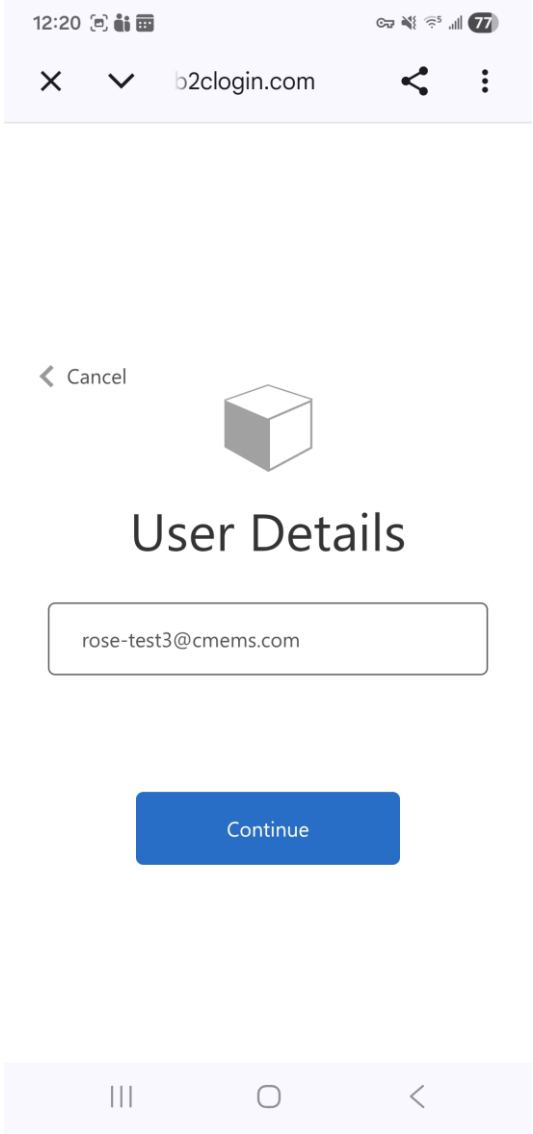


myCardioMEMS Patient Application Graphical User Interface Specification

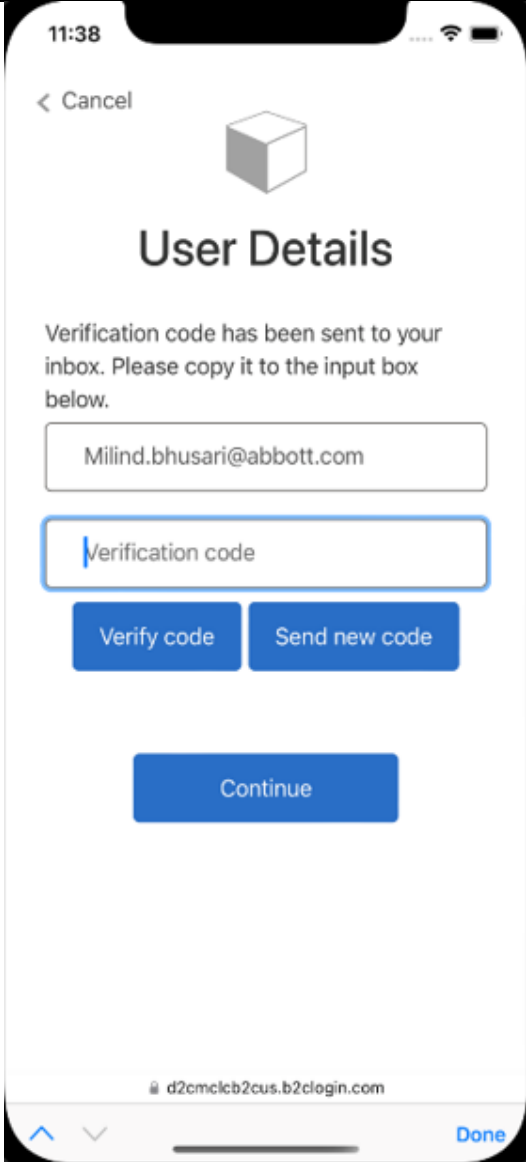
ID		English US Strings	Object Type	Project
	SEND VERIFICATION CODE ACTION TRIGGER COMPONENT VERIFICATION CODE INPUT COMPONENT (See Shared Elements: Components) ALPHABETICAL KEYBOARD COMPONENT (See Shared Elements: Components) NUMERICAL KEYBOARD COMPONENT (See Shared Elements: Components) VERIFY CODE ACTION TRIGGER COMPONENT SEND NEW VERIFICATION CODE ACTION TRIGGER COMPONENT CONTINUE ACTION TRIGGER COMPONENT (See Shared Elements: Components)			
PAGUI2071	8.12.1.23.3.1.2 USER DETAILS EMAIL ID PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2087	<Picture> USER DETAILS EMAIL ID PANEL (example: email input)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
PAGUI4593	<Picture> USER DETAILS EMAIL ID PANEL (example: email input confirmation)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI2710	<Picture> USER DETAILS EMAIL ID PANEL (example: code input)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI2072	8.12.1.23.3.1.3 USER DETAILS EMAIL ID PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2706	8.12.1.23.3.1.3.1 CANCEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2707	See Shared Elements: Components		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2711	8.12.1.23.3.1.3.2 USER DETAILS ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2712	See Shared Elements: Components		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2713	8.12.1.23.3.1.3.3 USER DETAILS HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2714	See Shared Elements: Components		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2632	8.12.1.23.3.1.3.4 REQUIRED FIELD ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4106	See Shared Elements: Components.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4845	8.12.1.23.3.1.3.5 INVALID EMAIL ERROR COMPONENT		Heading	CloudDRxRTR
PAGUI4847	<Picture> INVALID EMAIL ERROR COMPONENT Please enter a valid email address.		UI Design	CloudDRxRTR
PAGUI2673	8.12.1.23.3.1.3.6 CLAIM NOT VERIFIED ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2674	Component: CLAIM NOT VERIFIED ERROR COMPONENT State: Email/ Phone Display: Conditionally displayed in the USER DETAILS EMAIL ID PANEL when the software is unable to send a code to the input of the EMAIL ID INPUT COMPONENT or the PHONE ID INPUT COMPONENT. The CLAIM NOT VERIFIED ERROR COMPONENT is displayed in the email state when in the USER DETAILS EMAIL ID PANEL and the phone state when in the TWO FACTOR AUTHENTICATION PANEL. Label: <ots.config.claim_not_verified> - determined by a third party vendor Behavior: Text.	PString1075 Claim not verified: [Email Address] ots.config.claim_not_verified PString1084 E-mail address verified. You can now continue. ots.config.email_address_verified	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2675	<Picture> CLAIM NOT VERIFIED ERROR COMPONENT (Email, example) Claim not verified: [Email Address]		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2642	8.12.1.23.3.1.3.7 VERIFICATION INSTRUCTIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2643	Component: VERIFICATION INSTRUCTIONS LABEL COMPONENT State: Send/ Code/ Incorrect Code/ Too Many Attempts Display: Always displayed in the USER DETAILS EMAIL ID PANEL. The default state is send. Each state displays a different set of instructions. The VERIFICATION INSTRUCTIONS LABEL COMPONENT displays the code state when a verification code has succesfully been sent to the input of the EMAIL ID INPUT COMPONENT or the PHONE ID INPUT COMPONENT. The VERIFICATION INSTRUCTIONS LABEL COMPONENT displays the incorrect code state after the user selects the VERIFY CODE ACTION TRIGGER COMPONENT and the input in the VERIFICATION CODE INPUT COMPONENT does not match the code sent to the inputed email or phone number. The VERIFICATION INSTRUCTIONS LABEL COMPONENT is displayed in the too many attempts state when the user has failed to input the correct verification code three times. Label: determined by third party vendor <ots.config.verification_instructions_default>, <ots.config.verification_code_sent>, <ots.config.verification_code_incorrect>, <ots.config.verification_too_many_attempts> Behavior: Text.	PString1076 Verification is necessary. Please click Send button. ots.config.verification_instructions_default PString1077 Verification code has been sent to your inbox. Please copy it to the input box below. ots.config.verification_code_sent PString1078 That code is incorrect. Please try again. ots.config.verification_code_incorrect PString1079 You've made too many incorrect attempty. Please try again later. ots.config.verification_too_many_attempts	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2644	<Picture> VERIFICATION INSTRUCTIONS LABEL COMPONENT (default: Send, example) Verification is necessary. Please click Send button.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2657	<Picture> VERIFICATION INSTRUCTIONS LABEL COMPONENT (Code, example) Verification code has been sent to your inbox. Please copy it to the input box below.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2658	<Picture> VERIFICATION INSTRUCTIONS LABEL COMPONENT (Incorrect Code, example) <u>That code is incorrect. Please try again.</u>		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2659	<Picture> VERIFICATION INSTRUCTIONS LABEL COMPONENT (Too Many Attempts, example) You've made too many incorrect attempts. Please try again later.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2645	8.12.1.23.3.1.3.8 INFORMATION REQUIRED ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4107	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4846	Component: INVALID EMAIL ERROR COMPONENT State: Email/ Phone Display: Conditionally displayed in the USER DETAILS EMAIL ID PANEL when the user does not input a valid email address. Label: <msg.gb.emailInvalid> Behavior: Text.		UI Info	CloudDRxRTR
PAGUI2708	8.12.1.23.3.1.3.9 EMAIL ID INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2088	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2662	8.12.1.23.3.1.3.10 SEND VERIFICATION CODE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2663	Component: SEND VERIFICATION CODE ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS EMAIL ID PANEL until the user selects the SEND VERIFICATION CODE ACTION TRIGGER COMPONENT and a code is sent to input of the EMAIL ID INPUT COMPONENT. Label: <ots.config.send_verification_code> - determined by third party vendor Behavior: Text and trigger. Upon activation, the SEND VERIFICATION CODE ACTION TRIGGER COMPONENT sends a code to the input of the EMAIL ID INPUT COMPONENT.	PString1080 Send verification code ots.config_send_verification_code	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2664	<Picture> SEND VERIFICATION CODE ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2723	8.12.1.23.3.1.3.11 VERIFICATION CODE INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2724	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2725	8.12.1.23.3.1.3.12 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2726	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2727	8.12.1.23.3.1.3.13 NUMERIC KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2728	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2729	8.12.1.23.3.1.3.14 VERIFY CODE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2730	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4687	8.12.1.23.3.1.3.15 SEND NEW VERIFICATION CODE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4688	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI2733	8.12.1.23.3.1.3.16 CONTINUE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2734	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2080	8.12.1.23.3.2 USER DETAILS CONFIRMATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2084	8.12.1.23.3.2.1 USER DETAILS CONFIRMATION PANEL (DISPLAY)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2109	The USER DETAILS CONFIRMATION PANEL is displayed once the user selects the VERIFY CODE ACTION TRIGGER and verification is successful.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2108	The USER DETAILS CONFIRMATION PANEL consists of: CANCEL ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR

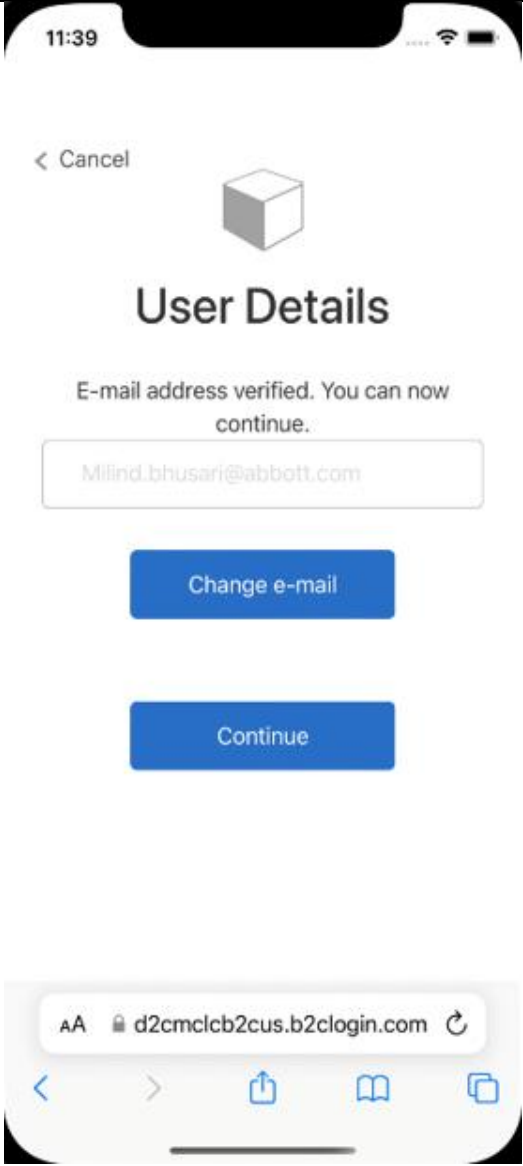
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	USER DETAILS ICON COMPONENT (See Shared Elements: Components) USER DETAILS HEADER COMPONENT (See Shared Elements: Components) CONFIRMATION INSTRUCTIONS LABEL COMPONENT EMAIL ID INPUT COMPONENT (See Shared Elements: Components) ALPHABETICAL KEYBOARD COMPONENT (See Shared Elements: Components) CHANGE EMAIL ACTION TRIGGER COMPONENT CONTINUE ACTION TRIGGER COMPONENT (See Shared Elements: Components)			CloudDRxPS
PAGUI2085	8.12.1.23.3.2.2 USER DETAILS CONFIRMATION PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2110	<Picture> USER DETAILS CONFIRMATION PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI2086	8.12.1.23.3.2.3 USER DETAILS CONFIRMATION PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2735	8.12.1.23.3.2.3.1 CANCEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2736	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2737	8.12.1.23.3.2.3.2 USER DETAILS ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2738	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2739	8.12.1.23.3.2.3.3 USER DETAILS HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2740	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2741	8.12.1.23.3.2.3.4 CONFIRMATION INSTRUCTIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2113	Component: CONFIRMATION INSTRUCTIONS LABEL COMPONENT State: N/A Display: Always displayed in the USER DETAILS CONFIRMATION PANEL Label: <ots.config.email_address_verified> determined by third party vendor Behavior: Text.	PString1084 E-mail address verified. You can now continue. ots.config.email_address_verified	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4644	<Picture> CONFIRMATION INSTRUCTIONS LABEL COMPONENT E-mail address verified. You can now continue.		UI Design	CloudDRxRTR
PAGUI2743	8.12.1.23.3.2.3.5 EMAIL ID INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2744	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4642	8.12.1.23.3.2.3.6 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRxRTR
PAGUI4643	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI2755	8.12.1.23.3.2.3.7 CHANGE EMAIL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2112	Component: CHANGE EMAIL ACTION TRIGGER State: N/A		UI Info	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Display: Always displayed in the USER DETAILS CONFIRMATION PANEL Label: N/A - determined by third party vendor Behavior: Text and trigger. Upon activation, the CHANGE EMAIL ACTION TRIGGER displays the USER DETAILS ID PANEL.			CloudDRxPS
PAGUI2756	<Picture> CHANGE EMAIL ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2757	8.12.1.23.3.2.3.8 CONTINUE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2758	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1117	8.12.1.23.3.3 USER DETAILS NEW PASSWORD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3141	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2354	8.12.1.24 FREEFORM MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2355	8.12.1.24.1 FREEFORM MESSAGE CARD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2356	The FREEFORM MESSAGE CARD PANEL is conditionally displayed if a Freeform Message is received by the Patient Application at least once since the Earliest History Calendar Date. A unique FREEFORM MESSAGE CARD PANEL is displayed for every Freeform Message that is received by the Patient Application. Multiple FREEFORM MESSAGE CARD PANELS may be displayed in the same screen or panel. The FREEFORM MESSAGE CARD PANEL is in the active state by default. The FREEFORM MESSAGE CARD PANEL CARD PANEL changes state from active to complete when the user completes all required actions in the card. The elements in the FREEFORM MESSAGE CARD PANEL listed below change behavior along these two states.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2357	The FREEFORM MESSAGE CARD PANEL is an element of: ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2358	The FREEFORM MESSAGE CARD PANEL consists of: FREEFORM BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) FREEFORM CONTENT LABEL COMPONENT MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2359	8.12.1.24.2 FREEFORM MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI2360	<Picture> FREEFORM MESSAGE CARD PANEL (default: active, example) NEW CLINICIAN MESSAGE FRIDAY, MAR 10 Please reduce your salt intake for the next week while we monitor your changes in pressures. Mark as Read		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2361	<Picture> FREEFORM MESSAGE CARD PANEL (complete, example) NEW CLINICIAN MESSAGE FRIDAY, MAR 10 Please reduce your salt intake for the next week while we monitor your changes in pressures.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2362	8.12.1.24.3 FREEFORM MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2363	8.12.1.24.3.1 FREEFORM BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2364	Component: FREEFORM BANNER COMPONENT State: Active/ Complete Display: Always displayed in the FREEFORM MESSAGE CARD PANEL. The default state is active. In the active state, the FREEFORM BANNER COMPONENT displays the active icon and the active background color. In the complete state, the FREEFORM BANNER COMPONENT displays the complete icon and complete background color. Label: <label.gb.new_message>; Text is all capitalized.	PString179 New Clinician Message HCM label.gb.new_message	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

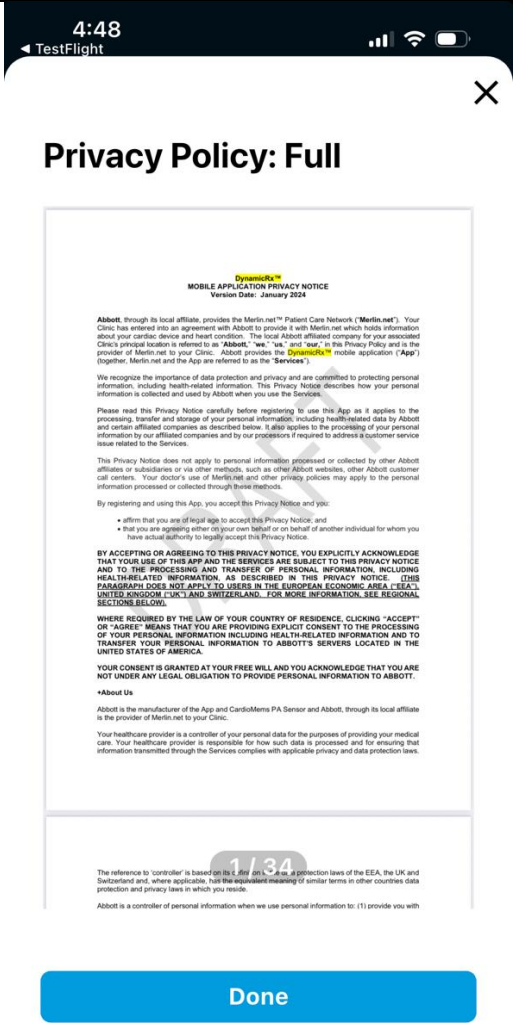
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Behavior: Icon and text. The FREEFORM BANNER COMPONENT changes state with the message card panel it is an element of.			
PAGUI2365	<Picture> FREEFORM BANNER COMPONENT (default: Active) 	PString179 New Clinician Message HCM label.gb.new_message	UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2366	<Picture> FREEFORM BANNER COMPONENT (Complete) 	PString179 New Clinician Message HCM label.gb.new_message	UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2367	8.12.1.24.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2368	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2369	8.12.1.24.3.3 FREEFORM CONTENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2370	Component: FREEFORM CONTENT LABEL COMPONENT State: N/A Display: Always displayed in the FREEFORM MESSAGE CARD PANEL. Label: N/A - content is determined by clinician input in the portal Behavior: Text. The FREEFORM CONTENT LABEL COMPONENT is determined by clinician input in the portal.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2371	<Picture> FREEFORM CONTENT LABEL COMPONENT (example) Please reduce your salt intake for the next week while we monitor your changes in pressures.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2375	8.12.1.24.3.4 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2376	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3925	8.12.1.25 FULL PRIVACY POLICY PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI3926	8.12.1.25.1 FULL PRIVACY POLICY PANEL (Display)		Heading	CloudDRxRTR CloudDRxPS

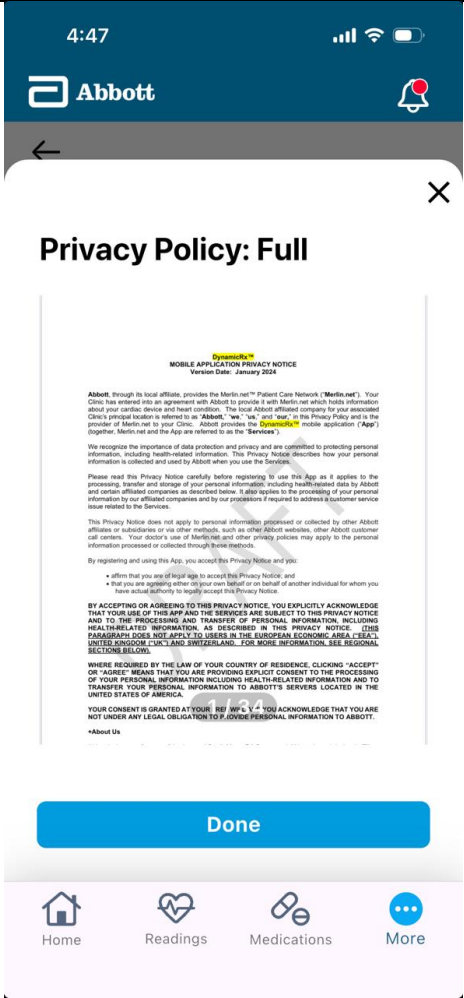
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3927	The FULL PRIVACY POLICY PANEL is displayed after the user selects the hyperlink within the SHORT PRIVACY POLICY DISPLAY COMPONENT. The FULL PRIVACY POLICY PANEL is also displayed after the user selects the PRIVACY POLICY ACTION TRIGGER COMPONENT in the ABOUT PANEL. From the POLICY SCREEN, the FULL PRIVACY POLICY PANEL covers the entire screen except the OS STATUS PANEL. From the ABOUT PANEL, the FULL PRIVACY POLICY PANEL covers the entire screen except the OS STATUS PANEL, GLOBAL HEADING PANEL, and GLOBAL NAVIGATION PANEL.		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3928	The FULL PRIVACY POLICY PANEL is an element of: POLICY SCREEN ABOUT PANEL		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3929	The FULL PRIVACY POLICY PANEL consists of: OS STATUS PANEL GLOBAL HEADING PANEL CLOSE PANEL ACTION TRIGGER COMPONENT (Shared Elements: Components) FULL PRIVACY POLICY HEADER COMPONENT FULL PRIVACY POLICY DISPLAY COMPONENT DONE ACTION TRIGGER COMPONENT (Shared Elements: Components) GLOBAL NAVIGATION PANEL (Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3930	8.12.1.25.2 FULL PRIVACY POLICY PANEL (Layout)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3931	<Picture> FULL PRIVACY POLICY PANEL (example: after the user selects the hyperlink within the SHORT PRIVACY POLICY DISPLAY COMPONENT)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4698	<Picture> FULL PRIVACY PANEL (example: ABOUT PANEL)			CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI3933	8.12.1.25.3 FULL PRIVACY POLICY PANEL (Elements)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4640	8.12.1.25.3.1 OS STATUS PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI4641	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS
PAGUI4706	8.12.1.25.3.1.1 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4707	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI3936	8.12.1.25.3.2 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3937	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI3938	8.12.1.25.3.3 FULL PRIVACY POLICY HEADER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3939	Component: FULL PRIVACY POLICY HEADER COMPONENT State: N/A Display: Always displayed in the FULL PRIVACY POLICY PANEL. Label: <label.gb.PP_full> Behavior: Text.	PString754 Privacy Policy: Full privacyPolicyFullTitle label.gb.PP_full	UI Info	CloudDRxRTR CloudDRxPS
PAGUI3940	<Picture> FULL PRIVACY POLICY HEADER COMPONENT Privacy Policy: Full		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3941	8.12.1.25.3.4 FULL PRIVACY POLICY DISPLAY COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3942	Component: FULL PRIVACY POLICY DISPLAY COMPONENT State: N/A Display: Always displayed in the FULL PRIVACY POLICY PANEL. Label: <other.pdf.Privacy_Policy_Full> - the Patient Application displays the PDF document located with the attached website link Behavior: Text. The content for the FULL PRIVACY POLICY DISPLAY COMPONENT is managed through content management in Merlin.net to enable updates as needed. The PDF document can be multiple pages with the numbr of pages displayed on the bottom and allows scrolling.	PString753 Display a PDF of the applicable document. pdfPrivacyPolicyFull other.pdf.Privacy_Policy_Fu ll	UI Info	CloudDRxRTR CloudDRxPS
PAGUI3943	<Picture> FULL PRIVACY POLICY DISPLAY COMPONENT (example)		UI Design	CloudDRxRTR CloudDRxPS



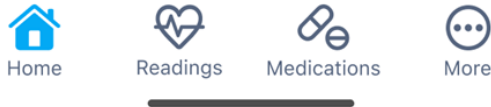
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Lorem ipsum dolor sit amet, consectetur adipiscing elit. Sed enim nunc, euismod vel porta tincidunt, condimentum eget risus. Nulla sodales, leo sit amet sagittis lobortis, sem risus rutrum nulla, sodales iaculis mi tellus pellentesque nulla. Curabitur sit amet felis et orci scelerisque semper. Pellentesque in dui viverra, cursus urna non, fermentum eros. Sed sed lectus vestibulum, semper augue non, auctor lorem. Aliquam ex nisl, finibus non aliquam sit amet, scelerisque venenatis dui. Quisque at sapien accumsan, finibus dolor vel, gravida magna. Nulla eleifend nibh ligula, quis sollicitudin magna dignissim at. Nunc imperdiet nunc vitae quam volutpat, at rhoncus eros dignissim. Ut convallis dui tortor, quis luctus erat egestas sed. Donec et scelerisque dolor, et posuere neque. Sed vitae justo ultrices, lobortis ipsum quis, interdum tortor.			
PAGUI3948	8.12.1.25.3.5 DONE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3949	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI4704	8.12.1.25.3.6 GLOBAL NAVIGATION PANEL		Heading	CloudDRxRTR
PAGUI4705	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI127	8.12.1.26 GLOBAL HEADING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI128	8.12.1.26.1 GLOBAL HEADING PANEL (DISPLAY)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI327	The GLOBAL HEADING PANEL displays secondary application functionality navigation.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI329	The GLOBAL HEADING PANEL is composed of: ABBOTT LOGO COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	NOTIFICATIONS ACTION TRIGGER COMPONENT			CloudDRxPS
PAGUI129	8.12.1.26.2 GLOBAL HEADING PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI717	<Picture> GLOBAL HEADING PANEL (without NOTIFICATIONS ACTION TRIGGER COMPONENT, appears in the POLICY WELCOME PANEL) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI718	<Picture> GLOBAL HEADING PANEL (with NOTIFICATIONS ACTION TRIGGER COMPONENT, no notifications) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI130	8.12.1.26.3 GLOBAL HEADING PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI330	8.12.1.26.3.1 LOGO COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI331	Component: LOGO COMPONENT State: N/A Display: Always displayed in the GLOBAL HEADING PANEL. Label: N/A Behavior: Icon.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI719	<Picture> LOGO COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI345	8.12.1.26.3.2 NOTIFICATIONS ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI346	Component: NOTIFICATIONS ACTION TRIGGER COMPONENT State: No Notifications/ Notifications Display: Conditionally displayed in the GLOBAL HEADING PANEL. NOTIFICATIONS ACTION TRIGGER COMPONENT is not displayed in the POLICY WELCOME SCREEN. The default state is no notifications. The state will change to notifications when there are active notifications in the NOTIFICATIONS SCREEN and an alert icon will be added to the display. Label: N/A Behavior: Icon and trigger. Upon activation, the NOTIFICATIONS ACTION TRIGGER COMPONENT will display the NOTIFICATIONS SCREEN.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI731	<Picture> NOTIFICATIONS ACTION TRIGGER COMPONENT (default: no notifications) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI732	<Picture> NOTIFICATIONS ACTION TRIGGER COMPONENT (notifications) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI131	8.12.1.27 GLOBAL NAVIGATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI132	8.12.1.27.1 GLOBAL NAVIGATION PANEL (DISPLAY)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI334	The GLOBAL NAVIGATION PANEL displays icon triggers that allow the user to navigate through the primary application screens.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4188	The GLOBAL NAVIGATION PANEL is composed of: HOME ACTION TRIGGER COMPONENT READINGS ACTION TRIGGER COMPONENT MEDICATIONS ACTION TRIGGER COMPONENT MORE ACTION TRIGGER COMPONENT		UI Info	CloudDRxRTR CloudDRxPS
PAGUI133	8.12.1.27.2 GLOBAL NAVIGATION PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI721	<Picture> GLOBAL NAVIGATIO PANEL (default: HOME SCREEN) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI134	8.12.1.27.3 GLOBAL NAVIGATION PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI337	8.12.1.27.3.1 HOME ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI338	Component: HOME ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Always displayed in the GLOBAL NAVIGATION PANEL. The default state is enabled. The state will be disabled if the user navigates to the HOME SCREEN and the color of the icon will change. Label: <label.gb.home> Home Behavior: Icon and trigger. Upon triggering, the HOME ACTION TRIGGER COMPONENT will change color to the disabled state and display the HOME SCREEN.	PString1047 Home home_label label.gb.home	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI724	<Picture> HOME ACTION TRIGGER COMPONENT (default: enabled)		UI Design	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

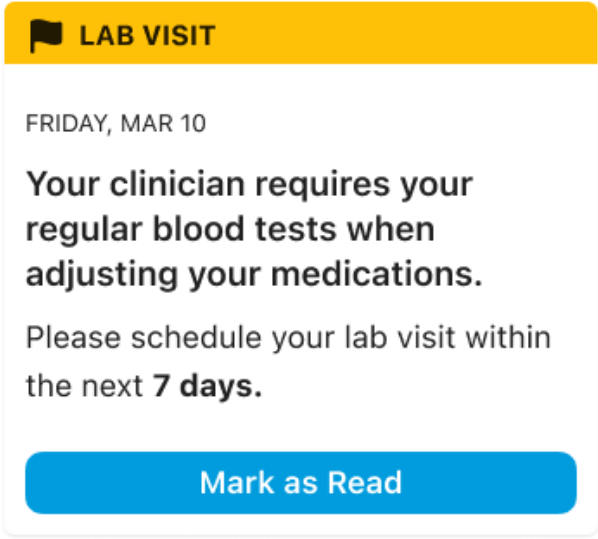
ID		English US Strings	Object Type	Project
	 Home			CloudDRxPS
PAGUI722	<Picture> HOME ACTION TRIGGER COMPONENT (disabled)  Home		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI339	8.12.1.27.3.2 READINGS ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI340	Component: READINGS ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Always displayed in the GLOBAL NAVIGATION PANEL. The default state is enabled. The state will be disabled if the READINGS SCREEN is displayed and the color of the icon will change. Label: <label.gb.readings_header> Readings Behavior: Icon and trigger. Upon being triggered, the READINGS ACTION TRIGGER COMPONENT will change color to the disabled component and display the READINGS SCREEN.	PString427 Readings readings label.gb.readings_header	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI725	<Picture> READINGS ACTION TRIGGER COMPONENT (default: enabled)  Readings		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI726	<Picture> READINGS ACTION TRIGGER COMPONENT (disabled)  Readings		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI341	8.12.1.27.3.3 MEDICATIONS ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4184	Component: MEDICATIONS ACTION TRIGGER COMPONENT State: Enabled/ Disabled/ Medication Change/ Elevation Change Display: Always displayed in the GLOBAL NAVIGATION PANEL. The default state is enabled. The state will be disabled if the MEDICATIONS SCREEN is displayed and the color of the icon will change. The state will change to medication change if there is an active DRX MEDICATION MESSAGE CARD PANEL or MEDICATION MESSAGE CARD PANEL in the NOTIFICATIONS SCREEN. The state will change to elevation change if there is a ELEVATION CHANGE MESSAGE CARD PANEL in the HOME SCREEN or NOTIFICATIONS SCREEN. Label: <label.gb.medications> Medications Behavior: Icon and trigger. Upon being triggered, the MEDICATIONS ACTION TRIGGER COMPONENT will change color to the disabled state and conditionally display the MEDICATIONS SCREEN when there are no active DRX MEDICATION MESSAGE CARD PANELS, MEDICATION MESSAGE CARD PANELS, or ELEVATION CHANGE MESSAGE CARD PANELS. Upon triggering in the medication change state, the MEDICATIONS ACTION TRIGGER COMPONENT will display the MEDICATION NOTIFICATION DIALOG PANEL. Upon triggering in the elevation change state, the MEDICATIONS ACTION TRIGGER COMPONENT displays the ELEVATION NOTIFICATION DIALOG PANEL. If there is both a medication intervention and an elevation change, the state is medication change change.	PString152 Medications medications label.gb.medications	UI Info	CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI727	<Picture> MEDICATIONS ACTION TRIGGER COMPONENT (default: enabled)  Medications		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI728	<Picture> MEDICATIONS ACTION TRIGGER COMPONENT (disabled)  Medications		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3420	8.12.1.27.3.3.1 MEDICATION NOTIFICATION DIALOG PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4201	See Shared Elements: Panels.		Info	CloudDRxRTR CloudDRxPS
PAGUI4462	8.12.1.27.3.3.2 ELEVATION NOTIFICATION DIALOG PANEL		Heading	CloudDRxRTR
PAGUI4464	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI332	8.12.1.27.3.5 MORE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI333	Component: MORE ACTION TRIGGER COMPONENT State: Enabled/Disabled Display: Always displayed in the GLOBAL NAVIGATION PANEL. The default state is enabled. The state will be disabled if the MORE SCREEN is displayed and the color of the icon will change. Label: <label.gb.more> More Behavior: Icon and trigger. Upon triggering, the MORE ACTION TRIGGER COMPONENT changes color to the disabled state and displays the MORE SCREEN.	PString1049 More more_label label.gb.more	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI558	<Picture> MORE ACTION TRIGGER COMPONENT (default: enabled)  More		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4222	<Picture> MORE ACTION TRIGGER COMPONENT (disabled)  More		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2329	8.12.1.28 LAB REQUEST MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2330	8.12.1.28.1 LAB REQUEST MESSAGE CARD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2331	The LAB REQUEST MESSAGE CARD PANEL is conditionally displayed if a Lab Request Message is received by the Patient Application at least once since the Earliest History Calendar Date. A unique LAB REQUEST MESSAGE CARD PANEL is		UI Info	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	displayed for every Lab Request Message that is received by the Patient Application. Multiple LAB REQUEST MESSAGE CARD PANELS may be displayed in the same screen or panel. The LAB REQUEST MESSAGE CARD PANEL is in the active state by default. The LAB REQUEST MESSAGE CARD PANEL CARD PANEL changes state from active to complete when the user completes all required actions in the card. The elements in the LAB REQUEST MESSAGE CARD PANEL listed below change behavior along these two states.			CloudDRxPS
PAGUI2332	The LAB REQUEST MESSAGE CARD PANEL is an element of: ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2333	The LAB REQUEST MESSAGE CARD PANEL consists of: LAB REQUEST BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) LAB REQUEST CONTENT LABEL COMPONENT LAB REQUEST INSTRUCTIONS LABEL COMPONENT MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2334	8.12.1.28.2 LAB REQUEST MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2335	<Picture> LAB REQUEST MESSAGE CARD PANEL (default: active, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2336	<Picture> LAB REQUEST MESSAGE CARD PANEL (complete, example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 LAB VISIT FRIDAY, MAR 10 Your clinician requires your regular blood tests when adjusting your medications. Please schedule your lab visit within the next 7 days.			
PAGUI2337	8.12.1.28.3 LAB REQUEST MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2338	8.12.1.28.3.1 LAB REQUEST BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2339	Component: LAB REQUEST BANNER COMPONENT State: Active/ Complete Display: Always displayed in the LAB REQUEST MESSAGE CARD PANEL. The default state is active. In the active state, the LAB REQUEST BANNER COMPONENT displays the active icon and the active background color. In the complete state, the LAB REQUEST BANNER COMPONENT displays the complete icon and complete background color. Label: <label.gb.lab_visit>; Text is all capitalized. Behavior: Icon and text. The LAB REQUEST BANNER COMPONENT changes state with the message card panel it is an element of.	PString71 Lab Visit HLR label.gb.lab_visit	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2340	<Picture> LAB REQUEST BANNER COMPONENT (default: Active) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2341	<Picture> LAB REQUEST BANNER COMPONENT (Complete) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2342	8.12.1.28.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2343	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2344	8.12.1.28.3.3 LAB REQUEST CONTENT LABEL COMPONENT		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



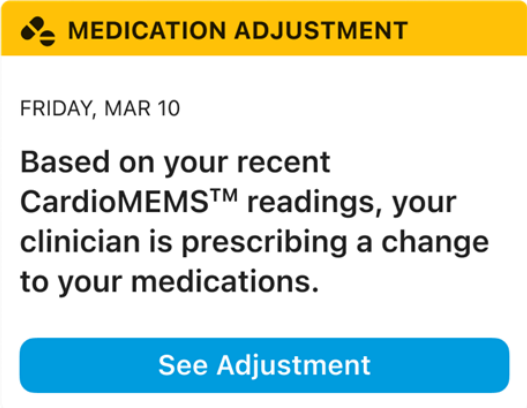
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI2345	Component: LAB REQUEST CONTENT LABEL COMPONENT State: N/A Display: Always displayed in the LAB REQUEST MESSAGE CARD PANEL. Label: <notif.all.lab_visit> Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2346	<Picture> LAB REQUEST CONTENT LABEL COMPONENT Your clinician requires your regular blood tests when adjusting your medications.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4431	8.12.1.28.3.4 LAB REQUEST INSTRUCTIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4432	Component: LAB REQUEST INSTRUCTIONS LABEL COMPONENT State: N/A Display: Always displayed in the LAB REQUEST MESSAGE CARD PANEL. Label: <notif.all.lab_visit_instruction> Behavior: Text.	PString745 Please schedule your lab visit within the next 7 days. HLRBody2 HLRBody notif.all.lab_visit_instruction	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4433	<Picture> LAB REQUEST INSTRUCTIONS LABEL COMPONENT Please schedule your lab visit within the next 7 days .		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2349	8.12.1.28.3.5 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2350	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2169	8.12.1.29 MEDICATION MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2170	8.12.1.29.1 MEDICATION MESSAGE CARD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2171	The MEDICATION MESSAGE CARD PANEL is conditionally displayed if a Medication Message is received by the Patient Application at least once since the Earliest History Calendar Date. A unique MEDICATION MESSAGE CARD PANEL is displayed for every Medication Message that is received by the Patient Application. Only the most recent MEDICATION MESSAGE CARD PANEL in the active state is displayed in the HOME SCREEN. Multiple MEDICATION MESSAGE CARD PANELS		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	may be displayed in the same screen or panel. The MEDICATION MESSAGE CARD PANEL is in the active state by default. The MEDICATION MESSAGE CARD PANEL changes state from active to complete when the user completes all required actions in the card. The MEDICATION MESSAGE CARD PANEL changes state from active to out of date if a more recent MEDICATION MESSAGE CARD PANEL is received before the user is able to complete the required actions, or if the clinician acknowledges the message on the patient's behalf. The elements in the MEDICATION MESSAGE CARD PANEL listed below change behavior along these four states.			
PAGUI2172	The MEDICATION MESSAGE CARD PANEL is an element of: HOME SCREEN ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2173	The MEDICATION MESSAGE CARD PANEL consists of: MEDICATION ADJUSTMENT BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) MEDICATION MESSAGE CONTENT LABEL COMPONENT MEDICATION UPDATE CARD PANEL (See Shared Elements: Components) SEE ADJUSTMENT ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2174	8.12.1.29.2 MEDICATION MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2175	<Picture> MEDICATION MESSAGE CARD PANEL (default: active, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2176	<Picture> MEDICATION MESSAGE CARD PANEL (complete, example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 MEDICATION ADJUSTMENT FRIDAY, MAR 10 Based on your recent CardioMEMS™ readings, your clinician is prescribing a change to your medications. See Adjustment			
PAGUI2241	8.12.1.29.3 MEDICATION MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2245	8.12.1.29.3.1 MEDICATION ADJUSTMENT BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2246	Component: MEDICATION ADJUSTMENT BANNER COMPONENT State: Active/ Complete/ Out of Date / Pre-acknowledged Display: Always displayed in the MEDICATION MESSAGE CARD PANEL. The default state is active. In the active state, the MEDICATION ADJUSTMENT BANNER COMPONENT displays the active icon and the active background color. In the complete state, the MEDICATION ADJUSTMENT BANNER COMPONENT displays the complete icon and complete background color. In the out of date and pre-acknowledged states, the MEDICATION ADJUSTMENT BANNER COMPONENT displays the out of date icon and the complete background color. Label: <label.gb.medication_adj>; Text is all capitalized. Behavior: Icon and text. The MEDICATION ADJUSTMENT BANNER COMPONENT changes state with the message card panel it is an element of the MEDICATION MESSAGE CARD PANEL.	PString264 Medication Adjustment notificationTitleMEDADJ label.gb.medication_adj	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2247	<Picture> MEDICATION ADJUSTMENT BANNER COMPONENT (default: Active) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2248	<Picture> MEDICATION ADJUSTMENT BANNER COMPONENT (Complete) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2249	<Picture> MEDICATION ADJUSTMENT BANNER COMPONENT (Out of Date) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2250	8.12.1.29.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

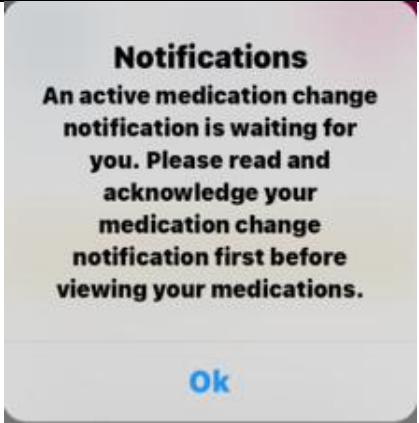
ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI2251	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2252	8.12.1.29.3.3 MEDICATION MESSAGE CONTENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4138	Component: MEDICATION MESSAGE CONTENT LABEL COMPONENT State: Active/ Out of Date/ Complete/ Pre-Acknowledged Display: Conditionally displayed in the MEDICATION MESSAGE CARD PANEL.The default state is Active Label: <notif.manual.med_change> Behavior: Text. In the Active and Out of Date and pre-acknowledged states, the MEDICATION MESSAGE CONTENT LABEL COMPONENT displays <notif.manual.med_change>. The MEDICATION MESSAGE CONTENT LABEL COMPONENT is displayed in the Out of Date (also known as Pre-Acknowledged) state when the patient's clinician has acknowledged the medication change before sending the notification. This scenario could arise when a clinician reviews a medication change directly with the patient and is sending the update to the application.	PString1008 Based on your recent CardioMEMS™ readings, your clinician is prescribing a change to your medications. medicationAdjustmentBody 1 notif.manual.med_change PString1068 Your clinician recommends an adjustment. medicationAdjustmentBody 1 notif.med_change	UI Info	CloudDRxRTR CloudDRxPS
PAGUI2254	<Picture> MEDICATION MESSAGE CONTENT LABEL COMPONENT Based on your recent CardioMEMS™ readings, your clinician is prescribing a change to your medications.	PString1008 Based on your recent CardioMEMS™ readings, your clinician is prescribing a change to your medications. medicationAdjustmentBody 1 notif.manual.med_change PString1064 You had {26} high CardioMEMS™ readings in a row. Your clinician recommends an adjustment. medicationAdjustmentBody 1	UI Design	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		notif.dtp.pa-on.med_change PString1068 Your clinician recommends an adjustment. medicationAdjustmentBody 1 notif.med_change		
PAGUI4440	8.12.1.29.3.4 MEDICATION UPDATE CARD PANEL		Heading	CloudDRxRTR
PAGUI4441	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI2259	8.12.1.29.3.5 SEE ADJUSTMENT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2260	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4189	8.12.1.30 MEDICATION NOTIFICATION DIALOG PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI4190	8.12.1.30.1 MEDICATION NOTIFICATION DIALOG PANEL (Display)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4191	The MEDICATION NOTIFICATION DIALOG PANEL is displayed after the user selects the MEDICATIONS ACTION TRIGGER COMPONENT in the medication change state. The MEDICATION NOTIFICATION DIALOG PANEL is displayed over the screen, does not displace any elements, and any part of the screen that is not covered by the MEDICATION NOTIFICATION DIALOG PANEL is greyed out.			CloudDRxRTR CloudDRxPS
PAGUI4192	The MEDICATION NOTIFICATION DIALOG PANEL is composed of: MEDICATION NOTIFICATION LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4193	8.12.1.30.2 MEDICATION NOTIFICATION DIALOG PANEL (Layout)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4194	<Picture> MEDICATION NOTIFICATION DIALOG PANEL		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

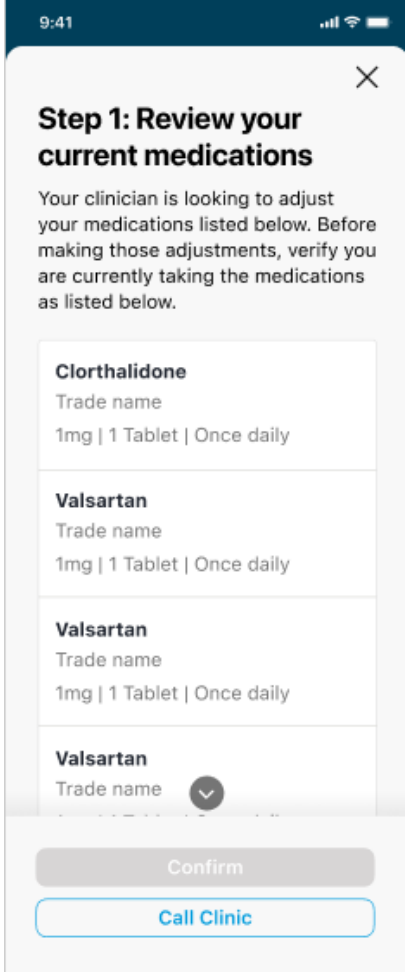
ID		English US Strings	Object Type	Project
				
PAGUI4195	8.12.1.30.3 MEDICATION NOTIFICATION DIALOG PANEL (Elements)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4196	8.12.1.30.3.1 MEDICATION NOTIFICATION LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4197	Component: MEDICATION NOTIFICATION LABEL COMPONENT State: N/A Display: Always displayed as a component of the MEDICATION NOTIFICATION DIALOG PANEL. Label: <msg.gb.medication_change> Behavior: Text.	PString666 You have a message about your medications waiting for you. Please read your message before viewing the Medication Screen. medAckPendingAlert msg.gb.medication_change	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4198	<Picture> MEDICATION NOTIFICATION LABEL COMPONENT 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4199	8.12.1.30.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4200	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI2462	8.12.1.31 MEDICATION REMEDIATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI2463	8.12.1.31.1 MEDICATION REMEDIATION PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2466	The MEDICATION REMEDIATION PANEL is displayed as a workflow after the user selects the SEE ADJUSTMENT ACTION TRIGGER COMPONENT when there is a medication update. The MEDICATION REMEDIATION PANEL is displayed over the entire screen		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2467	The MEDICATION REMEDIATION PANEL consists of: CONFIRM CURRENT MEDICATIONS PANEL MEDICATION UPDATE REMEDIATION PANEL MEDICATION UPDATE COMPLETE PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2464	8.12.1.31.2 MEDICATION REMEDIATION PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR
PAGUI2468	Not applicable - the MEDICATION REMEDIATION PANEL is a series of panels; layout details can be seen in MEDICATION REMEDIATION PANEL (Elements)		UI Info	CloudDRx 1.0 CloudDRxRTR
PAGUI2465	8.12.1.31.3 MEDICATION REMEDIATION PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR
PAGUI2521	8.12.1.31.3.1 CONFIRM CURRENT MEDICATIONS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2522	8.12.1.31.3.1.1 CONFIRM CURRENT MEDICATIONS PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2523	The CONFIRM CURRENT MEDICATIONS PANEL is conditionally displayed as the first panel in the MEDICATION REMEDIATION PANEL. This panel is displayed if there are medications that the patient is currently taking that are changing. The CONFIRM CURRENT MEDICATIONS PANEL is displayed over the entire screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2524	The CONFIRM CURRENT MEDICATIONS PANEL consists of: CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) REVIEW CURRENT MEDICATIONS HEADER COMPONENT CONFIRM MEDICATIONS INSTRUCTIONS LABEL COMPONENT CURRENT MEDICATION INFORMATION LABEL COMPONENT SCROLL DOWN ICON COMPONENT (See Shared Elements: Components) CURRENT MEDICATIONS WARNING BOX CONFIRM ACTION TRIGGER COMPONENT (See Shared Elements: Components) CALL CLINIC ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2525	8.12.1.31.3.1.2 CONFIRM CURRENT MEDICATIONS PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2526	<Picture> CONFIRM CURRENT MEDICATIONS PANEL (example: top of page, more medications to be displayed with scroll)		UI Design	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI4883	<Picture> CONFIRM CURRENT MEDICATIONS PANEL (example: bottom of page displayed with CURRENT MEDICATIONS WARNING BOX)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	9:41 × Step 1: Review your current medications Your clinician is looking to adjust your medications listed below. Before making those adjustments, verify you are currently taking the medications as listed below. Clorthalidone Trade name 1mg 1 Tablet Once daily Valsartan Trade name 1mg 1 Tablet Once daily Valsartan Trade name 1mg 1 Tablet Once daily Valsartan Trade name 1mg 1 Tablet Once daily ⚠️ If you are not taking one or more of these medications as listed, please contact your clinic. Confirm Call Clinic			
PAGUI2528	8.12.1.31.3.1.3 CONFIRM CURRENT MEDICATIONS PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4231	8.12.1.31.3.1.3.2 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4232	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI4390	8.12.1.31.3.1.3.3 REVIEW CURRENT MEDICATIONS HEADER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4391	Component: REVIEW CURRENT MEDICATIONS HEADER COMPONENT State: N/A Display: Always displayed in the CONFIRM CURRENT MEDICATIONS PANEL. Label: <label.gb.current_med_reconciliation> Behavior: Text.	PString1065 Review your current medications reviewYourCurrentMedications label.gb.review_medications PString1090 Step 1: Review your current medications reviewCurrentMeds label.gb.current_med_reconciliation	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4854	<Picture> REVIEW CURRENT MEDICATIONS HEADER COMPONENT Step 1: Review your current medications		UI Design	CloudDRxRTR
PAGUI2534	8.12.1.31.3.1.3.5 CONFIRM MEDICATIONS INSTRUCTIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2535	Component: CONFIRM MEDICATIONS INSTRUCTIONS LABEL COMPONENT State: N/A Display: Always displayed in the CONFIRM CURRENT MEDICATIONS PANEL. Label: <label.gb.currentMedsMessage> Behavior: Text.	PString1055 Confirm you are taking the following medications. If you didn't take one or more medications, please contact your clinic. label.gb.med_reconciliation PString1091 Your clinician is looking to adjust your medications listed below. Before making those adjustments, verify you are currently taking the medications as listed below. currentMedsMessage	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		label.gb.currentMedsMessage		
PAGUI4855	<Picture> CONFIRM MEDICATIONS INSTRUCTIONS LABEL COMPONENT STRING NEEDS CHANGE Your clinician is looking to adjust your medications listed below. Before making those adjustments, verify you are currently taking the medications as listed below.		UI Design	CloudDRxRTR
PAGUI2539	8.12.1.31.3.1.3.6 CURRENT MEDICATION INFORMATION LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2540	Component: CURRENT MEDICATION INFORMATION LABEL COMPONENT State: N/A Display: Always displayed in the CONFIRM CURRENT MEDICATIONS PANEL. Multiple CURRENT MEDICATION INFORMATION LABEL COMPONENT may be displayed in the CONFIRM CURRENT MEDICATIONS PANEL simultaneously. Display the current medication name, trade name, dosage amount, medication format, and number of times the medication is taken a day. Label: N/A - the label is determined by the medication information from the backend Behavior: Scrollable text and values.	PString140 Once daily qd data.Code.Med_Frequency_Cd.Code_Desc.QD PString141 Twice daily bid data.Code.Med_Frequency_Cd.Code_Desc.BID PString142 3 times daily tid data.Code.Med_Frequency_Cd.Code_Desc.TID PString143 4 times daily qid data.Code.Med_Frequency_Cd.Code_Desc.QID PString146 Every night at bedtime qhs data.Code.Med_Frequency_Cd.Code_Desc.QHS PString147 Daily before noon qam	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		data.Code.Med_Frequency _Cd.Code_Desc.QAM PString148 Daily after noon qpm data.Code.Med_Frequency _Cd.Code_Desc.QPM PString154 Tablet tablet data.Code.Dosage_Form_U nit_Cd.Code_Desc.TABLET PString155 Capsule capsule data.Code.Dosage_Form_U nit_Cd.Code_Desc.CAPSULE PString156 IV iv data.Code.Dosage_Form_U nit_Cd.Code_Desc.IV PString157 Injection injection data.Code.Dosage_Form_U nit_Cd.Code_Desc.INJECTIO N PString158 Patch patch data.Code.Dosage_Form_U nit_Cd.Code_Desc.PATCH PString159 Other other data.Code.Dosage_Form_U nit_Cd.Code_Desc.OTHER PString161 Spray spray		

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

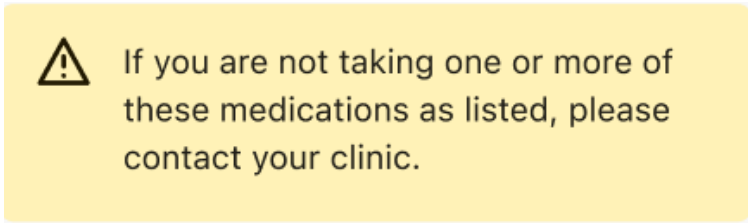


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		data.Code.Dosage_Form_U nit_Cd.Code_Desc.SPRAY		
PAGUI2541	<Picture> CURRENT MEDICATION INFORMATION LABEL COMPONENT (example, single medication) Rivaroxaban Zebeta 30mg 4 Tablets Once daily		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4439	<Picture> CURRENT MEDICATION INFORMATION LABEL COMPONENT (example, multiple medications) Clorthalidone Trade name 30mg 1 Tablet Once daily Furosemide Trade name 30mg 1 Tablet Once daily Telmisartan Zebeta 30mg 2 Tablets Once daily Rivaroxaban Zebeta 30mg 4 Tablets Once daily		UI Design	CloudDRxRTR
PAGUI2542	8.12.1.31.3.1.3.7 SCROLL DOWN ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2543	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4856	8.12.1.31.3.1.3.8 CURRENT MEDICATIONS WARNING BOX		Heading	CloudDRxRTR
PAGUI4857	Component: CURRENT MEDICATIONS WARNING BOX State: N/A Display: Conditionally displayed in the CONFIRM CURRENT MEDICATIONS PANEL when user scrolls to the bottom of the screen.	PString1092 If you are not taking one or more of these medications	UI Info	CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

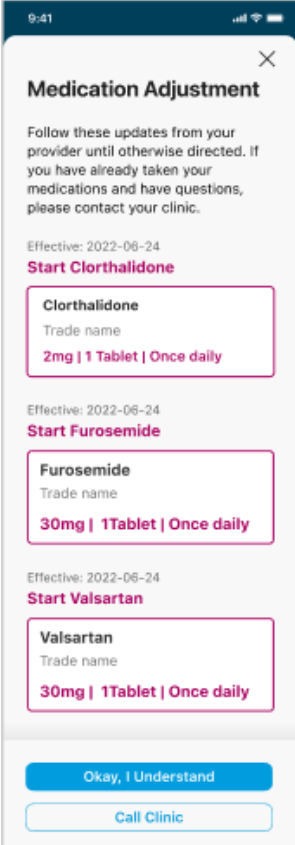
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Label: <label.gb.currentMeds_warning> Behavior: Text and icon.	as listed, please contact your clinic. currentMedsWarning label.gb.currentMeds_warning		
PAGUI4859	<Picture> CURRENT MEDICATIONS WARNING BOX 		UI Design	CloudDRxRTR
PAGUI2544	8.12.1.31.3.1.3.9 CONFIRM ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2547	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4356	8.12.1.31.3.1.3.10 CALL CLINIC ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4357	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI2550	8.12.1.31.3.4 MEDICATION UPDATE REMEDIATION PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI2551	8.12.1.31.3.4.1 MEDICATION UPDATE REMEDIATION PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2552	The MEDICATION UPDATE REMEDIATION PANEL is displayed after the CONFIRM CURRENT MEDICATIONS PANEL as the second panel in the MEDICATION REMEDIATION PANEL. The MEDICATION UPDATE REMEDIATION PANEL is also displayed directly from the MEDICATION MESSAGE CARD PANEL or DRX MEDICATION MESSAGE CARD PANEL when there is no current medications to confirm (i.e. a new medication is being added). The MEDICATION UPDATE REMEDIATION PANEL is displayed over the entire screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2553	The MEDICATION UPDATE REMEDIATION PANEL consists of: BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) MEDICATION UPDATE HEADER COMPONENT MEDICATION RECOMMENDATION LABEL COMPONENT MEDICATION UPDATE CARD PANEL SCROLL DOWN ICON COMPONENT (See Shared Elements: Components) UNDERSTAND ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	CALL CLINIC ACTION TRIGGER COMPONENT (See Shared Elements: Components)			
PAGUI2554	8.12.1.31.3.4.2 MEDICATION UPDATE REMEDIATION PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2555	<Picture> MEDICATION UPDATE REMEDIATION PANEL (example: top of page, need to scroll to view all medication adjustments) 9:41 ← × Step 2 : Medication Adjustment Follow these updates from your provider until otherwise directed. If you have already taken your medications and have questions, please contact your clinic. Effective: 2022-06-24 Increase Clorthalidone Clorthalidone Trade name 4mg Tablet Once daily 2mg 1 Tablet Once daily Effective: 2022-06-24 Start Furosemide Furosemide Trade name 30mg 1Tablet Once daily Effective: 2022-06-24 Stop Valsartan ▼ Okay, I Understand Call Clinic		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
PAGUI4885	<Picture> MEDICATION UPDATE REMEDIATION PANEL (example: only new medication added, no scrolling needed) 		UI Design	CloudDRxRTR
PAGUI2575	<Picture> MEDICATION UPDATE REMEDIATION PANEL (example: bottom of page displayed with all medication adjustments viewed)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
	9:41 ← × Step 2 : Medication Adjustment Follow these updates from your provider until otherwise directed. If you have already taken your medications and have questions, please contact your clinic. Effective: 2022-06-24 Increase Clorthalidone Clorthalidone Trade name 1mg Tablet Once daily 2mg 1 Tablet Once daily Effective: 2022-06-24 Start Furosemide Furosemide Trade name 30mg 1Tablet Once daily Effective: 2022-06-24 Stop Valsartan Valsartan Trade name 1mg 1 Tablet Once daily Okay, I Understand Call Clinic			
PAGUI2556	8.12.1.31.3.4.3 MEDICATION UPDATE REMEDIATION PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2557	8.12.1.31.3.4.3.1 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2558	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4233	8.12.1.31.3.4.3.2 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4234	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI2559	8.12.1.31.3.4.3.3 MEDICATION UPDATE HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2560	Component: MEDICATION UPDATE HEADER COMPONENT State: Medications to Confirm/ New Medication Display: Always displayed in the MEDICATION UPDATE REMEDIATION PANEL. If there are medications to confirm before adjusting, the display will include the step number. If there are no medications to confirm, the display will not include the step number. Label: <label.gb.medication_adj>, <label.gb.medication_adj_noConfirm> Behavior: Text.	PString264 Medication Adjustment notificationTitleMEDADJ label.gb.medication_adj PString1094 Step 2: Medication Adjustment label.gb.medication_adj_noConfirm	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2561	<Picture> MEDICATION UPDATE HEADER COMPONENT (example: New medication) Medication Adjustment		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4886	<Picture> MEDICATION UPDATE HEADER COMPONENT (example: no medications to confirm) Step 2 : Medication Adjustment		UI Design	CloudDRxRTR
PAGUI2562	8.12.1.31.3.4.3.4 MEDICATION RECOMMENDATION LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2563	Component: MEDICATION RECOMMENDATION LABEL COMPONENT State: N/A Display: Conditionally displayed in the MEDICATION UPDATE REMEDIATION PANEL when the panel is scrolled to the top. Label: <label.gb.provider_recomendations> Behavior: Text.	PString601 Follow these updates from your provider until otherwise directed. Refer to the Medications List for your most up to date medications. If you have	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		already taken your medications and have questions, please contact your clinic. medAdjustmentSuccessDescription label.gb.provider_recommendations		
PAGUI2564	<Picture> MEDICATION RECOMMENDATION LABEL COMPONENT Follow these updates from your provider until otherwise directed. Refer to the Medications List for your most up to date medications. If you have already taken your medications and have questions, please contact your clinic.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4435	8.12.1.31.3.4.3.5 MEDICATION UPDATE CARD PANEL		Heading	CloudDRxRTR
PAGUI4434	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI4590	8.12.1.31.3.4.3.6 SCROLL DOWN ICON COMPONENT		Heading	CloudDRxRTR
PAGUI4591	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI2570	8.12.1.31.3.4.3.7 UNDERSTAND ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2571	Component: UNDERSTAND ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Always displayed in the MEDICATION UPDATE REMEDIATION PANEL. The default state is disabled. Label: <label.gb.ok_understand> Behavior: Text. Upon activation in the disabled state, the UNDERSTAND ACTION TRIGGER COMPONENT does nothing. The UNDERSTAND ACTION TRIGGER COMPONENT changes from the disabled state to the enabled state when the bottom MEDICATION UPDATE CARD PANEL is displayed. Upon activation in the enabled state, the UNDERSTAND ACTION TRIGGER COMPONENT displays the MEDICATION UPDATE COMPLETE PANEL.	PString603 Okay, I Understand medAdjustmentSuccessButton1 label.gb.ok_understand	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2572	<Picture> UNDERSTAND ACTION TRIGGER COMPONENT (default: Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2573	<Picture> UNDERSTAND ACTION TRIGGER COMPONENT (Enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2600	8.12.1.31.3.4.3.8 CALL CLINIC ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI2601	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4236	8.12.1.31.3.5 MEDICATION UPDATE COMPLETE PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI4237	8.12.1.31.3.5.1 MEDICATION UPDATE COMPLETE PANEL (Display)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4238	The MEDICATION UPDATE COMPLETE PANEL is displayed after the MEDICATION UPDATE REMEDIATION PANEL as the third panel in the MEDICATION REMEDIATION PANEL. The MEDICATION UPDATE COMPLETE PANEL is displayed over the entire screen.		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4239	The MEDICATION UPDATE COMPLETE PANEL consists of: CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) MEDICATION UPDATE COMPLETE HEADER COMPONENT MEDICATION UPDATE COMPLETE LABEL COMPONENT CHECKMARK ICON COMPONENT MEDICATIONS ACTION TRIGGER COMPONENT DONE ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4240	8.12.1.31.3.5.2 MEDICATION UPDATE COMPLETE PANEL (Layout)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4241	<Picture> MEDICATION UPDATE COMPLETE PANEL		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4244	8.12.1.31.3.5.3 MEDICATION UPDATE COMPLETE PANEL (Elements)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4247	8.12.1.31.3.5.3.2 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI4248	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI4249	8.12.1.31.3.5.3.3 MEDICATION UPDATE COMPLETE HEADER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4250	Component: MEDICATION UPDATE COMPLETE HEADER COMPONENT State: N/A Display: Always displayed in the MEDICATION UPDATE COMPLETE PANEL. Label: <label.gb.med_adj_confirm> Behavior: Text.	PString1058 Your medication adjustment has been confirmed label.gb.med_adj_confirm	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4251	<Picture> MEDICATION UPDATE COMPLETE HEADER COMPONENT Your medication adjustment has been confirmed		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4252	8.12.1.31.3.5.3.4 MEDICATION UPDATE COMPLETE LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4253	Component: MEDICATION UPDATE COMPLETE LABEL COMPONENT State: N/A Display: Always displayed in the MEDICATION UPDATE COMPLETE PANEL. Label: <label.gb.med_adj_complete> Behavior: Text.	PString1059 Medication adjustment complete. Your medication list is now updated. label.gb.med_adj_complete	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4254	<Picture> MEDICATION RECOMMENDATION LABEL COMPONENT Medication adjustment complete. Your medication list is now updated.		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4288	8.12.1.31.3.5.3.5 CHECKMARK ICON COMPONENT		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4289	Component: CHECKMARK ICON COMPONENT State: N/A Display: Always displayed in the MEDICATION UPDATE COMPLETE PANEL. Label: N/A Behavior: Icon.		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4290	<Picture> CHECKMARK ICON COMPONENT		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4281	8.12.1.31.3.5.3.6 MEDICATIONS ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4282	Component: MEDICATIONS ACTION TRIGGER COMPONENT State: Enabled Display: Always displayed in the MEDICATION UPDATE COMPLETE PANEL. Label: <label.gb.go_med> Behavior: Text. Upon activation the MEDICATIONS ACTION TRIGGER COMPONENT displays the MEDICATIONS PANEL in the MEDICATIONS SCREEN.	PString1060 Go to Medications label.gb.go_med	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4287	<Picture> MEDICATIONS ACTION TRIGGER COMPONENT 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4285	8.12.1.31.3.5.3.7 DONE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4286	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI2576	8.12.1.32 MEDICATION UPDATE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2577	8.12.1.32.1 MEDICATION UPDATE CARD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2580	The MEDICATION UPDATE CARD PANEL is displayed as an element of the MEDICATION UPDATE REMEDIATION PANEL. A unique MEDICATION UPDATE CARD PANEL is displayed for every unique medication that is being updated in the MEDICATION REMEDIATION PANEL. The MEDICATION UPDATE CARD PANEL is conditionally displayed as an element of the DRX MEDICATION MESSAGE CARD and MEDICATION MESSAGE CARD when the state is complete. A unique MEDICATION UPDATE CARD PANEL is displayed for every unique medication was updated.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2588	The MEDICATION UPDATE CARD PANEL consists of: EFFECTIVE DATE LABEL COMPONENT MEDICATION UPDATE SUMMARY LABEL COMPONENT MEDICATION UPDATE INFORMATION LABEL COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2578	8.12.1.32.2 MEDICATION UPDATE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI2581	<Picture> MEDICATION UPDATE CARD PANEL (example, active) Effective: 2022-06-24 Increase Clorthalidone Clorthalidone Trade name 1mg Tablet Once daily 2mg 1 Tablet Once daily		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4436	<Picture> MEDICATION UPDATE CARD PANEL (example, complete) Effective: 2022-03-10 Increase Furosemide Furosemide Lasix 1mg Tablet Once daily 2mg 1 Tablet Once daily		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2579	8.12.1.32.3 MEDICATION UPDATE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2582	8.12.1.32.3.1 EFFECTIVE DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2583	Component: EFFECTIVE DATE LABEL COMPONENT State: N/A Display: Always displayed in the MEDICATION UPDATE CARD PANEL. The EFFECTIVE DATE LABEL COMPONENT displays the date the patient is recommended to enact the medication update. Label: <label.gb.med_update_effective_on> Behavior: Text and value.	PString171 Effective: {1} effective label.gb.med_update_effec tive_on	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2584	<Picture> EFFECTIVE DATE LABEL COMPONENT (example)		UI Design	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Effective: 2022-06-24			CloudDRxRTR CloudDRxPS
PAGUI2585	8.12.1.32.3.2 MEDICATION UPDATE SUMMARY LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2586	Component: MEDICATION UPDATE SUMMARY LABEL COMPONENT State: Active/Complete Display: Always displayed in the MEDICATION UPDATE CARD PANEL. Label: <label.gb.start>, <label.gb.stop>, <label.gb.change>, <label.gb.increase>, <label.gb.decrease>, <label.gb.hold>, <label.gb.discontinue> Behavior: Text and value. The UPDATE SUMMARY LABEL COMPONENT changes based on what kind of medication update is occurring and the name of the effected medication.	PString117 Decrease decrease label.gb.decrease PString118 Increase increase label.gb.increase PString125 Hold hold label.gb.hold PString127 Start start label.gb.start PString128 Change change label.gb.change PString699 Stop stop label.gb.stop	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2587	<Picture> MEDICATION UPDATE SUMMARY LABEL COMPONENT (example: change, active) Change Furosemide		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2589	<Picture> MEDICATION UPDATE SUMMARY LABEL COMPONENT (example: start, active) Start Clorthalidone		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2590	<Picture>MEDICATION UPDATE SUMMARY LABEL COMPONENT (example: increase, active) Increase Clorthalidone		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4438	<Picture> MEDICATION UPDATE SUMMARY LABEL COMPONENT (example: increase, complete) Increase Furosemide		UI Design	CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2591	8.12.1.32.3.3 MEDICATION UPDATE INFORMATION LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2592	Component: MEDICATION UPDATE INFORMATION LABEL COMPONENT State: Active/Complete Display: Always displayed in the MEDICATION UPDATE CARD PANEL. The MEDICATION UPDATE INFORMATION LABEL COMPONENT displays the medication generic name, brand name, original dosage, original format, original frequency, new dosage, new format, and new frequency. The original dosage, format, and frequency are crossed out to indicate that they are being changed. In the active state, the box around the information and the new medication information are displayed in an active color. In the complete state, they are displayed in the complete color. Label: N/A - the label is pulled from medication information in the backend. Behavior: Text and values.	PString140 Once daily qd data.Code.Med_Frequency _Cd.Code_Desc.QD PString141 Twice daily bid data.Code.Med_Frequency _Cd.Code_Desc.BID PString142 3 times daily tid data.Code.Med_Frequency _Cd.Code_Desc.TID PString143 4 times daily qid data.Code.Med_Frequency _Cd.Code_Desc.QID PString146 Every night at bedtime qhs data.Code.Med_Frequency _Cd.Code_Desc.QHS PString147 Daily before noon qam data.Code.Med_Frequency _Cd.Code_Desc.QAM PString148 Daily after noon qpm data.Code.Med_Frequency _Cd.Code_Desc.QPM PString154 Tablet tablet	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

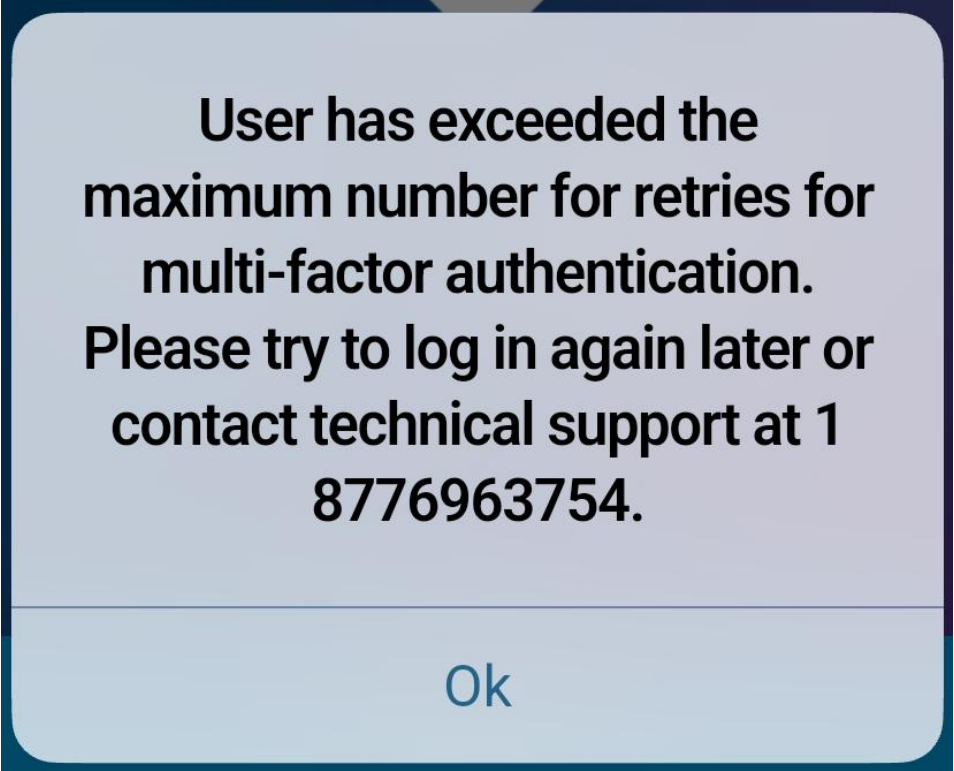
ID		English US Strings	Object Type	Project
		data.Code.Dosage_Form_Unit_Cd.Code_Desc.TABLET PString155 Capsule capsule data.Code.Dosage_Form_Unit_Cd.Code_Desc.CAPSULE PString156 IV iv data.Code.Dosage_Form_Unit_Cd.Code_Desc.IV PString157 Injection injection data.Code.Dosage_Form_Unit_Cd.Code_Desc.INJECTION PString158 Patch patch data.Code.Dosage_Form_Unit_Cd.Code_Desc.PATCH PString159 Other other data.Code.Dosage_Form_Unit_Cd.Code_Desc.OTHER PString161 Spray spray data.Code.Dosage_Form_Unit_Cd.Code_Desc.SPRAY		
PAGUI2593	<Picture> MEDICATION UPDATE INFORMATION LABEL COMPONENT (example, active) Clorthalidone Trade name 1mg Tablet Once daily 2mg 1 Tablet Once daily		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4437	<Picture> MEDICATION UPDATE INFORMATION LABEL COMPONENT (example, complete) Furosemide Lasix 1mg Tablet Once daily 2mg 1 Tablet Once daily		UI Design	CloudDRxRTR
PAGUI3625	8.12.1.33 MFA MAXIMUM ATTEMPTS DIALOG PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3626	8.12.1.33.1 MFA MAXIMUM ATTEMPTS DIALOG PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3627	The MFA MAXIMUM ATTEMPTS DIALOG PANEL is displayed when the user attempts to verify the MFA code in the TWO FACTOR AUTHENTION PANEL too many times (three times incorrectly).		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3628	The MFA MAXIMUM ATTEMPTS DIALOG PANEL consists of: MFA MAXIMUM ATTEMPTS INSTRUCTIONS COMPONENT OKAY ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3629	8.12.1.33.2 MFA MAXIMUM ATTEMPTS DIALOG PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3630	<Picture> MFA MAXIMUM ATTEMPTS DIALOG PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

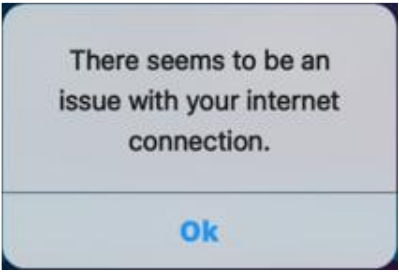
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI3631	8.12.1.33.3 MFA MAXIMUM ATTEMPTS DIALOG PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3632	8.12.1.33.3.1 MFA MAXIMUM ATTEMPTS INSTRUCTIONS COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3633	Component: MFA MAXIMUM ATTEMPTS INSTRUCTIONS COMPONENT State: N/A Display: Always displayed as a component of the MFA MAXIMUM ATTEMPTS DIALOG PANEL. Label: <msg.gb.mfa_max_attempts> Behavior: Text.	PString687 User has exceeded the maximum number for retries for multi-factor authentication. Please try to login again later. loginMaximumNumberOfRetries msg.gb.mfa_max_attempts PString1026	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

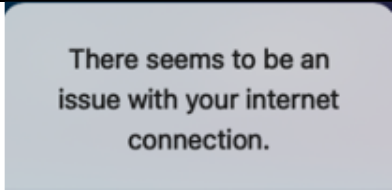
ID		English US Strings	Object Type	Project
		User has exceeded the maximum number for retries for multi-factor authentication. Please try to log in again later or contact technical support at {16}. loginMaximumNumberOfRetries msg.gb.mfa_max_attempts		
PAGUI3634	<Picture> MFA MAXIMUM ATTEMPTS INSTRUCTIONS COMPONENT User has exceeded the maximum number for retries for multi-factor authentication. Please try to log in again later or contact technical support at 1 8776963754.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3635	8.12.1.33.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3636	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3537	8.12.1.35 NO INTERNET CONNECTION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3538	8.12.1.35.1 NO INTERNET CONNECTION PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3539	The NO INTERNET CONNECTION PANEL is displayed when the Patient Application is unable to connect to the internet while the user is actively using the Patient Application. The NO INTERNET CONNECTION PANEL is displayed over the screen, does		UI Info	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	not displace any elements, and any part of the screen that is not covered by the NO INTERNET CONNECTION PANEL is greyed out.			
PAGUI3552	The NO INTERNET CONNECTION PANEL is an element of: STARTUP SCREEN POLICY SCREEN REGISTRATION SCREEN HOME SCREEN READINGS SCREEN MEDICATIONS SCREEN JOURNAL SCREEN NOTIFICATIONS SCREEN MORE SCREEN		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3540	The NO INTERNET CONNECTION PANEL is composed of: NO INTERNET CONNECTION LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3541	8.12.1.35.2 NO INTERNET CONNECTION PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3542	<Picture> NO INTERNET CONNECTION PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3543	8.12.1.35.3 NO INTERNET CONNECTION PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3544	8.12.1.35.3.1 NO INTERNET CONNECTION LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3545	Component: NO INTERNET CONNECTION LABEL COMPONENT State: N/A Display: Always displayed as a component of the NO INTERNET CONNECTION PANEL. Label: <msg.gb.no_internet> Behavior: Text.	PString676 There seems to be an issue with your internet connection. noInternetConnection msg.gb.no_internet	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3546	<Picture> NO INTERNET CONNECTION LABEL COMPONENT		UI Design	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3550	8.12.1.35.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3551	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI123	8.12.1.36 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI124	8.12.1.36.1 OS STATUS PANEL (DISPLAY)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI347	The OS STATUS PANEL displays an overview of key device information as the application is running.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI348	The OS STATUS PANEL is an element of: STARTUP SCREEN HOME SCREEN READINGS SCREEN MEDICATIONS SCREEN NOTIFICATIONS SCREEN MORE SCREEN POLICY SCREEN PROFILE PANEL SETTINGS PANEL EDUCATION PANEL ABOUT HEART FAILURE PANEL ABOUT HEART FAILURE SYSTEM PANEL LEARNING PANEL IN-APP TRAINING PANEL SUPPORT PANEL FAQ PANEL CONTACT PANEL ABOUT PANEL SHORT PRIVACY PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



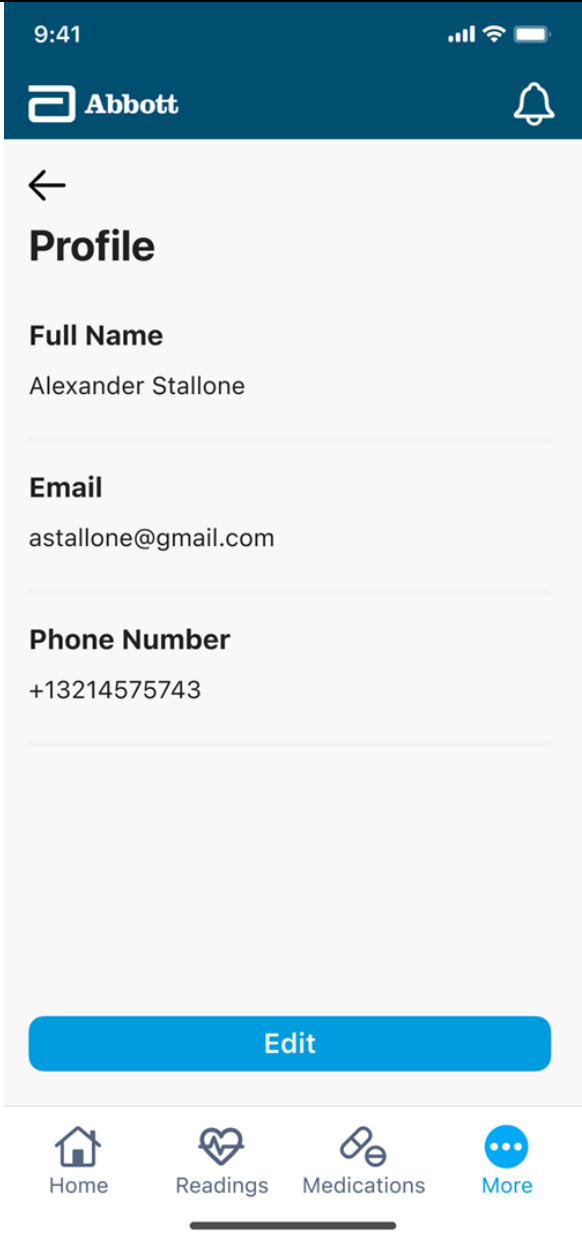
myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
	CONTACT PANEL USER MANUAL WEB VIEW			
PAGUI349	The OS STATUS PANEL is composed of: OS header elements are determined by SOUP provider. Panel components may vary based on language and OS level settings.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI125	8.12.1.36.2 OS STATUS PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI350	<Picture> OS STATUS PANEL (Apple device example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI351	<Picture> OS STATUS PANEL (Android device example) TBD		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI126	8.12.1.36.3 OS STATUS PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI352	OS STATUS PANEL elements are determined by the SOUP provider and the layout only serves as an example of what the panel may look like.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI580	8.12.1.41 PROFILE PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI934	8.12.1.41.1 PROFILE PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI937	The PROFILE PANEL is displayed after the user navigates to it from the MORE SCREEN.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI938	The PROFILE PANEL consists of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) PROFILE HEADER COMPONENT FULL NAME LABEL COMPONENT (See Shared Elements: Components) FULL NAME DISPLAY COMPONENT EMAIL LABEL COMPONENT (See Shared Elements: Components) EMAIL DISPLAY COMPONENT PHONE LABEL COMPONENT (See Shared Elements: Components) PHONE DISPLAY COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	EDIT ACCOUNT INFO ACTION TRIGGER COMPONENT EDIT ACCOUNT INFO PANEL			
PAGUI935	8.12.1.41.2 PROFILE PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI939	<Picture> PROFILE PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI936	8.12.1.41.3 PROFILE PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1632	8.12.1.41.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1633	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4679	8.12.1.41.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4680	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI941	8.12.1.41.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI942	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI945	8.12.1.41.3.4 PROFILE HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI946	Component: PROFILE HEADER COMPONENT State: N/A Display: Always displayed in the PROFILE PANEL. Label: <label.gb.profile> Behavior: Text.	PString205 Profile profile label.gb.profile	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI947	<Picture> PROFILE HEADER COMPONENT Profile		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI961	8.12.1.41.3.5 FULL NAME LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI962	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI963	8.12.1.41.3.6 FULL NAME DISPLAY COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI964	Component: FULL NAME DISPLAY COMPONENT State: N/A Display: Always displayed in the PROFILE PANEL. Label: N/A - content determined by content manager and user input Behavior: Text. The FULL NAME DISPLAY COMPONENT displays the name associated with the user.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI965	<Picture> FULL NAME DISPLAY COMPONENT (example)		UI Design	CloudDRx 1.0

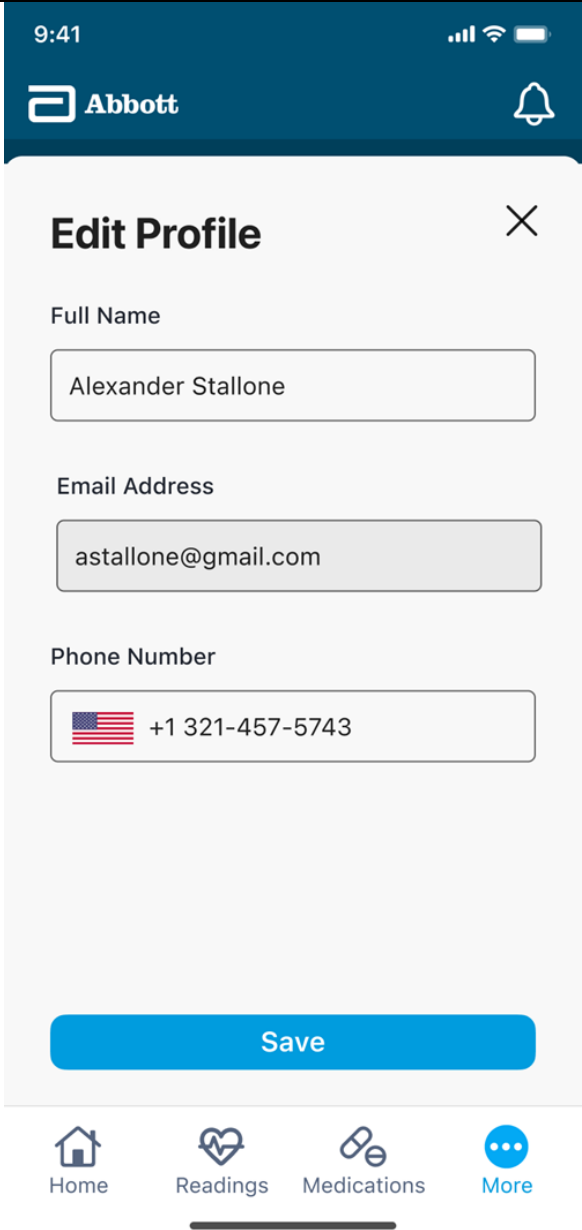
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

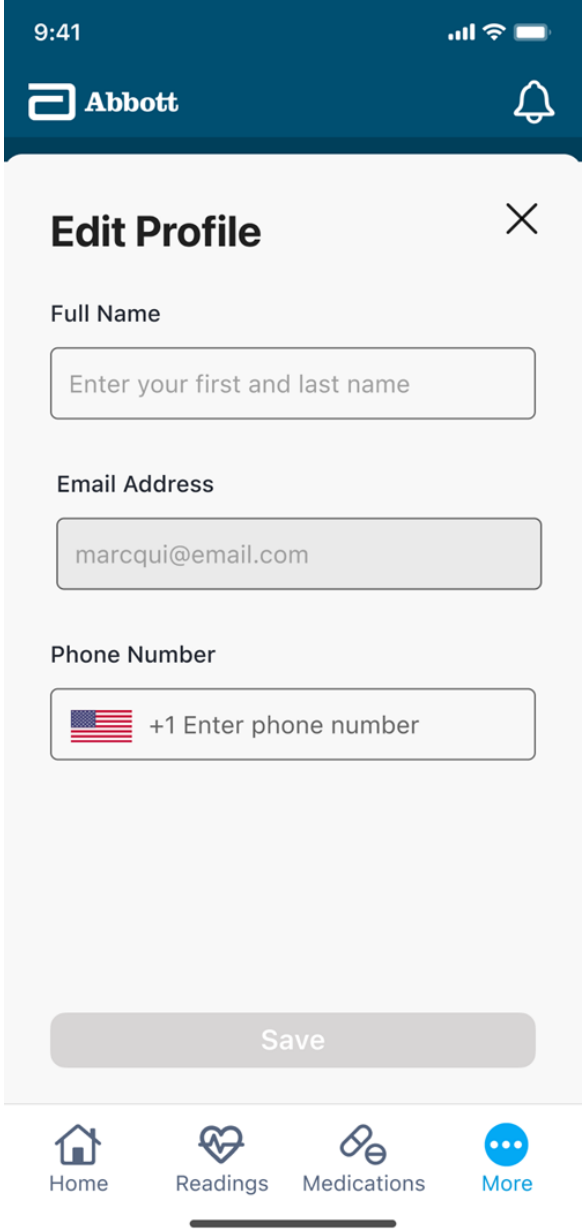
ID		English US Strings	Object Type	Project
	Alexandar Stallone			CloudDRxRTR CloudDRxPS
PAGUI966	8.12.1.41.3.7 EMAIL LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI967	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI968	8.12.1.41.3.8 EMAIL DISPLAY COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI969	Component: EMAIL DISPLAY COMPONENT State: N/A Display: Always displayed in the PROFILE PANEL. Label: N/A - content determined by content manager and user input Behavior: Text. The EMAIL DISPLAY COMPONENT displays the email associated with the user.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI970	<Picture> EMAIL DISPLAY COMPONENT (example) astallone@gmail.com		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI971	8.12.1.41.3.9 PHONE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI972	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI973	8.12.1.41.3.10 PHONE DISPLAY COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI974	Component: PHONE DISPLAY COMPONENT State: N/A Display: Always displayed in the PROFILE PANEL. Label: N/A - content determined by content manager and user input Behavior: Text. The PHONE DISPLAY COMPONENT displays the phone number associated with the user. If there is no number associated with the user, the display is left blank.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI975	<Picture> PHONE DISPLAY COMPONENT (example) (289) 534-5282		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI958	8.12.1.41.3.11 EDIT ACCOUNT INFO ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI959	Component: EDIT ACCOUNT INFO ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the PROFILE PANEL. Label: <label.gb.edit> Behavior: Text and trigger. Upon activation, display the EDIT ACCOUNT INFO PANEL.	PString613 Edit edit label.gb.edit	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI960	<Picture> EDIT ACCOUNT INFO ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI979	8.12.1.41.3.12 EDIT ACCOUNT INFO PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI980	8.12.1.41.3.12.1 EDIT ACCOUNT INFO PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI983	The EDIT ACCOUNT INFO PANEL is conditionally displayed in the PROFILE PANEL. The EDIT ACCOUNT INFO PANEL is displayed as an overlay over the entire screen except the OS STATUS PANEL and GLOBAL HEADING PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI984	The EDIT ACCOUNT INFO PANEL consists of: CLOSE PANEL TRIGGER COMPONENT (See Shared Elements: Components) EDIT ACCOUNT INFO HEADER COMPONENT FULL NAME LABEL COMPONENT (See Shared Elements: Components) FULL NAME INPUT COMPONENT (See Shared Elements: Components) EMAIL LABEL COMPONENT (See Shared Elements: Components) EMAIL INPUT COMPONENT (See Shared Elements: Components) PHONE LABEL COMPONENT (See Shared Elements: Components) PHONE INPUT COMPONENT (See Shared Elements: Components) ALPHABETICAL KEYBOARD COMPONENT (See Shared Elements: Components) NUMERIC KEYBOARD COMPONENT (See Shared Elements: Components) SAVE ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI981	8.12.1.41.3.12.2 EDIT ACCOUNT INFO PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4060	<Picture> EDIT ACCOUNT INFO PANEL (default: no changes)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
				
PAGUI4708	<Picture> EDIT ACCOUNT INFO PANEL (example: blank inputs)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI982	8.12.1.41.3.12.3 EDIT ACCOUNT INFO PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI986	8.12.1.41.3.12.3.1 CLOSE PANEL TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI987	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI988	8.12.1.41.3.12.3.2 EDIT ACCOUNT INFO HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1643	Component: EDIT ACCOUNT INFO HEADER COMPONENT State: N/A Display: Always displayed in the EDIT ACCOUNT INFO PANEL. Label: <label.gb.edit_profile> Behavior: Text.	PString614 Edit Profile editProfile label.gb.edit_profile	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1644	<Picture> EDIT ACCOUNT INFO HEADER COMPONENT Edit Profile		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI990	8.12.1.41.3.12.3.4 FULL NAME LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI991	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI992	8.12.1.41.3.12.3.5 FULL NAME INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI993	Component: FULL NAME INPUT COMPONENT State: No Input/ Input Display: Always displayed in the EDIT ACCOUNT INFO PANEL. The default status is no input. Label: <label.gb.enter_name>, <label.gb.enter_first_last> Behavior: Text and input box. In the no input state, display the default text. In the input state, display the input given from the ALPHABETICAL KEYBOARD COMPONENT or the currently saved name associated with the user.	PString615 Enter your first and last name nameValidation label.gb.enter_first_last	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4058	Component: FULL NAME INPUT COMPONENT State: No Input/ Input Display: Always displayed in the EDIT ACCOUNT INFO PANEL. The default status is no input. Label: <label.gb.enter_first_last>	PString615 Enter your first and last name nameValidation label.gb.enter_first_last	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Behavior: Text and input box. In the no input state, display the default text. In the input state, display the input given from the ALPHABETICAL KEYBOARD COMPONENT or the currently saved name associated with the user.			
PAGUI994	<Picture> FULL NAME INPUT COMPONENT (default: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI995	<Picture> FULL NAME INPUT COMPONENT (example: input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI996	8.12.1.41.3.12.3.6 EMAIL LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI997	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI998	8.12.1.41.3.12.3.7 EMAIL INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1657	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1002	8.12.1.41.3.12.3.8 PHONE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1003	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1004	8.12.1.41.3.12.3.9 PHONE INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1651	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1008	8.12.1.41.3.12.3.11 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1009	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

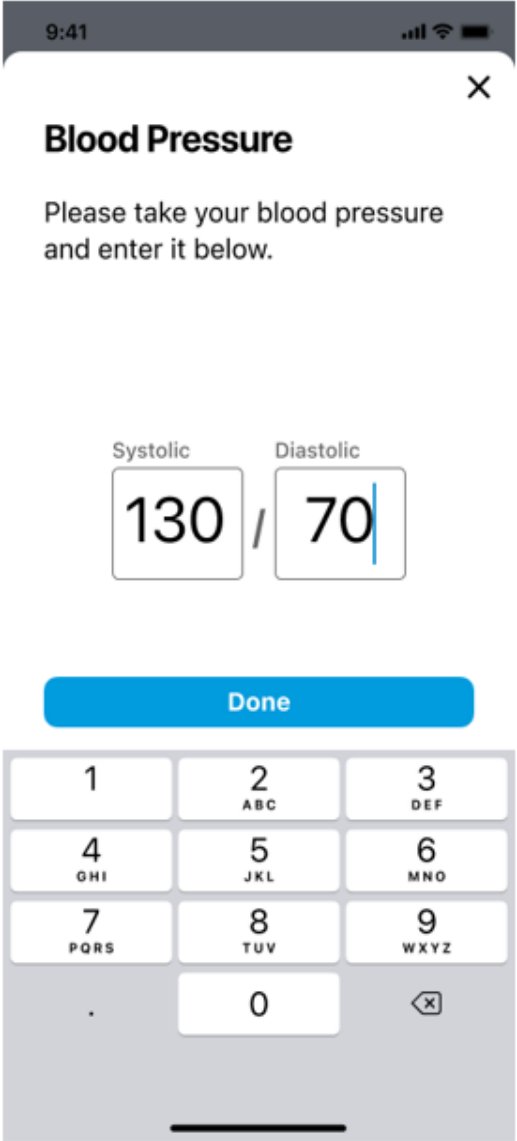
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4681	8.12.1.41.3.12.3.12 NUMERIC KEYBOARD COMPONENT		Heading	CloudDRxRTR
PAGUI4682	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI1056	8.12.1.41.3.12.3.13 SAVE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1057	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3210	8.12.1.42 RECORD BP READING PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3211	8.12.1.42.1 RECORD BP READING PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3212	The RECORD BP READING PANEL is conditionally displayed after the user selects the RECORD BP READING ACTION TRIGGER COMPONENT. The RECORD BP READING PANEL is displayed over the entire screen except the OS STATUS PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3235	The RECORD BP READING PANEL is an element of: HOME SCREEN READINGS SCREEN		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3213	The RECORD BP READING PANEL consists of: CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) RECORD BP READING HEADER COMPONENT RECORD BP READING INSTRUCTION LABEL COMPONENT RECORD SYSTOLIC READING INPUT COMPONENT FRACTION ICON COMPONENT RECORD DIASTOLIC READING INPUT COMPONENT LOW SYSTOLIC ERROR LABEL COMPONENT DONE ACTION TRIGGER COMPONENT (See Shared Elements: Components) NUMERIC KEYBOARD COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3214	8.12.1.42.2 RECORD BP READING PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3215	<Picture> RECORD BP READING PANEL (example: Input)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI3216	8.12.1.42.3 RECORD BP READING PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3217	8.12.1.42.3.1 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3218	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3219	8.12.1.42.3.2 RECORD BP READING HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3220	Component: RECORD BP READING HEADER COMPONENT State: N/A Display: Always displayed as a component of the RECORD BP READING PANEL. Label: <label.gb.blood_pressure> Behavior: Text.	PString616 Blood Pressure bloodPressure label.gb.blood_pressure	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3221	<Picture> RECORD BP READING HEADER COMPONENT Blood Pressure		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3222	8.12.1.42.3.3 RECORD BP READING INSTRUCTION COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3223	Component: RECORD BP READING INSTRUCTIOIN COMPONENT State: N/A Display: Always displayed as a component of the RECORD BP READING PANEL. Label: <label.gb.record_blood_pressure>, <label.bp_record>, <label.gb.bp_enter> Behavior: Text.	PString617 Please take your blood pressure and enter it below. recordBloodPressure label.gb.record_blood_pressure	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3224	<Picture> RECORD BP READING INSTRUCTION COMPONENT Please take your blood pressure and enter it below.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3225	8.12.1.42.3.4 RECORD SYSTOLIC READING INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3226	Component: RECORD SYSTOLIC READING INPUT COMPONENT State: Enabled/ Disabled Display: Always displayed as a component of the RECORD BP READING PANEL. The default state is disabled. Label: <label.gb.bp_systolic> Behavior: Input display. When enabled, display the input from the NUMERIC KEYBOARD PANEL and a cursor within the input box. When disabled, display an empty input box. Always display blood pressure name.	PString618 Systolic systolic label.gb.bp_systolic	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3227	<Picture> RECORD BP READING INPUT COMPONENT (default: Disabled)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

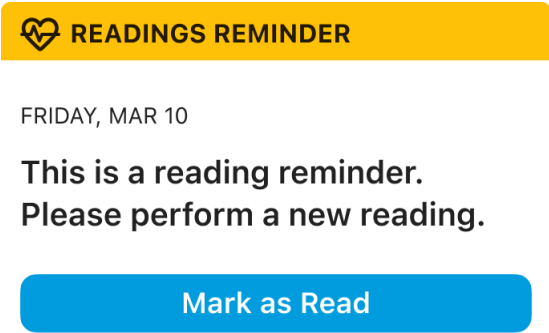
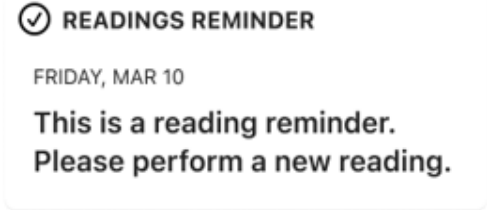
ID		English US Strings	Object Type	Project
				
PAGUI3228	<Picture> RECORD BP READING INPUT COMPONENT (example: Enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3236	8.12.1.42.3.5 FRACTION ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3237	Component: FRACTION ICON COMPONENT State: N/A Display: Always displayed as a component of the RECORD BP READING PANEL. Label: N/A Behavior: Icon.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3238	<Picture> FRACTION ICON COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3241	8.12.1.42.3.6 RECORD DIASTOLIC READING INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3242	Component: RECORD DIASTOLIC READING INPUT COMPONENT State: Enabled/ Disabled Display: Always displayed as a component of the RECORD BP READING PANEL. The default state is disabled. Label: <label.gb.bp_diastolic> Behavior: Input display. When enabled, display the input from the NUMERIC KEYBOARD PANEL and a cursor within the input box. When disabled, display an empty input box. Always display blood pressure name.	PString619 Diastolic diastolic label.gb.bp_diastolic	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3243	<Picture> RECORD BP READING INPUT COMPONENT (default: Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3244	<Picture> RECORD BP READING INPUT COMPONENT (example: Enabled)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Diastolic 70			
PAGUI3607	8.12.1.42.3.7 LOW SYSTOLIC ERROR LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3608	Component: LOW SYSTOLIC ERROR LABEL COMPONENT State: N/A Display: Conditionally displayed as a component of the RECORD BP READING PANEL after the user selects the DONE ACTION TRIGGER COMPONENT if the recorded systolic blood pressure is less than the recorded diastolic blood pressure. Label: <label.gb.low_systolic> Behavior: Text.	PString732 Systolic lower than Diastolic bpValidationError label.gb.low_systolic	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3609	<Picture> LOW SYSTOLIC ERROR LABEL COMPONENT Systolic lower than Diastolic		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3229	8.12.1.42.3.8 DONE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3230	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3233	8.12.1.42.3.10 NUMERIC KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3234	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2425	8.12.1.43 READING REMINDER MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2426	8.12.1.43.1 READING REMINDER MESSAGE CARD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2427	The READING REMINDER MESSAGE CARD PANEL is conditionally displayed if a Reading Reminder Message is received by the Patient Application at least once since the Earliest History Calendar Date. A unique READING REMINDER MESSAGE CARD PANEL is displayed for every Reading Reminder Message that is received by the Patient Application. Multiple READING REMINDER MESSAGE CARD PANELs may be displayed in the same screen or panel. The READING REMINDER MESSAGE CARD PANEL is in the active state by default. The READING REMINDER MESSAGE CARD PANEL CARD PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

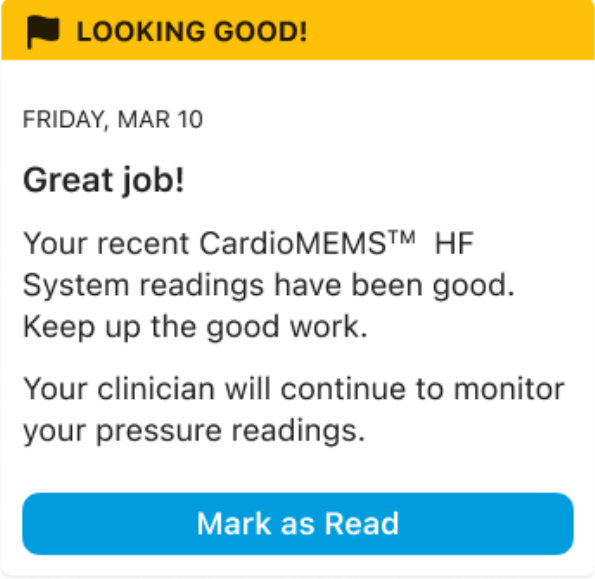
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	changes state from active to complete when the user completes all required actions in the card. The elements in the READING REMINDER MESSAGE CARD PANEL listed below change behavior along these two states.			
PAGUI2428	The READING REMINDER MESSAGE CARD PANEL is an element of: ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2429	The READING REMINDER MESSAGE CARD PANEL consists of: READING REMINDER BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) READING REMINDER INSTRUCTIONS LABEL COMPONENT MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2430	8.12.1.43.2 READING REMINDER MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2431	<Picture> READING REMINDER MESSAGE CARD PANEL (default: active, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2432	<Picture> READING REMINDER MESSAGE CARD PANEL (complete, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2433	8.12.1.43.3 READING REMINDER MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2434	8.12.1.43.3.1 READING REMINDER BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2435	Component: READING REMINDER BANNER COMPONENT State: Active/ Complete	PString74 Readings Reminder	UI Info	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Display: Always displayed in the READING REMINDER MESSAGE CARD PANEL. The default state is active. In the active state, the READING REMINDER BANNER COMPONENT displays the active icon and the active background color. In the complete state, the READING REMINDER BANNER COMPONENT displays the complete icon and complete background color. Label: <label.gb.readings_reminder>; Text is all capitalized. Behavior: Icon and text. The READING REMINDER BANNER COMPONENT changes state with the message card panel it is an element of.	HSN label.gb.readings_reminder		CloudDRxPS
PAGUI2436	<Picture> READING REMINDER BANNER COMPONENT (default: Active) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2437	<Picture> READING REMINDER BANNER COMPONENT (Complete) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2438	8.12.1.43.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2439	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2443	8.12.1.43.3.3 READING REMINDER INSTRUCTIONS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2444	Component: READING REMINDER INSTRUCTIONS LABEL COMPONENT State: N/A Display: Always displayed in the READING REMINDER MESSAGE CARD PANEL. Label: <notif.all.reading_reminder> Behavior: Text.	PString742 This is a reading reminder. Please perform a new reading. HSNBody1 notif.all.reading_reminder	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2445	<Picture> READING REMINDER INSTRUCTIONS LABEL COMPONENT This is a reading reminder. Please perform a new reading.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2446	8.12.1.43.3.4 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2447	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2377	8.12.1.44 READINGS REVIEWED MESSAGE CARD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2378	8.12.1.44.1 READINGS REVIEWED MESSAGE CARD PANEL (Dispaly)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2379	The READINGS REVIEWED MESSAGE CARD PANEL is conditionally displayed if a Readings Reviewed is received by the Patient Application at least once since the Earliest History Calendar Date. A unique READINGS REVIEWED MESSAGE CARD PANEL is displayed for every Readings Reviewed Message that is received by the Patient Application. Multiple READINGS REVIEWED MESSAGE CARD PANELS may be displayed in the same screen or panel. The READINGS REVIEWED MESSAGE CARD PANEL is in the active state by default. The READINGS REVIEWED MESSAGE CARD PANEL CARD PANEL changes state from active to complete when the user completes all required actions in the card. The elements in the READINGS REVIEWED MESSAGE CARD PANEL listed below change behavior between these two states.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2380	The READINGS REVIEWED MESSAGE CARD PANEL is an element of: ACTIVE NOTIFICATIONS PANEL COMPLETE NOTIFICATIONS PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2381	The READINGS REVIEWED MESSAGE CARD PANEL consists of: READINGS REVIEWED BANNER COMPONENT NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components) READINGS REVIEWED CONTENT HEADER LABEL COMPONENT READINGS REVIEWED CONTENT LABEL COMPONENT MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2382	8.12.1.44.2 READINGS REVIEWED MESSAGE CARD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2383	<Picture> READINGS REVIEWED MESSAGE CARD PANEL (default: active, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2384	<Picture> READINGS REVIEWED MESSAGE CARD PANEL (complete, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2385	8.12.1.44.3 READINGS REVIEWED MESSAGE CARD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2386	8.12.1.44.3.1 READINGS REVIEWED BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2387	Component: READINGS REVIEWED BANNER COMPONENT State: Active/ Complete Display: Always displayed in the READINGS REVIEWED MESSAGE CARD PANEL. The default state is active. In the active state, the READINGS REVIEWED BANNER COMPONENT displays the active icon and the active background color. In the complete state, the READINGS REVIEWED BANNER COMPONENT displays the complete icon and complete background color. Label: <label.gb.readings_reviewed>; Text is all capitalized. Behavior: Icon and text. The READINGS REVIEWED BANNER COMPONENT changes state with the message card panel it is an element of.	PString72 Looking Good! HDR label.gb.readings_reviewed	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2388	<Picture> READINGS REVIEWED BANNER COMPONENT (default: Active) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2389	<Picture> READINGS REVIEWED BANNER COMPONENT (Complete) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2390	8.12.1.44.3.2 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

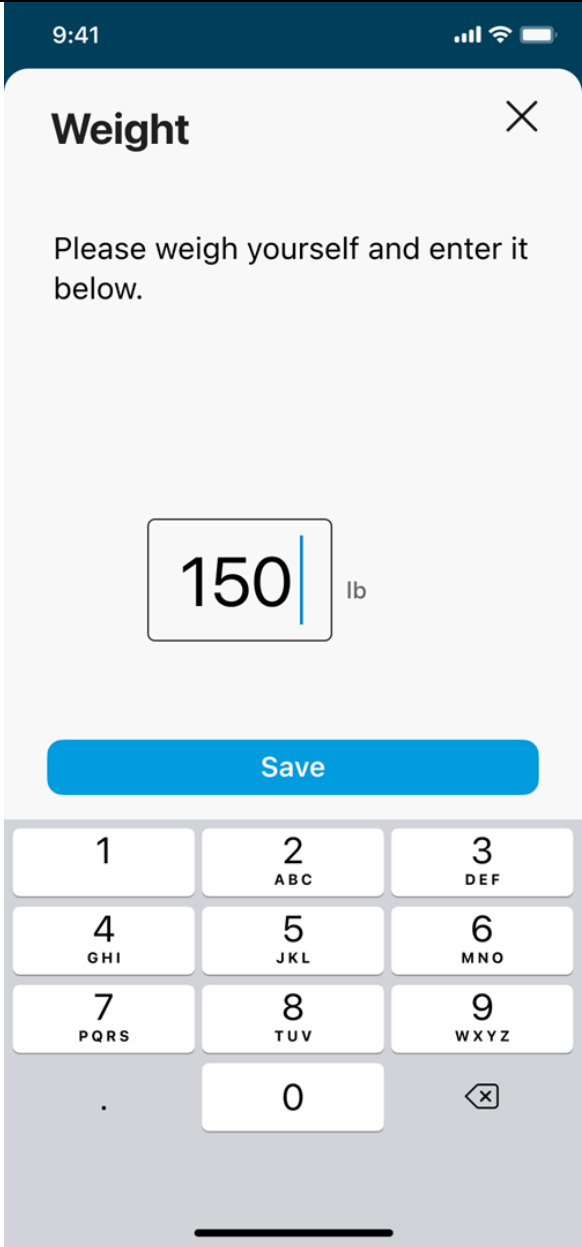
ID		English US Strings	Object Type	Project
PAGUI2391	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4423	8.12.1.44.3.3 READINGS REVIEWED CONTENT HEADER LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4424	Component: READINGS REVIEWED CONTENT HEADER LABEL COMPONENT State: N/A Display: Always displayed in the READINGS REVIEWED MESSAGE CARD PANEL. Label: <notif.standard.great_job> Behavior: Text.	PString740 Great job! HDRBody1 notif.standard.great_job	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4425	<Picture> READINGS REVIEWED CONTENT HEADER LABEL COMPONENT Great job!		UI Design	CloudDRxRTR
PAGUI2392	8.12.1.44.3.4 READINGS REVIEWED CONTENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2393	Component: READINGS REVIEWED CONTENT LABEL COMPONENT State: N/A Display: Always displayed in the READINGS REVIEWED MESSAGE CARD PANEL. Label: <notif.standard.looking_good_message> Behavior: Text.	PString747 Your recent CardioMEMS™ HF System readings have been good. Keep up the good work. We will continue to monitor your pressure readings. HDRBody2 notif.standard.looking_goo d_message PString1007 Your recent CardioMEMS™ HF System readings have been good. Keep up the good work. Your clinician will continue to monitor your pressure readings. HDRBody2 notif.standard.looking_goo d_message	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2394	<Picture> READINGS REVIEWED CONTENT LABEL COMPONENT		UI Design	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Your recent CardioMEMS™ HF System readings have been good. Keep up the good work. Your clinician will continue to monitor your pressure readings.			CloudDRxRTR CloudDRxPS
PAGUI2397	8.12.1.44.3.5 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2398	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI208	8.12.1.45 RECORD WEIGHT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI209	8.12.1.45.1 RECORD WEIGHT PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2896	The RECORD WEIGHT PANEL is conditionally displayed as an element of the MORE SCREEN after the user selects the RECORD WEIGHT ACTION TRIGGER COMPONENT. The RECORD WEIGHT PANEL is displayed over the entire screen except the OS STATUS PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI210	The RECORD WEIGHT PANEL consists of: CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) RECORD WEIGHT HEADER COMPONENT RECORD WEIGHT INSTRUCTION LABEL COMPONENT RECORD WEIGHT INPUT COMPONENT SAVE ACTION TRIGGER COMPONENT (See Shared Elements: Components) CONFIRM INPUT PANEL (See Shared Elements: Panels) NUMERIC KEYBOARD COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI211	8.12.1.45.2 RECORD WEIGHT PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2897	<Picture> RECORD WEIGHT PANEL (example: Input)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI212	8.12.1.45.3 RECORD WEIGHT PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

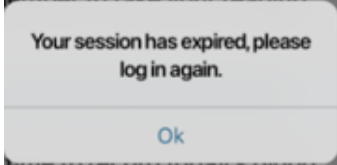
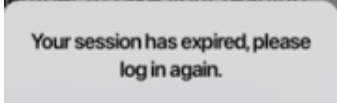
ID		English US Strings	Object Type	Project
PAGUI259	8.12.1.45.3.1 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI260	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI263	8.12.1.45.3.2 RECORD WEIGHT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI264	Component: RECORD WEIGHT HEADER COMPONENT State: N/A Display: Always displayed as a component of the RECORD WEIGHT PANEL. Label: <label.gb.weight> Behavior: Text.	PString521 Weight weight label.gb.weight	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2898	<Picture> RECORD WEIGHT HEADER COMPONENT Weight		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI265	8.12.1.45.3.3 RECORD WEIGHT INSTRUCTION COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI266	Component: RECORD WEIGHT INSTRUCTIOIN COMPONENT State: N/A Display: Always displayed as a component of the RECORD WEIGHT PANEL. Label: <label.gb.enter_weight> Behavior: Text.	PString435 Please weigh yourself and enter it below. weighYourselfAndEnterBelow label.gb.enter_weight	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2899	<Picture> RECORD WEIGHT INSTRUCTION COMPONENT Please weigh yourself and enter it below.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI267	8.12.1.45.3.4 RECORD WEIGHT INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI268	Component: RECORD WEIGHT INPUT COMPONENT State: Enabled/ Disabled Display: Always displayed as a component of the RECORD WEIGHT PANEL. The default state is disabled. Label: <date.Code.Weight_Unit_LB>, <data.Code.Med_Unit_Cd.Code_Desc.KG> Behavior: Input display. When enabled, display the input from the NUMERIC KEYBOARD PANEL and a curser within the input box. When disabled, display an empty input box. Always display input units.	PString342 kg kg data.Code.Med_Unit_Cd.Code_Desc.KG PString520 lb lb	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		data.Code.Weight_Unit_LB		
PAGUI2900	<Picture> RECORD WEIGHT INPUT COMPONENT (default: Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2901	<Picture> RECORD WEIGHT INPUT COMPONENT (example: Enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI271	8.12.1.45.3.5 SAVE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI272	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2923	8.12.1.45.3.6 CONFIRM INPUT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2924	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI269	8.12.1.45.3.7 NUMERIC KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI270	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3954	8.12.1.46 SESSION EXPIRED PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3955	8.12.1.46.1 SESSION EXPIRED PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3956	The SESSION EXPIRED PANEL is displayed when the user must login to the Patient Application again because their session has reached the time limit. The user is unable to use the Patient Application until they complete a new successful login. The SESSION EXPIRED PANEL is displayed over the screen, does not displace any elements, and any part of the screen that is not covered by the SESSION EXPIRED PANEL is greyed out.			CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3957	The SESSION EXPIRED PANEL is an element of:		UI Info	CloudDRx 1.0

myCardioMEMS Patient Application Graphical User Interface Specification

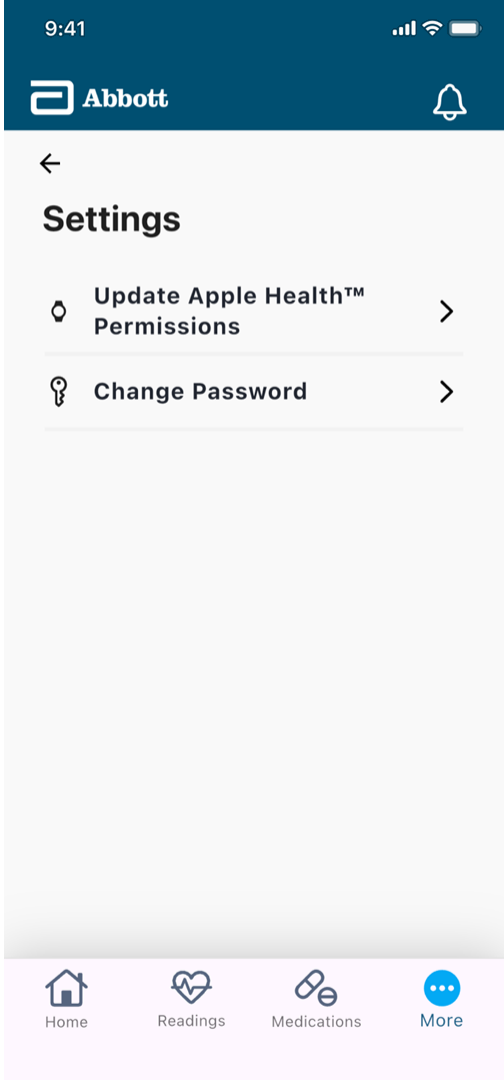
ID		English US Strings	Object Type	Project
	HOME SCREEN READINGS SCREEN MEDICATIONS SCREEN NOTIFICATIONS SCREEN MORE SCREEN			CloudDRxRTR CloudDRxPS
PAGUI3958	The SESSION EXPIRED PANEL is composed of: SESSION EXPIRED LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3959	8.12.1.46.2 SESSION EXPIRED PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3960	<Picture> SESSION EXPIRED PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3961	8.12.1.46.3 SESSION EXPIRED PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3962	8.12.1.46.3.1 SESSION EXPIRED LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3963	Component: SESSION EXPIRED LABEL COMPONENT State: N/A Display: Always displayed as a component of the SESSION EXPIRED PANEL. Label: <msg.gb.session_expired> Behavior: Text.	PString763 Your session has expired; please log in again. sessionExpiryAlertTitle msg.gb.session_expired	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3964	<Picture> SESSION EXPIRED LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3965	8.12.1.46.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3966	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI581	8.12.1.48 SETTINGS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

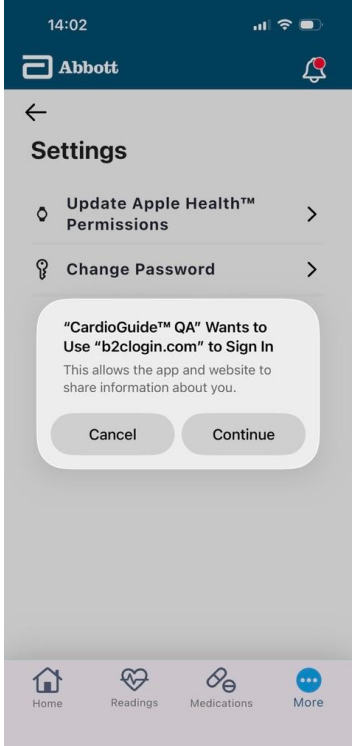
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1020	8.12.1.48.1 SETTINGS PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1024	The SETTINGS PANEL is displayed after the user navigates to it from the MORE SCREEN. The SETTINGS PANEL is displayed over the entire screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1025	The SETTINGS PANEL is composed of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) SETTINGS HEADER COMPONENT CONNECT HEALTH APP ACTION TRIGGER COMPONENT CONNECT HEALTH APP PANEL BP ERROR PANEL SUCCESSFUL SYNC PANEL PASSWORD ACTION TRIGGER COMPONENT SIGN IN POPUP DIALOG PANEL (See Shared Elements: Panels) USER DETAILS NEW PASSWORD PANEL SUCCESSFUL PASSWORD CHANGE PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1021	8.12.1.48.2 SETTINGS PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1023	<Picture> SETTINGS PANEL (example: Apple)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4628	<Picture> SETTINGS PANEL (example: PASSWORD ACTION TRIGGER COMPONENT triggered)		UI Design	CloudDRxRTR

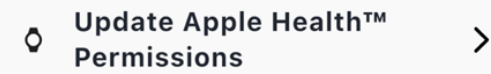


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1022	8.12.1.48.3 SETTINGS PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1660	8.12.1.48.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1661	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4606	8.12.1.48.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4607	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI1026	8.12.1.48.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1027	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1030	8.12.1.48.3.4 SETTINGS HEADER COMPONENT		Heading	CloudDRx 1.0

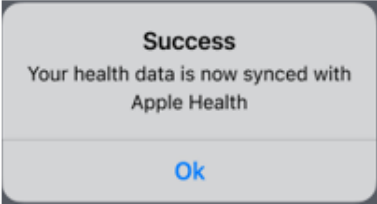
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1031	Component: SETTINGS HEADER COMPONENT State: N/A Display: Always displayed in the SETTINGS PANEL. Label: <label.gb.settings> Behavior: Text.	PString497 Settings settingsHome_title label.gb.settings	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1032	<Picture> SETTINGS HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1033	8.12.1.48.3.5 CONNECT HEALTH APP ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1034	Component: CONNECT HEALTH APP ACTION TRIGGER COMPONENT State: Apple/Android and Conected/ Disconnected Display: Always displayed in the SETTINGS PANEL. The state is determind by the native operating system. When in the Apple state, there is only one state. When in the Android state, the default state is disconnected. Label: <label.gb.update_apple_permissions>, <label.gb.enable_google_permissions>, <label.gb.revoke_google_permissions> Behavior: Icon, text, and trigger. Upon activation, in the Apple state and the Android Disconnected state, the CONNECT HEALTH APP ACTION TRIGGER COMPONENT displays the CONNECT HEALTH APP PANEL. Upon activation in the Android Connected state, the CONNECT HEALTH APP ACTION TRIGGER COMPONENT displays the PERMISSIONS REVOKED PANEL. The Android state changes from Disconnected to Connected when the user conects the applicable health app.	PString620 Update {19} Permissions updateAppleHealthPermissi ons label.gb.update_apple_per missions PString621 Enable {20} Permissions enableGoogleFitPermission s label.gb.enable_google_per missions PString622 Revoke {20} Permissions revokeGoogleFitPermission s label.gb.revoke_google_per missions	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1035	<Picture> CONNECT HEALTH APP ACTION TRIGGER COMPONENT (example: Apple) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1340	8.12.1.48.3.6 CONNECT HEALTH APP PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1663	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

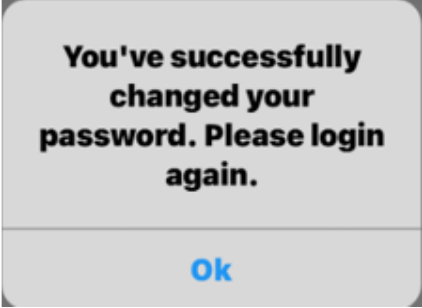
ID		English US Strings	Object Type	Project
PAGUI1664	8.12.1.48.3.7 BP ERROR PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1665	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1666	8.12.1.48.3.8 SUCCESSFUL SYNC PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1667	8.12.1.48.3.8.1 SUCCESSFUL SYNC PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1670	The SUCCESSFUL SYNC PANEL is conditionally displayed in the SETTINGS PANEL after the user successfully syncs with the Apple Health App for the first time. The SUCCESSFUL SYNC PANEL is displayed over the SETTINGS PANEL, does not displace any other elements, and any part of the screen that is not covered by the panel is greyed out.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1671	The SUCCESSFUL SYNC PANEL consists of: SUCCESSFUL SYNC LABEL COMPONENT OKAY ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1668	8.12.1.48.3.8.2 SUCCESSFUL SYNC PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1672	<Picture> SUCCESSFUL SYNC PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1669	8.12.1.48.3.8.3 SUCCESSFUL SYNC PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1673	8.12.1.48.3.8.3.1 SUCCESSFUL SYNC LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1674	Component: SUCCESSFUL SYNC LABEL COMPONENT State: N/A Display: Always displayed in the SUCCESSFUL SYNC PANEL. Label: N/A - determined by third party health application Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1675	<Picture> SUCCESSFUL SYNC LABEL COMPONENT		UI Design	CloudDRx 1.0



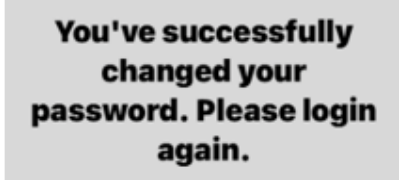
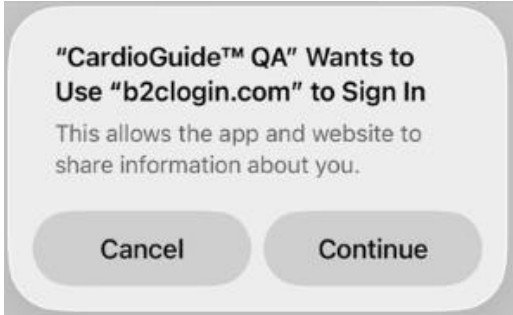
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Success Your health data is now synced with Apple Health			CloudDRxRTR CloudDRxPS
PAGUI1686	8.12.1.48.3.8.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1687	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1047	8.12.1.48.3.10 PASSWORD ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1048	Component: PASSWORD ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the SETTINGS PANEL. Label: <label.gb.password> Behavior: Icon, text, and trigger. Upon activation, the PASSWORD ACTION TRIGGER COMPONENT displays the SIGN IN POPUP DISPLAY PANEL.	PString218 Change Password changePasswordButton label.gb.change_password	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1049	<Picture> PASSWORD ACTION TRIGGER COMPONENT  Change Password >		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4632	8.12.1.48.3.11 SIGN IN POPUP DIALOG PANEL		Heading	CloudDRxRTR
PAGUI4633	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI3610	8.12.1.48.3.12 USER DETAILS NEW PASSWORD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3611	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3612	8.12.1.48.3.13 SUCCESSFUL PASSWORD CHANGE PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3613	8.12.1.48.3.13.1 SUCCESSFUL PASSWORD CHANGE PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3616	The USER DETAILS NEW PASSWORD PANEL is displayed as an element of the SETTINGS PANEL after the user has successfully changed their password in the USER DETAILS NEW PASSWORD PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3617	The SUCCESSFUL PASSWORD CHANGE PANEL consists of:		UI Info	CloudDRx 1.0

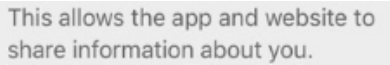
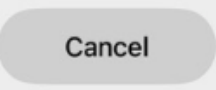
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	SUCCESSFUL CHANGE INSTRUCTIONS COMPONENT OKAY ACTION TRIGGER COMPONENT			CloudDRxRTR CloudDRxPS
PAGUI3614	8.12.1.48.3.13.2 SUCCESSFUL PASSWORD CHANGE PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3618	<Picture> SUCCESSFUL PASSWORD CHANGE PANEL 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3615	8.12.1.48.3.13.3 SUCCESSFUL PASSWORD CHANGE PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3619	8.12.1.48.3.13.3.1 SUCCESSFUL CHANGE INSTRUCTIONS COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3620	Component: SUCCESSFUL CHANGE INSTRUCTIONS COMPONENT State: N/A Display: Always displayed in the SUCCESSFUL PASSWORD CHANGE PANEL. Label: <msg.gb.password_changed> Behavior: Text.	PString698 You've successfully reset your password. You might want to make a note of your new password and store it in a safe place. successfulChangePassword AlertTitle msg.gb.password_changed PString1018 You've successfully changed your password. Please login again. successfulChangePassword AlertTitle msg.gb.password_changed	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3621	<Picture> SUCCESSFUL CHANGE INSTRUCTIONS COMPONENT		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

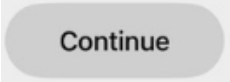
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI3622	8.12.1.48.3.13.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3623	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4608	8.12.1.49 SIGN IN POPUP DIALOG PANEL		Heading	CloudDRxRTR
PAGUI4609	Conditionally displayed when PASSWORD ACTION TRIGGER COMPONENT or LOGIN ACTION TRIGGER COMPONENT are activated. The SIGN IN POPUP DIALOG COMPONENT is displayed on top of the background screen which is greyed out. No elements of the background screen are displaced.		UI Info	CloudDRxRTR
PAGUI4611	The SIGN IN POPUP DIALOG PANEL consists of: SIGN IN POPUP HEADER COMPONENT SIGN IN POPUP LABEL COMPONENT CANCEL POPUP TRIGGER COMPONENT CONTINUE POPUP TRIGGER COMPONENT		UI Info	CloudDRxRTR
PAGUI4613	8.12.1.49.2 SIGN IN POPUP DIALOG PANEL (Layout)		Heading	CloudDRxRTR
PAGUI4610	<Picture> SIGN IN POPUP DIALOG PANEL 		UI Design	CloudDRxRTR
PAGUI4614	<Picture> SIGN IN POPUP DIALOG PANEL (example: displayed by triggering LOGIN ACTION TRIGGER COMPONENT on the STARTUP SCREEN)		UI Design	CloudDRxRTR

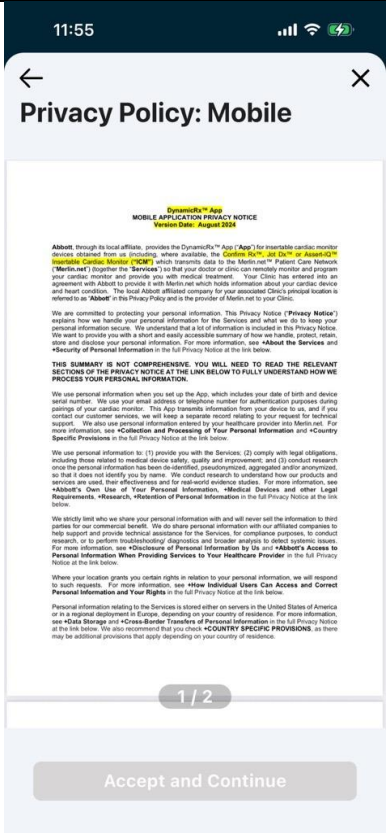
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4615	8.12.1.49.3 SIGN IN POPUP DIALOG PANEL (Elements)		Heading	CloudDRxRTR
PAGUI4616	8.12.1.49.3.1 SIGN IN POPUP HEADER COMPONENT		Heading	CloudDRxRTR
PAGUI4617	Component: SIGN IN POPUP HEADER COMPONENT State: N/A Display: Always displayed in the SIGN IN POPUP DIALOG PANEL Label: <ots.control.third_party_sign_in> Behavior: Text.	PString946 "{0}" Wants to Use "{36}" to Sign In ots.control.third_party_sign_in	UI Info	CloudDRxRTR
PAGUI4618	<Picture> SIGN IN POPUP HEADER COMPONENT (example)  "CardioGuide™ QA" Wants to Use "b2clogin.com" to Sign In		UI Design	CloudDRxRTR
PAGUI4619	8.12.1.49.3.2 SIGN IN POPUP LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI4620	Component: SIGN IN POPUP LABEL COMPONENT State: N/A Display: Always displayed in the SIGN IN POPUP DIALOG PANEL Label: <ots.control.third_party_info> Behavior: Text.	PString947 This allows the app and website to share information about you. ots.control.third_party_info	UI Info	CloudDRxRTR
PAGUI4621	<Picture> SIGN IN POPUP LABEL COMPONENT 		UI Design	CloudDRxRTR
PAGUI4622	8.12.1.49.3.3 CANCEL POPUP TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4623	Component: CANCEL POPUP TRIGGER COMPONENT State: N/A Display: Always displayed in the SIGN IN POPUP DIALOG PANEL Label: <ots.control.cancel> Behavior: Upon activation, the CANCEL POPUP TRIGGER COMPONENT, user returns to the previously displayed screen.	PString948 Cancel ots.control.cancel	UI Info	CloudDRxRTR
PAGUI4624	<Picture> CANCEL POPUP TRIGGER COMPONENT 		UI Design	CloudDRxRTR
PAGUI4625	8.12.1.49.3.4 CONTINUE POPUP TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4626	Component: CONTINUE POPUP TRIGGER COMPONENT State: N/A Display: Always displayed in the SIGN IN POPUP DIALOG PANEL Label: <ots.control.continue>	PString949 Continue ots.control.continue	UI Info	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Behavior: Upon activation in the STARTUP SCREEN, the CONTINUE POPUP TRIGGER COMPONENT proceeds to USER AUTHENTICATION PANEL. Upon activation in the SETTINGS PANEL, the CONTINUE POPUP TRIGGER COMPONENT proceeds to the USER DETAILS NEW PASSWORD PANEL.			
PAGUI4627	<Picture> CONTINUE POPUP TRIGGER COMPONENT 		UI Design	CloudDRxRTR
PAGUI4202	8.12.1.50 SHORT PRIVACY PANEL		Heading	CloudDRxRTR
PAGUI3901	8.12.1.50.1 SHORT PRIVACY POLICY PANEL (Display)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3902	The SHORT PRIVACY POLICY PANEL is displayed in the POLICY SCREEN after the user selects the CONSENT ACTION TRIGGER COMPONENT on the END USER LICENSE AGREEMENT PANEL, or after the user selects the LOGIN ACTION TRIGGER COMPONENT if there is an update to the Privacy Policy.		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3904	The SHORT PRIVACY POLICY PANEL consists of: OS STATUS PANEL BACK ACTION TRIGGER COMPONENT (Shared Elements: Components) CLOSE PANEL ACTION TRIGGER COMPONENT (Shared Elements: Components) SHORT PRIVACY POLICY HEADER COMPONENT SHORT PRIVACY POLICY DISPLAY COMPONENT ACCEPT ACTION TRIGGER COMPONENT		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3905	8.12.1.50.2 SHORT PRIVACY POLICY PANEL (Layout)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3906	<Picture> PRIVACY POLICY PANEL (example: top of the document)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI3907	<Picture> SHORT PRIVACY POLICY PANEL (example: bottom of document)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI3908	8.12.1.50.3 SHORT PRIVACY POLICY PANEL (Elements)		Heading	CloudDRxRTR CloudDRxPS
PAGUI4638	8.12.1.50.3.1 OS STATUS PANEL		Heading	CloudDRxRTR
PAGUI4639	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI3909	8.12.1.50.3.2 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3910	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI3911	8.12.1.50.3.3 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3912	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI3913	8.12.1.50.3.4 SHORT PRIVACY POLICY HEADER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3914	Component: SHORT PRIVACY POLICY HEADER COMPONENT	PString610	UI Info	CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	State: N/A Display: Always displayed in the SHORT PRIVACY POLICY PANEL. Label: <label.gb.PP_mobile> Behavior: Text.	Privacy Policy privacyPolicyTitle label.gb.PP PString798 Privacy Policy: Mobile privacyPolicyMobileTitle label.gb.PP_mobile		CloudDRxPS
PAGUI3915	<Picture> SHORT PRIVACY POLICY HEADER COMPONENT Privacy Policy: Mobile		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3916	8.12.1.50.3.5 SHORT PRIVACY POLICY DISPLAY COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3917	Component: SHORT PRIVACY POLICY DISPLAY COMPONENT State: N/A Display: Always displayed in the SHORT PRIVACY POLICY PANEL. Label: <other.pdf.Privacy_Policy> - the Patient Application displays the PDF document located with the attached website link Behavior: Text and trigger. The content for the SHORT PRIVACY POLICY DISPLAY COMPONENT is managed through content management in Merlin.net to enable updates as needed. The PDF document can be multiple pages with the numbr of pages displayed on the bottom and allows scrolling. If the user selects the hyperlink within the content the FULL PRIVACY POLICY PANEL is displayed.	PString703 Display a PDF of the applicable document. pdfPrivacyPolicy other.pdf.Privacy_Policy PString944 Display a PDF of the applicable document. pdfPrivacyPolicy other.pdf.Privacy_Policy	UI Info	CloudDRxRTR CloudDRxPS
PAGUI3918	<Picture> PRIVACY POLICY LABEL COMPONENT (example)		UI Design	CloudDRxRTR CloudDRxPS

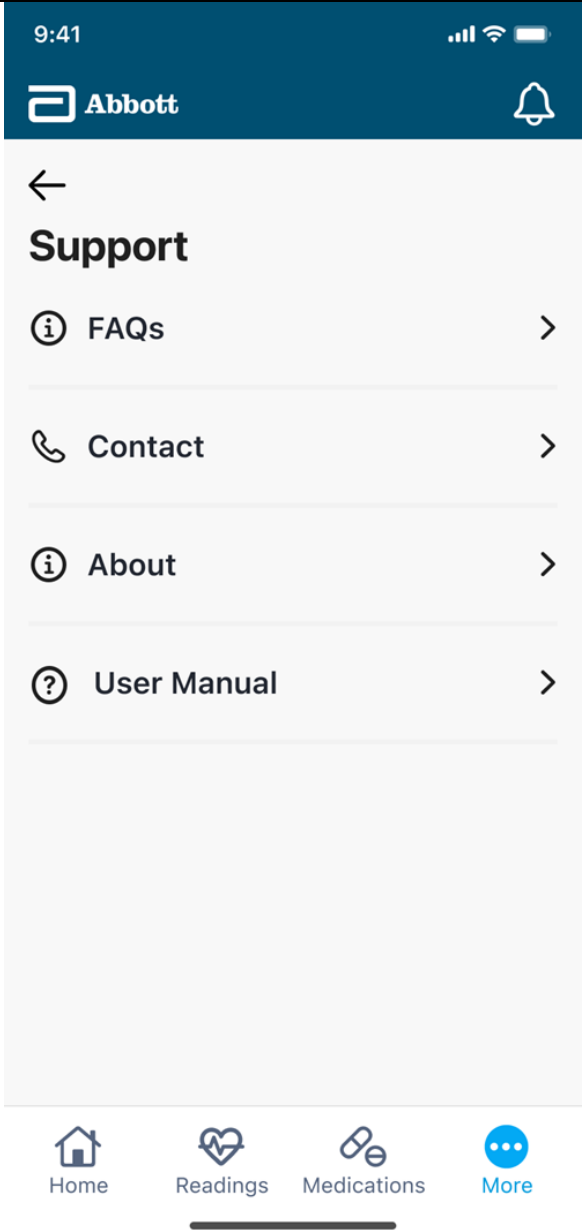
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Lorem ipsum dolor sit amet, consectetur adipiscing elit. Sed enim nunc, euismod vel porta tincidunt, condimentum eget risus. Nulla sodales, leo sit amet sagittis lobortis, sem risus rutrum nulla, sodales iaculis mi tellus pellentesque nulla. Curabitur sit amet felis et orci scelerisque semper. Pellentesque in dui viverra, cursus urna non, fermentum eros. Sed sed lectus vestibulum, semper augue non, auctor lorem. Aliquam ex nisl, finibus non aliquam sit amet, scelerisque venenatis dui. Quisque at sapien accumsan, finibus dolor vel, gravida magna. Nulla eleifend nibh ligula, quis sollicitudin magna dignissim at. Nunc imperdiet nunc vitae quam volutpat, at rhoncus eros dignissim. Ut convallis dui tortor, quis luctus erat egestas sed. Donec et scelerisque dolor, et posuere neque. Sed vitae justo ultrices, lobortis ipsum quis, interdum tortor.			
PAGUI4029	8.12.1.50.3.6 ACCEPT ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4024	Component: ACCEPT ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Always displayed in the SHORT PRIVACY POLICY AGREEMENT PANEL. Default state is disabled. Label: <label.gb.accept_continue> Behavior: Text and trigger. The ACCEPT ACTION TRIGGER COMPONENT is enabled when the user scrolls to the bottom of the SHORT PRIVACY POLICY DISPLAY COMPONENT. Upon activation of the ACCEPT ACTION TRIGGER COMPONENT in enabled state the HOME SCREEN is displayed if the user is reviewing the most up to date Privacy Policy. If the user is registering for the first time the PATIENT IDENTIFICATION PANEL is displayed.	PString678 Accept and Continue acceptAndContinue label.gb.accept_continue	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4025	<Picture> ACCEPT ACTION TRIGGER COMPONENT (default: disabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4026	<Picture> ACCEPT ACTION TRIGGER COMPONENT (enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI582	8.12.1.51 SUPPORT PANEL		Heading	CloudDRx 1.0

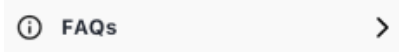


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI784	8.12.1.51.1 SUPPORT PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI787	The SUPPORT PANEL is displayed after the user navigates to it from the MORE SCREEN. The SUPPORT PANEL is displayed over the entire screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI788	The SUPPORT PANEL is composed of: OS STATUS PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) SUPPORT HEADER COMPONENT FAQ ACTION TRIGGER COMPONENT FAQ PANEL CONTACT ACTION TRIGGER COMPONENT CONTACT PANEL ABOUT ACTION TRIGGER COMPONENT ABOUT PANEL USER MANUAL ACTION TRIGGER COMPONENT USER MANUAL WEB VIEW		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI785	8.12.1.51.2 SUPPORT PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI789	<Picture> SUPPORT PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
				
PAGUI786	8.12.1.51.3 SUPPORT PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1804	8.12.1.51.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1805	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI790	8.12.1.51.3.2 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI791	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI794	8.12.1.51.3.3 SUPPORT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI795	Component: SUPPORT HEADER COMPONENT State: N/A Display: Always displayed in the SUPPORT PANEL. Label: <label.gb.support> Behavior: Text.	PString585 Support support label.gb.support	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI796	<Picture> SUPPORT HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI797	8.12.1.51.3.4 FAQ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI798	Component: FAQ ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the SUPPORT PANEL. Label: <label.gb.faqs> Behavior: Icon, text, and trigger. Upon activation, the FAQ ACTION TRIGGER COMPONENT displays the FAQ PANEL.	PString200 FAQs faqs label.gb.faqs	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI799	<Picture> FAQ ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI800	8.12.1.51.3.5 FAQ PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI801	8.12.1.51.3.5.1 FAQ PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
PAGUI853	The FAQ PANEL is displayed in the SUPPORT PANEL after the user selects the FAQ ACTION TRIGGER COMPONENT. The FAQ PANEL is displayed over the whole screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI854	The FAQ PANEL is composed of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) FAQ HEADER COMPONENT FAQ OVERVIEW LABEL COMPONENT FAQ EMERGENCY LABEL COMPONENT FAQ TOPIC JUMP ACTION TRIGGER COMPONENT FAQ TOPIC PANEL(See Shared Elements: Panels) BACK TO TOP ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI802	8.12.1.51.3.5.2 FAQ PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI804	<Picture> FAQ PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

[illegible]



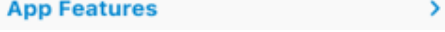
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1370	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4630	8.12.1.51.3.5.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4631	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI918	8.12.1.51.3.5.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI919	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI920	8.12.1.51.3.5.3.4 FAQ HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI921	Component: FAQ HEADER COMPONENT State: N/A Display: Always displayed in the FAQ PANEL. Label: <label.gb.faqs> Behavior: Text.	PString200 FAQs faqs label.gb.faqs	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI922	<Picture> FAQ HEADER COMPONENT FAQs		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI923	8.12.1.51.3.5.3.5 FAQ OVERVIEW LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI924	Component: FAQ OVERVIEW LABEL COMPONENT State: N/A Display: Always displayed in the FAQ PANEL. Label: N/A - determined by content manager Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI925	<Picture> FAQ OVERVIEW LABEL COMPONENT	PString802 The {33} for {21} is for patients with the {22} CardioMEMS™ PA pressure sensor and their caregivers. For simplicity “{21} for {22}” has been shortened to “{0}” for use throughout the	UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	The CardioGuide™ HF Mobile Application for CardioMEMS™ HF System is for patients with the Abbott CardioMEMS™ PA pressure sensor and their caregivers. For simplicity, "CardioGuide™ HF Mobile Application for CardioMEMS™ HF System" has been shortened to "CardioGuide™ HF App" for use throughout the following Frequently Asked Questions (FAQs).	following Frequently Asked Questions (FAQs). faqScreenSubHeadline label.faq_overview PString804 CardioGuide™ HF Mobile Application param.33.app_brand_long		
PAGUI1891	8.12.1.51.3.5.3.6 FAQ EMERGENCY LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1892	Component: FAQ EMERGENCY LABEL COMPONENT State: N/A Display: Always displayed in the FAQ PANEL. Label: <label.gb.emergency_instructions> Behavior: Text.	PString624 In a clinical emergency, do not use the app; call your local emergency number. faqWarningText label.gb.emergency_instructions PString1016 In an emergency, call your local emergency number. faqWarningText label.gb.emergency_instructions	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1893	<Picture> FAQ EMERGENCY LABEL COMPONENT In an emergency, call your local emergency number.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI926	8.12.1.51.3.5.3.7 FAQ TOPIC JUMP ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI927	Component: FAQ TOPIC JUMP ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the FAQ PANEL. A FAQ TOPIC JUMP ACTION TRIGGER COMPONENT is displayed for each FAQ TOPIC PANEL. The FAQ TOPIC JUMP ACTION TRIGGER COMPONENT displays the title of the general topic for the associated FAQ TOPIC PANEL. Label: N/A - determined by content manager Behavior: Text and trigger. Upon activation the FAQ TOPIC JUMP ACTION TRIGGER COMPONENT scrolls the FAQ PANEL down to the associated FAQ TOPIC PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

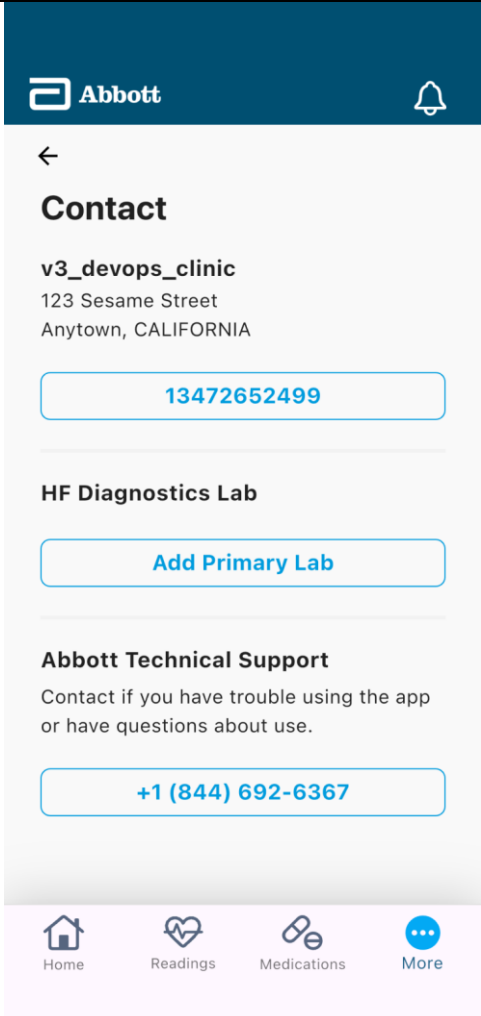
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI928	<Picture> SHOW FAQ INFORMATION ACTION TRIGGER EXPANDED COMPONENT (example: About the myCardioMEMS App) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1373	<Picture> SHOW FAQ INFORMATION ACTION TRIGGER EXPANDED COMPONENT (example: Getting Started) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1374	<Picture> SHOW FAQ INFORMATION ACTION TRIGGER EXPANDED COMPONENT (example: App Features) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1375	<Picture> SHOW FAQ INFORMATION ACTION TRIGGER EXPANDED COMPONENT (example: Troubleshooting) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1894	8.12.1.51.3.5.3.8 FAQ TOPIC PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1895	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1896	8.12.1.51.3.5.3.9 BACK TO TOP ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1897	Component: BACK TO TOP ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the FAQ PANEL. Label: N/A Behavior: Trigger. Upon activation, the BACK TO TOP ACTION TRIGGER COMPONENT brings the user to the top of the scrollable panel.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1899	<Picture> BACK TO TOP ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI805	8.12.1.51.3.6 CONTACT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI806	Component: CONTACT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the SUPPORT PANEL. Label: <label.gb.contact> Behavior: Icon, text, and trigger. Upon activation, the CONTACT ACTION TRIGGER COMPONENT displays the CONTACT PANEL.	PString201 Contact contact label.gb.contact	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

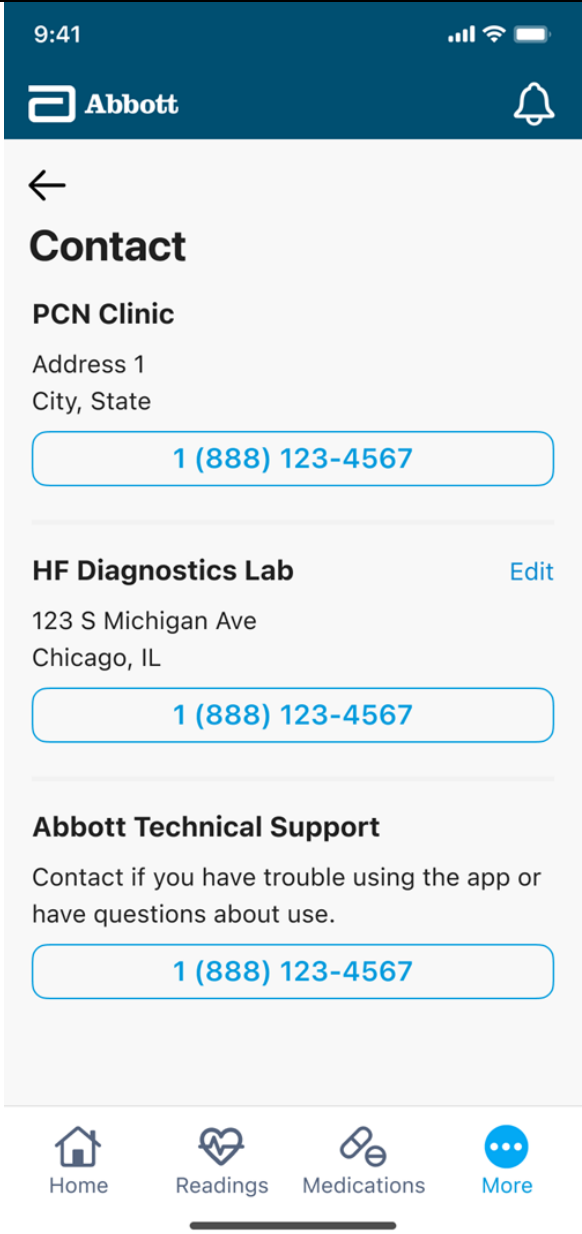


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI807	<Picture> CONTACT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI808	8.12.1.51.3.7 CONTACT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI809	8.12.1.51.3.7.1 CONTACT PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI813	The CONTACT PANEL is displayed in the SUPPORT PANEL after the user selects the CONTACT ACTION TRIGGER COMPONENT. The CONTACT PANEL is displayed over the whole screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI814	The CONTACT PANEL consists of: OS STATUS PANEL (Shared Elements: Panels) GLOBAL HEADING PANEL (Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) CONTACT HEADER COMPONENT PCN CLINIC HEADER COMPONENT PCN CLINIC ADDRESS LABEL COMPONENT CALL PCN CLINIC ACTION TRIGGER COMPONENT HF DIAGNOSTICS LAB HEADER COMPONENT EDIT PRIMARY LAB ACTION TRIGGER COMPONENT HF DIAGNOSTICS LAB ADDRESS LABEL COMPONENT ADD LAB ACTION TRIGGER COMPONENT PRIMARY LAB PANEL CALL HF DIAGNOSTICS LAB ACTION TRIGGER COMPONENT TECHNICAL SUPPORT HEADER COMPONENT TECHNICAL SUPPORT LABEL COMPONENT CALL ABBOTT TECHNICAL SUPPORT ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI810	8.12.1.51.3.7.2 CONTACT PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI811	<Picture> CONTACT PANEL (example: no HF Diagnostics Lab)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
				
PAGUI1906	<Picture> CONTACT PANEL (example: generic)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI812	8.12.1.51.3.7.3 CONTACT PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

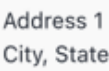


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1909	8.12.1.51.3.7.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1910	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4647	8.12.1.51.3.7.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4648	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI815	8.12.1.51.3.7.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI816	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI817	8.12.1.51.3.7.3.4 CONTACT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI818	Component: CONTACT HEADER COMPONENT State: N/A Display: Always displayed in the CONTACT PANEL. Label: <label.gb.contact> Behavior: Text.	PString201 Contact contact label.gb.contact	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI819	<Picture> CONTACT HEADER COMPONENT Contact		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI820	8.12.1.51.3.7.3.5 PCN CLINIC HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI821	Component: PCN CLINIC HEADER COMPONENT State: N/A Display: Always displayed in the CONTACT PANEL. The PCN CLINIC HEADER COMPONENT displays the clinic name associated with the patient in the backend. Label: N/A - determined by clinic information in the portal Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI822	<Picture> PCN CLINIC HEADER COMPONENT (example: generic) PCN Clinic		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI823	8.12.1.51.3.7.3.6 PCN CLINIC ADDRESS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI824	Component: PCN CLINIC ADDRESS LABEL COMPONENT State: N/A Display: Always displayed in the CONTACT PANEL. The PCN CLINIC ADDRESS LABEL COMPONENT displays the clinic address associated with the patient in the backend. Label: N/A - determined by clinic information in the portal Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI825	<Picture> PCN CLINIC ADDRESS LABEL COMPONENT (example: generic) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI826	8.12.1.51.3.7.3.7 CALL PCN CLINIC ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI827	Component: CALL PCN CLINIC ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the CONTACT PANEL. The CALL PCN CLINIC ACTION TRIGGER COMPONENT displays the telephone number for the clinic associated with the patient in the backend. Label: N/A - determined by clinic information in the portal Behavior: Text and trigger. Upon activation, the CALL PCN CLINIC ACTION TRIGGER COMPONENT opens the phone on the user's device and dials the displayed telephone number.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI828	<Picture> CALL PCN CLINIC ACTION TRIGGER COMPONENT (example: generic) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI829	8.12.1.51.3.7.3.8 HF DIAGNOSTICS LAB HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI830	Component: HF DIAGNOSTICS LAB HEADER COMPONENT State: No Input/ Input Display: Always displayed in the CONTACT PANEL. Default state is no input. Label: <label.gb.hf_diagnostics_lab> Behavior: Text. In the no input state, the HF DIAGNOSTICS LAB HEADER COMPONENT displays a generic label. In the input state, the HF DIAGNOSTICS LABE HEADER COMPONENT displays the most recent input from the user in the PRIMARY LAB PANEL.	PString625 HF Diagnostics Lab primaryLab label.gb.hf_diagnostics_lab	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI831	<Picture> HF DIAGNOSTICS LAB HEADER COMPONENT (example: generic, no input state) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1377	8.12.1.51.3.7.3.9 EDIT PRIMARY LAB ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1378	Component: EDIT PRIMARY LAB ACTION TRIGGER COMPONENT	PString613	UI Info	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	State: N/A Display: Conditionally displayed in the CONTACT PANEL.The EDIT PRIMARY LAB ACTION TRIGGER COMPONENT is displayed after the user has successfully saved information in the PRIMARY LAB PANEL. Label: <label.gb.edit> Behavior: Text and trigger. Upon activation, the EDIT PRIMARY LAB ACTION TRIGGER COMPONENT displays the PRIMARY LAB PANEL.	Edit edit label.gb.edit		CloudDRxRTR CloudDRxPS
PAGUI1379	<Picture> PRIMARY LAB ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI832	8.12.1.51.3.7.3.10 HF DIAGNOSTICS LAB ADDRESS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI833	Component: HF DIAGNOSTICS LAB ADDRESS LABEL COMPONENT State: N/A Display: Conditionally displayed in the CONTACT PANEL. The HF DIAGNOSTICS LAB ADDRESS LABEL COMPONENT is displayed after the user has successfully saved information in the PRIMARY LAB PANEL. Label: N/A - determined by user input in the PRIMARY LAB PANEL Behavior: Text. The HF DIAGNOSTIC LAB ADDRESS LABEL COMPONENT displays inputted information from the PRIMARY LAB PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI834	<Picture> HF DIAGNOSTICS LAB ADDRESS LABEL COMPONENT (example: generic) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1917	8.12.1.51.3.7.3.11 ADD LAB ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1918	Component: ADD LAB ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the CONTACT PANEL. The ADD LAB ACTION TRIGGER COMPONENT is displayed until the user has successfully saved information in the PRIMARY LAB PANEL. Label: <label.gb.add_lab> Behavior: Text and trigger. Upon activation, the ADD LAB ACTION TRIGGER COMPONENT displays the PRIMARY LAB PANEL.	PString626 Add Primary Lab addPrimaryLab label.gb.add_lab	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1919	<Picture> ADD LAB ACTION TRIGGER COMPONET 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1380	8.12.1.51.3.7.3.12 PRIMARY LAB PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1381	8.12.1.51.3.7.3.12.1 PRIMARY LAB PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1384	The PRIMARY LAB PANEL is displayed in the CONTACT PANEL after the user selects the PRIMARY LAB ACTION TRIGGER COMPONENT. The PRIMARY LAB PANEL is displayed over the whole screen except for the OS STATUS PANEL and the GLOBAL HEADING PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1385	The PRIMARY LAB PANEL is composed of: CLOSE PANEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) PRIMARY LAB HEADER COMPONENT LAB NAME LABEL COMPONENT LAB NAME INPUT COMPONENT LAB NAME ERROR COMPONENT STREET ADDRESS LABEL COMPONENT STREET ADDRESS INPUT COMPONENT CITY LABEL COMPONENT CITY INPUT COMPONENT STATE LABEL COMPONENT STATE INPUT COMPONENT PHONE LABEL COMPONENT PHONE INPUT COMPONENT SAVE ACTION TRIGGER COMPONENT ALPHABETICAL KEYBOARD COMPONENT NUMERIC KEYBOARD COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1382	8.12.1.51.3.7.3.12.2 PRIMARY LAB PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1386	<Picture> PRIMARY LAB PANEL (example: generic)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

	9:41 Abbott Edit Primary Lab Name HF Diagnostics Lab Street Address 123 S Michigan Ave City Chicago State/Province IL Phone Number  +1 Enter phone number Save  Home  Readings  Medications  More				
--	--	--	--	--	--

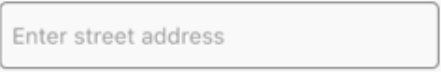
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1383	8.12.1.51.3.7.3.12.3 PRIMARY LAB PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1387	8.12.1.51.3.7.3.12.3.1 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1388	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1389	8.12.1.51.3.7.3.12.3.2 PRIMARY LAB HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1390	Component: PRIMARY LABL HEADER COMPONENT State: Add/ Edit Display: Always displayed in the PRIMARY LAB PANEL. The default state is add. Label: <label.gb.add_lab>, <label.gb.edit_lab> Behavior: Text. The PRIMARY LABEL HEADER COMPONENT changes state from add to edit after the user succefully saves input in the PRIMARY LAB PANEL AT LEAST ONCE.	PString626 Add Primary Lab addPrimaryLab label.gb.add_lab PString654 Edit Primary Lab editPrimaryLab label.gb.edit_lab	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1911	<Picture> PRIMARY LAB HEADER COMPONENT (default: add) Add Primary Lab		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1391	<Picture> PRIMARY LAB HEADER COMPONENT (example: edit) Edit Primary Lab		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1392	8.12.1.51.3.7.3.12.3.3 LAB NAME LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1393	Component: LAB NAME LABEL COMPONENT State: N/A Display: Always displayed in the PRIMARY LAB PANEL. Label: <label.gb.name> Behavior: Text.	PString210 Name name label.gb.name	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1394	<Picture> LAB NAME LABEL COMPONENT Name		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1395	8.12.1.51.3.7.3.12.3.4 LAB NAME INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1396	Component: LAB NAME INPUT COMPONENT State: No Input/ Input Display: Always displayed in the PRIMARY LAB PANEL. The default state is no input. Label: <label.gb.enter_lab_name> Behavior: Text and input box. In the no input state, display the default text. In the input state, display the input given from the ALPHABETICAL KEYBOARD COMPONENT or the currently saved lab name.	PString628 Enter lab name enterName label.gb.enter_lab_name	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1397	<Picture> LAB NAME INPUT COMPONENT (default: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1912	<Picture> LAB NAME INPUT COMPONENT (input, example: generic) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3697	8.12.1.51.3.7.3.12.3.5 LAB NAME ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3698	Component: LAB NAME ERROR COMPONENT State: N/A Display: Conditionally displayed in the PRIMARY LAB PANEL when the user selects away from the LAB NAME INPUT COMPONENT and there is no input. Label: <label.gb.please_enter_lab_name> Behavior: Text.	PString697 Please enter lab name pleaseEnterName label.gb.please_enter_lab_name	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3699	<Picture> LAB NAME ERROR COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1398	8.12.1.51.3.7.3.12.3.6 STREET ADDRESS LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1399	Component: STREET ADDRESS LABEL COMPONENT State: N/A Display: Always displayed in the PRIMARY LAB PANEL. Label: <label.gb.street_address> Behavior: Text.	PString629 Street Address streetAddress label.gb.street_address	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1400	<Picture> STREET ADDRESS LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1401	8.12.1.51.3.7.3.12.3.7 STREET ADDRESS INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1402	Component: STREET ADDRESS INPUT COMPONENT State: No Input/ Input Display: Always displayed in the PRIMARY LAB PANEL. The default state is no input. Label: <label.gb.enter_street_address> Behavior: Text and input box. In the no input state, display the default text. In the input state, display the input given from the ALPHABETICAL KEYBOARD COMPONENT or the currently saved street address.	PString630 Enter street address enterStreetAddress label.gb.enter_street_addre ss	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1403	<Picture> STREET ADDRESS INPUT COMPONENT (default: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1913	<Picture> STREET ADDRESS INPUT COMPONENT (input, example: generic) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1404	8.12.1.51.3.7.3.12.3.8 CITY LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1405	Component: CITY LABEL COMPONENT State: N/A Display: Always displayed in the PRIMARY LAB PANEL. Label: <label.gb.city> Behavior: Text.	PString213 City city label.gb.city	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1406	<Picture> CITY LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1407	8.12.1.51.3.7.3.12.3.9 CITY INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1408	Component: CITY INPUT COMPONENT State: No Input/ Input Display: Always displayed in the PRIMARY LAB PANEL. The default state is no input. Label: <label.gb.enter_city> Behavior: Text and input box. In the no input state, display the default text. In the input state, display the input given from the ALPHABETICAL KEYBOARD COMPONENT or the currently saved city.	PString631 Enter city enterCity label.gb.enter_city	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1409	<Picture> CITY INPUT COMPONENT (default: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1914	<Picture> CITY INPUT COMPONENT (input, example: generic)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1410	8.12.1.51.3.7.3.12.3.10 STATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1411	Component: STATE LABEL COMPONENT State: N/A Display: Always displayed in the PRIMARY LAB PANEL. Label: <label.gb.state_province> Behavior: Text.	PString214 State/Province stateProvince label.gb.state_province	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1412	<Picture> STATE LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1413	8.12.1.51.3.7.3.12.3.11 STATE INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1414	Component: STATE INPUT COMPONENT State: No Input/ Input Display: Always displayed in the PRIMARY LAB PANEL. The default state is no input. Label: <label.gb.enter_province> Behavior: Text and input box. In the no input state, display the default text. In the input state, display the input given from the ALPHABETICAL KEYBOARD COMPONENT or the currently saved state or province.	PString632 Enter province enterStateOrProvince label.gb.enter_province	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1415	<Picture> STATE INPUT COMPONENT (default: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1915	<Picture> STATE INPUT COMPONENT (input, example: generic) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1416	8.12.1.51.3.7.3.12.3.12 PHONE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4848	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI1419	8.12.1.51.3.7.3.12.3.13 PHONE INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4849	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI1422	8.12.1.51.3.7.3.12.3.14 SAVE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

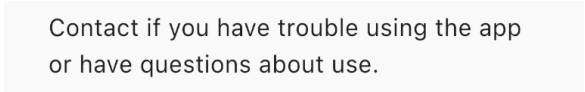


myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1423	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1424	8.12.1.51.3.7.3.12.3.15 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1425	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4645	8.12.1.51.3.7.3.12.3.16 NUMERIC KEYBOARD COMPONENT		Heading	CloudDRxRTR
PAGUI4646	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI835	8.12.1.51.3.7.3.13 CALL HF DIAGNOSTICS LAB ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI836	Component: CALL HF DIAGNOSTICS LAB ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the CONTACT PANEL. The CALL HF DIAGNOSTICS LAB ACTION TRIGGER COMPONENT displays the telephone number for the HF diagnostics laboratory associated with the patient in the backend. Label: N/A - determined by user input in the PRIMARY LAB PANEL Behavior: Text and trigger. Upon activation, the CALL HF DIAGNOSTICS LAB ACTION TRIGGER COMPONENT opens the phone on the user's device and dials the displayed telephone number.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI837	<Picture> CALL HF DIAGNOSTICS LAB ACTION TRIGGER COMPONENT (example: generic) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI838	8.12.1.51.3.7.3.14 TECHNICAL SUPPORT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI839	Component: TECHNICAL SUPPORT HEADER COMPONENT State: N/A Display: Always displayed in the CONTACT PANEL. Label: <label.tech_support>, <label.gb.technical_support> Behavior: Text.	PString247 Tech Support techSupport label.gb.tech_support PString634 {22} Technical Support abbottTechnicalSupport label.gb.technical_support	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI840	<Picture> HF DIAGNOSTICS LAB HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI841	8.12.1.51.3.7.3.15 TECHNICAL SUPPORT LABEL COMPONENT		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI842	Component: TECHNICAL SUPPORT LABEL COMPONENT State: N/A Display: Always displayed in the CONTACT PANEL. Label: <label.gb.technical_support_instructions> Behavior: Text.	PString635 Contact if you have trouble using the app or have questions about use. abbottTechnicalSupportSub Heading label.gb.technical_support_instructions	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI843	<Picture> TECHNICAL SUPPORT LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI847	8.12.1.51.3.7.3.16 CALL ABBOTT TECHNICAL SUPPORT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI848	Component: CALL ABBOTT TECHNICAL SUPPORT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the CONTACT PANEL. The CALL HF DIAGNOSTICS LAB ACTION TRIGGER COMPONENT displays the telephone number for the Abbott technical support team associated with the patient in the backend. Label: N/A - determined by content manager Behavior: Text and trigger. Upon activation, the CALL ABOUT TECHNICAL SUPPORT ACTION TRIGGER COMPONENT opens the phone on the user's device and dials the displayed telephone number.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI849	<Picture> CALL ABBOTT TECHNICAL SUPPORT ACTION TRIGGER COMPONENT (example: generic) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1806	8.12.1.51.3.8 ABOUT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1807	Component: ABOUT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the SUPPORT PANEL. Label: <label.gb.about> Behavior: Icon, text, and trigger. Upon activation, the ABOUT ACTION TRIGGER COMPONENT displays the ABOUT PANEL.	PString296 About about label.gb.about	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1808	<Picture> ABOUT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1809	8.12.1.51.3.9 ABOUT PANEL		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1810	8.12.1.51.3.9.1 ABOUT PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1813	The ABOUT PANEL is displayed in the SUPPORT PANEL after the user selects the ABOUT ACTION TRIGGER COMPONENT. The ABOUT PANEL is displayed over the whole screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1814	The ABOUT PANEL is composed of: OS STATUS PANEL (See Shared Elements: Panels) GLOBAL HEADING PANEL (See Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (See Shared Elements: Components) ABOUT HEADER COMPONENT APP MODEL LABEL COMPONENT APP VERSION LABEL COMPONENT GTIN LABEL COMPONENT PRIVACY POLICY ACTION TRIGGER COMPONENT FULL PRIVACY POLICY PANEL END USER LICENSE AGREEMENT ACTION TRIGGER COMPONENT END USER LICENSE AGREEMENT PANEL THIRD PARTY LICENSE ACTION TRIGGER COMPONENT THIRD PARTY LICENSE PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1811	8.12.1.51.3.9.2 ABOUT PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1815	<Picture> ABOUT PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
	9:41 Abbott About App Model Number: APP1003 App Software Version: 1.2.3 GTIN: 5415067027610 Privacy Policy End User License Agreement Third Party Licences Home Readings Medications More			
PAGUI1812	8.12.1.51.3.9.3 ABOUT PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

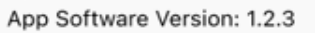


myCardioMEMS Patient Application Graphical User Interface Specification

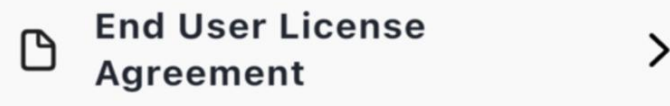
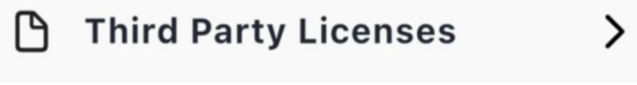
ID		English US Strings	Object Type	Project
PAGUI1816	8.12.1.51.3.9.3.1 OS STATUS PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1817	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4649	8.12.1.51.3.9.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR
PAGUI4699	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI1818	8.12.1.51.3.9.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1819	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1820	8.12.1.51.3.9.3.4 ABOUT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1821	Component: ABOUT HEADER COMPONENT State: N/A Display: Always displayed in the ABOUT PANEL Label: <label.gb.about> Behavior: Text.	PString296 About about label.gb.about	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1822	<Picture> ABOUT HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1823	8.12.1.51.3.9.3.5 APP MODEL LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1824	Component: APP MODEL LABEL COMPONENT State: N/A Display: Always displayed in the ABOUT PANEL. Label: <label.gb.model_number> Behavior: Text and value. The APP MODEL LABEL COMPONENT displays a label and the Patient Application model number for the device.	PString297 App Model Number: {2} modelNumber label.gb.model_number	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1825	<Picture> APP MODEL LABEL COMPONENT (example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1826	8.12.1.51.3.9.3.6 APP VERSION LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

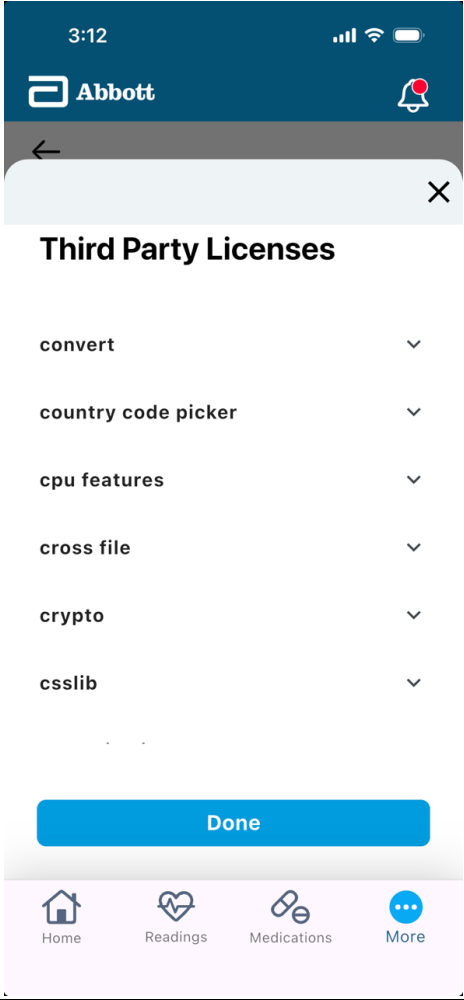
ID		English US Strings	Object Type	Project
PAGUI1827	Component: APP VERSION LABEL COMPONENT State: N/A Display: Always displayed in the ABOUT PANEL. Label: <label.gb.software_version> Behavior: Text and value: The APP VERSION LABEL COMPONENT displays a label and the Patient Application software version for the device.	PString298 App Software Version: {3} softwareVersion label.gb.software_version	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1828	<Picture> APP SOFTWARE LABEL COMPONENT (example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1829	8.12.1.51.3.9.3.7 GTIN LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1830	Component: GTIN LABEL COMPONENT State: N/A Display: Always displayed in the ABOUT PANEL. Label: <label.gb.gtin> Behavior: Text and value. The APP MODEL LABEL COMPONENT displays a label and the Global Trade Item Number for the Patient Application.	PString338 GTIN: {4} gtin label.gb.gtin	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1831	<Picture> GTIN LABEL COMPONENT (example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1870	8.12.1.51.3.9.3.8 PRIVACY POLICY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4054	Component: PRIVACY POLICY ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the ABOUT PANEL. Label: <label.gb.PP> Behavior: Icon, text, and trigger. Upon activation, the PRIVACY POLICY ACTION TRIGGER COMPONENT displays the FULL PRIVACY POLICY PANEL.	PString610 Privacy Policy privacyPolicyTitle label.gb.PP	UI Info	CloudDRxRTR CloudDRxPS
PAGUI1872	<Picture> PRIVACY POLICY ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1873	8.12.1.51.3.9.3.9 FULL PRIVACY POLICY PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1874	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1883	8.12.1.51.3.9.3.10 END USER LICENSE AGREEMENT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1884	Component: END USER LICENSE AGREEMENT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the ABOUT PANEL. Label: <label.gb.EULA> Behavior: Icon, text, and trigger. Upon activation, the END USER LICENSE AGREEMENT ACTION TRIGGER COMPONENT displays the END USER LICENSE AGREEMENT PANEL.	PString581 End User License Agreement termsOfUseTitle label.gb.EULA	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1885	<Picture> END USER LICENSE AGREEMENT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1886	8.12.1.51.3.9.3.11 END USER LICENSE AGREEMENT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1887	See Shared Elements: Panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3769	8.12.1.51.3.9.3.18 THIRD PARTY LICENSE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3770	Component: THIRD PARTY LICENSE ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the ABOUT PANEL. Label: <label.gb.third_party_license_title> Behavior: Icon, text, and trigger. Upon activation, the THIRD PARTY LICENSE ACTION TRIGGER COMPONENT displays the THIRD PARTY LICENSE PANEL.	PString715 Third Party Licenses thirdPartyLicenses label.gb.third_party_license_title	UI Info	CloudDRxRTR CloudDRxPS
PAGUI3771	<Picture> THIRD PARTY LICENSE ACTION TRIGGER COMPONENT 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3772	8.12.1.51.3.9.3.19 THIRD PARTY LICENSE PANEL		Heading	CloudDRxRTR CloudDRxPS
PAGUI3773	8.12.1.51.3.9.3.19.1 THIRD PARTY LICENSE PANEL (Display)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3776	The THIRD PARTY LICENSE PANEL is displayed after the user selects the THIRD PARTY LICENSE ACTION TRIGGER COMPONENT. The THIRD PARTY LICENSE PANEL displays a list of any third party software products used within the Patient Application. The THIRD PARTY LICENSE PANEL is displayed over the entire screen except the OS STATUS PANEL and GLOBAL HEADING PANEL. The THIRD PARTY LICENSE PANEL may be scrollable to allow the user to view all the contents of the panel. By default, it is displayed with THIRD PARTY DROPDOWN ACTION TRIGGER COMPONENTS closed.		UI Info	CloudDRxRTR CloudDRxPS

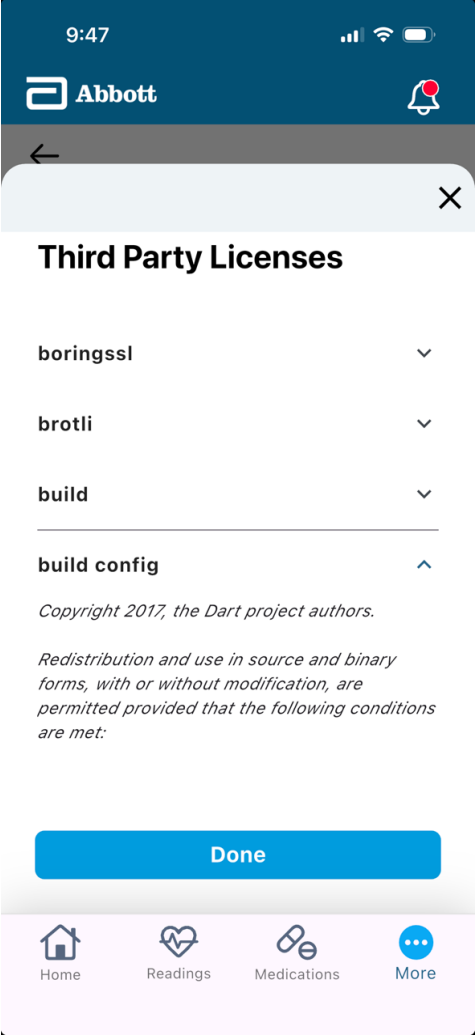
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3777	The THIRD PARTY LICENSE PANEL consists of: CLOSE PANEL ACTION TRIGGER COMPONENT (Shared Elements: Components) THIRD PARTY HEADER COMPONENT THIRD PARTY NAME COMPONENT THIRD PARTY DROPDOWN ACTION TRIGGER COMPONENT THIRD PARTY LIST ITEM DESCRIPTION COMPONENT DONE ACTION TRIGGER COMPONENT (Shared Elements: Components)		UI Info	CloudDRxRTR CloudDRxPS
PAGUI3775	8.12.1.51.3.9.3.19.2 THIRD PARTY LICENSE PANEL (Layout)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3778	<Picture> THIRD PARTY LICENSE PANEL (default: dropdowns closed) 		UI Design	CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4678	<Picture> THIRD PARTY LICENSE PANEL (example: one dropdown open) 		UI Design	CloudDRxRTR
PAGUI3774	8.12.1.51.3.9.3.19.3 THIRD PARTY LICENSE PANEL (Elements)		Heading	CloudDRxRTR CloudDRxPS
PAGUI3779	8.12.1.51.3.9.3.19.3.1 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3780	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI3781	8.12.1.51.3.9.3.19.3.2 THIRD PARTY HEADER COMPONENT		Heading	CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3782	Component: THIRD PARTY HEADER COMPONENT State: N/A Display: Always displayed in the THIRD PARTY LICENSE PANEL. Label: <label.gb.third_party_license_title> Behavior: Text.	PString715 Third Party Licenses thirdPartyLicenses label.gb.third_party_license_title	UI Info	CloudDRxRTR CloudDRxPS
PAGUI3783	<Picture> THIRD PARTY HEADER COMPONENT Third Party Licenses		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3787	8.12.1.51.3.9.3.19.3.4 THIRD PARTY NAME COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4670	Component: THIRD PARTY NAME COMPONENT State: N/A Display: Displayed in the THIRD PARTY LICENSE PANEL. Text is pulled from Application. Label: N/A Behavior: Text.		UI Info	CloudDRxRTR
PAGUI4671	<Picture> THIRD PARTY NAME COMPONENT convert		UI Design	CloudDRxRTR
PAGUI4672	8.12.1.51.3.9.3.19.3.5 THIRD PARTY DROPDOWN ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4674	Component: THIRD PARTY LIST ITEM DESCRIPTION COMPONENT State: Up/Down Display: THIRD PARTY DROPDOWN ACTION TRIGGER COMPONENT is displayed for each THIRD PARTY NAME COMPONENT. By default, it is displayed in the Down state. Label: N/A Behavior: Icon, and trigger. Upon activation in Down state, the THIRD PARTY DROPDOWN ACTION TRIGGER COMPONENT will change to Up state and displays the THIRD PARTY LIST ITEM DESCRIPTION COMPONENT specific to the THIRD PARTY NAME COMPONENT. Upon activation in Up state, the THIRD PARTY ACTION TRIGGER COMPONENT will change to Down state and hide the THIRD PARTY LIST ITEM DESCRIPTION COMPONENT specific to the THIRD PARTY NAME COMPONENT.		UI Info	CloudDRxRTR
PAGUI4673	<Picture> THIRD PARTY DROPDOWN ACTION TRIGGER COMPONENT (default: down) ▼		UI Design	CloudDRxRTR
PAGUI4675	8.12.1.51.3.9.3.19.3.6 THIRD PARTY LIST ITEM DESCRIPTION COMPONENT		Heading	CloudDRxRTR
PAGUI4676	Component: THIRD PARTY LIST ITEM DESCRIPTION COMPONENT State: Hidden/Visible Display: THIRD PARTY LIST ITEM DESCRIPTION COMPONENT is in Visible state and displayed for each THIRD PARTY NAME COMPONENT when THIRD PARTY ACTION TRIGGER COMPONENT is in the Up State. By default, it is in hidden state and not displayed because THIRD PARTY ACTION TRIGGER COMPONENT is in the Down state by default.		UI Info	CloudDRxRTR

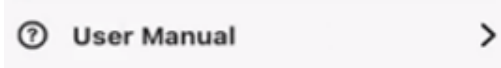


myCardioMEMS Patient Application Graphical User Interface Specification

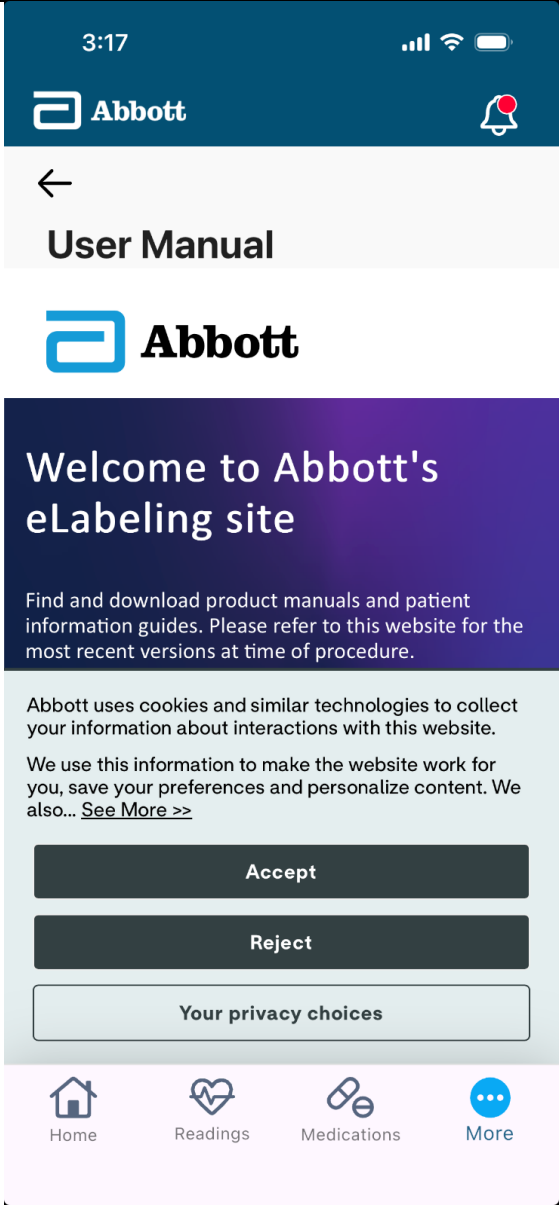
ID		English US Strings	Object Type	Project
	Label: N/A Behavior: Text. The THIRD PARTY LIST ITEM DESCRIPTION COMPONENT specific to the THIRD PARTY NAME COMPONENT will be in Visible state and displayed when the THIRD PARTY DROPDOWN ACTION TRIGGER COMPONENT is in the Up state. The THIRD PARTY LIST ITEM DESCRIPTION COMPONENT specific to the THIRD PARTY NAME COMPONENT will be in Hidden state and not displayed when the THIRD PARTY DROPDOWN ACTION TRIGGER COMPONENT is in the Down state.			
PAGUI4677	<Picture> THIRD PARTY LIST ITEM DESCRIPTION COMPONENT (example: Visible, THIRD PARTY ACTION TRIGGER COMPONENT: Up) <i>Copyright 2017, the Dart project authors.</i> <i>Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:</i> <i>* Redistributions of source code must retain the above copyright</i> <i>notice, this list of conditions and the following disclaimer. * Redistributions in binary form must reproduce the above</i> <i>copyright notice, this list of conditions and the following disclaimer in the documentation and/</i>		UI Design	CloudDRxRTR
PAGUI3790	8.12.1.51.3.9.3.19.3.7 DONE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3791	See Shared Elements: Components.		Info	CloudDRxRTR CloudDRxPS
PAGUI3495	8.12.1.51.3.10 USER MANUAL ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI3496	Component: USER MANUAL ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the ABOUT PANEL. Label: <label.gb.user_manual> Behavior: Icon, text, and trigger. Upon activation, the USER MANUAL ACTION TRIGGER COMPONENT displays the USER MANUAL WEB VIEW.	PString675 User Manual userManual label.gb.user_manual	UI Info	CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3497	<Picture> USER MANUAL ACTION TRIGGER COMPONENT 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3534	8.12.1.51.3.11 USER MANUAL WEB VIEW		Heading	CloudDRxRTR CloudDRxPS
PAGUI4652	8.12.1.51.3.11.1 USER MANUAL WEB VIEW (DISPLAY)		Heading	CloudDRxRTR
PAGUI3535	The USER MANUAL WEB VIEW shows the Abbott eLabeling site upon triggering the USER MANUAL ACTION TRIGGER COMPONENT.		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4655	The USER MANUAL WEB VIEW consists of: OS STATUS PANEL (Shared Elements: Panels) GLOBAL HEADING PANEL (Shared Elements: Panels) BACK ACTION TRIGGER COMPONENT (Shared Elements: Panels) USER MANUAL LABEL COMPONENT USER MANUAL WEB VIEW DISPLAY GLOBAL NAVIGATION PANEL (Shared Elements: Panels)		UI Info	CloudDRxRTR
PAGUI4653	8.12.1.51.3.11.2 USER MANUAL WEB VIEW (DISPLAY)		Heading	CloudDRxRTR
PAGUI3536	<Picture> USER MANUAL WEB VIEW (example: landing page)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

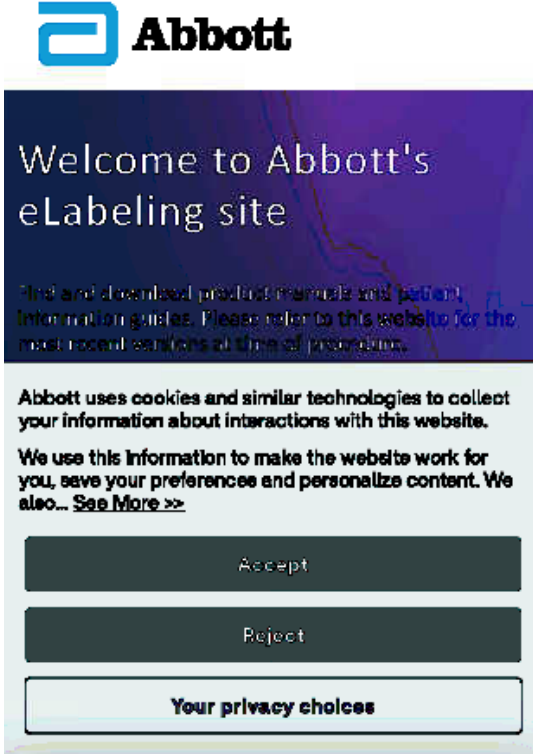
ID		English US Strings	Object Type	Project
				
PAGUI4654	8.12.1.51.3.11.3 USER MANUAL WEB VIEW (ELEMENTS)		Heading	CloudDRxRTR
PAGUI4656	8.12.1.51.3.11.3.1 OS STATUS PANEL		Heading	CloudDRxRTR
PAGUI4657	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI4658	8.12.1.51.3.11.3.2 GLOBAL HEADING PANEL		Heading	CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4659	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI4660	8.12.1.51.3.11.3.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4661	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI4664	8.12.1.51.3.11.3.4 USER MANUAL LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI4665	Component: USER MANUAL ACTION LABEL COMPONENT State: N/A Display: Always displayed in the USER MANUAL WEB VIEW. Label: <label.gb.user_manual> Behavior: Text.		UI Info	CloudDRxRTR
PAGUI4666	<Picture> USER MANUAL ACTION LABEL COMPONENT User Manual		UI Design	CloudDRxRTR
PAGUI4662	8.12.1.51.3.11.3.5 USER MANUAL WEB VIEW DISPLAY		Heading	CloudDRxRTR
PAGUI4663	Component: USER MANUAL WEB VIEW DISPLAY State: N/A Display: Always displayed in the USER MANUAL WEB VIEW Label: N/A Behavior: This is displaying the Abbott eLabeling site landing page: https://www.eifu.abbott/index.html .		UI Info	CloudDRxRTR
PAGUI4667	<Picture> USER MANUAL WEB VIEW DISPLAY (example)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

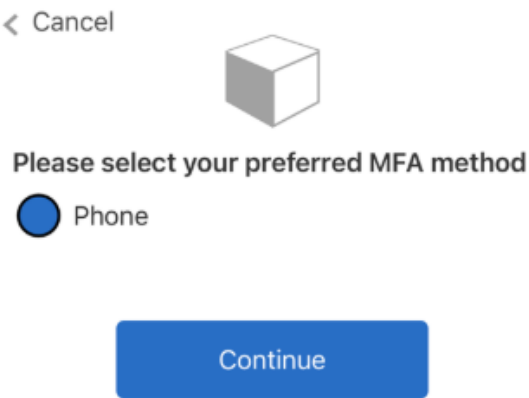
ID		English US Strings	Object Type	Project
				
PAGUI4668	8.12.1.51.3.11.3.6 GLOBAL NAVIGATION PANEL		Heading	CloudDRxRTR
PAGUI4669	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI2602	8.12.1.52 TWO FACTOR AUTHENTICATION PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2603	8.12.1.52.1 TWO FACTOR AUTHENTICATION PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2606	The TWO FACTOR AUTHENTICATION PANEL is displayed after the user selects the REGISTER ACTION TRIGGER COMPONENT or the LOGIN ACTION TRIGGER COMPONENT. The TWO FACTOR AUTHENTICATION PANEL is displayed over the entire screen except for the OS STATUS PANEL. The TWO FACTOR AUTHENTICATION PANEL is a series of panels. The USER DETAILS EMAIL ID PANEL is controlled by a third party and configurable for the Patient Application; the third party may change design details like strings, font, or icons. This section is for reference only.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2609	The TWO FACTOR AUTHENTICATION PANEL is an element of: USER AUTHENTICATION SCREEN PATIENT IDENTIFICATION PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2611	The TWO FACTOR AUTHENTICATION PANEL consists of: CONFIRM MFA METHOD PANEL SEND MFA CODE PANEL ENTER MFA CODE PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2604	8.12.1.52.2 TWO FACTOR AUTHENTICATION PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3306	N/A - TWO FACTOR AUTHENTICATION PANEL is a series of panels.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2605	8.12.1.52.3 TWO FACTOR AUTHENTICATION PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3307	8.12.1.52.3.1 CONFIRM MFA METHOD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3308	8.12.1.52.3.1.1 CONFIRM MFA METHOD PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3311	The CONFIRM MFA METHOD PANEL is displayed after the user selects the LOGIN ACTION TRIGGER COMPONENT. The CONFIRM MFA METHOD PANEL is displayed as a web browser over the entire screen except the OS STATUS PANEL. The CONFIRM MFA METHOD PANEL is controlled by a third party vendor and this section is for reference only.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3312	The CONFIRM MFA METHOD PANEL consists of: CANCEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) USER DETAILS ICON COMPONENT (See Shared Elements: Components) CONFIRM MFA INSTRUCTIONS COMPONENT MFA METHOD RADIO BUTTON COMPONENT CONTINUE ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3309	8.12.1.52.3.1.2 CONFIRM MFA METHOD PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3313	<Picture> CONFIRM MFA METHOD PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI3310	8.12.1.52.3.1.3 CONFIRM MFA METHOD PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2613	8.12.1.52.3.1.3.1 CANCEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2614	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2624	8.12.1.52.3.1.3.2 USER DETAILS ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2625	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3315	8.12.1.52.3.1.3.3 CONFIRM MFA INSTRUCTIONS COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3316	Component: CONFIRM MFA INSTRUCTIONS COMPONENT State: N/A Display: Always displayed in the CONFIRM MFA METHOD PANEL. Label: N/A - determined by a third party vendor Behavior: Text.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3317	<Picture> CONFIRM MFA INSTRUCTIONS COMPONENT Please select your preferred MFA method		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3318	8.12.1.52.3.1.3.4 MFA METHOD RADIO BUTTON COMPONENT		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3319	Component: MFA METHOD RADIO BUTTON COMPONENT State: Selected Display: Always displayed in the selected state in the CONFIRM MFA METHOD PANEL. Label: N/A - determined by a third party vendor Behavior: Trigger and text. Upon activation in the selected state, the MFA METHOD RADIO BUTTON COMPONENT does nothing. The MFA METHOD RADIO BUTTON COMPONENT is always Selected.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3320	<Picture> MFA METHOD RADIO BUTTON COMPONENT (default: Selected, example: Phone)  Phone		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4691	8.12.1.52.3.1.3.5 CONTINUE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4692	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI3324	8.12.1.52.3.2 SEND MFA CODE PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3325	8.12.1.52.3.2.1 SEND MFA CODE PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3326	The SEND MFA CODE PANEL is displayed after the user selects the CONTINUE ACTION TRIGGER COMPONENT in the CONFIRM MFA METHOD PANEL. The SEND MFA CODE PANEL is displayed as a web browser over the entire screen except the OS STATUS PANEL. The SEND MFA CODE PANEL is controlled by a third party vendor and this section is for reference only.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3327	The SEND MFA CODE PANEL consists of: CANCEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) USER DETAILS ICON COMPONENT (See Shared Elements: Components) MFA PHONE LABEL COMPONENT (See Shared Elements: Components) LOADING BAR ICON COMPONENT MFA PHONE ERROR COMPONENT SEND CODE ACTION TRIGGER COMPONENT CALL ME ACTION TRIGGER COMPONENT CONTINUE ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3328	8.12.1.52.3.2.2 SEND MFA CODE PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3329	<Picture> SEND MFA CODE PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	< Cancel  We have the following number on record for you. We can send a code via SMS or phone to authenticate you. ***-***-456 Send Code Call Me			
PAGUI3330	<Picture> SEND MFA CODE PANEL (example: Error Message)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

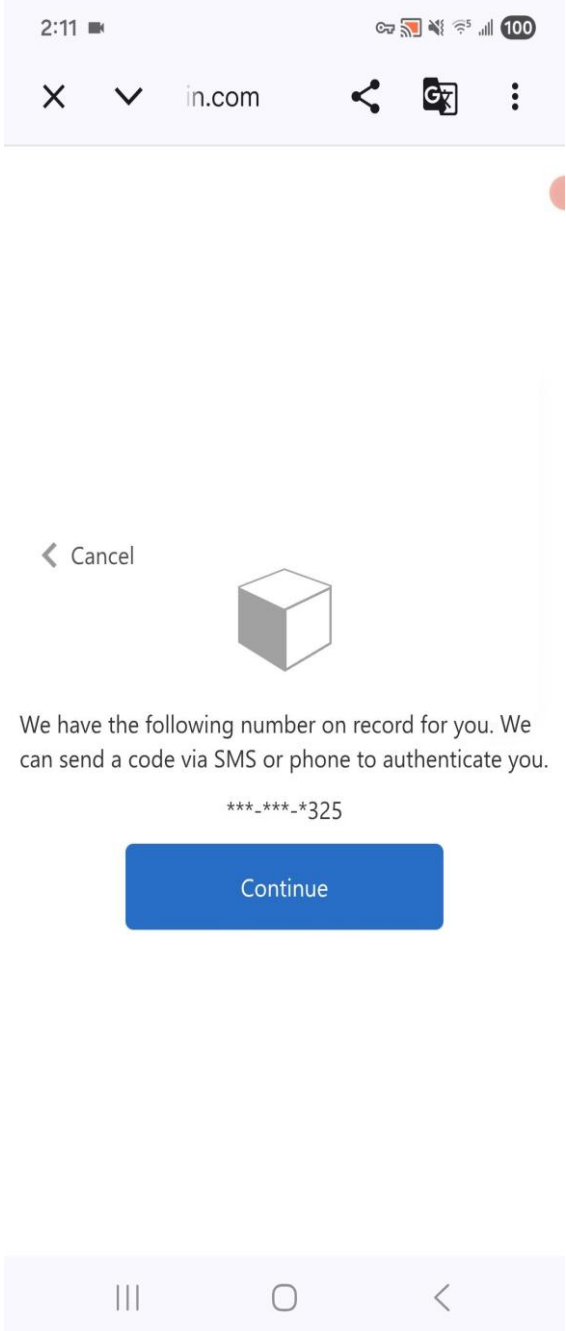
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	< Cancel  We have the following number on record for you. We can send a code via SMS or phone to authenticate you. ***-***-*456 The phone number you provided is busy or unavailable. Please check the number and try again. Send Code Call Me			
PAGUI3570	<Picture> SEND MFA CODE PANEL (example: Loading)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	< Cancel  We have the following number on record for you. We can send a code via SMS or phone to authenticate you. ***-***-818 * * * * Send Code Call Me			
PAGUI4890	<Picture> SEND MFA CODE PANEL (example: Android OS, Continue after Call Me)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3331	8.12.1.52.3.2.3 SEND MFA CODE PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3332	8.12.1.52.3.2.3.1 CANCEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3333	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3334	8.12.1.52.3.2.3.2 USER DETAILS ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3335	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4693	8.12.1.52.3.2.3.3 MFA PHONE LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI4694	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI3571	8.12.1.52.3.2.3.5 LOADING BAR ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3572	Component: LOADING BAR ICON COMPONENT State: N/A Display: Conditionally displayed in the SEND MFA CODE PANEL when the SEND MFA CODE PANEL is loading information. Loading information will commonly occur after the user selects the CALL ME ACTION TRIGGER COMPONENT while the panel is waiting for the call to complete, but it may also occur for other reasons. Label: N/A - determined by a third party vendor Behavior: Icon. The LOADING BAR ICON COMPONENT is an animated icon to demononstrate to the user that loading information is actively underway. The dots move from left to right and back.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3573	<Picture> LOADING BAR ICON COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3339	8.12.1.52.3.2.3.6 MFA PHONE ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3340	Component: MFA PHONE ERROR COMPONENT State: N/A Display: Conditionally displayed in the SEND MFA CODE PANEL when the user's phone number is not available to received a phone call. Label: N/A - determined by a third party vendor Behavior: Text		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3341	<Picture> MFA PHONE ERROR COMPONENT The phone number you provided is busy or unavailable. Please check the number and try again.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3342	8.12.1.52.3.2.3.7 SEND CODE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3343	Component: SEND CODE ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Always displayed in the SEND MFA CODE PANEL. The default state is enabled. Label: N/A - determined by a third party vendor Behavior: Text and trigger. Upon activation, the SEND CODE ACTION TRIGGER COMPONENT triggers the third party vendor to send a text message to the user's phone number and displays the ENTER CODE PANEL. The SEND CODE ACTION TRIGGER COMPONENT changes from enabled to disabled when the SEND MFA CODE PANEL is loading information.	PString974 Send Code ots.authentication.send_code	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3344	<Picture> SEND CODE ACTION TRIGGER COMPONENT (default: Enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3574	<Picture> SEND CODE ACTION TRIGGER COMPONENT (example: Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3345	8.12.1.52.3.2.3.8 CALL ME ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3346	Component: CALL ME ACTION TRIGGER COMPONENT State: Enabled/ Disabled Display: Always displayed in the SEND MFA CODE PANEL. The default state is enabled. Label: N/A - determined by a third party vendor Behavior: Text and trigger. Upon activation, the CALL ME ACTION TRIGGER COMPONENT triggers the third party vendor to call the user's phone number. If verification fails, display SEND MFA CODE PANEL. If verification is successful and the user is in the REGISTRATION SCREEN, conditionally display the CONNECT HEALTH APP PANEL (if supported by Mobile OS). If verification is successful and the user is in the USER AUTHENTICATION SCREEN, display either the HOME SCREEN (if supported by Mobile OS), or display the CONTINUE AFTER MFA CALL PANEL. The CALL ME ACTION TRIGGER COMPONENT changes from enabled to disabled when the SEND MFA CODE PANEL is loading information.	PString975 Call Me ots.authentication.call_me	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3347	<Picture> CALL ME ACTION TRIGGER COMPONENT (default: Enabled)		UI Design	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3575	<Picture> CALL ME ACTION TRIGGER COMPONENT (example: Disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4910	8.12.1.52.3.2.3.9 CONTINUE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4912	Component: CONTINUE ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the SEND MFA CODE PANEL, after user selects CALL ME ACTION TRIGGER COMPONENT. Display varies based on OS. Label: N/A - determined by a third party vendor Behavior: Text and Trigger. Upon activation after confirming login from the phone call, the CONTINUE ACTION TRIGGER COMPONENT displays the HOME SCREEN..		UI Info	CloudDRxRTR CloudDRxPS
PAGUI4913	<Picture> CONTINUE ACTION TRIGGER COMPONENT (default: Enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI3348	8.12.1.52.3.3 ENTER MFA CODE PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3349	8.12.1.52.3.3.1 ENTER MFA CODE PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3350	The ENTER MFA CODE PANEL is displayed after the user selects the SEND CODE ACTION TRIGGER COMPONENT. The ENTER MFA CODE PANEL is displayed as a web browser over the entire screen. The ENTER MFA CODE PANEL is controlled by a third party vendor and this section is for reference only.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3351	The ENTER MFA CODE PANEL consists of: CANCEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) USER DETAILS ICON COMPONENT (See Shared Elements: Components) MFA PHONE LABEL COMPONENT (See Shared Elements: Component) MFA INSTRUCTIONS ACTION TRIGGER COMPONENT VERIFICATION CODE ERROR LABEL COMPONENT MFA MAXIMUM ATTEMPTS DIALOG PANEL ENTER CODE RECEIVED ERROR COMPONENT VERIFICATION CODE INPUT COMPONENT (See Shared Elements: Components) VERIFY CODE ACTION TRIGGER COMPONENT (See Shared Elements: Components) ALPHABETICAL KEYBOARD COMPONENT (See Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3352	8.12.1.52.3.3.2 ENTER MFA CODE PANEL (Layout)		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3353	<Picture> ENTER MFA CODE PANEL < Cancel  We have the following number on record for you. We can send a code via SMS or phone to authenticate you. ***-***-*456 Enter your verification code below, or send a new code 546151 Verify Code		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3354	<Picture> ENTER MFA CODE PANEL (example: Error Message)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	< Cancel  We have the following number on record for you. We can send a code via SMS or phone to authenticate you. ***-***-456 Enter your verification code below, or send a new code The verification code you have entered does not match our records. Please try again, or request a new code. Verify Code			
PAGUI3355	8.12.1.52.3.3.3 ENTER MFA CODE PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3356	8.12.1.52.3.3.3.1 CANCEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3357	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3358	8.12.1.52.3.3.3.2 USER DETAILS ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI3359	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3369	8.12.1.52.3.3.3.3 MFA PHONE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3370	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3360	8.12.1.52.3.3.3.4 MFA INSTRUCTIONS ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3361	Component: MFA INSTRUCTIONS ACTION TRIGGER COMPONENT State: Black Text/ Blue Text Display: Always displayed in the ENTER MFA CODE PANEL. The MFA INSTRUCTIONS ACTION TRIGGER COMPONENT always displays both states. Label: N/A - determined by a third party vendor Behavior: Text and trigger. The first section of the text is always in the black text state. The second section of text is always in the blue text state. Upon activation in the blue text state, the MFA INSTRUCTIONS ACTION TRIGGER COMPONENT triggers the third party vendor to send another text message to the user's phone number.	PString976 Enter your verification code below, or send a new code ots.authentication.enter_or_send_new_verification_code	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3362	<Picture> MFA INSTRUCTIONS ACTION TRIGGER COMPONENT Enter your verification code below, or send a new code		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3363	8.12.1.52.3.3.3.5 VERIFICATION CODE ERROR LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3364	Component: VERIFICATION CODE ERROR LABEL COMPONENT State: N/A Display: Conditionally displayed in the ENTER MFA CODE PANEL when the user inputs an incorrect verification code. Label: N/A - determined by a third party vendor Behavior: Text.	PString978 The verification code you have entered does not match our records. Please try again, or request a new code. ots.authentication.verification_code_mismatch	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3365	<Picture> VERIFICATION CODE ERROR LABEL COMPONENT The verification code you have entered does not match our records. Please try again, or request a new code.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4695	8.12.1.52.3.3.3.7 ENTER CODE RECEIVED ERROR COMPONENT		Heading	CloudDRxRTR
PAGUI4696	Component: ENTER CODE RECEIVED ERROR COMPONENT State: N/A Display: Conditionally displayed in the ENTER MFA CODE PANEL when the user has entered a code and then deleted it. Label: N/A - determined by a third party vendor Behavior: Text.	PString1083 Please enter the verification code you received ots.authentication.enter_verification_code	UI Info	CloudDRxRTR
PAGUI4697	<Picture> ENTER CODE RECEIVED ERROR COMPONENT Please enter the verification code you received		UI Design	CloudDRxRTR
PAGUI3374	8.12.1.52.3.3.3.8 VERIFICATION CODE INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3375	See Shared Elements: Component		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3378	8.12.1.52.3.3.3.9 VERIFY CODE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3379	See Shared Elements: Components		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3382	8.12.1.52.3.3.3.10 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3383	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4702	8.12.1.52.3.3.3.11 MFA MAXIMUM ATTEMPTS DIALOG PANEL		Heading	CloudDRxRTR
PAGUI4703	See Shared Elements: Panels		UI Info	CloudDRxRTR
PAGUI3105	8.12.1.53 USER DETAILS NEW PASSWORD PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3106	8.12.1.53.1 USER DETAILS NEW PASSWORD PANEL (DISPLAY)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3107	The USER DETAILS NEW PASSWORD PANEL is displayed after the user selects the CONTINUE ACTION TRIGGER COMPONENT in the USER DETAILS CONFIRMATION PANEL. The user is able to set and confirm a new password out of the Patient Application.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3108	The USER DETAILS NEW PASSWORD PANEL consists of:		UI Info	CloudDRx 1.0

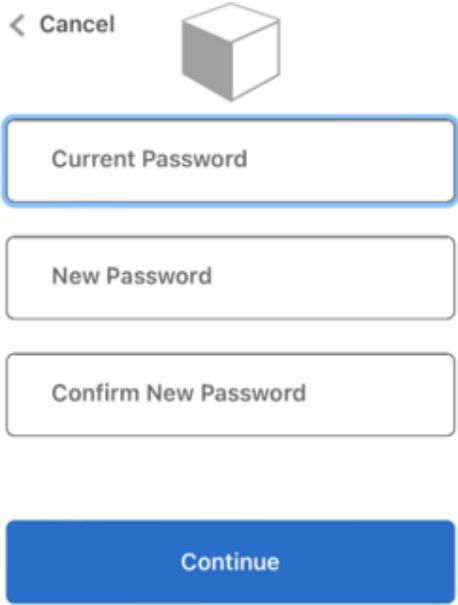
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	CANCEL ACTION TRIGGER COMPONENT (See Shared Elements: Components) USER DETAILS ICON COMPONENT (See Shared Elements: Components) REQUIRED FIELD ERROR COMPONENT (See Shared Elements: Components) INFORMATION REQUIRED ERROR COMPONENT (See Shared Elements: Components) USER DETAILS INVALID PASSWORD ERROR COMPONENT USER DETAILS DIFFERENT PASSWORD ERROR COMPONENT USER DETAILS PASSWORD MATCH ERROR COMPONENT USER DETAILS PASSWORD INCORRECT ERROR COMPONENT USER DETAILS OLD PASSWORD INPUT COMPONENT USER DETAILS NEW PASSWORD INPUT COMPONENT USER DETAILS CONFIRM PASSWORD INPUT COMPONENT ALPHABETICAL KEYBOARD COMPOENT (See Shared Elements: Components) CONTINUE ACTION TRIGGER COMPONENT (See Shared Elements: Components)			CloudDRxRTR CloudDRxPS
PAGUI3109	8.12.1.53.2 USER DETAILS NEW PASSWORD PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3110	<Picture> USER DETAILS NEW PASSORD PANEL (default: No Input, example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 The image shows a mobile application screen for changing a password. At the top left is a back arrow and the word 'Cancel'. Below this is a 3D cube icon. There are three text input fields stacked vertically: 'Current Password', 'New Password', and 'Confirm New Password'. At the bottom is a blue button labeled 'Continue'.			
PAGUI3111	<Picture> USER DETAILS NEW PASSORD PANEL (Input, example: error)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	< Cancel  A required field is missing. Please fill out all required fields and try again. This information is required. Current Password This information is required. New Password This information is required. Confirm New Password Continue			
PAGUI3112	8.12.1.53.3 USER DETAILS NEW PASSWORD PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3113	8.12.1.53.3.1 CANCEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3114	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3115	8.12.1.53.3.2 USER DETAILS ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3116	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

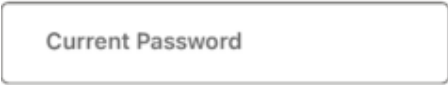
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI3117	8.12.1.53.3.3 USER DETAILS HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3118	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4111	8.12.1.53.3.4 REQUIRED FIELD ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4112	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3125	8.12.1.53.3.5 INFORMATION REQUIRED ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3126	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4113	8.12.1.53.3.6 USER DETAILS INVALID PASSWORD ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4114	Component: USER DETAILS INVALID PASSWORD ERROR COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS NEW PASSWORD PANEL when the input for the USER DETAILS OLD PASSWORD INPUT COMPONENT does not match the existing password. Label: N/A - determined by a third party vendor - <ots.config.invalid_credentials> Behavior: Text.	PString934 Invalid username or password. ots.config.invalid_credentials	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4115	<Picture> USER DETAILS INVALID PASSWORD ERROR COMPONENT Invalid username or password.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4116	8.12.1.53.3.7 USER DETAILS DIFFERENT PASSWORD ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4117	Component: USER DETAILS DIFFERENT PASSWORD ERROR COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS NEW PASSWORD PANEL when the input for the USER DETAILS NEW PASSWORD INPUT COMPONENT or the USER DETAILS CONFIRM PASSWORD INPUT COMPONENT is the same as the existing password. Label: N/A - determined by a third party vendor - <ots.config.different_password>	PString935 Please choose a new password that is different from your old password. UserMessageIfOldPassword Used	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

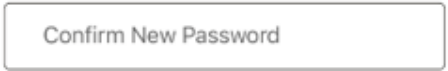


myCardioMEMS Patient Application Graphical User Interface Specification				
ID		English US Strings	Object Type	Project
	Behavior: Text.	ots.config.different_password		
PAGUI4118	<Picture> USER DETAILS DIFFERENT PASSWORD ERROR COMPONENT Please choose a new password that is different from your old password.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3119	8.12.1.53.3.8 USER DETAILS PASSWORD MATCH ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3120	Component: USER DETAILS PASSWORD MATCH ERROR COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS NEW PASSWORD PANEL when the input for the USER DETAILS NEW PASSWORD INPUT COMPONENT and the USER DETAILS CONFIRM PASSWORD INPUT COMONENT do not match. Label: N/A - determined by a third party vendor - <ots.config.password_mismatch> Behavior: Text.	PString936 The password entry fields do not match. Please enter the same password in both fields and try again. ots.config.password_mismatch	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3121	<Picture> USER DETAILS PASSWORD MATCH ERROR COMPONENT The password entry fields do not match. Please enter the same password in both fields and try again.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3122	8.12.1.53.3.9 USER DETAILS PASSWORD INCORRECT ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3123	Component: USER DETAILS PASSWORD INCORRECT COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS NEW PASSWORD PANEL when the input for the USER DETAILS NEW PASSWORD INPUT COMPONENT or the USER DETAILS CONFIRM PASSWORD INPUT COMONENT are filled out incorrectly. Label: N/A - determined by a thrid party vendor - <ots.config.password_rules> Behavior: Text.	PString937 The password entered is invalid. Please enter a valid password. A valid password must follow these rules: Must have between 8 and 23 characters -Must contain the following 4 character groups: Upper case letters (A,B,C, ... Z) Lower case letters (a, b, c, ... z) Numbers (0, 1, 2, ... 9) non-alphanumeric (~ ! # \$ ^ () - _ + { } [] ; , . ? ii) Must not include more than 2 consecutive characters (for example: 123 or abc)	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		ots.config.password_rules		
PAGUI3124	<Picture> USER DETAILS PASSWORD INCORRECT ERROR COMPONENT The password entered is invalid. Please enter a valid password. A valid password must follow these rules: Must have between 8 and 23 characters -Must contain the following 4 character groups: Upper case letters (A, B, C, ... Z) Lower case letters (a, b, c, ... z) Numbers (0, 1, 2, ... 9) non-alphanumeric (~ ! # \$ ^ () - _ + { } [] : ; , . ? ii) Must not include more than 2 consecutive characters (for example: 123 or abc)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4119	8.12.1.53.3.10 USER DETAILS OLD PASSWORD INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4120	Component: USER DETAILS NEW PASSWORD INPUT COMPONENT State: No Input/ Input/ Error Display: Always displayed in the USER DETAILS NEW PASSWORD PANEL. The default state is no input. Label: N/A - determined by a third party vendor Behavior: Text and input display. In the no input state, the USER DETAILS OLD PASSWORD INPUT COMPONENT displays label. In the input state, the USER DETAILS OLD PASSWORD INPUT COMPONENT displays the obscured input from the ALPHABETICAL KEYBOARD COMPONENT. The USER DETAILS OLD PASSWORD INPUT COMPONENT changes to the error state after the user selects the CONTINUE ACTION TRIGGER COMPONENT and validation is not successful. In the error state, the USER DETAILS OLD PASSWORD INPUT displays the obscured input from the ALPHABETICAL KEYBOARD COMPONENT and a different outline color.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4121	<Picture> USER DETAILS OLD PASSWORD INPUT COMPONENT (default: No Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4122	<Picture> USER DETAILS NEW PASSWORD INPUT COMPONENT (Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4123	<Picture> USER DETAILS NEW PASSWORD INPUT COMPONENT (Error, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3127	8.12.1.53.3.11 USER DETAILS NEW PASSWORD INPUT COMPONENT		Heading	CloudDRx 1.0

myCardioMEMS Patient Application Graphical User Interface Specification

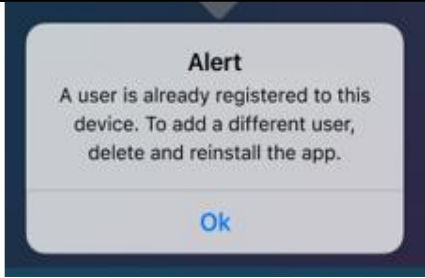
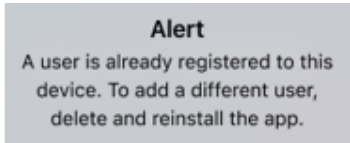
ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3128	Component: USER DETAILS NEW PASSWORD INPUT COMPONENT State: No Input/ Input/ Error Display: Always displayed in the USER DETAILS NEW PASSWORD PANEL. The default state is no input. Label: N/A - determined by a third party vendor Behavior: Text and input display. In the no input state, the USER DETAILS NEW PASSWORD INPUT COMPONENT displays label. In the input state, the USER DETAILS NEW PASSWORD INPUT COMPONENT displays the obscured input from the ALPHABETICAL KEYBOARD COMPONENT. The USER DETAILS NEW PASSWORD INPUT COMPONENT changes to the error state after the user selects the CONTINUE ACTION TRIGGER COMPONENT and validation is not successful. In the error state, the USER DETAILS NEW PASSWORD INPUT displays the obscured input from the ALPHABETICAL KEYBOARD COMPONENT and a different outline color.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3129	<Picture> USER DETAILS NEW PASSWORD INPUT COMPONENT (default: No Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3130	<Picture> USER DETAILS NEW PASSWORD INPUT COMPONENT (Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3131	<Picture> USER DETAILS NEW PASSWORD INPUT COMPONENT (Error, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3132	8.12.1.53.3.12 USER DETAILS CONFRIM PASSWORD INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3133	Component: USER DETAILS CONFRIM PASSWORD INPUT COMPONENT State: No Input/ Input/ Error Display: Always displayed in the USER DETAILS NEW PASSWORD PANEL. The default state is no input. Label: N/A - determined by a third party vendor Behavior: Text and input display. In the no input state, the USER DETAILS CONFRIM PASSWORD INPUT COMPONENT displays label. In the input state, the USER DETAILS CONFRIM PASSWORD INPUT COMPONENT displays the obscured input from the ALPHABETICAL KEYBOARD COMPONENT. The USER DETAILS CONFRIM PASSWORD INPUT COMPONENT changes to the error state after the user selects the CONTINUE ACTION TRIGGER COMPONENT and validation is not successful. In the error state, the USER DETAILS CONFRIM PASSWORD INPUT displays the obscured input from the ALPHABETICAL KEYBOARD COMPONENT and a different outline color.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3134	<Picture> USER DETAILS CONFIRM PASSWORD INPUT COMPONENT (default: No Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3135	<Picture> USER DETAILS CONFIRM PASSWORD INPUT COMPONENT (Input, example)		UI Design	CloudDRx 1.0



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI3136	<Picture> USER DETAILS CONFRIM PASSWORD INPUT COMPONENT (Error, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3137	8.12.1.53.3.13 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3138	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3139	8.12.1.53.3.14 CONTINUE ACTION TRIGGER		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3140	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1169	8.12.1.54 USER SAME DEVICE ERROR PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1170	8.12.1.54.1 USER SAME DEVICE ERROR PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1190	The USER SAME DEVICE ERROR PANEL is conditionally displayed as a component of the STARTUP SCREEN if a registration request fails validation. The USER SAME DEVICE ERROR PANEL is displayed over the STARTUP SCREEN and the background is greyed out. The USER SAME DEVICE ERROR PANEL is conditionally displayed as a component of the USER AUTHENTICATION SCREEN if a login request fails validation from the LOGIN ACTION TRIGGER COMPONENT		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1191	The USER SAME DEVICE ERROR PANEL consists of: USER SAME DEVICE ERROR LABEL COMPONENT OK ACTION TRIGGER COMPONENT (see Shared Elements: Components)		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1171	8.12.1.54.2 USER SAME DEVICE ERROR PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1192	<Picture> USER SAME DEVICE ERROR PANEL		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1172	8.12.1.54.3 USER SAME DEVICE ERROR PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1194	8.12.1.54.3.1 USER SAME DEVICE ERROR LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1198	Component: USER SAME DEVICE ERROR LABEL COMPONENT State: N/A Display: Always displayed in the USER SAME DEVICE ERROR PANEL. Label: <msg.gb.alert>, <msg.gb.same_device_error> Behavior: Text.	PString636 Alert alert msg.gb.alert PString921 A user is already registered to this device. To add a different user, delete and reinstall the app. alreadyRegistered msg.gb.same_device_error	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1199	<Picture> USER SAME DEVICE ERROR LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1195	8.12.1.54.3.2 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1200	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1735	8.12.1.55 VIDEO PLAYER PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1740	8.12.1.55.1 VIDEO PLAYER PANEL (DISPLAY)		Heading	CloudDRx 1.0

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1741	The VIDEO PLAYER PANEL displays a video and video controls over the entire screen.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1742	The VIDEO PLAYER PANEL is an element of: STARTUP SCREEN PAP READING INSTRUCTIONS PANEL ABOUT HEART FAILURE PANEL ABOUT HEART FAILURE SYSTEM PANEL IN-APP TRAINING PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1743	The VIDEO PLAYER PANEL is composed of: Video player elements are determined by SOUP provider. Panel components may vary based on language and OS level settings.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1744	8.12.1.55.2 VIDEO PLAYER PANEL (LAYOUT)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1745	<Picture> VIDEO PLAYER PANEL (example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

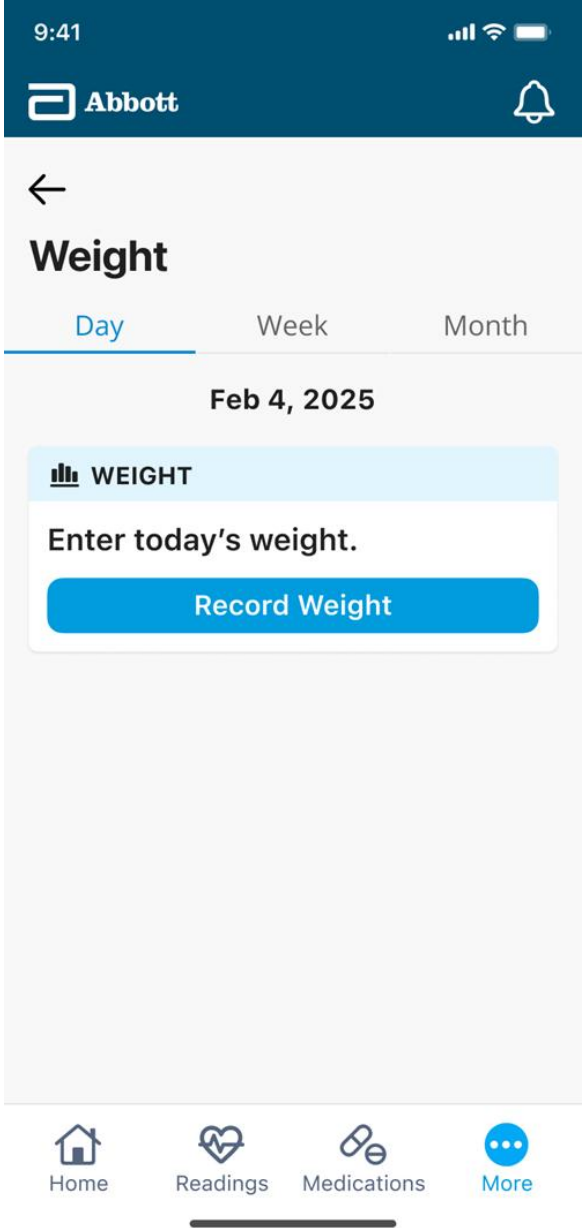
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI1747	8.12.1.55.3 VIDEO PLAYER PANEL (ELEMENTS)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1748	VIDEO PLAYER PANEL elements are determined by the SOUP provider and the layout only serves as an example of what the panel may look like.		UI Info	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI172	8.12.1.56 WEIGHT PANEL		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI173	8.12.1.56.1 WEIGHT PANEL (Display)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2851	The WEIGHT CARD PANEL is always displayed when the WEIGHT ACTION TRIGGER COMPONENT is activated.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI207	The WEIGHT CARD PANEL consists of: GLOBAL HEADING PANEL BACK ACTION TRIGGER COMPONENT WEIGHT HEADER COMPONENT INTERVAL TAB COMPONENT DATE LABEL COMPONENT DATE SELECTION TRIGGER COMPONENT HISTORICAL DATE SELECTION CALENDAR COMPONENT WEIGHT DISPLAY CARD GLOBAL NAVIGATION PANEL		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI174	8.12.1.56.2 WEIGHT PANEL (Layout)		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2891	<Picture> WEIGHT PANEL (example: Default: Day, No Input)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

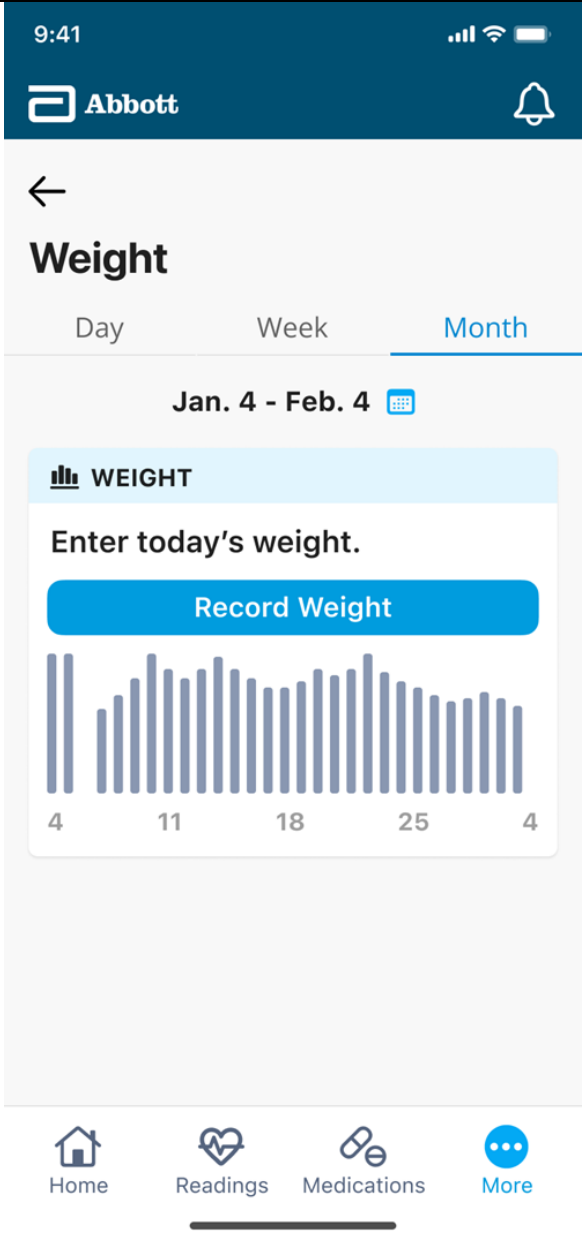
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4350	9:41 Abbott ← Weight Day Week Month Jan. 28 - Feb. 4 WEIGHT 150 152 149 149 154 149 28 29 30 1 2 3 4 Home Readings Medications More		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4349	<Picture> WEIGHT PANEL (example: Month, No Input)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

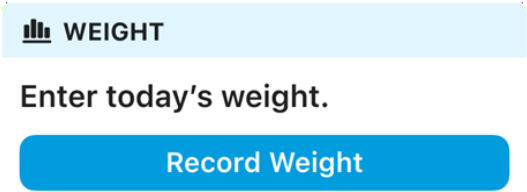
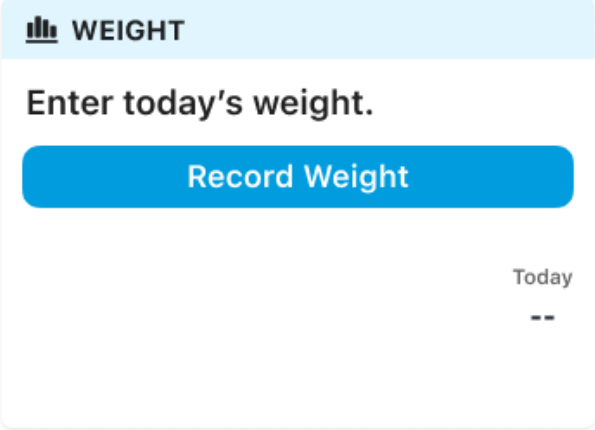
ID		English US Strings	Object Type	Project
				
PAGUI175	8.12.1.56.3 WEIGHT PANEL (Elements)		Heading	CloudDRx 1.0 CloudDRxRTR



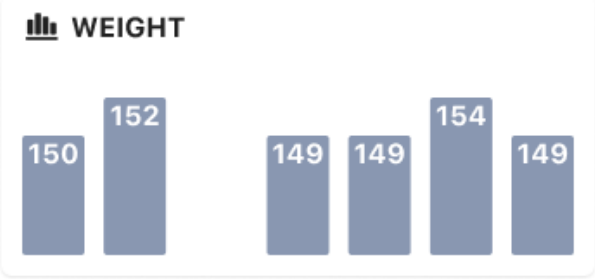
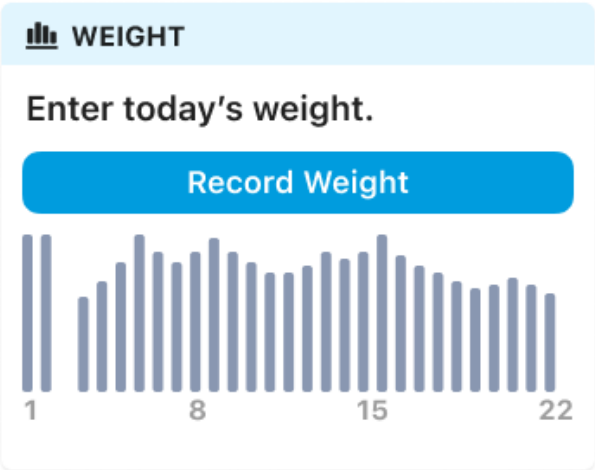
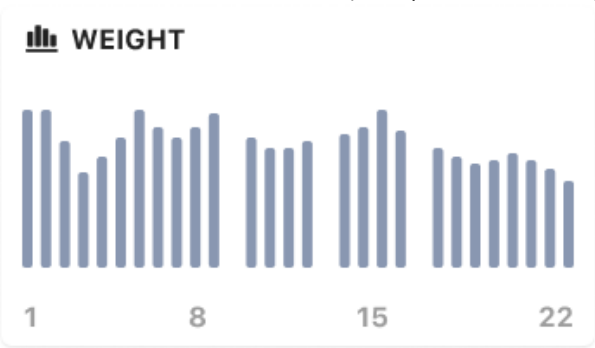
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI4539	8.12.1.56.3.1 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4540	See Shared Elements: Components.		Info	CloudDRxRTR
PAGUI237	8.12.1.56.3.2 WEIGHT HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4542	Component: WEIGHT HEADER COMPONENT State: N/A Display: Always displayed in the WEIGHT PANEL. Label: <label.gb.record_weight_reminder> Behavior: Text.	PString521 Weight weight label.gb.weight	UI Info	CloudDRxRTR
PAGUI4543	<Picture> WEIGHT HEADER COMPONENT 		UI Design	CloudDRxRTR
PAGUI154	8.12.1.56.3.3 INTERVAL TAB COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI155	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI156	8.12.1.56.3.4 DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI157	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI158	8.12.1.56.3.5 DATE SELECTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI159	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI160	8.12.1.56.3.6 HISTORICAL DATE SELECTION CALENDAR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI161	See Shared Elements: Components.		Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI241	8.12.1.56.3.7 WEIGHT DISPLAY CARD		Heading	CloudDRx 1.0

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI4553	8.12.1.56.3.7.1 WEIGHT DISPLAY CARD (Display)		Heading	CloudDRxRTR
PAGUI4554	The WEIGHT DISPLAY CARD is always displayed in the WEIGHT PANEL.		UI Info	CloudDRxRTR
PAGUI4556	The WEIGHT DISPLAY CARD consists of: WEIGHT BANNER COMPONENT WEIGHT REMINDER LABEL COMPONENT RECORD WEIGHT ACTION TRIGGER COMPONENT WEIGHT DISPLAY COMPONENT		UI Info	CloudDRxRTR
PAGUI2885	<Picture> WEIGHT DISPLAY CARD (default: Day, example: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4343	<Picture> WEIGHT DISPLAY CARD (default: Day, example: weight input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2886	<Picture> WEIGHT DISPLAY CARD (example: Week, example: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4348	<Picture> WEIGHT DISPLAY CARD (example: Week, Input)		UI Design	CloudDRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

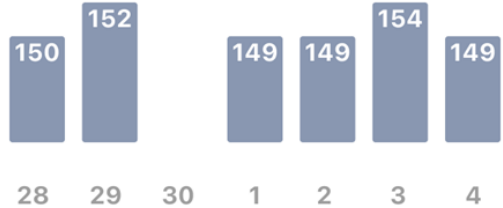
ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI2887	<Picture> WEIGHT DISPLAY CARD (example: Month, example: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4342	<Picture> WEIGHT DISPLAY CARD (example: Month, example:Weight entered) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4555	8.12.1.56.3.7.3 WEIGHT DISPLAY CARD (Elements)		Heading	CloudDRxRTR
PAGUI4583	8.12.1.56.3.7.3.1 WEIGHT BANNER COMPONENT		Heading	CloudDRxRTR
PAGUI238	Component: WEIGHT BANNER COMPONENT State: Enabled/Disabled	PString521 Weight weight	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Display: Always displayed in the WEIGHT CARD PANEL. If daily reading has been taken, display the banner text with white background as read-only (disabled). If daily reading has not been taken, banner header text with blue background and allow component activation upon click (enabled). The enabled state is the default state. Label: <label.gb.weight>; Text is all capitalized. Behavior: Icon and text.	label.gb.weight		
PAGUI4551	<Picture> WEIGHT BANNER COMPONENT (default: enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4557	<Picture> WEIGHT BANNER COMPONENT (disabled) 		UI Design	CloudDRxRTR
PAGUI239	8.12.1.56.3.7.3.2 WEIGHT REMINDER LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI240	Component: WEIGHT REMINDER LABEL COMPONENT State: N/A Display: Conditionally displayed in the WEIGHT CARD PANEL. If daily reading has been taken, remove display of WEIGHT REMINDER LABEL COMPONENT. If daily reading has not been taken, display WEIGHT REMINDER LABEL COMPONENT. Label: <label.gb.record_weight_reminder> Behavior: Text.	PString431 Enter today's weight. enterTodaysWeight label.gb.record_weight_re minder	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2889	<Picture> WEIGHT REMINDER LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI243	8.12.1.56.3.7.3.3 RECORD WEIGHT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI244	Component: RECORD WEIGHT ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the WEIGHT CARD PANEL. Default is to display the RECORD WEIGHT ACTION TRIGGER COMPONENT when daily reading is not taken. If daily reading has been taken, remove display of RECORD WEIGHT ACTION TRIGGER COMPONENT. Label: <labe.gb.record_weight> Behavior: Text and trigger. Upon activation, display the RECORD WEIGHT PANEL.	PString432 Record Weight recordWeight label.gb.record_weight	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2890	<Picture> RECORD WEIGHT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4558	8.12.1.56.3.7.3.4 WEIGHT DISPLAY COMPONENT		Heading	CloudDRxRTR
PAGUI4559	Component: WEIGHT DISPLAY COMPONENT State: Day/Week/Month, Reading taken/Reading not taken, kg/lb	PString342 kg kg	UI Info	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

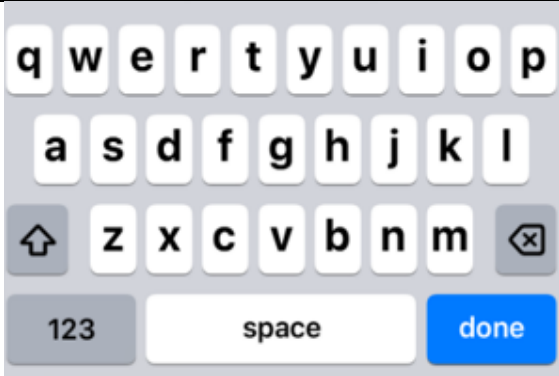
ID		English US Strings	Object Type	Project
	Display: Always displayed in the WEIGHT PANEL. Default state is Day. In the Day state, there is no display until a reading is taken and the current date's data is always displayed. In the Week and Month states, display the number of available readings per the interval specified in the INTERVAL SELECTION COMPONENT and DATE SELECTION COMPONENT. Label: <data.Code.Med_Unit_Cd.Code_Desc.KG> <data.Code.Weight_Unit_LB> Behavior: Graph or text. Day state: If no reading taken, nothing is displayed. If reading is taken, weight value is displayed with the weight units (kg/lb) depending on region user is registered on with their mobile device. Week: If no reading taken, graph is displayed with missing bar and dashes. If reading taken, bar is displayed with weight value. All available readings in the given period will be displayed. X-axis is labeled with dates. Month: If no reading taken, graph is displayed with missing bar. If reading taken, bar is displayed. All available readings in the given period will be displayed. X-axis is labeled with dates.	data.Code.Med_Unit_Cd.Code_Desc.KG PString520 lb lb data.Code.Weight_Unit_LB		
PAGUI4560	<Picture> WEIGHT DISPLAY COMPONENT (day, reading taken, lb) 		UI Design	CloudDRxRTR
PAGUI4575	<Picture> WEIGHT DISPLAY COMPONENT (day, reading taken, kg) 		UI Design	CloudDRxRTR
PAGUI4561	<Picture> WEIGHT DISPLAY COMPONENT (week, reading not taken) 		UI Design	CloudDRxRTR
PAGUI4562	<Picture> WEIGHT DISPLAY COMPONENT (week, reading taken) 		UI Design	CloudDRxRTR
PAGUI4563	<Picture> WEIGHT DISPLAY COMPONENT (month, reading not taken)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4581	8.12.1.56.3.8 GLOBAL NAVIGATION PANEL		Heading	CloudDRxRTR
PAGUI4582	See Shared Elements: Panels.		Info	CloudDRxRTR
PAGUI42	8.12.2 Components		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1535	8.12.2.1 ABBOTT LOGO COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1536	Component: ABBOTT LOGO COMPONENT State: N/A Display: Always displayed in the STARTUP SCREEN and GLOBAL HEADING PANEL. Label: N/A Behavior: Icon.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1538	<Picture> ABBOTT LOGO COMPONENT (example: STARTUP SCREEN) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1537	<Picture> ABBOTT LOGO COMPONENT (example: GLOBAL HEADING PANEL) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1058	8.12.2.2 ALPHABETICAL KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1059	Component: ALPHABETICAL KEYBOARD COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS EMAIL ID PANEL, USER DETAILS NEW PASSWORD PANEL, EDIT ACCOUNT INFO PANEL, and PRIMARY LAB PANEL. The ALPHABETICAL KEYBOARD COMPONENT is displayed after the user selects an input component that requires alphabetical input. The ALPHABETICAL KEYBOARD COMPONENT is controlled by the local OS of the user's device and is only utilized by the application. Label: N/A - determined by user device Operating System Behavior: Icon and trigger. Upon activation, the ALPHABETICAL KEYBOARD COMPONENT relays the user input to the input component display that was originally selected.		UI Info	CloudDRx 1.0 CloudDRxRTR
PAGUI1060	<Picture> ALPHABETICAL KEYBOARD COMPONENT (example)		UI Design	CloudDRx 1.0

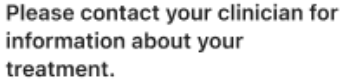
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxRTR CloudDRxPS
PAGUI1061	<Picture> ALPHABETICAL KEYBOARD COMPONENT (example: Android OS)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI850	8.12.2.3 BACK ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4052	Component: BACK ACTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the POLICY SCREEN, USER AUTHENTICATION SCREEN, SHORT PRIVACY POLICY PANEL, END USER LICENSE PANEL, DATA USE PANEL, FULL PRIVACY POLICY PANEL, PATIENT IDENTIFICATION PANEL, PROFILE PANEL, SETTINGS PANEL, EDUCATION PANEL, ABOUT HEART FAILURE PANEL, ABOUT HEART FAILURE SYSTEM PANEL, LEARNING PANEL, IN-APP TRAINING PANEL, SUPPORT PANEL, FAQ PANEL, CONTACT PANEL, ABOUT PANEL, MEDICATION UPDATE REMEDIATION PANEL, USER MANUAL WEB VIEW, and MEDICATION UPDATE COMPLETE PANEL Conditionally displayed in the CONNECT HEALTH APP PANEL; not displayed when the panel is in the REGISTRATION SCREEN. Label: N/A Behavior: Icon and trigger. Upon activation, the BACK ACTION TRIGGER COMPONENT displays the panel or screen that was displayed before the current panel was displayed. Upon activation in the PATIENT IDENTIFICATION PANEL or USER AUTHENTICATION SCREEN, the BACK ACTION TRIGGER COMPONENT displays the STARTUP SCREEN.		UI Info	CloudDRxRTR CloudDRxPS
PAGUI852	<Picture> BACK ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI324	8.12.2.6 CALL CLINIC ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4027	Component: CALL CLINIC ACTION TRIGGER COMPONENT State: Blue/ White Display: Conditionally displayed in the ELEVATION CHANGE MESSAGE CARD PANEL when the panel is in the active state. The CALL CLINIC ACTION TRIGGER COMPONENT is displayed in the white state.		UI Info	CloudDRxRTR CloudDRxPS



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Conditionally displayed in the CALL CLINIC DIRECTCALL MESSAGE CARD PANEL when the panel is in the active state. The CALL CLINIC ACTION TRIGGER COMPONENT is displayed in the blue state. Conditionally displayed in the DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL when the panel is in the active state. The CALL CLINIC ACTION TRIGGER COMPONENT is displayed in the blue state. Always displayed in the MEDICATION UPDATE REMEDIATION PANEL. The CALL CLINIC ACTION TRIGGER COMPONENT is displayed in the white state. Always displayed in the CONFIRM CURRENT MEDICATIONS PANEL. The CALL CLINIC ACTION TRIGGER COMPONENT is displayed in the white state. Label: <label.gb.call_clinic> Behavior: Text and trigger. Upon activation, the CALL CLINIC ACTION TRIGGER COMPONENT begins calling the number associated with the registered clinic for the user. The CALL CLINIC ACTION TRIGGER COMPONENT changes state based on where it is displayed.			
PAGUI2216	<Picture> CALL CLINIC ACTION TRIGGER COMPONENT (Blue) 	PString73 Call Clinic callClinic label.gb.call_clinic	UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2217	<Picture> CALL CLINIC ACTION TRIGGER COMPONENT (White) 	PString73 Call Clinic callClinic label.gb.call_clinic	UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4382	8.12.2.7 CALL CLINIC CONTENT LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4383	Component: CALL CLINIC CONTENT LABEL COMPONENT State: N/A Display: Always displayed in the CALL CLINIC DIRECTCALL MESSAGE CARD PANEL and DRX CALL CLINIC DIRECT CALL MESSAGE CARD PANEL. Label: <notif.manual.contact_clinic> Behavior: Text.	PString743 Please contact your clinician for information about your treatment. contactClinicBody notif.manual.contact_clinic	UI Info	CloudDRxRTR CloudDRxPS
PAGUI4384	<Picture> CALL CLINIC CONTENT LABEL COMPONENT 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2615	8.12.2.8 CANCEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2616	Component: CANCEL ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the USER DETAILS EMAIL ID PANEL, USER DETAILS CONFIRMATION PANEL, USER DETAILS NEW PASSWORD PANEL, CONFIRM MFA METHOD PANEL, SEND MFA CODE PANEL, and ENTER MFA CODE PANEL. Label: N/A - determined by a third party vendor Behavior: Text and trigger. Upon actication, returns to previous screen or panel.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Upon activation in the USER DETAILS EMAIL ID PANEL, the CANCEL ACTION TRIGGER COMPONENT displays the FORGOT PASSWORD PANEL. Upon activation in the TWO FACTOR AUTHENTICATION PANEL, the CANCEL ACTION TRIGGER COMPONENT displays the USER AUTHENTICATION SCREEN. Upon activation in the USER DETAILS NEW PASSWORD PANEL, the CANCEL ACTION TRIGGER COMPONENT displays the SETTINGS PANEL.			
PAGUI2617	<Picture> CANCEL ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI261	8.12.2.9 CLOSE PANEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4051	Component: CLOSE PANEL ACTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the LOGIN PANEL, PAP READING INSTRUCTIONS PANEL, RECORD BP READING PANEL, RECORD WEIGHT PANEL, CONFIRM INPUT PANEL, EDIT ACCOUNT INFO PANEL, SHORT PRIVACY POLICY PANEL, END USER LICENSE AGREEMENT PANEL, FULL PRIVACY POLICY PANEL, THIRD PARTY LICENSE PANEL, PATIENT IDENTIFICATION PANEL, PATIENT INFORMATION PANEL, PRIMARY LAB PANEL, CONNECT HEALTH APP PANEL, NOTIFICATIONS SCREEN, CONFIRM CURRENT MEDICATIONS PANEL, MEDICATION UPDATE REMEDIATION PANEL, and MEDICATION UPDATE COMPLETE PANEL. Label: N/A Behavior: Icon and trigger. Upon activation, the CLOSE PANEL ACTION TRIGGER COMPONENT displays the panel or screen that was previously displayed. Upon activation in the PATIENT IDENTIFICATION PANEL or the PATIENT INFORMATION PANEL, the CLOSE PANEL ACTION TRIGGER COMPONENT displays the STARTUP SCREEN. Upon activation in the SHORT PRIVACY POLICY PANEL, END USER LICENSE AGREEMENT PANEL, or DATA USE PANEL before a successful login or registration, the CLOSE PANEL ACTION TRIGGER COMPONENT displays the POLICY SCREEN. Upon activation in the CONFIRM CURRENT MEDICATIONS PANEL, the MEDICATION UPDATE REMEDIATION PANEL, and the MEDICATION UPDATE COMPLETE PANEL, the CLOSE PANEL ACTION TRIGGER COMPONENT displays the previous high level screen (HOME SCREEN or NOTIFICATIONS SCREEN).		UI Info	CloudDRxRTR CloudDRxPS
PAGUI720	<Picture> CLOSE PANEL ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2968	8.12.2.10 CONFIRM ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

Page 390 of 422



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	COMPONENT displays the SUCCESSFUL PASSWORD CHANGE PANEL. Upon activation in the TWO FACTOR AUTHENTICATION PANEL, the CONTINUE ACTION TRIGGER COMPONENT runs any required verification displays the HOME SCREEN.			
PAGUI2670	<Picture> CONTINUE ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI104	8.12.2.13 DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI105	Component: DATE LABEL COMPONENT State: Day/ Week/ Month Display: Always displayed in READINGS SCREEN and the WEIGHT SCREEN. Label: N/A - determined by third party vendor package Behavior: Text. The DATE LABEL COMPONENT labels the dates for the displayed data.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1614	<Picture> DATE LABEL COMPONENT (example: day) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1615	<Picture> DATE LABEL COMPONENT (example: week) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1616	<Picture> DATE LABEL COMPONENT (example: month) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI52	8.12.2.14 DATE SELECTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI53	Component: DATE SELECTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the READINGS SCREEN and the WEIGHT SCREEN. The DATE SELECTION TRIGGER COMPONENT is conditionally displayed when the INTERVAL TAB COMPONENT is in the week or month state. Label: N/A Behavior: Text and Trigger. Upon activation in the enabled state, the DATE SELECTION TRIGGER COMPONENT displays the HISTORICAL DATE SELECTION CALENDAR COMPONENT. The DATE SELECTION TRIGGER COMPONENT is disabled when the INTERVAL TAB COMPONENT displays an interval of 'Day'.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1617	<Picture> DATE SELECTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1880	8.12.2.15 DONE ACTION TRIGGER COMPONENT		Heading	CloudDRxRTR
PAGUI4050	Component: DONE ACTION TRIGGER COMPONENT State: blue/white Display: Always displayed in the THIRD PARTY LICENSE PANEL, RECORD BP READING PANEL, and, FULL PRIVACY POLICY PANEL. Conditionally displayed in the END USER LICENSE AGREEMENT PANEL in blue state. Always displayed in	PString204 Done done label.gb.done	UI Info	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	MEDICATION UPDATE COMPLETE PANEL in white state. Conditionally displayed in the END USER LICENSE AGREEMENT PANEL when opened from the ABOUT PANEL. Label: <label.gb.done> Behavior: Text and Trigger. Upon activation, the DONE ACTION TRIGGER COMPONENT displays the screen or panel that was previously displayed. Upon activation in the RECORD BP PANEL, the DONE ACTION TRIGGER COMPONENT triggers the required verification and displays the CONFIRM INPUT PANEL. Upon activation in the MEDICATION UPDATE COMPLETE PANEL, the DONE ACTION TRIGGER COMPONENT returns the user to the previously displayed high-level screen (i.e. HOME or NOTIFICATIONS).			
PAGUI1882	<Picture> DONE ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4178	8.12.2.16 DRX TREND MESSAGE LABEL COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4179	Component: DRX TREND MESSAGE LABEL COMPONENT State: N/A Display: Always displayed in the DRX MEDICATION MESSAGE CARD PANEL. Label: <notif.dtp.1_high>, <notif.dtp.many_high>, <notif.dtp.1_low>, <notif.dtp.many_low> Behavior: Text.	PString567 You have {26} {29} pressure readings in a row. html.gb.out_of_range_noInfo_bold PString734 We've noticed {26} high pressure reading. highPressureReadings notif.dtp.1_high PString735 We've noticed {27} low pressure reading. lowPressureReadings notif.dtp.1_low PString736 We've noticed {26} high pressure readings in a row. highPressureReadingsInARow notif.dtp.many_high PString737 We've noticed {27} low pressure readings in a row. lowPressureReadingsInARow	UI Info	CloudDRxRTR CloudDRxPS

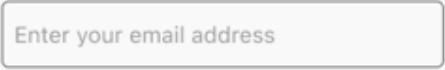
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		notif.dtp.many_low PString1003 You had {26} high CardioMEMS™ reading. highPressureReadings notif.dtp.1_high PString1004 You had {26} high CardioMEMS™ readings in a row. highPressureReadingsInARow notif.dtp.many_high PString1005 You had {27} low CardioMEMS™ reading. lowPressureReadings notif.dtp.1_low PString1006 You had {27} low CardioMEMS™ readings in a row. lowPressureReadingsInARow notif.dtp.many_low		
PAGUI4180	<Picture> DRX TREND MESSAGE LABEL COMPONENT TBD		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2965	8.12.2.17 EDIT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2966	Component: EDIT ACTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the CONFIRM INPUT PANEL. Label: <label.gb.edit> Behavior: Text and trigger. Upon activation, the EDIT ACTION TRIGGER COMPONENT displays the previously displayed input panel.	PString613 Edit edit label.gb.edit	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2967	<Picture> EDIT ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1013	8.12.2.18 EMAIL LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



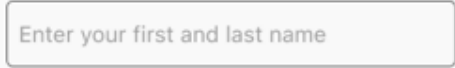
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI1014	Component: EMAIL LABEL COMPONENT State: N/A Display: Always displayed in PROFILE PANEL and EDIT ACCOUNT INFO PANEL. Label: <label.gb.email_address> Behavior: Text.	PString32 Email Address email label.gb.email_address	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1015	<Picture> EMAIL LABEL COMPONENT Email Address		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1653	8.12.2.19 EMAIL INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1654	Component: EMAIL INPUT COMPONENT State: No Input/ Input/ View Input Display: Always displayed in the PATIENT INFORMATION PANEL. The default status is no input. Always displayed in the EDIT ACCOUNT INFO PANEL in the View Input state. Label: <label.gb.enter_email_address> Behavior: Text and input box. - In the no input state, display the default text. - In the input state, display the input given from the ALPHABETICAL KEYBOARD COMPONENT or the currently saved email associated with the user. - In the View Input state, display the currently saved e-mail associated with the user.	PString642 Enter your email address emailFieldPlaceholder label.gb.enter_email_address	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1655	<Picture> EMAIL INPUT COMPONENT (default: no input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1656	<Picture> EMAIL INPUT COMPONENT (example: input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1658	<Picture> EMAIL INPUT COMPONENT (example: edit input) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2750	8.12.2.20 EMAIL ID INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2751	Component: EMAIL ID INPUT COMPONENT State: No Input/ Input/ Display Display: Always displayed in the USER DETAILS EMAIL ID PANEL and USER DETAILS CONFIRMATION PANEL. The default state is no input. The EMAIL ID INPUT COMPONENT is always in the display state when in the USER DETAILS CONFIRMATION PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

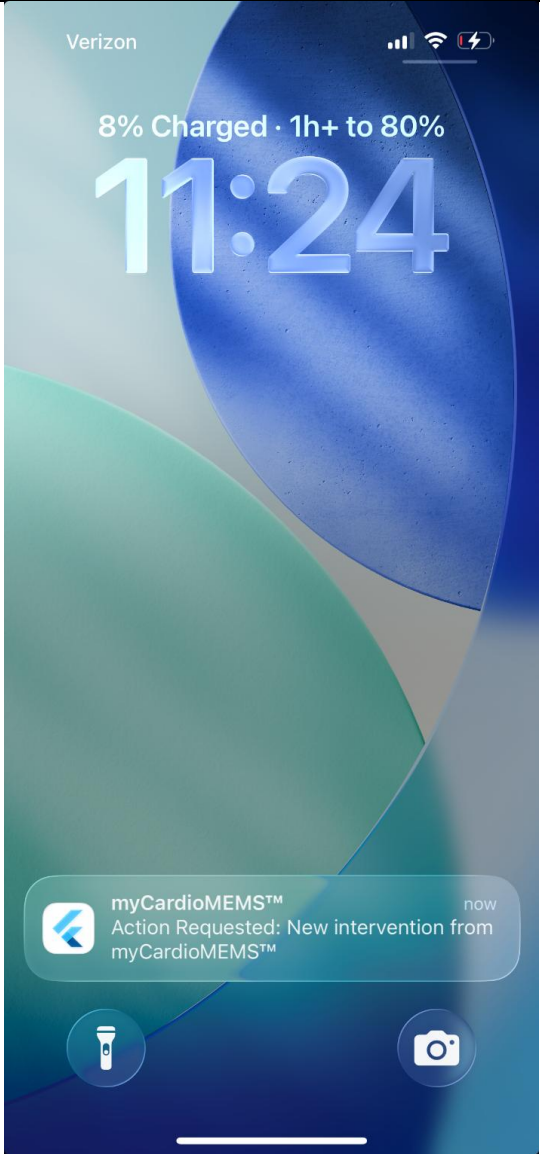
ID		English US Strings	Object Type	Project
	Label: N/A - determined by a third party vendor Behavior: Text and input display. In the no input state, the EMAIL ID INPUT COMPONENT displays label. In the input state, the EMAIL ID INPUT COMPONENT displays the input from the NUMERIC KEYBOARD COMPONENT and a cursor.			
PAGUI2752	<Picture> EMAIL ID INPUT COMPONENT (default: No Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2753	<Picture> EMAIL ID INPUT COMPONENT (Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2754	<Picture> EMAIL ID INPUT COMPONENT (Display, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1634	8.12.2.21 FULL NAME LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1635	Component: FULL NAME LABEL COMPONENT State: N/A Display: Always displayed in the PATIENT INFORMATION PANEL, PROFILE PANEL, and EDIT ACCOUNT INFO PANEL. Label: <label.gb.full_name> Behavior: Text.	PString643 Full Name fullName label.gb.full_name	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1636	<Picture> FULL NAME LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1637	8.12.2.22 FULL NAME INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1638	Component: FULL NAME INPUT COMPONENT State: Enabled/ Disabled Display: Always displayed in the PATIENT INFORMATION PANEL and EDIT ACCOUNT INFO PANEL. Label: <label.gb.enter_first_last> Behavior: Text and input display. When disabled, the FULL NAME INPUT COMPONENT displays label.When enabled, the FULL NAME INPUT COMPONENT displays the input from the ALPHABETIC KEYBOARD COMPONENT and a cursor. When the user deselects FULL NAME INPUT COMPONENT, the input is validated.	PString615 Enter your first and last name nameValidation label.gb.enter_first_last	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1639	<Picture> FULL NAME INPUT COMPONENT (default: disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1640	<Picture> FULL NAME INPUT COMPONENT (example: enabled)		UI Design	CloudDRx 1.0 CloudDRxRTR

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

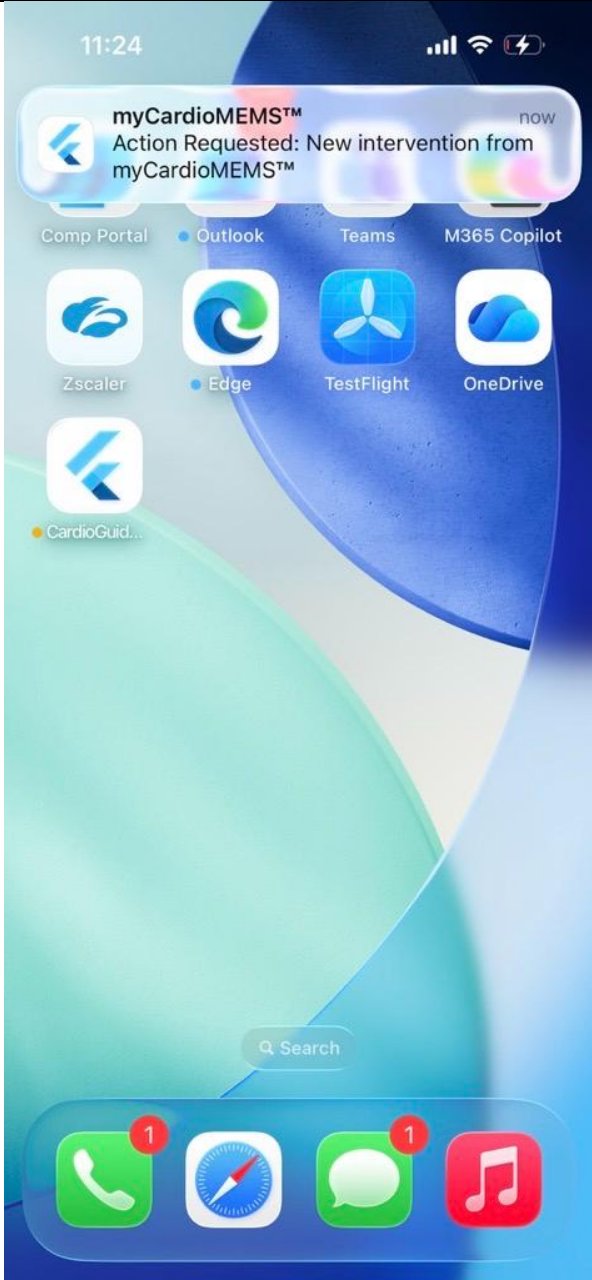
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Alexander Stallone			CloudDRxPS
PAGUI2448	8.12.2.23 HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2449	Component: HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT State: Lock Screen/ Banner and Apple/ Android Display: Conditionally displayed in the user's device lock screen when there is at least one new high priority message card. Conditionally displayed as a drop down banner in the user's screen when there is at least one new high priority message card. Label: <msg.gb.high_priority_notification> Behavior: Text. Depending on the user's device operating system, the HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT may allow the user to open the Patient Application. In the lock screen state, the HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT is displayed on the device lock screen until the user unlocks their device at least once. In the banner state, the HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT is displayed at the top of the device screen. The HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT changes from Apple to Android based on the user's device operating system.	PString644 {0} {n} Action Requested: New intervention from {0} msg.gb.high_priority_notifi cation	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4044	<Picture> HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT (default: Apple, Lock Screen)		UI Design	CloudDRxRTR CloudDRxPS

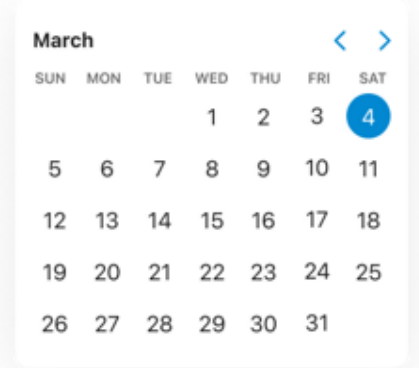
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4046	<Picture> HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT (default: Apple, Banner)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				
PAGUI4047	<Picture> HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT (Android, Lock Screen, light mode) TBD		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4048	<Picture> HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT (Android, Lock Screen, dark mode) TBD		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4049	<Picture> HIGH PRIORITY OUT OF APP NOTIFICATION COMPONENT (Android, Banner) TBD		UI Design	CloudDRxRTR CloudDRxPS
PAGUI106	8.12.2.24 HISTORICAL DATE SELECTION CALENDAR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI119	Component: HISTORICAL DATE SELECTION CALENDAR COMPONENT State: N/A Display: Conditionally displayed in the READINGS SCREEN and the WEIGHT SCREEN when the DATE SELECTION TRIGGER COMPONENT is activated. Displays a calendar object that can be used to specify the date range of readings to be viewed. USERS can select up to 180 days previous to the current date. If the INTERVAL TAB COMPONENT is set to 'Week' then the earliest 6 days of the 180 previous days is not selectable. If the INTERVAL TAB COMPONENT is set to 'Month' then the earliest 29 days of the 180 previous days is not selectable. Non-selectable dates are greyed out. Label: N/A - determined by third party vendor package Behavior: Calendar and Trigger. USER can select an available date from the calendar. Upon selection, the DATE LABEL COMPONENT will be updated to the newly selected date range and all available readings data for the specified date range will be displayed. USER can click anywhere outside of the calendar to close the calendar.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1618	<Picture> HISTORICAL DATE SELECTION CALENDAR COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4108	8.12.2.25 INFORMATION REQUIRED ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4109	Component: INFORMATION REQUIRED ERROR COMPONENT State: N/A Display: Conditionally displayed above the accosiated input component in the USER DETAILS EMAIL ID PANEL when the user has selected the SEND VERIFICATION CODE ACTION TRIGGER COMPONENT before completing all required input. Conditionally displayed above the associated input component in the USER DETAILS NEW PASSWORD PANEL when the user has selected the CONTINUE ACTION TRIGGER COMPONENT before completing all required input.	PString932 This information is required. ots.config.information_required	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Label: N/A - determined by third party vendor - <ots.config.information_required> Behavior: Text.			
PAGUI4110	<Picture> REQUIRED FIELD ERROR COMPONENT (example) This information is required.		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI49	8.12.2.26 INTERVAL TAB COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI50	Component: INTERVAL TAB COMPONENT State: Day/ Week/ Month Display: Always displayed in the READINGS SCREEN and WEIGHT PANEL. The INTERVAL TAB COMPONENT displays an interval of 'Week' as default in the READINGS SCREEN and displays an interval of 'Day' as a default in the WEIGHT PANEL. Label: label.gb.interval_tab Behavior: Trigger. Only one interval value may be selected at a time. Upon selection, the INTERVAL TAB COMPONENT displays data for the selected interval of time. Upon swiping between interval tabs, a refresh icon would be displayed prior to displaying data for the selected interval of time.	PString456 Day day label.gb.day PString691 Week week label.gb.week PString692 Month month label.gb.month	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1611	<Picture> INTERVAL TAB COMPONENT Day Week Month		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1612	<Picture> INTERVAL TAB COMPONENT (READINGS SCREEN default: week) Day Week Month		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1613	<Picture> INTERVAL TAB COMPONENT (example: month) Day Week Month		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2454	8.12.2.27 LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2455	Component: LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT State: Apple/ Android Display: Conditionally displayed in the user's device lock screen when there is at least one new low priority message card. Conditionally displayed as a drop down banner in the user's screen when there is at least one new low priority message card. Label: <msg.gb.low_priority_notification>	PString645 New notification from {0} msg.gb.low_priority_notification	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Behavior: Text. Depending on the user's device operating system, the LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT may allow the user to open the Patient Application. The LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT is displayed on the device lock screen until the user unlocks their device at least once. If the device is unlocked, the LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT is displayed at the top of the device screen. The LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT changes from Apple to Android based on the user's device operating system.			
PAGUI4040	<Picture> LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT (default: Apple, Lock Screen)		UI Design	CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

					
--	--	--	--	--	--

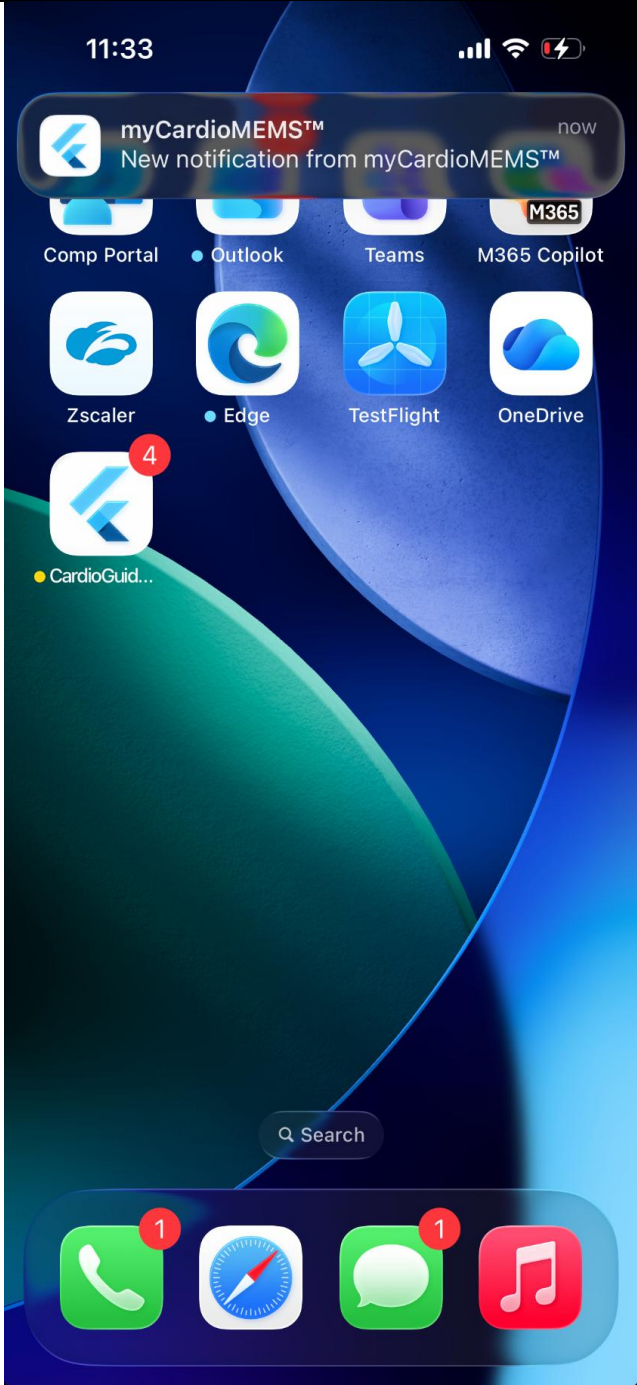
This sample user interface is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4850	<Picture> LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT (default: Apple, banner)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

				
PAGUI4042	<Picture> LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT (default: Android, Lock Screen, light mode)		UI Design	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	TBD			CloudDRxPS
PAGUI4043	<Picture> LOW PRIORITY OUT OF APP NOTIFICATION COMPONENT (default: Apple, Lock Screen, dark mode) TBD		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2212	8.12.2.28 MARK AS READ ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2213	Component: MARK AS READ ACTION TRIGGER COMPONENT State: Blue/ White Display: Conditionally displayed in ELEVATION CHANGE MESSAGE CARD PANEL, DRX FREEFORM MESSAGE CARD PANEL, ELEVATED PRESSURE MESSAGE CARD PANEL, LAB REQUEST MESSAGE CARD PANEL, FREEFORM MESSAGE CARD PANEL, READINGS REVIEWED MESSAGE CARD PANEL, PRESSURE READING MESSAGE CARD PANEL, and READING REMINDER MESSAGE CARD PANEL when the panel is in the active state. The MARK AS READ ACTION TRIGGER COMPONENT is displayed in the blue state. Conditionally displayed in the CALL CLINIC DIRECTCALL MESSAGE CARD PANEL and DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL, when the panel is in the active state. The MARK AS READ ACTION TRIGGER COMPONENT is displayed in the white state. Label: <label.gb.mark_as_read> Behavior: Text and trigger. Upon activation, the MARK AS READ ACTION TRIGGER COMPONENT changes the state of the associated message card to complete.	PString123 Mark as Read markAsRead label.gb.mark_as_read	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2214	<Picture> MARK AS READ ACTION TRIGGER COMPONENT (Blue) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2328	<Picture> MARK AS READ ACTION TRIGGER COMPONENT (White) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI498	8.12.2.29 MEDICATION LABEL ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI499	Component: MEDICATION LABEL ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in MEDICATIONS LIST PANEL and REMINDER GROUP LIST PANEL. Label: N/A - determined by clinician input in the portal Behavior: Text and trigger. Upon activation, the MEDICATION LABEL ACTION TRIGGER COMPONENT displays the ASSIGN REMINDER GROUP PANEL.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI500	<Picture> MEDICATION LABEL ACTION TRIGGER COPONENT (example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI389	8.12.2.30 MEDICATION LOG RADIO BUTTON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4039	Component: MEDICATION LOG RADIO BUTTON COMPONENT State: Selected/ Not Selected Display: Conditionally displayed in the MEDICATIONS LOG PANEL if medications log is enabled in the MEDICATIONS SETTINGS PANEL. Always displayed in the RECORD MEDICATIONS PANEL. Multiple MEDICATION LOG RADIO BUTTON COMPONENTS may be displayed simultaneously. When displayed in the RECORD MEDICATIONS PANEL, multiple MEDICATION LOG RADIO BUTTON COMPONENT are organized by reminder group. The default state is not selected. In the not selected state, display an empty check box. In the selected state, disply a filled check box. Label: <data.Code.Med_Frequency_Cd.Code_Desc.QD>, <data.Code.Med_Frequency_Cd.Code_Desc.BID>, <data.Code.Med_Frequency_Cd.Code_Desc.TID>, <data.Code.Med_Frequency_Cd.Code_Desc.QID>, <data.Code.Med_Frequency_Cd.Code_Desc.QHS>, <data.Code.Med_Frequency_Cd.Code_Desc.QAM>, <data.Code.Med_Frequency_Cd.Code_Desc.QPM>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.TABLET>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.CAPSULE>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.IV>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.INJECTION>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.PATCH>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.OTHER>, <data.Code.Dosage_Form_Unit_Cd.Code_Desc.SPRAY>, <data.Code.Med_Unit_Cd.Code_Desc.MG>, <data.Code.Med_Unit_Cd.Code_Desc.MCG>, <data.Code.Med_Unit_Cd.Code_Desc.MGKG>, <data.Code.Med_Unit_Cd.Code_Desc.MCGKGMIN>, <data.Code.Med_Unit_Cd.Code_Desc.ML>, <data.Code.Med_Unit_Cd.Code_Desc.IU>, <data.Code.Med_Unit_Cd.Code_Desc.IUKG>, <data.Code.Med_Unit_Cd.Code_Desc.MEQ>, <data.Code.Med_Unit_Cd.Code_Desc.MGML>, <data.Code.Med_Unit_Cd.Code_Desc.U>, <data.Code.Med_Unit_Cd.Code_Desc.MLMIN>, <data.Code.Med_Unit_Cd.Code_Desc.NGKGMIN>, <data.Code.Med_Unit_Cd.Code_Desc.MGKGMIN> Behavior: Text and trigger. Upon activation in the not selected state, the MEDICATION LOG RADIO BUTTON COMPONENT changes the radio button display and saves the selected medication as taken in the associated reminder group. Upon activation in the selected state, the MEDICATION LOG RADIO BUTTON COMPONENT changes the radio button display and saves the medication as not taken in the associated reminder group.	PString140 Once daily qd data.Code.Med_Frequency_Cd.Code_Desc.QD PString141 Twice daily bid data.Code.Med_Frequency_Cd.Code_Desc.BID PString142 3 times daily tid data.Code.Med_Frequency_Cd.Code_Desc.TID PString143 4 times daily qid data.Code.Med_Frequency_Cd.Code_Desc.QID PString154 Tablet tablet data.Code.Dosage_Form_Unit_Cd.Code_Desc.TABLET PString155 Capsule capsule data.Code.Dosage_Form_Unit_Cd.Code_Desc.CAPSULE PString156 IV iv data.Code.Dosage_Form_Unit_Cd.Code_Desc.IV PString157 Injection injection data.Code.Dosage_Form_Unit_Cd.Code_Desc.INJECTION N	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		PString158 Patch patch data.Code.Dosage_Form_Unit_Cd.Code_Desc.PATCH PString159 Other other data.Code.Dosage_Form_Unit_Cd.Code_Desc.OTHER PString161 Spray spray data.Code.Dosage_Form_Unit_Cd.Code_Desc.SPRAY		
PAGUI391	<Picture> MEDICATION LOG RADIO BUTTON COMPONENT (default: not selected, default: standard) Rivaroxaban 1mg Zebeta Tablet 2 times daily ☐		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI392	<Picture> MEDICATION LOG RADIO BUTTON COMPONENT (selected, default: standard) Rivaroxaban 1mg Zebeta Tablet 2 times daily ☒		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3336	8.12.2.32 MFA PHONE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3337	Component: MFA PHONE LABEL COMPONENT State: N/A Display: Always displayed in the SEND MFA CODE PANEL and the ENTER CODE PANEL. Label: N/A - determined by a third party vendor Behavior: Text. The text changes based on the last three digits of the user's phone number.	PString973 We have the following number on record for you. We can send a code via SMS or phone to authenticate you. ots.authentication.phone_number_message	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI3338	<Picture> MFA PHONE LABEL COMPONENT We have the following number on record for you. We can send a code via SMS or phone to authenticate you. ***-***-456		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1361	8.12.2.33 NEXT ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1362	Component: NEXT ACTION TRIGGER COMPONENT State: Disabled/ Enabled Display: Always displayed in the PATIENT IDENTIFICATION PANEL Label: <label.gb.next> Behavior: Text and trigger. Upon activation, in the disabled state, the NEXT ACTION TRIGGER COMPONENT does nothing. The NEXT ACTION TRIGGER COMPONENT is enabled after all required input for the displayed panel is complete. Upon activation, in the enabled state, the NEXT ACTION TRIGGER COMPONENT verifies the changes made to the displayed panel, saves the verified changes, and displays the next panel in the associated flow.	PString42 Next next label.gb.next	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1363	<Picture> NEXT ACTION TRIGGER COMPONENT (default: disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1364	<Picture> NEXT ACTION TRIGGER COMPONENT (enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2133	8.12.2.34 NOTIFICATION BANNER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4005	Component: NOTIFICATION BANNER COMPONENT State: Active/ Complete/ Out of Date Display: Always displayed in the ELEVATION CHANGE MESSAGE CARD PANEL. The default state is active. In the active state, the NOTIFICATION BANNER COMPONENT displays the active icon and the active background color. In the complete state, the NOTIFICATION BANNER COMPONENT displays the complete icon and complete background color. In the out of date state, the NOTIFICATION BANNER COMPONENT displays the active icon and the complete background color. Label: <label.gb.notification>; Text is all capitalized. Behavior: Icon and text. The NOTIFICATION BANNER COMPONENT changes state with the message card panel it is an element of.	PString646 Notification default label.gb.notification	UI Info	CloudDRxRTR CloudDRxPS
PAGUI2135	<Picture> NOTIFICATION BANNER COMPONENT (default: Active) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2136	<Picture> NOTIFICATION BANNER COMPONENT (Complete) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI2137	<Picture> NOTIFICATION BANNER COMPONENT (Out of Date) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI747	8.12.2.35 NOTIFICATION DATE LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI748	Component: NOTIFICATION DATE LABEL COMPONENT State: N/A Display: Always displayed in the DRX CALL CLINIC DIRECTCALL MESSAGE CARD PANEL, DRX FREEFORM MESSAGE CARD PANEL, MEDICATION MESSAGE CARD PANEL, DRX MEDICATION MESSAGE CARD PANEL, ELEVATED PRESSURE MESSAGE CARD PANEL, ELEVATION CHANGE MESSAGE CARD PANEL, CALL CLINIC DIRECTCALL MESSAGE CARD PANEL, LAB REQUEST MESSAGE CARD PANEL, FREEFORM MESSAGE CARD PANEL, READINGS REVIEWED MESSAGE CARD PANEL, PRESSURE READING MESSAGE CARD PANEL, READING REMINDER MESSAGE CARD PANEL, and MEDICATION CHANGE NOTIFICATION CARD PANEL. Label: N/A - determined by a third party vendor package Behavior: Text. The NOTIFICATION DATE LABEL COMPONENT displays the date the notification was sent to the applicaion.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI749	<Picture> NOTIFICATION DATE LABEL COMPONENT (example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1062	8.12.2.36 NUMERIC KEYBOARD COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1063	Component: NUMERIC KEYBOARD COMPONENT State: N/A Display: Conditionally displayed in the RECORD BP READING PANEL, RECORD WEIGHT PANEL, PRIMARY LAB PANEL, and EDIT ACCOUNT INFO PANEL. The NUMERIC KEYBOARD COMPONENT is displayed after the user selects an input component that requires numeric input. The NUMERIC KEYBOARD COMPONENT is controlled by the local OS of the user's device and is only utilized by the application. Label: N/A - determined by user device OS Behavior: Icon and trigger. Upon activation, the NUMERIC KEYBOARD COMPONENT relays the user input to the input component display that was originally selected.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1064	<Picture> NUMERIC KEYBOARD COMPONENT (example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI1065	<Picture> NUMERIC KEYBOARD COMPONENT (example: Android OS)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1442	8.12.2.37 OKAY ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1443	Component: OKAY ACTION TRIGGER COMPONENT State: N/A Display: Always displayed in the MFA MAXIMUM ATTEMPTS PANEL, NO INTERNET CONNECTION PANEL, NO BACKEND CONNECTION PANEL, USER SAME DEVICE ERROR PANEL, INVALID USER PANEL, FORGOT USER ID PANEL, ELEVATION NOTIFICATION DIALOG PANEL, MEDICATION NOTIFICATION DIALOG PANEL, REGISTRATION FAILED PANEL, PA PRESSURE ENABLED PANEL, PA PRESSURE DISABLED PANEL, MEDICATIONS ENABLED PANEL, MEDICATIONS DISABLED PANEL, UPDATE PERMISSIONS INSTRUCTIONS PANEL, BLOOD PRESSURE ERROR PANEL, SUCCESSFUL SYNC PANEL, SUCCESSFUL PASSWORD CHANGE PANEL, SESSION EXPIRED PANEL, EMAIL IN USE POPUP ERROR PANEL, and PERMISSIONS REVOKED PANEL. Label: N/A - determined by third party vendor Behavior: Text and trigger. Upon activation, the OKAY ACTION TRIGGER COMPONENT closes the panel that is being displayed and displays the panel that was previously displayed. Upon activation in the SUCCESSFUL PASSWORD CHANGE PANEL, the OKAY ACTION TRIGGER COMPONENT logs the user out of the Patient Application and displays the STARTUP SCREEN.	PString261 Ok ok label.gb.ok	UI Info	CloudDRx 1.0 CloudDRxRTR
PAGUI1444	<Picture> OKAY ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4476	8.12.2.38 OPTIMAL RANGE LABEL COMPONENT		Heading	CloudDRxRTR
PAGUI4328	Component: OPTIMAL RANGE LABEL COMPONENT State: Enabled, Disabled Display: Conditionally displayed in the CMEMS PAP TREND DISPLAY PANEL and CMEMS PAP READING CARD PANEL when PA Visibility is enabled. Label: <label.gb.optimal_range>, <label.gb.optimal_range_parameters>, <label.gb.mmhg>	PString452 mmHg mmHg label.gb.mmhg PString455	UI Info	CloudDRxRTR CloudDRxPS CloudDRxLMR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Behavior: Text. This component is displayed when enabled. This OPTIMAL RANGE LABEL COMPONENT is enabled when PA Visibility is enabled. This component is not displayed when disabled. This OPTIMAL RANGE LABEL COMPONENT is not displayed when PA Visibility is disabled.	Optimal Range: optimalRange label.gb.optimal_range PString1038 {17}-{18} label.gb.optimal_range_parameters		
PAGUI4418	<Picture> OPTIMAL RANGE LABEL COMPONENT Optimal range: 18-24mmHg		UI Design	CloudDRxRTR CloudDRxPS
PAGUI1690	8.12.2.39 PASSWORD LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1691	Component: PASSWORD LABEL COMPONENT State: N/A Display: Always displayed in the PATIENT INFORMATION PANEL Label: <label.gb.password> Behavior: Text.	PString33 Password passwordFieldLabel label.gb.password	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1692	<Picture> PASSWORD LABEL COMPONENT Password		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1700	8.12.2.40 PASSWORD REQUIREMENT LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1701	Component: PASSWORD REQUIREMENT LABEL COMPONENT State: Met/ Not Met Display: Always displayed in the PATIENT INFORMATION PANEL. The default state is not met. Label: <label.gb.password_characters>, <label.gb.password_cap_characters>, <label.gb.password_lower_characters>, <label.gb.password_special_characters> Behavior: Icon and text. The checkmark icons to the left of each requirement listed changes color as the specific requirement changes state.	PString471 At least 8 characters rule1 label.gb.password_characters PString472 One upper case letter rule2 label.gb.password_cap_characters PString473 One lower case letter rule3 label.gb.password_lower_characters PString474	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

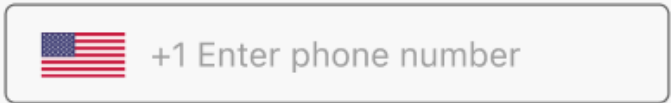
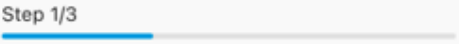
This confidential document is the property of Abbott and shall not be reproduced, distributed, disclosed or used without the express written consent of Abbott.



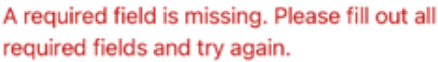
myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
		One number or symbol rule4 label.gb.password_special_characters		
PAGUI1702	<Picture> PASSWORD REQUIREMENT LABEL COMPONENT (default: not met) At least 8 characters One upper case letter One lower case letter One number or symbol		UI Design	CloudRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1703	<Picture> PASSWORD REQUIREMENT LABEL COMPONENT (example: partially met) At least 8 characters One upper case letter One lower case letter One number or symbol		UI Design	CloudRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1704	<Picture> PASSWORD REQUIREMENT LABEL COMPONENT (example: met) At least 8 characters One upper case letter One lower case letter One number or symbol		UI Design	CloudRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1016	8.12.2.41 PHONE LABEL COMPONENT		Heading	CloudRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1017	Component: PHONE LABEL COMPONENT State: N/A Display: Always displayed in PATIENT INFORMATION PANEL, PROFILE PANEL and EDIT ACCOUNT INFO PANEL. Label: <label.gb.phone_number> Behavior: Text.	PString215 Phone Number phoneNumber label.gb.phone_number	UI Info	CloudRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1018	<Picture> PHONE LABEL COMPONENT (display) Phone Number		UI Design	CloudRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1645	8.12.2.42 PHONE INPUT COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI1646	Component: PHONE INPUT COMPONENT State: Enabled/ Disabled Display: Always displayed in the PATIENT INFORMATION PANEL and EDIT ACCOUNT INFO PANEL. Label: <label.gb.enter_phone_number>	PString633 Enter phone number enterPhoneNumber label.gb.enter_phone_number	UI Info	CloudRx 1.0 CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Behavior: Text and input display. When disabled, the PHONE INPUT COMPONENT displays label.When enabled, the PHONE INPUT COMPONENT displays the input from the NUMERIC KEYBOARD COMPONENT and a cursor. When the user deselects PHONE INPUT COMPONENT, the input is validated.			
PAGUI4035	<Picture> PHONE INPUT COMPONENT 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4036	<Picture> PHONE INPUT COMPONENT (example: enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI1356	8.12.2.43 REGISTRATION PROGRESS BAR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1357	Component: REGISTRATION PROGRESS BAR COMPONENT State: Android/iOS Display: Always displayed in PATIENT IDENTIFICATION PANEL and PATIENT INFORMATION PANEL. Conditionally displayed in CONNECT HEALTH APP PANEL when in the REGISTRATION SCREEN. Label: <label.gb.step_progress> Behavior: Text and icon. The REGISTRATION PROGRESS BAR COMPONENT displays a bar that fills as the user progresses through the registration process. In the iOS state, there are 3 total steps. In the Android state, there are only 2 total steps.	PString650 Step {23}/{24} step label.gb.step_progress	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1358	<Picture> REGISTRATION PROGRESS BAR COMPONENT (example: PATIENT IDENTIFICATION PANEL, iOS) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1359	<Picture> REGISTRATION PROGRESS BAR COMPONENT (example: PATIENT INFORMATION PANEL, iOS) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1360	<Picture> REGISTRATION PROGRESS BAR COMPONENT (example: CONNECT HEALTH APP PANEL, iOS) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4475	<Picture> REGISTRATION PROGRESS BAR COMPONENT (example: CONNECT HEALTH APP PANEL, Android) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4474	<Picture> REGISTRATION PROGRESS BAR COMPONENT (example: CONNECT HEALTH APP PANEL, Android) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI4103	8.12.2.44 REQUIRED FIELD ERROR COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4104	Component: REQUIRED FIELD ERROR COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS EMAIL ID PANEL when the user has selected the SEND VERIFICATION CODE ACTION TRIGGER COMPONENT before completing all required input. Conditionally displayed in the USER DETAILS NEW PASSWORD PANEL under the USER DETAILS ICON COMPONENT when the user has selected the CONTINUE ACTION TRIGGER COMPONENT before completing all required input. Label: N/A - determined by third party vendor - <ots.config.field_missing> Behavior: Text.	PString933 A required field is missing. Please fill out all required fields and try again. ots.config.field_missing	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4105	<Picture> REQUIRED FIELD ERROR COMPONENT (example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1066	8.12.2.45 SAVE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1067	Component: SAVE ACTION TRIGGER COMPONENT State: Disabled/ Enabled Display: Always displayed in the RECORD BP READING PANEL, RECORD WEIGHT PANEL, EDIT ACCOUNT INFO PANEL, and PRIMARY LAB PANEL.The default state is disabled. Label: <label.gb.save> Behavior: Text and trigger. Upon activation, in the disabled state, the SAVE ACTION TRIGGER COMPONENT does nothing. Upon activation, in the enabled state, the SAVE ACTION TRIGGER COMPONENT verifies the changes made to the displayed panel, saves the verified changes, and closes the displayed panel. Upon activation in the enabled state in the RECORD WEIGHT PANEL, the SAVE ACTION TRIGGER COMPONENT displays the CONFIRM INPUT PANEL. The SAVE ACTION TRIGGER COMPONENT is enabled after all required input for the displayed panel is complete.	PString195 Save save label.gb.save	UI Info	CloudDRx 1.0 CloudDRxRTR
PAGUI1068	<Picture> SAVE ACTION TRIGGER COMPONENT (default: disabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1069	<Picture> SAVE ACTION TRIGGER COMPONENT (enabled) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI4352	8.12.2.46 SEE ADJUSTMENT TRIGGER COMPONENT		Heading	CloudDRxRTR CloudDRxPS
PAGUI4353	Component: SEE ADJUSTMENT TRIGGER COMPONENT State: Enabled/ Disbled Display: Conditionally displayed in the DRX MEDICATION MESSAGE CARD PANEL and the MEDICATION MESSAGE CARD PANEL when the panel is in the active state. The default state is enabled. The SEE ADJUSTMENT TRIGGER COMPONENT changes to the disabled state when it is displayed in a panel in the out of date or pre-acknowledged state.	PString1061 See Adjustment MEDADJ label.gb.see_adj	UI Info	CloudDRxRTR

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Label: <label.gb.see_adj> Behavior: Text and trigger. Upon activation in the enabled state, the SEE ADJUSTMENT TRIGGER COMPONENT displays the MEDICATION REMEDIATION PANEL.			
PAGUI4354	<Picture> SEE ADJUSTMENT ACTION TRIGGER COMPONENT (default: Enabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI4355	<Picture> SEE ADJUSTMENT ACTION TRIGGER COMPONENT (Disabled) 		UI Design	CloudDRxRTR CloudDRxPS
PAGUI2692	8.12.2.47 SEND NEW VERIFICATION CODE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2693	Component: SEND NEW VERIFICATION CODE ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS EMAIL ID PANEL after the user selects the SEND VERIFICATION CODE ACTION TRIGGER COMPONENT and a code is sent to input of the EMAIL ID INPUT COMPONENT if it is a registered e-mail address. Conditionally displayed in the TWO FACTOR AUTHENTICATION PANEL after the user selects the SEND VERIFICATION CODE ACTION TRIGGER COMPONENT and a code is sent to input of the PHONE ID INPUT COMPONENT. Label: <ots.config_send_new_code> - determined by third party vendor Behavior: Text and trigger. Upon activation, the SEND NEW VERIFICATION CODE ACTION TRIGGER COMPONENT sends a code to the input of the EMAIL ID INPUT COMPONENT or the PHONE ID INPUT COMPONENT.	PString1081 Send new code ots.config_send_new_code	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2694	<Picture> SEND NEW VERIFICATION CODE ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2515	8.12.2.49 SCROLL DOWN ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2516	Component: SCROLL DOWN ICON COMPONENT State: Visible/Invisible Display: Conditionally displayed in the visible state in the CONFIRM CURRENT MEDICATIONS PANEL, MEDICATIONS PANEL, and MEDICATION UPDATE REMEDIATION PANEL when the panel is not scrolled to display the last element. The SCROLL DOWN ICON COMPONENT is a translucent element that is displayed over other elements. Label: N/A Behavior: Icon. When triggered, navigates user to the bottom of the page and SCROLL DOWN ICON COMPONENT dissappears, switching to invisible state.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2517	<Picture> SCROLL DOWN ICON COMPONENT (example: visible) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
PAGUI758	8.12.2.50 TREATMENT COMPLETED LABEL COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI759	Component: TREATMENT COMPLETED LABEL COMPONENT State: N/A Display: Conditionally displayed in the PAP READING NOTIFICATION CARD PANEL and the MEDICATION CHANGE NOTIFICATION CARD PANEL. The TREATMENT COMPLETED LABEL COMPONENT is displayed when the panel it is displayed in is 'complete' Label: <label.gb.treatment_completed> Behavior: Text.	PString651 You completed this treatment. completedTreatment label.gb.treatment_completed PString1017 You confirmed this medication adjustment. completedTreatment label.gb.treatment_completed	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI760	<Picture> TREATMENT COMPLETED LABEL COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3527	8.12.2.51 UPDATE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3528	Component: UPDATE ACTION TRIGGER COMPONENT State: N/A Display: Always displayed as a component of the RECOMMENDED SOFTWARE UPDATE PANEL and the REQUIRED SOFTWARE UPDATE PANEL. Label: <msg.gb.update> Behavior: Text and trigger. Upon activation, the UPDATE ACTION TRIGGER COMPONENT displays the Patient Application in the local app store.	PString674 Update update label.gb.update	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3529	<Picture> UPDATE ACTION TRIGGER COMPONENT (example: RECOMMENDED SOFTWARE UPDATE PANEL) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI3533	<Picture> UPDATE ACTION TRIGGER COMPONENT (example: REQUIRED SOFTWARE UPDATE PANEL) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2618	8.12.2.52 USER DETAILS ICON COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2619	Component: USER DETAILS ICON COMPONENT State: N/A		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	Display: Always displayed in the USER DETAILS EMAIL ID PANEL, USER DETAILS CONFIRMATION PANEL, USER DETAILS NEW PASSWORD PANEL, CONFIRM MFA METHOD PANEL, SEND MFA CODE PANEL, and ENTER MFA CODE PANEL. Label: N/A - panel is controlled by a third party. Behavior: Icon			
PAGUI2620	<Picture> USER DETAILS ICON COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2621	8.12.2.53 USER DETAILS HEADER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2622	Component: USER DETAILS HEADER COMPONENT State: N/A Display: Always displayed in the USER DETAILS EMAIL ID PANEL, USER DETAILS CONFIRMATION PANEL, and USER DETAILS NEW PASSWORD PANEL. Label: <ots.config.user_details> - determined by third party vendor Behavior: Text.	PString1072 User Details ots.config.user_details	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2623	<Picture> USER DETAILS HEADER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2679	8.12.2.54 VERIFICATION CODE INPUT COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2680	Component: VERIFICATION CODE INPUT COMPONENT State: Input/ No Input Display: Conditionally displayed in the USER DETAILS EMAIL ID PANEL when the software has successfully sent a code to the input of the EMAIL ID INPUT COMPONENT. Always displayed in the ENTER MFA CODE PANEL. The default state is no input. Label: <ots.config.verification_code> - determined by third party vendor Behavior: Text and input box. In the no input state, display the default text. In the input state, display the input given from the NUMERIC KEYBOARD COMPONENT.	PString1073 Verification code ots.config.verification_code	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2681	<Picture> VERIFICATION CODE INPUT COMPONENT (default: No Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2682	<Picture> VERIFICATION CODE INPUT COMPONENT (Input, example) 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2689	8.12.2.55 VERIFY CODE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR



myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
				CloudDRxPS
PAGUI2690	Component: VERIFY CODE ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the USER DETAILS EMAIL ID PANEL after the software successfully sends a code to the input of the EMAIL ID INPUT COMPONENT. Always displayed in the ENTER MFA CODE PANEL. Label: <ots.config.verify_code> - determined by third party vendor Behavior: Text and trigger. Upon activation in the USER DETAILS EMAIL ID PANEL, the VERIFY CODE ACTION TRIGGER COMPONENT verifies the input of the VERIFICATION CODE INPUT COMPONENT and displays the CONNECT HEALTH APP PANEL. Upon activation in the ENTER MFA PANEL in the TWO FACTOR AUTHENTICATION PANEL, the VERIFY CODE ACTION TRIGGER COMPONENT verifies the input of the VERIFICATION CODE INPUT COMPONENT and displays the HOME SCREEN when the patient has logged in before. If patient has not logged in before, it will display the POLICY SCREEN. If verification is unsuccessful, the VERIFY CODE ACTION TRIGGER COMPONENT displays the MFA MAXIMUM ATTEMPTS DIALOG PANEL over the STARTUP SCREEN.	PString1074 Verify Code ots.config.verify_code	UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI2691	<Picture> VERIFY CODE ACTION TRIGGER COMPONENT 		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1729	8.12.2.56 VIDEO ARTICLE ACTION TRIGGER COMPONENT		Heading	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1730	Component: VIDEO ARTICLE ACTION TRIGGER COMPONENT State: N/A Display: Conditionally displayed in the PAP READING INSTRUCTIONS PANEL, ABOUT HEART FAILUE PANEL, ABOUT HEART FAILURE SYSTEM PANEL, and IN-APP TRAINING PANEL. The VIDEO ARTICLE ACTION TRIGGER COMPONENT is displayed when there is a video article available in the content management location for the particular panel. The number of VIDEO ARTICLE ACTION TRIGGER COMPONENTS displayed in a panel is determined by how many unique video articles there are stored in the content management location. The VIDEO ARTICLE ACTION TRIGGER COMPONENT displays a thumbnail and title for the associated video article. Label: N/A - determined by content manager Behavior: Trigger. Upon activation, the VIDEO ARTICLE ACTION TRIGGER COMPONENT displays the VIDEO PLAYER PANEL for the associated article.		UI Info	CloudDRx 1.0 CloudDRxRTR CloudDRxPS
PAGUI1731	<Picture> VIDEO ARTICLE ACTION TRIGGER COMPONENT (example)		UI Design	CloudDRx 1.0 CloudDRxRTR CloudDRxPS

myCardioMEMS Patient Application Graphical User Interface Specification

ID		English US Strings	Object Type	Project
	 Looking into the future: Beth's CardioMEMS™ HF story			



myCardioMEMS Patient Application Graphical User Interface Specification



myCardioMEMS Patient Application Graphical User Interface Specification



myCardioMEMS Patient Application Graphical User Interface Specification