

Patient Registration and Activation with CardioGuide and Merlin

Open Q3 Merlin portal [Merlin.net™ PCN](https://q3-hfcloud.products.abbott/apps/clinician/notifications)

- Username: q3_swvq3
- Password: Test!1234

Press the blue “Enroll Patient” button on the top right of the screen

Merlin.net™ PCN

https://q3-hfcloud.products.abbott/apps/clinician/notifications

Abbott Notifications All Patients Clinic Merlin.net™ Patient Care Network Heart Failure Portal Help SIGN OUT

Welcome, q3_test_1 q3_test_2

565 patients

565 Optimization

0 Maintenance

DYNAMIC TREATMENT PLAN

Dynamic Treatment Plan

139 enrolled

426 not enrolled

Filter by Clinician: Me

Search

Notifications: 317

FILTER BY: 317 All 231 High priority 86 Low priority

Patient / Clinician	Dynamic Treatment Plan Enabled	Notification / Date	Goal / Type	Last Measurement	Last Reading	PA Trend (Last 7 days)	Actions
ab_1, Verif DOB: 01-01-2001 No subscribing clinicians	Dynamic Treatment Plan	Patient has multiple notifications. To triage and clear, go to the Patient Page	16mmHg PA Diastolic	06-15-2026 PAP	24mmHg		
cg, TestAutomation_1606440 DOB: 01-01-2020 No subscribing clinicians		Heart Rate out of threshold/06-16-2026 One or more Suspect Readings/06-16-2026 Sensor Calibrated: Review PA Pressure Settings and Dynamic Treatment Plan/06-16-2026 1 reading above optimal range/06-16-2026	16mmHg PA Diastolic	06-16-2026 PAP	22mmHg		

Enter your patient’s first name, last name, and date of birth. Make the month and date of the birthday the same to avoid confusion when using POSTMAN. Press the blue “ENROLL” button

Abbott Notifications All Patients Clinic Merlin.net™ Patient Care Network Heart Failure Portal Help SIGN OUT

Enroll by Name and Date of Birth

First Name
* Sam14


Last Name
* RTR

Date of Birth
* 01-01-2000

* Required field

CANCEL ENROLL


You will now be directed to the patient profile section. Open the Sensors Portal backend [Sensors | CMHF Customer Care](#) and sign in with your Abbott account. Select the “Unclaimed” tab and copy one of the sensor serial numbers.

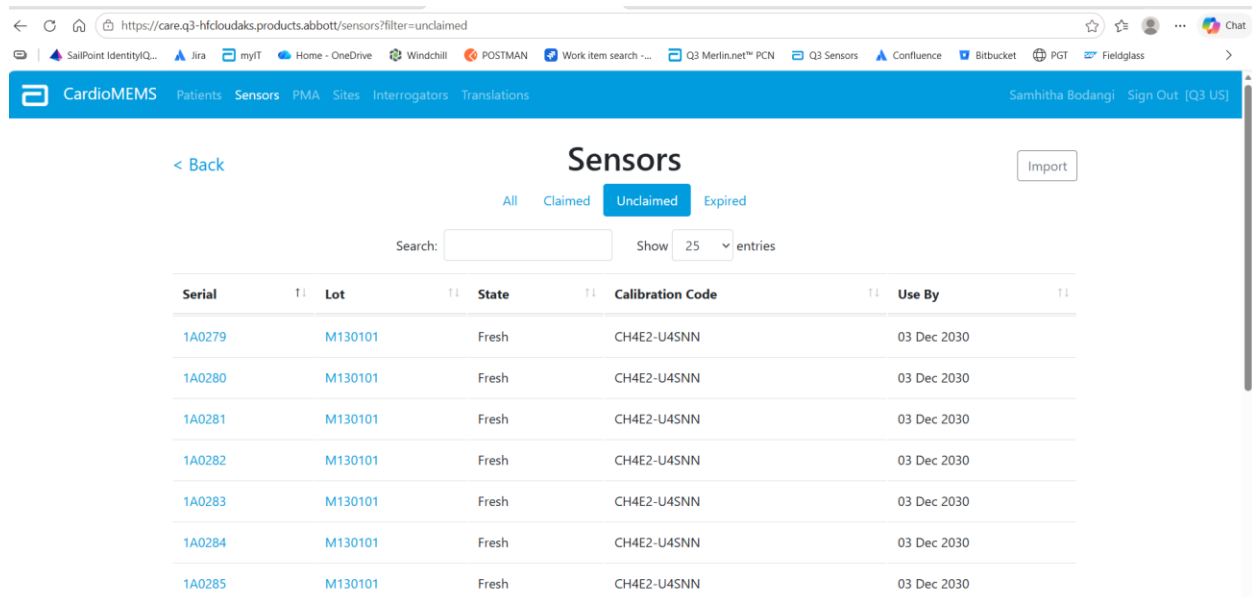

The screenshot shows the 'Sensors' page in the CardioMEMS portal. The 'Unclaimed' tab is selected. A table lists 7 unclaimed sensors with the following columns: Serial, Lot, State, Calibration Code, and Use By.

Serial	Lot	State	Calibration Code	Use By
1A0279	M130101	Fresh	CH4E2-U4SNN	03 Dec 2030
1A0280	M130101	Fresh	CH4E2-U4SNN	03 Dec 2030
1A0281	M130101	Fresh	CH4E2-U4SNN	03 Dec 2030
1A0282	M130101	Fresh	CH4E2-U4SNN	03 Dec 2030
1A0283	M130101	Fresh	CH4E2-U4SNN	03 Dec 2030
1A0284	M130101	Fresh	CH4E2-U4SNN	03 Dec 2030
1A0285	M130101	Fresh	CH4E2-U4SNN	03 Dec 2030

Go back to your patient profile page in Merlin and paste the sensor serial number in the “Sensor Serial #” box. Select “q3_test_2, q3_test_1” as the Testing Physician and press save.

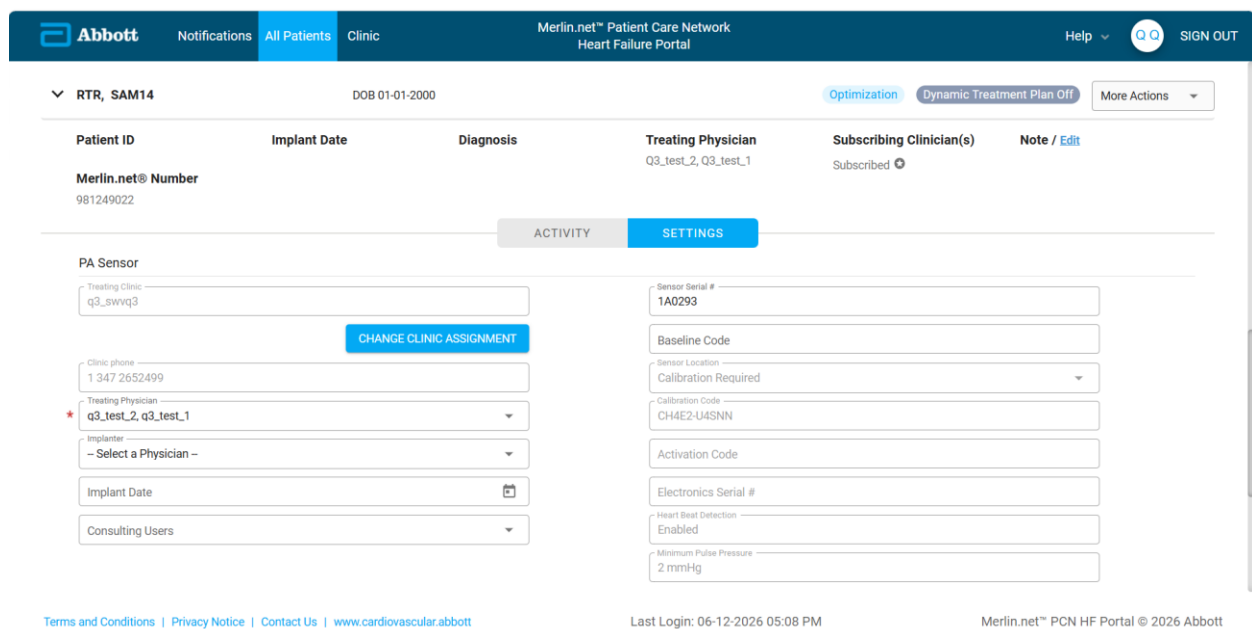

The screenshot shows the Merlin.net Patient Care Network interface. The patient profile for RTR, SAM14 (DOB 01-01-2000) is displayed. The 'Settings' tab is active, showing the 'PA Sensor' configuration. The 'Sensor Serial #' field is populated with '1A0293'. The 'Treating Physician' is set to 'q3_test_2, q3_test_1'. The 'Implant Date' is set to '1 347 2652499'. The 'Consulting Users' field is empty.

Patient Details:

- Patient ID: RTR, SAM14
- DOB: 01-01-2000
- Merlin.net® Number: 981249022

Sensor Settings (PA Sensor):


- Treating Clinic: q3_svvq3
- Baseline Code: [Empty]
- Sensor Location: Calibration Required
- Calibration Code: CH4E2-U4SNN
- Activation Code: [Empty]
- Electronics Serial #: [Empty]
- Heart Beat Detection: Enabled
- Minimum Pulse Pressure: 2 mmHg

In the Sensor backend, tap on the sensor serial number you pasted and it will open up the Sensor page. The patient name should be the name of your patient


The screenshot shows the CardioMEMS web application interface. At the top, there's a blue header with the CardioMEMS logo and navigation links: Patients, Sensors, PMA, Sites, Interrogators, and Translations. On the right of the header, it says 'Samhitha Bodangi' and 'Sign Out [Q3 US]'. The main content area has a title '< Sensors' and a large sensor ID '1A0293' with a 'Fresh' status. Below the ID, there are two columns of sensor data. The left column includes: Serial Number: 1A0293, Calibration Code: CH4E2-U4SNN, Baseline Code: ----, Activation Code: 4909, Patient: Sam14 RTR, and Implant Location. The right column includes: Current Search Pressure: ----, Set Search Pressure: ----, Min Signal Strength: 70.0, Heart Beat Detection: Enabled, Minimum Pulse Pressure: 2, Patm Baseline: ----, MFG SN: 2BE0010062, MFG Cal Code: 3347-077850-02836-0850, Lot #: M130101, MFG Date: 01 Jan 2023, Use by Date: 03 Dec 2030, Slope One: 77.85, Slope Two: -2.836, and F0 @ 860: 33.5. Below the sensor data, there are two sections: 'Recent Readings' and 'Calibrations', both stating 'No recent readings found.' and 'No calibrations found.' respectively.

Sensor Details for 1A0293	
Serial Number:	1A0293
Calibration Code:	CH4E2-U4SNN
Baseline Code:	----
Activation Code:	4909
Patient:	Sam14 RTR
Implant Location:	
Current Search Pressure:	----
Set Search Pressure:	----
Min Signal Strength:	70.0
Heart Beat Detection:	Enabled
Minimum Pulse Pressure:	2
Patm Baseline:	----
MFG SN:	2BE0010062
MFG Cal Code:	3347-077850-02836-0850
Lot #:	M130101
MFG Date:	01 Jan 2023
Use by Date:	03 Dec 2030
Slope One:	77.85
Slope Two:	-2.836
F0 @ 860:	33.5

Download the CardioGuide app on a phone. If on Android, you can download using apk file: [myCardioMEMS-q3-2.0.23.apk](#). The most recent build as of 6/16/2026 is 2.0.23. You can also download the app via google Play Store link on Android or Test Flight on IOS

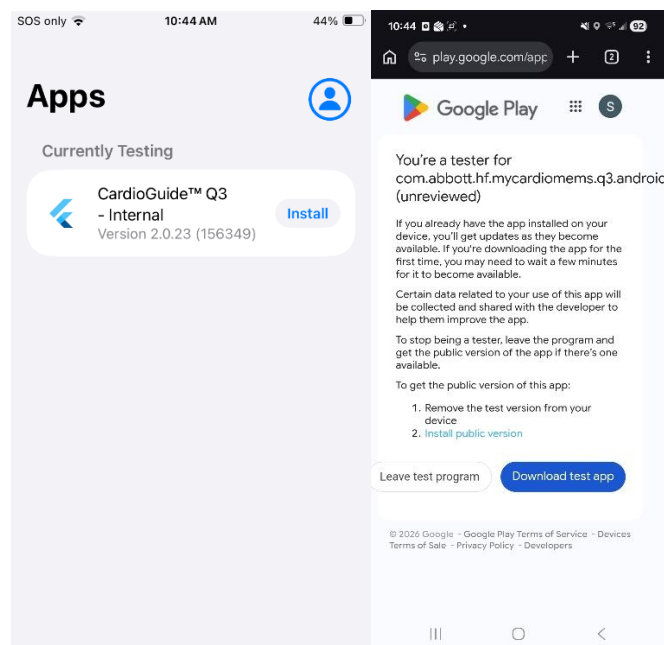


Once you install and open the app, you will be directed to the Startup Screen. Press on the “Register” button. Then press “Continue”. Scroll down to the end of the EULA document and press “Accept”. Scroll down to the end of the Privacy Policy document and press “Acknowledge”

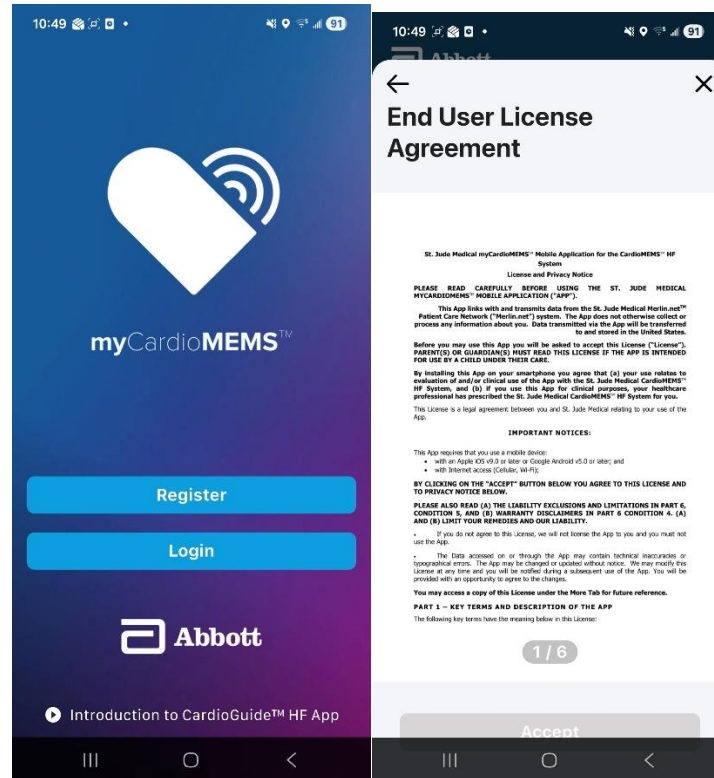


Enter the correct date of birth and Sensor Serial Number of you patient and press “Next.” Enter the full name (first and last), email address, phone number, and create a password. Then press “Register”.

The image displays two side-by-side screenshots of a mobile application interface for patient registration. The left screenshot, titled "Patient Registration" (Step: 1/2), shows input fields for "Date of birth" (01-01-2000) and "Sensor Serial Number" (1A0293). Below these is an informational note: "The sensor Serial Number can be found on the top left corner of the Patient ID card." This note is accompanied by an image of a Patient ID card with the Abbott logo and a red box highlighting the "SN" (Serial Number) field. A blue "Next" button is at the bottom. The right screenshot, titled "Patient Information" (Step: 2/2), shows input fields for "Email Address" (sam14@gmail.com) and "Phone Number" (+1 5086679568). It also includes a "Password" field with a strength indicator showing four checked criteria: "At least 8 characters", "One upper case letter", "One lower case letter", and "One number or symbol". A "Confirm Password" field is also present. Both screens have a back arrow and a close 'X' button at the top.

You are now directed to the user authentication screen. Enter the email address and the phone number for the patient you just created. You can toggle on the “Keep me logged in toggle”. This will ensure user stays logged in for 90 days, instead of automatically logging out after 24 hours.

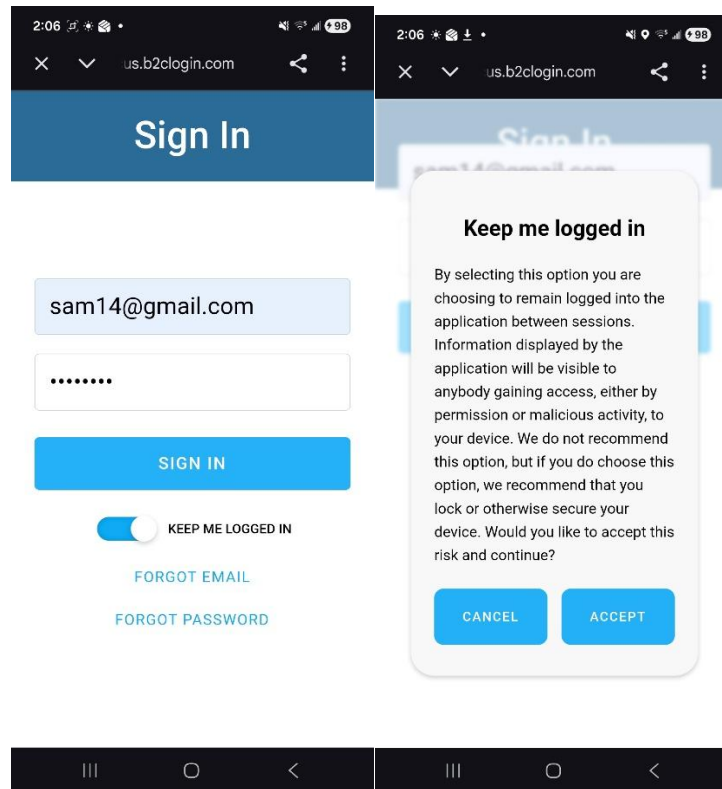

Press “Send Code” You should have gotten a code sent to the phone number associated with your patient. Enter the verification code and press “verify code”

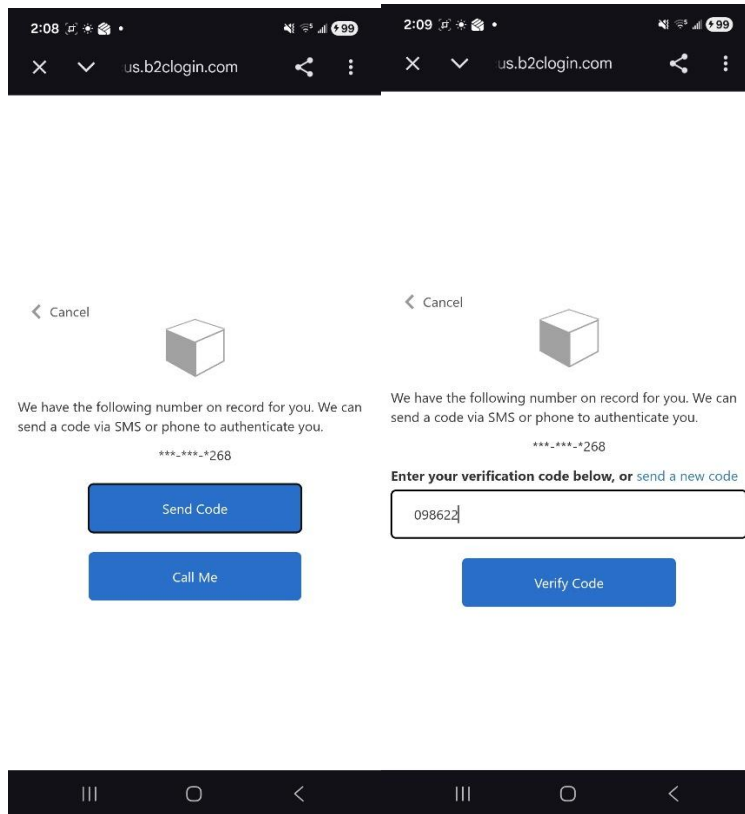

NOTE: on app builds prior to 2.0.23, you will be directed to the login screen after registering the patient. On 2.0.23, you will have to verify the phone number twice before landing on the home screen.

You should now be directed to the CardioGuide home screen. Open Merlin, and open the “CardioGuide HF App” section.” The “Patient” bar should say “Registered”. You can press “View Registered Application to see the name and email address on the app.

▼ RTR, SAM14 DOB 01-01-2000 Optimization Dynamic Treatment Plan Off More Actions ▼

Patient ID	Implant Date	Diagnosis	Treating Physician	Subscribing Clinician(s)	Note / Edit
Merlin.net® Number 981249022			Q3_test_2, Q3_test_1	Subscribed ⓘ	

ACTIVITY **SETTINGS**

★ Required field (view or enter details in Pacemaker/ICU section above) CANCEL SAVE

▼ CardioGuide™ HF App

Patient: Registered VIEW REGISTERED APPLICATION USER
 Medications managed on the mobile application: Off
 Send PA reading data to the mobile application: Off

> Goals & Treatment Plan
 > Medication History
 > Diagnosis
 > Pacemaker/ICU

Abbott Notifications All Patients Clinic Merlin.net™ Patient Care Network Heart Failure Portal Help SIGN OUT

▼ RTR, SAM14 DOB 01-01-2000 Optimization Dynamic Treatment Plan Off More Actions ▼

Patient ID	Implant Date	Diagnosis	Treating Physician	Subscribing Clinician(s)	Note / Edit
Merlin.net® Number 981249022			Q3_test_2, Q3_test_1	Subscribed ⓘ	

★ Required field CANCEL SAVE

▼ CardioGuide™ HF App

Patient: Registered
 Medications managed on the mobile application: Off
 Send PA reading data to the mobile application: Off

> Goals & Treatment Plan
 > Medication History
 > Diagnosis
 > Pacemaker/ICU

View registered application user ⓘ

CardioGuide™ HF App

The following user has registered the CardioGuide™ HF App for use with this patient

Name ▼	User Type	Medication Manager	Registered on	Email
Sam14 RTR	PAT	Yes	06-16-2026	sam14@gmail.com

Show per page 10 1 - 1 of 1 < >

CLOSE

Turn both the “Medications managed on the mobile application” and “Send PA reading data to the mobile application” On and press Confirm for both.

Goals & Treatment Plan

Set Goal

PA Diastolic

Optimization

TRENDS

PA Diastolic

PA Mean

12

High

20

16

12

CANCEL

SAVE

Right Heart Cath Implant Values

Last update:

RA Mean

mmHg

RV End Diastolic

mmHg

RV Systolic

mmHg

PA Systolic

mmHg

PCWP

mmHg

PA Diastolic

mmHg

DPG

n/a

mmHg

Cardiac Output

L/min

PA Mean

mmHg

PVR

n/a

Wood Units

[Terms and Conditions](#) | [Privacy Notice](#) | [Contact Us](#) | [www.cardiovascular.abbott](#)

Last Login: 06-18-2026 03:33 PM

Merlin.net™ PCN HF Portal © 2026 Abbott

In Merlin, scroll down to “Daily Medications” and press “Edit Daily Medications.” Press Confirm. Press the blue “+ Medication” button. Go to the “Formulary” tab and choose any medication, defining its strength, quantity, and frequency. You can add more medications using the “Save and Add Another” button and repeating the process of defining prescriptions. Then press “Save Medication”. Then Press “Save Daily Medications” and “Save”

Abbott

Notifications

All Patients

Clinic

Merlin.net™ Patient Care Network
Heart Failure Portal

Help

Q Q

SIGN OUT

ADD DAILY MEDICATION

Search Medications

FAVORITES

FORMULARY

Category	Class	Drug	Strength	Qty	Frequency
ADHD/Anti-Narcolepsy /Anti-Obesity/Anorexi...	ACE Inhibitors	Bempedoic Acid(Nexletol)	180MG TABS PO	0.5	BID (Twice a day)
Analgesic/Anti-Inflammatory/Migraine/Gout ...	Adenosine Triphosphate-Citrate Lyase (ACL) L...			1	QAM (Every day bef...
Cardiovascular Agents	Agents for Pheochromocytoma			1.5	QD (Every day)
Endocrine and Metabolic Agents	Alpha-Beta Blockers			2	QHS (Every night at ...
Hematological Agents	Anaphylaxis Therapy Agents			2.5	QID (Four times a da...
Nutritional Products	Angiotensin II Receptor Antagonists			3	QOD (Every other d...
				3.5	QPM (Every day afte...

Medication update will be sent to the patient application. Daily Medications will be displayed in the patient's application when medication is managed within the application. When dynamic treatment plan is ON, the medications are displayed when readings are within Optimal Range.

CANCEL

SAVE AND ADD ANOTHER

SAVE MEDICATION

Terms and Conditions

Privacy Notice

Contact Us

www.cardiovascular.abbott

Last Login: 06-18-2026 03:33 PM

Merlin.net™ PCN HF Portal © 2026 Abbott

▼ RTR, SAM14

DOB 01-01-2000

Optimization

Dynamic Treatment Plan Off

More Actions

Patient ID	Implant Date	Diagnosis	Treating Physician	Subscribing Clinician(s)	Note / Edit
Merlin.net® Number 981249022			Q3_test_2, Q3_test_1	Subscribed	

ACTIVITY

SETTINGS

• Bempedoic Acid(Nexletol) 180 MG, 1 TABS PO QD (Every day)

• Phenoxybenzamine HCl 10 MG, 1 CAPS PO QHS (Every night at bedtime)

+ MEDICATION

CANCEL

SAVE DAILY MEDICATIONS

INTERVENTIONS

1 READINGS OUT OF RANGE

3 READINGS OUT OF RANGE

7 READINGS OUT OF RANGE

High Range

over 20 mmhg

Low Range

under 12 mmhg

EDIT INTERVENTIONS

Terms and Conditions

Privacy Notice

Contact Us

www.cardiovascular.abbott

Last Login: 06-18-2026 03:33 PM

Merlin.net™ PCN HF Portal © 2026 Abbott

Scroll down and press “Edit interventions.” Select “+ Medication/Message” for any out of range category. You can titrate a daily medication, add a new one, or add a message.

▼ RTR, SAM14

DOB 01-01-2000

Optimization

Dynamic Treatment Plan Off

More Actions ▼

Patient ID	Implant Date	Diagnosis	Treating Physician	Subscribing Clinician(s)	Note / Edit
Merlin.net® Number 981249022			Q3_test_2, Q3_test_1	Subscribed ⓘ	

ACTIVITY

SETTINGS

EDIT DAILY MEDICATION

INTERVENTIONS

1 READINGS OUT OF RANGE

3 READINGS OUT OF RANGE

7 READINGS OUT OF RANGE

High Range
over 20 mmhg

+ MEDICATION/MESSAGE ▼

+ MEDICATION/MESSAGE ▼

+ MEDICATION/MESSAGE ▼

Low Range
under 12 mmhg

- Titrate Daily Medication
- Add New Medication
- Add Message ▶

+ MEDICATION/MESSAGE ▼

+ MEDICATION/MESSAGE ▼

CANCEL

SAVE INTERVENTIONS

Dynamic Treatment Plan: ⓘ Not eligible to enable Dynamic Treatment Plan

- 1) CardioGuide™ HF App: Registered
- 2) CardioGuide™ HF App Medication Manager: On
- 3) Daily Medications: Created