

Setup Notes

Tuesday, June 2, 2026 10:30 AM

Appium server is different from appium inspector

Install packages and softwares directly into C drive, not in Users

Always keep appium running in background when running automated test cases

- Run command *appium* in cmd or IntelliJ IDEA before starting session in Appium Inspector

Always work on you own branch first (make separate branch and make edits, do not work directly on the main branch)

Locator = used to identify specific element on webpage

X Path = XML path language

JSON To paste In Appium Inspector

```
{
  "platformName": "Android",
  "appium:platformVersion": "16",
  "appium:deviceName": "Samsung Galaxy A56",
  "appium:udid": "RZCY826JDLJ",
  "appium:automationName": "UiAutomator2",
  "appium:appPath": "C:\\work\\myCardioMEMS-q3-2.0.23.apk"
}
```

Mvn clean install -DskipTests

- Downloads all libraries
- Should say BUILD SUCCESS

Update branch with Global changes

Git Commands

Saturday, July 25, 2026 11:05 AM

Git Commands

Git branch

- Lists all local branches and highlights the one you are currently in

Git branch <name>

- Creates a new branch with that name
 - o Git branch SamhithaBranch

Git stash

- Saves local changes (does not save them to main branch, just the one on local computer)

Git stash drop

- Permanently deletes saved changes from history
- Do git stash first before doing git stash drop

Git pull

- Downloads changes from repo to current local branch

Git clone --branch CardioGuide2.0-Web-API_Integrations

<https://bitbucket.sjm.com/myc/cardioguideapp-testautomation.git>

Cd cardioguideapp-testautomation

Git checkout -b <name>

- Git checkout Allows you to switch between different branches
- Git checkout -b allows you to create a new branch and switch into it

Making PR

Git status

- Can see all your files that were modified. Red means they have not been staged/committed

```
Changes to be committed:
(use "git restore --staged <file>..." to unstage)
    deleted:    sensorgenerator/sensor import final.csv

Changes not staged for commit:
(use "git add <file>..." to update what will be committed)
(use "git restore <file>..." to discard changes in working directory)
    modified:   src/test/java/com/cardioguideapp/apppageactions/NotificationsPageActions.java
    modified:   src/test/java/com/cardioguideapp/tests/notification/DRXTMAPP5734Notification9_5Test.java
    modified:   src/test/resources/StringLabels.xlsx
    modified:   src/test/resources/androidmobilelocators/t1/LandingPage.properties
    modified:   src/test/resources/androidmobilelocators/v3/LandingPage.properties
    modified:   src/test/resources/framework.properties
    modified:   src/test/resources/iosmobilelocators/q3/LandingPage.properties
    modified:   src/test/resources/iosmobilelocators/t1/LandingPage.properties
    modified:   src/test/resources/iosmobilelocators/v3/LandingPage.properties
    modified:   src/test/resources/iosmobilelocators/v3/LandingPage.properties
    modified:   src/test/resources/suites/sanity.xml
    modified:   src/test/resources/testdata/DRXTMAPP5734Notification9_5Test.json

Untracked files:
(use "git add <file>..." to include in what will be committed)
    sensorgenerator/sensor import final.csv
    src/test/java/com/cardioguideapp/apppageactions/MoreOptionSettingsPageActions.java
    src/test/java/com/cardioguideapp/apppageobjects/MoreOptionSettingsPage.java
    src/test/java/com/cardioguideapp/demoTests/DebugHealthAccessTest.java
    src/test/java/com/cardioguideapp/mobilewrappers/MoreOptionSettingsPageActionsTest.java
    src/test/java/com/cardioguideapp/testdataconstants/DRXTMAPP415Registration13TestData.java
    src/test/java/com/cardioguideapp/tests/registration/DRXTMAPP415Registration13Test.java
    src/test/resources/iosmobilelocators/q3/moreOptionSettingsPage.properties
    src/test/resources/testdata/DRXTMAPP415Registration13Test.json
    src/test/resources/testdata/DebugHealthAccessTest.json
```

Git add <path of file>

- Adds files to stage

Git commit -m "message"

- Commits all the files you added

Git push

- Pushes your latest commit to your remote branch

Git rm --cached <path of file>

- Removes a file that has been pushed, but does not remove it from your local branch
- You need to do git commit and git push for your remote branch to be updated

Git push -u origin <branch name>

- Push to a branch for the first time, can use git push after

Tuesday, June 2, 2026 11:10 AM

Test case SignIn9, verifies workflow from Startup Screen to Patient Registration screen
Ends after pressing EULA and PP buttons

- Entering Bday and SN is after user accepts EULA and PP

In src>test>resources>testdata>DRXTMAPP16SignIn9Test.json

No information needed to change, but follow similar method to edit username and password of user to be logged into app

???

qa.name = Vikramraj Chandran
qa.emailId = vikramraj.chandran@abbott.com
qa.511Id = CHANDVX26
qa.appVersion = 2.0.18
qa.platform = android
qa.deviceModel = samsung s23
qa.osVersion = 15.0

Changed to

qa.name = Vikramraj Chandran
qa.emailId = vikramraj.chandran@abbott.com
qa.511Id = CHANDVX26
qa.appVersion = 2.0.20
qa.platform = android
qa.deviceModel = Samsung Galaxy S24
qa.osVersion = 16

Run test case

Run appium server (*appium*)
Start session in Appium inspector
Open new terminal in IntelliJ
Run correct command below (replace UIUD and path)

For Android :

MAC/Linux = mvn clean test -DplatformName=Android -Ddevice.udid=RZCXA1RK7FV -Dapp=/Users/vikramrajchandran/Documents/V3apk/myCardioMEMS-qa-2.0.14.apk

Windows = mvn clean test "-DplatformName=Android" "-Ddevice.udid=RFCY90ZGW6A" "-Dapp=C:\work\myCardioMEMS-qa.apk"

For IOS :

MAC/Linux = mvn clean test -DplatformName=iOS -Ddevice.udid=00008101-000E7D6E34EA001E

RESULTS

- App was downloaded
- It was opened
- All files access (allow permission was toggled on)
- Next went to startup screen
- Manual Verification required
 - o Character boldness cannot be automated
- Manual verification required again
 - o Trademark symbol rendered cannot be automated
- Opens Policy Welcome Panel
- Pressed Continue
- Goes to EULA
- Fails here
 - o Most likely xpath error

When opening src>test>java>com.cardioapp>tests>signin>DRXTMAP16SignIn9Test
Project JDK is not defined

On windows press Ctrl + Alt + Shift + S
Go to Project Settings>Project

Select the SDK (should be 22.0.2)
Press Apply then Ok

** Did not have app downloaded for first run, but not deleting it for 2nd run

Test Run Take 2 RESULTS

- Closed out old app
- Deleted app
- Reinstalled app
- Opened app and toggled allow permission
- Opens startup page
- Goes to Manual verification
- Goes to 2nd manual verification
- Goes to Policy screen
- Goes to EULA

As of June 2, 2026 the locators for EULA changed from Acknowledge to Consent and Continue
The 2.0.20 build still has Acknowledge button, so must change xpath for test to run correctly

Go to src>test>resources>androidmobilelocators>v3>eulaPage.properties
Change xpath based on Appium Inspector

Edited privacyPolicyPage.properties too

Goes to patient registration screen and the closes app
Says BUILD SUCCESS

Go to allure-report>index.html
Right click, open in browser

Ex Automating Test Case

Tuesday, June 2, 2026 2:57 PM

Ex: [Pull Request #248: Signin20Imr and Signin4Imr - Bitbucket](#)

Jira Ticket: [DRXTMAPP Tests - Kanban Board - Jira](#)

Testing to Run Finished Test Case

Trying to run it on 2.0.22

- Shuould use q3 env for this build

- Entering v3 email and password will not work

- Will run test case on 2.0.20

Original Index 0 patient sensorID = 1C4308

- Changed to 1A0191 (not being used so most likely not make difference)

- Cardioguideone20@gmail.com Is not valid email and pasword, change Index 0 (default user) to cardioguideone21@gmail.com

Manual intervention is required

Enter new patient and use manual OTP (have patient registered to personal phone num)

Copilot Test Case Tasks

Wednesday, June 3, 2026 10:29 AM

Gap analysis

Full blown test plan (automatic, semi automatic, manual)

Ask for rationale too

myCMEMSApp branch is development

Full Test Plan

Prompt 1 (Agent Mode):

Attached the mycmems-app folder and the requirements pdf

The folder mycmems-app is an app made in flutter and dart that helps clinicians manage patient care based on patient PA pressure readings. The CG_Requirements.pdf is a file that has all of the requirements and components that should be in the app and describes its functionality and design. Give me a document that is a comprehensive list of manual, automatic, and semi automatic test cases I should run to make sure all of the requirements are met in the app

Result

- Analyzed requirements using pdf parser
 - o Python -m pip install pypdf -q
- Parsed pdf and went through section headings and requirement IDs
- Identified primary requirement sections as Startup, Policy, Registration, Auth, Home, Readings, Medications, Notifications, More
- Created 4 documents
 - o Cg_requirements_extracted.txt = parsed requirement pdf as txt file
 - o Cg_requirements_high_level_headings.txt
 - o Cg_requirements_mined_summary.txt
 - o Cg_requirements_screen_headings_raw.txt
- Created a "testing" folder in mycmems-app>documentation
 - o Test cases list is in newly created CG_Requirements_Test_Catalog.md doc

Prompt 2:

Attached CG_Requirements_Test_Catalog.md

Make this catalog into a pdf and create separate tables to make it easier to read

Result

- Python -m pip install reportlab -q
- Created a CG_Requirements_Test_Catalog.pdf in mycmems-app>documentation>testing

Initial Observations from Test Plan:

- o Formatting is bad (will work in markdown file for now)

6.3 Widget/UI Automation

ID	Area	Candidate Test/Path	Pass Criteria
A-041	Welcome screen	test/widget_test/login/welcome_screen_test.dart	Enter CTAs and key UI states pass
A-042	Consent screens	test/widget_test/register/consent_screen_test.dart	On/Off/Back/next flow widgets pass
A-043	Home screen	test/widget_test/home/home_screen_test.dart	Test state render in expected states
A-044	Readings screen	test/widget_test/readings/screen/readings_screen_test.dart	Readings, chart, test, and BP area render correctly
A-045	Readings overlays	test/widget_test/readings/widgets/overlay_test.dart	Overlay test, error overlays pass
A-046	Med adjustment confirmation	test/widget_test/medications/medication_confirmation_test.dart	Confirmation UI and OK/Cancel buttons pass
A-047	Notifications screen	test/widget_test/notifications/notification_test.dart	Tab switching and board render pass
A-048	Settings screen	test/widget_test/account/settings_screen_test.dart	Settings test, actions render correctly
A-049	Contact/support/FAQ screens	test/widget_test/help/_screen_test.dart	Content and interaction widgets pass
A-050	Weight screen	test/widget_test/more/weight/weight_test.dart	Overlay, test, and entry state rendering pass

- Devices include both phones and tablets
- Separated into 18 different sections, about 6 sections each for manual, automatic, and semi automatic
- Language can be confusing (e.g. "data corruption")
- Has a lot of sign in and MFA requirements as manual

Prompt 3

Attached the cardioguideapp-testautomation folder

For testing we will only test on IOS iPhones and Android phones. The .md catalog has some confusing wording such as "Data corruption" and "Deterministic data". Make the language clearer and explicit as possible on what it is trying to verify. In the expected results section please reference the exact screen that it should land on if the test case passes or description if needed. The cardioguideapp-testautomation folder attached is the in-progress test automation we are working on. The src>test>tests folder is where we have the classes of automated test cases. In the test catalog, make the test cases automatic and only manual or semi-manual if necessary.

Results

- Still making many of the navigation requirements as manual
- Included column in Automated test case section that includes path to currently automated test cases in the cardioguideapp-testautomation folder

Prompt 4

Separate the list of test cases to run based on the main requirements of the app. Try to make the test cases automated and if necessary they can have manual verification within those cases. It is possible to automate most of the navigation of the app. Additionally, if test cases must be manual or semi manual give me a quick justification or rationale on why they cannot be automated.

Result

- Separated
- Still has a "Automatic test case" section where it shows path of classes
- Does not include any additional automated test cases, just the ones that are in the testautomation folder

TAKE 2: Using fresh chat

Prompt 1

Attached the mycmems-app folder and the requirements pdf

The folder mycmems-app is an app made in flutter and dart that helps clinicians manage patient care based on patient PA pressure readings. The CG_Requirements.pdf is a file that has all of the requirements and components that should be in the app and describes its functionality and design. Give me a document that is a comprehensive list of any manual, automatic, and semi automatic test cases I should run to make sure all of the requirements are met in the app. It is possible to automate the

navigation and verify the placement/presence of most components, buttons and images on the screens. Therefore, make test cases manual only if necessary, try to make most of them automatable. Use clear language and reference the exact names of screens or components that it is trying to verify. Separate the list of test cases to run based on the main requirements of the app.

***Define categories:

Result

- Separated test cases into 12 different section, each section has many automated cases, some semi and few manual
- Only lists the test case and ID, no expected result
 - o Has short description on what the test case is verifying
- References the screens by their all caps name

Prompt 2

Format the document into tables, similar to a traceability requirement matrix. Include the test case ID, a column of if the case is automatic, semi, or manual, a short description and what the test case is testing for and an expected result column

Prompt 3

Attached Drxtmapp4880E2eTest.java as Example Automated Test case

This is an example automated test case that has already been completed and what steps require manual intervention. Based on this edit the test case document

Results

- Says "matches requirement expectations" for the expected results columns

Prompt 4

Instead of "matches requirement expectations" be more specific on exactly what the expected result should show based on the app requirements

Prompt 5

Attached Drx_Patient_App_Test_Case_Summary.pdf

Use Drx_Patient_App_Test_Case_Summary.pdf as a template to create the test case document. Make sure to include all type of tests as following and more :

Boundary – Testing values that are outside the range of a specific software parameter and the software action that will occur when that value is selected or programmed, and the related function is exercised.

Edge Case – Testing under extreme, rare, or boundary conditions – situations that are unlikely but possible.

Use Case Testing – Testing the software based on real-world scenarios or workflows that reflect how end users (e.g., clinicians, patients, and/or technicians) will interact with the software to achieve a specific goal.

Fault Condition – Testing how software behaves when faults or failures occur, such as invalid inputs, communication errors, or unexpected user actions that is converse to conditions specified in the requirements.

Functional – Testing that ignores the internal mechanism or structure of a system or component and focuses on the output in response to selected inputs and execution.

Stress Testing – Testing for out-of-bounds conditions, testing system performance, or testing that exercises critical workflows under unusual or high load conditions.

Results

- Limited number of test cases

Prompt 6

Citation: [Best Practices for Writing Clear and Precise Test Cases](#)

The test case list is not as comprehensive as needed. Make sure to include all of these test cases as well:

Positive Test Cases: These test the system's expected behavior when valid inputs are provided.

Negative Test Cases: These test the system's behavior when invalid inputs are provided.

Boundary Test Cases: These test the edge limits of input fields.

Performance Test Cases: These test the system's performance under specific conditions like load, stress, or endurance.

Usability Test Cases: These ensure the software is user-friendly and meets UI/UX standards.

Security Test Cases: These test how well the application protects sensitive data and resists unauthorized access.

Prompt 7

Convert the document into a pdf document using the DRx_Patient_App_Test_Case_Summary.pdf as a template. Use the same color scheme and format.

Results: Still some formatting issues.

Prompt 8

Organize the test cases in the same way as DRx_Patient_App_Test_Case_Summary.pdf. For example, have a main header as the Test Design Type (such as Use Case testing, Boundary, Functional, etc). Then have such headings for each module within that main header. For the ID, use a similar naming convention where it references the Test Design Type and Module. It is not necessary to reference the execution type. When using the dark blue table headings, use white text color to make it easier to read

Prompt 9

Use white for the table row with dark blue background. For the Test Case Description, use similar language as the DRx_Patient_App_Test_Case_Summary.pdf test case descriptions. You do not need to include the "Requested QA Class" column or information

Results: Still not perfectly organized, will fix

Prompt 10

Attached is a folder of automated test cases that have already been completed. Cross check all cases with the comprehensive test plan you created. Tell me what cases have been successfully automated and what cases are still needed to finish.

Results: Just gave me a list of them in the chat, not detailed analytics

** Some test cases in document are combined: Verify that crash-free startup under no-network and slow-network conditions

Prompt 11

Similar to the last few pages of the DRx_Patient_App_Test_Case_Summary.pdf document, give me detailed information on what test cases have been run with a short description and what the test case was trying to verify. Perform a Gap analysis telling me how many test design types are missing.

Additionally, the test cases in your document have some that are combined (Verify that crash-free startup under no-network and slow-network conditions). I want all of the test cases to be split so they are only verifying one thing.


Testing Notes

Wednesday, June 17, 2026 12:34 PM

IP address changes based on WiFi, change SMS forwarder based on ipconfig results


Download the reports from the report-export folder

- This will be a zip folder with a html and pdf of both the allure and extent folder

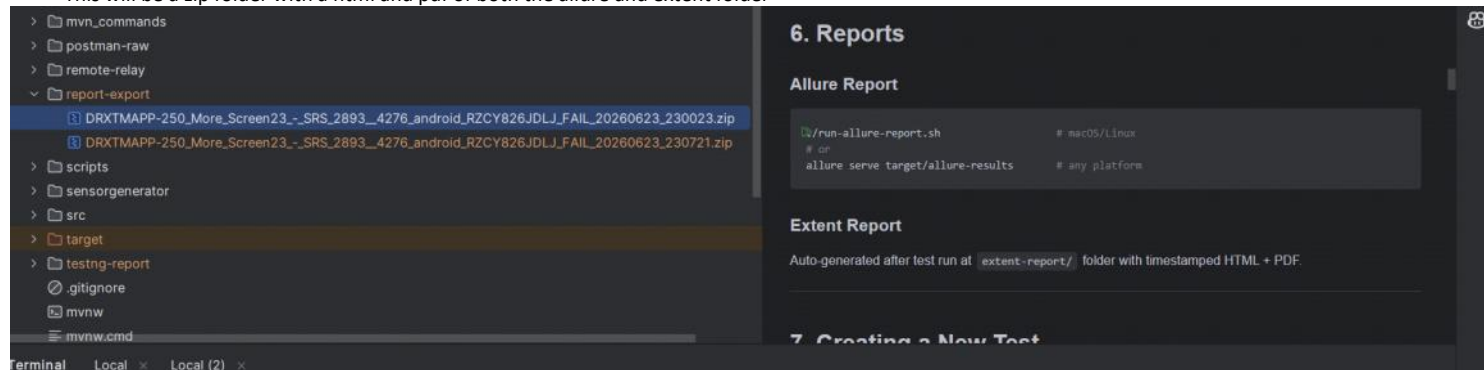


Errors in sanity.xml


Go to settings > build, execution, deployment > Maven > Importing
Make sure JDK Importer is the one installed on computer

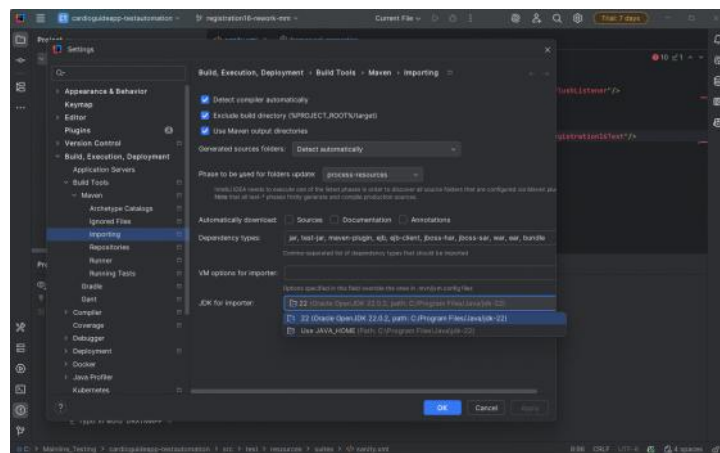


Go to File > Project Structure > Project
The SDK should be the one installed on computer

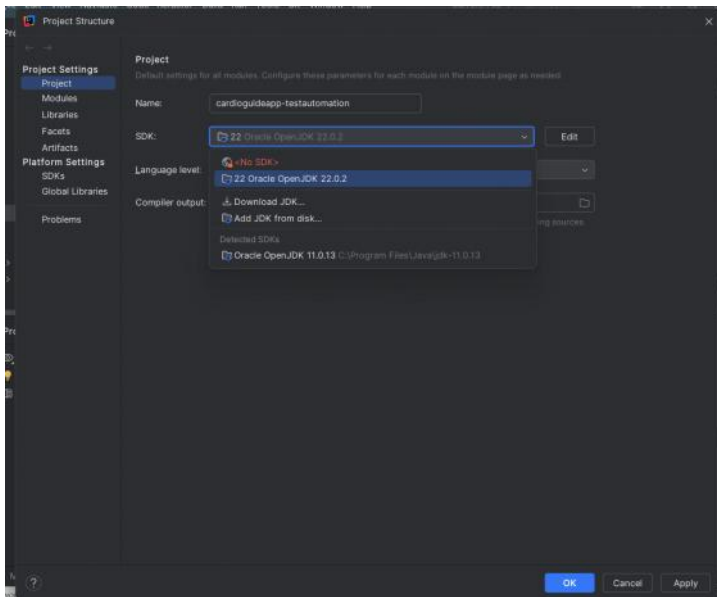

Go to Project Structure >Modules >+>Import Module
Import the pom.xml file from the cardioguide-testautomation folder

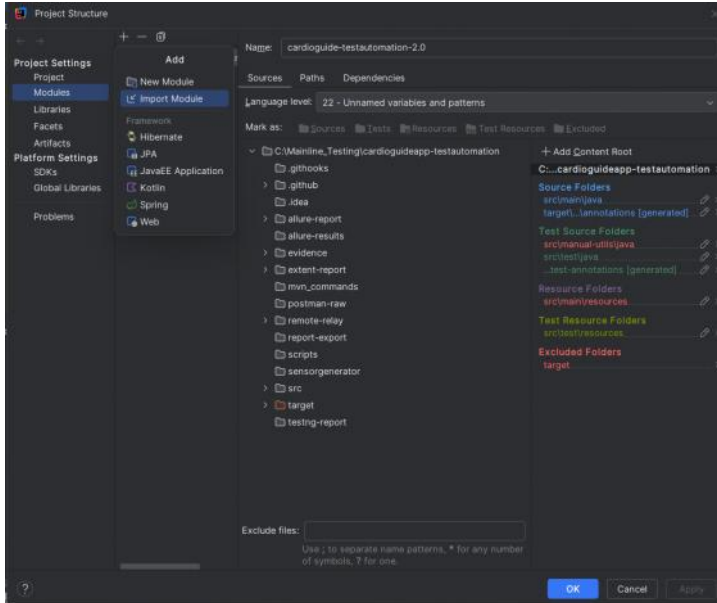

On the left side, open the "m". Press reload all maven projects

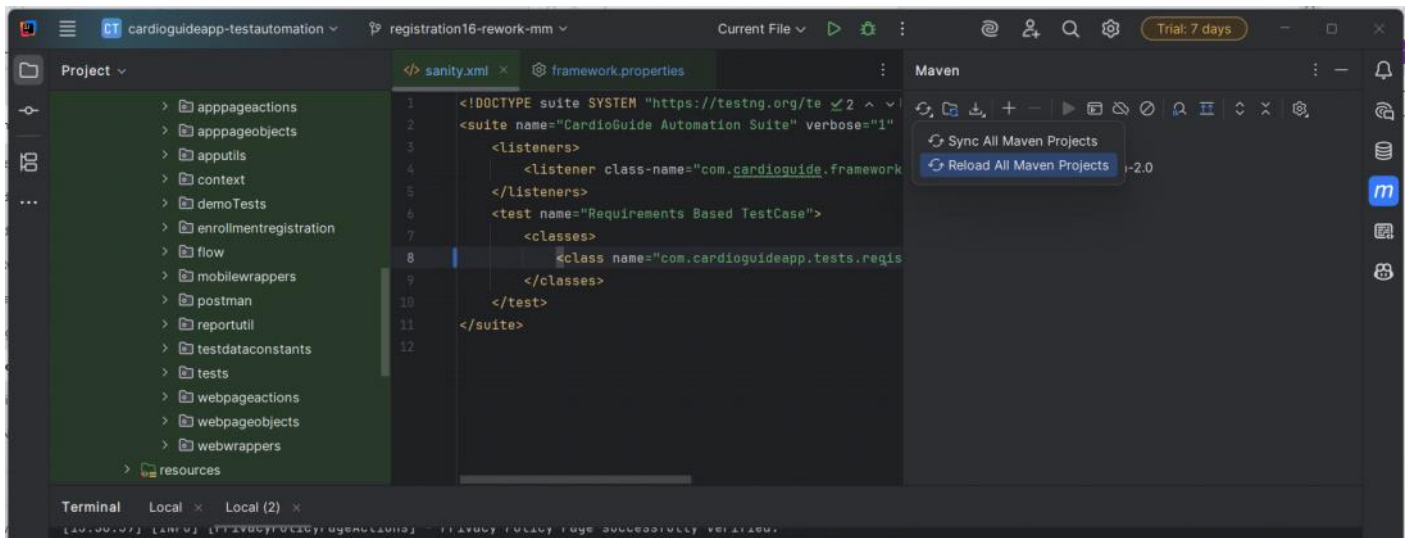

Demo Notes

Wednesday, June 24, 2026 8:11 AM

Tools that we install for automation

Tools should be installed for all users

Before execution

Pull the code from the branch

Ex, more screen 37

Git switch MoreScreen37Reworks

(rebase)

Git pull origin MoreScreen37Rework

git pull origin CardioGuide2.0-Web-API_Integrations

Stay on the branch and rebase it

Git stash

Git pull origin CardioGuide2.0-Web-API_Integrations

Close merge message tab

Mainline tester will do rebase after changes have been made

REBASE first to see if there are any merge conflicts

Keep a copy of the framework.properties

Git stash pop to recover changes

Execution.machine=windows

QA tester info

Download the report from report-export

Open the automation task ticket

- Go to the comments
- Attach the report

Clearly state what happened

Mainline testing completed

Execuated testcase on prisitine ____ setup
Pass or Faiol

Add comment to the PR also
Same as whats in the ticket

Copy the link from comment from automation ticket and past in PR

Framekwork: setting certain things

Macbook Automation Setup

Monday, June 29, 2026 12:43 PM

List IOS device udid: `xcrun xctrace list devices`

Make sure you use ipa when starting appium inspector

```
{  
  "platformName": "iOS",  
  "appium:platformVersion": "18.5",  
  "appium:deviceName": "iPhone SE",  
  "appium:automationName": "XCUITest",  
  "appium:udid": "00008110-001E5CD11EEA201E"  
}
```

When switching do different iPhones, make sure you trust the apple id in the VPN and Devices section of settings

Reworking test cases

Work with rahult o rework

Every weekday do some batch execution

Had to install appium driver install uiautomator2 and logger

Nvm path configuration

Mainline Testing

Tuesday, June 30, 2026 9:32 PM

			Reviewers	
	HKI10 new test auto, Onboarding flow, Xray cleanup, softFailure helper Car d oGui de2. 0- Web- Api _I ntegrations Jani, Nikhil - #361, last updated 1 hour ago			 ApprovedNeeds work  ApprovedNeeds work  ApprovedNeeds work
	Rework: DRXTMAPP-5611 Notification44 Car d oGui de2. 0- Web- Api _I ntegrations Kumar, Abhishek - #347, last updated 7 hours ago		Comments: 3	 ApprovedNeeds work  ApprovedNeeds work+3   
	Medication46 fail(stringLabelDuplicate) Car d oGui de2. 0- Web- Api _I ntegrations Kumar, Goutam - #341, last updated yesterday		Comments: 1	 ApprovedNeeds work  ApprovedNeeds work+3   
	ReworkNotification12 fails(Bug) Car d oGui de2. 0- Web- Api _I ntegrations Kumar, Goutam - #353, last updated yesterday		Comments: 2	 ApprovedNeeds work  ApprovedNeeds

				work+3   
	Notification9 1 LMR Fail(String Label Notification-header) Car d oGui de2 0- Web- Api_ ntegrations Kumar, Goutam - #360, last updated yesterday		Comment s: 1	 ApprovedNeeds work  ApprovedNeeds work+3   
	Get UserSync api integration Car d oGui de2 0- Web- Api_ ntegrations Chandran, Vikramraj - #350, last updated yesterday			 ApprovedNeeds work  ApprovedNeeds work  ApprovedNeeds work
	Notification17 4 - LMR Automation Car d oGui de2 0- Web- Api_ ntegrations Chandran, Vikramraj - #331, last updated yesterday		Comment s: 2	 ApprovedNeeds work  ApprovedNeeds work+3   
	DRXTMAPP-5699 Notification5 3 Car d oGui de2 0- Web- Api_ ntegrations Kumar, Abhishek - #335, last updated yesterday	Conflict	Comment s: 3	 ApprovedNeeds work  ApprovedNeeds

				work+3
	Notification17 3 - Windows exec issue fixed Car d oGui de2 0- Web- Api _I nt egr ations Chandran, Vikramraj - #346, last updated yesterday		Comments: 3	ApprovedNeeds work ApprovedNeeds work+3
	DRXDEVAPP-7566/More Screen25 - US Car d oGui de2 0- Web- Api _I nt egr ations Mohite, Pratik - #356, last updated 5 days ago	Conflict		ApprovedNeeds work ApprovedNeeds work ApprovedNeeds work
	Rework DRXTMAPP394 Registration10 Car d oGui de2 0- Web- Api _I nt egr ations O Priyanka, Priyanka - #334, last updated 6 days ago		Comments: 1	ApprovedNeeds work ApprovedNeeds work+2
	Priyanka/journal24 fix Car d oGui de2 0- Web- Api _I nt egr ations Patel, Priyanka - #348, last updated 6 days ago		Comments: 1	ApprovedNeeds work ApprovedNeeds work

				ApprovedNeeds work
	HKI9_1, HKI9_2 Peer execution rework Car d oGui de2 0- Web- Api_ Integrations Jani, Nikhil - #321, last updated 1 week ago		Comments: 4	 ApprovedNeeds work  ApprovedNeeds work  ApprovedNeeds work+2 •  • 
	DRXDEVAPP-9379/Reading29_2 Car d oGui de2 0- Web- Api_ Integrations Mohite, Pratik - #333, last updated on 18 Jun 2026	Conflict	Comments: 3	 ApprovedNeeds work  ApprovedNeeds work  ApprovedNeeds work
	Homescreen4 & Navigation4 (Pending Internal review)			

From <<https://bitbucket.sim.com/projects/MYC/repos/cardioguideapp-testautomation/pull-requests>>

Test Case Rework

Wednesday, July 15, 2026 12:44 PM

Signin3

-

Rework Notes

Monday, July 20, 2026 10:47 AM

Failing due to not pressing back button

- Registration9
- Signin26

Failed to click settings icon

- Nav13
- Journal9
- Registration31
- Registraion9

Client ID:

86A1D8829C494260B89DF497E5464EE1

Client Secret:

7377cf7ff35b39703a1021fdecfa976423736198928d06b9128f1ed820acc476

```
mvn clean test -DplatformName=Android -Ddevice.udid=f17437eb -  
Dapp=/Users/BODANSX/Downloads/29t1.apk -Dxray.enabled=true -Dxray.cloud.clientId=  
86A1D8829C494260B89DF497E5464EE1 -Dxray.cloud.clientSecret=  
7377cf7ff35b39703a1021fdecfa976423736198928d06b9128f1ed820acc476 -  
Djira.api.token=ATATT3xFfGF08lXeLMx7Vmf9U_LVX2zME1Dj5lfo1XI_zZTEth55FQAESMyErn9UMwb-  
HLUBAMU9y9ADcvc9JDz4CSK5dXThqtLIWid_  
1uphhm7lleAsvuWDMlexRoQAbw4F8Fy2wsm_aHXA5nn7YK_pOYYqojZvKXuA0brt3V53oHbkiZFvtE8=  
013D8E57 -Dxray.test.plan.key=DRXTMAPP-6225 -Dxray.execution.include.datetime=false -  
Dxray.execution.include.counter=false -Dxray.execution.prefix="LMR - Dry Run 3" -  
Dxray.project.key=DRXTMAPP
```


Cardio Guide App Notes

Monday, May 11, 2026 4:31 PM

Making patient

- Merlin
 - o Same as HF Cloud
 - Name, email, own phone #, sensor serial
 - o Sensor serial = 1st column of unclaimed from Sensor|CMHF Customer care backend
 - o Date of birth = 01-01-2000, Merlin is month, day year, but POSTMAN is day, month year
 - o For goals & treatment plan to be used, turn on the two buttons in CardioGuide HF App tab in Merlin
 - Activate patient using POSTMAN (environment v3)
 - o If patient is inactive, will show up in "Scheduled patient" tab and not the "Active patient" tab
 - After patient is registered in App, THEN dynamic treatment plan can be turned on in the Merlin portal
 - o Allows to call back patients and data (API)
 - "For Making Patient Active" folder
 - 5 files in total
 - ◆ Files 1-4 for making patient active for the first time
 - ◇ Send first file
 - ◇ Edit patient name, phone number and sensor serial tags in "Scripts>Pre-Request"
 - ◆ File 5 is like pillow, send readings to app after patient is active
 - ◇ Reading can be seen in "Activity" section of Merlin
 - ◇ If want to send reading to app, change in "Body"
 - ▶ Only need to change first 3 avg systolic, avg diastolic tags in main script
 - Make sure it is same as the long list of accepted numbers
 - ◆ 35 will pose "SUSPECT" reading as not in accepted numbers
 - ◆ SUSPECT means extremely high, won't even go app, outside of doctor determined range in Merlin account
 - Clinic is v3 and implanter is test
- Password: "Test!123" for email and patient account
- Can push messages for every 1, 3, and 7 high readings
 - o Such as Call clinic and take medication dosage

SAVE AFTER EVERY EDIT

Software

- Merlin = backend for patients
- CardioGuide App connects to pillow
 - o Test automation in Java
- Aquitter (ish) gives x paths for buttons which are used in test automation
 - o Assign x paths to string names
- POSTMAN
 - o Like API to do callbacks

5/13 Notes Patient Registration and Bug Verification

Thursday, May 14, 2026 9:05 AM

Official release build = 2.0.20

Should install app via google Play

- Past google play link in QR code generator and scan it
- Download it from Google play

Not updated (only 2.0.15 in Gplay),

- Copy apk from official folder
- Transfer apk file to mobile app and install by importing apk to files app and install on android phone
 - o Can scroll to confirm right version

Patient Registration

- First enroll patient in merlin portal
 - o 2 ways to register patient,
 - o Enroll patient button
 - DOB 01-01-1990
 - Treating Physician
 - Logging into Merlin using premade clinician
 - ◆ V3_test_last as physician
 - Fill out sensor Serial #
 - ◆ Choose unclaimed in backend (SWVT lot)
 - ◆ Enter in Merlin, pause for website to validate # with other patients
 - ◇ Website checks with other patients to check if other patients have same serial #, if not, box should be blue
 - o Patient is enrolled, but not active (scheduled patients tab)
 - Register patient in mobile app or POSTMAN

Register Patient w/ app

- Once button pressed, nav to Policy screen
- First opening screen is called startup screen
 - o Then got to EULA screen and Privacy screen
- Patient registration screen
 - o In merlin portal
 - Add DOB and serial # to app, press next
 - o Patient Infor
 - Fill up details
 - Should be same as testing, but can be different
 - Enter own phone #
 - PW Test!123
 - o Patient is registered to app, now login with email and PW

If try logging in with another patient, it won't allow as device can only be registered by 1 patient
Need to uninstall app and reinstall to register new patient

No policy screen when logging in, already created patient so skipped registration flow

Home Screen

Systematic BP

- Record BP for patient
- See same card on readings on global navigation panel
 - o 3 Panels, day, week, month
 - Day is today
 - Week is entire week
 - Month is entire month, data not recorded for entire month -- value is present

DTP in Merlin

- Add medications in merlin, after adding goal PA Pressure
- Turn on Medication managed in App in Merlin
- Use save and add another button to add more than 1 meds

In app, Medication adjustment notification if medication plan is changed in Merlin

IOS has extra button in More options (Update on Apple Health)

Bugs in Jira

- Filters, search work filters
- 3 spaces
 - o TestManagementApp = test cases
 - o SRS = requirements for project
 - o DevApp = dev builds (bug for app)
 - o DevPortal = Portal team (bugs related to portal)

App bugs (DevApp)

- Type = Bugs
- Status = Verify
 - o Bug has been already logged, dev team has checked bug, and now it needs to be verified
- Verify, view Workflow
 - o When people are working, status is reading and then in progress
 - o Once fixed, status is verify
 - o Once verified, assigned to Danny and we should be retesting the bug
- Labels, Not equal to Ready to Close

Example Bug

- Not been retested yet
- Revisit all the test description and read it all
- When understand bug, retest it in new build
- Perform steps that were mentioned to reproduce test
- Go to actual results and compare to expected results, bug is logged
- Make sure expected results are met
- Write up comment
 - o Environment details, what the build is, and what phone it is tested on
- Status is resolved or not
- Write expected and actual results
- Add evidence in comment s
- Once comment is added, add ready to close button

If bug is still unresolved

- Put up comment as status being unresolved
- Add second comment and reassign bug to original person, telling them bug is still not resolved

- Put it back to in progress

If want to log up bug,

- Log bug in Dev App space
- Press create
- Work type as Bug
- Add summary of what the bug is about
- Have pop up, this bug has been created
 - o Press it and land directly on ticket
 - o Enter description

Time of acquisition:2025-01-01T00:00:00Z

1. APK (Share point folder) : [GLB-HF-DRx-myCardioMEMS App - QA - All Documents](#)
2. Environment V3 <Merlin Portal> : <https://v3-hfcloud.products.abbott/apps/clinician/notifications>
Username : v3_test_pat
Pd : Test#1234
3. Backend V3 : [Sensors | CMHF Customer Care](#)
4. Bug Unresolved (Ex) : [\[DRXDEVAPP-7006\] Env V3-Notification Badge on CardioGuide App Icon Disappears After App Launch With...](#)
5. Bug Resolved (Ex): [\[DRXDEVAPP-7667\] Env V3 : Error message appears when tapping "Send a new code" with the verificatio...](#)

Defect Status: Resolved

The MFA Phone error message is displayed after 1 minute 40 seconds after user fails to complete Call Me MFA, thus satisfying the Expected result.

Details

Environment: Q3

Build Details (Android): Version: 2.0.28

Clinic Details: Username: q3_svvq3, Password: Test!1234

User Details: Sam15 RTR, Username: sam15@gmail.com, Test!123

Device Details: Samsung Galaxy A56, Android 16

Evidence:

5/14 More Options and Bug Testing

Tuesday, May 19, 2026 3:34 PM

More Options Tab

- Logged into App and on Home screen
- More Menu tab
 - o Profile option
 - Full name in Merlin Portal
 - Email Address
 - Phone number
 - Can edit only full name and phone number
- As of 2.0.22, can change country code
 - o If profile page updated
 - o Go to Merlin portal
 - View registered application user
 - Only find updates here, not on the main portal
 - o Support
 - FAQs, Contact, About, User Manual
 - FAQs have more info
 - Contact
 - Clinic contact details from merlin portal
 - Add primary lab,
 - o Just for user not in merlin
 - o About
 - Privacy Policy and EULA (info for user)
 - Has software version
 - GTIN - related to phone
 - App model number
 - o User manual
 - Guidelines, how to find product info
 - o Education
 - About HF
 - Videos about HF and CardioMEMS system
 - Learning (more info on everything)
 - o Weight
 - Recording weight for patient
 - 3 tables: day, week, month
 - Each tab has record weight button
 - For other days, use POSTMAN to enter weight
 - User calendar icon to see previous weeks and months weights trends
 - o Settings
 - Change password
 - For IOS, has update to apple health button there too
 - Whatever data in apple health app, it would be reflected here
 - Change password
 - Can change password, and it redirects you to login again

Uninstall on Android

- Settings, app info then uninstall

Loggin Bug in Jira

- Observation issues while testing (devApp space)
 - o Goes to dev team
 - o Create
 - Select space, devApp
 - Work type = bug
 - Status= open
 - Summary = what is the us
 - Env V3 - [what bug is]
 - No need for components
 - If know where the component is, can add that
 - Description = min thing
 - Environment in which the test was done
 - Build details
 - Test details
 - Clinic details (for merlin portal)
 - User details (what user was used to login)
 - Device details (android, IOS and update, and phone model)
 - Steps to reproduce
 - ◆ While doing what steps issue was identified
 - Actual result = what you observed
 - Expected result = what you should have seen
 - Evidence (screenshots or video)
 - ◆ Just result you observed, not every step
 - ◆ Buttin to add image and photo
 - ◆ Attach In description
 - ◆ When single attachment is added, then can add on plus signa nd directly add it
 - ◆ No attachment option if you just create a bug
 - ◆ Only added when added a screenshot/video from description
 - Assign is automatic, if known can put it
 - Severity option too, low high med
 - Labels
 - ◆ CG2.0 Defects
 - Affected version = Q#APP - 2.0.21
- Nees to create bug
- Space = drx_devapp
- Test cases is drs_testmanagementapp
 - o DRXTMAPP
- Bug retesting
 - o Reproduce it with steps mentioned
 - o In comment section
 - o Verify status
 - o Reproduce using states with new build with the new build
 - o In comment section
 - IF NTO SISSUES, AND IT WAS RESIVED, DO RESOLVED STATU
 - o TESTED ON, ON WHICH UPLOID AND TDATE UPDATE
 - Env tested
 - Device tested on
 - Merlin portal version as well
 - Help about, find merlin version

- Clinic details , and full name fir the patient, gmail, and password
- Add evidence that it is fixed ow
- Assign it to ourself, then add to comment and assign to danny again
- Label
 - o Ready to close label
 - o Danny knows bug is fixed

DOORS

- Requirements in CardioMEMS > Next Generation Products > Software
 - o App110 = CardioGuide App
 - Requirements
 - ☐ Have SRS and GUI Specification
 - ☐ SRS are the requirements
 - ☐ Test cases are created based on requirements
 - ☐ DRxRTR is a valid req for anyone NOT CLOUD
 - ☐ Sometimes req team modify reqs, then the description or ID would be modified
 - GUI Spec
 - ☐ Create test cases or update test cases based on GUI Spec
 - ☐ Name of screen
 - ☐ In test case step, provide name of screen
 - ☐ Be more specific based on GUI specific
 - ☐ Important for updating or creating test case
 - String export doc
 - ☐ Messages and strings on application

5/15 Notes POSTMAN

Friday, May 15, 2026 9:31 AM

POSTMAN

- Workspace = CG App Verification team
- Many Collections
 - o AK - IAPSS collection for V3
 - For making Patient Active folder
 - Register collection registers patient
 - Can directly log into app, no need for registering in the app
 - Copy sensor serial number
 - Run Register HES Copy
 - ◆ V3 environment
 - ◆ Copy same as Merlin portal
 - ◆ Just change environment and hit send
 - Create Patient HES Copy
 - ◆ Copy Sensor Serial #
 - ◆ Replace # in pre-request > scripts tab
 - ◆ Fill in patient details
 - ◇ Date = Year Month Date, in Merlin, its month day year
 - ◇ Change phone #
 - Sensor Calibration
 - ◆ Baseline code in Merlin
 - ◆ When run 3rd file, fill up baseline code
 - ◆ Pre-request, fill out sensor serial #
 - ◆ Send
 - Initial HES reading
 - ◆ Patient status still inactive
 - ◆ Change same patient details and serial #
 - ◆ Send
- In Merlin, don't need to do hard refresh, just tap readings and it should work
- IRL, postman not real
 - o User postman to ease testing
 - Patient still inactive, use 5th file to make active
- Now patient is active after sending 5th file

DTP

- If patient is having high reading or low reading
- Clinician can send medicine
- Ex
 - o Goal = 17
 - o Low = 13 systolic
 - o High = 21 diastolic
 - o User = 15, normal reading
 - o User = 24, high reading
 - o User = 11, low reading
- Readings are for 1 week period
- Save interventions in merlin
- Then can turn on dynamic treatment plan
 - o Clinician is going to use Merlin

- Clinician can turn on DTP for user
- Home Reading POSTMAN file
 - Go to Body
 - Change data section
 - Avg systolic = 22
 - Avg mean = 21
 - Avg diastolic = 20
 - Patient is in Active Patients tab in Merlin now
- Merlin time zone = California time
 - If clinic is in different time zone
 - Clinic = 5 AM
 - Triggered reading at 6:30 AM
 - Reading will not get posted
 - Time postings must be before clinic time
- Cannot revert back once send, once reading file is generated cannot revert back
- Suspect = high reading
 - When hit suspect, not trigger DTP
 - If trigger any other readings around it, always trigger suspect
 - Use another patient
- 2 DTP cards are triggered because of the 7 out of range readings

Installing App

- Extract folder
- Folders are environments
- Testing on QA (v3)
- Connect phone and copy apk folder and past it in the phone from computer
- On phone just tap it and it should install

Dashed on the medication screen

- Some have dashes, some don't
- Company from the medicine

Helpful Notes

- Always create patient first
- And enter it in pre request > scripts
- Send first reading to make patient active
- Body script is where edit info abt patient reading first
- Make the time a long time ago, just make a baseline
- Enter
 - Avg systolic
 - Avg mean
 - Avg diastolic
- File 3 is sensor calibration
 - Create patient
 - Calibration code
 - Then first one reading to make patient active
- Suspect = if system thinks data doesn't make sense
- Out of range = sensical value has been read, but it is out of range

5/19 Notes Final Features

Tuesday, May 19, 2026 9:31 AM

Non-Drx Interventions

- Go to Merlin, select specific patient
- More actions button, add interventions
 - o There are 6 options (radio buttons)
 - Reading reminder
 - Lab visit (lab requires blood test, schedule lab visit within next 7 days)
 - Treatment plan reminder (due to BP changes)
 - Can send custom message
 - Notification badge says "New clinician message"
 - o These are all non drx interventions
 - Click on one, press send
 - Notification should be sent to the app
 - Can be marked as read
 - ◆ See in complete tab in app
 - ◇ Complete screen shows read cards
 - See different texts on different notification tabs
- Daily medications sent to app are also considered non drx intervention

Drx intervention are only the DTP notifications (for high and low readings)

- Can titrate daily medication (hold = stop medicine)
 - o Select daily medication, select quantity and frequency
- Can add message to call clinic or custom message
 - o These would be considered drx intervention now
 - Only more actions messages would be considered non drx
 - Daily medications will be non drx

Out of Home Elevation

- Out of home elevation is change in atmospheric pressure
 - o This would not send an intervention
 - o Would send a notification to the app, but no intervention
- Go to postman
- Go to Send home readings (5th file)
- In body tags, AtmosphericPressure
 - o By default, value is 753
 - o Out of home gets generated when 35 is added to it (788+)
 - o Irrespective of if the reading is within the valid range or not, out of home notification will be generated
- In previous builds
 - o Label of notification badge is "Out of home elevation"
 - o Now it is "Away from home elevation"
 - o Even text is different
 - o Later on no badge will be on the home screen, everything will be on the notification screen
- Go to Merlin portal, go to activity tab
- Go to readings
- Out of home elevation status is there

7 Out of Range Readings

- Set optimal ranges
- Interventions for 7 high readings (3 total things)
 - o Adding a new medicine
 - o Asking to hold daily medicine
 - o Adding message to call clinics
- DTP plan must be on
- Go to postman
 - o Refresh postman every time you change patients
 - Sometimes POSTMAN overwrites, refresh before working on new patient
- Use 5th file to send high readings
 - o Enter patient details in scripts (pre-req)
- Log into app
- See week (May 13-May 19)
- Trigger 1 low or 1 high
- Create intervention in merlin
- Create readings for May 18, 17, or any of them
 - o No specific dates needed
- For every test, create new patient (for every new intervention case)
 - Create 6 patients
 - Test it individually for each intervention
 - o Once you trigger a intervention
 - o Cannot trigger another one for the same patient on the same day
 - o Triggering right now for 7 high
 - o Tomorrow need to trigger other readings, then you can get another intervention
 - o One day, only 1 intervention can be triggered
 - o Tomorrow you can trigger something different for this particular patient
 - o Cannot just keep triggering interventions
- Keep same values, change the date
 - o Have to be in a row
 - o Mentioned 3 interventions
 - Sends 3 separate notification badges
 - Text is different for each drx notification
- Ask user to confirm that they are currently taking old medication before accepting that they stop medicine
 - o Now earlier medication is not in the medicine list tab
- In Merlin portal
 - o Can see intervention cards for each DTP intervention
 - 3 cards, 1 for each intervention

Helpful Notes

- Don't create interventions for all 6 cases
- Only create one for the particular one you are testing
- If have 1 low and 3 low interventions
 - o Send 1 low readings few days ago
 - No intervention was triggered
 - Even 3 low readings, no intervention was triggered
 - This is because the past intervention is generated for TODAY
 - ◆ Notification was sent TODAY, even if reading was taken in past
- 1 intervention for 1 day
 - o 24 hour waiting period

Sending Weight and BP


- Get B2C token API file, just send it

- Got to Login page
 - o Send to get patient ID
 - Edit email address for patient
- Post vital(weight, fluid) page to send weight
- Post BP Reading page to send BP reading
 - o Mention patient ID
 - o Mention specific weight
 - o takenAtUtc tag
 - Epoch converter: <https://www.epochconverter.com/>
 - Mention date to get time stamp and a Unix value and paste into this tag
 - This is a millisecond value

Testing Notes

Monday, May 18, 2026 11:46 AM

- Nothing in Anabolic Steroid for both Q3 and V3

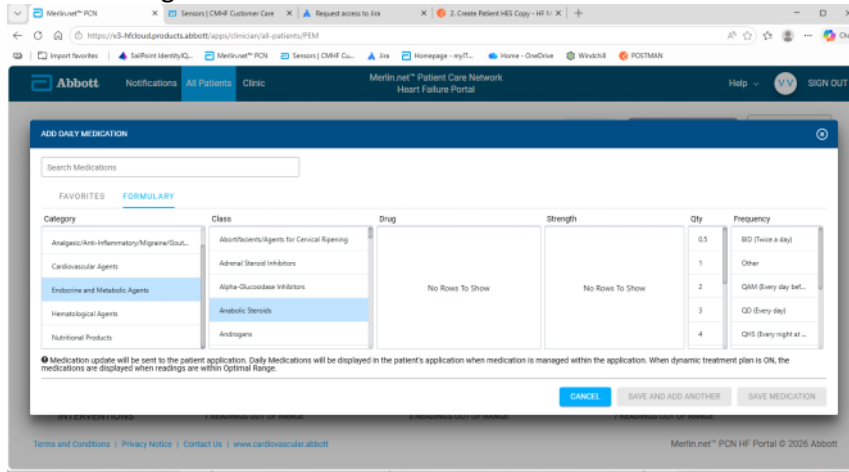

Once interventions are added, can delete them by pressing individual intervention when blue selected
Can remove them by pressing cancel if not saved for first time

Workspace = CG App Verification Team, not HF Mobile apps
AK - IAPI "For Making Patient Active"

In POSTAM Time of Acquisition tab, Year, Month, Day

PA Mean is the default indicator for high or low readings

Mean = $(1/3)(\text{Systolic}) + (2/3)(\text{Diastolic})$
[Estimation of Mean Pulmonary Artery Pressure - CHEST](#)

Must refresh merlin to see reading in app

For sending BP and Weight from POSTMA < Patient id is patientLoginSrcId tag

If the date on the phone is modified (not set to automatic) app does not login

Weight and BP readings will be updated in app once deleted and reinstalled app

If PA reading sent to app is disabled in Merlin, the CardioMEMS reading tag is greyed out AFTER reading is taken

Both are highlighted blue before reading is taken

In the readings tab, the week and month are just checks and x

For goals and treatment, optimization phase is after sensor is implanted
Clinician is trying to determine optimal metrics

Maintenance is after optimization phase

App very often says "Internet issue", assuming due to Asus 5G wifi
Sometimes multiple times on top of each other

This is a badge made from the developers, not sent by the android phone
**Don't think the app is sending phone notification, just a badge they created
Assuming by the style of badges on app vs android
Even in requirements page, have only 1 style version of badges

Interventions are not displayed on home screen, but only show up on notification section

Can only edit the HF primary lab phone number and information, not the clinic or tech support number

The call buttons do not call directly, only populate the phone number in the phone call app
User still needs to press the call button

Once change the time frame of the month, and press the day or week button, still only shows present time frame.

Readings data is kept for only past 6 months

When send atmospheric pressure reading less than 718 app also sends away from home elevation badge
Says, change in elevation still

Maximum number of times for MFA is 3 tries
Get notification

When DTP is turned off (option to keep meds but no interventions, past interventions are still shown on app

Must have the medication manager in app turned on for DTP to be turned on
When meds manager off, the page is blank, has title and text of what page is for

If editing strength of existing medicine, makes you confirm taking medicine before mark as read
In medication adjustment, red title says "Increase" or "decrease"

When registering patient, sensor SN and Bday must match patient. Name can be diff from Merlin
But in POSTMAN, name must be same as in app

The app keeps user logged in for 24 hrs, or 90 days with "Keep me logged in"

2.0.21

For Q3. Merlin portal and Sensors backend portal is different

Scroll bar and warning is present in medication page

In POSTMAN, use "Q3 Working" folder to make patient active

2.0.22

Must use Q3 environment from now on

Are able to change the country code for phone number

For Q3, PA diastolic is used as main metric

In POSTMAN Q3 working folder, file 5, the time of acquisition is 1 day off. It is one day ahead, if reading

for June 1 is entered, in Merlin, shows as May 31

**If both medication adjustment and message is used as intervention, the message can only be seen in the notification tab

Only medication adjustment is shown on the home screen

"High Priority" notification says Action Requested

If logged in first time on fresh app, have to accept EULA and privacy policy. When logged out, app is still registered to that user. When logging back in, do not need to accept again

Install APK from computer to phone:

- Once adb is installed, connect phone to computer
- Cd into the folder with the apk folder
 - o Run adb install myCardioMEMS-q3-2.0.25-157623.apk

Observations

Thursday, May 21, 2026 9:56 AM

Would change icons on the More > Support screen

- Make FAQ question mark and user manual as book or document icon

Keep intervention notification on home screen

In Notifications > Complete, more organization. After many notifications, can't remove them so could get very long (no scroll bar present)

Even if there is a less pressure sent, it still says "elevation", could be confusing?

Change in "altitude"?

Forgot email button only has "ok" button, there is no text telling user what to do

No FORGOT USER ID INSTRUCTIONS COMPONENT

Can only refresh on the home screen, not on the readings tab

Away from home elevation notification could also interfere with intervention notifications

Videos don't work on both 2.0.20 or 2.0.21 builds

Or 2.0.22, but at least there is close button

How would user know if reading was a suspect reading?

For writing test cases, how to address the different starting screens and pop ups for Android and IOS

~~6/5/2026~~


~~Noticing that after logging out and signing in, does not send OTP on the first time on both IOS and android.~~

Why is it called an "Elevated Pressure Message Card Panel" for the Treatment Plan reminder?

The ELEVATED PRESSURE MESSAGE CARD PANEL consists of:
TREATMENT PLAN BANNER COMPONENT
NOTIFICATION DATE LABEL COMPONENT (See Shared Elements: Components)
PRESSURE CHANGE CONTENT LABEL COMPONENT
PRESSURE CHANGE INSTRUCTIONS LABEL COMPONENT
MARK AS READ ACTION TRIGGER COMPONENT (See Shared Elements: Components)

8.12.1.17.2 ELEVATED PRESSURE MESSAGE CARD PANEL (Layout)

<Picture> ELEVATED PRESSURE MESSAGE CARD PANEL (default: active, example)

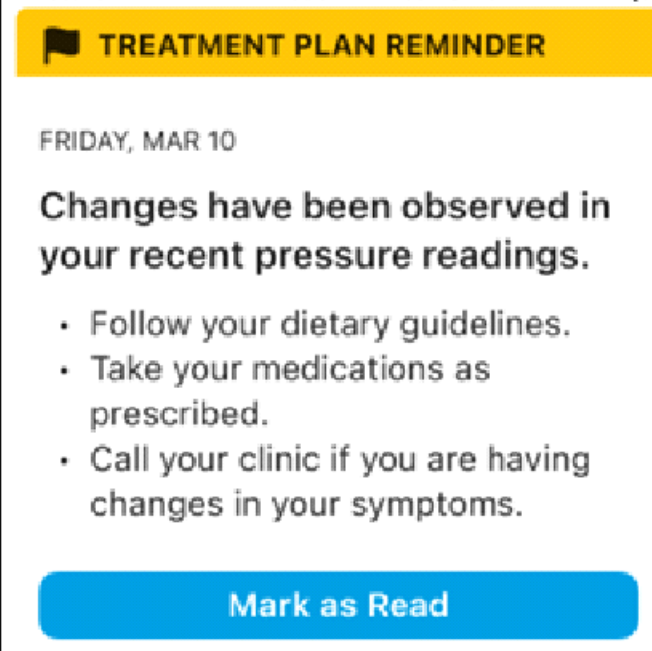

Can approve reading from Merlin to approve suspect reading
Shows up as approved on Sensors backend
Does not show up on app

2.0.26 Ad Hoc Testing

Change password issue

Requirements Notes

Wednesday, May 20, 2026 1:00 PM

DOORS Modular structure

- Strings -> GUI Spec -> Software Requirements
 - o GUI = graphical UI
 - Refers to the GUI of the DOORS program?
 - Or the app?
 - Different from UI requirement
 - UI = what system allows user to do
 - GUI = how it looks
 - Describes UI (wireframe and descriptions of all widgets and buttons)
 - And has details for messages shown to user
 - Screen = full page/display
 - Panel = collection of GUI components (taskbar)
 - Component = button, label, etc
 - Labeled in ALL CAPS
 - GUI screen View
 - Display = what is being shown to user
 - ◆ Lists everything visible on screen
 - ◇ Bulleted list ig
 - Layout = complete wireframe (picture ig)
 - Element = components
 - GUI Panel
 - Same as GUI screen view structure
 - o Strings
 - List of all strings
 - MVP = English translation

Customer Case Study

- Were confused about the medication list page
 - o Added warning sign and a scroll bar
 - NOTE: As of 2.0.20, I don't see warning sign
 - There is scroll button, but only when the notification is first sent to the user
 - ◆ No scroll button in complete tab or in medicine tab
 - o Confused about the confirmation of medications (before they stop them)
 - Used to have "Step 1", "step 2" labels on confirmation and acknowledgement pages, I do not see them now'
- Were confused about the out of home elevation
 - o Were not reading all of the text
 - Shortened it for new releases
 - Made message clearer (general wording)
 - o Were confused on what it meant
 - Changed from "out of home" to "away from home"
 - o The CardioMEMS readings tab just says "Remember to take reading", not "Reading is under review"

{0} = CardioGuide HF™

Requirements Document

- There is Required and Recommended software updates

- User does not have option to skip required update
- In document, setup goes name -> display -> UI info -> component list -> picture
- Keep me signed in became keep me logged in
 - Text is shorter too. There is no text telling
 - user that this is not recommended option
- Line graph is more sharp, not curvy

In DOORS

- CardioMEMS > Next Generation Products > Software > Configuration
 - File for all Strings
- CardioMEMS > Next Generation Products > Software > Requirements > GUI Spec file
 - Panels, components for app
- CardioMEMS > Next Generation Products > Software > Requirements > SRS file
 - Requirements for App

Building Notes

Friday, May 22, 2026 4:46 PM

App uses FCM for background notifications (when phone is locked)
Flutter_local_notification package for local notifications

Azure application used for analytics

Custom package with headless CMS for content management
- CMS = publish content easily

Flutter = framework, Dart = language

MVVM

- All user interaction occurs within View
 - o In charge of detecting user input
 - Forwards to ViewModel via Data Binding
 - Data Binding = connects UI to underlying data
- ViewModel implements commands that are used by View to bind data to
 - o Holds states of the UI
- Model is business logic and app data
 - o Business logic = rules/algorithms that controls data


POSTMAN CLI

Monday, June 8, 2026 11:10 AM

postman login --with-api-key PMAK-6a23393a28180a0001a82455-f61161159a59a6b975f93aca763d5554e4

| | |
|------------------------------|---|
| Postman Collect UID | 52864576-d9174d04-011e-4d5f-aca9-0001b00910e6 |
| Postman CLI environment ID | 31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85 |
| Register HES ID | 52864576-d3c5eacc-6783-4863-b76b-5f9a7f24b244 |
| Create Patient HES ID | 52864576-50086e72-309a-421d-98da-39a1e244a5e7 |
| Create Sensor Calibration ID | 52864576-ebb454ce-2f57-41cd-950a-8a5253994009 |
| Create Hospital Reading ID | 52864576-06b12575-0c09-4f79-9fc6-ce543bf81b70 |
| Create Home Reading ID | 52864576-47436003-a58d-4abf-8276-701281a3ceed |
| Br_hes_uuid | 33fab413-3f1b-4524-bf70-a03407985964 |
| Hes_token | zpVrV6XuAdPRmYxx2n2TFMsGMYt73wrr |

**** Br_hes_uuid and Hes_token could be different double check in POSTMAN browser**

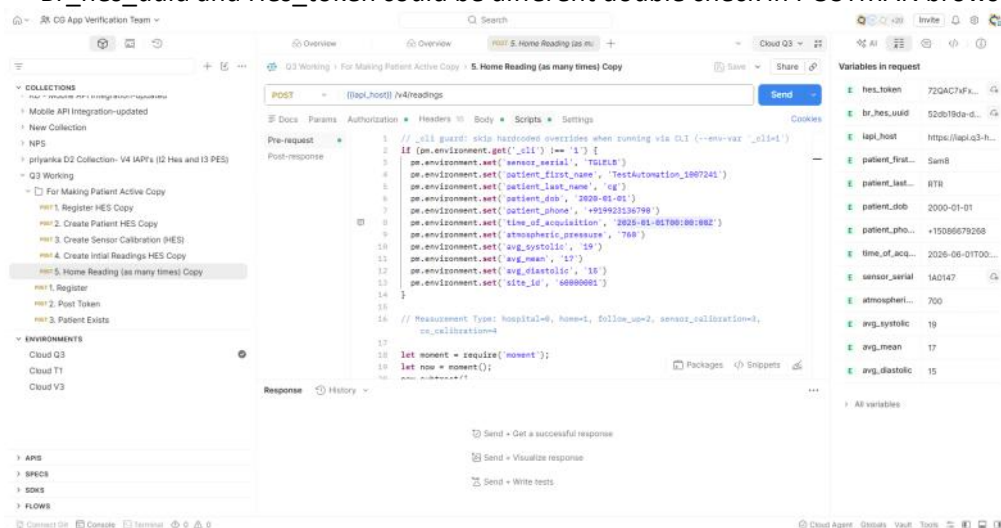

"Register HES"

postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "Register HES ID" --timeout-request 30000 --verbose

"Create Patient HES"

postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "CreatePatientHES ID" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "patient_first_name=EDIT_FIRST" --env-var "patient_last_name=EDIT_LAST" --env-var "patient_dob=EDIT_DOB" --env-var "patient_phone=EDIT_PHONE" --verbose

"Sensor Calibration"

postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "SensorCalibration ID" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "calibration_date=EDIT_CALDATE" --verbose

"Hospital Reading"

```
postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "HospitalReading ID" --
timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var
"hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var
"patient_first_name=EDIT_FIRST" --env-var "patient_last_name=EDIT_LAST" --env-var "patient_dob=EDIT_DOB" --
env-var "patient_phone=EDIT_PHONE" --env-var "time_of_acquisition=EDIT_CALDATE" --env-var
"atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --
verbose
```

"Home Reading"

```
postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "HomeReading ID" --timeout-
request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var
"hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var
"patient_first_name=EDIT_FIRST" --env-var "patient_last_name=EDIT_LAST" --env-var "patient_dob=EDIT_DOB" --
env-var "patient_phone=EDIT_PHONE" --env-var "time_of_acquisition=EDIT_CALDATE" --env-var
"atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --
verbose
```

```
mvn clean test -DplatformName=Android -Ddevice.udid=RZCX721WSBZ -Dapp=C:\Users\rajanax7
\Downloads\myCardioMEMS-q3-2.0.26-159262.apk -Dxray.enabled=true -Dxray.cloud.clientId=
57138CA16A6C40A8B85DDAC748FA2B72 -
Dxray.cloud.clientSecret=d271a2feaa3b03ee3b20e85e4340ecf0a6853d128cb7f0094e119fa15d78
1982 -Djira.api.token=ATATT3xFfGF0ffl81a1N9ZS-3McfTrrm4XsUXWotUafUrzd36c-
Rff66vtKkV7UK3_KxcIP_axXxlZtlc2J0pDseK4yJtkK7T1-
pXCdVAYVYI1r7RmG1S0HzmDwKFZe6RFEG4RA_hKBd_LLZN7Bm3JpbS4bVKpBm7ED9nllK
SLGICBgli1M6nIU=1F603FF4 -Dxray.test.plan.key=DRXTMAPP-5648 -
Dxray.execution.include.datetime=false -Dxray.execution.include.counter=false -
Dxray.execution.prefix="Trail Execution" -Dxray.project.key=DRXTMAPP
```

```
ATATT3xFfGF0KtzvH4y3VrrAU1MM8YINFAeP65qTU_st_3zy3B-
y5gtsmS48ny1BoIhoOL_VQOSss_VNcfVEiq6ssHyXy3KcV1pDqBLkJJ1V5bpCAD-3YguGmY6F9iWzsPUXb6xK
kr0FuxGsMs9dJk1bsPug9M8A3NIDlq-ikED0_fGoVVM1HQY=4F1D0395
```

May 31 2027 token

```
postman collection run "31616098-93f098bb-9a42-4cd5-8cf9-a6ffaddb51bf" -e "26913137-bb5246aa-
bec8-4fa5-a71d-765ea978c0bb" -i "52864576-f84a31f6-acca-4afc-bed3-03809766b495" -i
"52864576-3d04b979-13fa-4da9-ba7f-b9600896318e" --timeout-request 30000 --env-var "_cli=1" --env-
var "email=sam36@gmail.com" --env-var "password=Test!123" --env-var "user_group=user" --verbose
```

```
postman collection run '52864576-d9174d04-011e-4d5f-aca9-0001b00910e6' -e
'31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85' -i '52864576-f84a31f6-acca-4afc-
bed3-03809766b495' -i '52864576-3d04b979-13fa-4da9-ba7f-b9600896318e' --timeout-request
```

```
30000 --env-var '_cli=1' --env-var 'email=sam35@gmail.com' --env-var 'password=Test!123' --env-var  
'user_group=user' --verbose
```


POSTMAN CLI Filled In

Wednesday, June 10, 2026 11:38 AM

postman login --with-api-key PMAK-6a23393a28180a0001a82455-f61161159a59a6b975f93aca763d5554e4
Prasuna api key PMAK-6a56a2beb6a9b3000120face-ad762148f31fce6baa8398972fcc6d97a7

| | |
|------------------------------|---|
| Postman Collect UID | 52864576-d9174d04-011e-4d5f-aca9-0001b00910e6 |
| Postman CLI environment ID | 31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85 |
| Register HES ID | 52864576-d3c5eacc-6783-4863-b76b-5f9a7f24b244 |
| Create Patient HES ID | 52864576-50086e72-309a-421d-98da-39a1e244a5e7 |
| Create Sensor Calibration ID | 52864576-ebb454ce-2f57-41cd-950a-8a5253994009 |
| Create Hospital Reading ID | 52864576-06b12575-0c09-4f79-9fc6-ce543bf81b70 |
| Create Home Reading ID | 52864576-47436003-a58d-4abf-8276-701281a3ceed |
| Br_hes_uuid | 52db19da-d402-4adf-8a83-fc9f71e31a2b |
| Hes_token | 72QAC7xFcxsWs8rncWZJR3qxdWQzwr |

**** Br_hes_uuid and Hes_token could be different double check in POSTMAN browser**

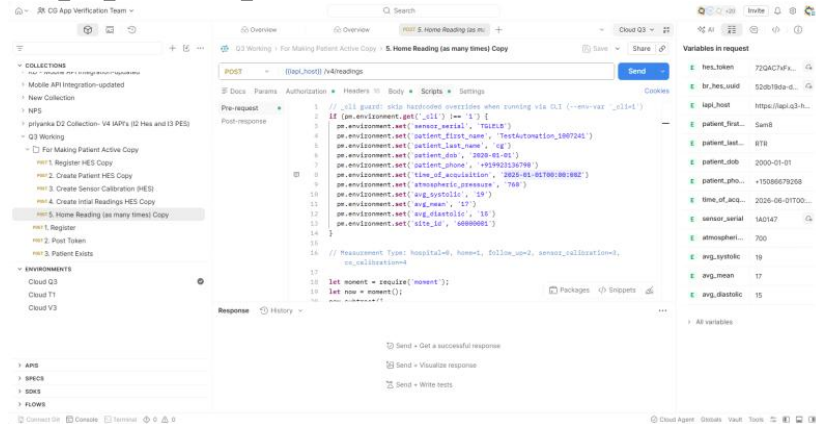

"Register HES"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-d3c5eacc-6783-4863-b76b-5f9a7f24b244" --timeout-request 30000 --verbose

"Create Patient HES"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-50086e72-309a-421d-98da-39a1e244a5e7" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "patient_first_name=EDIT_FIRST" --env-var "patient_last_name=RTR" --env-var "patient_dob=2000-01-01" --env-var "patient_phone=+15086679268" --verbose

"Sensor Calibration"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-ebb454ce-2f57-41cd-950a-8a5253994009" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "calibration_date=2025-01-01T00:00:00Z" --verbose

"Hospital Reading"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-06b12575-0c09-4f79-9fc6-ce543bf81b70" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "patient_first_name=EDIT_FIRST" --env-var "patient_last_name=RTR" --env-var "patient_dob=2000-01-01" --env-var "patient_phone=+15086679268" --env-var "time_of_acquisition=2025-01-01T00:00:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --verbose

"Home Reading"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-47436003-a58d-4abf-8276-701281a3ceed" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "patient_first_name=EDIT_FIRST" --env-var "patient_last_name=RTR" --env-var "patient_dob=2000-01-01" --env-var "patient_phone=+15086679268" --env-var "time_of_acquisition=2026-01-01T00:00:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --

verbose

30 days ago is the 5th july

```
postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-47436003-a58d-4abf-8276-701281a3ceed" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=7b0eb185-945e-4666-a119-65c7248d4c15" --env-var "hes_token=QqwXsqg2piyPbRHEYrmvNJV3gZtqpwr" --env-var "sensor_serial=1A3941" --env-var "patient_first_name=Sam47" --env-var "patient_last_name=RTR" --env-var "patient_dob=2020-01-01" --env-var "patient_phone="+15086679268" --env-var "time_of_acquisition=2026-07-05T10:20:00-06:00" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=16" --verbose
```

```
postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-86b8eeaa-9e22-47f2-a54f-e09cdbf69975" --timeout-request 30000 --env-var "_cli=1" --env-var "sensor_serial=TM4MGU" --env-var "patient_dob=2020-01-01" --env-var "patient_name=CGPatient Pat05175333" --env-var "email=cgta514552@gmail.com" --env-var "password=Test!123" --env-var "phone="+13399704117" --verbose
```

```
postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-86b8eeaa-9e22-47f2-a54f-e09cdbf69975" --timeout-request 30000 --env-var "_cli=1" --env-var "sensor_serial=1A4242" --env-var "patient_dob=2020-01-01" --env-var "patient_name=Nav6 Android" --env-var "email=nav6android@gmail.com" --env-var "password=Test!123" --env-var "phone="+17813230226" --verbose
```

POSTMAN CLI Filled In 2

Wednesday, June 10, 2026 11:38 AM

postman login --with-api-key PMAK-6a23393a28180a0001a82455-f61161159a59a6b975f93aca763d5554e4
Prasuna api key PMAK-6a56a2beb6a9b3000120face-ad762148f31fce6baa8398972fcc6d97a7

| | |
|------------------------------|---|
| Postman Collect UID | 52864576-d9174d04-011e-4d5f-aca9-0001b00910e6 |
| Postman CLI environment ID | 31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85 |
| Register HES ID | 52864576-d3c5eacc-6783-4863-b76b-5f9a7f24b244 |
| Create Patient HES ID | 52864576-50086e72-309a-421d-98da-39a1e244a5e7 |
| Create Sensor Calibration ID | 52864576-ebb454ce-2f57-41cd-950a-8a5253994009 |
| Create Hospital Reading ID | 52864576-06b12575-0c09-4f79-9fc6-ce543bf81b70 |
| Create Home Reading ID | 52864576-47436003-a58d-4abf-8276-701281a3ceed |
| Br_hes_uuid | 52db19da-d402-4adf-8a83-fc9f71e31a2b |
| Hes_tokem | 72QAC7xFxcS8rncWZJR3qxdWQzwrr |

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i
"52864576-86b8eeaa-9e22-47f2-a54f-e09cddf69975" --timeout-request 30000 --env-var "_cli=1" --env-var "sensor_serial=1A4300" --env-var
"patient_dob=2020-01-01" --env-var "patient_name=Homescreen14 Android" --env-var "email=homescreen14@gmail.com" --env-var "password=Test!
123" --env-var "phone=+17813230226" --verbose

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-98ab1f66-82cb-4953-927f-a79f7f9bb30b?action=share&creator=54964642&ctx=documentation&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/collection/32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316?action=share&creator=54964642&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-52d7d30d-6920-44c7-8a93-53e71cfc1e7d?action=share&creator=54964642&ctx=documentation&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-1c16b8d7-4448-43d0-a74c-5a6173b9e1b2?action=share&creator=54964642&ctx=documentation&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-1c16b8d7-4448-43d0-a74c-5a6173b9e1b2?action=share&creator=54964642&ctx=documentation&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

Collection id 32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316
Postman cli enviromne tid 41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f
Create patient hjes id 32592187-98ab1f66-82cb-4953-927f-a79f7f9bb30b

Command 1

postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i
"32592187-52d7d30d-6920-44c7-8a93-53e71cfc1e7d" --timeout-request 30000 --verbose

Command 2

postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i
"32592187-1c16b8d7-4448-43d0-a74c-5a6173b9e1b2" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=6c8ae382-a3cf-43db-9e97-
d9c3b805ec4a" --env-var "hes_token=fqdKQnu1EphybwfTf3ZkjRPqJWo9Uwrr" --env-var "sensor_serial=1A7999" --env-var
"patient_first_name=Sam" --env-var "patient_last_name=RTR" --env-var "patient_dob=2020-01-01" --env-var "patient_phone=+15086679268" --env-
var "time_of-acquisition=2026-07-03T10:20:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-

var "avg_diastolic=15" --verbose


<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/52864576-50086e72-309a-421d-98da-39a1e244a5e7?action=share&creator=54964642&ctx=documentation&active-environment=31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/collection/52864576-d9174d04-011e-4d5f-aca9-0001b00910e6?action=share&creator=54964642&active-environment=31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/collection/52864576-d9174d04-011e-4d5f-aca9-0001b00910e6?action=share&creator=54964642&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/collection/32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f16?action=share&creator=54964642&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

**** Br_hes_uuid and Hes_token could be different double check in POSTMAN browser**

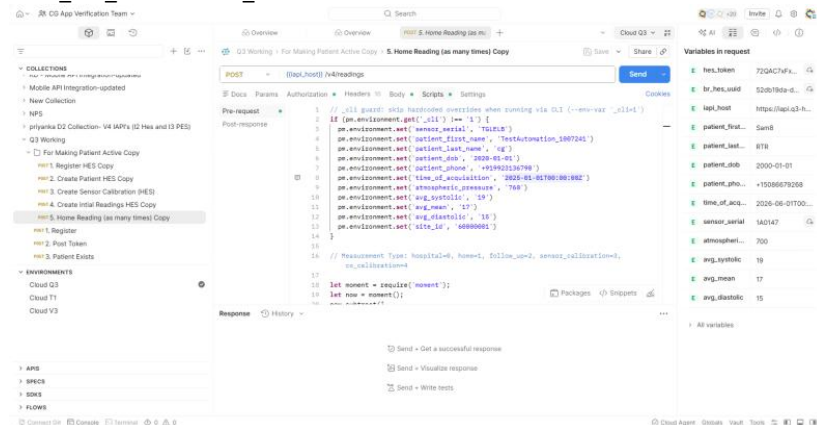

"Register HES"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-d3c5eacc-6783-4863-b76b-5f9a7f24b244" --timeout-request 30000 --verbose

"Create Patient HES"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-50086e72-309a-421d-98da-39a1e244a5e7" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=57f8ca55-2b80-4883-851d-7b9437f52026" --env-var "hes_token=TKy2HRKABZeyrZqQVwon7oHDSj53grr" --env-var "sensor_serial=1A4300" --env-var "patient_first_name=Homescreen14" --env-var "patient_last_name=Android" --env-var "patient_dob=2020-01-01" --env-var "patient_phone="+17813230226" --verbose

"Sensor Calibration"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-ebb454ce-2f57-41cd-950a-8a5253994009" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=57f8ca55-2b80-4883-851d-7b9437f52026" --env-var "hes_token=TKy2HRKABZeyrZqQVwon7oHDSj53grr" --env-var "sensor_serial=1A4300" --env-var "calibration_date=2026-01-01T10:20:00Z" --verbose

"Hospital Reading"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-06b12575-0c09-4f79-9fc6-ce543bf81b70" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=57f8ca55-2b80-4883-851d-7b9437f52026" --env-var "hes_token=TKy2HRKABZeyrZqQVwon7oHDSj53grr" --env-var "sensor_serial=1A4300" --env-var "patient_first_name=Homescreen14" --env-var "patient_last_name=Android" --env-var "patient_dob=2020-01-01" --env-var "patient_phone="+17813230226" --env-var "time_of_acquisition=2026-07-01T10:20:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --verbose

"Home Reading"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-47436003-a58d-4abf-8276-701281a3ceed" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=

```
57f8ca55-2b80-4883-851d-7b9437f52026" --env-var "hes_token=TKy2HRKABZeyzrZqQVwon7oHDSj53grr" --env-var "sensor_serial=1A4300" --env-var "patient_first_name=Homescreen14" --env-var "patient_last_name=Android" --env-var "patient_dob=2020-01-01" --env-var "patient_phone="+17813230226" --env-var "time_of_acquisition=2026-08-07T10:20:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=23" --env-var "avg_mean=22" --env-var "avg_diastolic=21" --verbose
```

30 days ago is the 5th july

```
postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-47436003-a58d-4abf-8276-701281a3ceed" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=7b0eb185-945e-4666-a119-65c7248d4c15" --env-var "hes_token=QqwXsqg2piyPbRHEYrmvNJV3gZtqpwr" --env-var "sensor_serial=1A3941" --env-var "patient_first_name=Sam47" --env-var "patient_last_name=RTR" --env-var "patient_dob=2020-01-01" --env-var "patient_phone="+15086679268" --env-var "time_of_acquisition=2026-07-05T10:20:00-06:00" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=16" --verbose
```

T1 POSTMAN CLI Filled In

Wednesday, June 10, 2026 11:38 AM

postman login --with-api-key PMAK-6a23393a28180a0001a82455-f61161159a59a6b975f93aca763d5554e4
Prasuna api key PMAK-6a56a2beb6a9b3000120face-ad762148f31fce6baa8398972fcc6d97a7

| | |
|------------------------------|---|
| Postman Collect UID | 32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316 |
| Postman CLI environment ID | 41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f |
| Register HES ID | 32592187-52d7d30d-6920-44c7-8a93-53e71cfc1e7d |
| Create Patient HES ID | 32592187-98ab1f66-82cb-4953-927f-a79f7f9bb30b |
| Create Sensor Calibration ID | 32592187-52dafb95-e31f-4053-aa4b-f07904d90f33 |
| Create Hospital Reading ID | 32592187-255f16a7-a4a5-49fd-8dab-b3343524793a |
| Create Home Reading ID | 32592187-1c16b8d7-4448-43d0-a74c-5a6173b9e1b2 |

32592187-906e0b7c-ae6b-4541-8c3f-4f7d21f64e80

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-30c8e0fd-bbcf-4016-812b-fb5cd96bd382?action=share&creator=54964642&ctx=documentation&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-fbb546e9-5be9-4017-9958-ee5da03fbf13?action=share&creator=54964642&ctx=documentation&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

B2C: 32592187-30c8e0fd-bbcf-4016-812b-fb5cd96bd382
Fetch contacts 32592187-fbb546e9-5be9-4017-9958-ee5da03fbf13

```
postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i  
"32592187-906e0b7c-ae6b-4541-8c3f-4f7d21f64e80" --timeout-request 30000 --env-var "_cli=1" --env-var "sensor_serial=1B0130" --env-var  
"patient_dob=2020-01-01" --env-var "patient_name=Nav6 Androidt1" --env-var "email=nav6t1@gmail.com" --env-var "password=Test!123" --env-var  
"phone=+17813230226" --verbose
```

```
postman collection run '32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316' -e '41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f' -i  
'32592187-30c8e0fd-bbcf-4016-812b-fb5cd96bd382' -i '32592187-fbb546e9-5be9-4017-9958-ee5da03fbf13' --timeout-request 30000 --env-var '_cli=  
1' --env-var 'subscription_key=66a7084ab0e44f74a049d7a8e052e32c' --env-var 'email=support5android@gmail.com' --env-var 'password=Test!123' --  
env-var 'user_group=user' --verbose
```

Fetch Contacts:

```
postman collection run '32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316' -e '41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f' -i  
'32592187-30c8e0fd-bbcf-4016-812b-fb5cd96bd382' -i '32592187-fbb546e9-5be9-4017-9958-ee5da03fbf13' --timeout-request 30000 --env-var '_cli=  
1' --env-var 'email=support5android@gmail.com' --env-var 'password=Test!123' --env-var "user_group=user" --verbose
```

"Register HES"

```
postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i  
"32592187-52d7d30d-6920-44c7-8a93-53e71cfc1e7d" --timeout-request 30000 --verbose
```

"Create Patient HES"

```
postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i  
"32592187-98ab1f66-82cb-4953-927f-a79f7f9bb30b" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=496cc4e7-d8c9-4ed1-  
ad53-9567962fed89" --env-var "hes_token=57K4pjpqLKFxpckQW3NNYwcsX66a3Arr" --env-var "sensor_serial=1B1380" --env-var  
"patient_first_name=Sam80" --env-var "patient_last_name=RTR" --env-var "patient_dob=2020-01-01" --env-var "patient_phone=+15086679268" --  
verbose
```

"Sensor Calibration"

```
postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i  
"32592187-52dafb95-e31f-4053-aa4b-f07904d90f33" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=496cc4e7-d8c9-4ed1-
```

```
ad53-9567962fed89" --env-var "hes_token=57K4pjppqLKFXpckQW3NNYwcsX66a3Arr" --env-var "sensor_serial=1B1380" --env-var "calibration_date=2026-01-01T10:20:00Z" --verbose
```

"Hospital Reading"

```
postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i  
"32592187-255f16a7-a4a5-49fd-8dab-b3343524793a" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=496cc4e7-d8c9-4ed1-  
ad53-9567962fed89" --env-var "hes_token=57K4pjppqLKFXpckQW3NNYwcsX66a3Arr" --env-var "sensor_serial=1B1380" --env-var  
"patient_first_name=Sam80" --env-var "patient_last_name=RTR" --env-var "patient_dob=2020-01-01" --env-var "patient_phone="+15086679268" --  
env-var "time_of_acquisition=2026-07-01T00:00:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --  
env-var "avg_diastolic=15" --verbose
```

"Home Reading"

```
postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i  
"32592187-1c16b8d7-4448-43d0-a74c-5a6173b9e1b2" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=496cc4e7-d8c9-4ed1-  
ad53-9567962fed89" --env-var "hes_token=57K4pjppqLKFXpckQW3NNYwcsX66a3Arr" --env-var "sensor_serial=1B1380" --env-var  
"patient_first_name=Sam80" --env-var "patient_last_name=RTR" --env-var "patient_dob=2020-01-01" --env-var "patient_phone="+15086679268" --  
env-var "time_of_acquisition=2026-07-02T10:20:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --  
env-var "avg_diastolic=15" --env-var "measurement_type=1" --env-var "site_id=600000001" --verbose
```

5/13

Tuesday, May 19, 2026 3:36 PM

Should Installnapp via google play

Link to qr c ode

Aste goolge play link in qr code and scan it

Download it from google play

Android apk install

02.0.20

(official build)

On ggolge blu, 2.0.15 on aop (dependency not updated on goggle pkay)

Copy apkl from official

OA folder, -qa-apk file, transfer to mobile app and isntall by importing apk to files aop and install on android aopop

Scroll and see version 20 to confirm right version

RGITRTION

In order to register

ENROLL PATEITN

In merlin portal

2 process

Env couldb eb q3 later, v3 for now

Backend = v3 sensors portal

Enrrol pateitn button

Sam1 Bodangi1

DOB 02=02-1990

Treating physican

Loging tinot merlin portan usign premade clinican

Select v3_test_last physical

Fill out Sesnor Seiral #

Unclaimed in backend pick any sesnor (choose SWVT Lot sesnors)

Pause tfor wesbite to validate with other pateitnes to check if other

patients have same serial #, if not, then it should be blue

Patient is enrolled, but to active (schedules patient's tab)

Enrolled patient via merlin portal

Register it on mobile application

1, register via mobile app or via POSTMAN (either option ok)

AP WAY:

Register button

- Once button pressed, nav to Policy screen, 1st screen is called Startup screen
 - o EULA screen
 - o Scroll and accept, then same w/ privacy
- Patient registration screen
- Now go to merlin portal
 - o Add DOB and serial # to app
 - Next
 - o Patient Info
 - Fill up details
 - Should be same but for testing, can be
 - Enter own #
 - PW Test123
 - o Patient is registered to app, now login with email and PW

Try to login with another patient

Give patient created previously after already logging in, it won't allow

Need to uninstall application and reinstall to register new patient

No policy screen when logging in already created patient, do not see policy (skipped registration flow)

Home Screen

Systematic BP

Record blood pressure for patient

See same card on reading in global navigation panel

3 panels in reading screen, day week and month

Day is today

Week is entire week data plan is --

Month is entire month, data not recorded for entire month -- value is not present

DTP

in Merlin

Add medications in merlin, after assign goal PAP

Turn on Medication managed in App in Merlin

Save and add another to add more than 1 meds

Medication adjustment = notification is medication plan changed

If app not deleted only able to login back to original patient

Skipped policy screen and landed on home screen, already accepted it before

iOS has extra update on Apple Health that Android doesn't have

BUGS

Jira

Filters, search work filters

3 spaces, Test Management App = test cases

SRS = res for project

DevApp = dev builds

DevPortal = portal team

If bug related to portal, DevPortal

If bug for App, do DevApp

App Bugs (DevApp)

Type = Bugs

Status Verify

BUG HAS BEEN already logged, dev team has checked bug and now it needs to be verified

Verify, View workflow

Once bug is logged, open

When people working, reading, and in progress

Once it's fixed, go to verify

Once verified, assigned to Danny and we should be retesting the bug

Labels, Ready to Close Not equal to

Ex bug

Not been retested

If I lock it up,

Revisit all the test description

Reading all of it

When understood the bug, retested in the newest build

Perform steps that were mentioned to reproduce

Go to actual results and compare to expected results, bug is logged

Make sure expected results are met

Write up a comment

Env details, that the env is what the build is and what phone it is tested on

Status is resolved or not

Write expected and actual results

Add evidence in comment

Once comment is added add a Ready to Close label

If bug is still unresolved, \

Put up comment as status being unresolved

Add second and reassign the bug to aomsoen, telling them the bug is still not resolved

And put it back to back in progress, not verify and change assignee

If want to log up a bug,

Log bug in DevApp space

That team work on bug

Press Create,

Work type as Bug

Add summary of what the bug is about

Pop up, this bug has been created, press it, lands on the ticket, [P100045S056C.pdf](#) Summary has been mentioned

Need to enter the description

Nanda, Sourav joined the conversation.

Chat has been turned on for this meeting.

Bodangi, Samhitha and 2 others were invited to the meeting.

Nanda, Sourav named the meeting Samhitha KT Session.

10:29 AM Meeting started

1. APK (Share point folder) : [GLB-HF-DRx-myCardioMEMS App - QA - All Documents](#)
2. Environment V3 <Merlin Portal> : <https://v3-hfcloud.products.abbott/apps/clinician/notifications>
Username : v3_test_pat
Pd : Test#1234
3. Backend V3 : [Sensors | CMHF Customer Care](#)

4. Bug Unresolved (Ex) : [\[DRXDEVAPP-7006\] Env V3-Notification Badge on CardioGuide App Icon Disappears After App Launch With...](#)
5. Bug Resolved (Ex): [\[DRXDEVAPP-7667\] Env V3 : Error message appears when tapping "Send a new code" with the verificatio...](#)

5/14

Tuesday, May 19, 2026 4:51 PM

Logged into app

Home screen

Go thru more menu

Profile option

- Full name in merlin portal
- Email address
- Phone number
- | | |
|----------|--------------------------------------|
| Can edit | Only edit full name and phone number |
|----------|--------------------------------------|
- Don't unlock country code option
 - o Please enter valid number
 - o Updates in profile page
 - o Go to merlin portal
 - o View registered application user
 - Only find updates one here not in main portal

Support

- FAQs, contact abt, use manual

FAQs

- More info

Contact

- Clinic contact details from merlin portal
- Add primary lab,
 - o Just for user not in merlin

About

-- ant privacy policy

EULA

- Info for the user

Jas software version

GTIN- related to phone,

App model number = phone

User manual

- Guidelines and model, how to find product info

Education

- Abt hf
 - o Videos abt hf and cardiomech systems
 - o Learning
 - More info on everything

Weight

Recording weight for patient

- 3 tabs, day, week, month
 - o Each had record weight button

For other days use POSTMAN to enter days,

Use calendar icon to see previous weeks and months weights trends

Seeting

- Change password
- For ios, update to apple health button there too
- Whatever data in apple health app, it would be reflected here
- Change password
 - o Can change password, and it redirects you to login again

Install

- Settings, app info then uninstall

Jira

- Observation issues while testing (devApp space)
 - o Goes to dev team
 - o Create
 - Select space, devApp
 - Work type = bug
 - Status = open
 - Summary = what is the issue
 - o Env V3 - [what bug is]
 - No need for components
 - o If know where the component is, can add that
 - Description = min thing
 - o Environment in which the test was done
 - o Build details
 - o Test details
 - o Clinic details (for merlin portal)
 - o User details (what user was used to login)
 - o Device details (android, iOS and update, and phone model)
 - o Steps to reproduce
 - ◆ While doing what steps issue was identified
 - o Actual result = what you observed
 - o Expected result = what you should have seen
 - o Evidence (screenshots or video)
 - ◆ Just result you observed, not every step
 - ◆ Button to add image and photo
 - ◆ Attach in description
 - ◆ When single attachment is added, then can add on plus sign and directly add it
 - ◆ No attachment option if you just create a bug

- ◆ Only added when added a screenshot/video from description
 - Assign is automatic, if known can put it
 - Severity option too, low hig med
 - Lagbels
 - ◆ CG2.0 Defects
 - Affected version = Q#APP - 2.0.21
- Nees to create bug
- Spoace = drx_devapp
- Test cases is drs_testmanagementapp
 - DRXTMAPP
- Bug retesting
 - Reproduce it with steps mentioned
 - In comment section
 - Verify status
 - Reproduce usign stes with new build with the new build
 - In comment section
 - IF NTO SISSUES, AND IT WAS RESIVED, DO RESOLVED STATU
 - TESTED ON, ON WHICH UPLOID AND TDATE UPDATE
 - Env tested
 - Device tested on
 - Merlin portal version as well
 - Help about, find merlin version
- Clinic details , and full name fir the pateitn, gmal, and password
- Add evidence tht it is fixed ow
- Assign it to ourself, then add to comment and adsgn to to danny again
- Label
 - Ready to close label
 - Danny knows bug is fixed

| DOORS | Database |
|---|----------|
| Propj reqs are in CardioMEMS
Nex gen produces
Software
App110
Requirements
Have optios
2 things
SRS anf GUI specififac
SRS are the requiirememnts
For test cases, createbased on the requirements

DRxRTR, then it is a valid req for anyone NOT CLOUD
Sometimes req team modify reqs, then the description or ID would be
modifies

GUI spec
Create test cases or update test cases
Based on GUI spec | |

| | |
|--|--|
| <p>Name of sceen
In test case step, provide name of the screen
Be more specific based on GUI speciifc
Important for updating or creating test case</p> <p>String exort doc
Messages on the application</p> | |
|--|--|

5/15

Tuesday, May 19, 2026 4:51 PM

Use POSTMAN

CG APP Verification team Workspace

Many collections

Formal V3

For maing pateitn active

Register collection registers pateitn

Directly login into the app, no need for registering in the app

Copy sensor serial number

Run Register HES Copy

V3 envioemnent

Copuy same as Merlin portal enviornment

Just send number 1 (Register HES Copy)

Just change env and hit send

2nd file Create Patient HES Copy

- Copy sensor serial #
- Replace sensor serial # in Scripts tabl
- Fil in Pateint details

Year month day

Chanfge phonw

Pre-requist

1. Sensor calibration

Base;line code

When run 3rd file, it ewill fill up baseline code

Got to pre-request, fill out sensor rseial #

And press send

1. Initial readings HES Copy

- a. Ptentit istatus, inactive
- b. Readings in
- c. Also change same patient details and seiral number

Press send

Fdont need to do ahrd refresh, just tap readigns again and it should work

In irl, postamn is not real

Wqhen send it, merlin portal

Use postman to ease our use

Ned to stil make it active
Go to page 5, to make pateint active

Patient is active

Merlin folllows mm-dd-yy

AK - IAPSs collection

DTP

If pateitn is habng high reading, or low reading
When clinician sends mediicine,

Goal:17

Low:13 systolic

Hgh 21 diastolic

User = 15, normal readiung

24 = 24 high reading

11 = 11 low reading

Readingsd are for 1 week period \

Save interventions

Then can turn on dynamioc treatment plan

Merlin

- Cliniciam is going to use Merlin
- Clinician can turn on dtp for user

5th file Hime reading

Go to body

Merlin time zome

Users and location

California time

Pateint is in active pateints now

If clinic is in diff time zone

Clinic = 5 am

triggered reading 6 :30am

Will not get posted

Time postings must be before clinic time

Need to change data section

Avg diastolic = 22

Avg mean = 21

Systolic = 20

Cannot revert back once send, once reading file is generated cannot revert back

Suspect = high erasing

When hit suspect, not trigger STP

Website if trigger any other readings around it, always trigger suspect

Use another patient

The two dtp cards are triggered because of the 7 out of range readings

Extract folder

Folders are environments

Testing on QA (v3)

Connect phone and copy aok folder and paste it into the phone from computer

On phone just tap on it and it should install

Dashes on the medication screen

- Some have dashes, some don't
- Company from the medicine

Always do create patient first

And enter it in pre request > scripts

Send first reading to make patient active

Body script is where edit info about patient reading first

Make the time a long time ago, just make a baseline

Avg systolic

Avg mean

Avg diastolic

File 3 is sensor calibration

Create patient
Calibration code
Then first one reading to make patient active

Suspect = if system thinks data doesn't make sense
Out of range = sensor value has been read, but it is out of range

5/19

Tuesday, May 19, 2026 3:42 PM

None drx and out of home elevation

Patient ad_62

Logged in into the app

NON DRX INTERVENTIONS

For non drx app

Go to merlin

Got to specific patenmt

More actions

- Add interventions

6 optiojs (radio button)

Non drx intervention

Click on one

And press send

And a notification should be sent to the app

Can mark as read

Can see it in complete tab

Reading reviewed

- Can see diff texts on diff tabs

Complete screen will show read cards

Can send readign reminder

Can send lab visit (requires blood test, scehdule lab visit within next 7 days)

Treatemt plan reminderd due to bp changes

Can send a custo message

- New clinician message

Daily medications we send to app

- Non drx ta

Drx interventiosn are for DTP (high and low readigns)

Titrate daily medicne (hold = stop medicine)

Can select daily medications, select quantity and frequency

Can add message to call clinic or custom messge

- These would be considered drx intervention now

Only more actions message would be considered non drx
Saily medications will be non drx

Generate OUT OF HOME ELEVSATOON

Go to postman

Send home readings (5th file)

Out of home elevetaion - change in atmospheric pressure
Would not send intervention

Will send a notification card to app, but no intervention

5th file (honme reading)

Im body tags

Atmospjeric pressure

- Bu default wau is 753

Out of home gets generatored when add 35 to it (788+)

Reading is within valid range or out of range (irrespective of this), out of home notification will be generated

Prevous buuikds

Label is out of home elevations

Now is away from home elevation

Even text is different

Later on, nothing will be on the homen screen, everything will be on the notification screen

Go to merlin portal, go to activiti tab

Go to readings

Out of home elevation Status

7 OUT OF RANGE READINGS

Set optimal

Interventions, 7 high readingsd

Asking to hold daily medicin

Adding new medicine

And adding message to call clinci (3 things)

DTP plan must be on

Go to postman

Refresh postman every time you change patients

Sometimes overwrites, refresh before workin on new patient

Send high readings
Use 5th file home readings file

Enter patient details in scripts (pre-req)

Log into app

See the week mentions (May 13-May 19)

Trigger 1 low or 1 high
Create intervention in merlin
Create reading for May 18, 17 or any
No specific dates needed

For every patient, create new patient
Once you trigger intervention
Can't trigger for same patient on same day

Triggering right now for 7 high
Tomorrow, need to trigger other readings, then only you can get

One day, only 1 intervention can be triggered

Tomorrow you can trigger something diff for this particular patient

Can't keep trigger

May 13 1 low
May 14, 1 high

Create 6 patients,
Test it individually for each intervention

Keep same values, change date

Have to be in a row
Mentioned 3 interventions

Sends 3 separate badges
Text is different for dx notifications

Asks user to confirm that they are currently taking old medication, before accepting that they stop
medication

Early medication is not in medication list tab

In merline portal

See intervention cards for each DTP intervention
3 cards, 1 for each intervention

Don't create interventions for all of them

Only for the particular one you are testing

Have 1 low and 3 ;pw

Have 1 low a few days ago

No intervention was triggered

Even 3 lpow

No intervention was generated

That is because the past intervention is generated for TODAY

Even if the reading was taken in the past

1 intervention one day

24 hour waiting period

Trigger

Get B2C token API file

Login screen

Send to get patient ID

Post vital (weight, fluid)

Mention patient ID

Mention eright

tekemAtUtc

- Separate
- Epoch converter
- Mention date to get time stamp and a Unix value and past that in
 - o This is milliseconds value

For BP reading, use Post Bp Reading tab

<https://www.epochconverter.com/>

5/20

Wednesday, May 20, 2026 1:30 PM

By Friday, 5/22, expected to very quickly set up app and go through any bugs if needed

6/11

Thursday, June 11, 2026 2:38 PM

Windows use mvn clean

Only diff bwt n is slashes and path are diff (backward)

Windows, only do android

2 way to do it

Mvn starts test, starts running whatever is in sanity.xml file

They will do

Copy and paste multiple classes, put in all the tests

Run one after another

Not bad, but

In command itself, can add flag -Dtest = (and list all tests
-Dtest=MFAlogin, Next

Separate with commas

Kranti, or someone
Where the actual list is

Once you have list

Start with

Run all 5 of them

Run all of them for both android and ios

Keep track of what test cases are passing and failing

Running 2 cases back to back

1 test passed one failed

Will generate report-export folder

Allure >html
How many tests ran and how many passes

Show all
See which pass and which failed

Be on cardioguide branch

Make sure that

Runnin it on ios an android

6/12

Friday, June 12, 2026 11:21 AM

Sign in tests work

Will first

Specific common user COMMON_USER

Whatever user is defined in framework files

Specific index

In case ppl want to use different suser and diff patients

Sign in and mfa create users and registration too

All other test cases should use common users

MobileTestFlow

WebTestFlow (for merlin)

landignPage is defined, can go to langingPageActions all have functions

At end of test, will report

Create reports (generate no matter pass or fail)

Every run allure report gets clear

Report = has html, used results data

Results = contins data used in report

Sensor generate

Import.csv gets changes every time it gets changed

Main = framework

Fdrivfermanager = used to load and install app

Basetest = template test, all test cases extend basetest

Commonmobilepageactions = common functions

Add functons to apppageactions

Based on what screen add to that file

Apiendpoint

Postman folder = postman cli

Webpageactions = merlin actions

Webwrappers = have to create wrapper afgter creating apppageaction

Functioncould be in preconditioon or normal test floiw

Originallly were separate

Adds wrapper for function, so it knows that it is for normal test case

Or add one in preconditions to know its for preconditions

Now preconditions are part of test case so the precondition wrapper could not be used anymore

Resources = config files, framework files

And has locators

Has xpaths for all buttons

Manual and Automated Verification

PR issues conflicts

Wolfridge

Document KT to prasuna

Document KT from Rahul

Contact ishani and rahul on how they document

Definitions

Wednesday, May 13, 2026 8:22 AM

High NT proBNP = HF

GDMT = DTP (Dynamic Treatment Plan)

preserved ejection fraction (HFpEF) and with reduced ejection fraction (HFrEF)

CardioMEMS = MicroElectroMechanical Systems

For NYHA Class II and III, mild/moderate systems

PA pressure = pressure from heart to lung

BP pressure = pressure from heart to body (cuff)

Pulmonary blood pressure is normally a lot lower than regular blood pressure. Normal pulmonary artery pressure is 11-20 mm Hg at rest.

From <<https://www.heart.org/en/health-topics/high-blood-pressure/the-facts-about-high-blood-pressure/pulmonary-hypertension-high-blood-pressure-in-the-heart-to-lung-system>>

Systolic = top

Diastolic = bottom

Normal <120 : <80

High BP 130-139 : 80-89

[Understanding Blood Pressure Readings | American Heart Association](#)

CM1200 = Patient system (HERO and Wolfridge)

The CardioMEMS HF system consists of

- implantable wireless sensor
- a patient and hospital electronics system (HERO + Pillow)
- patient database (Merlin)

From <<https://pmc.ncbi.nlm.nih.gov/articles/PMC8037296/>>

Take daily readings using Pillow

Readings transmit to Merlin Portal

New gen = updated and smaller reader, more phone-like

Hospital system = standing monitor

[P100045S056C.pdf](#)

DOORS = Dynamic Object-Oriented Requirements System

Requirements management application

Modular structure = how it organizes data

In hierarchy

GTIN = global trade item number (kind of like barcode)

Waveform = graph (signal v. time)

CardioMEMS works using LC circuit (inductor capacitor)

Voltage oscillates between inductor and capacitor

Capacitor resists changes in voltage (stores energy in electric field)

Inductor resists changes in current (stores energy in magnetic field)

Sensor remains in situ (in original place)

[CardioMEMS Heart Failure System: An Up-to-Date Review - PMC](#)

Suspect readings

- Why in web app if already in PES and HES
- Suspect readings can only be changed by user if due to technical issues
 - o Can change in Q3
 - Go to readings

HERO

Wednesday, May 27, 2026 10:23 AM

Testing multiple sensors and devices interferes with signal strength

- Signal Strength stayed constant once moved away

Metal close by can interfere with the sensor (gave negative pressure reading)

Best signal is when PA sensor is horizontal (parallel to Abbott logo)

Atmospheric pressure and PA pressure are directly related

CardioMEMS works at frequency between 30 and 37.5 MHz (extremely rare to go beyond this)

- Possible devices that could interfere are radios, walkie talkies, science equipment (anything on same band)

Error #8 is if the sensor is likely operating outside of range

https://www.bfarm.de/SharedDocs/Kundeninfos/EN/01/2023/03660-23_kundeninfo_en.pdf?blob=publicationFile

Search Pressure is the value the device starts on (benchmark)

Sensor should slowly come down or up to the actual pressure

- Hypothetically, changing search pressure should not change the mean

- The farther going away from 18, the time it takes to find signal could increase

ADT Setup

Thursday, May 28, 2026

12:33 PM

Aerodynamic tester

Simulates PA pressure

Range from 15 to 65 mmHg

If sensor already has a static baseline, the ADT pressure is added on top of the static pressure

Can recalibrate sensor using POSTMAN to change the baseline code

Changes the atmospheric pressure to change the PA pressure (directly related)

In IRL, use Blueridge to recalibrate the sensor

IRL, Search pressure will be different from patient to patient

Will follow the trend, in increments of 2 mmHg

Sourav KT 5/13/2026

Sunday, June 14, 2026 10:24 PM

KT received by Sourav on 5/13/2026

Official release build = 2.0.20

Should install app via google Play

- Past google play link in QR code generator and scan it
- Download it from Google play

Not updated (only 2.0.15 in Gplay),

- Copy apk from official folder
- Transfer apk file to mobile app and install by importing apk to files app and install on android phone
 - o Can scroll to confirm right version

Patient Registration

- First enroll patient in merlin portal
 - o 2 ways to register patient,
 - o Enroll patient button
 - DOB 01-01-1990
 - Treating Physician
 - Logging into Merlin using premade clinician
 - ◆ V3_test_last as physician
 - Fill out sensor Serial #
 - ◆ Choose unclaimed in backend (SWVT lot)
 - ◆ Enter in Merlin, pause for website to validate # with other patients
 - ◇ Website checks with other patients to check if other patients have same serial #, if not, box should be blue
 - o Patient is enrolled, but not active (scheduled patients tab)
 - Register patient in mobile app or POSTMAN

Register Patient w/ app

- Once button pressed, nav to Policy screen
- First opening screen is called startup screen
 - o Then got to EULA screen and Privacy screen
- Patient registration screen
 - o In merlin portal
 - Add DOB and serial # to app, press next
 - o Patient Infor
 - Fill up details
 - Should be same as testing, but can be different
 - Enter own phone #
 - PW Test!123
 - o Patient is registered to app, now login with email and PW

If try logging in with another patient, it won't allow as device can only be registered by 1 patient
Need to uninstall app and reinstall to register new patient

No policy screen when logging in, already created patient so skipped registration flow

Home Screen

Systematic BP

- Record BP for patient
- See same card on readings on global navigation panel
 - o 3 Panels, day, week, month
 - Day is today
 - Week is entire week
 - Month is entire month, data not recorded for entire month -- value is present

DTP in Merlin

- Add medications in merlin, after adding goal PA Pressure
- Turn on Medication managed in App in Merlin
- Use save and add another button to add more than 1 meds

In app, Medication adjustment notification if medication plan is changed in Merlin

IOS has extra button in More options (Update on Apple Health)

Bugs in Jira

- Filters, search work filters
- 3 spaces
 - o TestManagementApp = test cases
 - o SRS = requirements for project
 - o DevApp = dev builds (bug for app)
 - o DevPortal = Portal team (bugs related to portal)

App bugs (DevApp)

- Type = Bugs
- Status = Verify
 - o Bug has been already logged, dev team has checked bug, and now it needs to be verified
- Verify, view Workflow
 - o When people are working, status is reading and then in progress
 - o Once fixed, status is verify
 - o Once verified, assigned to Danny and we should be retesting the bug
- Labels, Not equal to Ready to Close

Example Bug

- Not been retested yet
- Revisit all the test description and read it all
- When understand bug, retest it in new build
- Perform steps that were mentioned to reproduce test
- Go to actual results and compare to expected results, bug is logged
- Make sure expected results are met
- Write up comment
 - o Environment details, what the build is, and what phone it is tested on
- Status is resolved or not
- Write expected and actual results
- Add evidence in comment s
- Once comment is added, add ready to close button

If bug is still unresolved

- Put up comment as status being unresolved
- Add second comment and reassign bug to original person, telling them bug is still not resolved
- Put it back to in progress

If want to log up bug,

- Log bug in Dev App space
- Press create
- Work type as Bug
- Add summary of what the bug is about
- Have pop up, this bug has been created
 - o Press it and land directly on ticket
 - o Enter description

Time of acquisition:2025-01-01T00:00:00Z

1. APK (Share point folder) : [GLB-HF-DRx-myCardioMEMS App - QA - All Documents](#)
2. Environment V3 <Merlin Portal> : <https://v3-hfcloud.products.abbott/apps/clinician/notifications>
Username : v3_test_pat
Pd : Test#1234
3. Backend V3 : [Sensors | CMHF Customer Care](#)
4. Bug Unresolved (Ex) : [\[DRXDEVAPP-7006\] Env V3-Notification Badge on CardioGuide App Icon Disappears After App Launch With...](#)
5. Bug Resolved (Ex): [\[DRXDEVAPP-7667\] Env V3 : Error message appears when tapping "Send a new code" with the verificatio...](#)

Ameya KT 5/14/2026

Sunday, June 14, 2026 10:26 PM

KT received by Ameya on 5/14/2026

More Options Tab

- Logged into App and on Home screen
- More Menu tab
 - o Profile option
 - Full name in Merlin Portal
 - Email Address
 - Phone number
 - Can edit only full name and phone number
 - As of 2.0.22, can change country code
 - o If profile page updated
 - o Go to Merlin portal
 - View registered application user
 - Only find updates here, not on the main portal
 - o Support
 - FAQs, Contact, About, User Manual
 - FAQs have more info
 - Contact
 - Clinic contact details from merlin portal
 - Add primary lab,
 - o Just for user not in merlin
 - o About
 - Privacy Policy and EULA (info for user)
 - Has software version
 - GTIN - related to phone
 - App model number
 - o User manual
 - Guidelines, how to find product info
 - o Education
 - About HF
 - Videos about HF and CardioMEMS system
 - Learning (more info on everything)
 - o Weight
 - Recording weight for patient
 - 3 tables: day, week, month
 - Each tab has record weight button
 - For other days, use POSTMAN to enter weight
 - User calendar icon to see previous weeks and months weights trends
 - o Settings
 - Change password
 - For IOS, has update to apple health button there too
 - Whatever data in apple health app, it would be reflected here
 - Change password
 - Can change password, and it redirects you to login again

Uninstall on Android

- Settings, app info then uninstall

Loggin Bug in Jira

- Observation issues while testing (devApp space)
 - o Goes to dev team
 - o Create
 - Select space, devApp
 - Work type = bug
 - Status= open
 - Summary = what is the us
 - Env V3 - [what bug is]
 - No need for components
 - If know where the component is, can add that
 - Description = min thing
 - Environment in which the test was done
 - Build details
 - Test details
 - Clinic details (for merlin portal)
 - User details (what user was used to login)
 - Device details (android, IOS and update, and phone model)
 - Steps to reproduce
 - ◆ While doing what steps issue was identified
 - Actual result = what you observed
 - Expected result = what you should have seen
 - Evidence (screenshots or video)
 - ◆ Just result you observed, not every step
 - ◆ Buttin to add image and photo
 - ◆ Attach In description
 - ◆ When single attachment is added, then can add on plus signa nd directly add it
 - ◆ No attachment option if you just create a bug
 - ◆ Only added when added a screenshot/video from description
 - Assign is automatic, if known can put it
 - Severity option too, low high med
 - Labels
 - ◆ CG2.0 Defects
 - Affected version = Q#APP - 2.0.21
- Nees to create bug
- Space = drx_devapp
- Test cases is drs_testmanagementapp
 - o DRXTMAPP
- Bug retesting
 - o Reproduce it with steps mentioned
 - o In comment section
 - o Verify status
 - o Reproduce using states with new build with the new build
 - o In comment section
 - IF NTO SISSUES, AND IT WAS RESIVED, DO RESOLVED STATU
 - o TESTED ON, ON WHICH UPLOID AND TDATE UPDATE
 - Env tested
 - Device tested on

- Merlin portal version as well
 - Help about, find merlin version
- Clinic details , and full name fir the patient, gmail, and password
- Add evidence that it is fixed ow
- Assign it to ourself, then add to comment and assign to danny again
- Label
 - Ready to close label
 - Danny knows bug is fixed

DOORS

- Requirements in CardioMEMS > Next Generation Products > Software
 - App110 = CardioGuide App
 - Requirements
 - Have SRS and GUI Specification
 - SRS are the requirements
 - Test cases are created based on requirements
 - DRxRTR is a valid req for anyone NOT CLOUD
 - Sometimes req team modify reqs, then the description or ID would be modified
 - GUI Spec
 - Create test cases or update test cases based on GUI Spec
 - Name of screen
 - In test case step, provide name of screen
 - Be more specific based on GUI specific
 - Important for updating or creating test case
 - String export doc
 - Messages and strings on application

Sourav KT 5/15/2026

Sunday, June 14, 2026 10:26 PM

KT received by Sourav on 5/15/2026

POSTMAN

- Workspace = CG App Verification team
- Many Collections
 - o AK - IAPSs collection for V3
 - For making Patient Active folder
 - Register collection registers patient
 - Can directly log into app, no need for registering in the app
 - Copy sensor serial number
 - Run Register HES Copy
 - ◆ V3 environment
 - ◆ Copy same as Merlin portal
 - ◆ Just change environment and hit send
 - Create Patient HES Copy
 - ◆ Copy Sensor Serial #
 - ◆ Replace # in pre-request > scripts tab
 - ◆ Fill in patient details
 - ◇ Date = Year Month Date, in Merlin, its month day year
 - ◇ Change phone #
 - Sensor Calibration
 - ◆ Baseline code in Merlin
 - ◆ When run 3rd file, fill up baseline code
 - ◆ Pre-request, fill out sensor serial #
 - ◆ Send
 - Initial HES reading
 - ◆ Patient status still inactive
 - ◆ Change same patient details and serial #
 - ◆ Send
- In Merlin, don't need to do hard refresh, just tap readings and it should work
- IRL, postman not real
 - o User postman to ease testing
 - Patient still inactive, use 5th file to make active
- Now patient is active after sending 5th file

DTP

- If patient is having high reading or low reading
- Clinician can send medicine
- Ex
 - o Goal = 17
 - o Low = 13 systolic
 - o High = 21 diastolic
 - o User = 15, normal reading
 - o User = 24, high reading
 - o User = 11, low reading
- Readings are for 1 week period
- Save interventions in merlin

- Then can turn on dynamic treatment plan
 - o Clinician is going to use Merlin
 - o Clinician can turn on DTP for user
- Home Reading POSTMAN file
 - o Go to Body
 - o Change data section
 - Avg systolic = 22
 - Avg mean = 21
 - Avg diastolic = 20
 - o Patient is in Active Patients tab in Merlin now
- Merlin time zone = California time
 - o If clinic is in different time zone
 - Clinic = 5 AM
 - Triggered reading at 6:30 AM
 - Reading will not get posted
 - Time postings must be before clinic time
- Cannot revert back once send, once reading file is generated cannot revert back
- Suspect = high reading
 - o When hit suspect, not trigger DTP
 - o If trigger any other readings around it, always trigger suspect
 - o Use another patient
- 2 DTP cards are triggered because of the 7 out of range readings

Installing App

- Extract folder
- Folders are environments
- Testing on QA (v3)
- Connect phone and copy apk folder and past it in the phone from computer
- On phone just tap it and it should install

Dashed on the medication screen

- Some have dashes, some don't
- Company from the medicine

Helpful Notes

- Always create patient first
- And enter it in pre request > scripts
- Send first reading to make patient active
- Body script is where edit info abt patient reading first
- Make the time a long time ago, just make a baseline
- Enter
 - o Avg systolic
 - o Avg mean
 - o Avg diastolic
- File 3 is sensor calibration
 - o Create patient
 - o Calibration code
 - o Then first one reading to make patient active
- Suspect = if system thinks data doesn't make sense
- Out of range = sensical value has been read, but it is out of range

Ameya KT 5/19/2026

Sunday, June 14, 2026 10:26 PM

KT received by Ameya on 5/19/2025

Non-Drx Interventions

- Go to Merlin, select specific patient
- More actions button, add interventions
 - o There are 6 options (radio buttons)
 - Reading reminder
 - Lab visit (lab requires blood test, schedule lab visit within next 7 days)
 - Treatment plan reminder (due to BP changes)
 - Can send custom message
 - Notification badge says "New clinician message"
 - o These are all non drx interventions
 - Click on one, press send
 - Notification should be sent to the app
 - Can be marked as read
 - ◆ See in complete tab in app
 - ◇ Complete screen shows read cards
 - See different texts on different notification tabs
- Daily medications sent to app are also considered non drx intervention

Drx intervention are only the DTP notifications (for high and low readings)

- Can titrate daily medication (hold = stop medicine)
 - o Select daily medication, select quantity and frequency
- Can add message to call clinic or custom message
 - o These would be considered drx intervention now
 - Only more actions messages would be considered non drx
 - Daily medications will be non drx

Out of Home Elevation

- Out of home elevation is change in atmospheric pressure
 - o This would not send an intervention
 - o Would send a notification to the app, but no intervention
- Go to postman
- Go to Send home readings (5th file)
- In body tags, AtmosphericPressure
 - o By default, value is 753
 - o Out of home gets generated when 35 is added to it (788+)
 - o Irrespective of if the reading is within the valid range or not, out of home notification will be generated
- In previous builds
 - o Label of notification badge is "Out of home elevation"
 - o Now it is "Away from home elevation"
 - o Even text is different
 - o Later on no badge will be on the home screen, everything will be on the notification screen
- Go to Merlin portal, go to activity tab
- Go to readings
- Out of home elevation status is there

7 Out of Range Readings

- Set optimal ranges
- Interventions for 7 high readings (3 total things)
 - o Adding a new medicine
 - o Asking to hold daily medicine
 - o Adding message to call clinics
- DTP plan must be on
- Go to postman
 - o Refresh postman every time you change patients
 - Sometimes POSTMAN overwrites, refresh before working on new patient
- Use 5th file to send high readings
 - o Enter patient details in scripts (pre-req)
- Log into app
- See week (May 13-May 19)
- Trigger 1 low or 1 high
- Create intervention in merlin
- Create readings for May 18, 17, or any of them
 - o No specific dates needed
- For every test, create new patient (for every new intervention case)
 - Create 6 patients
 - Test it individually for each intervention
 - o Once you trigger a intervention
 - o Cannot trigger another one for the same patient on the same day
 - o Triggering right now for 7 high
 - o Tomorrow need to trigger other readings, then you can get another intervention
 - o One day, only 1 intervention can be triggered
 - o Tomorrow you can trigger something different for this particular patient
 - o Cannot just keep triggering interventions
- Keep same values, change the date
 - o Have to be in a row
 - o Mentioned 3 interventions
 - Sends 3 separate notification badges
 - Text is different for each drx notification
- Ask user to confirm that they are currently taking old medication before accepting that they stop medicine
 - o Now earlier medication is not in the medicine list tab
- In Merlin portal
 - o Can see intervention cards for each DTP intervention
 - 3 cards, 1 for each intervention

Helpful Notes

- Don't create interventions for all 6 cases
- Only create one for the particular one you are testing
- If have 1 low and 3 low interventions
 - o Send 1 low readings few days ago
 - No intervention was triggered
 - Even 3 low readings, no intervention was triggered
 - This is because the past intervention is generated for TODAY
 - ◆ Notification was sent TODAY, even if reading was taken in past
- 1 intervention for 1 day
 - o 24 hour waiting period

Sending Weight and BP

- Get B2C token API file, just send it
- Got to Login page
 - o Send to get patient ID
 - Edit email address for patient
- Post vital(weight, fluid) page to send weight
- Post BP Reading page to send BP reading
 - o Mention patient ID
 - o Mention specific weight
 - o takenAtUtc tag
 - Epoch converter: <https://www.epochconverter.com/>
 - Mention date to get time stamp and a Unix value and paste into this tag
 - This is a millisecond value

Rahul KT

Friday, June 12, 2026 3:20 PM

KT received on 6/11/2026 and 6/12/2026

Running Test Cases

Windows and Mach have different commands. The main difference is the slashes are backwards and the file path are different.

When executing on Windows, can only run test cases on Android

Run test cases on cardioguide branch

Commands to run test cases:

For Android :

MAC/Linux = mvn clean test -DplatformName=Android -Ddevice.udid=RZCXA1RK7FV -

Dapp=/Users/vikramrajchandran/Documents/V3apk/myCardioMEMS-qat1.apk

Windows = mvn clean test "-DplatformName=Android" "-Ddevice.udid=e2ea92f8" "-


Dapp=C:/Users/CVIKRA/Downloads/V3apk/myCardioMEMS-qa-2.0.14.apk"

For IOS :

MAC/Linux = mvn clean test -DplatformName=iOS -Ddevice.udid=00008101-000E7D6E34EA001E

Two ways to run multiple test cases:

1. Mvn clean command starts executing test, running whatever is in sanity.xml file
 - a. Copy and paste classes and put different test case paths in each tag

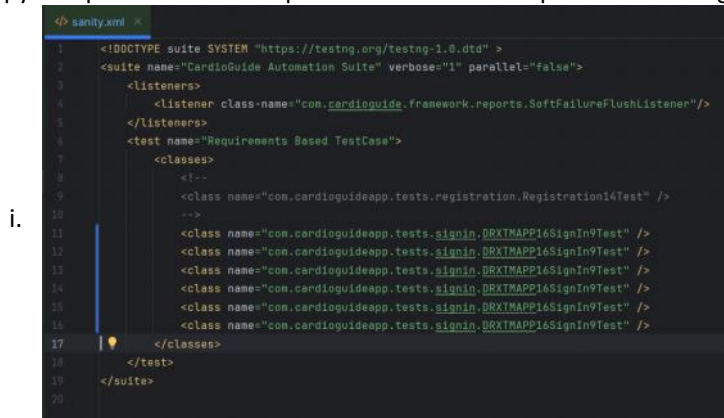

- b. Runs one after another
2. In command itself, can add flag -Dtest=
 - a. List all test cases, separated by comma (no path necessary, just name)
 - i. mvn clean test "-DplatformName=Android" "-Ddevice.udid=e2ea92f8" "-Dapp=C:/Users/CVIKRA/Downloads/V3apk/myCardioMEMS-qa-2.0.14.apk" "-Dtest=DRXTMAPP16SignIn9Test,DRXTMAPP16SignIn9Test"

Kranti or someone has actual list. Once you have actual list

- Run all of them at once
- Run all of them for both android and IOS (on mac)
- Keep track of what test cases are passing and failing

Find Report

- 1 test case passed one failed
- Report generated in report_export folder
 - o Allure > html file
 - Has how many tests ran and how many passed
 - Press show all to see which ones passed and which failed

Note from Danny on DSU 2 6/12/26

*Must perform 2 step verification on pristine code system. No edits should be necessary
Every time a run is performed, should restart and pull all the code again

Springboot setup

Guidelines

- Connect test plan to execution
- Create defects
 - o Marked as (old), not as good

Execution tracker

- Used by any person who is going to execute a test case
- All test cases here are good to run
- Not all 206 tc, some still need updates in tc steps
- Of entered here, its good to go and good to execute
- Type of test case column C
 - o Whether it is manual, automated, etc
 - o Make sure view is only used by us
 - Not by everyone else
- Once start execution
 - o Put execution type as manual (col D)
put status = in progress
 - o Do type in D
 - o Once pass for Android, then do in progress for iOS
 - o Then do Pass/fail once both executions done
 - o Do comments in column L
- Go to deviation sheet and add info there
- Automation
 - o Run through batch (ALL of them)
 - Do at night for auto FULL
 - o Unless they are auto-partial
 - o Only run tc that don't have tags such as in progress/in review
 - o Partially ones, do whenever
- Manual tc are already run
 - o No more manual tc
 - o In review/in progress automated tc
 - No manual tc to run, we can pick, up the auto in review/in [progress to run for manual
 - Put date when you pick them up in column D (execution type)
- Add bug in the bugs column if you logged anything
- Take screenshots and name them for manual exec
- Do on env t1
- Make sure postman cli is correct make it 1.31.3 t

Make tickets for android/iOs
Link them to test plan

Need to capture evidence and add comments
- save in sharepoint, don't need to link in jira

Continue execution even if it fails

Link tc in linked work item (even for automation)
Relates to ...

Annotations must be in present test
Add more info, not just copying expected result

Manual Testing

Thursday, June 4, 2026 10:52 AM

Requirement testing

- Have to test the requirement
- Test app based on the specification on the requirements
- Test cases were based on the requirements
 - o Check if they were satisfied or not

Ad hoc

- Not with any planned workflow
- Randomly take system
- Randomly play with system
 - o Break the system and find bug
- Do whatever you feel like
- Not following any requirements

Test Flight

- Use Test Flight to download ipa files to IOS phones
- No need for Macbook

- Notification Banner is not present on the HOME SCREEN as mentioned in the Release Notes for 2.0.26.

- For a new Patient when Daily Medication is added and saved on Merlin Portal and when user logs in to the CG application, the Medication Adjustment Card is not visible or present on the Active Notifications Screen.

Tried to add new Daily medication and checked if notification is received or not. No notification was received.

- Also, Non-DRx interventions like Call Clinic, Reading Reminder etc.. are not received or visible in Active Notifications Screen.

- For DRx interventions also, not receiving DRX Medication Adjustment card on Active Notifications Panel even after DTP Plan is triggered.

- Away from Home Elevation card is not triggered on the CG application but it's visible on Merlin Portal.

- The days sequence on the Weight Card is not in correct order for the Week and Month views

Ex Test Cases

Thursday, June 4, 2026 1:52 PM

DRXTMAPP-2337

[\[DRXTMAPP-2337\] Notification19 4-SRS 2867, 3967 - Jira](#)

Verifies that PA values are correctly displayed
Running test case on Q3, even though it says V3

Testing on Android 16 Samsung Galaxy A16

NOTE: In precondition DRXTMAPP 4874, Step 5 says to select any of the Unclaimed sensors
- Was told to only use SWVT sensors, but will follow test case steps for now

Use the new tab button when trying to view pre condition steps

In precondition 2337, does not say to open the app before step 3 or to allow the file access permission.
Does not say to press allow or don't allow for the notification panel

Downloaded app using Google Play Store QR code

[Manual Test Cases BODANSX.xlsx](#)

Two-Step Verification

Tuesday, June 16, 2026 1:58 PM

Cloning: git clone --branch CardioGuide2.0-Web-API_Integrations
<https://bitbucket.sjm.com/scm/myc/cardioguideapp-testautomation.git>

Download all libraries" mvn clean install -DskipTests

Need phone to be connected to wifi for springboot to work

Time Issue

Monday, August 3, 2026 2:08 PM

Deleted App, Restarted phone

1st Download

| | |
|-----------------------|----------|
| Open: | 1 11 34 |
| - Notification accept | 1 11 48 |
| To startup | 4 econds |

2nd download

| | |
|-----------------------|---------|
| Open: | 1 12 43 |
| - Notification accept | 1 12 55 |
| To startup | instant |

- 3rd downlpad

| | |
|-----------------------|---------|
| Open: | 1 13 32 |
| - Notification accept | 44 |
| To startup | instant |

12
24
Instant

54
10
17

41
54
Instant

18
31
Instant

57
10
instant

Dry Runs 2

Wednesday, August 5, 2026 11:39 AM

```
mvn clean test -DplatformName=Android -Ddevice.udid=<device udid> -Dapp=C:\Users\TIPIRLX\Downloads\myCardioMEMS-q3-2.0.26-159262.apk -Dxray.enabled=true -Dxray.cloud.clientId=<client ID> -Dxray.cloud.clientSecret=<client secret> -Djira.api.token=<JIRA Tolen> -Dxray.test.plan.key=DRXTMAPP-5648 -Dxray.execution.include.datetime=false -Dxray.execution.include.counter=false -Dxray.execution.prefix="LMR - Dry Run 1" -Dxray.project.key=DRXTMAPP
```

Client ID:

86A1D8829C494260B89DF497E5464EE1

Client Secret:

7377cf7ff35b39703a1021fdecfa976423736198928d06b9128f1ed820acc476

If green, ready for manual execution

Do not run automation if not complete incolumn E

1A4242

```
postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-86b8eeaa-9e22-47f2-a54f-e09cdbf69975" --timeout-request 30000 --env-var "_cli=1" --env-var "sensor_serial=1A4242" --env-var "patient_dob=2020-01-01" --env-var "patient_name=Nav6 Android" --env-var "email=nav6android@gmail.com" --env-var "password=Test!123" --env-var "phone="+17813230226" --verbose
```

Postman Collection UID:- 52864576-d9174d04-011e-4d5f-aca9-0001b00910e6

Postman CLI Environment ID: 31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85

Register API:- 52864576-86b8eeaa-9e22-47f2-a54f-e09cdbf69975

| | |
|------------------------------|---|
| Postman Collect UID | 52864576-d9174d04-011e-4d5f-aca9-0001b00910e6 |
| Postman CLI environment ID | 31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85 |
| Register HES ID | 52864576-d3c5eacc-6783-4863-b76b-5f9a7f24b244 |
| Create Patient HES ID | 52864576-50086e72-309a-421d-98da-39a1e244a5e7 |
| Create Sensor Calibration ID | 52864576-ebb454ce-2f57-41cd-950a-8a5253994009 |
| Create Hospital Reading ID | 52864576-06b12575-0c09-4f79-9fc6-ce543bf81b70 |
| Create Home Reading ID | 52864576-47436003-a58d-4abf-8276-701281a3ceed |
| Br_hes_uuid | 52db19da-d402-4adf-8a83-fc9f71e31a2b |
| Hes_tokem | 72QAC7xFxsWs8rncWZJR3qxwDWqzwrr |

Udid 74d93de0-b02c-4078-8204-980c409b6898

Token EpQBqEELrBytLfwfejuzzxzXyjYPBArr

"Create Patient HES"

```
postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-50086e72-309a-421d-98da-39a1e244a5e7" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=74d93de0-b02c-4078-8204-980c409b6898" --env-var "hes_token=EpQBqEELrBytLfwfejuzzxzXyjYPBArr" --env-var "sensor_serial=1A4242" --env-var
```

"patient_first_name=Nav6" --env-var "patient_last_name=Android" --env-var "patient_dob=2020-01-01" --env-var "patient_phone=+15086679268" --verbose

"Sensor Calibration"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-ebb454ce-2f57-41cd-950a-8a5253994009" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=74d93de0-b02c-4078-8204-980c409b6898" --env-var "hes_token=EpQBqEELrBytLfwfejuzzxzXyjYPBArr" --env-var "sensor_serial=1A4242" --env-var "calibration_date=2026-01-01T10:20:00Z" --verbose

"Hospital Reading"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-06b12575-0c09-4f79-9fc6-ce543bf81b70" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=74d93de0-b02c-4078-8204-980c409b6898" --env-var "hes_token=EpQBqEELrBytLfwfejuzzxzXyjYPBArr" --env-var "sensor_serial=1A4242" --env-var "patient_first_name=Nav6" --env-var "patient_last_name=Android" --env-var "patient_dob=2020-01-01" --env-var "patient_phone=+15086679268" --env-var "time_of_acquisition=2026-07-01T10:20:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --verbose

"Home Reading"

postman collection run "52864576-d9174d04-011e-4d5f-aca9-0001b00910e6" -e "31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85" -i "52864576-47436003-a58d-4abf-8276-701281a3ceed" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=74d93de0-b02c-4078-8204-980c409b6898" --env-var "hes_token=EpQBqEELrBytLfwfejuzzxzXyjYPBArr" --env-var "sensor_serial=1A4242" --env-var "patient_first_name=Nav6" --env-var "patient_last_name=Android" --env-var "patient_dob=2020-01-01" --env-var "patient_phone=+15086679268" --env-var "time_of_acquisition=2026-07-02T10:30:00Z" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --verbose

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-98ab1f66-82cb-4953-927f-a79f7f9bb30b?action=share&source=copy-link&creator=54964642>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-98ab1f66-82cb-4953-927f-a79f7f9bb30b?action=share&source=copy-link&creator=54964642>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/32592187-98ab1f66-82cb-4953-927f-a79f7f9bb30b?action=share&creator=54964642&ctx=documentation&active-environment=41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f>

<https://abbott-hf-software.postman.co/workspace/CG-App-Verification-Team~725d34a1-048c-40d6-bc30-7f558fc909b5/request/52864576-50086e72-309a-421d-98da-39a1e244a5e7?action=share&creator=54964642&ctx=documentation&active-environment=31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85>

mvn clean test -DplatformName=Android -Ddevice.udid=<Device udid> -Dapp=<App .apk path> -Dxray.enabled=true -Dxray.doud.dientId=<dient ID> -Dxray.doud.dientSecret=<dient secret> -Djira.api.token=<JIRA Token> -Dxray.test.plan.key=DRXTMAPP-5648 -Dxray.execution.indude.datetime=alse -Dxray.execution.indude.counter=alse -Dxray.execution.prefix="LMR - Dry Run 1" -Dxray.project.key=DRXTMAPP

From <https://abbott.sharepoint.com/:w:/r/sites/US-EPHF-CMEMSDynamicRx/_layouts/15/Doc.aspx?sourcedoc=%7BDC7EC4C7-8DFE-4AA5-896C-5BCEBB075A96%7D&file=XXXXXXXX%20CardioGuide%20Test%20Automation%20NPS%20Documentation.docx&wdLOR=c66B5CB53-30D7-4BF0-B5D2-DC8ADC618EA9&action=default&mobileredirect=true>

Monday, August 17, 2026 10:15 PM

[Characteristics of Pulmonary Artery Pressure Waveform for Differential Diagnosis of Chronic Pulmonary Thromboembolism and Primary Pulmonary Hypertension | JACC](#)

Overview

Wednesday, June 10, 2026 12:25 PM

Indeoenent

- Heart how it works
- Theory electrical and mech
- PA why speical
 - o Waveform
 - What it notch (dicotic
 - What you can see
 - Characteris, illness
 - Velocity
 - Derivatives
- June 30
- July 17

Definitions

Friday, June 12, 2026 10:38 AM

Heart beat = frequency of beats/minute
how often blood gets pushed out
Normal range = 60-100 bpm at rest

Stroke volume of blood ejected by heart with each beat (mL/beat)

Cardiac output = amount of blood pumped by heart per minute (Heart rate x stroke volume = mL/minute)

[Physiology, Cardiac Output - StatPearls - NCBI Bookshelf](#)
[Understanding Cardiac Output and What It Means](#)

Systolic = maximum pressure when heart beats and pumps blood OUT of body
Indicates amount of pressure heart exerts to pump blood throughout circulatory system
High values could mean heart is working harder than it should be
Sign of hypertension

Diastolic = minimum pressure when heart rests and refills with blood to prepare for next beat
Indicates resistance to blood flow (movement)
High value can indicate narrow blood vessels (low movement)

<https://www.youtube.com/watch?v=EH6e2JuxAho>
[Systolic vs. Diastolic Blood Pressure - GoodRx](#)

Heart disease = umbrella term for any condition effecting heart structure/function

[Know the Differences: Cardiovascular Disease, Heart Disease, Coronary Heart Disease](#)

Stroke = when blood flow to part of brain is cut off, or brain vessel in brain is blocked or burst

Heart Attack = AKA myocardial infarction
Blood flow to the heart is blocked, usually due to fat buildup or blood clot in coronary arteries

[Stroke vs. Heart Attack: The Difference | Capitol Cardiology Associates](#)

Heart Failure

Sunday, June 14, 2026 11:09 PM

Heart failure = heart muscle no longer able to pump oxygen rich blood to rest of body

Can affect right, left, or both

HRrEF = heart cannot contract very well (Systolic heart failure)

HF with REDUCED ejection fraction

HFpEF = heart is too stiff and does not full up with blood easily (diastolic heart failure)

HF with PRESERVED ejection fraction

[Heart failure: MedlinePlus Medical Encyclopedia](#)

Anatomy of Heart

Sunday, June 14, 2026 11:04 PM

Oxygen rich blood enters the heart from the lungs and goes out to the body

Lungs --> O2 rich blood --> heart --> Body

Oxygen poor blood enters the heart from the body and goes out to the lungs

Body --> O2 poor blood --> heart --> Lungs

[Heart failure: MedlinePlus Medical Encyclopedia](#)

Useful Info

Monday, June 29, 2026 9:51 AM

Windows PC name: APLTmuQyd2OJQL6

Macbook name: M9H7Y7QC4Y

GitHub username: samhitha-bodangi_ABT

Client ID:

86A1D8829C494260B89DF497E5464EE1

Client Secret:

7377cf7ff35b39703a1021fdecfa976423736198928d06b9128f1ed820acc476

mvn clean test -DplatformName=**Android** -Ddevice.udid=f17437eb -
Dapp=/Users/BODANSX/Downloads/29t1.apk -Dxray.enabled=true -Dxray.cloud.clientId=
86A1D8829C494260B89DF497E5464EE1 -Dxray.cloud.clientSecret=
7377cf7ff35b39703a1021fdecfa976423736198928d06b9128f1ed820acc476 -
Djira.api.token=ATATT3xFfGF08lXeLMx7Vmf9U_LVX2zME1Dj5lfo1XI_zZTEth55FQAESMyErn9UMwb-
HLUBAMU9y9ADcvc9JDz4CSK5dXThqtLIWid_
1uphbm7lleAsvuWDMlexRoQAbw4F8Fy2wsm_aHXA5nn7YK_pOYYqojZvKXuA0brt3V53oHbkiZFvtE8=
013D8E57 -Dxray.test.plan.key=DRXTMAPP-6225 -Dxray.execution.include.datetime=false -
Dxray.execution.include.counter=false -Dxray.execution.prefix="**LMR - Dry Run 3**" -
Dxray.project.key=DRXTMAPP

| | | |
|---------------|---|--------|
| DRXTMAPP-424 | Medication24 - SRS 3952, 3691, 3982, 3736
Verifies user medication acknowledgement flow the Medication Screen correctly displays all current medications in their confirmed date with valid details. | FAIL |
| | | FAIL |
| DRXTMAPP-301 | Medication9 - SRS 3715, 3784
Verifies medication acknowledgement flow for medications triggered due to 7 Low readings. | FAIL |
| DRXTMAPP-2057 | Medication65 - SRS 4252 | PASS |
| DRXTMAPP-2056 | Medication64 - SRS 3824 | PASS |
| DRXTMAPP-4098 | Medication9_2 - SRS 3715,3784 | PASS |
| DRXTMAPP-531 | Notification17_1 - SRS 3819, 3903, 3967, 4296, 3751 | FAIL |
| DRXTMAPP-559 | Notification8_1 - SRS 3910 | PASS |
| DRXTMAPP-425 | Navigation4 - SRS 2918 | PASS |
| DRXTMAPP-798 | Medication45 - SRS 4161, 4213, 3817 | Passed |
| DRXTMAPP-429 | Navigation5 - SRS 3730, 2918 | FAIL |
| DRXTMAPP-431 | Navigation6 - SRS 3730, 2918 | |
| DRXTMAPP-4875 | Medication68 - SRS 3787, 4417, 3914, 3818, 4620, 4621, 4609 | FAIL |

Links to Important Docs

Monday, July 20, 2026 11:05 AM

CG Android framework setup:

CG Mac framework setup: [CG Test Framework Setup Guide \(Macbook\).docx](#)

Empty jira bug

[\[DRXTMAPP-5776\] Env Q3 - No Error Message Displayed When User Fails the "Call Me" MFA - Jira](#)

Personal mainline testing doc: [https://abbott-](https://abbott-my.sharepoint.com/:w:/r/personal/samhitha_bodangi_abbott_com/Documents/KT%20Documents/Test_Case_Verification.docx?d=w8b9f8fc392ef498a9e167ee13d4ff11f&csf=1&web=1&e=nd9hUG)

[my.sharepoint.com/:w:/r/personal/samhitha_bodangi_abbott_com/Documents/KT%20Documents/Test_Case_Verification.docx?d=w8b9f8fc392ef498a9e167ee13d4ff11f&csf=1&web=1&e=nd9hUG](https://abbott-my.sharepoint.com/:w:/r/personal/samhitha_bodangi_abbott_com/Documents/KT%20Documents/Test_Case_Verification.docx?d=w8b9f8fc392ef498a9e167ee13d4ff11f&csf=1&web=1&e=nd9hUG)

```
// Step 2: Find requested key from configured key column.
// Fallback for StringLabels exports where key may be in either Column B(index 1) or Column C(index 2).
List<Row> matchingRows = new ArrayList<>();
LinkedHashSet<Integer> keyColumnCandidates = new LinkedHashSet<>();
keyColumnCandidates.add(keyColumn);
keyColumnCandidates.add(1);
keyColumnCandidates.add(2);
for (int keyColumnCandidate : keyColumnCandidates) {
    matchingRows = findRowsByExactKey(sheet, formatter, keyColumnCandidate, normalizedKey);
    if (!matchingRows.isEmpty()) {
        if (keyColumnCandidate != keyColumn) {
            LogUtil.warn("[StringLabel][WARN] Key " + normalizedKey + " not found using configured key
column index " + keyColumn + ". Using fallback key column index " + keyColumnCandidate + ".");
        }
        break;
    }
}
LogUtil.info("[StringLabel] Matching rows found for key " + normalizedKey + ": " + matchingRows.size());
```

```
// Step 2: Find requested key from Column B and select exact project row.
List<Row> matchingRows = findRowsByExactKey(sheet, formatter, keyColumn, normalizedKey);
LogUtil.info("[StringLabel] Matching rows found for key " + normalizedKey + ": " + matchingRows.size());
```

```
# Clinic-specific overrides (patched by ConfigReader.get("clinic", "US"))
merlinClinicUserName.us=ust1swvt1
merlinClinicPassword.us=Test!1234
merlinClinicUserName.nonus=krt1swvt1
merlinClinicPassword.nonus=Test!1234
```

T1/webConstants.properties

```
# Clinic-specific overrides (patched by ConfigReader.get("clinic", "US"))
merlinClinicUserName.us=ust1swvt1
merlinClinicPassword.us=Test!1234
merlinClinicUserName.nonus=krt1swvt1
merlinClinicPassword.nonus=Test!1234
```

```
public String getPatientProfileSettingsPageTreatingPhysicianOptionClinicAware() {
    // Environment + clinic aware treating physician option selection
    String environment = ConfigReader.get("environment", "q3").trim().toLowerCase();
    String clinic = ConfigReader.get("clinic", "US").trim().toLowerCase();
    LogUtil.debug("Getting treating physician locator for environment: " + environment + ", clinic: " + clinic);

    if ("nonus".equals(clinic) || "t1".equals(environment)) {
        return "//span[normalize-space()='Lee, Danny']";
    } else {
        return getPatientProfileSettingsPageTreatingPhysicianOption();
    }
}
```

Precon 4873, sst4ep 1

- Shoul db e "LOG In" (space missing)
- Pre con 582 step 3
 - o Add text to go to active or scheduled section
 - o Is step 6 necessary?

List of Tasks

Wednesday, August 12, 2026

1:14 PM

Mainline testing

Update notification tcs

Update weight8_2 tc

Update contact tc (name change)

Manual execution

Documentation

1. T1 Environment Values: Postman Collection UID:- 32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316
Postman CLI Environment ID: 41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f
Register API:- 32592187-906e0b7c-ae6b-4541-8c3f-4f7d21f64e80
2. Q3 Environment Values: Postman Collection UID:- 52864576-d9174d04-011e-4d5f-aca9-0001b00910e6
Postman CLI Environment ID: 31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85
Register API:- 52864576-86b8eeaa-9e22-47f2-a54f-e09cddf69975
1. postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i "32592187-906e0b7c-ae6b-4541-8c3f-4f7d21f64e80" --timeout-request 30000 --env-var "_cli=1" --env-var "sensor_serial=1B1128" --env-var "patient_dob=2000-01-01" --env-var "patient_name=Homescreen13 iOS" --env-var "email=homescreen13ios@gmail.com" --env-var "password=Test!123" --env-var "phone="+17813230226" --verbose

Android: <https://zhrnd.atlassian.net/browse/DRXTMAPP-6244>

Ios: <https://zhrnd.atlassian.net/browse/DRXTMAPP-6247>

Links to KT Documents

Monday, August 17, 2026 10:23 PM

Setting up PTC: [KT_PTC.docx](#)

CardioGuide Test Automation NPS Documentation and Test cases:
[XXXXXXX CardioGuide Test Automation NPS Documentation.docx](#)

Keyboard Commands

Tuesday, June 30, 2026 3:27 PM

right-click the file in Finder, press and hold the **Option (⌘)** key, and select "**Copy [File Name] as Pathname**". Then, press **Command (⌘) + V** to paste the path wherever you need it. [[1](#), [2](#), [3](#)]

IntelliJ file search (Windows) = Ctrl+Shift+N

IntelliJ file search (Mac) = Cmd+Shift+O

IntelliJ multiple line comment (Mac) = Cmd + Options + /

IntelliJ multiple line comment (Windows) = Ctrl + Shift + /

In Finder, do Shift + Cmd + G to paste and go to folder path

Screenshots

Windows: prt sc or Windows + Shift + S

Mac = Cmd + Shift + 3 (capture whole screen)

Cmd + Shift + 4 (capture certain part of screen)

Cmd + Shift + 5 (screen recording)

Final Checklist

Monday, August 17, 2026 10:11 PM

| Test Case | iOS | Android |
|-----------------|------|---------|
| Homescreen2 | Pass | Pass |
| More_Weight8_2 | Pass | Pass |
| More_Support2 | Pass | Pass |
| More_Support5 | Pass | Pass |
| Notification9_5 | | Pass |
| Notification9_6 | | Pass |
| Notification9_7 | | |
| Notification9_8 | | |
| Registration13 | | |

Cell Logs Analysis

Friday, June 19, 2026 1:14 PM

~~P260400055~~

~~P260400015~~

P260400054

~~P260300386~~

~~P260300406~~

P260300455

P260300093

P260400055

Friday, June 19, 2026 1:19 PM

2026-06-08
Sent via WiFi, not Cellular

Was able to see multiple carriers, but could not connect to any of them

```
INFO |2026-06-08 14:56:00.061|interrogatorServer[878]|file|(getTimestampHandler.cpp)|line|(65)|Msg|(Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR |2026-06-08 14:56:00.199|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO |2026-06-08 14:56:00.204|interrogatorServer[878]|file|(getTimestampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and ppp0 ip address: (null))
INFO |2026-06-08 14:56:00.456|interrogatorServer[878]|file|(getTimestampHandler.cpp)|line|(69)|Msg|(Kernel IP routing table
INFO |2026-06-08 14:56:00.459|interrogatorServer[878]|file|(getTimestampHandler.cpp)|line|(74)|Msg|(Received timeStamp message reply; -1)
INFO |2026-06-08 14:56:00.546|cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: T-Mobile status available)
INFO |2026-06-08 14:56:00.560|cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: AT&T status available)
ERROR |2026-06-08 14:56:00.625|interrogatorServer[878]|file|(getTimestampHandler.cpp)|line|(101)|Msg|(getTimestamp Failed)
INFO |2026-06-08 14:56:00.639|cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: 313 100 status available)
INFO |2026-06-08 14:56:00.685|cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: Verizon status available)
INFO |2026-06-08 14:56:00.826|cellServer[535]|file|(cellMgrState.c)|line|(394)|Msg|(Got carrier list!)
INFO |2026-06-08 14:56:00.891|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Switch Carriers)
INFO |2026-06-08 14:56:01.017|cellServer[535]|file|(cellMgrState.c)|line|(427)|Msg|(Attempt to connect to carrier /jacs_0/operator/310260)
ERROR |2026-06-08 14:56:01.674|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:01.753|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:01.799|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:01.848|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:01.901|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:01.903|cellServer[535]|file|(cellMgrState.c)|line|(433)|Msg|(Carrier switch failed!)
INFO |2026-06-08 14:56:01.927|cellServer[535]|file|(cellMgrState.c)|line|(175)|Msg|(Cell powered down)
INFO |2026-06-08 14:56:01.927|cellServer[535]|file|(cellMgrState.c)|line|(175)|Msg|(Cell powered down)
INFO |2026-06-08 14:56:01.943|cellServer[535]|file|(cellMgrState.c)|line|(427)|Msg|(Attempt to connect to carrier /jacs_0/operator/310410)
ERROR |2026-06-08 14:56:01.995|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310410
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:02.042|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310410
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:02.109|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310410
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:02.171|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310410
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:02.204|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310410
returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 14:56:02.227|cellServer[535]|file|(cellMgrState.c)|line|(433)|Msg|(Carrier switch failed!)
```

After 4 attempts to switch carriers, tried to find new carrier list. Now, can't even see any carriers

```
INFO |2026-06-08 14:56:03.118|interrogatorServer[878]|file|(getTimestampHandler.cpp)|line|(104)|Msg|(Cellular output ppp0: error fetching interface information: Device not found)
INFO |2026-06-08 14:56:03.129|interrogatorServer[878]|file|(getTimestampHandler.cpp)|line|(105)|Msg|(wlan0 output wlan0: flags=4098<BROADCAST,MULTICAST> mtu 1500)
INFO |2026-06-08 14:56:03.160|watchdogServer[1072]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (30030) is less than 15 mins)
INFO |2026-06-08 14:56:03.190|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Update Carrier List)
ALERT |2026-06-08 14:56:03.215|cellServer[535]|file|(cellMgrState.c)|line|(378)|Msg|(Default STATE_FAULT)
ALERT |2026-06-08 14:56:03.261|cellServer[535]|file|(cellMgrState.c)|line|(386)|Msg|(Query for new carrier list...)WARNING |2026-06-08 14:56:03.340
|cellServer[535]|file|(cmsgwCell.c)|line|(478)|Msg|(Failed to issue dbus call on ofonoDBus, returned -53, error: Method "Scan" with signature "" on interface "org.ofono.NetworkRegistration" doesn't
exist
ERROR |2026-06-08 14:56:03.396|cellServer[535]|file|(cellMgrState.c)|line|(398)|Msg|(Unable to get carrier list!)
INFO |2026-06-08 14:56:03.413|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: FAULT!)
```

Multiple failures to get carrier list

```
INFO |2026-06-08 14:57:38.750|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
WARNING |2026-06-08 14:57:38.845|dataTransferService[1062]|file|(cmsgwHttp.c)|line|(205)|Msg|(curl_easy_perform() failed: Timeout was reached)
WARNING |2026-06-08 14:57:38.847|dataTransferService[1062]|file|(inetServiceWorker.c)|line|(181)|Msg|(FAILURE: Get TimeStamp (-1))
ERROR |2026-06-08 14:57:38.850|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(591)|Msg|(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
ERROR |2026-06-08 14:57:38.851|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=wlan0)
INFO |2026-06-08 14:57:38.852|interrogatorServer[878]|file|(getTimestampHandler.cpp)|line|(65)|Msg|(Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR |2026-06-08 14:57:39.102|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
```

User Switches on WiFi, cellular trying to connect in the background, finds carriers but cannot connect to any of them

```
INFO |2026-06-08 14:59:18.427|watchdogServer[1072]|file|(watchdogSupport.h)|line|(342)|Msg|(Data received on shared memory wd_reset)
INFO |2026-06-08 14:59:18.429|watchdogServer[1072]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (0) is less than 15 mins)
INFO |2026-06-08 14:59:19.645|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/dialogs/DialogStackElement.qml - 380 - noConnectionDialog Connection Settings
ButtonClicked)
INFO |2026-06-08 14:59:19.651|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1458)|Msg|(GUI: Show dialog 'None')
INFO |2026-06-08 14:59:20.568|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SettingConnectionPageElement.qml - 107 - Wifi Menu Item Clicked)
INFO |2026-06-08 14:59:21.490|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SetupWifiPageElement.qml - 85 - Wifi ON/OFF Button Clicked : ON)
INFO |2026-06-08 14:59:21.506|CardioMEMS_GuiApp[896]|file|(NetworkInfoController.cpp)|line|(228)|Msg|(Wifi is on ? 0)
INFO |2026-06-08 14:59:21.511|CardioMEMS_GuiApp[896]|file|(NetworkInfoController.cpp)|line|(234)|Msg|(Powering on Wifi)
INFO |2026-06-08 14:59:21.555|CardioMEMS_GuiApp[896]|file|(NetworkInfoController.cpp)|line|(270)|Msg|(NetworkInfoController: initiating wifi scan)
INFO |2026-06-08 14:59:22.577|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1222)|Msg|(GuiController::enable Power off button)
INFO |2026-06-08 14:59:31.578|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SetupWifiPageElement.qml - 221 - Refresh Wifi called)
```

INFO |2026-06-08 14:59:31.831| CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(wifiOnOffRefreshTimer triggered)
INFO |2026-06-08 14:59:33.440| watchdogServer[1072]|file|(watchdogSupport.h)|line|(342)|Msg|(Data received on shared memory wd_reset)
INFO |2026-06-08 14:59:33.441| watchdogServer[1072]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (0) is less than 15 mins)
INFO |2026-06-08 14:59:36.370| CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SetupWifiPageElement.qml - 145 - Available Network Item Clicked : iPhone)
INFO |2026-06-08 14:59:42.800| CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SettingWifiJoinPageElement.qml - 173 - Password Show/Hide Button Clicked)
INFO |2026-06-08 14:59:48.475| watchdogServer[1072]|file|(watchdogSupport.h)|line|(342)|Msg|(Data received on shared memory wd_reset)
INFO |2026-06-08 14:59:48.495| watchdogServer[1072]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (0) is less than 15 mins)
INFO |2026-06-08 14:59:57.975| CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SettingWifiJoinPageElement.qml - 277 - Join Button Clicked)
INFO |2026-06-08 15:00:02.423| cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: T-Mobile status available)
INFO |2026-06-08 15:00:02.470| cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: AT&T status available)
INFO |2026-06-08 15:00:02.566| cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: Verizon status available)
INFO |2026-06-08 15:00:02.690| cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: 313 100 status available)
INFO |2026-06-08 15:00:02.793| cellServer[535]|file|(cellMgrState.c)|line|(394)|Msg|(Got carrier list!)
INFO |2026-06-08 15:00:02.847| cellServer[535]|file|(cellMgrState.c)|line|(175)|Msg|(Cell powered down)
INFO |2026-06-08 15:00:02.924| cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Switch Carriers)
INFO |2026-06-08 15:00:02.972| cellServer[535]|file|(cellMgrState.c)|line|(427)|Msg|(Attempt to connect to carrier /jacs_0/operator/310260)
INFO |2026-06-08 15:00:03.503| watchdogServer[1072]|file|(watchdogSupport.h)|line|(342)|Msg|(Data received on shared memory wd_reset)
INFO |2026-06-08 15:00:03.504| watchdogServer[1072]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (0) is less than 15 mins)
ERROR |2026-06-08 15:00:03.719| cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260 returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 15:00:03.787| cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260 returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 15:00:03.871| cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260 returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 15:00:03.996| cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260 returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 15:00:04.091| cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260 returned -53: Method "Register" with signature "" on interface "org.ofono.NetworkOperator" doesn't exist
ERROR |2026-06-08 15:00:04.092| cellServer[535]|file|(cellMgrState.c)|line|(433)|Msg|(Carrier switch failed!)

Device works on WiFi, received an IP address via WiFi

(configureHttpDefaults: configureHttpDefaults success.)

INFO |2026-06-08 15:00:49.728| dataTransferService[1062]|file|(cmemsvrIf.c)|line|(549)|Msg|(doGet:549: Attempting get of <https://papi.g1.cardiomemshf.com/v6/timestamp to file /tmp/getTimeStamp.rsp>)

INFO |2026-06-08 15:00:49.759| interrogatorServer[878]|file|(getTimeStamp.cpp)|line|(53)|Msg|([38;5;27;1mSending TimeStamp Msg to Backend [0m)
INFO |2026-06-08 15:00:51.665| cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-08 15:00:53.813| dataTransferService[1062]|file|(inetServiceWorker.c)|line|(189)|Msg|(SUCCESS: Get TimeStamp)
ERROR |2026-06-08 15:00:53.822| interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(591)|Msg|(cannot find ppp0 operstate file, /sys/class/net/ ppp0/operstate)
INFO |2026-06-08 15:00:53.829| interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(65)|Msg|(Wifi enabled? TRUE and wlan0 ip address: 172.20.10.4)
ERROR |2026-06-08 15:00:54.162| interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO |2026-06-08 15:00:54.174| interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and ppp0 ip address: (null))
INFO |2026-06-08 15:00:54.464| interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(69)|Msg|(Kernel IP routing table
INFO |2026-06-08 15:00:54.477| interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(74)|Msg|(Received timeStamp message reply; 200)
INFO |2026-06-08 15:00:54.483| interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(78)|Msg|(Get TimeStamp Successful)

Still tries to connect to cellular, but fails

INFO |2026-06-08 15:01:05.189| cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Switch Carriers)
INFO |2026-06-08 15:01:05.206| cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Update Carrier List)
ALERT |2026-06-08 15:01:05.231| cellServer[535]|file|(cellMgrState.c)|line|(378)|Msg|(Default STATE_FAULT)
ALERT |2026-06-08 15:01:05.248| cellServer[535]|file|(cellMgrState.c)|line|(386)|Msg|(Query for new carrier list...)
WARNING |2026-06-08 15:01:05.286| cellServer[535]|file|(cmsgwCell.c)|line|(478)|Msg|(Failed to issue dbus call on ofonoDbus, returned -53, error: Method "Scan" with signature "" on interface "org.ofono.NetworkRegistration" doesn't exist
ERROR |2026-06-08 15:01:05.329| cellServer[535]|file|(cellMgrState.c)|line|(398)|Msg|(Unable to get carrier list!)
INFO |2026-06-08 15:01:05.335| cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: FAULTI)
ALERT |2026-06-08 15:01:05.357| cellServer[535]|file|(cellMgrState.c)|line|(580)|Msg|(Setting to STATE_FAULT)
ALERT |2026-06-08 15:01:05.391| cellServer[535]|file|(cellMgrState.c)|line|(580)|Msg|(Setting to STATE_FAULT)
ALERT |2026-06-08 15:01:05.398| cellServer[535]|file|(cellMgrState.c)|line|(580)|Msg|(Setting to STATE_FAULT)

Was able to connect to cellular, but disconnects after

INFO | 2026-06-06 15:27:01.810 | interrogatorServer[900] | file | (networkSettings.c) | line | (111) | Msg | (networkConfigWiFi file /usr/cardiomems/var/config/uiConfig.json)
ERROR | 2026-06-06 15:27:01.901 | interrogatorServer[900] | file | (sendData.cpp) | line | (469) | Msg | (The wifi or cell is coming back as enabled in check connections function)
ERROR | 2026-06-06 15:27:01.971 | interrogatorServer[900] | file | (sendData.cpp) | line | (472) | Msg | (check connections return value: 1)
ERROR | 2026-06-06 15:27:02.059 | interrogatorServer[900] | file | (sendData.cpp) | line | (473) | Msg | (check connections wifi results: 0)
ERROR | 2026-06-06 15:27:02.060 | interrogatorServer[900] | file | (sendData.cpp) | line | (474) | Msg | (check connections cell results: 1)
INFO | 2026-06-06 15:27:02.144 | interrogatorServer[900] | file | (sendData.cpp) | line | (165) | Msg | (send geTimestamp for send readings)
INFO | 2026-06-06 15:27:02.227 | interrogatorServer[900] | file | (sendData.cpp) | line | (166) | Msg | (checked for connections and now will send gettimestamp)
INFO | 2026-06-06 15:27:02.229 | interrogatorServer[900] | file | (getTimeStamp.cpp) | line | (45) | Msg | ([38;5;27;1mGetTimeStampEv event: GetTimeStamp [0m)
INFO | 2026-06-06 15:27:02.331 | dataTransferService[1065] | file | (cmemsSrvf.c) | line | (682) | Msg | (configureHttpDefaults: configureHttpDefaults success.)

INFO | 2026-06-06 15:27:16.637 | CardioMEMS_GuiApp[904] | file | (GuiController.cpp) | line | (896) | Msg | (Turning on Wifi and cell)
INFO | 2026-06-06 15:27:16.638 | CardioMEMS_GuiApp[904] | file | (NetworkInfoController.cpp) | line | (228) | Msg | (Wifi is on ? 0)
INFO | 2026-06-06 15:27:16.641 | CardioMEMS_GuiApp[904] | file | (GuiController.cpp) | line | (2266) | Msg | (Turning on cell)
INFO | 2026-06-06 15:27:16.732 | CardioMEMS_GuiApp[904] | file | (cmswNetwork.c) | line | (1471) | Msg | (Success to issue dbus method call SetProperty on net.connman.Manager: /net/connman/technology/cellular)
INFO | 2026-06-06 15:27:24.921 | watchdogServer[1073] | file | (watchdogSupport.h) | line | (366) | Msg | (Inactivity time (15012) is less than 15 mins)
WARNING | 2026-06-06 15:27:32.430 | dataTransferService[1065] | file | (cmswHttp.c) | line | (205) | Msg | (curl_easy_perform() failed: Timeout was reached)
WARNING | 2026-06-06 15:27:32.436 | dataTransferService[1065] | file | (inetServiceWorker.c) | line | (181) | Msg | (FAILURE: Get TimeStamp (-1))
ERROR | 2026-06-06 15:27:32.476 | interrogatorServer[900] | file | (cmswNetwork.c) | line | (591) | Msg | (cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
ERROR | 2026-06-06 15:27:32.478 | interrogatorServer[900] | file | (cmswNetwork.c) | line | (531) | Msg | (Error with ioctl call, iface=wlan0)
INFO | 2026-06-06 15:27:32.529 | interrogatorServer[900] | file | (getTimeStampHandler.cpp) | line | (65) | Msg | (Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR | 2026-06-06 15:27:32.751 | interrogatorServer[900] | file | (cmswNetwork.c) | line | (531) | Msg | (Error with ioctl call, iface=ppp0)
INFO | 2026-06-06 15:27:32.764 | interrogatorServer[900] | file | (getTimeStampHandler.cpp) | line | (66) | Msg | (Cell enabled? FALSE and ppp0 ip address: (null))
INFO | 2026-06-06 15:27:32.994 | interrogatorServer[900] | file | (getTimeStampHandler.cpp) | line | (69) | Msg | (Kernel IP routing table
INFO | 2026-06-06 15:27:33.005 | interrogatorServer[900] | file | (getTimeStampHandler.cpp) | line | (74) | Msg | (Received timeStamp message reply; -1)
ERROR | 2026-06-06 15:27:33.007 | interrogatorServer[900] | file | (getTimeStampHandler.cpp) | line | (101) | Msg | (getTimeStamp Failed)
INFO | 2026-06-06 15:27:33.616 | interrogatorServer[900] | file | (getTimeStampHandler.cpp) | line | (104) | Msg | (Cellular output ppp0: error fetching interface information: Device not found

User turns on WiFi

INFO | 2026-06-06 15:29:19.573 | CardioMEMS_GuiApp[904] | file | (GuiController.cpp) | line | (1416) | Msg | (qrc:/dialogs/DialogStackElement.qml - 380 - noConnectionDialog Connection Settings ButtonClicked)
INFO | 2026-06-06 15:29:19.581 | CardioMEMS_GuiApp[904] | file | (GuiController.cpp) | line | (1458) | Msg | (GUI: Show dialog 'None')
INFO | 2026-06-06 15:29:25.098 | watchdogServer[1073] | file | (watchdogSupport.h) | line | (342) | Msg | (Data received on shared memory wd_reset)
INFO | 2026-06-06 15:29:25.099 | watchdogServer[1073] | file | (watchdogSupport.h) | line | (366) | Msg | (Inactivity time (0) is less than 15 mins)
INFO | 2026-06-06 15:29:25.900 | CardioMEMS_GuiApp[904] | file | (GuiController.cpp) | line | (1416) | Msg | (qrc:/pages/SettingConnectionPageElement.qml - 107 - Wifi Menu Item Clicked)
INFO | 2026-06-06 15:29:28.041 | CardioMEMS_GuiApp[904] | file | (GuiController.cpp) | line | (1416) | Msg | (qrc:/pages/SetupWifiPageElement.qml - 85 - Wifi ON/OFF Button Clicked : ON)
INFO | 2026-06-06 15:29:28.057 | CardioMEMS_GuiApp[904] | file | (NetworkInfoController.cpp) | line | (228) | Msg | (Wifi is on ? 0)
INFO | 2026-06-06 15:29:28.061 | CardioMEMS_GuiApp[904] | file | (NetworkInfoController.cpp) | line | (234) | Msg | (Powering on Wifi)
INFO | 2026-06-06 15:29:28.097 | CardioMEMS_GuiApp[904] | file | (NetworkInfoController.cpp) | line | (270) | Msg | (NetworkInfoController: initiating

wifi scan)

INFO | 2026-06-06 15:29:29.138 | CardioMEMS_GuiApp[904] | file |(GuiController.cpp)| line |(1222)| Msg |(GuiController::enable Power off button)
INFO | 2026-06-06 15:29:38.124 | CardioMEMS_GuiApp[904] | file |(GuiController.cpp)| line |(1416)| Msg |(qrc:/pages/SetupWifiPageElement.qml -
221 - Refresh Wifi called)
INFO | 2026-06-06 15:29:38.988 | CardioMEMS_GuiApp[904] | file |(GuiController.cpp)| line |(1416)| Msg |(wifiOnOffRefreshTimer triggered)
INFO | 2026-06-06 15:29:40.116 | watchdogServer[1073] | file |(watchdogSupport.h)| line |(342)| Msg |(Data received on shared memory wd_reset)

Cellular is still trying in the background. Found carrier, but unable to switch

INFO | 2026-06-06 15:30:15.639 | cellServer[539] | file |(cellMgrState.c)| line |(427)| Msg |(Attempt to connect to carrier /jacs_0/operator/310410)
ERROR | 2026-06-06 15:30:15.661 | cellServer[539] | file |(cmswCell.c)| line |(386)| Msg |(Failed to issue dbus method call Register on
org.ofono.NetworkOperator path /jacs_0/operator/310410 returned -53: Method "Register" with signature "" on interface
"org.ofono.NetworkOperator" doesn't exist
ERROR | 2026-06-06 15:30:15.687 | cellServer[539] | file |(cmswCell.c)| line |(386)| Msg |(Failed to issue dbus method call Register on
org.ofono.NetworkOperator path /jacs_0/operator/310410 returned -53: Method "Register" with signature "" on interface
"org.ofono.NetworkOperator" doesn't exist
ERROR | 2026-06-06 15:30:15.703 | cellServer[539] | file |(cmswCell.c)| line |(386)| Msg |(Failed to issue dbus method call Register on
org.ofono.NetworkOperator path /jacs_0/operator/310410 returned -53: Method "Register" with signature "" on interface
"org.ofono.NetworkOperator" doesn't exist
ERROR | 2026-06-06 15:30:15.716 | cellServer[539] | file |(cmswCell.c)| line |(386)| Msg |(Failed to issue dbus method call Register on
org.ofono.NetworkOperator path /jacs_0/operator/310410 returned -53: Method "Register" with signature "" on interface
"org.ofono.NetworkOperator" doesn't exist
ERROR | 2026-06-06 15:30:15.755 | cellServer[539] | file |(cmswCell.c)| line |(386)| Msg |(Failed to issue dbus method call Register on
org.ofono.NetworkOperator path /jacs_0/operator/310410 returned -53: Method "Register" with signature "" on interface
"org.ofono.NetworkOperator" doesn't exist
ERROR | 2026-06-06 15:30:15.760 | cellServer[539] | file |(cellMgrState.c)| line |(433)| Msg |(Carrier switch failed!)
INFO | 2026-06-06 15:30:15.762 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Check For Cell)
INFO | 2026-06-06 15:30:15.767 | cellServer[539] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-06 15:30:15.823 | cellServer[539] | file |(cellMgrState.c)| line |(311)| Msg |(SUCCESS: Cell Carrier: Service Name: cellular_
262017840798815_context1 State: idle)
INFO | 2026-06-06 15:30:15.828 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Switch Carriers)
INFO | 2026-06-06 15:30:15.834 | cellServer[539] | file |(cellMgrState.c)| line |(427)| Msg |(Attempt to connect to carrier /jacs_
0/operator/313100)
ERROR | 2026-06-06 15:30:15.847 | cellServer[539] | file |(cmswCell.c)| line |(386)| Msg |(Failed to issue dbus method call Register on
org.ofono.NetworkOperator path /jacs_0/operator/313100 returned -53: Method "Register" with signature "" on interface
"org.ofono.NetworkOperator" doesn't exist
ERROR | 2026-06-06 15:30:15.882 | cellServer[539] | file |(cmswCell.c)| line |(386)| Msg |(Failed to issue dbus method call Register on
org.ofono.NetworkOperator path /jacs_0/operator/313100 returned -53: Method "Register" with signature "" on interface
"org.ofono.NetworkOperator" doesn't exist

Again, found carriers, but can't connect

INFO | 2026-06-06 15:31:35.292 | cellServer[539] | file |(cmswCell.c)| line |(331)| Msg |(Carrier found: T-Mobile status available)
INFO | 2026-06-06 15:31:35.370 | cellServer[539] | file |(cmswCell.c)| line |(331)| Msg |(Carrier found: AT&T status available)
INFO | 2026-06-06 15:31:35.435 | cellServer[539] | file |(cmswCell.c)| line |(331)| Msg |(Carrier found: 313 100 status available)
INFO | 2026-06-06 15:31:35.539 | cellServer[539] | file |(cellMgrState.c)| line |(394)| Msg |(Got carrier list!)
INFO | 2026-06-06 15:31:35.638 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Cell Never Connected)
INFO | 2026-06-06 15:31:35.719 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Switch Carriers)
INFO | 2026-06-06 15:31:35.799 | cellServer[539] | file |(cellMgrState.c)| line |(427)| Msg |(Attempt to connect to carrier /jacs_0/operator/310260)
INFO | 2026-06-06 15:31:40.245 | watchdogServer[1073] | file |(watchdogSupport.h)| line |(342)| Msg |(Data received on shared memory wd_reset)
INFO | 2026-06-06 15:31:40.246 | watchdogServer[1073] | file |(watchdogSupport.h)| line |(366)| Msg |(Inactivity time (0) is less than 15 mins)
INFO | 2026-06-06 15:31:46.857 | CardioMEMS_GuiApp[904] | file |(NetworkDataModel.cpp)| line |(402)| Msg |(m_wifilsConnected changed to: 1)

But for some reason, says was able to connect

INFO | 2026-06-06 15:31:46.862 | CardioMEMS_GuiApp[904] | file |(GuiController.cpp)| line |(1222)| Msg |(GuiController::enable Power off button)
INFO | 2026-06-06 15:31:55.269 | watchdogServer[1073] | file |(watchdogSupport.h)| line |(366)| Msg |(Inactivity time (15022) is less than 15 mins)
INFO | 2026-06-06 15:31:58.011 | CardioMEMS_GuiApp[904] | file |(GuiController.cpp)| line |(1416)| Msg |(qrc:/elements/HeaderElement.qml -
118 - Connection return Button clicked)
INFO | 2026-06-06 15:31:58.026 | CardioMEMS_GuiApp[904] | file |(GuiController.cpp)| line |(1218)| Msg |(GuiController::disable Power off for
SettingsSavedReadings)
INFO | 2026-06-06 15:31:58.029 | interrogatorServer[900] | file |(sendData.cpp)| line |(82)| Msg |([38;5;27;1mStartSendData event: SendData
[0m)
INFO | 2026-06-06 15:31:58.030 | interrogatorServer[900] | file |(sendData.cpp)| line |(138)| Msg |(Either send data flow just started or response
received for other commands.)

INFO | 2026-06-06 15:31:58.131 | interrogatorServer[900] | file |(networkSettings.c) | line |(111) | Msg |(networkConfigWiFi file /usr/cardiomems/var/config/uiConfig.json)

INFO | 2026-06-06 15:31:58.179 | interrogatorServer[900] | file |(networkSettings.c) | line |(116) | Msg |(networkConfigWiFi Child Exist)

ERROR | 2026-06-06 15:31:58.206 | interrogatorServer[900] | file |(sendData.cpp) | line |(469) | Msg |(The wifi or cell is coming back as enabled in check connections function)

ERROR | 2026-06-06 15:31:58.256 | interrogatorServer[900] | file |(sendData.cpp) | line |(472) | Msg |(check connections return value: 1)

ERROR | 2026-06-06 15:31:58.305 | interrogatorServer[900] | file |(sendData.cpp) | line |(473) | Msg |(check connections wifi results: 1)

ERROR | 2026-06-06 15:31:58.324 | interrogatorServer[900] | file |(sendData.cpp) | line |(474) | Msg |(check connections cell results: 1)

INFO | 2026-06-06 15:31:58.467 | dataTransferService[1065] | file |(cmemsSrvIf.c) | line |(549) | Msg |(doGet:549: Attempting get of <https://papi.g1.cardiomemshf.com/v6/timestamp to file /tmp/getTimeStamp.rsp>)

INFO | 2026-06-06 15:31:58.528 | interrogatorServer[900] | file |(getTimeStamp.cpp) | line |(53) | Msg |([38;5;27;1mSending TimeStamp Msg to Backend [0m)

INFO | 2026-06-06 15:32:01.507 | dataTransferService[1065] | file |(inetServiceWorker.c) | line |(189) | Msg |(SUCCESS: Get TimeStamp)

ERROR | 2026-06-06 15:32:01.538 | interrogatorServer[900] | file |(cmswNetwork.c) | line |(591) | Msg |(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)

INFO | 2026-06-06 15:32:01.553 | interrogatorServer[900] | file |(getTimeStampHandler.cpp) | line |(65) | Msg |(Wifi enabled? TRUE and wlan0 ip address: 192.168.1.212)

ERROR | 2026-06-06 15:32:01.885 | interrogatorServer[900] | file |(cmswNetwork.c) | line |(531) | Msg |(Error with ioctl call, iface=ppp0)

INFO | 2026-06-06 15:32:01.907 | interrogatorServer[900] | file |(getTimeStampHandler.cpp) | line |(66) | Msg |(Cell enabled? FALSE and ppp0 ip address: (null))

INFO | 2026-06-06 15:32:02.143 | interrogatorServer[900] | file |(getTimeStampHandler.cpp) | line |(69) | Msg |(Kernel IP routing table

INFO | 2026-06-06 15:32:02.148 | interrogatorServer[900] | file |(getTimeStampHandler.cpp) | line |(74) | Msg |(Received timeStamp message reply; 200)

INFO | 2026-06-06 15:32:02.167 | interrogatorServer[900] | file |(getTimeStampHandler.cpp) | line |(78) | Msg |(Get TimeStamp Successful)

INFO | 2026-06-06 15:32:02.180 | interrogatorServer[900] | file |(timeSynchHandler.cpp) | line |(80) | Msg |(synchronize the time)

INFO | 2026-06-06 15:32:02.191 | interrogatorServer[900] | file |(timeSynchHandler.cpp) | line |(88) | Msg |(timestamp is 2026-06-06T15:32:04Z)

Failed to connect to cellular. Tries multiple times before "no connection screen" is shown

```
INFO |2026-06-05 15:03:53.503|CardioMEMS_GuiApp[918]|file|(GuiController.cpp)|line|(896)|Msg|(Turning on Wifi and cell)
INFO |2026-06-05 15:03:53.504|CardioMEMS_GuiApp[918]|file|(NetworkInfoController.cpp)|line|(228)|Msg|(Wifi is on ? 0)
INFO |2026-06-05 15:03:53.507|CardioMEMS_GuiApp[918]|file|(GuiController.cpp)|line|(2266)|Msg|(Turning on cell)
INFO |2026-06-05 15:03:53.595|CardioMEMS_GuiApp[918]|file|(cmsgwNetwork.c)|line|(1471)|Msg|(Success to issue dbus method call SetProperty on net.connman.Manager:
/net/connman/technology/cellular)
INFO |2026-06-05 15:04:04.364|watchdogServer[1088]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (15008) is less than 15 mins)
WARNING |2026-06-05 15:04:08.728|dataTransferService[1079]|file|(cmsgwHttp.c)|line|(205)|Msg|(curl_easy_perform() failed: Timeout was reached)
WARNING |2026-06-05 15:04:08.729|dataTransferService[1079]|file|(inetServiceWorker.c)|line|(181)|Msg|(FAILURE: Get TimeStamp (-1))
ERROR |2026-06-05 15:04:08.731|interrogatorServer[914]|file|(cmsgwNetwork.c)|line|(591)|Msg|(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
ERROR |2026-06-05 15:04:08.732|interrogatorServer[914]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=wlan0)
INFO |2026-06-05 15:04:08.733|interrogatorServer[914]|file|(getTimestampHandler.cpp)|line|(65)|Msg|(Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR |2026-06-05 15:04:08.919|interrogatorServer[914]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO |2026-06-05 15:04:08.925|interrogatorServer[914]|file|(getTimestampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and ppp0 ip address: (null))
INFO |2026-06-05 15:04:09.435|interrogatorServer[914]|file|(getTimestampHandler.cpp)|line|(69)|Msg|(Kernel IP routing table
INFO |2026-06-05 15:04:09.465|interrogatorServer[914]|file|(getTimestampHandler.cpp)|line|(74)|Msg|(Received timeStamp message reply; -1)
ERROR |2026-06-05 15:04:09.470|interrogatorServer[914]|file|(getTimestampHandler.cpp)|line|(101)|Msg|(getTimestamp Failed)
INFO |2026-06-05 15:04:09.948|interrogatorServer[914]|file|(getTimestampHandler.cpp)|line|(104)|Msg|(Cellular output ppp0: error fetching interface information: Device not found
```

User turns on WiFi

```
INFO |2026-06-05 15:06:08.467|CardioMEMS_GuiApp[918]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SettingConnectionPageElement.qml - 107 - Wifi Menu Item Clicked)
INFO |2026-06-05 15:06:15.878|CardioMEMS_GuiApp[918]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SetupWifiPageElement.qml - 85 - Wifi ON/OFF Button Clicked : ON)
INFO |2026-06-05 15:06:15.899|CardioMEMS_GuiApp[918]|file|(NetworkInfoController.cpp)|line|(228)|Msg|(Wifi is on ? 0)
INFO |2026-06-05 15:06:15.904|CardioMEMS_GuiApp[918]|file|(NetworkInfoController.cpp)|line|(234)|Msg|(Powering on Wifi)
INFO |2026-06-05 15:06:15.953|CardioMEMS_GuiApp[918]|file|(NetworkInfoController.cpp)|line|(270)|Msg|(NetworkInfoController: initiating wifi scan)
INFO |2026-06-05 15:06:16.379|CardioMEMS_GuiApp[918]|file|(NetworkInfoController.cpp)|line|(81)|Msg|(NetworkController: wifi scan complete)
INFO |2026-06-05 15:06:16.976|CardioMEMS_GuiApp[918]|file|(GuiController.cpp)|line|(1222)|Msg|(GuiController::enable Power off button)
```

But device was able to connect to WiFi, but immediately disconnects. Was able to find carriers, but not connect. User then refreshes Wifi

```
INFO |2026-06-05 15:06:53.608|cellServer[539]|file|(cellMgrState.c)|line|(311)|Msg|(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840797883
_context1 State: idle)
INFO |2026-06-05 15:06:53.617|cellServer[539]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Connect Cell Service)
INFO |2026-06-05 15:06:53.624|cellServer[539]|file|(cellMgrState.c)|line|(347)|Msg|(Attempt 0 of 7 to connect to cell service: cellular_262017840797883_context1)
WARNING |2026-06-05 15:06:53.660|cellServer[539]|file|(cmsgwNetwork.c)|line|(1428)|Msg|(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)
INFO |2026-06-05 15:06:54.303|cellServer[539]|file|(cellMgrState.c)|line|(347)|Msg|(Attempt 1 of 7 to connect to cell service: cellular_262017840797883_context1)
WARNING |2026-06-05 15:06:55.874|cellServer[539]|file|(cmsgwNetwork.c)|line|(1428)|Msg|(Failed to issue dbus method call Connect on net.connman.Service: Input/output error)
INFO |2026-06-05 15:06:55.876|cellServer[539]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Wait for Connman Cell Online)
INFO |2026-06-05 15:06:55.879|cellServer[539]|file|(cellMgrState.c)|line|(148)|Msg|(Network state=true)
INFO |2026-06-05 15:06:55.889|cellServer[539]|file|(cellMgrState.c)|line|(148)|Msg|(Network state=false)
INFO |2026-06-05 15:06:55.898|cellServer[539]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Online)
ALERT |2026-06-05 15:06:55.904|cellServer[539]|file|(cellMgrState.c)|line|(490)|Msg|(NOW ONLINE!)
ALERT |2026-06-05 15:06:55.911|cellServer[539]|file|(cellMgrState.c)|line|(495)|Msg|(Failedonlineattempts:0,max:3
INFO |2026-06-05 15:06:55.920|cellServer[539]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Offline)
INFO |2026-06-05 15:06:55.936|cellServer[539]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Wait For Offline)
```

Was able to send reading via WiFi

```
ERROR |2026-06-05 15:18:09.581|interrogatorServer[914]|file|(sendData.cpp)|line|(469)|Msg|(The wifi or cell is coming back as enabled in check connections function)
ERROR |2026-06-05 15:18:09.649|interrogatorServer[914]|file|(sendData.cpp)|line|(472)|Msg|(check connections return value: 1)
ERROR |2026-06-05 15:18:09.695|interrogatorServer[914]|file|(sendData.cpp)|line|(473)|Msg|(check connections wifi results: 1)
ERROR |2026-06-05 15:18:09.709|interrogatorServer[914]|file|(sendData.cpp)|line|(474)|Msg|(check connections cell results: 1)
INFO |2026-06-05 15:18:09.715|interrogatorServer[914]|file|(sendData.cpp)|line|(165)|Msg|(send getTimestamp for send readings )
INFO |2026-06-05 15:18:09.795|dataTransferService[1079]|file|(cmemsSrvIf.c)|line|(549)|Msg|(doGet:549: Attempting get of https://papi.g1.cardiomemshf.com/v6/timestamp to file
/tmp/getTimestamp.rsp)
INFO |2026-06-05 15:18:09.852|interrogatorServer[914]|file|(getTimestamp.cpp)|line|(53)|Msg|( [38;5;27;1mSending TimeStamp Msg to Backend [0m)
INFO |2026-06-05 15:18:12.786|dataTransferService[1079]|file|(inetServiceWorker.c)|line|(189)|Msg|(SUCCESS: Get TimeStamp)
ERROR |2026-06-05 15:18:12.800|interrogatorServer[914]|file|(cmsgwNetwork.c)|line|(591)|Msg|(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
INFO |2026-06-05 15:18:12.830|interrogatorServer[914]|file|(getTimestampHandler.cpp)|line|(65)|Msg|(Wifi enabled? TRUE and wlan0 ip address: 10.0.0.201)
ERROR |2026-06-05 15:18:13.131|interrogatorServer[914]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO |2026-06-05 15:18:13.132|interrogatorServer[914]|file|(getTimestampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and ppp0 ip address: (null))
```

P260300093

Friday, June 19, 2026 3:20 PM

```
ERROR |2026-06-03 19:27:38.538|CardioMEMS_GuiApp[896]|file|(cmsgwNetwork.c)|line|(1760)|Msg|(Process connmand is running)
INFO |2026-06-03 19:27:38.542|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(896)|Msg|(Turning on Wifi and cell)
INFO |2026-06-03 19:27:38.545|CardioMEMS_GuiApp[896]|file|(NetworkInfoController.cpp)|line|(228)|Msg|(Wifi is on ? 0)
INFO |2026-06-03 19:27:38.548|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(2266)|Msg|(Turning on cell)
INFO |2026-06-03 19:27:38.636|CardioMEMS_GuiApp[896]|file|(cmsgwNetwork.c)|line|(1471)|Msg|(Success to issue dbus method call SetProperty on
net.connman.Manager: /net/connman/technology/cellular)
INFO |2026-06-03 19:27:38.644|cellServer[535]|file|(cellMgrState.c)|line|(164)|Msg|(Cell powered up)
INFO |2026-06-03 19:27:38.645|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Start)
INFO |2026-06-03 19:27:38.666|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Check For Cell)
INFO |2026-06-03 19:27:38.669|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:27:41.697|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:27:43.456|watchdogServer[1068]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (15009) is less than 15 mins)
INFO |2026-06-03 19:27:44.741|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:27:47.769|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:27:50.802|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:27:53.835|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
WARNING |2026-06-03 19:27:54.087|dataTransferService[1059]|file|(cmsgwHttp.c)|line|(205)|Msg|(curl_easy_perform() failed: Timeout was reached)
WARNING |2026-06-03 19:27:54.088|dataTransferService[1059]|file|(inetServiceWorker.c)|line|(181)|Msg|(FAILURE: Get TimeStamp (-1))
ERROR |2026-06-03 19:27:54.090|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(591)|Msg|(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
ERROR |2026-06-03 19:27:54.091|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=wlan0)
INFO |2026-06-03 19:27:54.092|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(65)|Msg|(Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR |2026-06-03 19:27:54.260|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO |2026-06-03 19:27:54.261|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and ppp0 ip address: (null))

INFO |2026-06-03 19:27:59.902|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:28:02.962|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:28:05.983|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:28:06.030|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Switch Carriers)
INFO |2026-06-03 19:28:06.057|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Update Carrier List)
ALERT |2026-06-03 19:28:06.091|cellServer[535]|file|(cellMgrState.c)|line|(378)|Msg|(Default STATE_FAULT)
ALERT |2026-06-03 19:28:06.117|cellServer[535]|file|(cellMgrState.c)|line|(386)|Msg|(Query for new carrier list...)
INFO |2026-06-03 19:28:13.479|watchdogServer[1068]|file|(watchdogSupport.h)|line|(342)|Msg|(Data received on shared memory wd_reset)
INFO |2026-06-03 19:28:13.480|watchdogServer[1068]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (0) is less than 15 mins)
WARNING |2026-06-03 19:28:26.719|dataTransferService[1059]|file|(cmsgwHttp.c)|line|(205)|Msg|(curl_easy_perform() failed: Timeout was reached)
WARNING |2026-06-03 19:28:26.731|dataTransferService[1059]|file|(inetServiceWorker.c)|line|(181)|Msg|(FAILURE: Get TimeStamp (-1))
ERROR |2026-06-03 19:28:26.748|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(591)|Msg|(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
ERROR |2026-06-03 19:28:26.777|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=wlan0)
INFO |2026-06-03 19:28:26.799|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(65)|Msg|(Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR |2026-06-03 19:28:27.327|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO |2026-06-03 19:28:27.329|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and ppp0 ip address: (null))
INFO |2026-06-03 19:28:27.530|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(69)|Msg|(Kernel IP routing table)
INFO |2026-06-03 19:28:27.535|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(74)|Msg|(Received timeStamp message reply; -1)
ERROR |2026-06-03 19:28:27.546|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(101)|Msg|(getTimeStamp Failed)
INFO |2026-06-03 19:28:28.109|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(104)|Msg|(Cellular output ppp0: error fetching interface information:
Device not found)
INFO |2026-06-03 19:28:28.111|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(105)|Msg|(wlan0 output wlan0: flags=4098<BROADCAST,MULTICAST>
mtu 1500)

INFO |2026-06-03 19:28:58.526|watchdogServer[1068]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (45044) is less than 15 mins)
WARNING |2026-06-03 19:28:59.095|dataTransferService[1059]|file|(cmsgwHttp.c)|line|(205)|Msg|(curl_easy_perform() failed: Timeout was reached)
WARNING |2026-06-03 19:28:59.096|dataTransferService[1059]|file|(inetServiceWorker.c)|line|(181)|Msg|(FAILURE: Get TimeStamp (-1))
ERROR |2026-06-03 19:28:59.097|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(591)|Msg|(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
ERROR |2026-06-03 19:28:59.098|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=wlan0)
INFO |2026-06-03 19:28:59.099|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(65)|Msg|(Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR |2026-06-03 19:28:59.600|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO |2026-06-03 19:28:59.612|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and ppp0 ip address: (null))
INFO |2026-06-03 19:29:00.072|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(69)|Msg|(Kernel IP routing table)
INFO |2026-06-03 19:29:00.074|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(74)|Msg|(Received timeStamp message reply; -1)
ERROR |2026-06-03 19:29:00.080|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(101)|Msg|(getTimeStamp Failed)
INFO |2026-06-03 19:29:01.151|interrogatorServer[878]|file|(getTimeStampHandler.cpp)|line|(104)|Msg|(Cellular output ppp0: error fetching interface information:
Device not found)
```

INFO |2026-06-03 19:29:28.084|cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: T-Mobile status available)
 INFO |2026-06-03 19:29:28.086|cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: AT&T status available)
 INFO |2026-06-03 19:29:28.087|cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: Verizon status available)
 INFO |2026-06-03 19:29:28.088|cellServer[535]|file|(cmsgwCell.c)|line|(331)|Msg|(Carrier found: 313 100 status available)
 INFO |2026-06-03 19:29:28.089|cellServer[535]|file|(cellMgrState.c)|line|(394)|Msg|(Got carrier list!)
 INFO |2026-06-03 19:29:28.090|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Switch Carriers)
 INFO |2026-06-03 19:29:28.091|cellServer[535]|file|(cellMgrState.c)|line|(427)|Msg|(Attempt to connect to carrier /jacs_0/operator/310260)
 INFO |2026-06-03 19:29:28.555|watchdogServer[1068]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (75073) is less than 15 mins)
 ERROR |2026-06-03 19:29:28.950|cellServer[535]|file|(cmsgwCell.c)|line|(386)|Msg|(Failed to issue dbus method call Register on org.ofono.NetworkOperator path /jacs_0/operator/310260 returned -5: Modem is not registered to the network)
 INFO |2026-06-03 19:29:29.953|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Check Carrier Status)
 INFO |2026-06-03 19:29:29.957|cellServer[535]|file|(cellMgrState.c)|line|(452)|Msg|(Checking state for carrier: /jacs_0/operator/310260)
 INFO |2026-06-03 19:29:30.012|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Wait for Carrier Current)
 INFO |2026-06-03 19:29:30.138|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Carrier Current)
 INFO |2026-06-03 19:29:30.145|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Check For Cell)
 INFO |2026-06-03 19:29:30.190|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
 WARNING |2026-06-03 19:29:32.034|dataTransferService[1059]|file|(cmsgwHttp.c)|line|(205)|Msg|(curl_easy_perform() failed: Timeout was reached)

INFO |2026-06-03 19:29:32.527|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(2210)|Msg|(Power button enabled)

INFO |2026-06-03 19:29:32.531|CardioMEMS_GuiApp[896]|file|(AudioController.cpp)|line|(264)|Msg|(Qt Stop Audio play is called)
 INFO |2026-06-03 19:29:32.632|interrogatorServer[878]|file|(dataTransferError.cpp)|line|(98)|Msg|(Switching Machine to Daily)
 INFO |2026-06-03 19:29:32.965|CardioMEMS_GuiApp[896]|file|(AudioController.cpp)|line|(255)|Msg|(Qt Start Audio play is called)
 INFO |2026-06-03 19:29:33.285|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
 INFO |2026-06-03 19:29:33.299|cellServer[535]|file|(cellMgrState.c)|line|(311)|Msg|(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798086_context1 State: idle)
 INFO |2026-06-03 19:29:33.305|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Connect Cell Service)
 INFO |2026-06-03 19:29:33.309|cellServer[535]|file|(cellMgrState.c)|line|(347)|Msg|(Attempt 0 of 7 to connect to cell service: cellular_262017840798086_context1)
 WARNING |2026-06-03 19:29:33.325|cellServer[535]|file|(cmsgwNetwork.c)|line|(1428)|Msg|(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)
 INFO |2026-06-03 19:29:33.668|CardioMEMS_GuiApp[896]|file|(AudioController.cpp)|line|(264)|Msg|(Qt Stop Audio play is called)
 INFO |2026-06-03 19:29:33.765|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Wait for Connman Cell Online)
 INFO |2026-06-03 19:29:33.857|cellServer[535]|file|(cellMgrState.c)|line|(148)|Msg|(Network state=true)
 INFO |2026-06-03 19:29:33.859|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Online)
 ALERT |2026-06-03 19:29:33.868|cellServer[535]|file|(cellMgrState.c)|line|(490)|Msg|(NOW ONLINE!)

INFO |2026-06-03 19:29:43.567|watchdogServer[1068]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (0) is less than 15 mins)
 INFO |2026-06-03 19:29:43.781|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SettingConnectionPageElement.qml - 70 - Cellular Button Clicked : OFF)
 INFO |2026-06-03 19:29:44.154|cellServer[535]|file|(cellMgrState.c)|line|(148)|Msg|(Network state=false)
 INFO |2026-06-03 19:29:44.156|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Offline)
 INFO |2026-06-03 19:29:44.158|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Check For Cell)
 INFO |2026-06-03 19:29:44.182|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:29:44.351|cellServer[535]|file|(cellMgrState.c)|line|(311)|Msg|(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798086_context1 State: idle)
 INFO |2026-06-03 19:29:44.371|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Connect Cell Service)
 INFO |2026-06-03 19:29:44.395|cellServer[535]|file|(cellMgrState.c)|line|(347)|Msg|(Attempt 0 of 7 to connect to cell service: cellular_262017840798086_context1)
 WARNING |2026-06-03 19:29:44.424|cellServer[535]|file|(cmsgwNetwork.c)|line|(1428)|Msg|(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)
 INFO |2026-06-03 19:29:44.613|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(1471)|Msg|(Success to issue dbus method call SetProperty on net.connman.Manager: /net/connman/technology/cellular)
 WARNING |2026-06-03 19:29:44.890|cellServer[535]|file|(cmsgwNetwork.c)|line|(1428)|Msg|(Failed to issue dbus method call Connect on net.connman.Service: Operation aborted)
 INFO |2026-06-03 19:29:44.896|cellServer[535]|file|(cellMgrState.c)|line|(175)|Msg|(Cell powered down)
 INFO |2026-06-03 19:29:44.915|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Wait for Connman Cell Online)
 INFO |2026-06-03 19:29:44.930|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Cell not powered)
 INFO |2026-06-03 19:29:44.937|cellServer[535]|file|(cellMgrState.c)|line|(540)|Msg|(Cell not powered)
 INFO |2026-06-03 19:29:46.751|CardioMEMS_GuiApp[896]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SettingConnectionPageElement.qml - 70 - Cellular Button Clicked : ON)
 INFO |2026-06-03 19:29:46.938|interrogatorServer[878]|file|(cmsgwNetwork.c)|line|(1471)|Msg|(Success to issue dbus method call SetProperty on net.connman.Manager: /net/connman/technology/cellular)
 INFO |2026-06-03 19:29:46.951|cellServer[535]|file|(cellMgrState.c)|line|(164)|Msg|(Cell powered up)
 INFO |2026-06-03 19:29:46.954|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Start)
 INFO |2026-06-03 19:29:47.065|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Check For Cell)

INFO |2026-06-03 19:29:47.071|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
 INFO |2026-06-03 19:29:50.114|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO |2026-06-03 19:29:50.150|cellServer[535]|file|(cellMgrState.c)|line|(311)|Msg|(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798086_context1 State: idle)
 INFO |2026-06-03 19:29:50.156|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Connect Cell Service)
 INFO |2026-06-03 19:29:50.157|cellServer[535]|file|(cellMgrState.c)|line|(347)|Msg|(Attempt 0 of 7 to connect to cell service: cellular_262017840798086_context1)
 WARNING |2026-06-03 19:29:50.175|cellServer[535]|file|(cmwNetwork.c)|line|(1428)|Msg|(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)
 INFO |2026-06-03 19:29:50.857|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Wait for Connman Cell Online)
 INFO |2026-06-03 19:29:50.947|cellServer[535]|file|(cellMgrState.c)|line|(148)|Msg|(Network state=true)
 INFO |2026-06-03 19:29:50.949|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Online)
ALERT |2026-06-03 19:29:50.958|cellServer[535]|file|(cellMgrState.c)|line|(490)|Msg|(NOW ONLINE!)

ERROR |2026-06-03 19:29:53.523|interrogatorServer[878]|file|(sendData.cpp)|line|(469)|Msg|(The wifi or cell is coming back as enabled in check connections function)
 ERROR |2026-06-03 19:29:53.582|interrogatorServer[878]|file|(sendData.cpp)|line|(472)|Msg|(check connections return value: 1)
 ERROR |2026-06-03 19:29:53.650|interrogatorServer[878]|file|(sendData.cpp)|line|(473)|Msg|(check connections wifi results: 0)
 ERROR |2026-06-03 19:29:53.701|interrogatorServer[878]|file|(sendData.cpp)|line|(474)|Msg|(check connections cell results: 1)
 INFO |2026-06-03 19:29:53.777|interrogatorServer[878]|file|(sendData.cpp)|line|(165)|Msg|(send geTimestamp for send readings)
 INFO |2026-06-03 19:29:53.810|interrogatorServer[878]|file|(sendData.cpp)|line|(166)|Msg|(checked for connections and now will send gettimestamp)
 INFO |2026-06-03 19:29:53.891|interrogatorServer[878]|file|(getTimeStamp.cpp)|line|(45)|Msg|([38;5;27;1mGetTimeStampEv event: GetTimeStamp [0m)
 INFO |2026-06-03 19:29:53.926|dataTransferService[1059]|file|(cmemsSrvIf.c)|line|(682)|Msg|(configureHttpDefaults: configureHttpDefaults success.)
 INFO |2026-06-03 19:29:53.941|dataTransferService[1059]|file|(cmemsSrvIf.c)|line|(549)|Msg|(doGet:549: Attempting get of <https://papi.g1.cardiomemshf.com/v6/timestamp/to/file/tmp/getTimeStamp.rsp>)
 INFO |2026-06-03 19:29:53.960|interrogatorServer[878]|file|(getTimeStamp.cpp)|line|(53)|Msg|([38;5;27;1mSending TimeStamp Msg to Backend [0m)
 INFO |2026-06-03 19:29:58.574|watchdogServer[1068]|file|(watchdogSupport.h)|line|(342)|Msg|(Data received on shared memory wd_reset)
 INFO |2026-06-03 19:29:58.576|watchdogServer[1068]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (0) is less than 15 mins)
INFO |2026-06-03 19:29:59.251|dataTransferService[1059]|file|(inetServiceWorker.c)|line|(189)|Msg|(SUCCESS: Get TimeStamp)

P260300406

Monday, June 22, 2026 6:11 AM

```
INFO | 2026-06-04 15:26:19.513 | waveformServer[924] | file | (data-model.c) | line | (1101) | Msg | ( [38;5;202m Turning RF Tx OFF [0m)
INFO | 2026-06-04 15:26:19.552
| CardioMEMS_GuiApp[901] | file | (GuiController.cpp) | line | (1433) | Msg | (GUI: Show screen
'SendReading')
INFO | 2026-06-04 15:26:19.557
| CardioMEMS_GuiApp[901] | file | (AudioController.cpp) | line | (264) | Msg | (Qt Stop Audio play is called )
INFO | 2026-06-04 15:26:19.559
| CardioMEMS_GuiApp[901] | file | (AudioController.cpp) | line | (264) | Msg | (Qt Stop Audio play is called )
INFO | 2026-06-04 15:26:19.565
| CardioMEMS_GuiApp[901] | file | (AudioController.cpp) | line | (264) | Msg | (Qt Stop Audio play is called )
ERROR | 2026-06-04 15:26:19.705 | waveformServer[924] | file | (data-model.c) | line | (382) | Msg | (. is not a
file)
ERROR | 2026-06-04 15:26:19.793 | waveformServer[924] | file | (data-model.c) | line | (382) | Msg | (.. is not a
file)
INFO | 2026-06-04 15:26:19.831
| interrogatorServer[882] | file | (sendData.cpp) | line | (82) | Msg | ( [38;5;27;1mStartSendData event:
SendData [0m)
INFO | 2026-06-04 15:26:19.914 | interrogatorServer[882] | file | (sendData.cpp) | line | (138) | Msg | (Either
send data flow just started or response received for other commands.)
INFO | 2026-06-04 15:26:19.926
| CardioMEMS_GuiApp[901] | file | (AudioController.cpp) | line | (255) | Msg | (Qt Start Audio play is called )
INFO | 2026-06-04 15:26:20.088 | pressureServer[918] | file | (interrogator-
controller.c) | line | (443) | Msg | ([interrogatorControllerControlDevice]
INFO | 2026-06-04 15:26:20.099
| interrogatorServer[882] | file | (networkSettings.c) | line | (111) | Msg | (networkConfigWiFi file
/usr/cardiomems/var/config/uiConfig.json)
ERROR | 2026-06-04 15:26:20.168 | interrogatorServer[882] | file | (sendData.cpp) | line | (469) | Msg | (The
wifi or cell is coming back as enabled in check connections function)
ERROR | 2026-06-04 15:26:20.266 | interrogatorServer[882] | file | (sendData.cpp) | line | (472) | Msg | (check
connections return value: 1)
ERROR | 2026-06-04 15:26:20.310 | interrogatorServer[882] | file | (sendData.cpp) | line | (473) | Msg | (check
connections wifi results: 0)
ERROR | 2026-06-04 15:26:20.404
| interrogatorServer[882] | file | (sendData.cpp) | line | (474) | Msg | (check connections cell results: 1)
INFO | 2026-06-04 15:26:20.459 | interrogatorServer[882] | file | (sendData.cpp) | line | (165) | Msg | (send
getTimestamp for send readings )
INFO | 2026-06-04 15:26:20.574 | interrogatorServer[882] | file | (sendData.cpp) | line | (166) | Msg | (checked
for connections and now will send gettimestamp)
```

Fails to get cellular

```
INFO | 2026-06-04 15:26:35.373
| CardioMEMS_GuiApp[901] | file | (GuiController.cpp) | line | (2266) | Msg | (Turning on cell)
INFO | 2026-06-04 15:26:35.447
| CardioMEMS_GuiApp[901] | file | (cmswNetwork.c) | line | (1471) | Msg | (Success to issue dbus method call
```

SetProperty on net.connman.Manager: /net/connman/technology/cellular)
INFO | 2026-06-04 15:26:35.458 | cellServer[535] | file |(cellMgrState.c)| line |(164)| Msg |(Cell powered up)
INFO | 2026-06-04 15:26:35.459 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Start)
INFO | 2026-06-04 15:26:35.470 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Check For Cell)
INFO | 2026-06-04 15:26:35.471 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:26:38.483 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:26:41.280
| watchdogServer[1070] | file |(watchdogSupport.h)| line |(366)| Msg |(Inactivity time (15012) is less than 15 mins)
INFO | 2026-06-04 15:26:41.507 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:26:44.545 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:26:47.559 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:26:50.607 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
WARNING | 2026-06-04 15:26:50.836
| dataTransferService[1061] | file |(cmswHttp.c)| line |(205)| Msg |(curl_easy_perform() failed: Timeout was reached)
WARNING | 2026-06-04 15:26:50.838
| dataTransferService[1061] | file |(inetServiceWorker.c)| line |(181)| Msg |(FAILURE: Get TimeStamp (-1))
ERROR | 2026-06-04 15:26:50.840
| interrogatorServer[882] | file |(cmswNetwork.c)| line |(591)| Msg |(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
ERROR | 2026-06-04 15:26:50.842
| interrogatorServer[882] | file |(cmswNetwork.c)| line |(531)| Msg |(Error with ioctl call, iface=wlan0)
INFO | 2026-06-04 15:26:50.854
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(65)| Msg |(Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR | 2026-06-04 15:26:51.004
| interrogatorServer[882] | file |(cmswNetwork.c)| line |(531)| Msg |(Error with ioctl call, iface=ppp0)
INFO | 2026-06-04 15:26:51.025
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(66)| Msg |(Cell enabled? FALSE and ppp0 ip address: (null))

Fails again

INFO | 2026-06-04 15:26:53.622 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:26:56.288
| watchdogServer[1070] | file |(watchdogSupport.h)| line |(366)| Msg |(Inactivity time (30020) is less than 15 mins)
INFO | 2026-06-04 15:26:56.649 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:26:59.666 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:27:02.677 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell


```

service)
INFO | 2026-06-04 15:27:02.682 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated
to: Switch Carriers)
INFO | 2026-06-04 15:27:02.706 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated
to: Update Carrier List)
ALERT | 2026-06-04 15:27:02.736 | cellServer[535] | file |(cellMgrState.c)| line |(378)| Msg |(Default
STATE_FAULT)
ALERT | 2026-06-04 15:27:02.749 | cellServer[535] | file |(cellMgrState.c)| line |(386)| Msg |(Query for new
carrier list...)
INFO | 2026-06-04 15:27:11.296
| watchdogServer[1070] | file |(watchdogSupport.h)| line |(366)| Msg |(Inactivity time (45028) is less than
15 mins)
WARNING | 2026-06-04 15:27:22.882
| dataTransferService[1061] | file |(cmswHttp.c)| line |(205)| Msg |(curl_easy_perform() failed: Timeout
was reached)
WARNING | 2026-06-04 15:27:22.895
| dataTransferService[1061] | file |(inetServiceWorker.c)| line |(181)| Msg |(FAILURE: Get TimeStamp (-1))
ERROR | 2026-06-04 15:27:22.901
| interrogatorServer[882] | file |(cmswNetwork.c)| line |(591)| Msg |(cannot find ppp0 operstate file,
/sys/class/net/ppp0/operstate)
ERROR | 2026-06-04 15:27:22.903
| interrogatorServer[882] | file |(cmswNetwork.c)| line |(531)| Msg |(Error with ioctl call, iface=wlan0)
INFO | 2026-06-04 15:27:22.932
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(65)| Msg |(Wifi enabled? FALSE and
wlan0 ip address: (null))
ERROR | 2026-06-04 15:27:23.072
| interrogatorServer[882] | file |(cmswNetwork.c)| line |(531)| Msg |(Error with ioctl call, iface=ppp0)
INFO | 2026-06-04 15:27:23.085
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(66)| Msg |(Cell enabled? FALSE and
ppp0 ip address: (null))
INFO | 2026-06-04 15:27:23.393
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(69)| Msg |(Kernel IP routing table
INFO | 2026-06-04 15:27:23.405
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(74)| Msg |(Received timeStamp
message reply; -1)
ERROR | 2026-06-04 15:27:23.457
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(101)| Msg |(getTimeStamp Failed)
INFO | 2026-06-04 15:27:25.234
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(104)| Msg |(Cellular output ppp0: error
fetching interface information: Device not found

WARNING | 2026-06-04 15:27:57.456
| dataTransferService[1061] | file |(inetServiceWorker.c)| line |(181)| Msg |(FAILURE: Get TimeStamp (-1))
ERROR | 2026-06-04 15:27:57.459
| interrogatorServer[882] | file |(cmswNetwork.c)| line |(591)| Msg |(cannot find ppp0 operstate file,
/sys/class/net/ppp0/operstate)
ERROR | 2026-06-04 15:27:57.473
| interrogatorServer[882] | file |(cmswNetwork.c)| line |(531)| Msg |(Error with ioctl call, iface=wlan0)
INFO | 2026-06-04 15:27:57.520
| interrogatorServer[882] | file |(getTimeStampHandler.cpp)| line |(65)| Msg |(Wifi enabled? FALSE and
wlan0 ip address: (null))

```

ERROR | 2026-06-04 15:27:57.870
|interrogatorServer[882]|file|(cmswNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO | 2026-06-04 15:27:57.901
|interrogatorServer[882]|file|(getTimeStampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and
ppp0 ip address: (null))
INFO | 2026-06-04 15:27:58.079
|interrogatorServer[882]|file|(getTimeStampHandler.cpp)|line|(69)|Msg|(Kernel IP routing table
INFO | 2026-06-04 15:27:58.100
|interrogatorServer[882]|file|(getTimeStampHandler.cpp)|line|(74)|Msg|(Received timeStamp
message reply; -1)
ERROR | 2026-06-04 15:27:58.111
|interrogatorServer[882]|file|(getTimeStampHandler.cpp)|line|(101)|Msg|(getTimeStamp Failed)
INFO | 2026-06-04 15:27:58.650
|interrogatorServer[882]|file|(getTimeStampHandler.cpp)|line|(104)|Msg|(Cellular output ppp0: error
fetching interface information: Device not found

Found carriers but still failed

INFO | 2026-06-04 15:28:11.381
|watchdogServer[1070]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (105113) is less than
15 mins)
INFO | 2026-06-04 15:28:19.204|cellServer[535]|file|(cmswCell.c)|line|(331)|Msg|(Carrier found: T-
Mobile status available)
INFO | 2026-06-04 15:28:19.221|cellServer[535]|file|(cmswCell.c)|line|(331)|Msg|(Carrier found: AT&T
status available)
INFO | 2026-06-04 15:28:19.281|cellServer[535]|file|(cmswCell.c)|line|(331)|Msg|(Carrier found:
Verizon status available)
INFO | 2026-06-04 15:28:19.340|cellServer[535]|file|(cmswCell.c)|line|(331)|Msg|(Carrier found: 313
100 status available)
INFO | 2026-06-04 15:28:19.359|cellServer[535]|file|(cellMgrState.c)|line|(394)|Msg|(Got carrier list!)
INFO | 2026-06-04 15:28:19.379|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated
to: Switch Carriers)
INFO | 2026-06-04 15:28:19.414|cellServer[535]|file|(cellMgrState.c)|line|(427)|Msg|(Attempt to
connect to carrier /jacs_0/operator/310260)
INFO | 2026-06-04 15:28:21.571|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated
to: Check Carrier Status)
INFO | 2026-06-04 15:28:21.615|cellServer[535]|file|(cellMgrState.c)|line|(452)|Msg|(Checking state
for carrier: /jacs_0/operator/310260)
INFO | 2026-06-04 15:28:21.796|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated
to: Wait for Carrier Current)
INFO | 2026-06-04 15:28:21.858|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated
to: Carrier Current)
INFO | 2026-06-04 15:28:21.888|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated
to: Check For Cell)
INFO | 2026-06-04 15:28:21.954|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell
service)
INFO | 2026-06-04 15:28:25.167|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell
service)
INFO | 2026-06-04 15:28:26.401
|watchdogServer[1070]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (120121) is less than
15 mins)
INFO | 2026-06-04 15:28:28.196|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell
service)

WARNING | 2026-06-04 15:28:29.780
|dataTransferService[1061]|file|(cmswHttp.c)|line|(205)|Msg|(curl_easy_perform() failed: Timeout was reached)
WARNING | 2026-06-04 15:28:29.782
|dataTransferService[1061]|file|(inetServiceWorker.c)|line|(181)|Msg|(FAILURE: Get TimeStamp (-1))
ERROR | 2026-06-04 15:28:29.785
|interrogatorServer[882]|file|(cmswNetwork.c)|line|(591)|Msg|(cannot find ppp0 operstate file, /sys/class/net/ppp0/operstate)
ERROR | 2026-06-04 15:28:29.786
|interrogatorServer[882]|file|(cmswNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=wlan0)
INFO | 2026-06-04 15:28:29.799
|interrogatorServer[882]|file|(getTimeStampHandler.cpp)|line|(65)|Msg|(Wifi enabled? FALSE and wlan0 ip address: (null))
ERROR | 2026-06-04 15:28:30.319
|interrogatorServer[882]|file|(cmswNetwork.c)|line|(531)|Msg|(Error with ioctl call, iface=ppp0)
INFO | 2026-06-04 15:28:30.325
|interrogatorServer[882]|file|(getTimeStampHandler.cpp)|line|(66)|Msg|(Cell enabled? FALSE and ppp0 ip address: (null))

Finally connects

ERROR | 2026-06-04 15:28:30.545
|interrogatorServer[882]|file|(getTimeStampHandler.cpp)|line|(139)|Msg|(getTimeStamp failed and will not retry further)
INFO | 2026-06-04 15:28:30.556
|interrogatorServer[882]|file|(sendData.cpp)|line|(82)|Msg|([38;5;27;1mStartSendData event: SendData [0m)
INFO | 2026-06-04 15:28:30.561
|interrogatorServer[882]|file|(sendData.cpp)|line|(116)|Msg|(getTimeStampRetry exhausted, sending message to GUI to display connection error)
INFO | 2026-06-04 15:28:30.573
|CardioMEMS_GuiApp[901]|file|(GuiController.cpp)|line|(1458)|Msg|(GUI: Show dialog 'NoConnection')
INFO | 2026-06-04 15:28:30.583
|CardioMEMS_GuiApp[901]|file|(GuiController.cpp)|line|(2210)|Msg|(Power button enabled)
INFO | 2026-06-04 15:28:30.588
|CardioMEMS_GuiApp[901]|file|(AudioController.cpp)|line|(264)|Msg|(Qt Stop Audio play is called)
INFO | 2026-06-04 15:28:30.688
|interrogatorServer[882]|file|(dataTransferError.cpp)|line|(98)|Msg|(Switching Machine to Daily)
INFO | 2026-06-04 15:28:30.981
|CardioMEMS_GuiApp[901]|file|(AudioController.cpp)|line|(255)|Msg|(Qt Start Audio play is called)
INFO | 2026-06-04 15:28:31.206|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)
INFO | 2026-06-04 15:28:31.686
|CardioMEMS_GuiApp[901]|file|(AudioController.cpp)|line|(264)|Msg|(Qt Stop Audio play is called)
INFO | 2026-06-04 15:28:31.990
|CardioMEMS_GuiApp[901]|file|(AudioController.cpp)|line|(255)|Msg|(Qt Start Audio play is called)
INFO | 2026-06-04 15:28:33.089
|CardioMEMS_GuiApp[901]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/dialogs/DialogStackElement.qml - 380 - noConnectionDialog Connection Settings ButtonClicked)
INFO | 2026-06-04 15:28:33.095
|CardioMEMS_GuiApp[901]|file|(GuiController.cpp)|line|(1458)|Msg|(GUI: Show dialog 'None')
INFO | 2026-06-04 15:28:34.195

|CardioMEMS_GuiApp[901]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SettingConnectionPageElement.qml - 70 - Cellular Button Clicked : OFF)

INFO |2026-06-04 15:28:34.285|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)

INFO |2026-06-04 15:28:34.576|interrogatorServer[882]|file|(cmswNetwork.c)|line|(1471)|Msg|(Success to issue dbus method call SetProperty on net.connman.Manager: /net/connman/technology/cellular)

INFO |2026-06-04 15:28:34.654|cellServer[535]|file|(cellMgrState.c)|line|(175)|Msg|(Cell powered down)

INFO |2026-06-04 15:28:34.655|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Cell not powered)

INFO |2026-06-04 15:28:34.657|cellServer[535]|file|(cellMgrState.c)|line|(540)|Msg|(Cell not powered)

INFO |2026-06-04 15:28:37.017|CardioMEMS_GuiApp[901]|file|(GuiController.cpp)|line|(1416)|Msg|(qrc:/pages/SettingConnectionPageElement.qml - 70 - Cellular Button Clicked : ON)

INFO |2026-06-04 15:28:37.348|interrogatorServer[882]|file|(cmswNetwork.c)|line|(1471)|Msg|(Success to issue dbus method call SetProperty on net.connman.Manager: /net/connman/technology/cellular)

INFO |2026-06-04 15:28:37.408|cellServer[535]|file|(cellMgrState.c)|line|(164)|Msg|(Cell powered up)

INFO |2026-06-04 15:28:37.414|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Start)

INFO |2026-06-04 15:28:37.548|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Check For Cell)

INFO |2026-06-04 15:28:37.554|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)

INFO |2026-06-04 15:28:40.651|cellServer[535]|file|(cellMgrState.c)|line|(304)|Msg|(Query for cell service)

INFO |2026-06-04 15:28:40.661|cellServer[535]|file|(cellMgrState.c)|line|(311)|Msg|(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798032_context1 State: idle)

INFO |2026-06-04 15:28:40.685|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Connect Cell Service)

INFO |2026-06-04 15:28:40.690|cellServer[535]|file|(cellMgrState.c)|line|(347)|Msg|(Attempt 0 of 7 to connect to cell service: cellular_262017840798032_context1)

WARNING |2026-06-04 15:28:40.703|cellServer[535]|file|(cmswNetwork.c)|line|(1428)|Msg|(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)

INFO |2026-06-04 15:28:41.386|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Wait for Connman Cell Online)

INFO |2026-06-04 15:28:41.440|watchdogServer[1070]|file|(watchdogSupport.h)|line|(342)|Msg|(Data received on shared memory wd_reset)

INFO |2026-06-04 15:28:41.446|watchdogServer[1070]|file|(watchdogSupport.h)|line|(366)|Msg|(Inactivity time (0) is less than 15 mins)

INFO |2026-06-04 15:28:41.594|cellServer[535]|file|(cellMgrState.c)|line|(148)|Msg|(Network state=true)

INFO |2026-06-04 15:28:41.599|cellServer[535]|file|(cellMgrState.c)|line|(207)|Msg|(State updated to: Online)

ALERT |2026-06-04 15:28:41.623|cellServer[535]|file|(cellMgrState.c)|line|(490)|Msg|(NOW ONLINE!)

ERROR | 2026-06-04 15:28:43.075 | interrogatorServer[882] | file |(sendData.cpp)| line |(469)| Msg |(The wifi or cell is coming back as enabled in check connections function)
ERROR | 2026-06-04 15:28:43.106 | interrogatorServer[882] | file |(sendData.cpp)| line |(472)| Msg |(check connections return value: 1)
ERROR | 2026-06-04 15:28:43.133 | interrogatorServer[882] | file |(sendData.cpp)| line |(473)| Msg |(check connections wifi results: 0)
ERROR | 2026-06-04 15:28:43.167 | interrogatorServer[882] | file |(sendData.cpp)| line |(474)| Msg |(check connections cell results: 1)
INFO | 2026-06-04 15:28:43.185 | interrogatorServer[882] | file |(sendData.cpp)| line |(165)| Msg |(send geTimestamp for send readings)
INFO | 2026-06-04 15:28:43.199 | interrogatorServer[882] | file |(sendData.cpp)| line |(166)| Msg |(checked for connections and now will send gettimestamp)
INFO | 2026-06-04 15:28:43.242 | interrogatorServer[882] | file |(getTimeStamp.cpp)| line |(45)| Msg |([38;5;27;1mGetTimeStampEv event: GetTimeStamp [0m)
INFO | 2026-06-04 15:28:43.249 | dataTransferService[1061] | file |(cmemsSrvlf.c)| line |(682)| Msg |(configureHttpDefaults: configureHttpDefaults success.)
INFO | 2026-06-04 15:28:43.250 | dataTransferService[1061] | file |(cmemsSrvlf.c)| line |(549)| Msg |(doGet:549: Attempting get of <https://papi.g1.cardiomemshf.com/v6/timestamp to file /tmp/getTimeStamp.rsp>)
INFO | 2026-06-04 15:28:43.316 | interrogatorServer[882] | file |(getTimeStamp.cpp)| line |(53)| Msg |([38;5;27;1mSending TimeStamp Msg to Backend [0m)
INFO | 2026-06-04 15:28:50.492 | dataTransferService[1061] | file |(inetServiceWorker.c)| line |(189)| Msg |(SUCCESS: Get TimeStamp)

Device connects to cellular

WARNING | 2026-06-04 15:17:06.808 | cellServer[539] | file |(cmswCell.c)| line |(612)| Msg |(Failed sd_bus_message_enter_container: Unknown error -6)
INFO | 2026-06-04 15:17:06.821 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Cell not powered)
INFO | 2026-06-04 15:17:06.826 | cellServer[539] | file |(cellMgrState.c)| line |(540)| Msg |(Cell not powered)
WARNING | 2026-06-04 15:17:06.877 | cellServer[539] | file |(cmswNetwork.c)| line |(1354)| Msg |(Failed to issue dbus method call GetTechnologies on net.connman.Manager: The name net.connman was not provided by any .service files)
INFO | 2026-06-04 15:17:06.878 | cellServer[539] | file |(cellMgrState.c)| line |(550)| Msg |(Cell initialized to not powered)

INFO | 2026-06-04 15:17:06.879 | cellServer[539] | file |(cellMgrState.c)| line |(540)| Msg |(Cell not powered)

INFO | 2026-06-04 15:17:08.395 | cellServer[539] | file |(cellMgrState.c)| line |(164)| Msg |(Cell powered up)

INFO | 2026-06-04 15:17:08.396 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Start)

INFO | 2026-06-04 15:17:08.428 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Check For Cell)

INFO | 2026-06-04 15:17:08.429 | cellServer[539] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)

INFO | 2026-06-04 15:17:11.515 | cellServer[539] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)

INFO | 2026-06-04 15:17:11.526 | cellServer[539] | file |(cellMgrState.c)| line |(311)| Msg |(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798032_context1 State: idle)

INFO | 2026-06-04 15:17:11.604 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Connect Cell Service)

INFO | 2026-06-04 15:17:11.607 | cellServer[539] | file |(cellMgrState.c)| line |(347)| Msg |(Attempt 0 of 7 to connect to cell service: cellular_262017840798032_context1)

WARNING | 2026-06-04 15:17:11.623 | cellServer[539] | file |(cmswNetwork.c)| line |(1428)| Msg |(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)

INFO | 2026-06-04 15:17:12.379 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Wait for Connman Cell Online)

INFO | 2026-06-04 15:17:12.439 | cellServer[539] | file |(cellMgrState.c)| line |(148)| Msg |(Network state=true)

INFO | 2026-06-04 15:17:12.440 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Online)

ALERT | 2026-06-04 15:17:12.442 | cellServer[539] | file |(cellMgrState.c)| line |(490)| Msg |(NOW ONLINE!)

Disconnects then reconnects

INFO | 2026-06-04 15:17:39.262 | cellServer[539] | file |(cellMgrState.c)| line |(148)| Msg |(Network state=false)

INFO | 2026-06-04 15:17:39.263 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Offline)

INFO | 2026-06-04 15:17:39.264 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Check For Cell)

INFO | 2026-06-04 15:17:39.265 | cellServer[539] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)

INFO | 2026-06-04 15:17:39.851
| interrogatorServer[850] | file |(tempmonitor.cpp)| line |(107)| Msg |(DOES NOT EXIST: /usr/cardiomems/var/tmp/armBarometer)

INFO | 2026-06-04 15:17:40.852
| interrogatorServer[850] | file |(tempmonitor.cpp)| line |(107)| Msg |(DOES NOT EXIST: /usr/cardiomems/var/tmp/armBarometer)

INFO | 2026-06-04 15:17:41.854
| interrogatorServer[850] | file |(tempmonitor.cpp)| line |(107)| Msg |(DOES NOT EXIST: /usr/cardiomems/var/tmp/armBarometer)

INFO | 2026-06-04 15:17:42.526 | cellServer[539] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)

INFO | 2026-06-04 15:17:42.532 | cellServer[539] | file |(cellMgrState.c)| line |(311)| Msg |(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798032_context1)

State: idle)
INFO | 2026-06-04 15:17:42.534 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Connect Cell Service)
INFO | 2026-06-04 15:17:42.535 | cellServer[539] | file |(cellMgrState.c)| line |(347)| Msg |(Attempt 0 of 7 to connect to cell service: cellular_262017840798032_context1)
WARNING | 2026-06-04 15:17:42.541 | cellServer[539] | file |(cmswNetwork.c)| line |(1428)| Msg |(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)
INFO | 2026-06-04 15:17:42.789 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Wait for Connman Cell Online)
INFO | 2026-06-04 15:17:42.828 | cellServer[539] | file |(cellMgrState.c)| line |(148)| Msg |(Network state=true)
INFO | 2026-06-04 15:17:42.829 | cellServer[539] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Online)
ALERT | 2026-06-04 15:17:42.831 | cellServer[539] | file |(cellMgrState.c)| line |(490)| Msg |(NOW ONLINE!)

WARNING | 2026-06-04 15:19:46.526 | cellServer[535] | file |(cmswCell.c)| line |(612)| Msg |(Failed sd_bus_message_enter_container: Unknown error -6)
INFO | 2026-06-04 15:19:46.550 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Cell not powered)
INFO | 2026-06-04 15:19:46.551 | cellServer[535] | file |(cellMgrState.c)| line |(540)| Msg |(Cell not powered)
WARNING | 2026-06-04 15:19:46.559 | cellServer[535] | file |(cmswNetwork.c)| line |(1354)| Msg |(Failed to issue dbus method call GetTechnologies on net.connman.Manager: The name net.connman was not provided by any .service files)
INFO | 2026-06-04 15:19:46.560 | cellServer[535] | file |(cellMgrState.c)| line |(550)| Msg |(Cell initialized to not powered)
INFO | 2026-06-04 15:19:46.597 | cellServer[535] | file |(cellMgrState.c)| line |(540)| Msg |(Cell not powered)
INFO | 2026-06-04 15:19:48.141 | cellServer[535] | file |(cellMgrState.c)| line |(164)| Msg |(Cell powered up)
INFO | 2026-06-04 15:19:48.142 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Start)
INFO | 2026-06-04 15:19:48.181 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Check For Cell)
INFO | 2026-06-04 15:19:48.181 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:19:51.207 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:19:51.213 | cellServer[535] | file |(cellMgrState.c)| line |(311)| Msg |(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798032_context1 State: idle)
INFO | 2026-06-04 15:19:51.216 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Connect Cell Service)
INFO | 2026-06-04 15:19:51.264 | cellServer[535] | file |(cellMgrState.c)| line |(347)| Msg |(Attempt 0 of 7 to connect to cell service: cellular_262017840798032_context1)
WARNING | 2026-06-04 15:19:51.275 | cellServer[535] | file |(cmswNetwork.c)| line |(1428)| Msg |(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)
INFO | 2026-06-04 15:19:51.861 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Wait for Connman Cell Online)
INFO | 2026-06-04 15:19:51.883 | cellServer[535] | file |(cellMgrState.c)| line |(148)| Msg |(Network

state=true)
INFO | 2026-06-04 15:19:51.884 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Online)
ALERT | 2026-06-04 15:19:51.885 | cellServer[535] | file |(cellMgrState.c)| line |(490)| Msg |(NOW ONLINE!)

Disabled cell service, but trying to connect anywat

INFO | 2026-06-04 15:22:34.821
| CardioMEMS_GuiApp[901] | file |(GuiController.cpp)| line |(1416)| Msg |(qrc:/pages/StartPageElement.qml - 95: Start Button clicked)
INFO | 2026-06-04 15:22:34.832
| CardioMEMS_GuiApp[901] | file |(AudioController.cpp)| line |(264)| Msg |(Qt Stop Audio play is called)
INFO | 2026-06-04 15:22:34.841
| CardioMEMS_GuiApp[901] | file |(GuiController.cpp)| line |(1433)| Msg |(GUI: Show screen 'SignalAquisition')
INFO | 2026-06-04 15:22:34.848
| CardioMEMS_GuiApp[901] | file |(AudioController.cpp)| line |(264)| Msg |(Qt Stop Audio play is called)
INFO | 2026-06-04 15:22:34.851
| CardioMEMS_GuiApp[901] | file |(GuiController.cpp)| line |(820)| Msg |(Turning off Wifi)
INFO | 2026-06-04 15:22:34.852
| CardioMEMS_GuiApp[901] | file |(NetworkInfoController.cpp)| line |(228)| Msg |(Wifi is on ? 0)
INFO | 2026-06-04 15:22:34.861
| CardioMEMS_GuiApp[901] | file |(GuiController.cpp)| line |(2050)| Msg |(Turning off cell and Wi-Fi since reading is in progress)
ERROR | 2026-06-04 15:22:35.028
| CardioMEMS_GuiApp[901] | file |(cmswNetwork.c)| line |(1760)| Msg |(Process connmand is running
INFO | 2026-06-04 15:22:35.031
| CardioMEMS_GuiApp[901] | file |(NetworkInfoController.cpp)| line |(228)| Msg |(Wifi is on ? 0)
INFO | 2026-06-04 15:22:35.033
| CardioMEMS_GuiApp[901] | file |(GuiController.cpp)| line |(2242)| Msg |(Initialize cell network module)
**INFO | 2026-06-04 15:22:35.036
| CardioMEMS_GuiApp[901] | file |(GuiController.cpp)| line |(2246)| Msg |(Disabling cell network)**
INFO | 2026-06-04 15:22:35.159 | cellServer[535] | file |(cellMgrState.c)| line |(148)| Msg |(Network state=false)
INFO | 2026-06-04 15:22:35.160 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Offline)
INFO | 2026-06-04 15:22:35.161 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Check For Cell)
INFO | 2026-06-04 15:22:35.162 | cellServer[535] | file |(cellMgrState.c)| line |(304)| Msg |(Query for cell service)
INFO | 2026-06-04 15:22:35.197
| CardioMEMS_GuiApp[901] | file |(AudioController.cpp)| line |(255)| Msg |(Qt Start Audio play is called)
INFO | 2026-06-04 15:22:35.327 | cellServer[535] | file |(cellMgrState.c)| line |(311)| Msg |(SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798032_context1 State: idle)
INFO | 2026-06-04 15:22:35.328 | cellServer[535] | file |(cellMgrState.c)| line |(207)| Msg |(State updated to: Connect Cell Service)
INFO | 2026-06-04 15:22:35.329 | cellServer[535] | file |(cellMgrState.c)| line |(347)| Msg |(Attempt 0 of 7 to connect to cell service: cellular_262017840798032_context1)
WARNING | 2026-06-04 15:22:35.403 | cellServer[535] | file |(cmswNetwork.c)| line |(1428)| Msg |(Failed to issue dbus method call Disconnect on net.connman.Service: Not connected)

P260300455

Monday, June 22, 2026 7:15 AM

Comments

Monday, June 22, 2026 9:35 AM

| |
|---|
| carrier switch happened and it took
2026-05-08 17:41:14.135 = Cell ON
2026-05-08 17:42:49.767 = Got Cell IP Address |
| 2026-06-06 11:15:04.226 = Cell ON
2026-06-06 11:15:13.332 = Got Cell IP Address |
| 2026-06-05 08:58:21.857 = Cell ON
2026-06-05 08:58:31.047 = Got Cell IP Address |
| 2026-06-05 15:03:34.371 = Cell OFF and Wi-Fi OFF
2026-06-05 15:03:37.461 = Rx OFF
2026-06-05 15:03:53.503 = Turn ON Wi-Fi and Cell after taking a reading
2026-06-05 15:03:53.507 = Cell ON
2026-06-05 15:05:47.327 = NoConnection
2026-06-05 15:06:55.904 = Cell ONLINE, but no IP address |
| 2026-06-04 15:25:44.565 =Cell powered down
2026-06-04 15:26:19.513 = RF OFF - Reading completed
2026-06-04 15:26:20.168 = Wi-Fi / Cell turn on command sent
2026-06-04 15:26:20.678 = 1st Time stamp sent
2026-06-04 15:27:02.682 = Carrier Switch
2026-06-04 15:28:50.556 = Got ppp0 IP address
Sent reading using Cell = Took lot of time because of Carieer switch |

| |
|------------|
| P260400055 |
| P260400015 |
| P260400054 |
| P260300386 |
| P260300406 |
| P260300455 |
| P260300455 |
| P260300093 |

| | |
|------------|---|
| P260400055 | INFO 2026-06-08 14:56:00.826 cellServer[535] file (cellMgrState.c) line (394) Msg (Got carrier list!)
ERROR 2026-06-08 14:56:02.227 cellServer[535] file (cellMgrState.c) line (433) Msg (Carrier switch failed!) = Unable to carrier switch
ERROR 2026-06-08 14:56:03.396 cellServer[535] file (cellMgrState.c) line (398) Msg (Unable to get carrier list!) = Unable to get Carrier li st
WARNING 2026-06-08 14:57:38.847 dataTransferService[1062] file (inetServiceWorker.c) line (181) Msg (FAILURE: Get TimeStamp (-1)) =
TimeStamp failure on cellular
INFO 2026-06-08 14:59:21.490
 CardioMEMS_GuiApp[896] file (GuiController.cpp) line (1416) Msg (qrc:/pages/SetupWifiPageElement.qml - 85 - Wifi ON/OFF Button
Clicked : ON) = Turn Wi-Fi on
INFO 2026-06-08 15:00:02.793 cellServer[535] file (cellMgrState.c) line (394) Msg (Got carrier list!) = Cellular trying to connect in
background
ERROR 2026-06-08 15:00:04.092 cellServer[535] file (cellMgrState.c) line (433) Msg (Carrier switch failed!)
INFO 2026-06-08 15:00:53.829 interrogatorServer[878] file (getTimeStampHandler.cpp) line (65) Msg (Wifi enabled? TRUE and wlan0 ip
address: 172.20.10.4)
INFO 2026-06-08 15:00:54.174 interrogatorServer[878] file (getTimeStampHandler.cpp) line (66) Msg (Cell enabled? FALSE and ppp0 ip
address: (null)) = No IP address from cellular
INFO 2026-06-08 15:00:54.483 interrogatorServer[878] file (getTimeStampHandler.cpp) line (78) Msg (Get TimeStamp Successful) = Got
TimeStamp from Wi-Fi |
| P260400015 | INFO 2026-06-06 15:27:16.641 CardioMEMS_GuiApp[904] file (GuiController.cpp) line (2266) Msg (Turning on cell)
WARNING 2026-06-06 15:27:32.436 dataTransferService[1065] file (inetServiceWorker.c) line (181) Msg (FAILURE: Get TimeStamp (-1))
INFO 2026-06-06 15:27:32.764 interrogatorServer[900] file (getTimeStampHandler.cpp) line (66) Msg (Cell enabled? FALSE and ppp0 ip
address: (null))
INFO 2026-06-06 15:27:33.616 interrogatorServer[900] file (getTimeStampHandler.cpp) line (104) Msg (Cellular output ppp0: error
fetching interface information: Device not found
INFO 2026-06-06 15:29:10.740 CardioMEMS_GuiApp[904] file (GuiController.cpp) line (1458) Msg (GUI: Show dialog 'NoConnection')
INFO 2026-06-06 15:29:28.041
 CardioMEMS_GuiApp[904] file (GuiController.cpp) line (1416) Msg (qrc:/pages/SetupWifiPageElement.qml - 85 - Wifi ON/OFF Button
Clicked : ON)
INFO 2026-06-06 15:30:15.823 cellServer[539] file (cellMgrState.c) line (311) Msg (SUCCESS: Cell Carrier: Service Name: cellular_
262017840798815_context1 State: idle)
ERROR 2026-06-06 15:30:15.960 cellServer[539] file (cellMgrState.c) line (433) Msg (Carrier switch failed!)
INFO 2026-06-06 15:32:01.553 interrogatorServer[900] file (getTimeStampHandler.cpp) line (65) Msg (Wifi enabled? TRUE and wlan0 ip
address: 192.168.1.212)
ERROR 2026-06-06 15:32:01.885 interrogatorServer[900] file (cmswNetwork.c) line (531) Msg (Error with ioctl call, iface=ppp0)
INFO 2026-06-06 15:32:01.907 interrogatorServer[900] file (getTimeStampHandler.cpp) line (66) Msg (Cell enabled? FALSE and ppp0 ip
address: (null))
INFO 2026-06-06 15:32:01.507 dataTransferService[1065] file (inetServiceWorker.c) line (189) Msg (SUCCESS: Get TimeStamp) |

| | |
|------------|--|
| P260400054 | INFO 2026-06-05 19:01:32.336 cellServer[537] file (cellMgrState.c) line (394) Msg (Got carrier list!)
INFO 2026-06-05 19:01:35.534 cellServer[537] file (cellMgrState.c) line (311) Msg (SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840795179_context1 State: idle)
INFO 2026-06-05 19:01:35.546 cellServer[537] file (cellMgrState.c) line (207) Msg (State updated to: Connect Cell Service)
ALERT 2026-06-05 19:01:36.923 cellServer[537] file (cellMgrState.c) line (490) Msg (NOW ONLINE!)
WARNING 2026-06-05 19:01:43.055 dataTransferService[1061] file (inetServiceWorker.c) line (181) Msg (FAILURE: Get TimeStamp (-1))
INFO 2026-06-05 19:01:43.131 interrogatorServer[900] file (getTimeStampHandler.cpp) line (66) Msg (Cell enabled? TRUE and ppp0 ip address: 10.155.39.26)
INFO 2026-06-05 19:01:49.040 CardioMEMS_GuiApp[904] file (GuiController.cpp) line (1416) Msg (qrc:/pages/SettingConnectionPageElement.qml - 70 - Cellular Button Clicked : OFF)
INFO 2026-06-05 19:01:50.853 cellServer[537] file (cellMgrState.c) line (175) Msg (Cell powered down)
INFO 2026-06-05 19:01:52.238 CardioMEMS_GuiApp[904] file (GuiController.cpp) line (1416) Msg (qrc:/pages/SettingConnectionPageElement.qml - 70 - Cellular Button Clicked : ON)
INFO 2026-06-05 19:01:52.572 cellServer[537] file (cellMgrState.c) line (164) Msg (Cell powered up)
INFO 2026-06-05 19:01:55.732 cellServer[537] file (cellMgrState.c) line (311) Msg (SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840795179_context1 State: idle)
INFO 2026-06-05 19:01:56.633 cellServer[537] file (cellMgrState.c) line (207) Msg (State updated to: Online)
INFO 2026-06-05 19:02:04.364 dataTransferService[1061] file (inetServiceWorker.c) line (189) Msg (SUCCESS: Get TimeStamp) |
| P260300386 | INFO 2026-06-05 15:04:08.925 interrogatorServer[914] file (getTimeStampHandler.cpp) line (66) Msg (Cell enabled? FALSE and ppp0 ip address: (null))
INFO 2026-06-05 15:04:09.948 interrogatorServer[914] file (getTimeStampHandler.cpp) line (104) Msg (Cellular output ppp0: error fetching interface information: Device not found)
INFO 2026-06-05 15:06:15.878 CardioMEMS_GuiApp[918] file (GuiController.cpp) line (1416) Msg (qrc:/pages/SetupWifiPageElement.qml - 85 - Wifi ON/OFF Button Clicked : ON)
INFO 2026-06-05 15:06:53.608 cellServer[539] file (cellMgrState.c) line (311) Msg (SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840797883_context1 State: idle)
INFO 2026-06-05 15:06:54.303 cellServer[539] file (cellMgrState.c) line (347) Msg (Attempt 1 of 7 to connect to cell service: cellular_262017840797883_context1)
WARNING 2026-06-05 15:06:55.874 cellServer[539] file (cmswNetwork.c) line (1428) Msg (Failed to issue dbus method call Connect on net.connman.Service: Input/output error)
ALERT 2026-06-05 15:06:55.904 cellServer[539] file (cellMgrState.c) line (490) Msg (NOW ONLINE!)
INFO 2026-06-05 15:06:55.920 cellServer[539] file (cellMgrState.c) line (207) Msg (State updated to: Offline)
ERROR 2026-06-05 15:18:09.709 interrogatorServer[914] file (sendData.cpp) line (474) Msg (check connections cell results: 1)
INFO 2026-06-05 15:18:12.786 dataTransferService[1079] file (inetServiceWorker.c) line (189) Msg (SUCCESS: Get TimeStamp)

Remarks: Device connects to cellular, be immediately disconnects, within less than a second. |
| P260300406 | ERROR 2026-06-04 15:26:20.404 interrogatorServer[882] file (sendData.cpp) line (474) Msg (check connections cell results: 1)
INFO 2026-06-04 15:26:35.458 cellServer[535] file (cellMgrState.c) line (164) Msg (Cell powered up)
WARNING 2026-06-04 15:26:50.838 dataTransferService[1061] file (inetServiceWorker.c) line (181) Msg (FAILURE: Get TimeStamp (-1))
INFO 2026-06-04 15:26:51.025 interrogatorServer[882] file (getTimeStampHandler.cpp) line (66) Msg (Cell enabled? FALSE and ppp0 ip address: (null))
INFO 2026-06-04 15:27:02.682 cellServer[535] file (cellMgrState.c) line (207) Msg (State updated to: Switch Carriers)
INFO 2026-06-04 15:27:02.706 cellServer[535] file (cellMgrState.c) line (207) Msg (State updated to: Update Carrier List)
ALERT 2026-06-04 15:27:02.736 cellServer[535] file (cellMgrState.c) line (378) Msg (Default STATE_FAULT)
WARNING 2026-06-04 15:27:22.895 dataTransferService[1061] file (inetServiceWorker.c) line (181) Msg (FAILURE: Get TimeStamp (-1))
INFO 2026-06-04 15:27:23.085 interrogatorServer[882] file (getTimeStampHandler.cpp) line (66) Msg (Cell enabled? FALSE and ppp0 ip address: (null))
WARNING 2026-06-04 15:27:57.456 dataTransferService[1061] file (inetServiceWorker.c) line (181) Msg (FAILURE: Get TimeStamp (-1))
INFO 2026-06-04 15:28:19.359 cellServer[535] file (cellMgrState.c) line (394) Msg (Got carrier list!)
WARNING 2026-06-04 15:28:29.780 dataTransferService[1061] file (cmswHttp.c) line (205) Msg (curl_easy_perform() failed: Timeout was reached)
WARNING 2026-06-04 15:28:29.782 dataTransferService[1061] file (inetServiceWorker.c) line (181) Msg (FAILURE: Get TimeStamp (-1))
INFO 2026-06-04 15:28:40.661 cellServer[535] file (cellMgrState.c) line (311) Msg (SUCCESS: Cell Carrier: Telekom.de (T-Mobile Telekom.de) Service Name: cellular_262017840798032_context1 State: idle)
ALERT 2026-06-04 15:28:41.623 cellServer[535] file (cellMgrState.c) line (490) Msg (NOW ONLINE!) |
| P260300455 | |
| P260300093 | |

P260300443

cellular_262017840796356_context1

P260300435

Service Name: cellular_262017840796353_context1

P260300386

Service Name: cellular_262017840797883_context1

HERO Test Case

Thursday, May 14, 2026 9:06 AM

Connect Device

- Waveform Cord and pins
 - o Put 2 pin thing on end of 6 pin array
 - o GND goes on black pins, orange mid, yellow
 - o Black pins on J36 blank
- 12 V DC Supply
- Ethernet and USB Cable
- Modem
- Put waveform generator cable after hero is boot up

Setup

- Ttermpro
 - o Serial (usually COM3)
 - o Setup > Serial > Speed = 115200 > New Setting
- Waveform generator
 - o Need to install national instruments labview and runtime engine
 - Download runtime engine first then reboot then the shorter name one
 - o Extract pressure waveform generator and it should run
 - o Install driver for waveform cable (COM waveform titled folder)
 - The waveform cable should be found in the PORTS section, not USB Section
 - In device manager, rename COM port for waveform cable as COM 15 or above

Waveform generator replaces patient readings

To optimize device, using low power processor

- Need to turn off cellular and wifi while taking reading to optimize device

Not visible at end of execution, but wifi bars were

No cellular bars

Wifi

Same test case was executed for both dt and tmobile sims

Did not find any issues

Confusion is

No cellular

If it comes

If it is visible after 20 sec or if it is visible

Cellular Experiment 6/18

Thursday, June 18, 2026 10:55 AM

MSSI

- Software for testing purposes
- Has all the individual parts working like real board (cellular, wifi, etc)
- We own the linux and app software
 - o Application software has 3 parts
 - Data acquired from sensor and does an algorithm
 - If score is less than 40, reading is approved, otherwise suspect
 - ◆ Doctor can view waveform and ignore reading or ignore suspect classification

Experiment Context

- (Expected) In hospital, need to do setup of the device
 - o In merlin, once sensor is implanted, take the HERO and sensor
 - Once powered up for first time, its in setup mode (pairing mode)
 - At time of setup, device uses cellular
- (Actual) In some hospitals, when tak

When `1st reading is taking

Wifi and cellular bith off

Cellular is not cming up

Only happens 1st taime taking readings

Sees no connection (no cellular)

Turn on turn off celular

(sometimes work)

Or try wifi

lonly for 1st readings, not for 2nd ones

Common and ophone (lincx libraies)

Put some debug code on it

3 pARTS In debug code

And check logs

Should be persisitent after multiple
readigns

Change modem to new one

In hospital modems are currently DT

Device Information

| | |
|------------|--------------|
| SN | P25050017 |
| Env | V3 |
| IP Address | 172.20.10.14 |

Linux Commands

Login at root

Vi /usr/sbin/conman-systemd-wrapper

Login at root

Vi /usr/sbin/conman-systemd-wrapper

Use -d for debug in logs

Ofono

Vi /usr/sbin/ofonod

Vi Lib/systemd/system/ofono.service

Another file for persistent logs

/etc/systemd/journald.conf

Make storage-persistent

Enable debug conman logs:

Vi /usr/sbin/conman-systemd-wrapper

Change the lin

/usr/bin/conmand -n -d \$EXTRA_PARAM

Vi Lib/systemd/system/ofono.service

Exec-Start=/usr/

mkdir /usr/cardiomems/var/data_to_sned/cellogs

Logically tell 2 directories are equal (logs keep getting delated each time, make persisten)

Ln -s var/log/journal/ /usr/cardiomems/var/data_to_sned/cellogs

Changes service file to do ts

systemctl daemon-reload

systemctl restart conman

Make ofono logs consistent

Systemctl edit ofono

Ls /usr/cardiomems/var/data_to_sned/cellogs/journal/...

Started creating logs

Vire epersitne logs

Jounralctl -u ofono -b -l (show past logs)

827 0a6 (never delete folder)

Experiment

1. System resert
2. Took (old modem) reading
3. Before and after it connect
4. After cpmnecting, saw 1 bar

New modem

Old mdoem we did setup

New modem

Tried to take rading

2. New modem, tried setup
 - a. No connection
3. Setup with old modem
showed mwelcome screen

Tried to take reading with new modem

In btwn start and actual readin (4 min wait)

Were able to reproduce

Connected briefly for a second, lost it

No connection

Saved reading

Taking 2nd reading (delaying, waited 4 min) before rading (new mdoem)

- Connecting for lot longer
- Was able to connect eventually

eMMC file system check (starts) boot up one

After

Scoring: PASSED

RF Tx off when reading is taken

18 devices has 2 iss

1. Cell did not connect at all (4 in total)
2. Took time to connect (12 total)
 - a. Stopped with no connection screen

Sometimes, couldn't connect, (no connection)

Goes inside,

Have option to connect to wifi or toggle cellular

When toggle, connect cellular (if not, connect to wifi)

After completing read

Connectivity check results 1 (connected to cellular)

Sent timestamp to back end

Turn on wifi and cell

On device

Wifi and cell would be connect based on previous settings

ONLINE means connected

Suddenly

Network state = true to network state = false

Started carrier switch

Cannot connect backend, no network

If gets IP address, it will send

Table blank = [no ip address

timestamp = -1 (didn't work)

Send later = saved reading,

Want timestamp of tx off

For 2nd reading expecting to see timestamp

Ppp0 = only cellular

Windchill

Friday, May 15, 2026 3:55 PM

Change content

Actions > replace Content

Upload redline first

Replace content

Then clean one

Process > author

History

Right click, mark as redline

Change order > Change activity > personal

Change activity = what was changed

Change Order doesn't see changes

Created redline locally, make clear version locally

Api sync quart

Old api was rudimentray dump

One fuz was creating specific apis that would work

No confirmation yet that app is using new API sync

If we officially test bugs

Can track in sperate bug (sync testing)

Not any official testing yet