

CardioGuide and Merlin Basics

Monday, August 17, 2026 10:59 PM


Installing the App

On Android, for formal testing, you should install app via google Play

- Paste google play link in QR code generator and scan it
 - o Link:
- Download it from Google play

On Android, if google Play store link is not updated, you can use the APK file

- [App Binaries To QA](#)
- Download the apk to local computer
- In your command prompt, cd into the folder with the apk
- Run **adb install <apk file name>** to download the apk file onto the phone

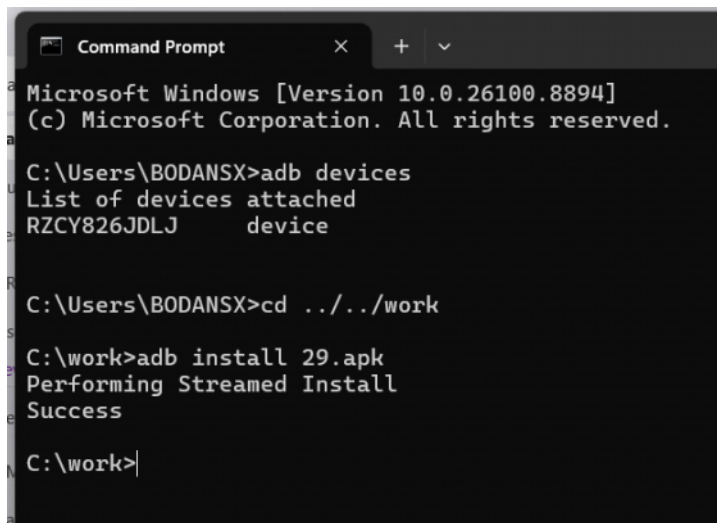

```
Microsoft Windows [Version 10.0.26100.8894]
(c) Microsoft Corporation. All rights reserved.

C:\Users\BODANSX>adb devices
List of devices attached
RZCY826JDLJ    device

C:\Users\BODANSX>cd ../../work

C:\work>adb install 29.apk
Performing Streamed Install
Success

C:\work>
```

On iOS, should install app using Test Flight

- Use Test Flight to download ipa files to IOS phones
- No need for Macbook

Patient Registration

- First enroll a patient in the merlin portal
 - o Q3 Merlin Portal: q3-hfcloud.products.abbott/apps/clinician/notifications
 - Clinic Username: q3_svwq3
 - Clinic Password: Test!1234
 - o T1 Merlin Portal: t1-hfcloud.products.abbott/apps/clinician/notifications
 - US Clinic Username: ust1swvt1
 - Korean Clinic Username: krt1swvt1
 - Clinic Password: Test#1234
 - o Press the Enroll patient button
 - Enter Patient first and last name
 - DOB 01-01-2020
 - Once on the Patient profile settings,

- Select Treating Physician (You are logging into Merlin using premade clinician)
 - ◆ Lee, Danny as physician
- Fill out sensor Serial #
 - ◆ Choose unclaimed in backend (SWVT lot)
 - ◇ Q3 Sensors: [Sign In | CMHF Customer Care](#)
 - ◇ T1 Sensors: [Sign In | CMHF Customer Care](#)
 - ◆ Enter copied sensor in Merlin, pause for website to validate # with other patients
 - ◇ The website checks with other patients to check if other patients have same serial #, if not, box should be blue
- Patient is enrolled, but not active, so it will still be scheduled patients tab
 - You can register a patient in mobile app or POSTMAN

Register Patient w/ app

- Once register button pressed, navigate to Policy screen
- First opening screen is called startup screen
 - Then go to EULA screen and Privacy screen
- Patient registration screen
 - In merlin portal
 - Add DOB and serial # to app, press next
 - Patient Information
 - Fill up details
 - Should be same as testing, but can be different
 - Enter own phone #
 - PW Test!123
 - Patient is registered to app, now login with email and PW

If try logging in with another patient, it won't allow as device can only be registered by 1 patient
Need to uninstall app and reinstall to register new patient

No policy screen when logging in, already created patient so skipped registration flow

Home Screen

Systematic BP

- Record BP for patient
- See same card on readings on global navigation panel
 - 3 Panels, day, week, month
 - Day is today
 - Week is entire week
 - Month is entire month, data not recorded for entire month -- value is present

DTP in Merlin

- Add medications in merlin, after adding goal PA Pressure
- Turn on Medication managed in App in Merlin
- Use save and add another button to add more than 1 meds

In app, Medication adjustment notification if medication plan is changed in Merlin

IOS has extra button in More options (Update on Apple Health)

More Options Tab

- Logged into App and on Home screen

- More Menu tab
 - Profile option
 - Full name in Merlin Portal
 - Email Address
 - Phone number
 - Can edit only full name and phone number
- As of 2.0.22, can change country code
 - If profile page updated
 - Go to Merlin portal
 - View registered application user
 - Only find updates here, not on the main portal
 - Support
 - FAQs, Contact, About, User Manual
 - FAQs have more info
 - Contact
 - Clinic contact details from merlin portal
 - Add primary lab,
 - Just for user not in merlin
 - About
 - Privacy Policy and EULA (info for user)
 - Has software version
 - GTIN - related to phone
 - App model number
 - User manual
 - Guidelines, how to find product info
 - Education
 - About HF
 - Videos about HF and CardioMEMS system
 - Learning (more info on everything)
 - Weight
 - Recording weight for patient
 - 3 tables: day, week, month
 - Each tab has record weight button
 - For other days, use POSTMAN to enter weight
 - User calendar icon to see previous weeks and months weights trends
 - Settings
 - Change password
 - For IOS, has update to apple health button there too
 - Whatever data in apple health app, it would be reflected here
 - Change password
 - Can change password, and it redirects you to login again

POSTMAN

Monday, August 17, 2026 11:06 PM

POSTMAN

- Workspace = CG App Verification team
- Many Collections
 - o AK - IAPSS collection for V3
 - o Q3 Working for Q3
 - o T1 Official for T1
 - For making Patient Active folder
 - Register collection registers patient
 - Can directly log into app, no need for registering in the app
 - Copy sensor serial number
 - Run Register HES Copy
 - ◆ V3 environment
 - ◆ Copy same as Merlin portal
 - ◆ Just change environment and hit send
 - Create Patient HES Copy
 - ◆ Copy Sensor Serial #
 - ◆ Replace # in pre-request > scripts tab
 - ◆ Fill in patient details
 - ◇ Date = Year Month Date, in Merlin, its month day year
 - ◇ Change phone #
 - Sensor Calibration
 - ◆ Baseline code in Merlin
 - ◆ When run 3rd file, fill up baseline code
 - ◆ Pre-request, fill out sensor serial #
 - ◆ Send
 - Initial HES reading
 - ◆ Patient status still inactive
 - ◆ Change same patient details and serial #
 - ◆ Send
- In Merlin, don't need to do hard refresh, just tap readings and it should work
- IRL, postman not real
 - o User postman to ease testing
 - Patient still inactive, use 5th file to make active
- Now patient is active after sending 5th file

DTP

- If patient is having high reading or low reading
- Clinician can send medicine
- Ex
 - o Goal = 17
 - o Low = 13 systolic
 - o High = 21 diastolic
 - o User = 15, normal reading
 - o User = 24, high reading
 - o User = 11, low reading
- Readings are for 1 week period
- Save interventions in merlin

- Once interventions are added, you can delete them by pressing individual intervention when blue selected
- Then can turn on dynamic treatment plan
 - Clinician is going to use Merlin
 - Clinician can turn on DTP for user
- Home Reading POSTMAN file
 - Go to Body
 - Change data section
 - Avg systolic = 22
 - Avg mean = 21
 - Avg diastolic = 20
 - Patient is in Active Patients tab in Merlin now
- Merlin time zone = California time
 - If clinic is in different time zone
 - Clinic = 5 AM
 - Triggered reading at 6:30 AM
 - Reading will not get posted
 - Time postings must be before clinic time
- Cannot revert back once send, once reading file is generated cannot revert back
- Suspect = high reading
 - When hit suspect, not trigger DTP
 - If trigger any other readings around it, always trigger suspect
 - Use another patient
- 2 DTP cards are triggered because of the 7 out of range readings

Dashed on the medication screen

- Some have dashes, some don't
- Company from the medicine

Sending Weight and BP

- Get B2C token API file, just send it
- Got to Login page
 - Send to get patient ID (patientLoginSrcId tag)
 - Edit email address for patient
- Post vital(weight, fluid) page to send weight
- Post BP Reading page to send BP reading
 - Mention patient ID
 - Mention specific weight
 - takenAtUtc tag
 - Epoch converter: <https://www.epochconverter.com/>
 - Mention date to get time stamp and a Unix value and paste into this tag
 - This is a millisecond value

Helpful Notes

- Always create patient first
- And enter it in pre request > scripts
- Send first reading to make patient active
- Body script is where edit info abt patient reading first
- Make the time a long time ago, just make a baseline
- Enter
 - Avg systolic

- Avg mean
 - Avg diastolic
- File 3 is sensor calibration
 - Create patient
 - Calibration code
 - Then first one reading to make patient active
- Suspect = if system thinks data doesn't make sense
 - Use this formula to avoid suspect readings: $\text{Mean} = (1/3)(\text{Systolic}) + (2/3)(\text{Diastolic})$
- Out of range = sensical value has been read, but it is out of range

POSTMAN CLI

Monday, August 17, 2026 11:07 PM

Generate POSTMAN API Key by going to Settings Icon > Account Settings > API keys

Authenticate your account:

postman login --with-api-key PMAK-6a23393a28180a0001a82455-f61161159a59a6b975f93aca763d5554e4

	Q3 Environment	T1 Environment
Postman Collect UID	52864576-d9174d04-011e-4d5f-aca9-0001b00910e6	32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316
Postman CLI environment ID	31616098-5883cfa9-36bd-46f2-9b69-f6443b2cef85	41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f
Register HES ID	52864576-d3c5eacc-6783-4863-b76b-5f9a7f24b244	32592187-52d7d30d-6920-44c7-8a93-53e71cfc1e7d
Create Patient HES ID	52864576-50086e72-309a-421d-98da-39a1e244a5e7	32592187-98ab1f66-82cb-4953-927f-a79f7f9bb30b
Create Sensor Calibration ID	52864576-ebb454ce-2f57-41cd-950a-8a5253994009	32592187-52dafb95-e31f-4053-aa4b-f07904d90f33
Create Hospital Reading ID	52864576-06b12575-0c09-4f79-9fc6-ce543bf81b70	32592187-255f16a7-a4a5-49fd-8dab-b3343524793a
Create Home Reading ID	52864576-47436003-a58d-4abf-8276-701281a3ceed	32592187-1c16b8d7-4448-43d0-a74c-5a6173b9e1b2
Get B2C ID		32592187-30c8e0fd-bbcb-4016-812b-fb5cd96bd382
Fetch Contacts ID		32592187-fbb546e9-5be9-4017-9958-ee5da03fbf13

"Register HES"

postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "Register HES ID" --timeout-request 30000 --verbose

In the output of Register HES command, you will see this output. The value inside uuid should be used for the remaining commands in place of PASTE_UUID_HERE. The value inside token should be used in place of PASTE_TOKEN_HERE

```
□ For Making Patient Active
└ 1. Register HES Copy

POST https://iapi.t1-hfcloudaks.cardiomemshf.com/v4/interrogators/?current_software=2022.1&uuid=3bdbdd7b-87f8-420e-b240-fb57decd595d&serial_number=HES_925680
201 Created * 734 ms time * 434 B ↑ 1.77 KB ↓ size * 9 ↑ 23 ↓ headers * 0 cookies
↓ application/xml * text * xml * utf8 * 277 B
<?xml version="1.0" encoding="UTF-8"?>
<interrogator>
  <id type="integer">238105</id>
  <uuid>3bdbdd7b-87f8-420e-b240-fb57decd595d</uuid>
  <serial-number>HES_925680</serial-number>
  <current-software>2022.1</current-software>
  <expected-software>2022.1</expected-software>
  <synced-at/>
  <token>9rYzknw4Nf27C6dWvBuY6qZtPPoxrArr</token>
  <environment>T1</environment>
</interrogator>
```

"Create Patient HES"

postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "CreatePatientHES ID" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "patient_first_name=EDIT_FIRST" --env-var "patient_last_name=EDIT_LAST" --env-var "patient_dob=EDIT_DOB" --env-var "patient_phone=EDIT_PHONE" --verbose

"Sensor Calibration"

postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "SensorCalibration ID" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "calibration_date=EDIT_CALDATE" --verbose

"Hospital Reading"

postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "HospitalReading ID" --timeout-request 30000 --env-var "_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --env-var "patient_first_name=EDIT_FIRST" --env-var "patient_last_name=EDIT_LAST" --env-var "patient_dob=EDIT_DOB" --env-var "patient_phone=EDIT_PHONE" --env-var "time_of_acquisition=EDIT_CALDATE" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --verbose

"Home Reading"

```
postman collection run "Postman Collection UID" -e "Postman CLI Environment ID" -i "HomeReading ID" --timeout-request 30000 --env-var  
"_cli=1" --env-var "br_hes_uuid=PASTE_UUID_HERE" --env-var "hes_token=PASTE_TOKEN_HERE" --env-var "sensor_serial=EDIT_SENSOR" --  
env-var "patient_first_name=EDIT_FIRST" --env-var "patient_last_name=EDIT_LAST" --env-var "patient_dob=EDIT_DOB" --env-var  
"patient_phone=EDIT_PHONE" --env-var "time_of_acquisition=EDIT_CALDATE" --env-var "atmospheric_pressure=760" --env-var "avg_systolic=  
19" --env-var "avg_mean=17" --env-var "avg_diastolic=15" --verbose
```

"Get B2C + Fetch Contacts"

```
postman collection run "32592187-45073f8f-b2b9-4df0-aad4-a1e6a423f316" -e "41856171-d9a206ed-e2b2-431f-89b0-0f04da6b571f" -i  
"32592187-906e0b7c-ae6b-4541-8c3f-4f7d21f64e80" --timeout-request 30000 --env-var "_cli=1" --env-var "sensor_serial=1B0130" --env-var  
"patient_dob=2020-01-01" --env-var "patient_name=Nav6 Androidt1" --env-var "email=nav6t1@gmail.com" --env-var "password=Test!123" --  
env-var "phone="+17813230226" --verbose
```

DRx and DTP

Monday, August 17, 2026 11:07 PM

Non-Drx Interventions

- Go to Merlin, select specific patient
- More actions button, add interventions
 - o There are 6 options (radio buttons)
 - Reading reminder
 - Lab visit (lab requires blood test, schedule lab visit within next 7 days)
 - Treatment plan reminder (due to BP changes)
 - Can send custom message
 - Notification badge says "New clinician message"
 - o These are all non drx interventions
 - Click on one, press send
 - Notification should be sent to the app
 - Can be marked as read
 - ◆ See in complete tab in app
 - ◇ Complete screen shows read cards
 - See different texts on different notification tabs
- Daily medications sent to app are also considered non drx intervention

Drx intervention are only the DTP notifications (for high and low readings)

- Can titrate daily medication (hold = stop medicine)
 - o Select daily medication, select quantity and frequency
- Can add message to call clinic or custom message
 - o These would be considered drx intervention now
 - Only more actions messages would be considered non drx
 - Daily medications will be non drx

Out of Home Elevation

- Out of home elevation is change in atmospheric pressure
 - o This would not send an intervention
 - o Would send a notification to the app, but no intervention
- Go to postman
- Go to Send home readings (5th file)
- In body tags, AtmosphericPressure
 - o By default, value is 753
 - o Out of home gets generated when 35 is added to it (788+)
 - o Irrespective of if the reading is within the valid range or not, out of home notification will be generated
- In previous builds
 - o Label of notification badge is "Out of home elevation"
 - o Now it is "Away from home elevation"
 - o Even text is different
 - o Later on no badge will be on the home screen, everything will be on the notification screen
- Go to Merlin portal, go to activity tab
- Go to readings
- Out of home elevation status is there

7 Out of Range Readings

- Set optimal ranges
- Interventions for 7 high readings (3 total things)
 - o Adding a new medicine
 - o Asking to hold daily medicine
 - o Adding message to call clinics
- DTP plan must be on
- Go to postman
 - o Refresh postman every time you change patients
 - Sometimes POSTMAN overwrites, refresh before working on new patient
- Use 5th file to send high readings
 - o Enter patient details in scripts (pre-req)
- Log into app
- See week (May 13-May 19)
- Trigger 1 low or 1 high
- Create intervention in merlin
- Create readings for May 18, 17, or any of them
 - o No specific dates needed
- For every test, create new patient (for every new intervention case)
 - Create 6 patients
 - Test it individually for each intervention
 - o Once you trigger a intervention
 - o Cannot trigger another one for the same patient on the same day
 - o Triggering right now for 7 high
 - o Tomorrow need to trigger other readings, then you can get another intervention
 - o One day, only 1 intervention can be triggered
 - o Tomorrow you can trigger something different for this particular patient
 - o Cannot just keep triggering interventions
- Keep same values, change the date
 - o Have to be in a row
 - o Mentioned 3 interventions
 - Sends 3 separate notification badges
 - Text is different for each drx notification
- Ask user to confirm that they are currently taking old medication before accepting that they stop medicine
 - o Now earlier medication is not in the medicine list tab
- In Merlin portal
 - o Can see intervention cards for each DTP intervention
 - 3 cards, 1 for each intervention

Helpful Notes

- Don't create interventions for all 6 cases
- Only create one for the particular one you are testing
- If have 1 low and 3 low interventions
 - o Send 1 low readings few days ago
 - No intervention was triggered
 - Even 3 low readings, no intervention was triggered
 - This is because the past intervention is generated for TODAY
 - ◆ Notification was sent TODAY, even if reading was taken in past
- 1 intervention for 1 day
 - o 24 hour waiting period

Setup Notes

Tuesday, June 2, 2026 10:30 AM

Appium server is different from appium inspector

Install packages and softwares directly into C drive, not in Users

Always keep appium running in background when running automated test cases

- Run command *appium* in cmd or IntelliJ IDEA before starting session in Appium Inspector

Always work on you own branch first (make separate branch and make edits, do not work directly on the main branch)

Locator = used to identify specific element on webpage

X Path = XML path language

JSON To paste In Appium Inspector

```
{  
  "platformName": "Android",  
  "appium:platformVersion": "16",  
  "appium:deviceName": "Samsung Galaxy A56",  
  "appium:udid": "RZCY826JDLJ",  
  "appium:automationName": "UiAutomator2",  
  "appium:appPath": "C:\\work\\myCardioMEMS-q3-2.0.23.apk"  
}
```

Mvn clean install -DskipTests

- Downloads all libraries
- Should say BUILD SUCCESS

Update branch with Global changes


Troubleshooting

Wednesday, June 17, 2026 12:34 PM


If you have issues with Springboot and automated OTP, remember that IP address changes based on WiFi. Change SMS forwarder IP address based on ipconfig or ifconfig results

Errors in sanity.xml


Go to settings > build, execution, deployment > Maven > Importing
Make sure JDK Importer is the one installed on computer

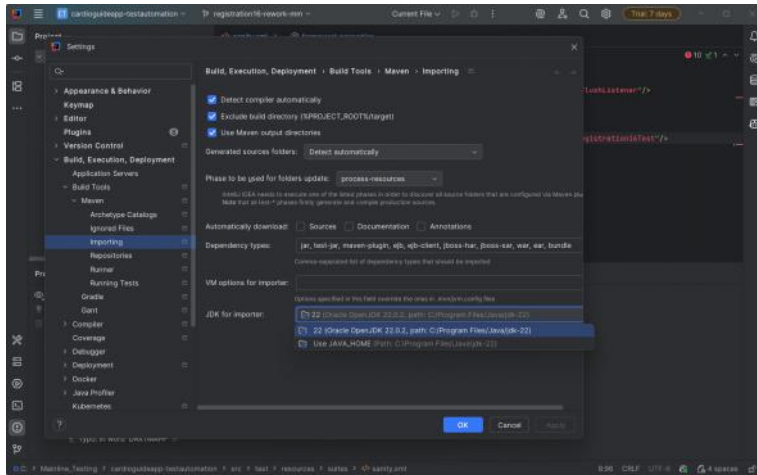


Go to File > Project Structure > Project
The SDK should be the one installed on computer

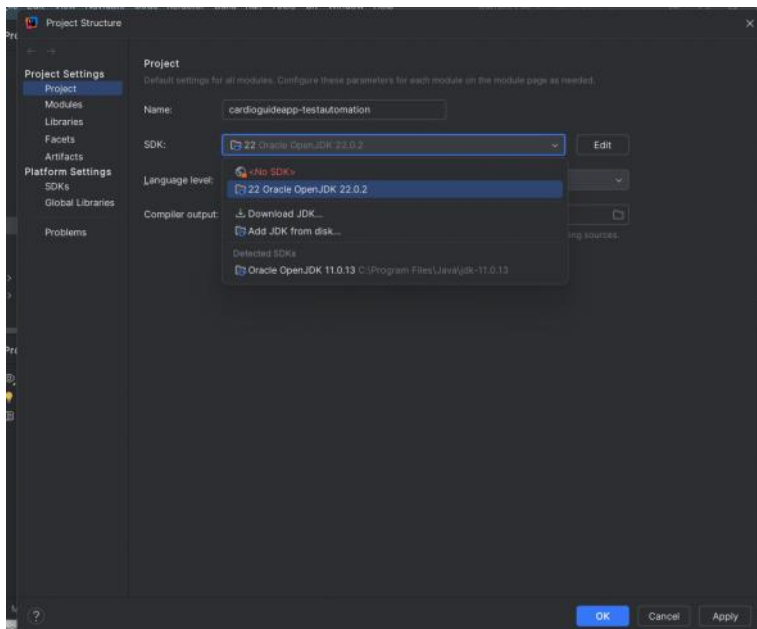

Go to Project Structure > Modules > +> Import Module
Import the pom.xml file from the cardioguide-testautomation folder

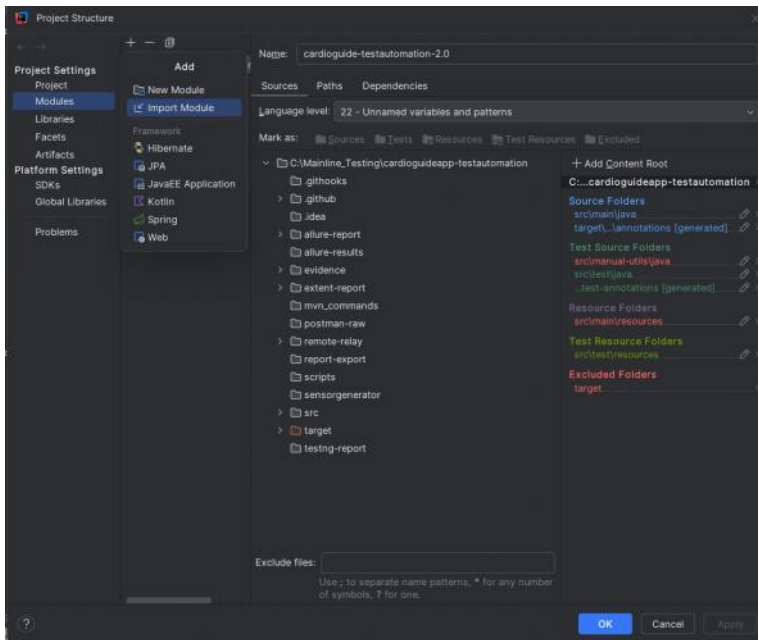

On the left side, open the "m". Press reload all maven projects

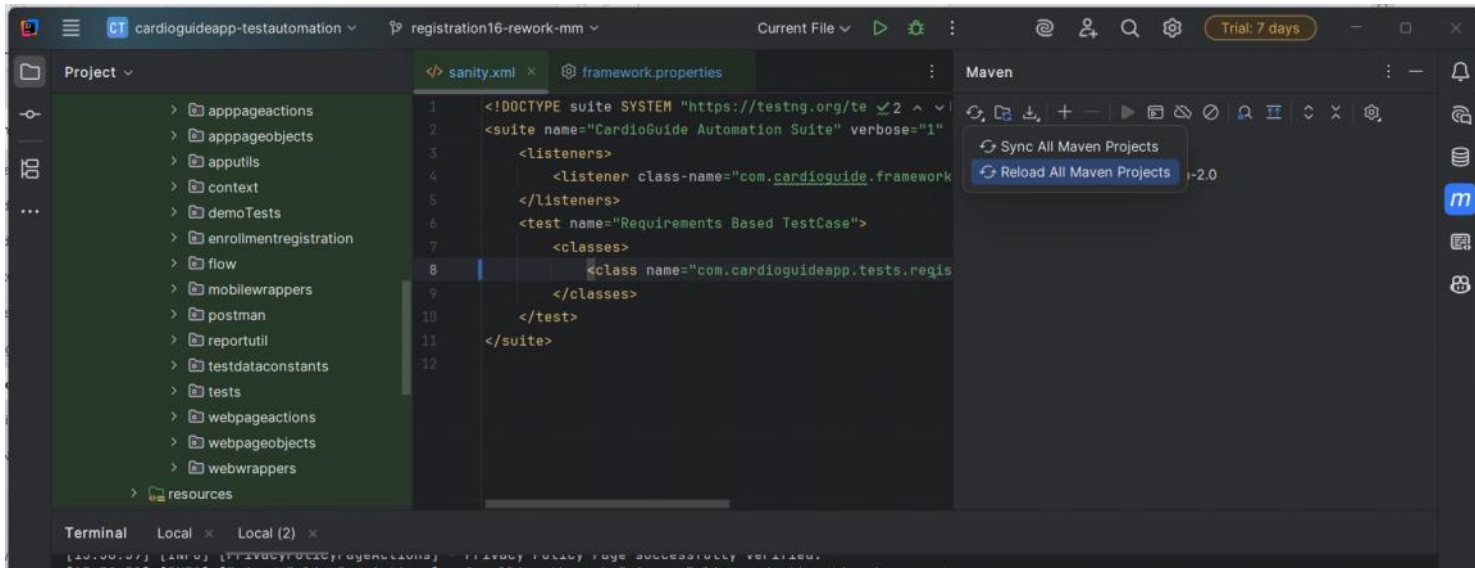

Manual Testing

Monday, August 17, 2026 10:26 PM

Requirement testing

- When you are testing the requirements/specifications of the application
- Test cases are made based on the requirements, and goal is to check if they were satisfied or not

Ad Hoc Testing

- Not with any planned workflow
- Randomly take system
- Randomly play with system
 - o Break the system and find bug
- Do whatever you feel like
- Not following any requirements

HERO Setup

Monday, August 17, 2026 10:58 PM

Connect Device

- Waveform Cord and pins
 - o Put 2 pin thing on end of 6 pin array
 - o GND goes on black pins, orange mid, yellow
 - o Black pins on J36 blank
- 12 V DC Supply
- Ethernet and USB Cable
- Modem
- Put waveform generator cable after hero is boot up

Setup

- Ttermpro
 - o Serial (usually COM3)
 - o Setup > Serial > Speed = 115200 > New Setting
- Waveform generator
 - o Need to install national instruments labview and runtime engine
 - Download runtime engine first then reboot then the shorter name one
 - o Extract pressure waveform generator and it should run
 - o Install driver for waveform cable (COM waveform titled folder)
 - The waveform cable should be found in the PORTS section, not USB Section
 - In device manager, rename COM port for waveform cable as COM 15 or above

Waveform generator replaces patient readings

Automation Execution

Monday, August 17, 2026 11:32 PM

Xray Client ID: 86A1D8829C494260B89DF497E5464EE1

Xray Client Secret ID: 7377cf7ff35b39703a1021fdecfa976423736198928d06b9128f1ed820acc476

```
mvn clean test -DplatformName=Android -Ddevice.udid=<Device udid> -Dapp=<App .apk path> -  
Dxray.enabled=true -Dxray.cloud.clientId=<client ID> -Dxray.cloud.clientSecret=<client secret> -  
Djira.api.token=<JIRA Token> -Dxray.test.plan.key=DRXTMAPP-5648 -  
Dxray.execution.include.datetime=false -Dxray.execution.include.counter=false -  
Dxray.execution.prefix="LMR - Dry Run 3" -Dxray.project.key=DRXTMAPP
```

```
mvn clean test -DplatformName=Android -Ddevice.udid=f17437eb -  
Dapp=/Users/BODANSX/Downloads/29t1.apk -Dxray.enabled=true -Dxray.cloud.clientId=  
86A1D8829C494260B89DF497E5464EE1 -Dxray.cloud.clientSecret=  
7377cf7ff35b39703a1021fdecfa976423736198928d06b9128f1ed820acc476 -  
Djira.api.token=ATATT3xFfGF08IXeLMx7Vmf9U_LVX2zME1Dj5If01XI_zZTEth55FQAESMyErn9UMwb-  
HLUBAMU9y9ADcvc9JDz4CSK5dXThqtLIWid_  
1uphhm7IleAsvuWDMlexRoQAbw4F8Fy2wsm_aHXA5nn7YK_pOYYqojZvKXuA0brt3V53oHbkiZFvtE8=  
013D8E57 -Dxray.test.plan.key=DRXTMAPP-6225 -Dxray.execution.include.datetime=false -  
Dxray.execution.include.counter=false -Dxray.execution.prefix="LMR - Dry Run 3" -  
Dxray.project.key=DRXTMAPP
```

Get JIRA API Token from: <https://id.atlassian.com/manage-profile/security/api-tokens>

Bug Testing

Monday, August 17, 2026 10:59 PM

Bugs in Jira

- Filters, search work filters
- 3 spaces
 - o TestManagementApp = test cases
 - o SRS = requirements for project
 - o DevApp = dev builds (bug for app)
 - o DevPortal = Portal team (bugs related to portal)

App bugs (DevApp)

- Type = Bugs
- Status = Verify
 - o Bug has been already logged, dev team has checked bug, and now it needs to be verified
- Verify, view Workflow
 - o When people are working, status is reading and then in progress
 - o Once fixed, status is verify
 - o Once verified, assigned to Danny and we should be retesting the bug
- Labels, Not equal to Ready to Close

Example Bug

- Not been retested yet
- Revisit all the test description and read it all
- When understand bug, retest it in new build
- Perform steps that were mentioned to reproduce test
- Go to actual results and compare to expected results, bug is logged
- Make sure expected results are met
- Write up comment
 - o Environment details, what the build is, and what phone it is tested on
- Status is resolved or not
- Write expected and actual results
- Add evidence in comment s
- Once comment is added, add ready to close button

If bug is still unresolved

- Put up comment as status being unresolved
- Add second comment and reassign bug to original person, telling them bug is still not resolved
- Put it back to in progress

If want to log up bug,

- Log bug in Dev App space
- Press create
- Work type as Bug
- Add summary of what the bug is about
- Have pop up, this bug has been created
 - o Press it and land directly on ticket
 - o Enter description

Example Bug Tested (Resolved): [\[DRXDEVAPP-9612\] Env Q3 - Broken String Header in Third Party Licenses Panel - Jira](#)

Example Bug Tested (Un-resolved):

[\[DRXDEVAPP-7006\] Env V3-Notification Badge on CardioGuide App Icon Disappears After App Launch With...](#)

Example Comment:

Defect Status: Resolved

The MFA Phone error message is displayed after 1 minute 40 seconds after user fails to complete Call Me MFA, thus satisfying the Expected result.

Details

Environment: Q3

Build Details (Android): Version: 2.0.28

Clinic Details: Username: q3_swvq3, Password: Test!1234

User Details: Sam15 RTR, Username: [sam15@gmail.com](#), Test!123

Device Details: Samsung Galaxy A56, Android 16

Evidence:

Logging a Defect

Tuesday, August 18, 2026 12:10 AM

Loggin Bug in Jira

- Observation issues while testing (devApp space)
 - o Goes to dev team
 - o Create
 - Select space, devApp
 - Work type = bug
 - Status= open
 - Summary = what is the us
 - Env V3 - [what bug is]
 - No need for components
 - If know where the component is, can add that
 - Description = min thing
 - Environment in which the test was done
 - Build details
 - Test details
 - Clinic details (for merlin portal)
 - User details (what user was used to login)
 - Device details (android, IOS and update, and phone model)
 - Steps to reproduce
 - ◆ While doing what steps issue was identified
 - Actual result = what you observed
 - Expected result = what you should have seen
 - Evidence (screenshots or video)
 - ◆ Just result you observed, not every step
 - ◆ Buttin to add image and photo
 - ◆ Attach In description
 - ◆ When single attachment is added, then can add on plus signa nd directly add it
 - ◆ No attachment option if you just create a bug
 - ◆ Only added when added a screenshot/video from description
 - Assign is automatic, if known can put it
 - Severity option too, low high med
 - Labels
 - ◆ CG2.0 Defects
 - Affected version = Q#APP - 2.0.21
- Nees to create bug
- Space = drx_devapp
- Test cases is drs_testmanagementapp
 - o DRXTMAPP
- Bug retesting
 - o Reproduce it with steps mentioned
 - o In comment section
 - o Verify status
 - o Reproduce using states with new build with the new build
 - o In comment section
 - IF NTO SISSUES, AND IT WAS RESIVED, DO RESOLVED STATU

- TESTED ON, ON WHICH UPLOID AND TDATE UPDATE
 - Env tested
 - Device tested on
 - Merlin portal version as well
 - Help about, find merlin version
- Clinic details , and full name fir the patient, gmail, and password
- Add evidence that it is fixed ow
- Assign it to ourself, then add to comment and assign to danny again
- Label
 - Ready to close label
 - Danny knows bug is fixed

Example Bug Logged: [\[DRXDEVAPP-9368\] Env Q3 - Red Notification Badge on CardioMEMS Application Icon do not reset the count even after the Notification Cards are acknowledged. - Jira](#)

DOORS


Tuesday, August 18, 2026 12:55 AM

DOORS

- DOORS = Dynamic Object-Oriented Requirements System
- Requirements in CardioMEMS > Next Generation Products > Software
 - o App110 = CardioGuide App
 - Requirements
 - ☐ Have SRS and GUI Specification
 - ☐ SRS are the requirements
 - ☐ Test cases are created based on requirements
 - ☐ DRxRTR is a valid req for anyone NOT CLOUD
 - ☐ Sometimes req team modify reqs, then the description or ID would be modified
 - GUI Spec
 - ☐ Create test cases or update test cases based on GUI Spec
 - ☐ Name of screen
 - ☐ In test case step, provide name of screen
 - ☐ Be more specific based on GUI specific
 - ☐ Important for updating or creating test case
 - String export doc
 - ☐ Messages and strings on application

In DOORS

- CardioMEMS > Next Generation Products > Software > Configuration
 - o File for all Strings
- CardioMEMS > Next Generation Products > Software > Requirements > GUI Spec file
 - o Panels, components for app
- CardioMEMS > Next Generation Products > Software > Requirements > SRS file
 - o Requirements for App

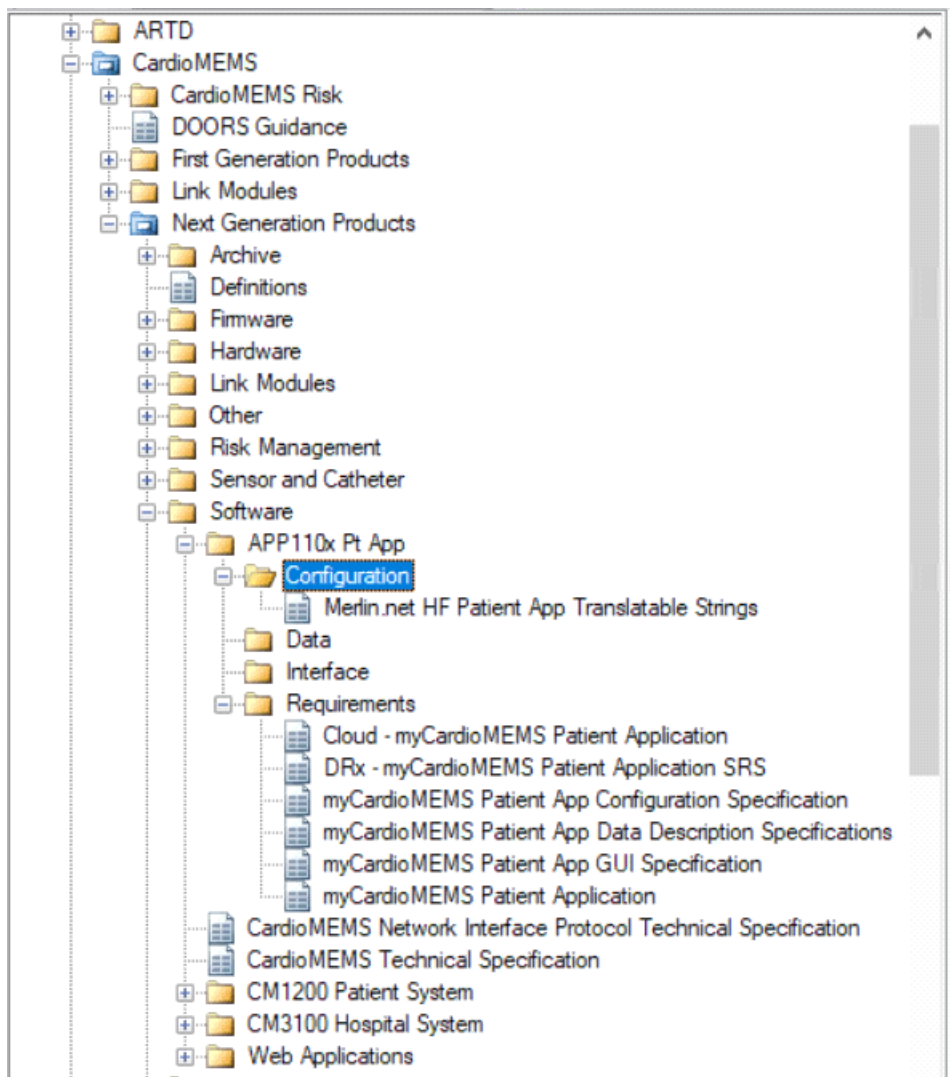

Links to KT Documents

Monday, August 17, 2026 10:27 PM

HERO and CMEMS Onboarding Document: [HERO On Boarding Document - Wolf Ridge SW - Meddev Confluence](#)

CardioGuide Automation

CG Android framework setup: [CG Test Framework Setup Guide \(Android\).docx](#)

CG Mac framework setup: [CG Test Framework Setup Guide \(Macbook\).docx](#)

Mainline Testing/Running Test Cases: [Mainline Testing Document.docx](#)

Springboot Setup for Automated OTP: [Springboot Setup Document.docx](#)

CardioGuide Test Automation NPS (Non Product Software) Documentation and Test cases:
[XXXXXXX CardioGuide Test Automation NPS Documentation.docx](#)

Miscellaneous

Setting up PTC: [KT_PTC.docx](#)

Setting up GitHub Copilot with IntelliJ: [KT Github copilot.docx](#)

Using Device Farm: [testfarm tool-automationtesting.docx](#)

Setting up Waveform Generator for HERO: [90662299 CardioMEMS Waveform Generator SW Executable Cover Sheet \(1\).docx](#)

Running Formal Manual Executions: [RTR LV Execution Guidelines.docx](#)

Using Merlin, Patient Registration

All things CardioGuide (Unfinished): [Merlin Patient Registration.docx](#)

Getting Started with CardioGuide: [CardioGuideIntroAndSetup.docx](#)

Keyboard Commands

Monday, August 17, 2026 10:27 PM

1. Macbook Specific Comands

- a. Copy file as path: Right click file in Finder, press and hold the Option key, and select Copy [File Name] as Pathname
- b. Go to folder path: Shift + Cmd + G, and paste the path of the file
- c. On Xcode test build with Cmd + B

2. IntelliJ commands

- a. File search on Windows: Ctrl + Shift + N
- b. Multiple Line comment on Windows: Ctrl + Shift + /
- c. File search on Mac: Cmd + Shift + O
- d. Multiple Line comment on Mac: Cmd + Options + /

3. Screenshots

- a. Windows: prt sc (F10) or Windows + Shift + S
- b. Mac:
 - i. Capture whole screen: Cmd + Shift + 3
 - ii. Capture certain part of screen: Cmd + Shift + 4
 - iii. Screen recording: Cmd + Shift 5

Git Commands

Monday, August 17, 2026 10:34 PM

Git branch

- Lists all local branches and highlights the one you are currently in

Git branch <name>

- Creates a new branch with that name
 - o Git branch SamhithaBranch

Git stash

- Saves local changes (does not save them to main branch, just the one on local computer)

Git stash drop

- Permanently deletes saved changes from history
- You must do git stash first before doing git stash drop

Git pull

- Downloads changes from repo to current local branch

Git pull origin <branch name>

- Downloads changes from a different branch on top of the branch you are currently on

Git clone --branch CardioGuide2.0-Web-API_Integrations

<https://bitbucket.sjm.com/myc/cardioguideapp-testautomation.git>

- Clones a local copy of the CG test automation framework

Git checkout -b <name>

- Git checkout Allows you to switch between different branches
- Git checkout -b allows you to create a new branch and switch into it

Making a PR

Git status

- Can see all your files that were modified
- Red colored files means they have not been staged/committed
- Green colored files means they have been committed

```
Changes to be committed:
(use "git restore --staged <file>..." to unstage)
    deleted:    sensorgenerator/sensor import final.csv

Changes not staged for commit:
(use "git add <file>..." to update what will be committed)
(use "git restore <file>..." to discard changes in working directory)
    modified:   src/test/java/com/cardioguideapp/apppageactions/NotificationsPageActions.java
    modified:   src/test/java/com/cardioguideapp/tests/notification/DRXTMAPP5734Notification9_5Test.java
    modified:   src/test/resources/StringLabels.xlsx
    modified:   src/test/resources/androidmobilelocators/t1/LandingPage.properties
    modified:   src/test/resources/androidmobilelocators/v3/LandingPage.properties
    modified:   src/test/resources/framework.properties
    modified:   src/test/resources/iosmobilelocators/q3/LandingPage.properties
    modified:   src/test/resources/iosmobilelocators/t1/LandingPage.properties
    modified:   src/test/resources/iosmobilelocators/v3/LandingPage.properties
    modified:   src/test/resources/suites/sanity.xml
    modified:   src/test/resources/testdata/DRXTMAPP5734Notification9_5Test.json

Untracked files:
(use "git add <file>..." to include in what will be committed)
    sensorgenerator/sensor import final.csv
    src/test/java/com/cardioguideapp/apppageactions/MoreOptionSettingsPageActions.java
    src/test/java/com/cardioguideapp/apppageobjects/MoreOptionSettingsPage.java
    src/test/java/com/cardioguideapp/demoTests/DebugHealthAccessTest.java
    src/test/java/com/cardioguideapp/mobilewrappers/MoreOptionSettingsPageActionsTest.java
    src/test/java/com/cardioguideapp/testdataconstants/DRXTMAPP415Registration13TestData.java
    src/test/java/com/cardioguideapp/tests/registration/DRXTMAPP415Registration13Test.java
    src/test/resources/iosmobilelocators/q3/moreOptionSettingsPage.properties
    src/test/resources/testdata/DRXTMAPP415Registration13Test.json
    src/test/resources/testdata/DebugHealthAccessTest.json
```

Git add <path of file>

- Adds files to stage

Git commit -m "message"

- Commits all the files you added

Git push

- Pushes your latest commit to your remote branch

Git rm --cached <path of file>

- Removes a file that has been pushed, but does not remove it from your local branch
- You need to do git commit and git push for your remote branch to be updated

Git reset --hard origin/<branch name>

- Resets local branch to the same as the remote branch

Git push -u origin <branch name>

- Push to a branch for the first time, can use git push after

Git restore <file>

- Discards changes to a file in local branch
- If file is committed but not pushed, must do git restore <file> twice

Git checkout CardioGuide2.0-Web-API_Integrations -- <file path>

- Replaces the file of the local branch with the file from the CardioGuide branch